


Vitafit


Smart Medical Grade Floor Scale

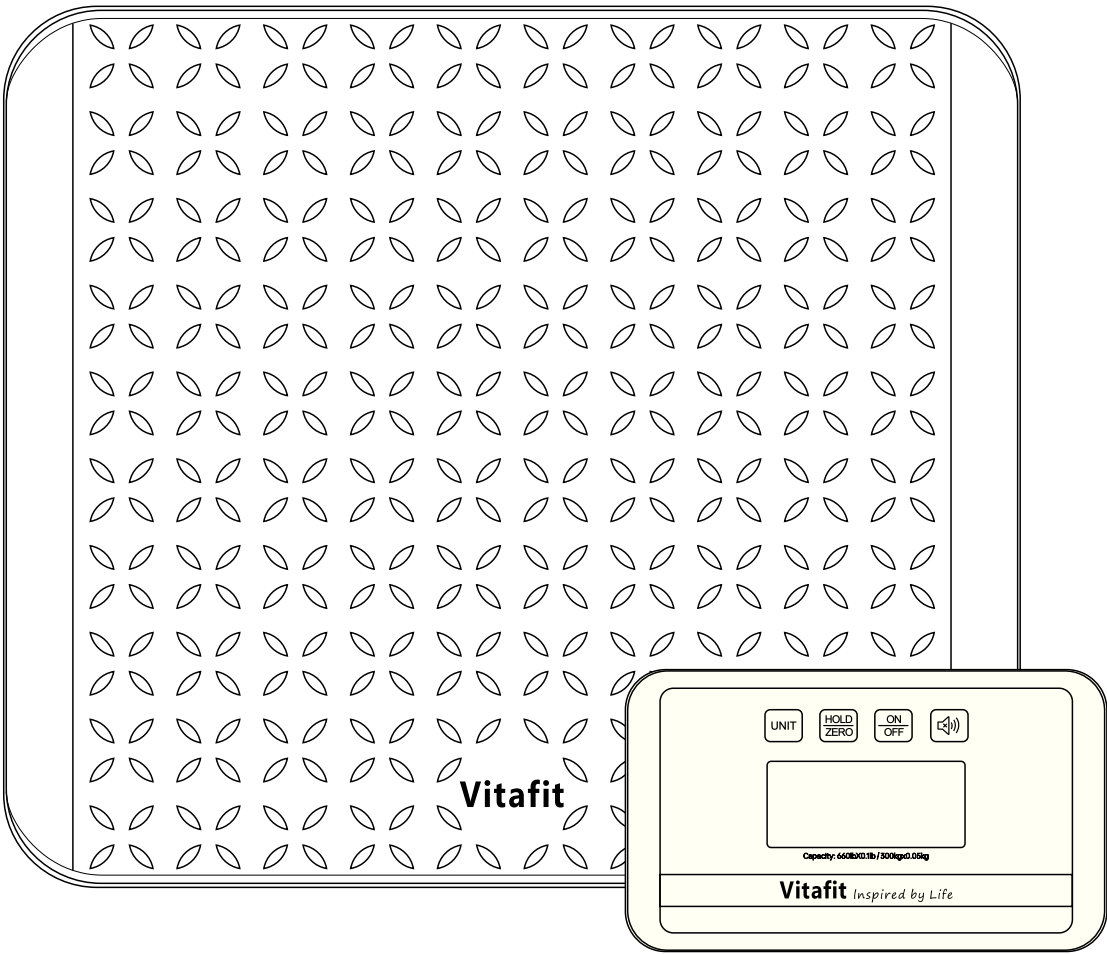


Model:VT739

User Manual

Please read this User Manual carefully before using the device

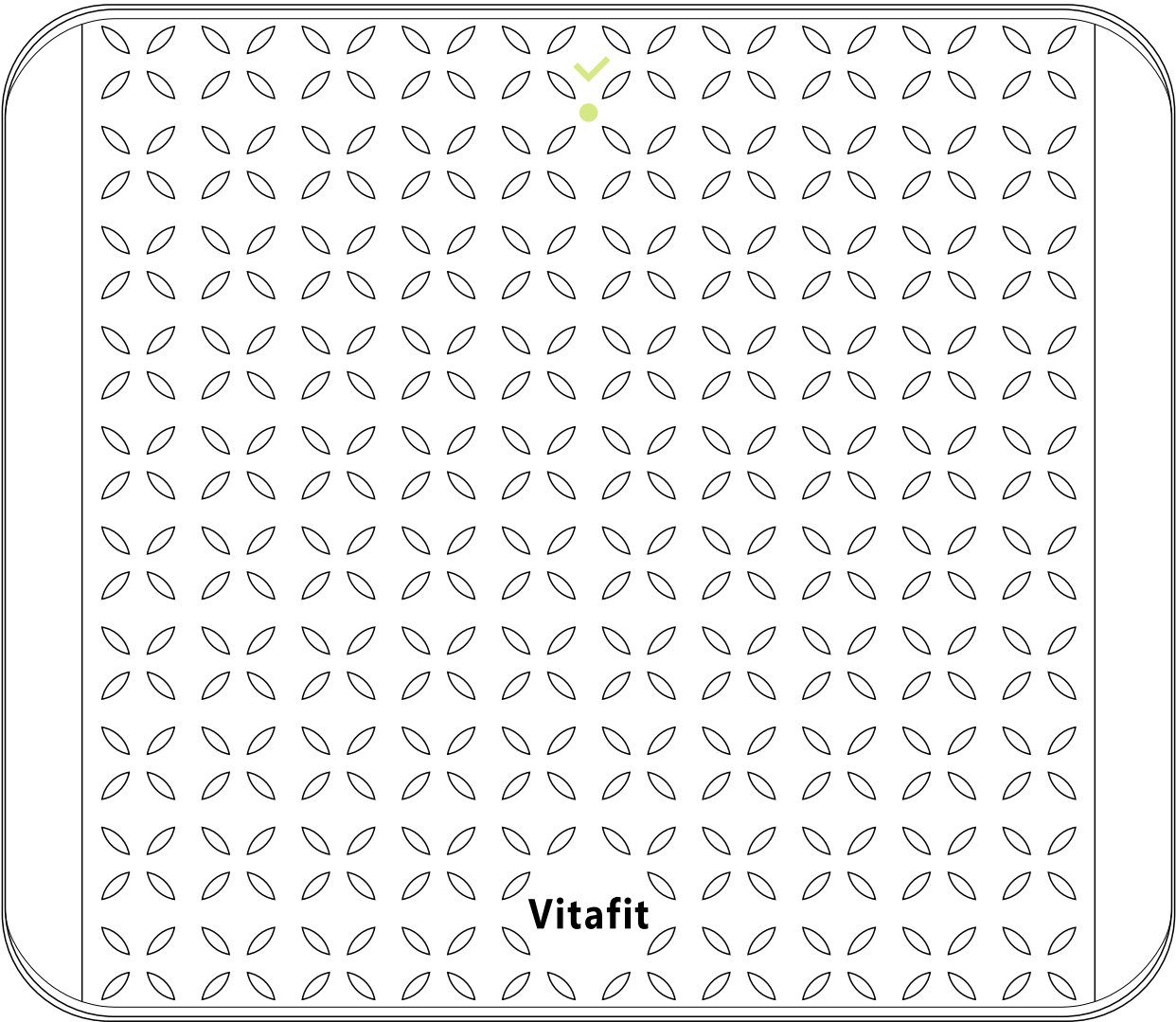

▲ Scale Base

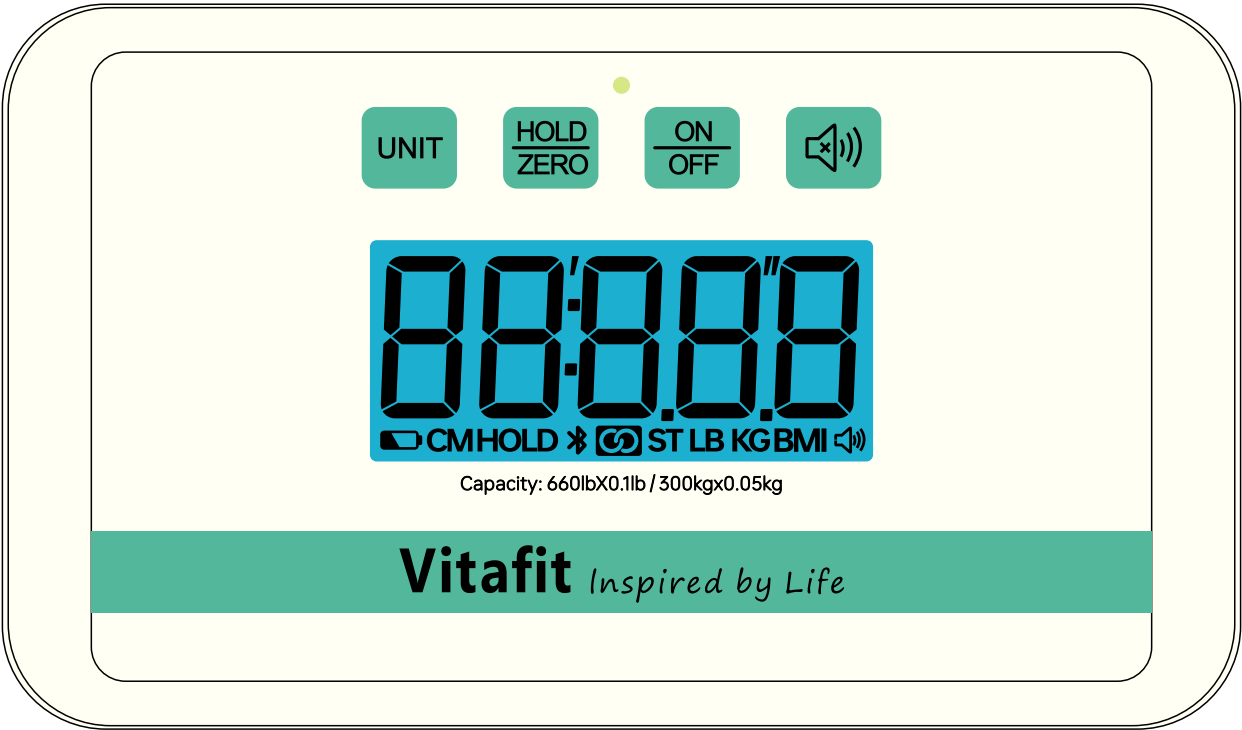

▲ Remote Controller

Note: Please ensure that the scale is placed on a hard and flat surface before power on.
The green icons " ✓ " and " ● " in the picture are the light indication diagrams.

Core Indicator Function

Indicator	Location	Description
	Scale base	When both icon of scale base flashing alternately, it indicates that the scale base and remote controller is connecting.
	Scale base and remote controller	When the scale base and remote controller flashing in green, it indicates that the scale is ready for weighing. All in all, flashing for weighing, no flashing no weighing.
	Scale base	When the scale base keep flashing, it indicates that weighing in progress. When keep lighting, it indicates that weighing completed.
	Scale base and remote controller	When red icon flashing on scale base , it indicates that the scale base is low battery; When red icon flashing on remote controller , it indicates that the remote controller is low battery. Please replace it in time to avoid the risk of leakage caused by over-discharge.

Brief Operation Procedures


Note: All operations are controlled by the remote controller button.

Function	Button	Description
Power on		1. Press the button  to power on the scale base and remote controller
Connecting	NA	1. After power on, when the icon  and  flashing on screen of remote controller alternately, it indicates that scale base and remote controller is connecting; When the "0.0" displayed and the  icon flashing , it means the connection is successful, and ready for weighing. (TIPS: When green icon  flashing, It's ready for weighing; When green icon  off, it's not for weighing.)
Weighing	NA	1. Step on the scale base to weighing, stand still until get stable weight and the result displayed on screen.
HOLD ZERO		1.Press button  to hold the weight, and HOLD icon will appear on the screen. 2. Under the hold function, the weighing function is disabled, meantime the  icon on remote controller is off. If without any operation, the hold value will remain on the screen until the scale is automatically turned off. 3. Press the button  again to clear the weight, then you can step on the scale to weigh again.
Unit conversion		1. Press the button  to convert unit among kg, lb and st. 2. Remark: When the scale connected with app, please convert unit via app, and the unit button function of remote controller is invalid.
Sound		1. Press the button  to turn on or turn off the sound function. Remark: The  icon is light on the screen when sound function on, when sound function off, this icon is also off.
Power off		1. Press the button  around 3 second to power off the scale. 2. When the screen turn off, it means power off completed.
Automatic Power off	NA	1. When without any operation around 60 seconds, scale will power off.


IOS QUICK START GUIDE




1. Search" Vitafit" via Apple Store.

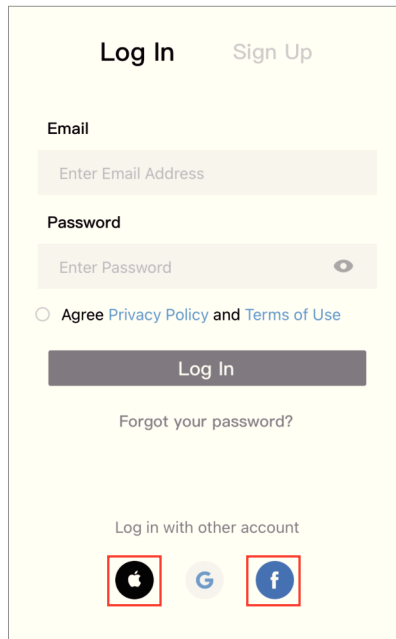


2. Convenient login by Apple or Facebook account.

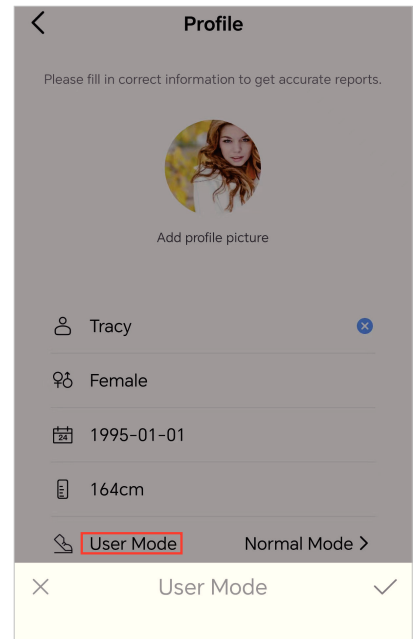


3. Choose suitable Mode. Note: Weight Only Mode specialized for pregnant or having a pacemaker, etc.

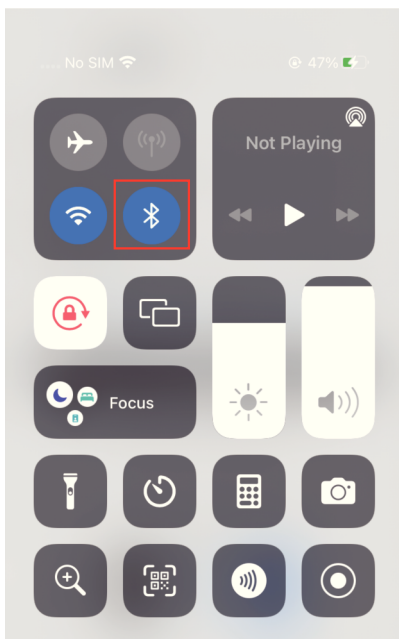

4. Turn on Bluetooth on your phone.

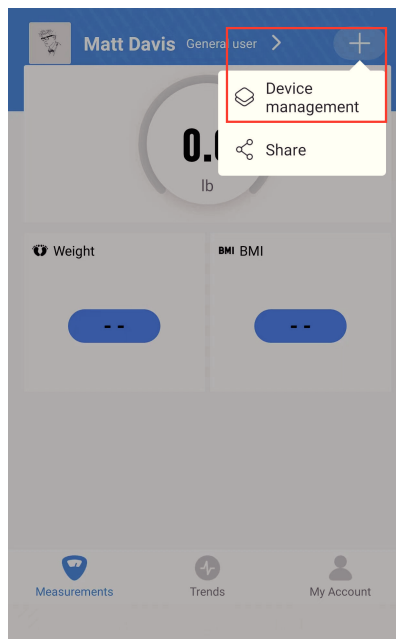

5. Pairing with Vitafit Scale.

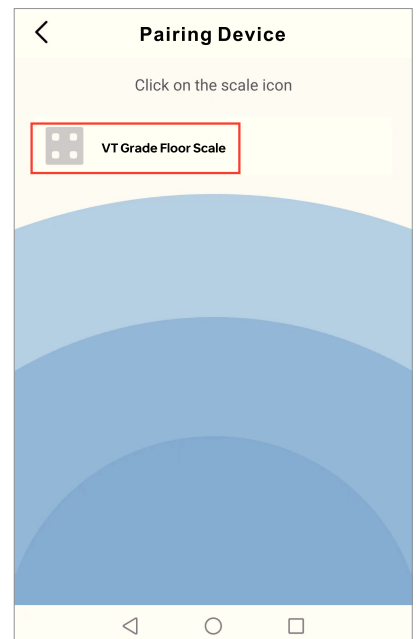


Notes:

1. For body fat measurement, make sure your feet bare and dry when measuring.
2. For excellent experience, recommend to reassemble batteries for the first beginning use.
3. UI will be different depended on scale type: Smart Weight, Smart Body Fat, VT Grade Floor Scale.


IOS QUICK START GUIDE




6. Turn on Vitafit App and step on for measuring.

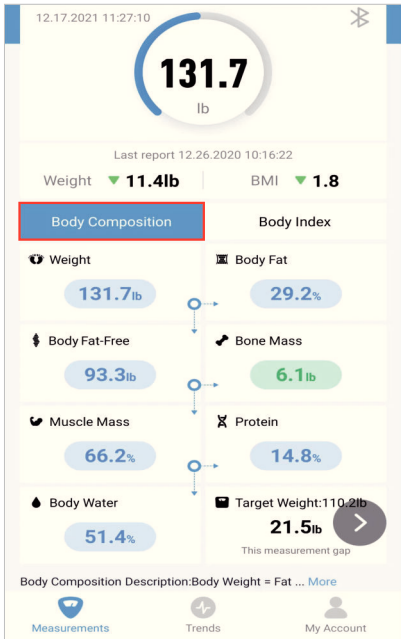



7. Measurement data.

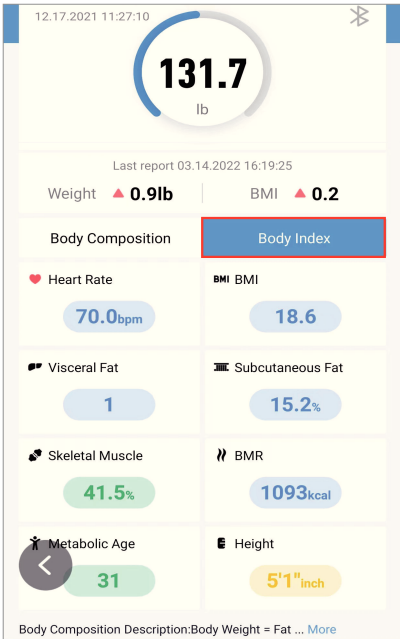

8. Set up for sync with Apple Health.

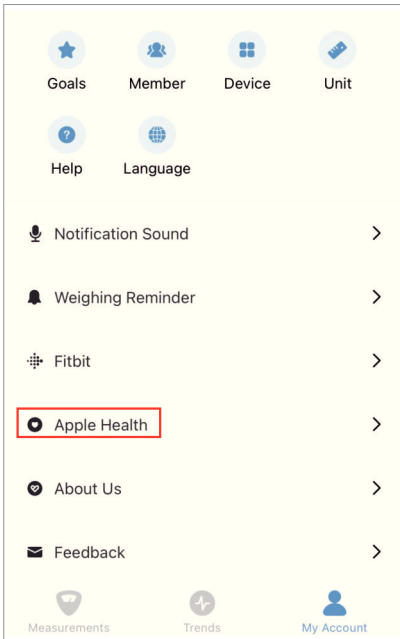


9. View measurement data in Apple Health after sync.

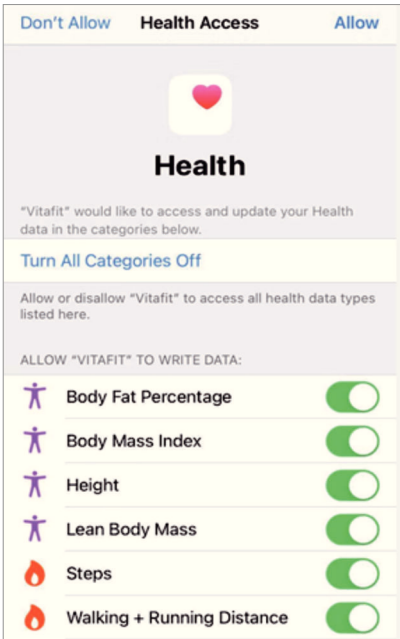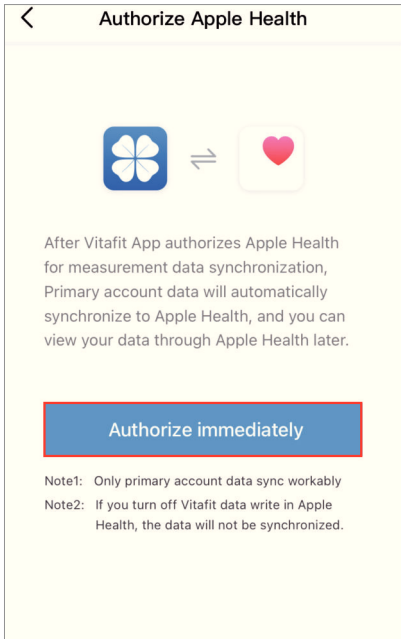


Notes:

1. For body fat measurement, make sure your feet bare and dry when measuring.
2. For excellent experience, recommend to reassemble batteries for the first beginning use.
3. UI will be different depended on scale type: Smart Weight, Smart Body Fat, VT Grade Floor Scale.


ANDROID QUICK START GUIDE

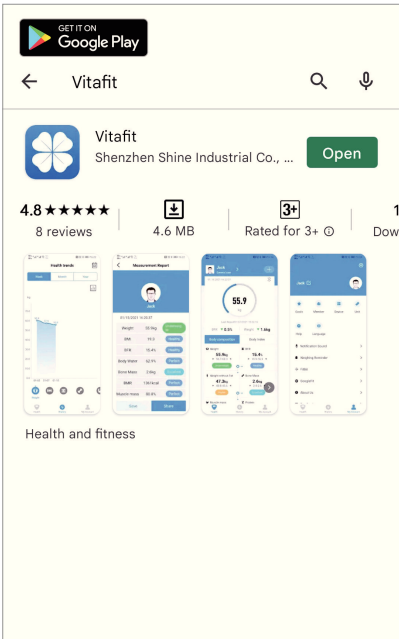


1. Search” Vitafit” via Google Play Store.

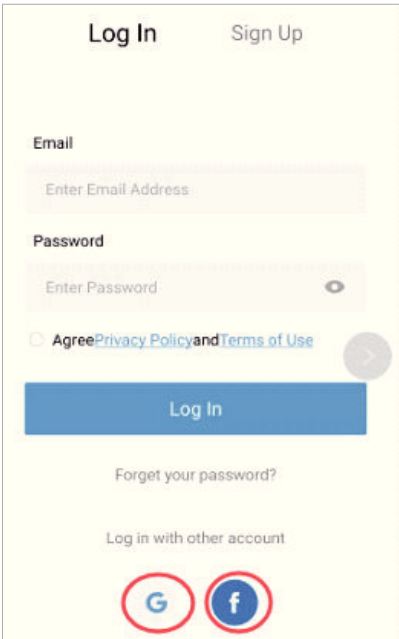


2. Convenient login by Google or Facebook account.

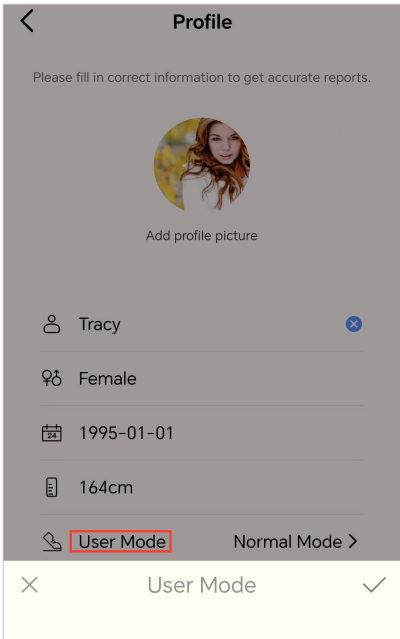


3. Choose suitable Mode. Note: Weight Only Mode specialized for pregnant or having a pacemaker, etc.

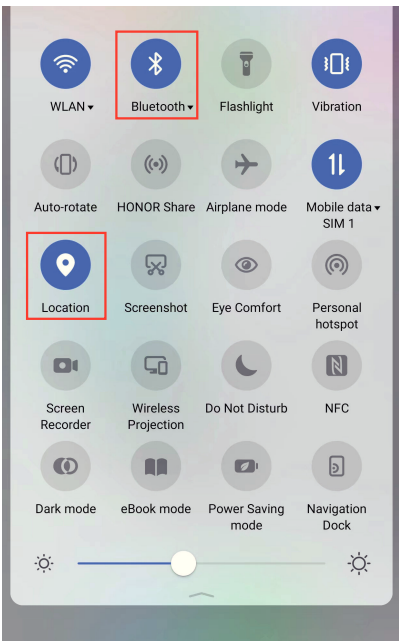

4. Turn on Bluetooth and Location on your phone, **OR ELSE THE APP DON'T GET RESULTS ACCORDING TO ANDROID REGULATIONS.**

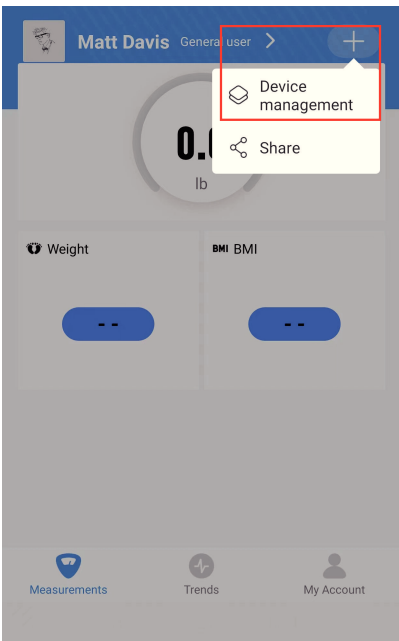


5. Pairing with Vitafit Scale.

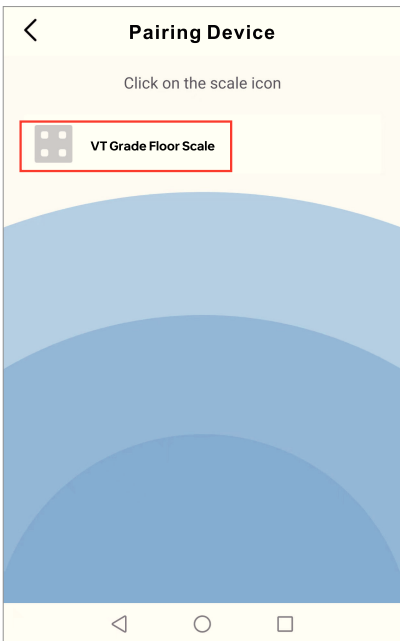


Notes:

1. For body fat measurement, make sure your feet bare and dry when measuring.
2. For excellent experience, recommend to reassemble batteries for the first beginning use.
3. UI will be different depended on scale type: Smart Weight, Smart Body Fat, VT Grade Floor Scale.



ANDROID QUICK START GUIDE

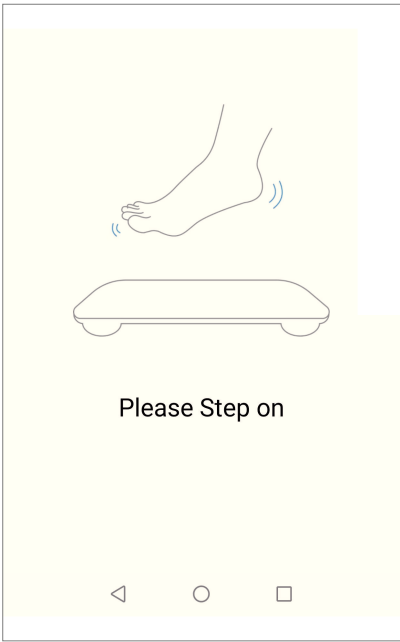


6. Turn on Vitafit App and step on for measuring.

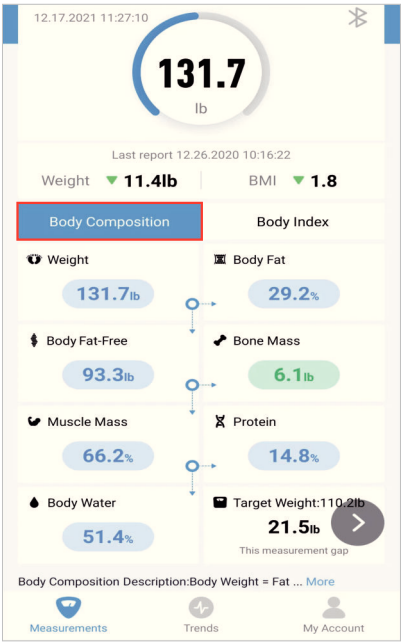


7. Measurement data.

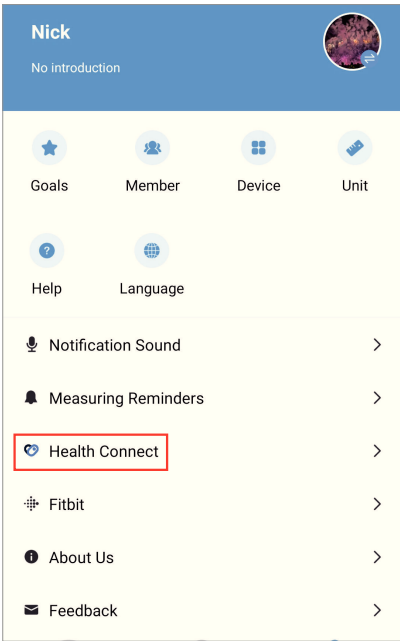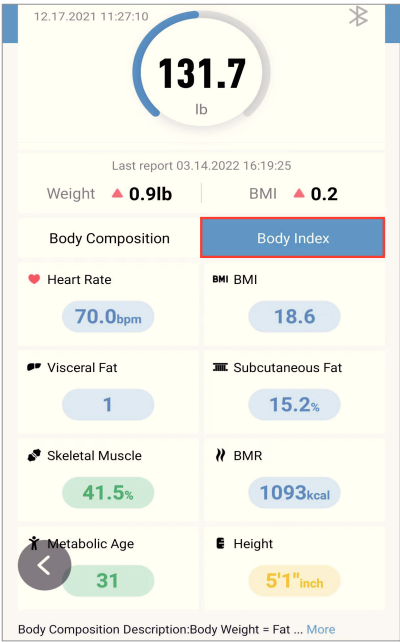

8. Set up for sync with Google Fit.

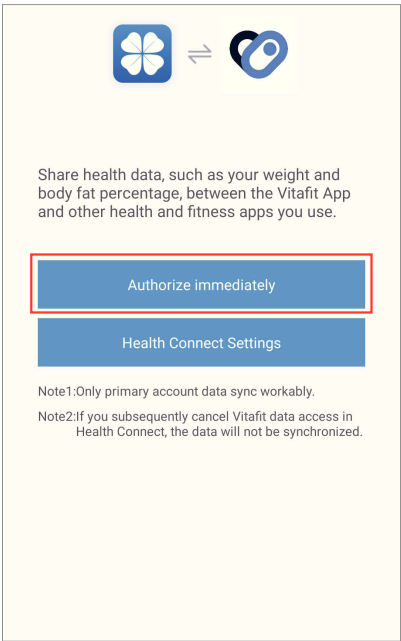

9. View measurement data in Google Fit after sync.

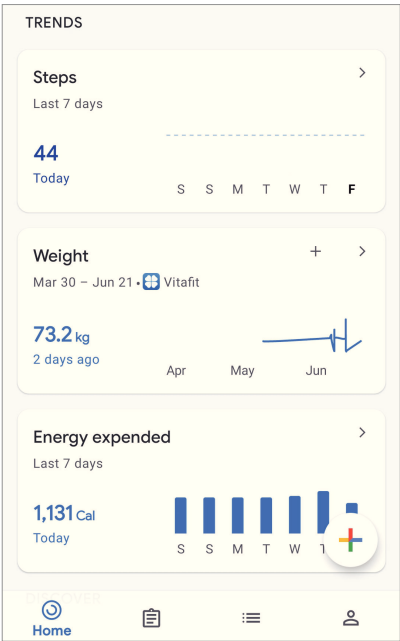

Notes:

1. For body fat measurement, make sure your feet bare and dry when measuring.
2. For excellent experience, recommend to reassemble batteries for the first beginning use.
3. UI will be different depended on scale type: Smart Weight, Smart Body Fat, VT Grade Floor Scale.

CONTENTS

01-06 English

07-12 Deutsch

13-18 Français

19-24 Español

25-30 Italiano

Safety Cautions

1. This scale is not intended for commercial use.
2. This scale should not be used to diagnose or treat any medical condition. Consult your primary care physician before making changes to your diet, exercise plans, or physical activities.
3. If you are pregnant, have a pacemaker, or other internal device, you should not heavily rely on the data provided by this product. Please note that all data should be used as a reference and is not substitute for medical advice from a licensed healthcare professional.
4. Check the scale before use. Do not use the scale if it is damaged in any way.
5. The platform is slippery when wet. Make sure the platform and your feet are dry before stepping on for measurement.

Use and Care

1. To weigh safely and accurately, stand with your feet shoulder-width apart while on the scale.
2. DO NOT strike or drop the scale.
3. DO NOT stand on the edge of the scale or jump on it.
4. DO NOT stand on the scale platform for a prolonged period of time or when not in operation.
5. DO NOT store or operate near heat sources or in humid environments.
6. Handle the scale with care, as it is a precision instrument.
7. Clean the scale with a clean and dry cloth. Avoid getting water or any other liquid inside the scale.
8. DO NOT use chemical or abrasive cleaning agents to clean the scale
9. Keep the scale in a cool and dry environment.
10. Always keep the scale in a horizontal position, do not store vertically.
11. If the scale fails to turn on, check whether the batteries are installed correctly. Replace the batteries if necessary.
12. If an error appears on the display or if the scale does not turn off automatically, remove the batteries for over 30 seconds, then reinstall them.
13. DO NOT overload the scale (maximum weight capacity:300kg/660lb/47st).
14. DO NOT mix new and used batteries or different battery brands to power the scale.

Product Introduction

Thank you for choosing VitaFit Smart body weight scale. This scale is your personal health assistant that could track and record your weight via VitaFit app. We sincerely hope that you enjoy our VitaFit scale.

Note: All operations are controlled by the remote controller button.


Specifications

Weight Capacity:	Max 300kg/660lb/47st
Weight unit:	kg/lb/st
Graduation:	0.05kg/0.1lb/0.1st
Dimensions of Scale Base:	381x330x38mm
Dimensions of Remote Controller:	180x105x45.3mm
Batteries:	4 X AA for Scale Base, 3 X AA for Remote Controller(Pre-Installed)


Operation

Weighing without VitaFit App

1. Remove the isolator tab from the battery compartment.
2. Place the scale base on a hard and flat floor. Carpeted or uneven floors can negatively affect the accuracy. For Remote Controller, you can set it up on the table or keep it handheld.
3. Press the on/off button of remote controller to turn on the scale
When the "0.0" displayed on the screen and the "●" icon is flashing, it indicates that the connection between scale base and remote controller is successful. Now it's ready for weighing.

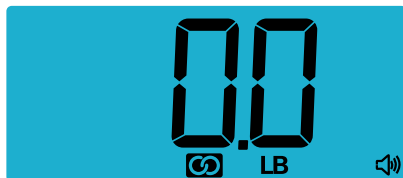

4. Step on the scale to start weighing, When “✔” icon keep on scale base flashing, it indicates that weighing process, please stand still unit get stable weighing result.

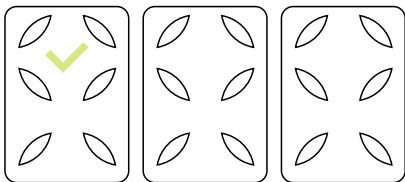

5. After the weight result is stable, the LCD screen will display the weight , and the buzzer will emit a BI sound when the sound is turned on. At this time, you can press the HOLD button to keep the current result. If need to weigh again, please press the HOLD button again to clear the result, and then step on the scale to take measurement.



Weighing with VitaFit App


The VitaFit Bluetooth Scale can connect with your mobile phone or Ipad via Bluetooth . The Body Weight and BMI result data will be showed at the App named “VitaFit”, which is easy to download from IOS App Store or Google Play.

Download App and Register ID

Scan the QR code or Download “ VitaFit” via IOS App Store or Google Play.





Register Your own ID, Then add personal data and confirm OK.

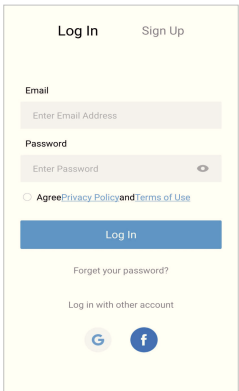

Pairing the VitaFit Scale with your smartphone

Turn on Bluetooth® on your smartphone.
NOTE: For Android users, Location and Bluetooth® must be turned on for Android 6.0 and later.
OR ELSE THE APP DON'T GET RESULTS ACCORDING TO ANDROID REGULATIONS.

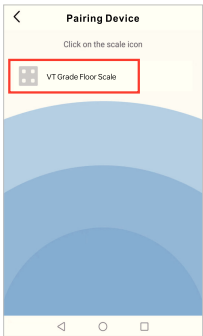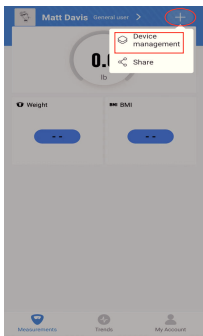

1. Remove the insulator strip from battery compartment
2. Place the scale on a hard flat floor surface, Turn on the remote controller.
3. Open VitaFit app, and tap “+” on the top right corner of the “**Device Management**” page to pair .
4. When the app appears " Please Step on", please step on the scale base with one foot to activate it, then the device name will be appeared on the app.