Enjoy!

Please visit the PicoCreate website for more products operation and maintenance tutorials.

Note: The UI interface is only for reference, and shall be subject to the actual UI.

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
Quick Start Guide

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
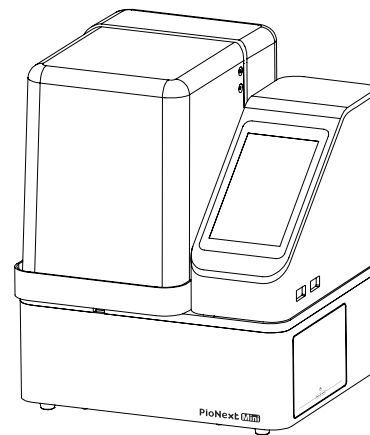
Usage scenarios


This product is designed for indoor use only

Please refer to this guide to complete the installation and printing of the printer.
*Safety Notice: Do not power on the printer before installation is complete.


Component Introduction

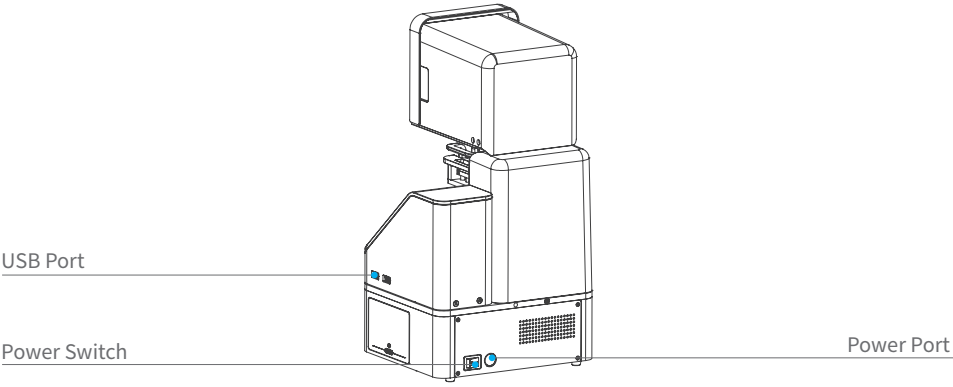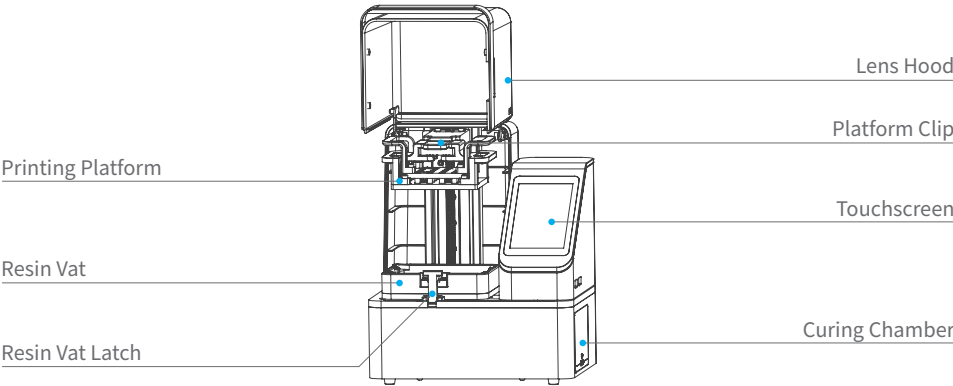List of Accessories

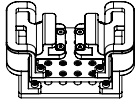Printing Platform
×1

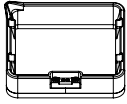Resin Vat
×2


Resin Vat Dust Cover
×2


Wash Mini
×1


U Disk
×1

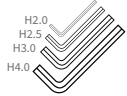
Hex Keys Set
×1

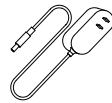
PioNext Mini
Power Adaptor ×1

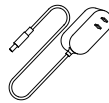
Wash Mini
Power Adaptor ×1

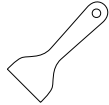Plastic Blade
×1

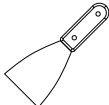Metal Blade
×1

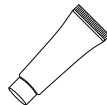Lubricant
×1

Paintbrush
×1

Disposable Latex Gloves
×1


Filters
×4

Qualified Certificate
×1


Quick Start
Guide ×1


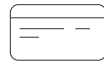
After-Sales
Service Card ×1

Quick Start Operation Schematic


① Unpack the Wash Mini, power adapter, and other accessories from the box. Then, remove the top foam.

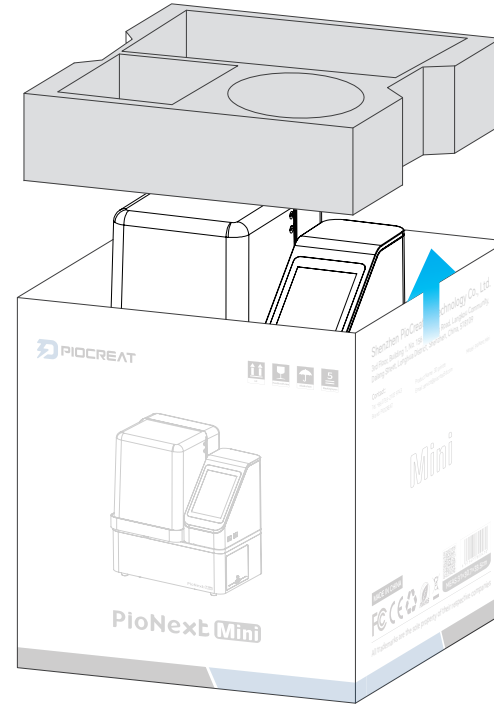② Connect power to PioNext Mini and switch on to start.

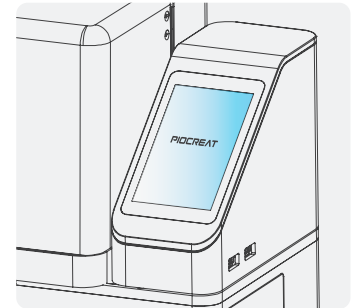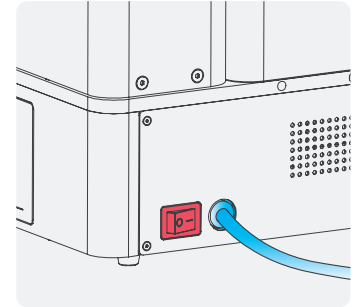First Printing

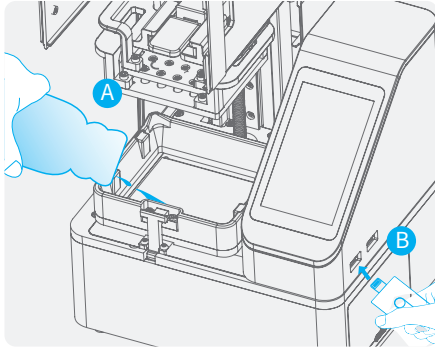
- A. Fill the resin vat.
B. Insert the included USB drive.

1


2

3


4

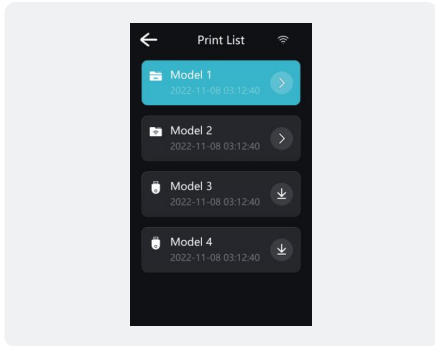Select the pre-sliced file stored in the device.

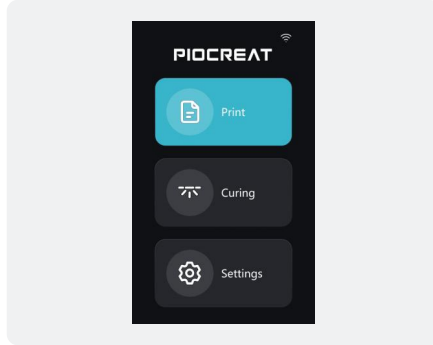
Tap "Print" on the home page.

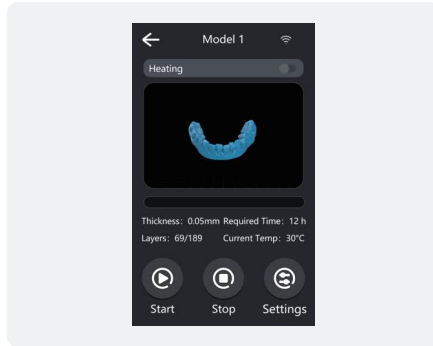Click "Start" to start the first print.

First Curing

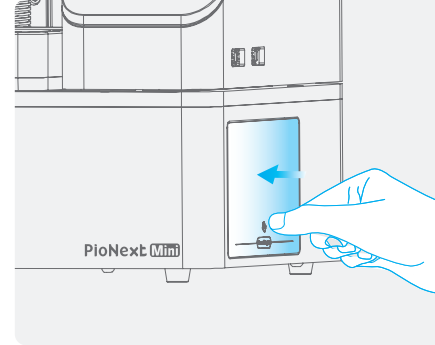
Press it inward to eject the curing chamber, then remove it.

1


2

3


4

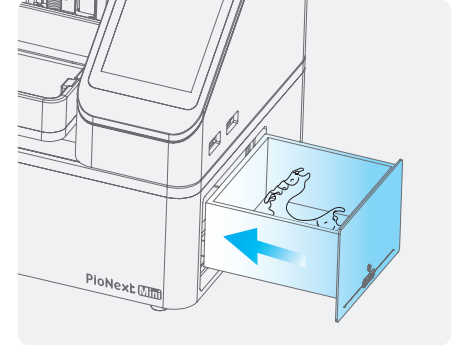Place the model inside, and reinstall the curing chamber into the machine.

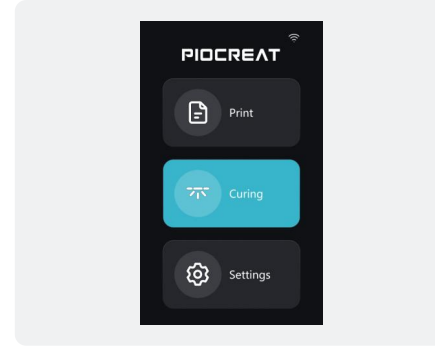
On the main interface, tap "Curing".

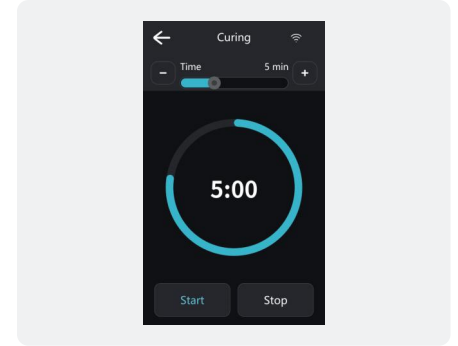
Set the curing time, tap "Start" to begin the first curing process.

Ultrasonic Cleaning

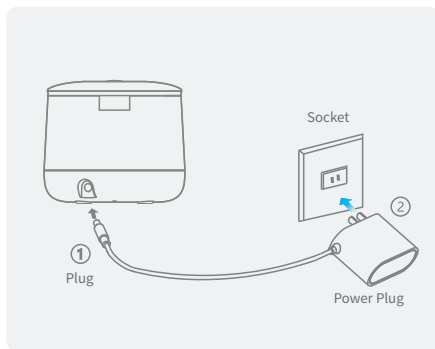With the power cord connected, the device is in standby with its digital display reading **888**.

1

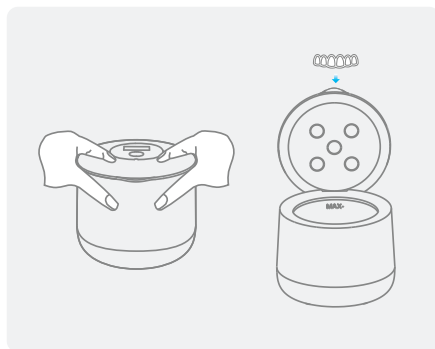
Place the machine on a level surface and open the lid. Then, gently place the dental model into the cleaning tank.

2

Fill the tank with clean water or alcohol, being sure not to exceed the "MAX" water level line.


3

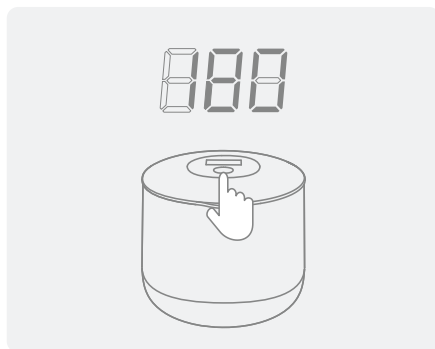
Press the button to initiate the ultrasonic function, which will operate for 3 minutes.

4


Remove Model & Replace Resin Vat

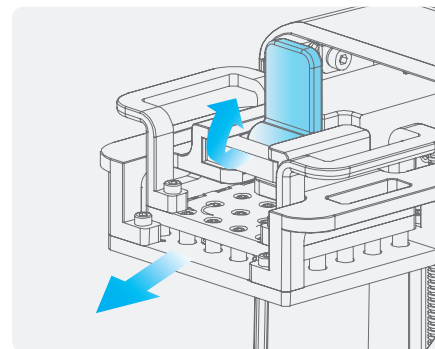
Lift the platform latch up to release and remove the printing platform.

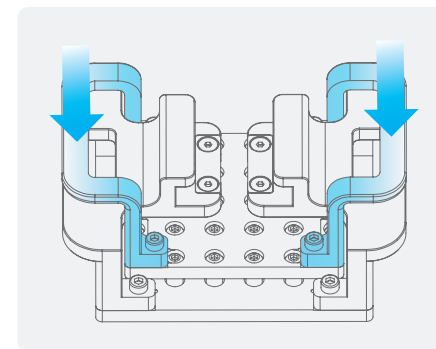Press down on the printing platform with both hands to detach the model from it.

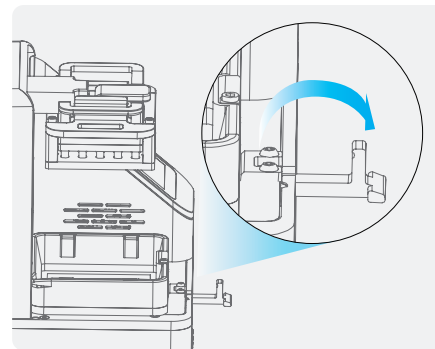Pull the resin vat latch outward to release and remove the vat.

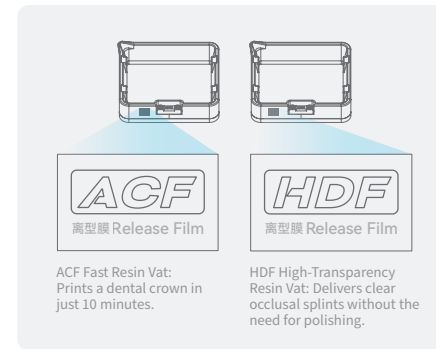
ACF Fast Resin Vat:
Prints a dental crown in
just 10 minutes.

HDF High-Transparency
Resin Vat: Delivers clear
occlusal splints without the
need for polishing.

When installing the resin vat,
note the type of pre-installed release liner.

Connect Device to Network & Install Software

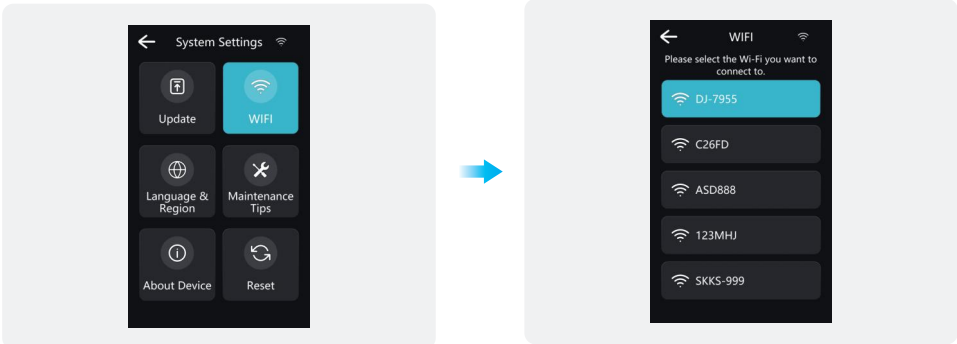In the main interface of the device, go to “Settings” → “System Settings” → “Wi-Fi”, select the Wi-Fi you want to connect to, and network the device.

Piocreat BOX

Download the installation package from www.piocreat3d.com
Read the installation package from a random USB flash drive

Tips: This product supports CHITUBOX Dental. Contact customer service to receive a free 3-month CHITUBOX Dental slicing software license.

Technical Specifications

PioNext Mini		
Body Parameters	Printing Size	80mm×70mm×80mm
	Machine Size	200mm×170mm×241mm
	Package Size	310mm×307mm×340mm
	N.W.	4.86kg
	G.W.	6.8kg
Printing Parameters	Printing Technology	LCD
	Light Source Wavelength	385nm、405nm
	Resolution	9K (8520*4320px)
	Pixel Size	18μm
	Layer Thickness	0.01-0.1mm
	Printing Speed	Max.60mm/h
	Touch Screen	3.97-inch Capacitive Color Touch Screen
System Parameters	Slicing Software	Piocreat BOX
	Data Transmission	USB/WIFI
	Input Voltage	100-240V~, 50/60Hz
Supported Consumables	Power	120W
	Ceramic Crown Resin / Surgical Guide Resin / Crown & Bridge Resin / Temporary Restoration Resin / Denture Base Resin / Removable Die Model Resin / Ortho Model Resin / Ortho Model 2.0 Resin and other official authorized resin	

IC ID Requirement

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference;
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

The device meets the exemption from the routine evaluation limits in section 2.5 of RSS 102 and compliance with RSS-102 RF exposure, users can obtain Canadian information on RF exposure and compliance.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. This equipment should be installed and operated with a minimum distance of 20 centimeters between the radiator and your body.

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FCC Requirement

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference;
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator & your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

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