


Re:Balans®

Quick Start

1. Re:Balans preparation

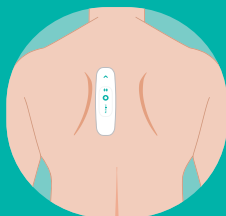


- Remove device from pouch. Trim hair, if needed.
- Clean the skin and let it dry.
- Place the device 2-4 cm / 0.8 - 1.6" left or right of spine.

2. Apply the device



- Remove release liner #1

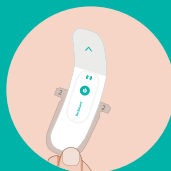

- Apply the first half with arrow up




- Remove release liner #2



- Apply device firmly

3. Activate the device

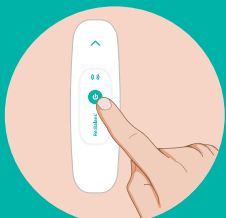

Press and hold the  button until green light. If orange light consult the user manual.

Safety Information

Intended Purpose

The Re:Balans is intended for use under the direction of a physician, for the non-invasive monitoring of patients with fluid management-related health conditions. The device measures thoracic bioimpedance in patients to assist the physician in fluid management assessment. The device is indicated for adult patients with fluid overload such as end-stage renal disease, and patients at risk of dehydration. The device does not generate any real-time alarms for consideration by the user at home or by a healthcare professional. Data from the device should be used in the context of all clinical data to make determinations of a patient's fluid status.

Contraindications

The device should not be used by persons with:

- Known allergies to hydrogel or acrylic tapes
- Non-intact skin
- Implantable pulse generators
- Scheduled MRI or diathermy
- Pregnancy



Storage Conditions:



Rx only

IP54

FC

FCC-ID: 2A6Z8-PCT0010

Model Number: PCT0010A



Quantity: 1



A Mode Sensors product



Mode Sensors AS
Sluppenvegen 6
7037, Trondheim, Norway
support@modesensors.com
www.modesensors.com