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


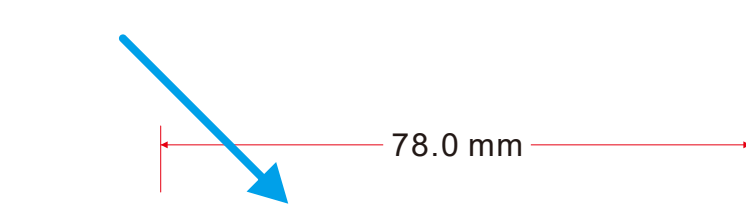
中文/English/日本語/Español/Français/Italiano

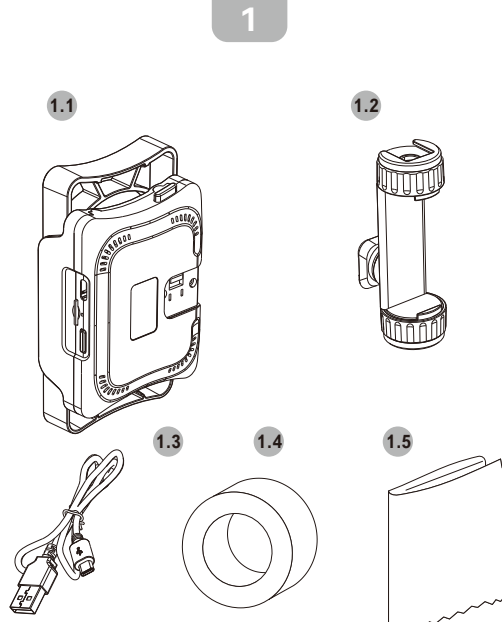
2

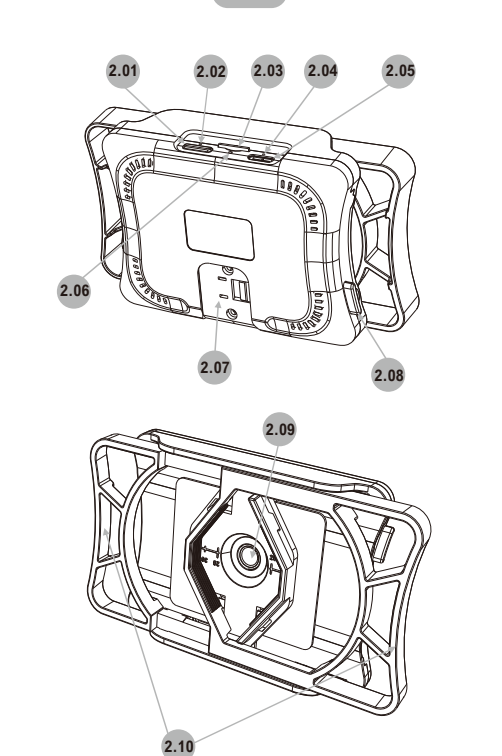3

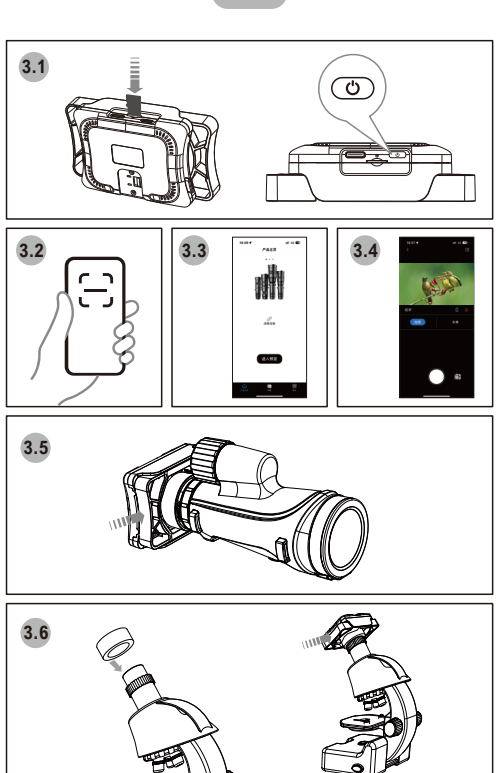4

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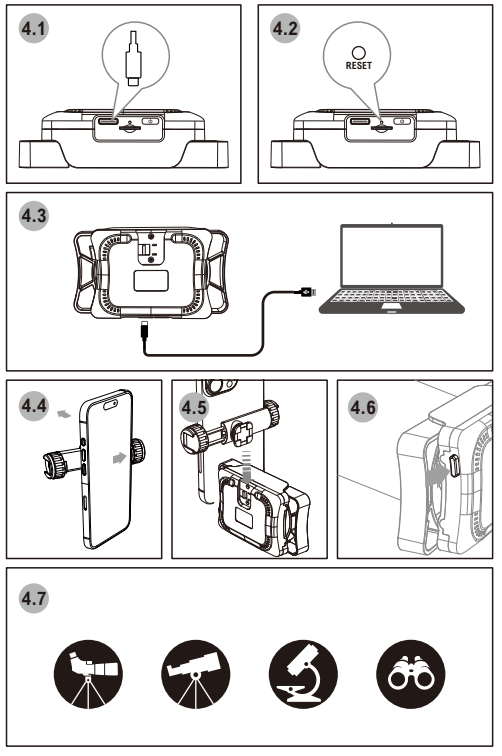English

Note: This section is an English instruction content using the illustration in the first column on the left as an example.

1

Product Accessories List

- 1.1 Electronic Eyepiece
- 1.2 Smartphone Stretch Clip
- 1.3 USB Type-C Charging Cable
- 1.4 EVA Ring (for connecting and fixing small-sized eyepieces)
- 1.5 Cleaning Cloth

2

Introduction to the Electronic Eyepiece

- 2.01 Power and WiFi Switch
- 2.02 WiFi Indicator Light
- 2.03 TF Card Slot
- 2.04 Charging/Working Status Indicator Light
- 2.05 USB Type C Interface
- 2.06 Reset Button
- 2.07 Smartphone Clip Slot
- 2.08 Smartphone Clip Release Switch
- 2.09 Lens
- 2.10 Clip Left and Right Pressure Switches

3

Eyepiece Usage Instructions

- 3.1 Insert the TF card and press the button for 3 seconds to turn on the device. (Note: The memory card needs to be formatted for the first use, and the APP will remind you with a pop-up window.) 3.2 Scan the QR code in the manual with your smartphone to download the APVision software. 3.3 After installation, open the APP and click "Connect Device" to go to the device search/WiFi hotspot interface. Please find the WiFi name: Apexel-OC001-xxxxx, default password: 12345678. Note: Double-click the button to turn WiFi on/off. 3.4 Electronic Eyepiece Indicator Light Description Blue light steady on - Device is powered on Blue light flashing - Device is low on battery Yellow light flashing - WiFi is waiting to connect Yellow light steady on - WiFi is connected 3.5 When the WiFi connection is successful, it will jump to the real-time preview screen of the APP, and operations such as shooting, recording, viewing, and downloading are all completed within the APP. 3.6 Connecting to a Telescope First, place the EVA ring over the eyepiece clip without releasing and attach it to the telescope's eyepiece end to start shooting and using. 3.7 Connecting to a Microscope First, place the EVA ring over the microscope's eyepiece end, then press both sides of the eyepiece clip without releasing and attach it to the microscope's eyepiece end to start shooting and using. (Note: Using the EVA ring increases the outer diameter size to make the clip more secure. Some microscope eyepiece ends with a diameter larger than 025MM do not need to use the EVA ring.)

4

Other Instructions

- 4.1 Charging Connect the power adapter, with Type-C for power supply. The red light will stay on during charging, and it takes about 2-3 hours to fully charge. The green light will stay on when fully charged. The battery life is approximately 150 minutes. 4.2 Reset Use a card removal pin to briefly press the reset button. The device will shut down, and all indicator lights will go out. 4.3 Data Transfer Insert the USB Type-C end of the data cable into the port, and the other end into the computer's USB port for data transfer and reading. At this point, briefly press the electronic eyepiece switch button to cycle through the selection modes: "Real Time" > "Still" > "Movie". 4.4 *Note: If you need to mount your smartphone on the clip for observation, please first place your smartphone into the stretch clip. 4.5 Then, insert the smartphone along with the stretch clip into the electronic eyepiece end. 4.6 Press the unlock button and slide up to remove the smartphone and stretch clip. 4.7 The electronic eyepiece is a cutting-edge device that seamlessly links your smartphone to a range of optical instruments, from microscopes to telescopes. It enables direct viewing and recording of observations on your phone, simplifying the capture of stunning images and videos. Plus, it facilitates easy sharing of these incredible moments, enhancing the user experience and making observation more engaging and fun. 4.8 If there is no operation within 3 minutes in the powered-on state, the electronic eyepiece will automatically shut down. Please turn it on again if you need to use it. 4.9 Thank you very much for choosing our smartphone electronic eyepiece. Please read the manual carefully before using this lens, and we suggest you keep it well for future reference.

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FCC Caution:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.