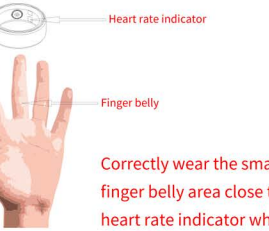
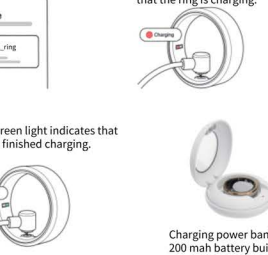
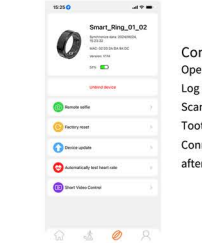
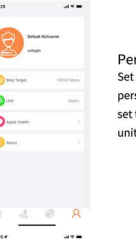
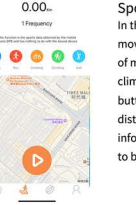
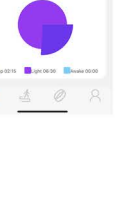
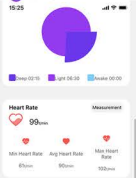
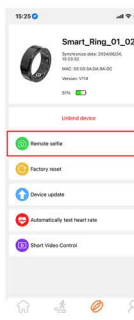
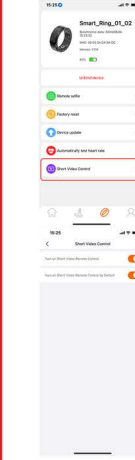
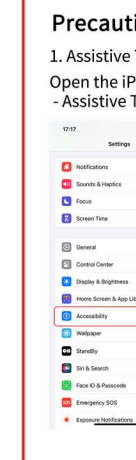
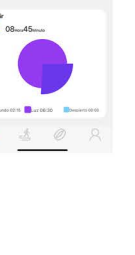
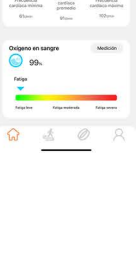
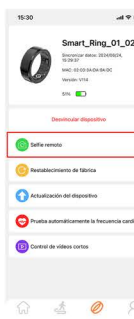


Instruction manual Smart ring JRing	Charge activation When using it for the first time, please charge and activate the device. If factory settings are restored on the app, please charge and reactivate. 1.Ring charging: Use the matching charging cable to connect the ring or charging power bank box 2. After the power is turned on, the ring indicator light is always red. 3.Indicate light status: Charging Red light is always on. Charging completed Green light is always on 4.Charging time: 60 minutes 5.Battery Volume check: Check the battery Volume at the device connection on the main page of the APP 6.Charging compartment: If the box charger is out of power, you need to plug in the Type-C cable to charge. The red light turns off to indicate that the box charger has completed charging.  The optional package comes with a 200mah charging box	Please note that the Ring Service is not intended to diagnose, treat or prevent any disease or medical. The information and data in the condition service are for reference only and are not a substitute for health professionals or physician services. If you have any questions about a medical condition, it is recommended that you consult a professional physician.  Heart rate indicator Finger belly Correctly wear the smart ring schematic, finger belly area close to heart rate indicator when wearing	1. Charge the ring for 60 minutes before using it for the first time. 2. Search on Google Play or App Store JRing or scan the QR code below to download the APP  JRing 3.Select your ring in the app and pair it. 4.A flashing red light indicates that the ring is charging. 5.A flashing green light indicates that the ring has finished charging.  Charging power bank box 200 mah battery built-in  JRing Connecting ring Open JRing APP and register Log in to the APP account and find the device. Scan the device to find the corresponding blue Tooth name Smart_Ring, click Connect, it will be displayed Ring information after the connection is successful. 	Personal settings Set up related information in the personal information interface, you can set the number of step targets, units and other information.  Sports In the mobile app interface, find the movement interface, you can choose a variety of motion mode,including: walking,running, riding, climbing, Golf, etc, click on the motion button to start recording the distance,pace,duration, and heart rate information of the motion.(This function needs to be hit GPS on mobile phone) 	Motion data On the homepage of the mobile APP, the number of steps,distance,movement duration, and calorie data recorded on the day can be displayed.You can also click to view the sports data of the nearly week or one month. The number of target steps can be set on the APP personal information interface.  Sleep Smart ring can show the total sleep last night and the length of light sleep and deep sleep. You can also click to view the sleep data for nearly a week or a month. 	Heart rate On the homepage of the mobile app,select the heart rate,click to start the in dicat light is turned on, and the measurement is completed for 30 to 60 seconds.  Blood oxygen On the home page of the mobile APP, select blood oxygen, click to start measurement, the indicator will light up green, start testing for 30-60 seconds to complete the measurement, and the mobile APP will display the current blood oxygen, blood oxygen and fatigue value. 	Open the remote selfie interface of the app and shake to take a photo (Remote photography can only be used in JRing App)  	Short video remote control 1. After the mobile phone is connected to the device and the device Bluetooth, enter the device interface and click the "Short Video Control" option. 2. After entering the short video control interface, turn on the two control switches on the right. (Turn on switch 1 to implement the short video control function, then turn on switch 2. Each time you connect, switch 1 is turned on by default.) 3. After turning on the switch, you can control the short video software (tiktok, Youtube, etc.) remotely (switch to the next video) and turn the page of the phone by shaking the ring 	Precautions for using Apple mobile phones 1. Assistive Touch Settings Open the iPhone's Settings - Accessibility - Touch - Assistive Touch Settings - Turn on 2. After entering the short video control interface, turn on the two control switches on the right. (Turn on switch 1 to implement the short video control function, then turn on switch 2. Each time you connect, switch 1 is turned on by default.) 3. After turning on the switch, you can control the short video software (tiktok, Youtube, etc.) remotely (switch to the next video) and turn the page of the phone by shaking the ring 	Common problem 1. Causes of abnormal heart rate and blood oxygen data measured by ring 1) It is generally caused by poor contact between the sensor of the ring and the human body when measuring. 2) When measuring, please note that the sensor is in full contact with the finger. 2. Reasons for abnormal sleep data 1) Sleep monitoring is to simulate people's self-perceived state of falling asleep and waking up time, neededWear it correctly. 2) Sleep too late or sleep before wearing, there may be errors. 3) Sleep data is not monitored during the day, and the default sleep monitoring is from 9:00 PM to 11:00 noon the next day. More frequently asked questions please see the App[My —User Manual]	Notice 1) Do not charge in humid and watery environments 2) Please clean the magnetic charging hole on the back of the ring with a clean flannel regularly to ensure that the magnetic charging hole of the ring and the magnetic charger fit together to ensure normal charging. This product is not equipped with a power adapter. To ensure the safety of your family and property, When charging, you can use the computer USB interface or choose a power adapter whose output does not exceed 5V1A. Please purchase power adapters through regular channels and avoid using low-quality or shoddy power adapters to avoid bursting or catching fire. Notice 1. The measurement results of this product are for reference only. Do not be used as a medical use and basis for the guidance of the doctor. This measurement result can not use for self-diagnosis and treatment. 2. The waterproof level of this product is 5ATM waterproof and deep waterproof. 3. The right to modify the content of this manual without any notifications will not be notified separately. Some functions are different in the corresponding software version, which is normal.
MANUAL DE INSTRUCCIONES ANILLO INTELIGENTE JRing	Activación de carga Cuando lo utilice por primera vez, cargue y active el dispositivo. Si se restablecen los ajustes de fábrica en la aplicación, cárguela y vuelva a activarla. 1. Carga del anillo: utiliza el cable de carga correspondiente para conectar el anillo. 2. Después de encender la alimentación, la luz indicadora del anillo siempre está roja. 3.Indicar el estado de la luz: Cargando La luz roja siempre está encendida. Carga completada La luz verde siempre está encendida 4.Tiempo de carga: 60 minutos 5.Verificación del volumen de la batería: verifique el volumen de la batería en el dispositivo conexión en la página principal de la APP 6. Compartimento de carga: si el cargador de la caja no tiene energía, debe enchufarlo en el cable tipo C para cargar. La luz roja se apaga para indicar que el El cargador de caja ha completado la carga.  El paquete opcional viene con una caja de carga de 200 mah.	Tenga en cuenta que el Servicio Ring no está destinado a diagnosticar,tratar o prevenir cualquier enfermedad o enfermedad médica.La información y datos en el estado del servicio son para Solo como referencia y no sustituyen a los productos de salud.profesionales o servicios médicos. Si tiene alguna pregunta sobre una condición médica, se recomienda consultar a un médico profesional.  Indicador de frecuencia cardíaca Ventre dedo Llevar correctamente el esquema del anillo inteligente,área del vientre del dedo cerca de indicador de frecuencia cardíaca cuando se usa	1. Carga el anillo durante 60 minutos. antes de usarlo por primera vez. 2. Busque en Google Play o App Store JRing o escanea el código QR a continuación para descargar la APLICACIÓN  JRing 3.Selecciona tu anillo en la aplicación y emparejela. 4. Una luz roja intermitente indica que el anillo se está cargando. 5. Una luz verde intermitente indica que el anillo ha terminado de cargarse.  Caja de carga del banco de energía Batería de 200 mah incorporada  JRing Anillo de conexión Abra la APLICACIÓN JRing y regístrese Inicie sesión en la cuenta de la aplicación y busque el dispositivo.Escanee el dispositivo para encontrar el azul correspondiente Nombre del diente Smart_Ring, haga clic Conéctese, se mostrará la información del timbre. después de que la conexión sea exitosa. 	Configuraciones personales Configurar información relacionada en el interfaz de información personal, puede establecer el número de objetivos de pasos, unidades y otra información.  Deportes En la interfaz de la aplicación móvil, busque el Interfaz de movimiento, puedes elegir una variedad del modo de movimiento, que incluye: caminar, correr, montar a caballo,escalar, Golf, etc, haga clic en el movimiento botón para comenzar a grabar el distancia, ritmo, duración y frecuencia cardíaca información del movimiento. (Esta función necesita ser golpeado GPS en el teléfono móvil) 	Datos de movimiento En la página de inicio de la aplicación móvil, el número de pasos, distancia, movimiento duración y datos de calorías registrados Se puede mostrar el día. Puedes Haz clic también para ver los datos deportivos de la Casi una semana o un mes. El número El número de pasos objetivo se puede configurar en la aplicación.interfaz de información personal.  Dormir El anillo inteligente puede mostrar el sueño total de la noche anterior y la duración del sueño ligero y del sueño profundo.También puedes hacer clic para ver los datos del sueño. durante casi una semana o un mes. 	Ritmo cardíaco En la página de inicio de la aplicación móvil, seleccione el ritmo cardíaco, haga clic para iniciar la indicación La luz se enciende y la medición es completado durante 30 a 60 segundos.  Oxígeno en sangre En la página de inicio de la aplicación móvil, seleccione oxígeno en sangre, haga clic para iniciar la medición,el indicador se iluminará en verde, comience a probar durante 30-60 segundos para completar la medición,y la aplicación móvil mostrará la información actualoxígeno en sangre, oxígeno en sangre y valor de fatiga. 	Abra la interfaz de selfie remota de la aplicación y agitar para tomar una foto (La fotografía remota solo se puede utilizar en la aplicación JRing)  	Vídeo corto de control remoto 1. Después de que el teléfono móvil esté conectado al dispositivo y al Bluetooth del dispositivo, ingrese a la interfaz del dispositivo y haga clic en "Control de video corto" para seleccionar artículo. 2. Después de ingresar a la interfaz de control de video breve, encienda los dos interruptores de control a la derecha. (Encienda el interruptor 1 para realizar la función de control de video corto, luego encienda el interruptor 2 y encienda el interruptor 1 de forma predeterminada para cada conexión) 3. Después de encender el interruptor, agitando el anillo, puedes controlar el software de video corto (tiktok, Youtube, etc.) a través del aire (cambiar el video a la siguiente canción) y controlar la página al encender el teléfono móvil. 	Precauciones al utilizar teléfonos móviles Apple 1. Configuración táctil de asistencia Abre los ajustes de tu iPhone - Accesibilidad - Toque - Ajustes de Assistive Touch - Encender 	Problema comun 1. Causas de frecuencia cardíaca anormal y datos de oxígeno en sangre medidos por anillo 1) Generalmente es causado por un mal contacto entre el sensor del anillo y el cuerpo humano al medir. 2) Al medir, tenga en cuenta que el sensor está en contacto total con el dedo. 2. Razones de los datos de sueño anormales 1) La monitorización del sueño sirve para simular el estado de caída autopercibido por las personas.Hora de dormir y despertarse, necesaria Úselo correctamente. 2) Duerma demasiado tarde o duerma antes de usarlo, puede haber errores. 3) Los datos del sueño no se controlan durante el día y los datos predeterminados El seguimiento del sueño es desde las 21:00 hasta las 11:00 del mediodía del día siguiente Para preguntas más frecuentes, consulte la aplicación [Mi —Manual de usuario]	Aviso 1) No cargar en ambientes húmedos y acuosos 2) Limpie el orificio de carga magnético en la parte posterior del anillo con un Limpe la franja con regularidad para asegurarse de que el orificio de carga magnético del anillo y el cargador magnético encajan para garantizar la carga normal. Este producto no está equipado con una fuente de alimentación. adaptador. Para garantizar la seguridad de su familia y su propiedad, al cargar, Puede utilizar la interfaz USB del ordenador o elegir un adaptador de corriente cuya salida no supere los 5V1A. Por favor compre adaptadores de corriente a través de canales regulares y evite el uso de baja calidad o Adaptadores de corriente de mala calidad para evitar estallar o incendiarse. Aviso 1. Los resultados de las mediciones de este producto son solo de referencia. No debe utilizarse como uso médico ni como base para la orientación del doctor. Este resultado de medición no se puede utilizar para el autodiagnóstico y el tratamiento. 2.El nivel de impermeabilidad de este producto es 5ATM Resistente al agua y al agua profunda. 3.El derecho a modificar el contenido de este manual sin ningún Las notificaciones no se notificarán por separado. Alguno Las funciones son diferentes en el versión de software correspondiente, lo cual es normal.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction