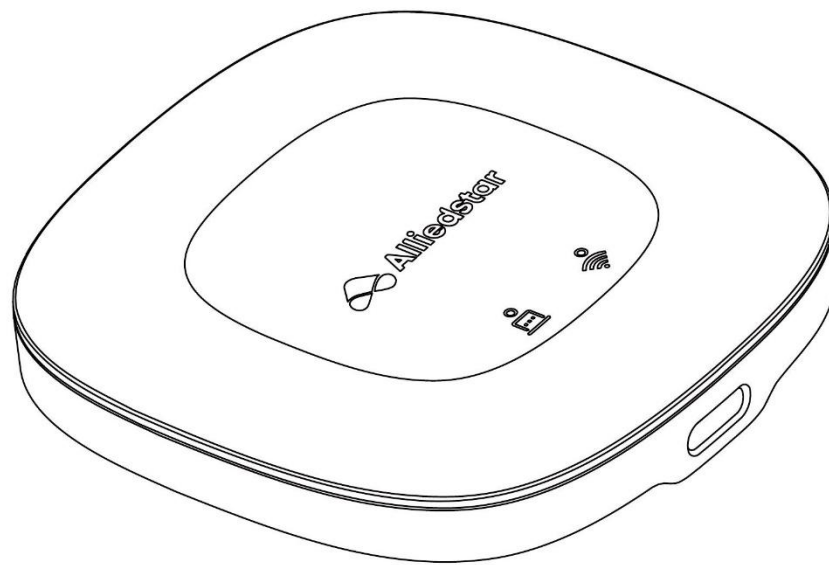
Intraoral Scanner Wireless Adapter

Instruction Manual

All rights reserved.

Copyright © 2026 Alliedstar.

The content of this instruction is subject to change.

Alliedstar and its subsidiaries shall not be liable for any indirect damages arising from the preparation, presentation, or use of the information contained in this instruction, including any errors. No part of this instruction can be reproduced without the permission of Alliedstar.

Trademarks

All trademarks and registered trademarks are the property of their respective owners.



Lifetime: 6 years

Contents

EN/English	4
Package Contents	4
Setup Guide.....	5
LED Indicator Status	6
Troubleshooting.....	7
Environmental Requirements	7
Regulatory Information	7
Labels and Symbols	13
ES/ Español.....	15
Contenido del paquete	15
Guía de instalación	16
Estado de los indicadores.....	17
Solución de problemas.....	17
Requisitos ambientales	17
Información normativa.....	18
Etiquetas y símbolos	22
PT/Português	24
Conteúdo da Embalagem	24
Guia de Instalação	25
Status do Indicador	26
Solução de Problemas.....	26
Requisitos Ambientais	26
Informações Regulatórias	27
Rótulos e Símbolos	30

EN/English


Package Contents

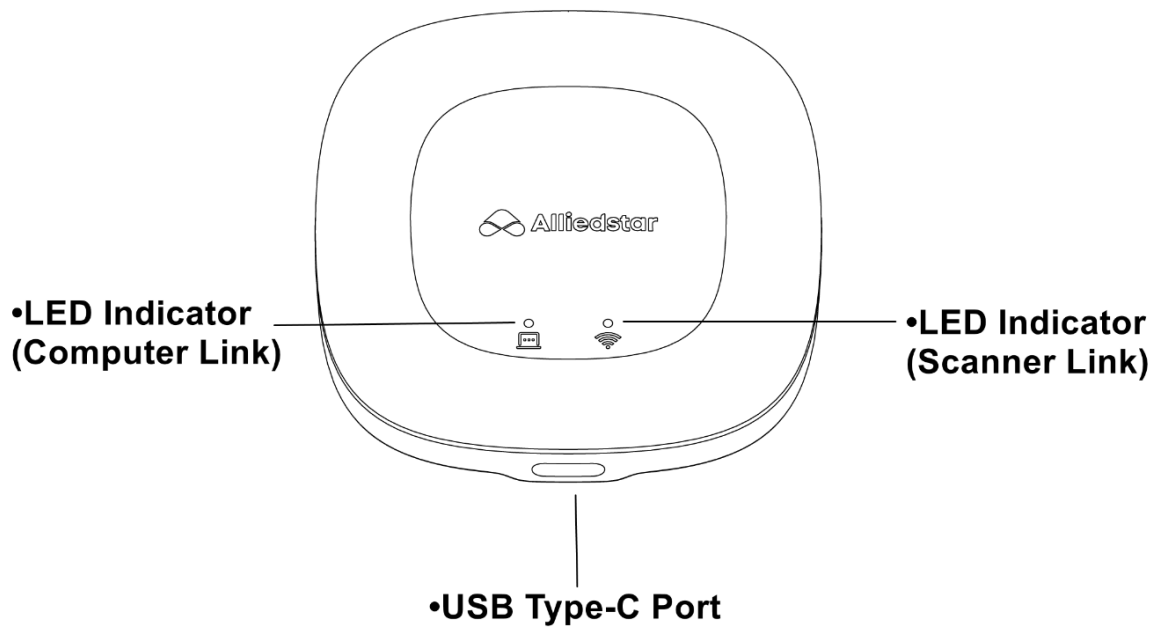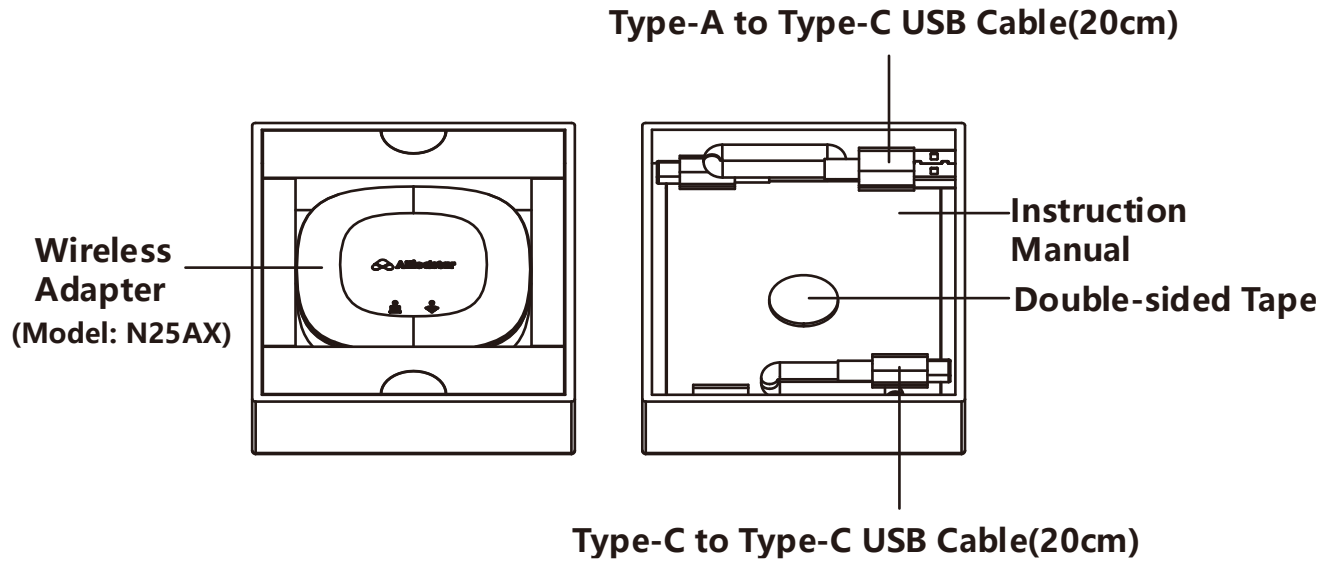Setup Guide

Part A: Physical Setup


1. Place the wireless adapter on a stable, flat surface near your computer.
2. (Optional) For vertical mounting, stick the provided double-sided tape to the back of the adapter and affix it to a clean, vertical surface.
3. Connect the included USB cable from the adapter to an available USB port on your computer.

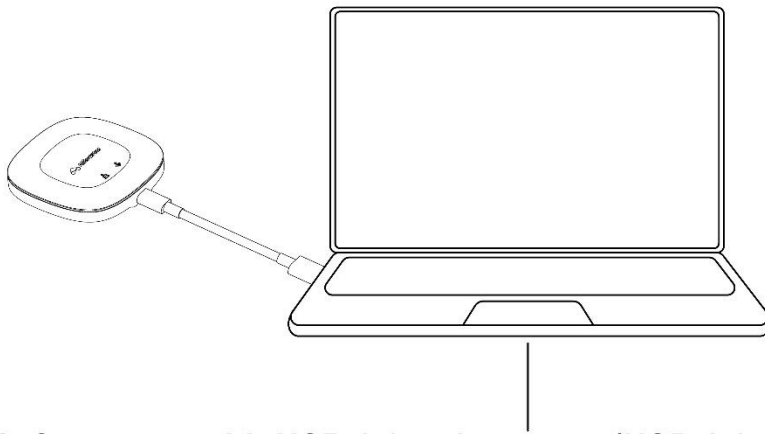
1. **Computer with USB 2.0 or later port (USB 3.0, 3.1, etc. supported)**
2. **Windows 10 (build 18362 or later) or Windows 11, 64-bit**

Part B: Software Setup

4. Ensure ScanPro software version 1.0.40.0 or later is installed.
5. Power on your computer and your intraoral scanner.
6. Open the ScanPro software.
7. From the scanner list, click **Connect** next to your scanner model.

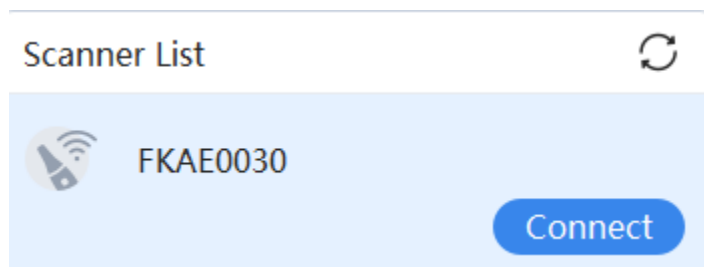8. Once connected, the button will display **Disconnect** and the battery status will appear. You can now begin scanning.

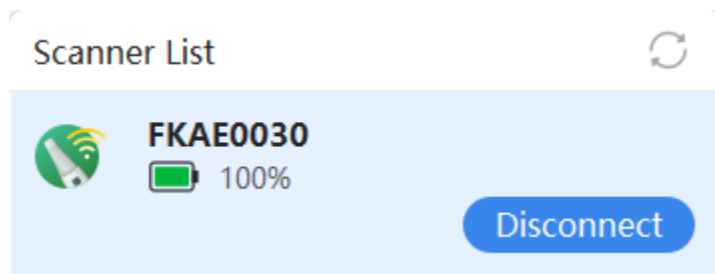Part C: Wireless Communication & Network Security

- **Local Connectivity Statement:** Alliedstar Intraoral Scanner and its dedicated Wireless Adapter utilize standard Wi-Fi protocols for localized peer-to-peer communication. This connection is strictly dedicated to the transmission of control commands and captured oral image data between the hardware components.
- **Network Isolation:** The system operates within a closed-loop wireless environment. It does not require, nor is it capable of, connecting to the public internet (WAN). All data processing and 3D model reconstruction are performed locally on the computer. This architecture ensures high-speed data stability and maximum network security by remaining completely isolated from external networks.

Note: Users cannot manually configure this device to access the public internet, ensuring the integrity of the clinical workflow.

LED Indicator Status

Status			Meaning
Starting Up	Flashing Green	Flashing Green	Device is starting
Ready to Connect	Solid Green	Breathing Green	Ready to connect
Connected	Solid Green	Solid Green	Connected to scanner and PC
Updating Firmware	Flashing Orange	–	Firmware is being updated
Hardware Error	–	Flashing Orange	Hardware error – restart device

Troubleshooting

Problem	Solution
No connection	1. Try another USB port. 2. Restart your computer.
Cannot find scanner Wi-Fi	1. Make sure the scanner is powered on (status light breathing). 2. Keep the scanner within range. 3. Disconnect and reconnect the wireless adapter.
Connection unstable or slow	1. Move the adapter and scanner closer. 2. Remove large metal objects or interference sources between them. 3. Change the Wi-Fi channel on the scanner if possible.
Indicator not Solid Green	See " <i>Indicator Status</i> " section for details.

For technical support, visit <https://support.allied-star.com/>

Environmental Requirements

Condition	Temperature	Relative Humidity	Atmospheric Pressure
Operating	15°C ~ 30°C	10% RH ~ 65% RH	70 kPa ~ 106 kPa
Transportation	-10 °C ~ 50°C	10 % RH ~ 95% RH	60 kPa ~ 106 kPa
Storage	-10°C ~ 50°C	10 % RH ~ 85% RH	70 kPa ~ 106 kPa

Regulatory Information

FCC compliance information statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception,

which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/ TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1) This device may not cause harmful interference.
- 2) This device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user's authority to operate the equipment.

5.9 and 6GHz Wi-Fi Devices

FCC regulations restrict the operation of this device to indoor use only. This device cannot be used to provide connection between separate buildings or structures. The operation of this is prohibited on oil platforms, cars, trains, boats, and aircraft, except the device is permitted in large aircraft while flying above 10,000 feet. Operation of transmitters in the 5.915-7.125 GHz band is prohibited for control of or communication with unmanned aircraft systems.

FCC RF Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for

an uncontrolled environment. End users must follow the specific operating instructions for satisfying RF exposure compliance. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. This equipment has been SAR-evaluated for use in hand. SAR measurements are based on a 20 cm spacing from the body and that compliance is achieved at that distance.

CE Mark Warning

This is a class B product. In a domestic environment, this product may cause radio interference, in which case the user may be required to take adequate measures.

OPERATING FREQUENCY (the maximum transmitted power)

5150MHz-5250MHz(23dBm)

5735MHz-5835MHz(25mW)

5925MHz-6425MHz(23dBm)

EU declaration of conformity

Alliedstar hereby declares that the device is in compliance with the essential requirements and other relevant provisions of directives 2014/53/EU, 2011/65/EU and (EU)2015/863.

Alliedstar hereby declares that the device is in compliance with the essential requirements and other relevant provisions of the Radio Equipment Regulations 2017.

RF Exposure Information

This device meets the EU requirements (2014/53/EU Article 3.1a) on the limitation of exposure of the general public to electromagnetic fields by way of health protection.

This device has been tested and meets the ICNIRP exposure guidelines and the European Standard EN 62209-2. SAR is measured with this device at a separation of 20 cm to the body, while transmitting at the highest certified output power level in all frequency bands of this device. Carry this device at least 20 cm away from your body to ensure exposure levels remain at or below the as-tested levels.

National restrictions

Attention: This device may only be used indoors in all EU member states and EFTA countries.

	AT	BE	BG	CH	CY	CZ	DE	DK
	EE	EL	ES	FI	FR	HR	HU	IE
	IS	IT	LI	LT	LU	LV	MT	NL
	NO	PL	PT	RO	SE	SI	SK	UK(NI)

Attention: This device may only be used indoors in Great Britain.

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This product must be powered by a device whose output meets the requirements of ES1 and PS1.

ISED Warning

This device complies with Innovation, Science and Economic Development Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device. This Class B digital apparatus complies with Canadian ICES-003.

IC: 29544-N25AX

The device for operation in the band 5150-5250 MHz is only for indoor use to reduce the potential for harmful interference to co-channel mobile satellite systems;

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. This equipment should be installed and operated with a minimum distance of 20 centimeters between the radiator and your body.

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




Labels and Symbols

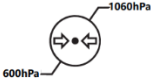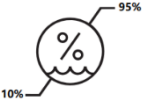	Manufacturer
	Date of Manufacture
	Model number
	Serial number
	Direct current
	Separate collection for waste electric and electronic equipment (WEEE) is required.
	Atmospheric pressure limitation
	Fragile, handle with care
	Humidity limitation
	Keep away from rain
	Stacking limit by number

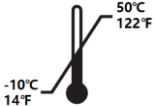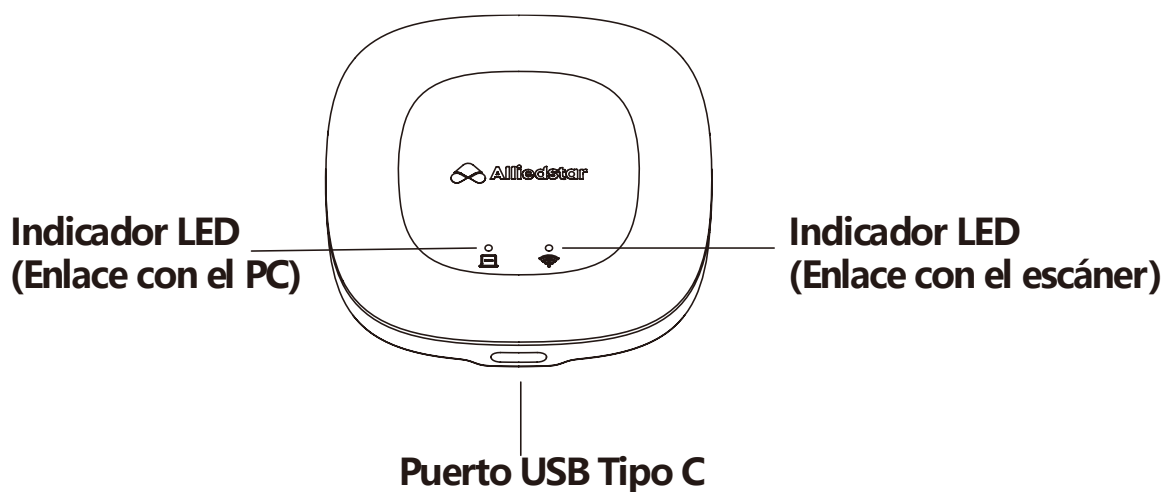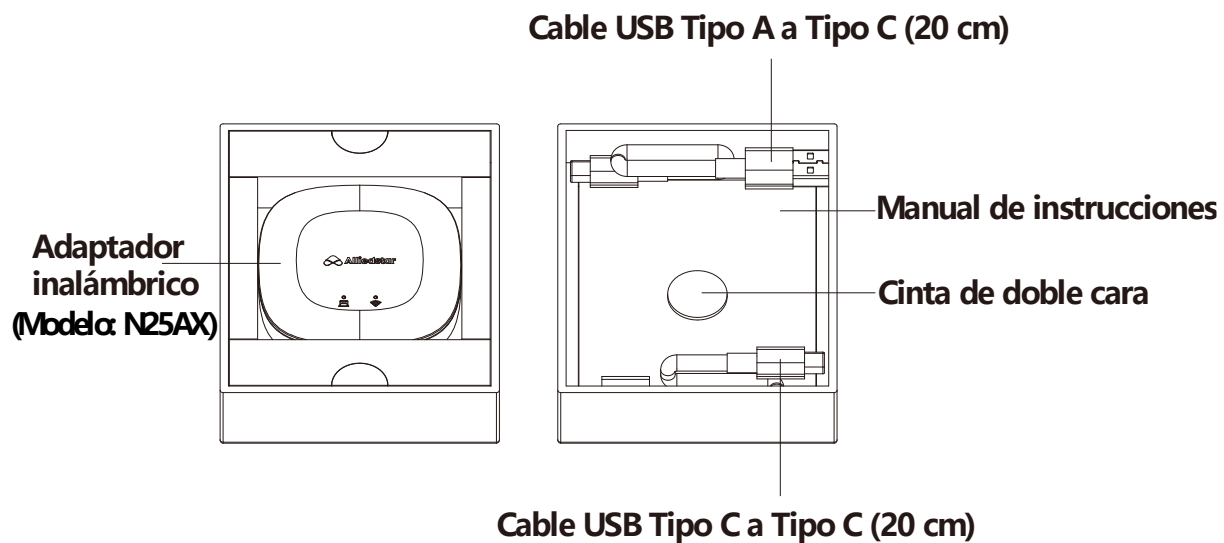	Temperature limits
	Consult instructions for use. Follow the link to the eIFU: ifu.allied-star.com .
	Catalogue number

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Alliedstar Medical Equipment Co., Ltd.