
Smart Wireless Station

User Manual

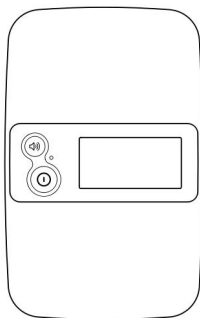Contact us: service@getbabytone.com

Website: www.getbabytone.com

Contents

User Manual.....	English 1-9
Benutzerhandbuch.....	Deutsch 10-21
Manuale d'uso.....	Italiano 22-32
Manual de usuario.....	Español 33-43
Manuel de l'utilisateur.....	Français 44-54

User Manual

Contents

1 Instructions for Safe Sleep	1
2 Warnings and Cautions	1
3 Overview	2
3.1 Station	2
3.2 Charging	2
4 Pairing the Station with the Main Unit	3
5 Monitoring with the Station	4
6 Stop Measuring	5
7 Cleaning	6
8 Troubleshooting	6
9 Specifications	7
10 Symbols	7
11 FCC Statement	8

1 Instructions for Safe Sleep

- Room-share instead of co-sleeping.
- Place your baby alone and on their back to sleep.
- Use a firm sleep surface. Keep soft objects out of crib.
- Use a crib that meets safety standards. Do not use cradle for sleep.
- Don't use Baby Sleep Monitor as an excuse for unsafe sleep.

2 Warnings and Cautions

- Do not use this device during MRI examination.
- Do not store the device in the following locations: locations in which the device is exposed to direct sunlight, high temperatures or levels of moisture, or heavy contamination; locations near to sources of water or fire; or locations that are subject to strong electromagnetic influences.
- Do not use the device in a combustible environment.
- Never submerge the device in water or other liquids.
- Do not clean the device with acetone or other volatile solutions.
- Do not drop this device or subject it to strong impact.
- The device is provided non-sterile.
- Do not place this device in pressure vessels or gas sterilization device.
- Do not dismantle the device, as this could cause damage or malfunctions or impede the operation of the device.
- Do not open the device cover without authorization. The cover should only be opened by a qualified service personnel.

3 Overview

3.1 Station

The station, coupled with the Babytone Sleep Sock (Baby Sleep Monitor), is used by parents to track the baby's real-time oxygen level and pulse rate. It provides reminders when they are out of the target range.

1. Display screen


2. Power button

- Press for 2 seconds to power on/off.

- Press to turn off the screen.

3. Volume button

- Press to adjust the volume.
- Press for 2 seconds to check the serial number of the station and the paired device.

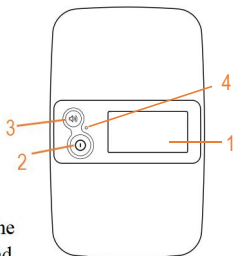Note:

Press either button to wake up the screen or turn off the reminder.

4. Indicator


- Green Slow Blinking: Battery $\geq 5\%$
- Orange Slow Blinking: Battery $< 5\%$
- Orange On: Charging
- Green On: Fully charged

3.2 Charging


- Charge the battery before use.
- Connect the station to computer or USB charging adapter using the provided cable.

4 Pairing the Station with the Main Unit


- a) Turn on both the station and the main unit.

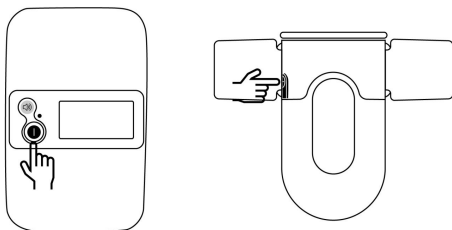- b) Hold the power button and the volume button of the station simultaneously to enable Bluetooth and search for the main unit.

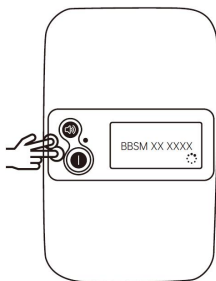- c) Press the volume button of the station to pair it with the main unit.
- d) When “√” pops up on the screen, the station and main unit are successfully paired.

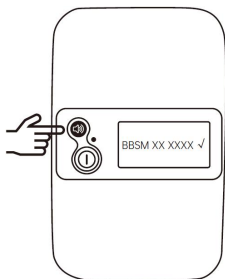
Note: You can check the serial number on the back of the main unit.

5 Monitoring with the Station


Press the power button for 2 seconds to turn on the station.

When the device is normally used for measurement, the screen of the station will display as follows:

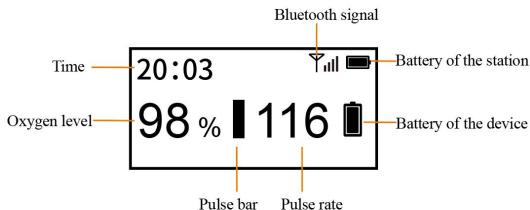
When the following situations occur, please check the device and station as soon as possible.

- The sensor cable of the device falls off.

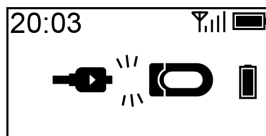
- The station is disconnected from the device.

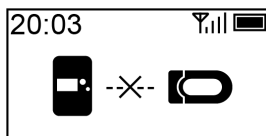
6 Stop Measuring

Unplug the sensor cable, the countdown on the screen will begin.

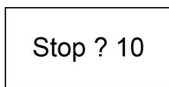Screen of the device

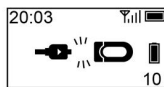Screen of the station

- During the countdown, if you plug the sensor cable again, the record will be resumed.
- After the countdown, the data will be ready for sync.

7 Cleaning

Use a soft cloth moistened with water or alcohol to clean the device surface.

8 Troubleshooting

Problem	Possible Cause	Possible Solution
No readings on screen of the station.	The sensor is not positioned and secured properly	Ensure that the device is worn correctly.
	Baby is moving frequently.	It will recover when baby stays still.
	Foot is too cold.	Make the foot warmer.
The station does not turn on or no response.	Battery may be low.	Charge battery and try again.
	Device might be damaged.	Please contact the seller or the after-service email.
The Station is disconnected from the device.	The device is off.	Turn on the device.
	The device is too far away.	Move the device closer.
	The Bluetooth connector failed to pair with the device.	Hold the power button and the volume button of the station at the same time to search for the device. Then, Press the volume button to pair it with the device.

9 Specifications

Weight	103g
Size	90mm×56mm×15 mm
Battery	3.8Vd.c., Li-ion rechargeable
Charge time	About 4 hours
Battery life	720 hours for typical use
Wireless	Bluetooth 4.0 BLE

10 Symbols

Symbol	Description
	Manufacturer
	Date of manufacture
	Serial number
	Indicates that the product should not be discarded as unsorted waste but must be sent to separate collection facilities for recovery and recycling.
	Refer to instruction manual
	MR unsafe.

IP22	Indicates that the product is protected against solid foreign objects of 12,5 mm Ø and greater; and protected against vertically falling water drops when enclosure tilted up to 15° .
	CE marking
	Authorized representative in the European community
	UKCA marking
	UK Responsible Person
	Indicates that the product complies with the applicable FCC requirements
	Non-ionizing radiation
	Indicates that the marked item or its material is part of a recovery or recycling process.
	Our products and packaging can be recycled, don't throw them away! Find where to drop them off on the www.quefairedemesdechets.fr site (Only applicable for French market).

11 FCC Statement

FCC Warning:

FCC ID: 2AD XK-1660

Any Changes or modifications not expressly approved by the

party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

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Shenzhen Viatom Technology Co., Ltd.
Prolinx GmbH
Prolinx Global LTD