


LITEC



Super Sky Lounge Massage Chair LC7800S+/LC8800

FINEVER

User's Guide

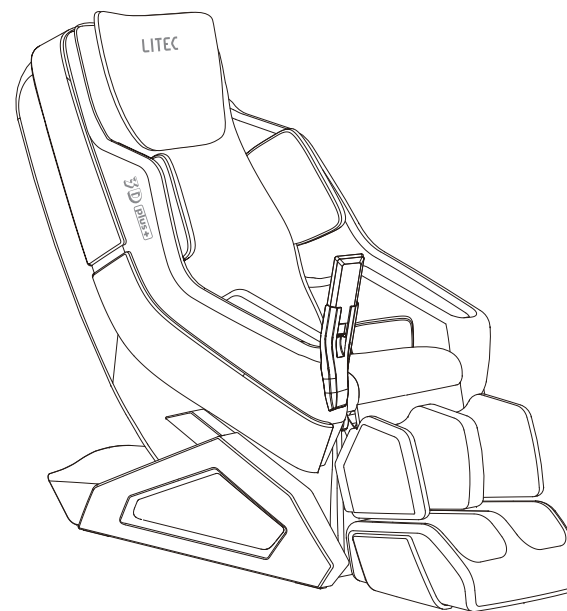

© 2015 Anhui Jiugong Wellness Industry Co., Ltd. All rights reserved.
LITEC, The LITEC Logo are trademarks of Shanghai Jiugong Industry Co., Ltd.,
registered in The People's Republic of China and other countries.



M_LC7800S+ -12 2015

Printed in China

English

MENU

Safety Precautions	1-4
Major Features And Functions	5
Components And Descriptions	6-9
Before Use	9-12
Adjust The Chair	12-14
Start Massage	15-17
Display of pictures	18-19
Music Play	20-21
End Of The Use	21
Troubleshooting	22
Regular Maintenance For Safety	23
Care And Storage	23
Specifications	24

Safety Precautions

- To clearly indicate the degrees of danger and damage, symbols “warning” and “caution” should shown below to highlight important safety instructions.
- Please adhere to all precautions.
- The Operating Manual contains important information about how to use the chair safely and correctly, avoiding possible damage to property or injury to persons using the chair and others.

 WARNING	This action could result in serious injury or death.
 CAUTION	This action could result in minor personal injury, or damage to property.

[Marginal Data]

	Prohibited (The symbol on the left means “Do not disassemble.”)
	Required Action (The symbol on the left means “Unplug the chair”)

- Remember to keep this Operating Manual in an easily accessible place.

 WARNING	
If you find any abnormal condition, please stop operation immediately. Electric shock or fire may result in if the chair is operated under such circumstances. Please turn off the main power switch on the chair and unplug the plug from the wall outlet. Contact the vendor for inspection or repairing.	
	● If there is a threat of lightning, abnormal noise (any abnormality), turn off the main power switch. For the safety, contact the vendor for inspection. Do not try to repair any malfunction of the chair.
	● Any moisture into the chair ● Any obstacle into the chair ● Any drooping down, tipping over or broken on the out housing
 CAUTION	
	● If you feel uncomfortable during operation, please stop using immediately. Consult with your doctors. ● Use the product only when it is placed on a flat surface. ● Before operation, make sure the switch is in order. Or it will result in injury and accident.

 WARNING	
 PROHIBITED	The following persons must not use this product. Persons who have thrombosis, serve aneurism, acute varicosity, or skin infection, including inflammation of the hypodermis.
	Persons with any of the following conditions should seek medical advice before using. 1. People who have electronic medical devices easily susceptible to electromagnetic interference. 2. People with malignant tumors. 3. People with serious heart disease 4. Woman who is in menstrual period or pregnant 5. Sensory disturbance caused by peripheral circulatory disorder due to diabetes 6. People with sensory disturbance 7. Persons being treated for personal injury 8. High fever (38° C or higher), such as acute inflammation, high blood pressure or chills 9. Persons with sickness or under bed rest orders 10. Osteoporosis, spinal fractures and acute pain 11. Inflammation due to sprains and muscle strain. It is not recommended for those who are in uncomfortable, sickness.

Keep massage to no more than 15 minutes in length. Do not continue to massage an area-specific of any particular of your body. Please stop for at least 30 minutes before reuse. Continue using for extend periods of time can result in over-stimulation on muscles and nervous.

 WARNING	
 FORBIDDEN WET HAND	Do not use the chair when your body is wet. It could result in electric shock or malfunction
 GUARD AGAINST DAMP	Do not spill water or other liquid on the chair. Do not place it in the humid environment like bathroom. It may result in electric shock or malfunction
 REQUIRED ACTION	Make sure no one especially pets is in the way when operation.
 DO NOT DISASSEMBLE	Do not try to disassembly or repair the chair, power cord and controller. It may result in fire and malfunction.
 PROHIBITED	Do not destroy the power cord. The power cord should not be broken or damaged. Don't bend, stretch, twist or wrap the power cord. Don't hang up objects on the power cord. The broken cord will result in electric shock, short circuit or fire.



WARNING

 PROHIBITED	This product is not intended for use by children and persons with reduced mental capabilities. Keep children and persons with reduced mental capabilities away from the chair during operation.
	Do not use the product with anything other than an AC 220 –240V power source. Do not wrap the power cord around the main body or it will result in electric shock, short circuit and fire.
	Do not let children sit and play on the product or it will result in injury and malfunction.
	During the operation, do not put your fingers into the gap between the housing and cushion or it will in injury and accident.
	Do not use the chair to massage head or stomach, or it may result in nausea or injury.
	Do not use the chair unclothed. It may harm your skin.
	Should the fabric cover be broken or removed, please do not use the product.
	This appliance is not intended for use by persons (including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance
	The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

If you feel uncomfortable during use, please stop the product immediately and consult with your doctors.



CAUTION

 REQUIRED ACTION	In the event of power failure, switch off the Main Switch and disconnect the power plug to prevent any sudden operation if the power is restored unexpectedly.
 FORBIDDEN WET HAND	Do not touch outlet, plug and power cord with wet hands.
	- After use or in the event of power failure, please disconnect the power plug to prevent electric leakage, electric shock and accident - In the event of cleaning the product, please disconnect the power plug to prevent electric shock and fire. - When disconnecting the power plug from the socket, be sure to hold the plug instead of pull over the power cord. Or it will result in electric shock, short circuit and fire. Please check power plug and power cord as schedule.
 PROHIBITED	- Do not use the product immediately after drinking and dinner. Or it will result in uncomfortable. - Do not hang up heavy objects with power cord and control wire. Do not drop the controller or it will result in electric shock and fire - Do not use the product near the fire or heater device. Do not smoke while you are using. Do not place the product on the floor where there is underground heating or other heating device. Or it will result in fire.



CAUTION

 PROHIBITED	- Do not pull or crash the fabric cover with sharp implement, like knife, or it will result in injury and malfunction. - Do not unduly twist stretch the power cord or it will result in fire. - Should any broken on the power cord and power plug, or any loosen on the outlet, please do not use the product. Or it will result in electric shock, fire and accident. - During the operation, do not put your fingers into the gap between the housing and cushion or it will in injury and accident. - Do not force overloading on each function course or it will result in malfunction and injury. - Do not use the product when you fall in asleep. Use with unconsciousness will result in injury and uncomfortable. - Do not unduly pull over the product from the wooden floor, carpet or other place which is easily scratched. - Please carry and place the product carefully to avoid scratching the floor.



CAUTION

 PROHIBITED	This product is designed for home use and not for commercial purposes.
	Do not operate under blanket. Excessive heating will occur and cause fire and accident.
	Make sure the well connection between the plug and the outlet. Wipe dust away from them. Or it will result in electric shock, short circuit or fire.
	Do not stand or tread on the product. Do not force overloading on the product. Or it will result in injury, malfunction and accident.
	Do not use the product near fire resource or flammable gas. Or it will result in product distortion, malfunction and fire accident.
	Do not use the product where aerosol spray products are being used or where oxygen is being administered. Or it will result in product distortion, malfunction and fire accident.
	Do not clean or spray the product with abrasive cleaning detergent, alcohol or insecticide. Or it will result in product distortion, electric shock and fire accident.
	Please keep children and persons with reduced mental capabilities away from the running product. Or it will result in injury.
	This product is not intended for use by persons with reduced mental capabilities or under bed rest orders. It may result in personal uncomfortable.

Major Features And Functions

1 SENSOR FIT:

The system will detect the user's body figure smartly, before the massage begins.

2 3D CORE:

Based on the user's back curve, the 3D function will conduct the massage, as if done by human.

3 ZERO GRAVITY :

Imitating the zero gravity condition in the space capsule, it would help the user relax totally.

4 HIP SHAPE:

Through three groups of airbags at both sides of hip and below hip, it can help tighten up and rectify hip, to achieve a dimensional hip shaping effect perfectly.

5 HUMANE TOUCH SCREEN:

The highly user-friendly humane interface system can lead the user to enjoy thousands of massage skills.

6 AUTO PROGRAM:

Various programs have been designed by Japanese massage masters, which would provide the most genuine massage skills.

7 MAGUNES:

MAGUNES is developed based on rich experience and performance, with strong power and excellent muting effect.

8 AEROSLIM:

Flexible airbags is an aero system with strong function and extensive moving range, the massage chair with this technology realizes the wide-range potent traction of human body, so that it can help users keep slim.

9 CARBON FIBER HEATING:



The most advanced carbon fiber heating film in the world adopted, which is safe and highly-effective, making sure that heating function is provided for user's feet and legs.


10 SHIATSU FOR HAND:

It helps promoting blood circulation for hands and front arms, under natural sitting habit.

Components And Descriptions

● Body

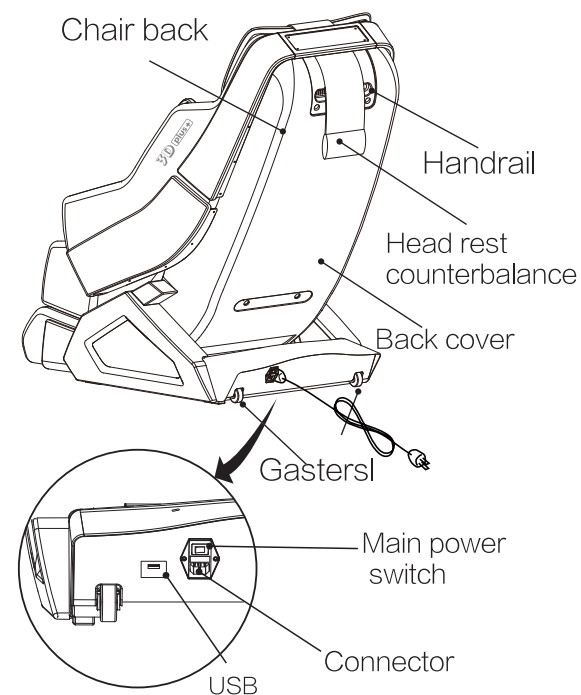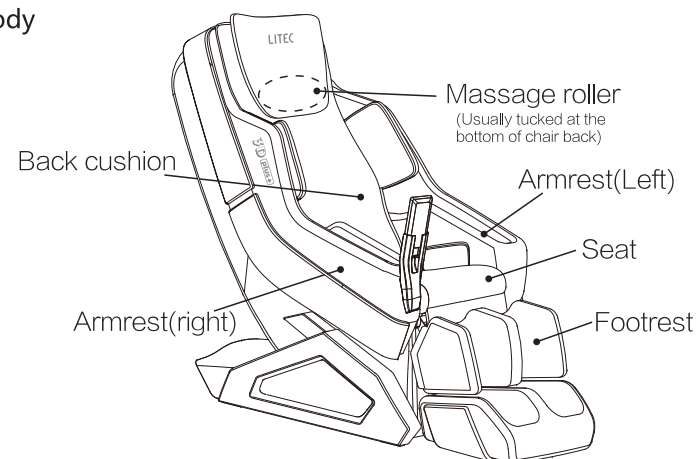

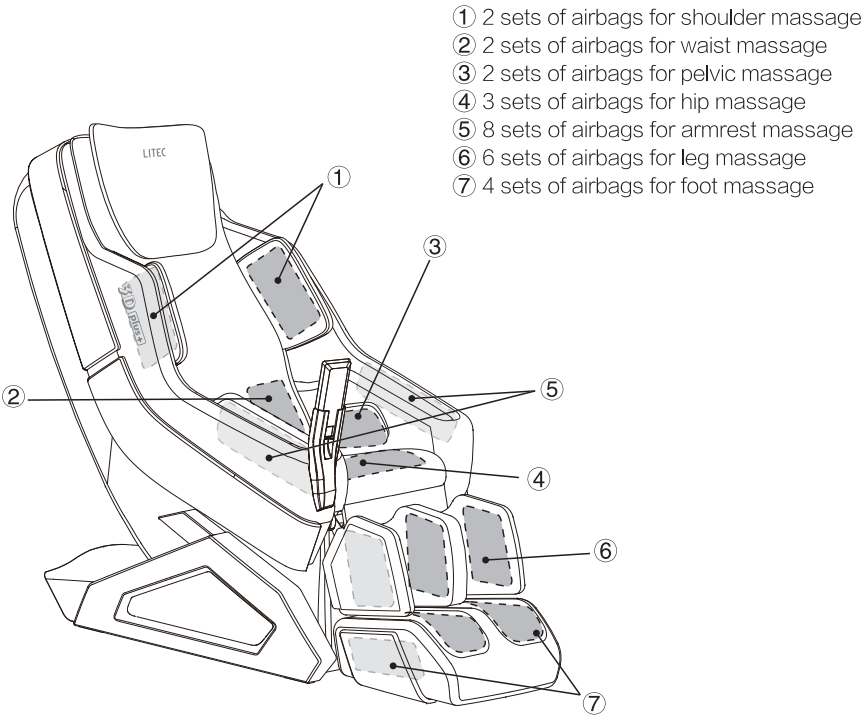


●Power cord and outlet

One end of the power cord will be plugged into the outlet of armchair, and the other end will be plugged into the mains jacket.
Attention 1. Make sure the main power switch is closed before it is plugged into the outlet.
Attention: 2. In case of disrepair of the power cord, the arbitrary dismantling or change is forbidden, and it should be repaired or changed by the dealer or qualified professionals.

●Main power switch

Open the main power switch, the armchair moves into the standby mode. Before the start button of the controller is pressed, only the seat angle and the leg part angle are adjustable.

●Controller

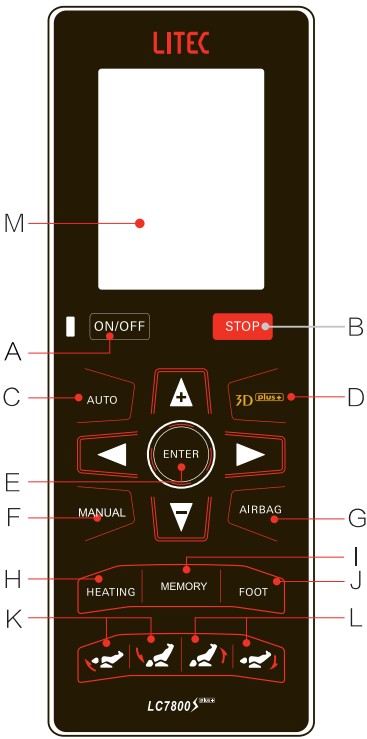

A	ON/OFF	Enter or exit the standby mode. when the main power switch is pressed, the red indicator light will be on, and it will be off after the switch is switched off.
B	STOP	The button should be pressed in case of emergency or in case that you feel the abnormal operation of the armchair. All the airbags shall be deflated and other moving parts shall stop operating.
C	AUTO	There are 6 optional automatic massage process. Press any button and the corresponding automatic massage process can be operated. In the meantime, the massage strength,3D smart function or overlaying hot compressing is adjustable.
D	3D PLUS	By pressing the button, you can enter the interface where you can adjust the strength of 3D on the massager.
E	NAVIGATION AND ENTER	By pressing the button, you can enter or browse the menu, move the cursor, make choice or confirmation, adjust 3D smart function, massage speed, width and strength, set the massage range and adjust the strength and range of the optional procedure.
F	MANUAL	You can control and choose manual roller procedures, by pressing this button.
G	AIRBAG	You can control and choose manual airbag procedures, by pressing this button.
H	HEATING	You can start or stop heating function, by pressing this button.
I	MEMORY	You can save and customize the manual mode according to your preference.
J	FOOT	You can control and choose manual pedicure procedures, by pressing this button.
K	LEG ANGLE ADJUSTMENT	By pressing ascend or descend button, you can adjust the leg angle.
L	SEAT ANGLE ADJUSTMENT	By pressing ascend or descend button, you can adjust the seat angle.
M	LIQUID CRISTAL DISPLAY	The display of operating menu and operating conditions.

●Audio wire and audio socket

Use the audio connecting wire of the accessory to connect the music player and the audio socket, so as to play the music by operating the music player. if you want to use the Bluetooth to connect, then this audio wire is not needed.

Warning	- Do not splash other liquid or water onto the armchair or wire controller. - Do not operate with wet hands.
Caution	- Do not sit on the wire controller or drop, trample it; do not pull the wire of the controller by force. - Do not press the LCD or touch the panel with sharp or hard articles, otherwise it may cause damages.

Before Use


- Assemble the chair according to the Assembly Guide. The product is packaged with the chair back lying flatwise. Before use, the chair back should be raised.
- When using the armchair, do not increase or decrease the room temperature all of sudden. The rapid change of temperature may lead to the rusting of metal parts, and the operation of the armchair may also be disrupted by the short circuit or other faults caused by the moisture of parts of the electrical appliances.

About the placement venue


Please ensure that there is about 50 cm space in front of and behind the armchair, to prevent the armchair touches the wall or other articles when it is inclined.

Caution	- Do not use it in the bathroom, sauna house or other places with high humidity. - Place the machine on a flat and slippery-free floor for operation. - Do not use it in places with sundries and dust or with high ambient temperature (above 40°C). - Do not place the machine in places with direct sunlight, or near heating equipments with high temperature. - Do not let cigarette and other open fire-inducing items approach the armchair.
----------------	---

- Lay a cushion under the armchair to prevent damages to the floor and for reduction of operation noise. (Never lay an electric carpet)
- Due to the fact that the armchair is heavy, make sure that the floor will not be damaged when placing it on the floor.
- Please put the armchair to places 1 meter away from the TV, radio and other audio and visual machine, otherwise the interruption noise may occur.
- The armchair is off limit to users with weight over 100 kilograms.

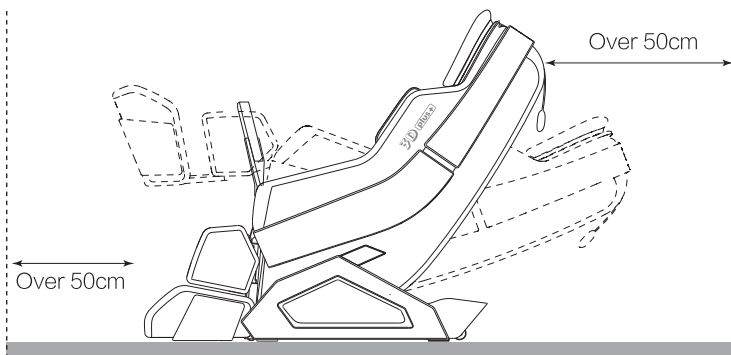


Confirm the articles inside the packaging case

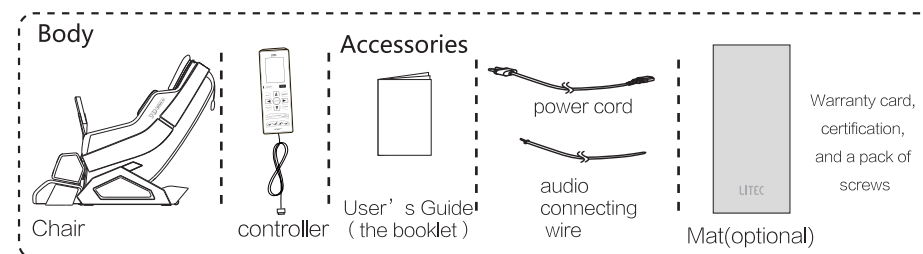


How to move

When moving the product, close the main power switch, and unplug the power plug.

When using the gasterl, straighten the chair back, grasp the two sides of leg parts to lift the upper part, and move the armchair by moving the gastersl.

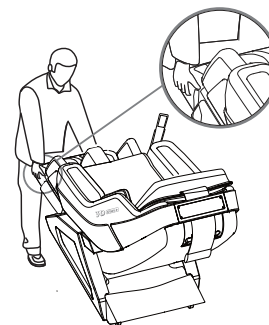

When lifting and moving the machine, lay flatwise the chair back, catch the handrail of the chair back to lift it so as to move the armchair.

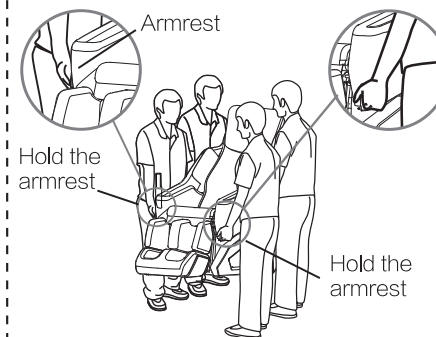

In order to prevent damages to the floor, please move the armchair after laying cushions on the floor..

In case of the falling of the machine, this may damage the floor or cause personnel injury, so please hold fast to the armrest and leg parts.

Caution	● The armchair cannot be moved if the passage is less than 80 cm wide. ● Contact the dealer in case you cannot move the armchair normally. ● Due to the heavy weight of the armchair, make sure that it will not damage your waist or other parts of your body when moving the machine. ● The remote controller is connected to the armchair, so do not drop it when moving the machine. ● Do not move the armchair with person loaded. ● Please move the armchair after unplugging the outlet. Otherwise it may cause damages to the power plug, leading to fire disaster or electric shock.
----------------	--

How to connect the power source

- ① Plug the power cord to the switch outlet on the back of the armchair and plug the other end to the power outlet.
- ② Turn on the main power switch at the back of the chair.




※In case that the massage rollers gets stuck, turn on the main power switch and press the start button to return the armchair to the inception condition.

⚠ Warning	● Do not plug the power source to an external transformer.
⚠ Caution	● Before plugging the power source to an outlet, make sure the main power source is closed. ● Please fully plug the power plug. ● Before use, carefully read the specifications, and set the armchair in light of the requirement to ensure the normal operation.

Connect power cord and wire controller

In case of first use of the armchair or reuse of it after a long time, please examine soundness of the power cord, plug, the armchair outlet, main power switch, and the wire controller, and then connect the power cord to the outlet, connect the wire controller to the chair.

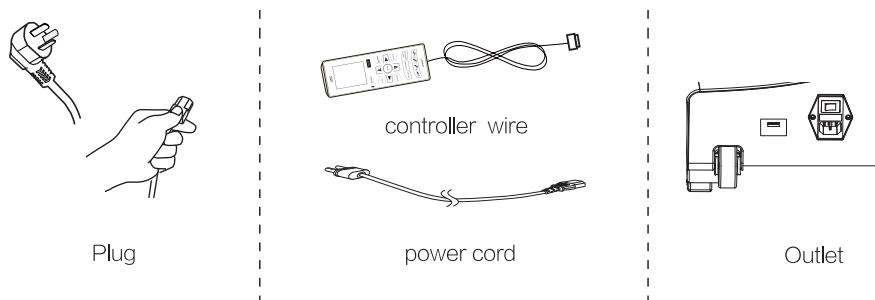


Confirm the material


Before using the armchair, confirm whether the material is damaged or not according to the following steps:

- ① Unzip the back cushion

Turn over this way

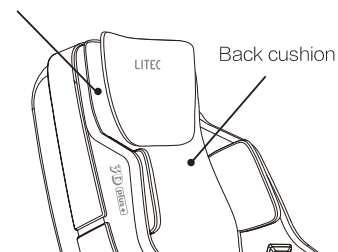

- ② Confirm whether the material of the armchair is damaged or not by putting your hands into the space between the back cushion and the chair back.

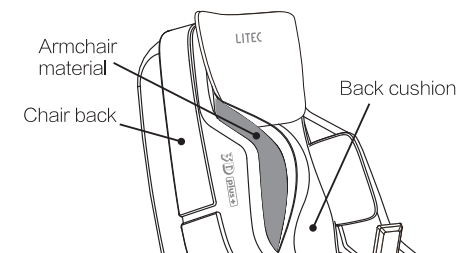

Caution

Before using the chair, confirm that whether there is damage in the material of the armchair or of other parts. In case of damage, even the smallest tearing, in the chair, please stop using it immediately, turn off the power switch, unplug the chair, and contact the customer service center. Never use the armchair with dismantled cushion or damaged material, otherwise there will be serious safety hazard.

Adjust The Chair

Sitting in the chair

Under general circumstances, the back rollers are laid folded in the left and right side of the chair back, and at the same time, the footrest are descended to the lowest position. Please ensure the position of the massage rollers and the footrest for the sake of safety. If the massage rollers or footrest are laid otherwise, restart the chair. After all parts return to the initial condition, sit on the armchair slowly and put your foot into the footrest.


After use of the chair and before standing up from the chair, remove your feet off the footrest first.



Caution

- After the footrest are raised, do not sit on the chair again, otherwise great pressure may be caused to the footrest, leading to the falling, damage or other accident of the armchair.
- Raise the chair before you use it, enjoy the massage from light level, and then adjust the angle of your lying position in the light of your personal demand.
- Never stand on the footrest, otherwise the footrest or even the whole chair will be damaged.

Adjust the angel of the chair

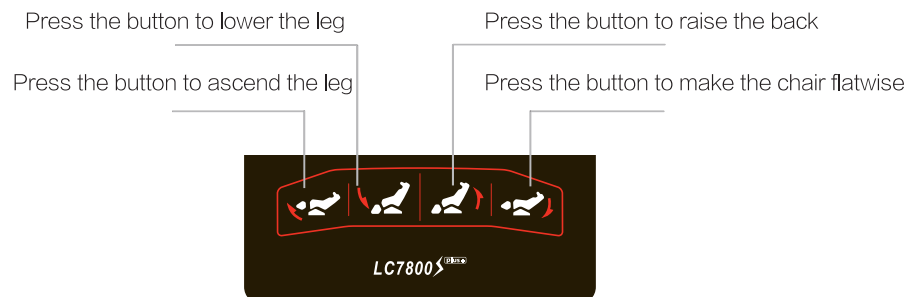

Long time-pressing function

Pressing the angle adjustment button for over 3 seconds with the indicator lights flashing, the armchair will adjust angles automatically until the button or reverse button is pressed again.


Attention: when the main power switch at the back of the armchair is on, the angle is adjustable even the on-off key of the wire controller is turned off. However, the angle cannot be adjusted in case of emergency stop.

Warning

Ensure there is no person (especially children) or pets around when you adjust the angle of the armchair. Moreover, ensure at least 50 cm space in front of and behind the armchair, and there is no obstruction in the trajectory.

Adjust the footrest

If you feel that your legs cannot be fully stretched, you may properly stretch your leg to lengthen the footrest, and the lengthening range for the chair is 160 mm. When the armchair is returned, you only need to retrieve your legs, or use the force of the floor produced when legs are descended.

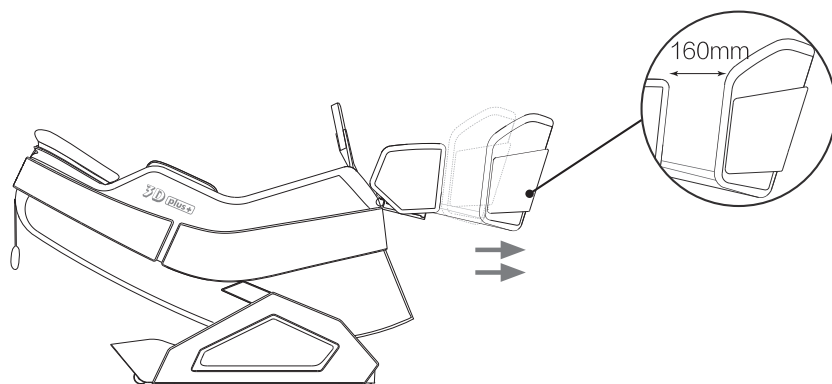


Warning


Ensure there is no person (especially children) or pets around when you descend the footrest of the armchair. Moreover, please do not let your hands or other articles under your footrest.

Adjust the headrest

- The position of the headrest is only adjustable upside and downside. It keeps balance depending on the counterbalance attached to the headrest on the back.
- Put the headrest on the right position before you use the chair. Incorrect position may cause dysfunction or discomfort.

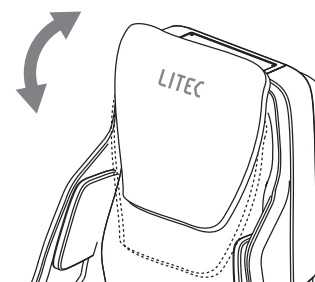

- If you want to have stronger neck massage, please lift the headrest over to the top and fully release rollers 3D travel effect.


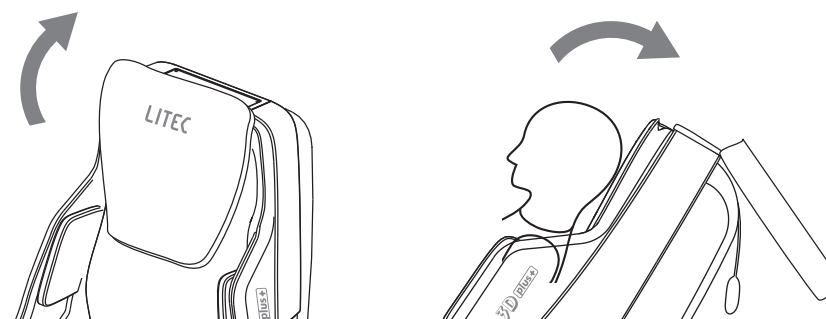


CAUTION

When the headrest is lifted, the back pressure will be stronger. It is suggested to massage less than 15 minutes every time. There is a counter balancer inside the headrest. Please do not play with it. Or it will result in injury. Please keep children away from the counter balancer

Start Massage


1. Connect the main power, and turn on the switch

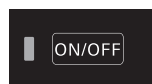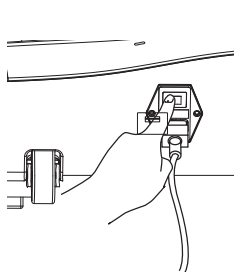

Switch, emergency stop.

2. Turn on the switch, and enter the interface of AUTO procedures


The screen would twinkle, to guide users to choose AUTO message procedures.

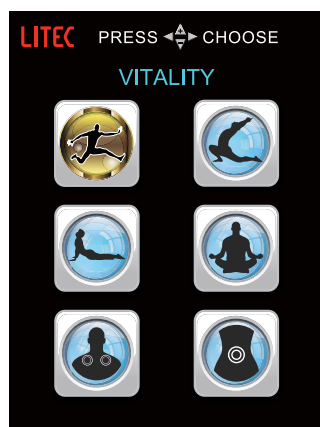

Press other manual buttons, to choose manual massage.




Press Customize button, to enter the MEMORY interface.

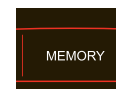

3. Navigation and enter button

By pressing the button, you can enter or browse the menu, move the cursor, make choice or confirm the Auto procedure, manual procedure, route, skills, strength, and range.

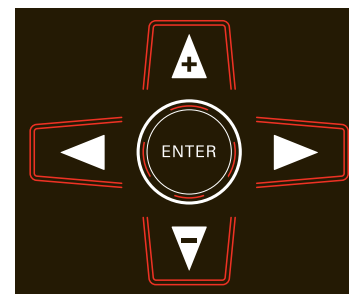

4. By pressing 3D button and Navigation buttons, you can adjust the 3D strength during the massage.



5. STOP

The massage chair will stop operating, after you press this STOP button. Meanwhile, for the safety of the user, all mechanical movements should be locked, and all airbags shall be deflated.



Safety Notice:

If any problem occurs during the massage, the chair will get into an emergency stop state automatically.

For the safety of the user, all mechanical movements should be locked, and all airbags shall be deflated.

Press the ON/OFF button, after you confirm that nothing is abnormal nor dangerous, and the massage chair will get back to its initial state.



Caution:

The STOP button is only meant for the protection of the user, under emergency state. Generally, please do not use this button as the switch to change the procedures, remove settings, or get back to its initial state.

When operating the automatic mode, the roller massage, air massage, and foot massage should be processed in due procedure, and the users cannot change them arbitrarily. The heat compressing is not set in advance, where start or dose is chosen according to the user's demand.

The opening and closing of heat compressing will not impact the counting of the massage procedure. When the massage procedure is completed, the heat compressing will be closed as well.

The operating methods, body detecting, and screen display would be similar for all the 6 AUTO massage procedures.

When operating the automatic mode, the speed the roller massage and foot massage proceeds in due procedure, and the user can not change arbitrarily. The strength of 3D air massage is set as the middle level in advance, which can be adjusted according to the demand of the users.

During the operation of automatic mode, the 3D strength, roller strength, airbag strength, and pedicure strength are shown below the screen.

Press any button, including the AUTO button, to get back to the main menu.


Display of pictures

If chosen, it will appear as colored;
if not, it will not appear or just appear as white.


Massage
main menu




Body
detecting

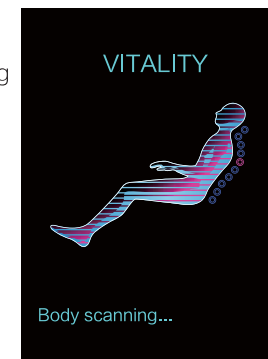


Roller
massage

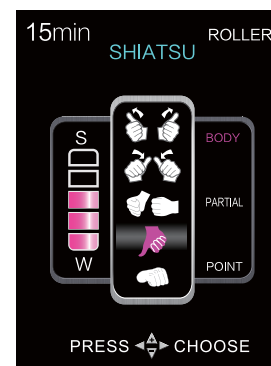


Route

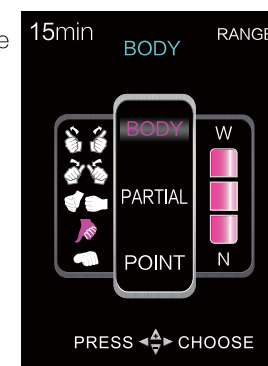


Roller
strength

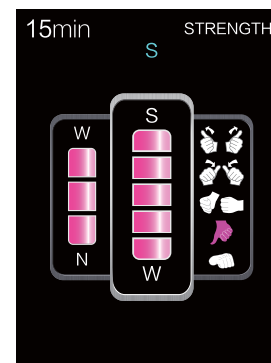


Roller
width

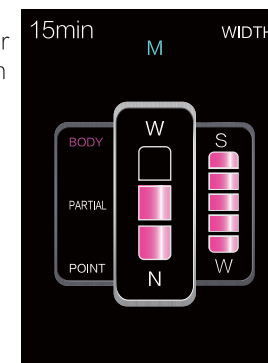


Airbags

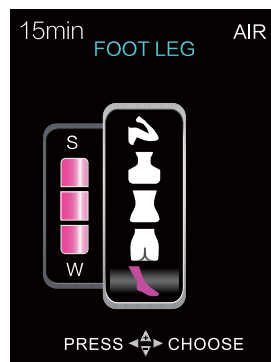


Airbags strength

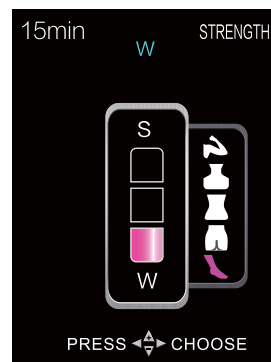


Customize

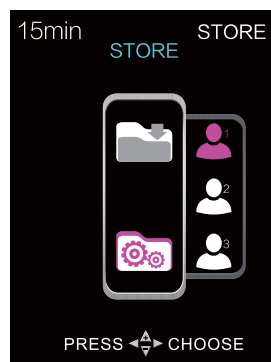

RUN

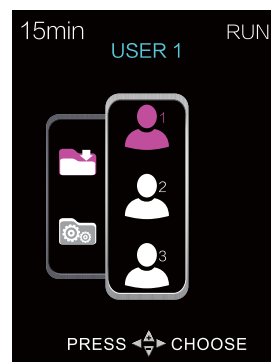

Manual 3D

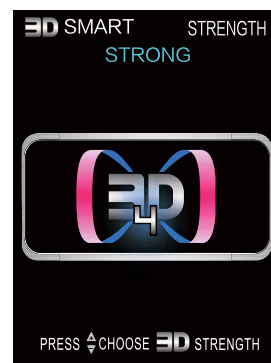


AUTO 3D



Music Play

There is a built-in speaker in the machine, the user may connect the music player and the chair by Bluetooth or audio connecting wire(bundled), then operate the music the music player and listen to the music from the speaker. Beautiful music can be played during massage to cultivate artistic atmosphere, combined with various massage skills, to promote the massage effect hugely.

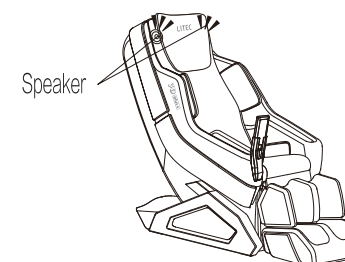

Attention

While playing the music, turn down the volume of the player, and turn up gradually after you adapt yourself to it in case it may cause shock or irritation

This device complies with Part 15 of the FCC Rules / Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes :

- (1) l'appareil ne doit pas produire de brouillage, et
- (2) l'utilisateur de l'appareil doit accepter tout brouillage radio électrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Music Play

Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power

(e.i.r.p.) is not more than that necessary for successful communication.

Conformément à la réglementation d'Industrie Canada, le présent émetteur radio peut fonctionner avec une antenne d'un type et d'un gain maximal (ou inférieur) approuvé pour l'émetteur par Industrie Canada. Dans le but de réduire les risques de brouillage radioélectrique à l'intention des autres utilisateurs, il faut choisir le type d'antenne et son gain de sorte que la puissance isotrope rayonnée équivalente (p.i.r.e.) ne dépasse pas l'intensité nécessaire à l'établissement d'une communication satisfaisante.



Attention

End Of The Use

Raise the seat completely, and descend the legs fully. Turn off the chair by pressing the ON/OFF button. After the chair goes back to the initial condition (the airbags get deflated completely, the massage rollers move to the top, the foot massage stops, and the hot compressing closes), the foot part should be removed from the leg part. After it stands steady on the floor, stand up from the chair. Do not stand on the footrest to prevent falling down or damaging the armchair.)

Turn off the main power switch, and unplug the power plug from the socket.




Attention

After you use the chair, turn off the main power switch, and unplug the power plug from the socket.
→Otherwise, accident or injuries may be caused by children or pets.
After you use the chair, be sure to raise back completely and lower the footrest completely.
After you use the chair, please clean it and place it in dry places.
When unplugging the power cord, please do not pull it directly; instead, you should unplug it pinching the power plug.
Please be sure to put the power cord away in prevention of tripping over.

Troubleshooting

If there is still a problem after trying the solutions given in this troubleshooting guide, stop using the chair immediately to prevent accident or injury, and unplug the power plug from the wall outlet. Contact the vendor or your nearest service center for assistance.

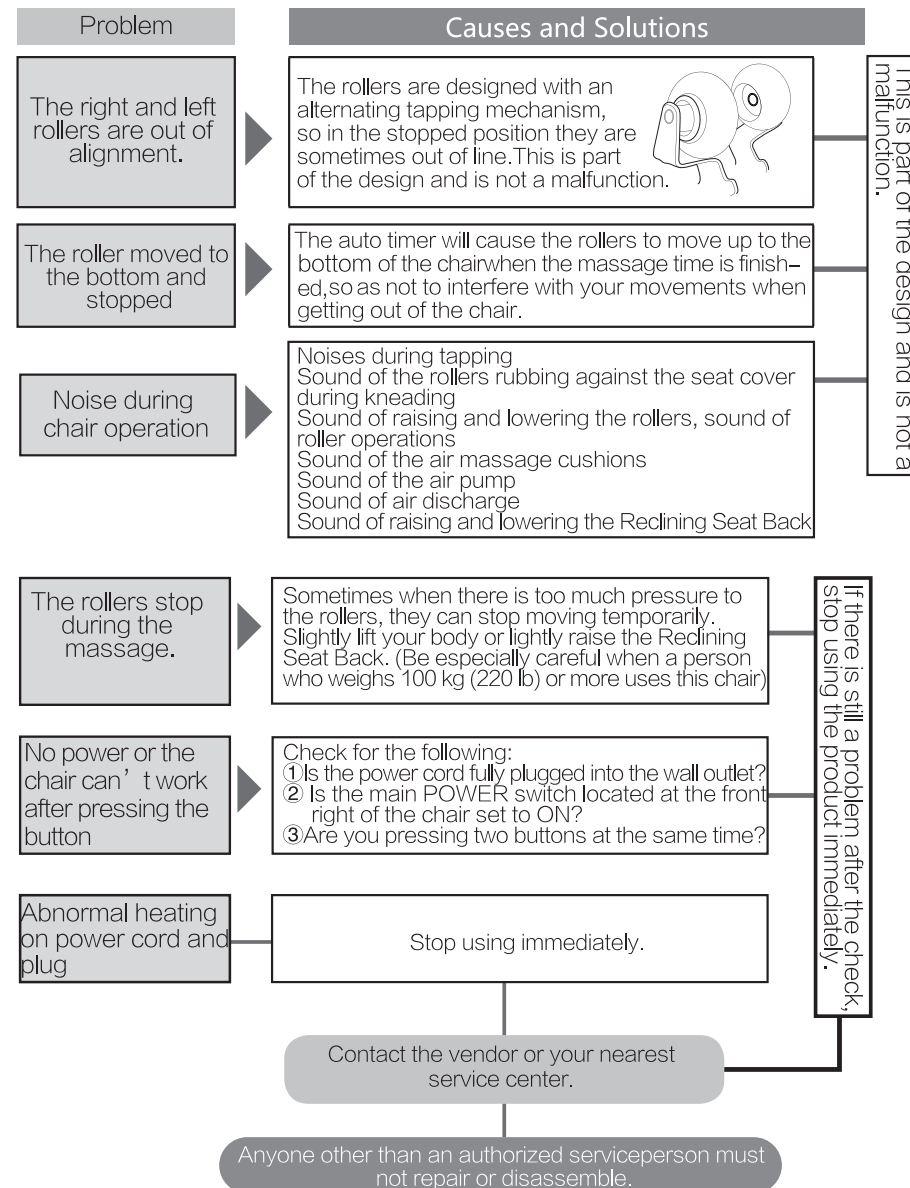

Regular Maintenance For Safety

Regular Maintenance Will Preserve the Life of Your Chair

Have you experienced any of these problems?

- The chair does not start with the Main POWER switch set to ON.
- The power cord has been cracked or scratched.
- The power cord or plug becomes abnormally hot.
- The chair produces a burnt smell, abnormal sound or vibration.
- The auto timer function does not work in some cases.
- Other failures.

If you find anything listed on the left, stop using and unplug the power cord from the outlet to prevent failure or an accident. Be sure to contact the vendor or our service center for maintenance.

To use safely

The chair fabric (roller fabric) is a consumable part, and its service life is 3 years (conditions of use: two 15-minute sessions a day, every day). For safe use, periodic inspection by the vendor is recommended before the end of its service life.

*Durable years vary depending on conditions of use.

* After about five years of use, a professional inspection is recommended even if there are no symptoms or problems. For inspection fees, contact the vendor.

Care And Storage

Storage

- After wiping off any dust or dirt, store your chair in a place that is relatively dry.
- If the chair is left unused for a while, put a cover on it that will protect it from dust and dirt.
- Pets can sometimes chew on power cords, so try to keep pets away from the chair even when it is being stored.

⚠ Caution ● Do not install the chair where it can come into contact with direct heat, such as direct sunlight or a radiator.

Care and Maintenance

⚠ Caution ● When cleaning the chair, be sure to unplug the power plug from the wall outlet first. Do not plug or unplug the power plug with wet hands.

[Caring for the Back Pad and Covers]

Soak a piece of soft cloth in lukewarm water in which a small amount of mild detergent has been added. Then thoroughly wring it out before cleaning off any spots or grime.

⚠ Caution ● The fabric color can change if the chair is in contact with vinyl covers or wall coverings for a prolonged period of time.

[Caring for Other Plastic Parts]

Soak a piece of soft cloth in lukewarm water in which a small amount of mild detergent has been added, thoroughly wring it out, and clean off any spots or grime. Then wipe with a dry cloth so that no detergent residue remains.

⚠ Caution ● Do not use alcohol-based or other household cleaning products or abrasives because they can scratch or crack the finish or cause color to fade.


[Caring for the Remote Control]

Wipe dirt off the remote control with a dry cloth.

⚠ Caution ● Never use a moistened or wet cloth. May result in a mechanical failure.

Specifications

Model:	LC7800S+/LC8800
Voltage:	AC220-240V~ 50/60HZ
Power:	130W
Time:	15 mins
Size:	width*length*height about 78cm*110(180)cm*80(130)cm Numbers in the brackets are when the chair reclines the most
Weight:	95KG
Outside material:	synthetic leather, fabric
Supplied items:	Power cord, audio cord, instruction book, mat(optional) warranty card, certification, and a pack of screws
Manufacturer:	Anhui Jiugong Wellness Industry Co., Ltd

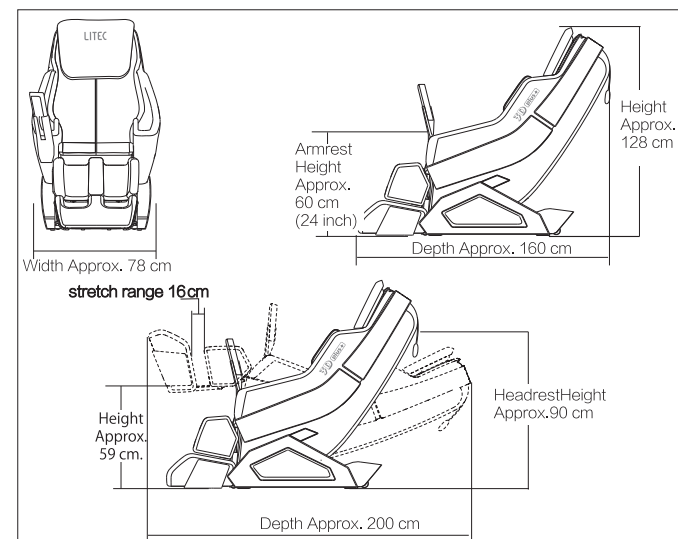

Correct Disposal of this product



This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.