Live In The Know

37degree

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How to wear

Install APP & wristband charging


Install APP (Figure 1) : scan the QR code at the bottom of the packaging to install APP.

(Figure 1)

Wristband charging (Figure 2) :

1. When 37° Journey wristband starts charging, please aim the metal contact point on the back of the tracker at the point in charging adapter, then the tracker and the charger will be automatically buckled;
2. Please plug the USB port of the charger into the computer or any other USB interface for charging; when charging, the charger' s green light is flashing;
3. The charging time is about 2 hours, and green light is always on when it is fully charged.

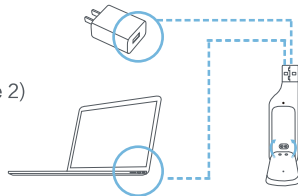(Figure 2)

How to wear

Install APP & wristband binding


Assemble the wristband (Figure 3) : please put the tracker into the wrist strap.

(Figure 3)

Wristband binding (Figure 3) : when keeping the wristband in the state of charging, switch the mobile phone Bluetooth on. Search the wristband according to the APP tips and refresh the APP for viewing the result until the completion of the binding.


(Figure 4)

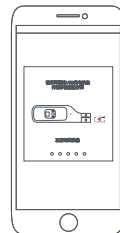How to wear

The best way to measure

Figure 5 The wristband is best worn about one inch from the ulnar styloid of left hand.

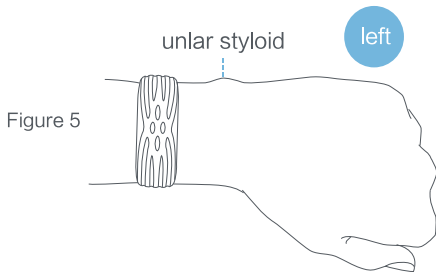Specification

Basic parameters

	Journey
	Weight : 24g
	Full length : 238 mm
	Adjustable length : 133 — 209 mm
	Waterproof level : IP67
	Battery type : Lithium polymer battery
	Battery capacity : 157 mAh
	Connected by : Bluetooth 4.1
	Working temperature : -10° C — 45° C
	System requirement : System for Android 4.3 and above version & IOS7 and above

Specification

Indicator lights:

Green light ---- Fully charged


Green flash ---- Charging

Indigo speed flash ---- Wristband paired successfully

Yellow flash ---- Wristband turning off or entering flight mode

Yellow light ---- Wristband turning on or exiting flight mode

Red flash ---- Please open the app and connect the wristband to synchronize data.

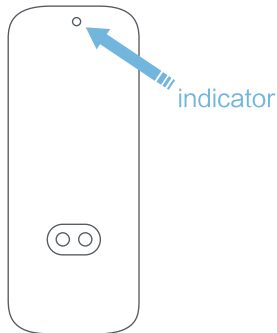Specification

Names and content of manufacturing substances

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Remarks

○ Means the toxic substance content in the materials used in the part conforms to SJ/T 11363–2006 regulations.

× Means the toxic substance content in at least one of the materials used in the part is higher than the limit stated in SJ/T 11363–2006.

This product conforms to the RoHS Directive of the E.U.

For those with sensitive skin and in hot and humid environments, do not wear the bracelet tightly for extended periods of time. Otherwise, rashes and itches may occur due to the skin's lack of opportunity to breathe.

Product declaration

This product measures biorhythms by monitoring ambulatory blood pressure changes to help the user live a healthier life. It cannot be used as a basis for any medical diagnosis.

The 37° Wristband is not categorized as medical apparatus. The 37° Wristband, its related applications, and information thereby obtained shall not be used in medical diagnosis, treatment, or relieving of any disease or any conditions. To prevent injuries from excessive exercise or user's self-diagnosis, please consult a licensed health professional before making changes to your exercise, sleep, nutrition, or medication.

Note: Ambulatory blood pressure is not the same concept as clinical blood pressure used in medical settings

Warranty and after-sales service

Please visit our official website www.37body.com for instructions, solutions to common problems and quick start guide.

Product warranty

The “limited warranty” starts from the date of purchase and expires 1 year later. After it expires, you can't continue to enjoy the after-sales service. In the case where a store/retailer (importer or distributor) provides you with prolonged warranty as a part of the sales, 37° will not take any responsibility. For more relevant information, please contact the store/retailer. This limited warranty will automatically expire without notice and be unenforceable in the following cases (such cases are left to the discretion of 37°):

1. The product is opened or any part of it is changed by the customer;
2. The product is repaired or dismantled by non-37° workers;
3. The serial number is wiped out, changed or made difficult to recognize by any means.
4. The product is purchased from dealers unauthorized by 37° (please contact 37° for the list of authorized dealers).

Guarantee of repair, replacement and refund

37° wristband's promise of after-sales service:

Warranty and after-sales service

We promise to offer service of repair, replacement and refund after sales according to “Law of the People’s Republic of China on Protection of Consumer Rights and Interests” and “Law of the People’s Republic of China on Product Quality” . Within the warranty period, you can have the product repaired, replaced or returned by virtue of your receipt or proof of purchase according to this warranty.

1. Within the warranty period, you can have the product repaired, replaced or returned by virtue of your receipt or proof of purchase according to this warranty;
2. Within 7 days from the date of purchase, you can enjoy free services such as return and replacement if 37° wristband has any performance problem;
3. Within 15 days from the date of purchase, you can enjoy replacement and other free services if 37° wristband has any performance problem;
4. Within 1 year from the date of purchase, you can choose to have 37° wristband repaired, replaced or others after getting it tested in 37° after-sales service center, by virtue of your receipt or proof of purchase, if 37° wristband has any performance problem.

Product warranty restrictions

In situations as follow, which are not in accordance with the warranty conditions, customers can’t apply for any warranty service stipulated in this limited warranty or other documents:

Warranty and after-sales service

1. Damage due to abuse, misuse or carelessness (including bending or dropping the product, not limited to the aforesaid situation);
2. Damage due to other use that is not as expected or recommended;
3. Damage due to the refit of 37° product, including damage to the normal function of the product by exposing the product in the harmful environment, damage to the outer layer of the product or visible cracks in the surface, etc;
4. Damage due to the use of any product, accessory, software or service that are not made or provided by 37° (specific situations are left to the discretion of 37°).

Warranty and after-sales service

Contact 37°

If you'd like more detailed information about warranty, please ring 37° hotline.

Hotline: 400-995-3773

Refund and replace

1. Purchase from a retailer;
2. Purchase from our website

Purchase from a retailer

If you purchase a 37° product from a retailer, you must contact the retailer for replacing or refunding.

Purchase from our website

You can get full refund on the product purchased from 37° website within 7 days unconditionally. If you're not completely satisfied with the product purchased from our website for any reason, please contact us for refunding. To receive the full refund (exclusive of shipping fees), please make sure that the product is sent back unused with its original package.

Note: All sales return must be authorized by 37° beforehand.

For more information about sales return, please ring 37° hotline 400-995-3773.

Warning

37degree

FCC Regulations:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.`

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiated radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.

Warning

- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

37degree

FCC Note:

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user 's authority to operate the equipment.

RF Exposure Information

This device meets the government' s requirements for exposure to radio waves.

This device is designed and manufactured not to exceed the emission limits for exposure to radio frequency (RF)

Warning

energy set by the Federal Communications Commission of the U.S. Government.
This device complies with FCC radiation exposure limits set forth for an uncontrolled environment.

Assessment of compliance of the product with the requirements relating to the Radio and Telecommunication Terminal Equipment Directive (EC Directive 1999/5/EC) was performed by Telefication BV (Notified Body No.0560),

CE 0560

37degree