





Quick Start Guide

Wearable Breast Pump S2 Pro

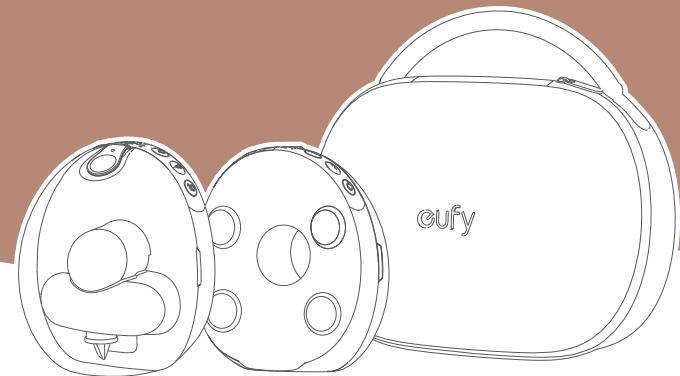

User Manual/Video/Support

Breast Pump Overview

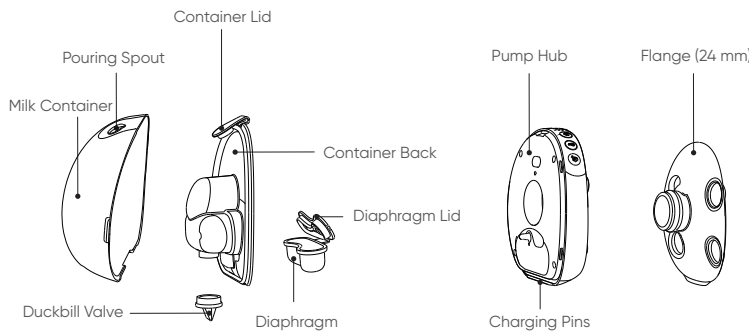

Accessories also include:



Charging Case



Magnetic Charging Cable



Stain-proof Mat




Extra Duckbill Valves




Sizing Inserts (17/19/21 mm)



USB-C Charging Cable

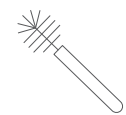

Cleaning Brush

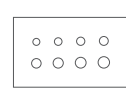


Nipple Measurement Card


Charging

Fully charge the pump hub before initial use; it won't operate while charging.

- Option 1: 2-in-1 Magnetic Charging Cable

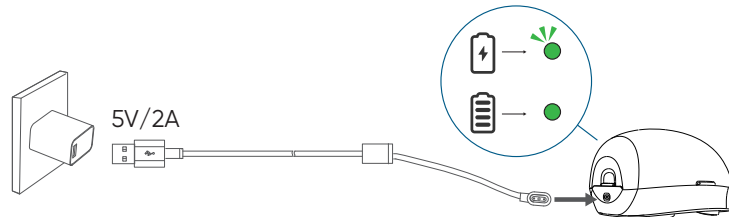

- Option 2: Charging Case

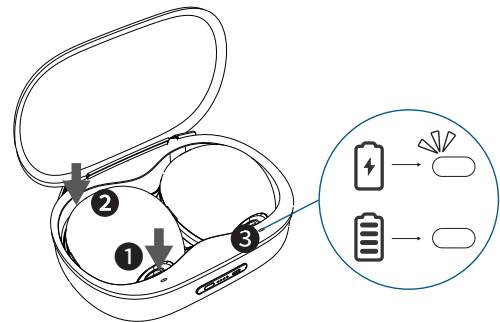


Step 1: Align the pump hub's charging pins with the port inside the case until you hear they click.

Step 2: Let go of the pumps to allow the entire pump to be placed in the slot.


Step 3: Ensure the indicator light flashes white, indicating the pumps are charging properly.

Why does my charging case or pump hub not boot up?
If you encounter a pump hub or charging case that does not work, try the following steps:
1. Press and hold the power button on the pump hub for 12 seconds to reset and see if it turns on.
2. Be sure to use a 5V/2A or higher current rating adapter.
3. Adjust or swap the position of the breast pumps to ensure proper contact.
4. Use the included magnetic charging cable for direct charging.


Disassembling

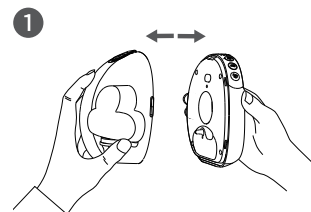


Remove the container from the pump hub.

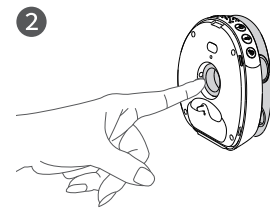


Press to release the flange from the pump hub.

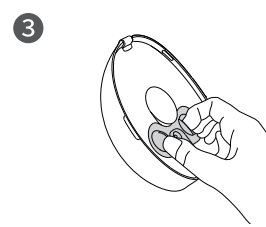


Pull out the diaphragm.

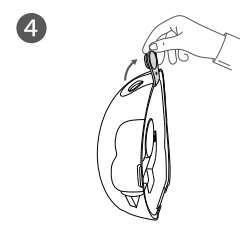

Open the container lid.

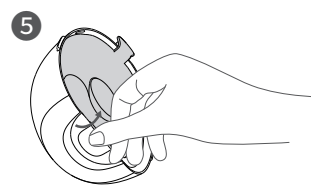

Remove the container back.

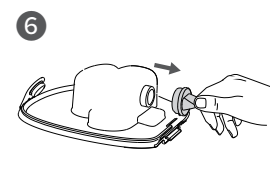

Remove the duckbill valve.

Cleaning and Sterilizing

- Before Initial Use:** Sterilize all washable parts by boiling for up to 5 minutes, stirring constantly to prevent the parts from sticking to the container and melting. Alternatively, steam-sterilize for a maximum of 5 minutes.
- Daily Cleaning:** Rinse all washable parts with warm, soapy water, then air dry. Sterilize if needed, using the same methods as the initial cleaning.
- Lifespan of washable parts:** The estimated life of the washable parts is 1-3 months, and actual life will vary with use and sterilization frequency. Replace them if you notice a decrease in suction and performance (particularly the duckbill valve), or if you have stored them too long when planning to pump for your next child.


- Do not use a microwave to sterilize.



- All washable parts (**except pump hub**) can be sterilized as follows.

- The pump hub is NOT washable. Always keep the pump hub dry and clean after each pumping session. If milk is spilled on the pump hub, wipe it clean with a dry cloth to ensure proper charging.**

Dishwasher

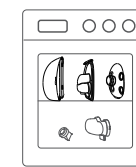

Cooking Pot



Baby Bottle Sterilizer




Anker Innovations Limited. All rights reserved, registered in the United States and other countries.

Anker Innovations Limited
Unit 56, 8th Floor, Tower 2, Admiralty Centre, 18 Harcourt Road, Hong Kong

Model: T600P
5100500xxxx_V1

Assembling

1

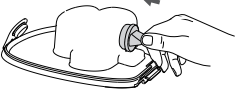


Attach the duckbill valve to the small hole inside the container. Flick it to ensure it does not fall out.

2




Press the container back to the container firmly.

3

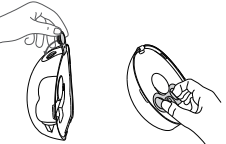


Fit the diaphragm into the contoured position at the back of the container, then press the lid firmly to ensure a secure seal.

4

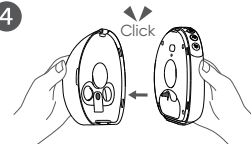


Align the peg on the pump hub with the hole in the diaphragm, and press the assembled container onto the hub until it clicks.

5

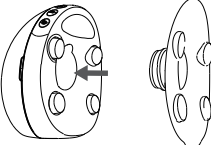

Press the flange onto the pump


6

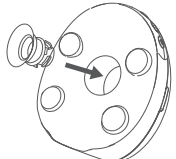

Use the Nipple Measurement Card to select an appropriate sizing insert, and place it in the flange tunnel.


Wearing

1

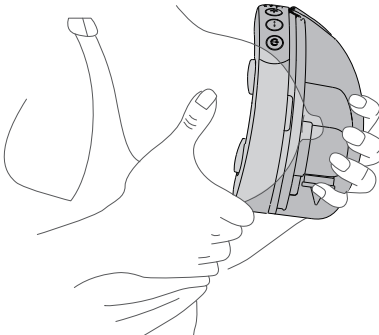

Slightly lean forward to make sure your nipple is centered in the flange tunnel. The breast pump should be firmly against your breast.

2

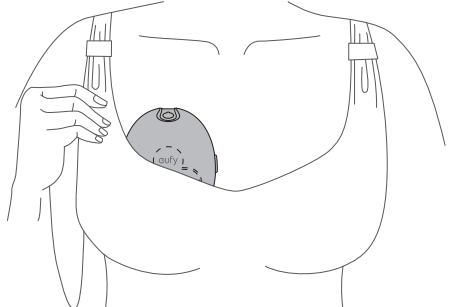

Use a full-coverage nursing bra with a flap to create a tight seal between your breasts and the flange.

 **Two criteria will help you determine if the breast pump is positioned correctly:**

- Whether you feel a noticeable suction that can pull out milk;
- Whether it is comfortable to wear.

If the wearing position doesn't meet either of these criteria, you need to adjust the position of the pump or change to a more supportive nursing bra.

Using



Button Controls		
	Press and hold for 3 seconds	Turn Breast Pump On / Off • The breast pump will turn off automatically after pumping for 30 minutes or pausing for 10 minutes.
	Press	Pause / Continue
	Press	Switch between Stimulation Mode and Expression Mode • Stimulation Mode: In this mode, the breast pump begins with a quick and short sucking pattern to get your milk to start flowing. The stimulation mode will last 2 minutes, then it will automatically switch to expression mode. There are 7 levels in stimulation mode. • Expression Mode: If breast milk comes out early, switch to expression mode. In this mode, the breast pump begins with a slow sucking pattern for milk expression, sucking more deeply and more slowly. There are 10 levels in expression mode.
	Press and hold for 2 seconds	Turn Nipple Tunnel Light On / Off
	Press	Turn Massage mode On / Off. There are 3 levels in massage mode.
	Press	Turn Heating On / Off • Adjust heating level via the app. • Automatically stops after 10 minutes of heating.

	Press	Adjust the suction/heating/massage level to the most comfortable level. Increase / Decrease the Suction/Heating/Massage Level
Any Button	Press	Check the remaining battery level.

Light & Display Guide



	Flashing White	Charging in Process
	Steady White	Fully Charged
	Steady Orange	Battery Level <20%
	Flashing Orange	Battery Level <10%
	Steady White	In Stimulation / Expression / Massage Mode
	Flashing White	The corresponding function is selected to switch stimulation / expression / massage levels.
	Steady Orange	In Heating Mode
	Flashing Orange	The function is selected to switch heating levels.
	Steady White	Levels
	Steady White	Pumping time
	Left / Right Indicator Flashing Alternately	Left / Right Side Pairing in the App
	Left / Right Indicator White	Left / Right Side Paired in the App

 **Can't feel suction or have zero or very little output?**

- Inspect the duckbill valve and diaphragm for damage, including cracks, chips, holes, or tears. If a part is damaged, discontinue use and purchase a replacement.
- Check that the assembly is correct and that the five pieces (duckbill valve, container back, flange tunnel, diaphragm, and pump hub) are perfectly sealed.
- Adjust the wearing position, make sure the nipple is at the center of the flange tunnel, the pump is sealed to your breast, and you can feel noticeable suction.
- Experiment with different flange or insert sizes.
- Try different suction and rhythms.

It takes a bit of time for your body to adapt to a new pump, so don't be discouraged if the results aren't good at the very beginning.


Transferring Milk

1



Pause or turn off the breast pump. Lean forward when taking out the pump to prevent leakage.


2

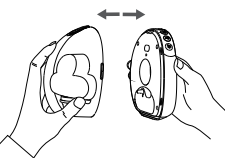

Remove the container from the pump hub.

 Do not remove any parts from the container before pouring out the milk to avoid leakage.

3

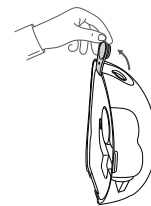

Open the container lid.

4



Tilt the container with the pouring spout facing down. Remove the diaphragm if needed to collect any trapped milk under it.

Downloading the App

Download and install the eufy app, then sign in or create an account. Follow in-app instructions to add this device and complete the setup.



Customer Service

 12-month limited warranty on pump hub
3-month limited warranty on washable parts

 (US/CA) +1 (800) 988 7973

 support@eufy.com (US/CA)
support.au.eufy.com (AU)

 @EufyOfficial @EufyOfficial @EufyOfficial

 Australian and New Zealand customers, please scan the QR code for the warranty information.