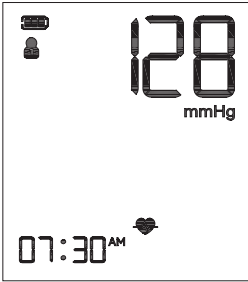
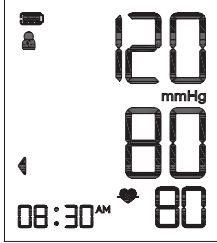
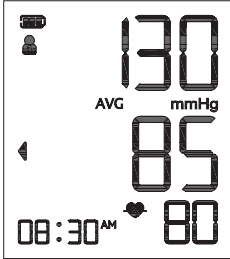
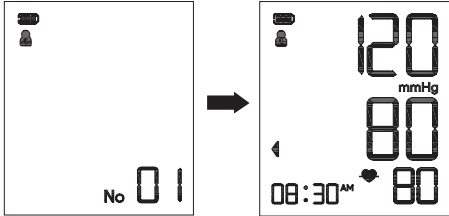
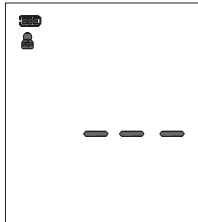


Unit Operation 23 3. Testing After cuff inflation, air will slowly subside as indicated by the corresponding cuff pressure value. A flashing "♥" will appear simultaneously on screen signaling heart beat detection.  Note: Remain relaxed during testing. Avoid speaking or moving body parts.	Unit Operation 24 4. Result Display The screen will display measurements for systolic and diastolic blood pressure with voice broadcast. An indicator representing the current measurement will appear next to the corresponding WHO Classification.  Note: Refer to Page 36~37 for detail WHO Blood Pressure Classification Information.
Unit Operation 25 Irregular Heartbeat Indicator If the monitor detects an irregular heart rhythm two or more times during the measuring process, the Irregular Heartbeat Symbol "♥" appears on screen along with measurement results. Irregular heartbeat rhythm is defined as rhythm that is either 25% slower or faster than the average rhythm detected while measuring systolic blood pressure and diastolic blood pressure. Consult your physician if the Irregular Heartbeat Symbol "♥" frequently appears with your test results. 5. Deleting/Storing Test Results User may delete their current test result due to unfavorable testing conditions or for any other reason. To delete the last test result, press the "SET" button after result is displayed. If result is not deleted, it will automatically store by date within the previously configured Memory Group. Note: Be sure the appropriate Memory Group selection is made prior to testing. If the number of tests surpasses the allotted 60 memories per group, the most recent tests will appear first, thus eliminating the oldest readings. Power Off The "OFF" button can be pressed to turn off the unit in any mode. The unit can turn off the power itself about 3 minutes no operation in any mode.	Unit Operation 26 Safety Precaution: If pressure in cuff becomes too extreme while testing, press the "OFF" button to turn power off. The cuff pressure will rapidly dissipate once the unit is off. Memory Check and Last 3 Tests Average With power off, press the "ON" button to activate screen display. After the unit performs a self diagnosis, the screen will display the average test results from the last 3 readings of the last group used. The "AVG" symbol will appear along with the corresponding WHO Blood Pressure Indicator. To check the average results from other groups, select the desired group first prior to activating the "ON" button in the off position (See "Select Memory Group" on page 16). 
Unit Operation 27 Press the "ON" button again, you may check past test results. Upon activating test results, you can press the "ON" button to scroll through all test results stored in memory. The LCD will display the last memory as NO: 01 reading. ("ON" button means turn forward.)  Note: Past test results will only be displayed from the most recently used memory group. To check past test results in other memory groups, you must Press corresponding button or select the desired group and then turn monitor off. (See Select Memory Group on Page 15.)	Unit Operation 28 Memory Deletion Memory for a selected group may be deleted while in Memory Check mode. Press and hold the "ON" button 3 second for approximately 3 seconds to delete all memory records from the selected group with voice broadcast "Memory Clear". And then transfer into testing mode. Press the "OFF" button to turn the unit off.  Note: Memory cannot be recovered once it has been deleted.