
XFT-2017

MiniAndy

CPR Manikin

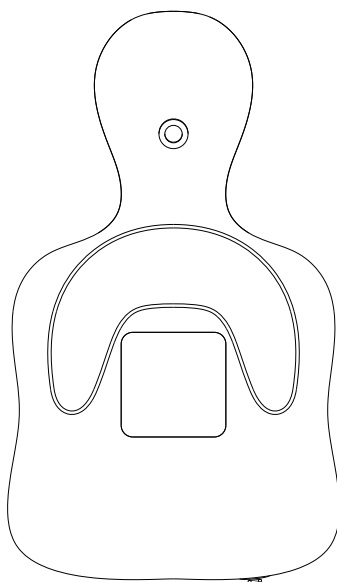Instruction Manual

FCC Notice:

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference. (2) This device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

WARNING: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

RF Exposure

The equipment has been evaluated to meet general RF exposure requirement. The device can be used in portable condition without restriction.

EN	Instruction Manual	01
CN	使用说明书	08

List of Symbols

	Date of manufacture
	Manufacturer
	Serial Number
	CE Marking
	The product complies with the applicable FCC requirements
	Fragile
	Keep upward
	Keep dry
	Don't roll while moving
	This product should not be disposed of together with other domestic waste, and the disposal of electronic waste should follow local regulations or government regulations
	Temperature limitation
	Humidity limitation
	Atmospheric pressure limitation

Product Overview

Intended Use

The XFT-2017 is an inflatable CPR training manikin designed for chest compression and mouth-to-mouth ventilation practice. It features real-time feedback on compression depth and rate via Bluetooth connection with the XCPR App on a smartphone or tablet.


Ideal for personal, institutional, and professional CPR training, the device supports accurate, data-driven CPR instruction and performance evaluation.

Package Contents

	CPR Manikin	1
	Compression Module	1
	Quick Start Guide	1
	Air Pump	(Optional)
	AAA Battery	(Optional)

CPR Manikin

User Interface

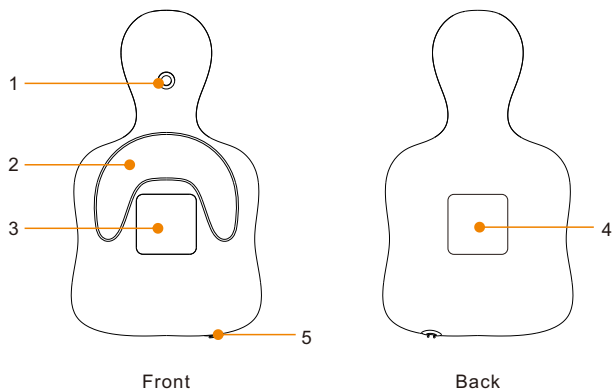
View of the CPR Manikin


1	Manikin mouth	The contact point for mouth-to-mouth resuscitation training.
2	Lung bag	Simulates chest rise during rescue breathing.
3	Compression Area	Target area for performing chest compressions.
4	Compression Module Area	Dedicated slot for installing/replacing Compression Module.
5	Inflation/ Deflation Port	Valve port for adjusting internal pressure with an air pump.

Compression Module

User Interface

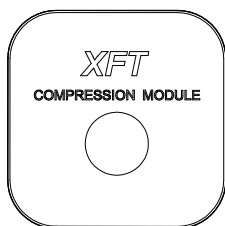Top

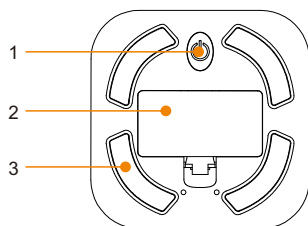Bottom

View of the Compression Module


1	Power On/Off	Turns the compression module on and off.
2	Battery Compartment	For installing and replacing batteries.
3	Stabilizing Feet	Secures and stabilizes the module position.

Operation Guide


Instructions

1. Open the compression module's battery cover, insert batteries, press the power button to turn it on.

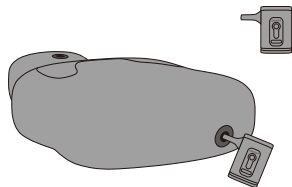2. Inflate the manikin using an air pump. Note: Avoid over-inflation.

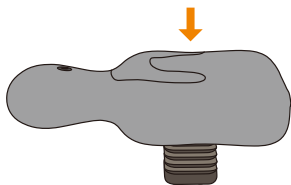
3. Position the manikin securely on the compression module.

4. Proceed directly to CPR training if no App is required.

5. For app-guided training: Install the XCPR App, pair the manikin with the App, then start CPR training.

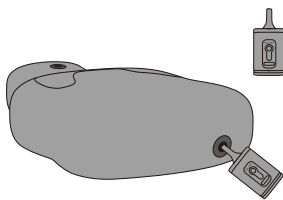6. Post-training: Deflate the manikin using an air pump, then store all components.

Maintenance

- Avoid Harsh Chemicals. Do not use alcohol, bleach, strong acids/alkalis, or corrosive cleaners, as they may damage the material.
- Do not use steel wool, hard brushes, or sharp objects to avoid scratching the silicone/PVC surface.
- Clean with mild soapy water (pH-neutral detergent) and a soft cloth. For stubborn stains, dampen with water and gently scrub.
- For disinfection, use diluted sodium hypochlorite (1:50 ratio) or <70% isopropyl alcohol wipes. Limit contact time to 1 minute, then wipe with water immediately. Air-dry naturally after cleaning.
- Avoid direct sunlight or high heat. Store away from heat sources and sharp objects.

Specification

Power: 3V (2×AAA batteries)

Manikin Size: 605×376×0.75 mm(±5mm, Uninflated state)

Compression Module Size: 90×90×120 mm (±5mm)

Shutdown current: ≤ 20μA

Working Current: ≤300mA

Storage Temperature.: -20°C~55°C

Operating Temperature.: 5°C~40°C

Humidity: ≤80% (non-condensing)

Pressure: 86~106 kPa

Working Condition

Temperature: 5°C~40°C

Humidity: ≤ 80% (Non-condensing)

Atmospheric pressure: 86kPa~106kPa

Storage Condition

Temperature: -20°C~55°C

Humidity: ≤ 93% (Non- condensing)

Atmospheric pressure: 70kPa~106kPa

After-Sales Service

Warranty Card

Product Name: _____ Model No.: _____

Purchase Date: _____ Product Serial No.: _____

Buyer's Information: _____

Distributor's Information: _____

Manufacturer: Shenzhen XFT Medical Limited

Add: Room 203, Building 1, Biomedicine Innovations Industrial Park,
#14 Jinhui Road, Pingshan District, Shenzhen, China

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