


Begin Your Session

- Customize the temperature and dose to your desired preferences
- Keep track of
 - How many draws you've had
 - How much you've inhaled per draw
 - The total amount inhaled
 - The average temperature of your session



Record How You Feel

- Rate how much you enjoyed the session
- Select any mental or physical side effects you experienced during or after your session
- Track how effective your session was at relieving the symptoms of your medical condition

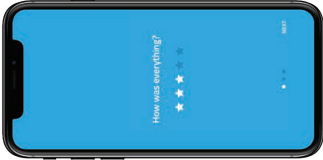

Cleaning Mouthpiece

- Step 1:**
- Detach the mouthpiece by pressing upwards on the front edge



- Step 2:**
- Remove the extended pathway insert



- Step 3:**
- Soak the mouthpiece and extended pathway / insert in soapy water or isopropyl alcohol for 15 minutes. Do not use harsh chemicals or heat (more than 15 minutes could affect exterior finish)

- Step 4:**
- After soaking for 15 minutes, use the scraping tool to remove any debris / residue from the mouthpiece and extended vapor pathway

Cleaning Heating Chamber

- Step 1:**
- Before cleaning, make sure heating chamber temperature indicator is GREEN before proceeding
 - If temperature indicator is OFF, please turn on RYAH. Wait until the indicator turns GREEN and ensure temperature indicator is GREEN before proceeding



Caution: RED temperature indicator means chamber is hot!

- Step 2:**
- Pull out the bottom chamber lid



- Step 3:**
- Detach mouthpiece by pressing upwards on the front edge



- Step 4:**
- Gently insert the pipe cleaners included in the RYAH cleaning kit through the bottom chamber and mouthpiece openings to remove any residue within the heating chamber.



WHAT'S IN THE BOX

- RYAH device
- Cleaning kit
- Manual