

SOBERLINK 7.0

Device Quick Start Guide



Table of Contents


Device Diagram.....	1
Accessories.....	2
Download the App.....	3
Pair the Device.....	3
Preparing to Submit a Test.....	4
Submitting a Test.....	5
Client Text Notifications.....	6
Device Operating Information.....	7
Troubleshooting.....	8
Intended Use.....	9
Device Recalibration.....	9
Disclaimers.....	9

Device Diagram



Accessories

Device Charger

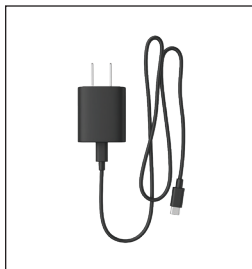

Insert one end of the cable into the charging base and the other end into the port on the bottom of the Soberlink Device.

Plug the charging base into an electrical outlet.

Mouthpiece



Insert the small, flanged end of the Mouthpiece into the Device.

Note: Always use a Soberlink Mouthpiece when submitting a test.

Device Case



To protect the Device when not in use, keep it stored in the Device Case.

The Mouthpiece can be stored inside the Device Case.

Download the App

1. Download the **Soberlink** app.





2. Activate your device by calling Soberlink at 714.975.7200

If you purchased Soberlink through your Treatment Provider, contact them to activate.

Pair New Device

Make sure Bluetooth is enabled on the iOS or Android device. Open the Soberlink App and follow the screen prompts.

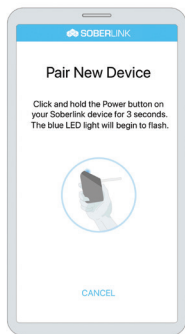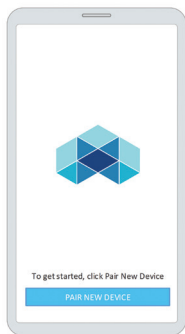

Preparing to Submit a Test

This document, and any other provided training material, should be fully understood by the user before submitting a test.

Before activating the Device and submitting a test, the Client is required to **set up a two-factor authenticated User account**. To activate call Soberlink at 714-975-7200 or contact your Treatment Provider.

Wait 20 minutes after eating, drinking, or smoking before submitting a test.

Failure to wait may result in inaccurate test results.

Rinsing out the mouth with water does not act as a substitute for the 20-minute waiting period.



⚠ To ensure proper device function keep it fully charged. To avoid damage or injury caused by the product only use accessories provided by the manufacturer.

Find an area free from alcohol, cleaners, solvent vapors, and smoke before submitting a test.

Use the Device within the operating temperature range (32°F -122°F).

Breathe normally before submitting a test.

Remove sunglasses, hats, or other objects covering the face.


Submitting a Test

Open the Soberlink App. The App should remain open for the duration of testing.

Press the Power button to turn on the Soberlink Device.


Follow the prompts on the App, and wait for the green light to flash before submitting a test.

When prompted to blow, look directly into the camera, breathe in deeply and deliver a firm, continuous breath for 4 seconds. The Device will click when the test is complete.

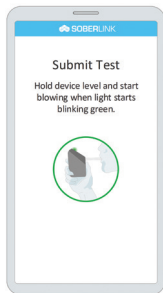

Note: Do not hold the mouthpiece or obstruct the camera.

The test results will take approximately 30 seconds to display on the App and send to the Sober Sky Web Portal. Once the test has successfully sent, close the App.

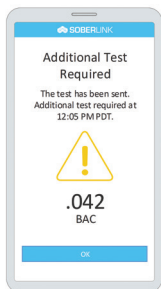

Note: After a positive BAC (Breath Alcohol Concentration) is received, wait 15 minutes before submitting another test.

Client Text Notifications

Once a schedule is enabled on the Sober Sky Web Portal, Soberlink's automated system will send text message reminders directly to the Client's phone.

Note: During the Client's setup of their two-factor authenticated User account, they will consent to receive SMS alerts & notifications to be Opted-In to notifications. These preferences can be managed in their User settings.

ALERTS AND NOTIFICATIONS

Email Opted-In	Edit
SMS Opted-In	Edit

49382



REMINDER TEXT: This is your 8:00 AM PST test reminder.

Reminder texts are sent before and during the test window.

49382



MISSED TEST: You have missed 1 or more tests. Send a test now.

Missed notification after test window closes

49382



RETEST REQUIRED: A test is required at 8:18 AM PST.

Notification of retesting requirement

Device Operating Information

Technical Specifications

Dimensions	4.75" x 2.50" x 1.25"
Weight	5.35 oz. (151.66 grams)
Sensor Technology	Deep lung, Dart fuel-cell sensor
Battery Type	Lithium-ion rechargeable battery, holds 10+ days
Power Supply	USB-C charger, Output: 5V 2A
Breath Sample Duration	4 seconds
Operating Temperature	32°F - 122°F (0°C - 50°C)
Relative Humidity	80% non-condensing up to 87 degrees F (31°C), decreasing linearly to 50% non-condensing at 104 degrees F (40 C°)
Min. Ambient Air Pressure	0.74 atm
Pollution Degree	2
Sensor Accuracy	+/- .005 BAC
Detection Range	.000 - .400% BAC
Calibration	Every 1,500 tests
Camera	Color, low light compensation
Test Data	BAC, time of test, and photo
Components	Breathalyzer, mouthpieces, charger, case
Manufactured	QSR 820, ISO 9001, ISO 13485, AS9100
Warranty	1 year limited

The user does not need to perform any tests at start up, the Device performs a self-test when it is powered on and prior to a test submission. The self-test is meant to verify the proper operation of the test storage memory, alcohol fuel cell, and imager. Failure of any of these components indicates that the Device is unable to properly record test data. In this case, an error code is displayed to the user and testing will not proceed.

Safety

There are no anticipated risks or discomforts utilizing the Device, however, if the Device is used in rapid succession, dizziness may occur. For correct and effective use of the Device, always consult this manual to familiarize yourself with the nature of the potential hazard and actions to be taken to avoid them. The Device should not be used to determine if it is safe for the user to drive. To ensure proper Device function, keep it fully charged.

The radiated output power of the device is below the FCC/ISED radio frequency exposure limits when hand-held during operation or when used at 2 cm touch distance to the body.

Troubleshooting

Predictable Use Problems

- Failure to use correct power supply
- Failure to re-charge the Device
- Failure to use provided Soberlink mouthpiece
- Failure to follow proper test preparation
- Failure to provide a sufficient breath sample (i.e. 4 continuous seconds)
- Failure to wait 20 minutes after eating, drinking or smoking before submitting a breath sample
- Failure to comply with the instructions prompted by the App
- Failure to wait 15 minutes between submitting breath samples
- Failure to activate the Device with the Sober Sky Web Portal

The table below lists a range of Device errors and is provided to help eliminate confusion and prevent downtime by supplying troubleshooting steps.

FOR FURTHER ASSISTANCE, OR TO RETURN A MALFUNCTIONING DEVICE, CONTACT SOBERLINK HEALTHCARE AT SUPPORT@SOBERLINK.COM OR 714-975-7200.

Error Code	Solution
4000 - 10999	Wait approximately 30 minutes and retest in an area with adequate cellular coverage. If error persists, disable cellular service and attempt to submit a test via Wi-Fi connection. Alternatively, pair with different iOS or Android device and attempt to submit a test. Return the Device if error persists.
11000 - 11999	The Device is unable to perform a test. Follow the prompts on the App screen and retest. Return the Device if error persists.
12000 - 12999	Wait for the Device to power off, then power on the Device and retest. Return the Device if error persists.
17000 - 17999	Ensure Bluetooth is enabled on the iOS or Android device and that the Soberlink Device is properly paired. If error persists re-pair the Bluetooth connection. If error persists, contact Soberlink Healthcare.

Maintenance and Warranty

If any equipment is used in a manner not specified by manufacturer, the protection provided by the equipment may be impaired.

DO NOT use alcohol or other chemical cleaners to clean the Device, only water and a soft cloth.

DO NOT submerge the Device in water.

Repairs of the Device may be performed only by Soberlink Healthcare or an authorized service technician. Only original parts may be used in the Device. For Warranty information, review the Soberlink General Terms and Conditions found at www.soberlink.com/general-tc.

Intended Use

The Soberlink 7.0 Device is intended to quantitatively measure alcohol in human breath. Measurements obtained by this Device are used in the diagnosis of alcohol intoxication. For Prescription use and OTC use.

DEVICE RECALIBRATION IS REQUIRED AFTER 1,500 TESTS. FAILURE TO DO SO MAY CAUSE INACCURATE TEST RESULTS.

Users will receive text message notifications when the Device is nearing recalibration.

The Device will need to be returned to the manufacturer for recalibration.

Contact Soberlink Healthcare at support@soberlink.com, or 714-975-7200 to request a Return Authorization number.

Disclaimers

- Do not use this Device to determine if it is safe for the user, or anyone else, to operate machinery or a motor vehicle.
- Do not drink and drive. Always use a non-drinking designated driver after consuming alcohol.
- Do not attempt to open the Device, or perform any maintenance on the Device.
- Do not blow any food, drink, or smoke into the Device.
- Do not drop the Device.
- Do not share mouthpieces.

NOTICE: This device complies with Part 15 of the FCC Rules [and contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS standard[s]].

Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

L'émetteur/récepteur exempt de licence contenu dans le présent appareil est conforme aux CNR d'Innovation, Sciences et Développement économique Canada applicables aux appareils radio exempts de licence.

L'exploitation est autorisée aux deux conditions suivantes:

1. L'appareil ne doit pas produire de brouillage, et
2. L'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

Changes or modification made to this equipment not performed by Soberlink may void the FCC authorization to operate this equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help

Contact Us

If you have any questions, please contact us at support@soberlink.com or 844-975-7200.

Soberlink Healthcare LLC
16787 Beach Blvd. #211
Huntington Beach, CA 92647

Copyright 2026 Soberlink. All rights reserved.

QSG-SL7-26-002 | PN755-00009 S7