



Solve the problem of ear cleaning for the entire family

Open bebird Note3

IOS
Search for "bebird" in App Store

Android
Search for "bebird" in Google Play




Scan QR code to download APP



Enter APP
Select bebird Note3

Open the cover
Power on automatically

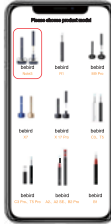



Open the cover
Power on automatically

Blue lighting
flashing

Connect Product WIFI

Connect product
WIFI "bebird-Note3-xxxx"

Back to bebird APP
Click on "Start" to use it

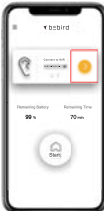Click on the tweezers sign to
click between tweezers'
opening and tweezers clamping

Miniature robotic arm
Open the visual tweezers
Clamp ear wax





Accessory replacement

Cover the tweezers with the ear spoon
Fasten the ear spoon

Cover the tweezers with
the ear spoon
Fasten the ear spoon

The ear spoon
must be assembled in
place before using
to prevent falling off.

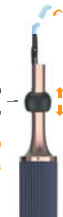

Ear canal protection ring
Move up and down to
limit insertion depth

Please put the ear canal
protection ring on the rod
when you use it lying on
your back

Storage base

Spin the cover of the base to open it
Swap accessories

Spin
Open the storage base



Accessory

Classic Q-elastic Ear Spoon x 4
Smooth use with no harm to the ear
A set separated for each family member



Microdots Cover x 1
Cover the ear plating rod in the base cover
It turns into an electronic magnifier



Plug the ear plating rod into the base cover
It turns into an electronic magnifier

Special Cleaning Ear
Spoon x 4
Clean the earwax
Deep dig into the inner ear for hidden corners



Ear Canal Protection Plug x 1
Avoid accidental injury to the
vulnerable ear canal



Curved-necked Annular
Ear Spoon x 4
Easy to scraping
Professional in removing stubborn earwax



Alcohol Cotton Swab x 5
Softening Cotton Swab x 5
Clean the ear canal
Softens stubborn earwax



Be careful not to let the
cotton swab touch the ear
drum and the eardrum

Mini Brush x 1
Clean hair, tweezers and ear scoops



Type-C cable x 1
Fast and convenient charging



Hazardous substance

Part name	Hazardous material					
	(Pb)	(Hg)	(Cd)	(Cr(VI))	(PBB)	(PBDE)
Printed Circuit Components	X	O	O	O	O	O
Charging Base	X	O	O	O	O	O
Case	O	O	O	O	O	O
Plastic Rubber Accessory	O	O	O	O	O	O

This form is formulated in accordance with SJ/T 11364
 O : represent the hazardous material is below 0.01(0.0025%~0.01)
 X : represent the hazardous material is over 0.01(0.0025%~0.01)

Troubleshooting

Phenomenon	Solution
Start up and no response	Please try several more times according to the correct boot-up method. Please use the attached charging base for charging. The charging time is about 1h.
Camera does not turn on	Shut down and restart the ear cleaning rod.
APP fails quickly	Delete the existing APP and re-scan the QR code for installation.
The power indicator of the ear cleaning rod flashes (Max light flashing rapidly)	Charge the ear cleaning rod with the charging base from the original packing.
WiFi connection is unstable, sometimes it is disconnected in seconds	Restart ear cleaning rod, reconnect bebird APP.
The product failed to display the real-time image after successfully connecting the bebird APP	Disconnect the power from WiFi, restart the product and reconnect the WiFi.

If the above operation will still be to solve the problem, please contact Customer Service for further assistance.

Attention

- When cleaning lens, please use it carefully with professional optical cotton lens.
- Make sure the ear spoon is well installed before using, to prevent the ear spoon from falling off.
- Please don't use it while you are moving, and protect from being hit by others.
- It is forbidden for children to use this product done in order to avoid accidental harm.
- Please do not put the product in liquid to avoid damage.
- Product built-in rechargeable lithium battery, so it must be charged once a month if it is not used frequently.
- Avoid product exposure to sunshine directly, especially the accessories, in order to avoid softening.
- During using, the temperature of the product will rise slightly (Data from bebird lab, the maximum temperature of the machine doesn't exceed about 45°C). Please rest assured to use.

Product parameter

Brand	bebird
Product Name	bebird smart visual ear cleaning rod
Model	bebird nose3
Net weight	128g
Network standard	IEEE 802.11b/g/n
Antenna	Wire antenna
Working frequency	2.4GHz
Image transmission rate	20fps
Image sensor	CMOS
Working temperature	+10 ~ 50℃
Battery	300mAh Lithium battery
Battery life	About 70 minutes
Charging time	1h
Input parameters	5V/220-6A
Working voltage	4.2V
Rated power	1.5W
Lens diameter	3.7mm
Optimum focal length	1.5 ~ 2mm
Pixel	1000W
Gyroscope	6 axis

bebird smart visual ear cleaning rod
 Model: bebird Nose3
 210-153770
 FCC ID:2AUJ3-BEBIRDNOSE3

Warranty card

User information	
User name:	Telephone number:
Purchase date:	
Seller:	Telephone number:

Solve the problem of ear cleaning for the entire family
 NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installations.
 This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.
 However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 -Reorient or relocate the receiving antenna.
 -Increase the separation between the equipment and receiver.
 -Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 -Consult the dealer or an experienced radio/TV technician for help.
 Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
 This device complies with Part 15 of the FCC Rules.
 Operation is subject to the following two conditions:
 (1) this device may not cause harmful interference, and
 (2) this device must accept any interference received, including interference that may cause undesired operation.

Solve the problem of ear cleaning for the family using the following steps:

1. **State the problem** has been tested and tested to comply with the limits of the FCC's Part 15 Class B digital device limits. FCC Rules. These limits are designed to provide reasonable protection against harmful interference. This equipment has been tested and found to comply with these limits. This equipment generates uses and can radiate radio frequency energy and, if not installed and used properly, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, you can try the following measures:

- (1) Reorient or relocate the receiving antenna.
- (2) Increase the separation between the equipment and receiver.
- (3) Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- (4) Consult the dealer or an experienced radio/TV technician for help.

2. **Identify the cause of the problem** The manufacturer is not responsible for compliance with civil law's authority to generate interference.

3. **Identify the solution** This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.