SNORE CIRCLE

Manufacturer: VFLY Electronics Co., Ltd.
Address: Unit 2103-2104, Changhong Sci-tech Building, Science
and Technology Park, Nanshan District, Shenzhen, China
Service email: kf@vfly.cn Service hotline: 400-8058192

Smart Electronic Muscle Stimulator
YA4600 User Manual

Hereby , we VVFLY Electronics Co.,Ltd.
declares that this device is in compliance with the essential
requirements and other relevant provisions of Directive 2014/53/EU.
This product can be used across EU member states. Adapter shall be
installed near the equipment and shall be easily accessible. The
operating temperature of the device between -10°C to 50°C.

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
Contents

1	Smart Electronic Muscle Stimulator (YA4600)	05
2	Powering On/Off	06
3	Battery and Power	07
4	Downloading YZFTools	08
5	Device Binding	09
6	Data Synchronization	09
7	Usage Guidelines	10
8	Personalized App Settings	10
9	Wearing Instructions	11
10	FAQs	13
11	Precautions	15
12	Declaration of Hazardous Substances or Elements	17
13	After-sales Service	18
14	FCC	18

1.Smart Electronic Muscle Stimulator(YA4600)

As an intelligent wearable electronic device, this smart throat massager is designed to help you enjoy quality sleep and embrace a healthier lifestyle. We hope our product assists you in achieving better, healthier sleep!

Product parameters

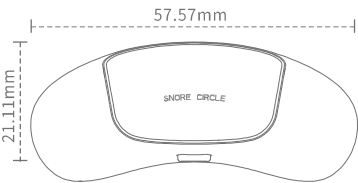Product name: Smart Electronic Muscle Stimulator
Model: YA4600
Weight: 12.6g
Modulation wave frequency:15KHz±10%
Run time: About 15 hours

Brand: Snore Circle
Dimensions: 57.57*21.11*14.78mm
Intermediate frequency: 3-8kHz+10%
Battery capacity:3.8V/130mAh

Product components: Main component x1, Charging Dockx1, Wearing Guide x1, Charging Cable x1, Electrode Pad x10, Product instructionx1

Recommended for: Those without contraindications, mild to moderate apnea ,that is AHI <30;


Unsuitable users: People with the following conditions: pacemakers, implantable defibrillators or any other implanted metal or electronic devices; acute suppurative inflammation, a tendency to bleed, or malignant tumors in the submandibular region; severe cardiovascular or cerebrovascular problems; ventilatory disorders caused by nasal disease; an allergy to electrodes or currents; skin ulcers on the site of adhesion; pregnant women or those in labor; and children under 16.

2.Powering On/Off

Power on

To turn the device on, simply press the power key for 3 seconds.

Once powered on, the power indicator will turn on for 10seconds and the Bluetooth indicator will flash for 30seconds before automatically turning off.

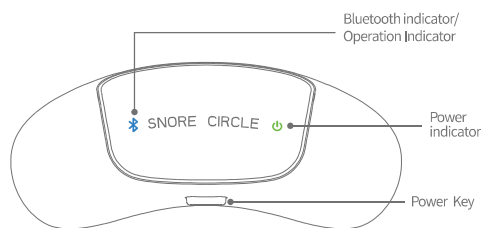Power off

To turn the device off, simply press the power key for 3 seconds.


The power and Bluetooth indicators will go off once the device is successfully powered off.

- 1) The power indicator will be on for less than 10 seconds after slightly press power key during or after power-on. In order to save the power, it will go off automatically.
- 2) If you want to confirm whether the device is working, slightly press the power key to check whether the power indicator is on. It indicates that the device is in working state when it is on.

3.Battery and Power

- 1) Checking the battery level Once the device is powered on or the power key is pressed, you can determine the battery level by checking the color of the power indicator. A yellow light indicates the device battery is low and requires charging.

A yellow light indicates the device has sufficient power and does not require charging.

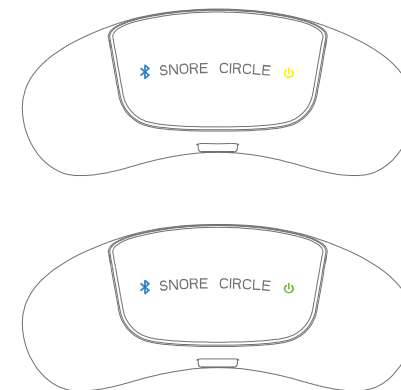To check the battery level in detail, connect the device to the "YZFTools" App. After binding the device to the "YZFTools" App, navigate to [Basic Settings] in the device control interface to view the battery status.

(2) This device uses a rechargeable lithium polymer battery.

(3) It takes two hours to fully charge the battery.

(4) To charge the device, insert the USB charging cable into the charging dock and place the device onto the dock. A flashing power light indicates the device is charging, while a steady light indicates charging is complete.

4. Downloading YZFTools

How to download:

(1) iPhone : Search YZFTools in the App Store to download and install the app.

(2) Android : Search YZFTools in Google Play to download and install the app.

(3) Scan the QR code below.

It is recommended that the device shall be paired with the app.

(1) Only iPhones with iOS 9.0 or above or phones with Android 6.0 or above are supported.

(2) Some Android phones require users to enable the location permission for the app. The Bluetooth permission must be enabled on your phone to search and pair a device.

(3) The device only needs to be paired once, unless it is later unpaired.

5. Device Binding

(1) Enable Bluetooth on your phone and open the "YZFTools" App.

(2) When the device is powered on (or shortly after), briefly press the power button, the Bluetooth indicator will flash for 30 seconds, during which binding is possible. After 30 seconds, the Bluetooth indicator turns off, and the device's Bluetooth will be disabled. Briefly press the power button again to re-enable Bluetooth.

(3) Place the device close to your phone, then tap [Add Device] on the home screen or the [+] button in the upper-right corner to initiate device binding. If the Bluetooth indicator changes from flashing to a steady light, it indicates the connection is successful.

(4) If the device is not detected, please ensure Bluetooth remains enabled, keep the device near the phone, and retry the search, following the App's instructions for device binding.

6. Data Synchronization

If the device is not paired with or disconnects from the phone during use, recorded data will be temporarily stored in the device. After use, we recommend syncing the data to the App for complete statistics and trend analysis.

Method 1: Quick Sync

(1) Open the "YZFTools" App and go to the home screen.

(2) Place the App close to the device and gently press the power button.

(3) Tap the device card to access the device control interface. A synchronization progress window will appear. Please wait patiently until the process completes, after which the window will automatically close.

Method 2: Manual Sync

(1) On the App home screen, tap the device card (e.g., YA4600).

- (2) Select [Basic Settings], then tap [Sync Data] button.
- (3) Wait for the progress window to complete (avoid switching apps or locking the screen).

7. Usage Guidelines

- (1) After powering on, the device must connect to the phone via Bluetooth. Enable Bluetooth, open the "YZFTools" App, and tap the paired device card to enter the sleep monitoring interface. During use, it is recommended to keep the device and phone within ≤ 1 meter, with no obstructions or interference between them.
- (2) Wear the device under the jaw, ensuring the electrodes maintain close contact with the skin. If the intervention intensity feels insufficient, slightly adjust the electrodes toward the Adam's apple (avoid forcibly increasing intervention intensity to prevent discomfort).
- (3) The App must remain active in the foreground (sleep monitoring screen) during use. Screen-off operation is allowed (reduces phone power consumption without affecting connectivity).
- (4) If the App runs in the background, the system may forcefully terminate the process, switching the device to algorithm compatible with offline mode (AI real-time analysis is unavailable in offline mode).
- (5) If the Bluetooth connection is interrupted during use, the device will automatically activate the algorithm compatible with offline mode. To restore the online AI recognition algorithm, simply tap the corresponding device card on the App's home screen, and the device will immediately switch to the online AI algorithm mode.

8. Personalized App Settings

After device binding, users may customize the device's intervention intensity, intervention delay based on individual needs for optimal comfort.

- (1) Ensure Bluetooth is enabled and the device is powered on.

- (2) Open the "YZFTools" App and place it near the device.
- (3) On the App's home screen, tap the device card to enter the device control interface.
- (4) Briefly press the device's power button, the Bluetooth indicator will flash.
- (5) Once the device interface shows "Connected", tap [Basic Settings] button to adjust default settings.
- (6) Operation indicator: The device's operation indicator is enabled by default. If the light affects your sleep, you may manually turn off this feature.
- (7) The device offers three intervention intensity levels: Low (gentle adaptation) as the default setting, Medium (balanced efficacy), and High (rapid intervention), recommended for users with greater tolerance. We advise first-time users to begin with the Low setting and gradually adjust to higher intensities.
- (8) Intervention delay: For users who are easily disturbed during the falling asleep phase, you may adjust the intervention delay parameter based on personal sleep habits.

9. Wearing Instructions

- (1) It is recommended to shave any beard or hair under the chin and clean the site of adhesion with soapy water before use and keep the area clean and dry to ensure maximum adhesion. Avoid overshaving as this may cause skin irritation. Please refer to the image below for the correct positioning:

- (2) Tear off the protective film on the electrode pad and ensure the groove is nearest your throat. lift your chin slightly and use both hands to stick the electrode pad to the correct position under your chin. This helps to prevent poor adhesion of the electrode pad as a result of movement to the soft neck tissue.
- (3) To ensure the pad is firmly adhered to the skin, press the entire electrode pad firmly with your palm for about 30 seconds after application to ensure that it is in full contact with the skin.

(4) Power on the device and attach it to the electrode pad, as shown below:

Note: If you are prone to sweating during sleep or find the electrode pad falls off easily, you may use medical tape that meets the requirements of local policies and regulations to hold the device in place.

Note: Scan the QR code below for videos.

10. FAQs

(1) The device wakes me up before I properly fall asleep. What should I do?

It is recommended to pair the device with the YZFTools app and use the Intervention Delay to solve the problem.

(2) It hurts slightly when I use the device. What should I do?

It is recommended to pair your device with the YZFTools app and reduce the Maximum intervention intensity in Settings. You can also experience the Maximum intervention intensity when pairing and adjust the Maximum intervention intensity according to the highest level you can tolerate. For detailed steps, please view Section 8 on personalized settings.

(3) Now that the device can work without connecting to the phone, why use?

The device has two operation modes: one is based on AI algorithm via mobile hashrate to identify suspicious breathing events as well as snoring; the other is based on matching algorithm via chip hashrate inside the device for identification. The AI algorithm is higher in recognition rate and accuracy for suspicious breathing events compared to the offline model. Meanwhile, the AI algorithm can be constantly upgraded.

(4) Does the app work on all phones?

The app supports phones with iOS 9.0 or above or Android 6.0 or above. Currently, 99% of the phones on the market can be used.

(5) Do I need to sync data every day?

If the device is not connected to the phone or is disconnected while in use, the statistics will be temporarily stored in the device. When the device is in operation, it is recommended to synchronize the data to the app for the complete statistic data as well as the trends.

(6)will the power indicator stay on after the device is powered on?

No.When the device is powered on, the power indicator will turn on for about 10 secondes ,after which it will automatically turn off, Gently press the power key to turn it on again for another 10seconds.

(7) How do I know if the device is working?

Gently press the power key, the device is working if the power and Bluetooth indicators turn on. If the do not, please power on or charge the device. The operation indicator(the Bluetooth indicator) will light when the device intervenes.

(8)Is it possible to turn off the flashing light when the device is in use?

Yes, it is. when snoring is detected during use, the device will automatically intervene during the coming respiratory cycle. The light on the device will flash during intervention. After pairing you device with the YZFTools app,you can go to Settings to disable the Operation indicator.

(9) How does the device monitors breathing?

The device is worn under the chin to collect data on skin movements when the user breathes. It will identify whether the user is snoring or there is apnea. Therefore, it is necessary to stick the electrode pads firmly and wear the device well.

(10) How to diagnose whether I have apnea and its level in mild or moderate?

To get the most accurate diagnosis, please consult with doctor for sleep monitoring for a report to know your accurate AHI condition. Or use our breathing monitoring products for initial screening. To order, visit our official website: www.snorecircle.com/shop. In addition, there are about half of the people who snore loudly suffering from sleep apnea. In such cases, snoring is a sign indicating the airway is getting blocked and narrowed, normally with the symptom of drowsy during the day, and awakened by

suffocation during the sleep at night, If you can't figure out whether you snore or not, maybe you can ask the person sleeping next to you, He or she will tell you about the snoring condition and whether it is accompanied by apnea.

(11) Is the strength level of the device intervention constant?

When the device is intervening, it will start from the minimum level, and gradually increase if the intervention is unsuccessful, until it works or the Maximum intervention intensity is reached. If no more suspicious breathing events for a while, the intervention level will be automatically reduced.

11.Precautions

(1) Use and store the device in a normal temperature range (-10°C to +50°C).

(2) Keep the device away from sources of fire or high temperatures (e.g. electric heaters, microwaves, ovens, etc.).

(3) Keep the device away from high impact or strong vibration to avoid damage.

(4) Do not use chemicals, detergents, etc. to clean the device. Keep the device dry; do not wash it.

(5) Do not attempt to remove or replace the built-in lithium polymer battery on your own so as to avoid damage.

(6) Do not throw the device into fire to prevent the possibility of an explosion.

(7) This device is not a toy. Please keep it out of the reach of children to prevent accidental swallowing or other dangerous incidents.

(8) The electrode pads contained in this product are disposable. Do not reuse them.

(9) After use of this product, the electrode pads should be processed in accordance with the local laws and regulations, and must not be disposed of at will.

(10) Discontinue use immediately if you experience any discomfort during use of this product.

(11) This product must not be used in conjunction with other medical device products.

(12) Please do not attach the product to the eyes, the carotid sinus or elsewhere near the cardiac shadow area.

12. Declaration of Hazardous Substances or Elements

In accordance with the relevant regulations and requirements of the Administrative Measures on the Control of Pollution Caused by Electronic Information Products of the People's Republic of China, the following table includes the component names of the electronic product and the toxic and hazardous substances contained within.

Component Name	Toxic and Hazardous Substances and Their Concentration					
	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent Chromium (Cr6+)	IPolybrominated Biphenyls (PBB)	Palybsomiated Diphenyl Ethers (PBDE)
Housing	○	○	○	○	○	○
Keys	○	○	○	○	○	○
User Manual	○	○	○	○	○	○
Standard Battery	○	○	○	○	○	○
Electrodes	○	○	○	○	○	○
Package box	○	○	○	○	○	○
Metallic Components	○	○	○	○	○	○

1. "O" indicates that the concentration of the hazardous substance in all homogeneous materials in the component is below the limit requirement of the SJ/T11363-2006 standard.

2. "x" indicates that the content of the toxic and hazardous substance in at least one homogeneous material of the component exceeds the limit requirement of the SJ/T11363-2006 standard.

Environmental Friendly Use Period

ComponentSmart	EnvironmentaEFriendly Use Period
Electronic Throat Massager	

*The above use period is based on normal operating conditions (temperature/humidity)
Product Standard:Q/YZFD 008-2021

13.After-sales Service

- (1) This product is covered under a one-year warranty from the date of purchase.
- (2) Artificial damage is not covered by warranty.
- (3) It is recommended to maintain a two-week trial period, with the first week as the adjustment period.
- (4) The electrode pad is disposable, it needs to replace with a new piece per night. Order at VVFLY official website www.snorecircle.com/shop.

FCC Caution

§15.19 Labelling requirements. This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, includinginterference that may cause undesired operation.

§15.21 Information to user.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

§15.105 Information to the user.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment of and on, the user is encouraged to try to correct the interference by one or more of the following measueres:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

*RF warning for Portable device:

The device has been evaluated to meet general RF exposure requirement.
The device can be used in portable exposure condition without restriction.

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