
WUUK

WUUK Smart Sock

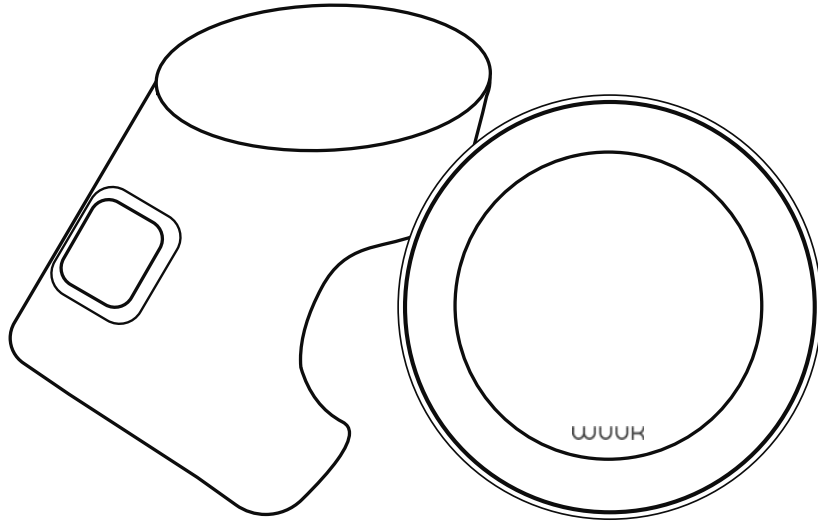WUUK

English	02
Deutsch	17
Español	29
Français	41
Italiano	53

TABLE OF CONTENTS


WARNINGS / IMPORTANT SAFETY INSTRUCTIONS	03
WHAT'S INCLUDED	04
PRODUCT OVERVIEW	05
HOW DOES THE PRODUCT WORK	07
SETTING UP THE SYSTEM	08
CUSTOMER SERVICE	13
SAFETY STATEMENT	14

WARNINGS / IMPORTANT SAFETY INSTRUCTIONS


Remember:

- This product is meant for use with healthy children up to 18 months of age. If your children have a pre-existing medical condition, consult your doctor before use.
- This product is an aid, it is not a substitute or should not be substituted for care and oversight of an adult or consultation with medical professionals.
- This product is designed and intended for tracking a baby's sleep patterns, sleep training, and general wellness purposes. False positive or false negative readings about your baby's heart rate are a potential risk of this product.
- Do not use this product as a diagnostic tool.
- Do not use this product as an excuse for unsafe sleep habits.


WHAT'S INCLUDED

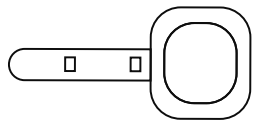
Sock Sensor

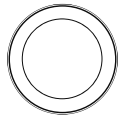Smart Sock
Base Station

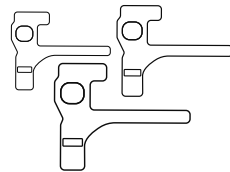Smart Sock
Wrap * 3 Sizes

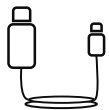Base Station
Type-C Power
Cable (1M)


Power
Adapter (5V 1.5A)

User Manual

PRODUCT OVERVIEW

Smart Sock

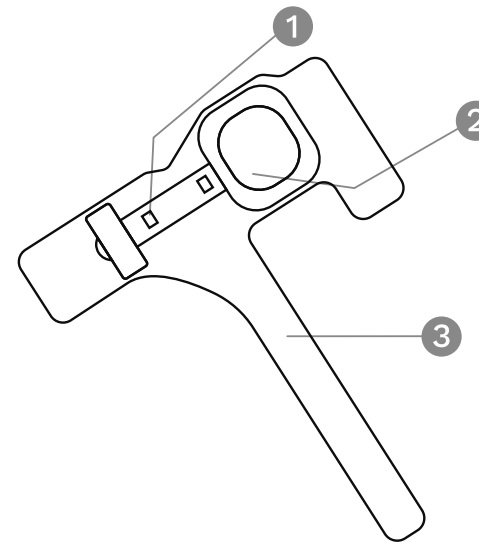
1 Indicator light

2 Sock Sensor


3 Sock Wrap

Base Station

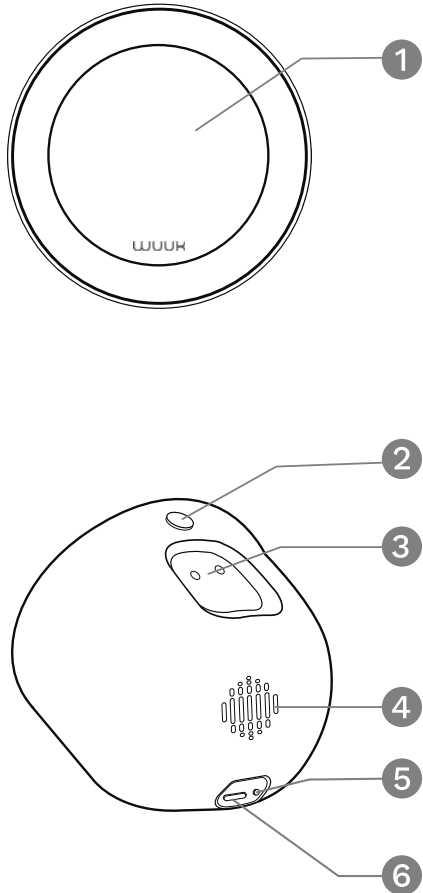- 1 Screen
- 2 Snooze Button
- 3 Magnetic Charging Pin
- 4 Speaker
- 5 Sync Button
- 6 Power Port

HOW DOES THE PRODUCT WORK

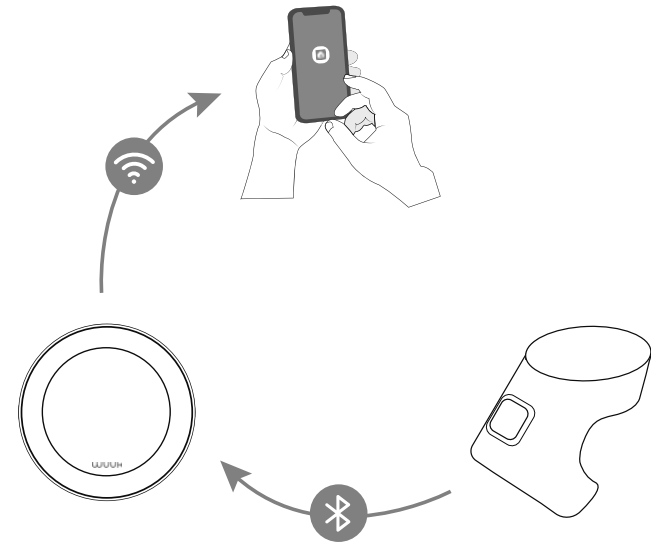
The Smart Sock sensor collects data and transmits it to the base station via Bluetooth and the base station will send information to WUUK Baby server via Wi-Fi.

The Sock Base Station connects to the WUUK Baby App to provide a complete sleep improvement solution for both the baby and the entire family.


SETTING UP THE SYSTEM

Power On the Sock Base Station

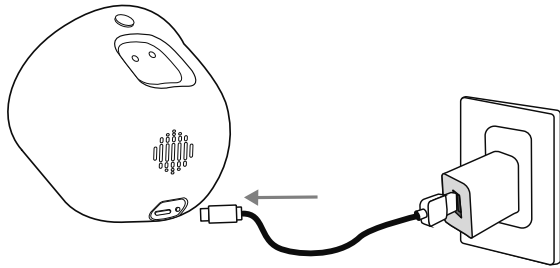Connect the Sock Base Station to a power outlet

Charge the Smart Sock

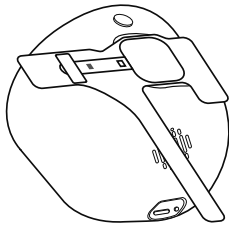Place the Smart Sock on the Sock Base. The sensor's indicator light will flash red.

Connect the Sock Base Station to the WUUK Baby App

Download the WUUK Baby app from the App Store (iOS devices) or Google Play Store (Android devices).


Go to the app homepage. Tap Add Device, and add the WUUK Smart Sock to your WUUK Baby app.


Pick a size that fits well

Not too tight, not too loose, with enough room for your baby's skin to breathe. This ensures comfort and improves the accuracy of the sensors.

Align and attach the sensor

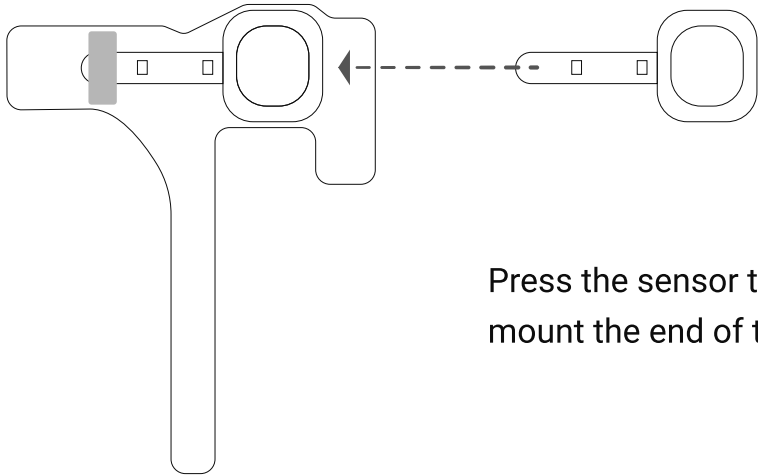Press the sensor to the sock and mount the end of the sensor.


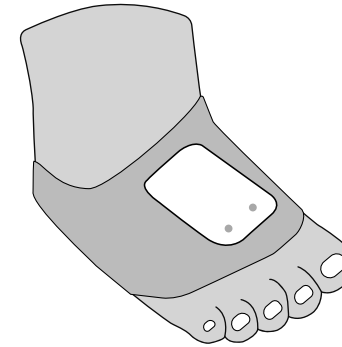
Right Foot

Dome on top of foot

Left Foot


Dome on bottom of foot

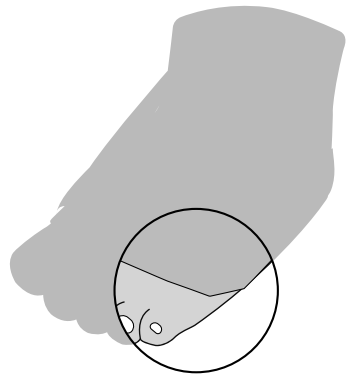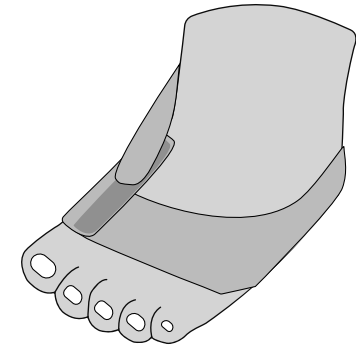
Interchangeable Sock

The sock is reversible and can be worn on either foot. For proper placement on either foot, place the notch on the outside of the foot behind the pinky toe.


Fit the Sock

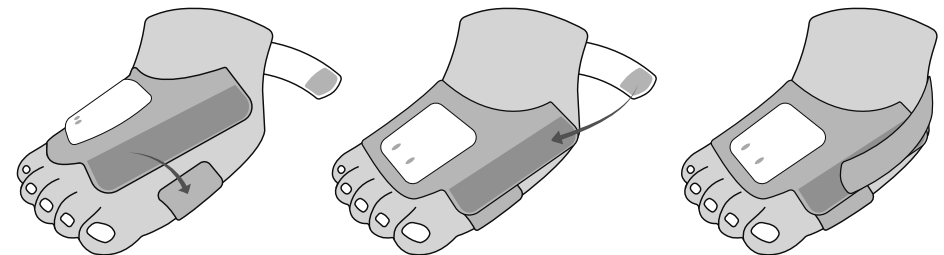Understand the Snooze Button

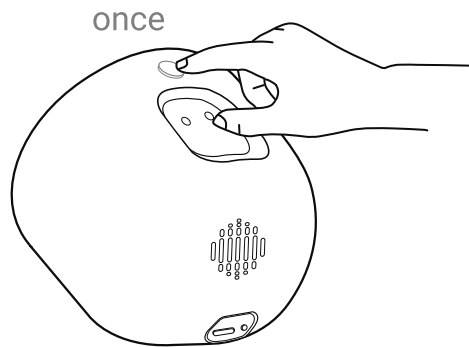Snooze the Alarms for 10 minutes: Press the Snooze Button once.

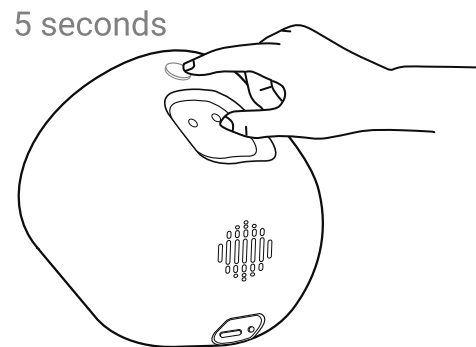Turn off the Screen: Press and hold the Snooze Button for 5 seconds.

CUSTOMER SERVICE

For troubleshooting, check the WUUK BABY app > Account > Help Center for more references and live chat with Customer Support agents if needed.

Warranty: 12-month Limited Warranty

Call Us: United States +1 (415) 881-5715

Email Us: hello@wuuklabs.com

WUUK Labs Corp. All rights reserved.

WUUK and WUUK Logo are trademarks of WUUK Labs Corp., registered in the United States and other countries. All other trademarks are the property of their respective owners.

SAFETY STATEMENT

FCC Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator and your body.

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This equipment complies with Canada radiation exposure limits set forth for an uncontrolled environment.

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This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:
This device may not cause interference.
This device must accept any interference, including interference that may cause undesired operation of the device.

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**Product Model: BAB10
IC: 29239-BAB10**

**PMN: WUUK Smart Sock Station
HVIN: BAB10**

Declaration of Conformity

- When charging, please place the device in an environment that has a normal room temperature and good ventilation.
- It is recommended to charge the device in an environment with a temperature that ranges from 5°C to 25°C.
- Please ensure to use only the charger offered by the manufacturer. Using an unauthorized charger may cause danger and violate the authorization of the device and the guarantee article.
- This charger is for indoor use only.

IC RF Statement

When using the station, maintain a distance of 20 cm from the body to ensure compliance with RF exposure requirements.

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