USER INSTRUCTIONS

Widex SoundConnect

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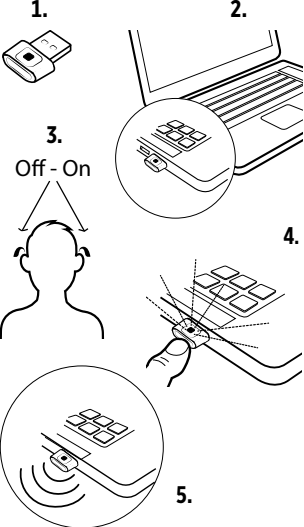
EU Directives
Directive 2014/53/EU
WSAUD A/S hereby declares that this SoundConnect model is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU. The device contains radio transmitters operating at: 2.4 GHz, 0.9 mW EIRP. A copy of the Declaration of Conformity can be found at: <https://global.widex.com/doc>

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ENGLISH
Read these instructions carefully before you start using the device.
Intended use
The device is intended to transmit audio signals from a Personal Computer directly to supported hearing aids.
NOTE:
The device uses a USB A connector. If your computer is without an USB A configuration, you may use a converter.
HOW TO USE
Connect the device to your computer. (Illu. 1 and 2)
Restart your hearing aids. The hearing aids will then be in pairing mode for 3 minutes. (Illu. 3)
Press and hold the button on the device until the light indicator is flashing blue, then release. (Illu. 4) The SoundConnect is now in pairing mode.

Information regarding disposal
Do not dispose of hearing aids, hearing aid accessories and batteries along with ordinary household waste. Hearing aids, batteries and hearing aid accessories should be disposed of at sites intended for electrical waste and electronic equipment, or given to your hearing care professional for safe disposal.

Símbolo	Título/Descripción
	Manufacturer: The product is produced by the manufacturer whose name and address are stated next to the symbol. If appropriate, the date of manufacture may also be stated.
	Catalogue number: The product's catalogue (item) number.
	Consult instructions for use The user instructions contain important cautionary information (warnings/precautions) and must be read before using the product.
	Warning: Text marked with a warning symbol must be read before using the product.

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Make sure the SoundConnect device and your hearing aids are within 1 m range of each other. The pairing sequence will start automatically. The light indicator will glow constant green, when the pairing is successfully completed and if sound is played from the computer the audio will be streamed to your hearing aids. (Illu. 5)

The SoundConnect can only be paired to one set of hearing aids. To pair a new set of hearing aids, you must perform the same steps and the new pairing will overwrite the existing one.

NOTE
After pairing and you unplug the SoundConnect and reconnect the device, it will keep the pairing.

RANGE
Line of site up to 8 meters
Outside in open air up to 5 meters

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MAINTENANCE
Clean your device using a soft, dry cloth. Do not use any liquids or chemicals and do not immerse in liquids. When the device is not in use keep it in a dry location out of reach of children and pets. Never try to open or repair the device yourself. This may only be done by authorised technicians.

Símbolo	Título/Descripción
	WECE mark: "Not for general waste". When a product is to be discarded, it must be sent to a designated collection point for recycling and recovering to prevent the risk of harm to the environment or human health as a result of the presence of hazardous substances.
	CE mark The product is in conformity with the requirements set out in European CE marking directives.
	RCM mark The product complies with electrical safety, EMC and radio spectrum regulatory requirements for products supplied to the Australian or New Zealand market.
	Medical device Indication that the device is a medical device.

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Problem	Possibile causa	Solution
The device does not work	1. Not paired 2. Out of range	1. Pair with hearing aids 2. Move into range
No sound	1. Not paired 2. Device is out of the hearing aids' range 3. Audio from PC is muted 4. Another audio output device is selected on the computer	1. Repeat the pairing process 2. Move into range 3. Unmute PC audio, and try to turn up the volume 4. Windows: Click the loudspeaker icon in the taskbar, and choose the Widex SoundConnect as output device
The device does not respond, and the light indicator is solid orange	Device is in a service mode intended for the hearing care professional	Press the button to return to regular use
Device does not respond, and the light indicator is flashing red	Device error	Contact your hearing care professional
Audio from "Microsoft Teams" is coming out of another audio device, even while the SoundConnect is the main audio output device on the computer	On Windows operation systems "Microsoft Teams" has its own audio settings that overrules the ones for the operating system	Go into settings in "Microsoft Teams", choose "Devices" and set the "Speaker" option to the SoundConnect

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