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
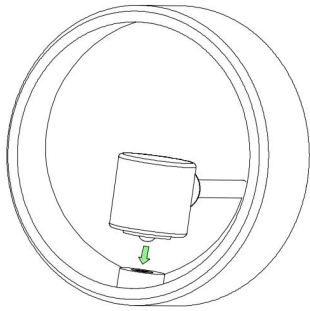
Smart Ring User Manual

Charge and power on

Connect the ring device to the magnetic charging cable. When the red light comes on and the green light flashes, turn on the device.

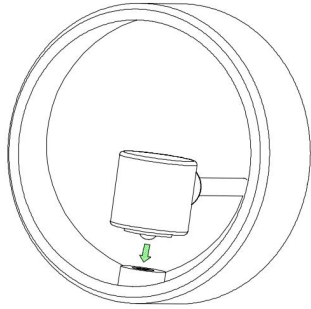Notes

Before using and operating the device, please read and observe the following precautions to ensure the best performance of the equipment and avoid dangerous and illegal conditions.

In daily use, please choose the right size for your finger, avoid wearing too loose or too tight. Please keep the contact part of the ring dry, and keep the sensor and heart rate blood oxygen lamp clean to avoid dirt affecting the test data. If you have any symptoms such as redness and swelling during wearing, please stop using it immediately and consult your doctor.

Security


Please ensure that the charging cable and the device are kept in a dry environment during the charging process. Refrain from touching the charging cable with wet hands. Maintain the cleanliness of the surface of the ring magnetic suction charging port to prevent poor contact at the charging port, which may lead to charging failure. Avoid exposing the charging cable to water and other liquids. Do not disassemble or modify the battery, insert foreign objects into it, or immerse the battery in water or other liquids for an extended period, as this may result in battery leakage, overheating, fire, or explosion. Do not drop, squeeze, or puncture the ring battery. Avoid subjecting the battery to external pressure that could cause an internal short - circuit or overheating of the battery.

1. Charging

Connect the ring device to the charging device. Ensure that the ring charging cable port aligns with the ring charging point, and adjust the alignment until the ring charging indicator light (red light) illuminates.


Note:

- After the device is fully charged, the red light of the charging indicator goes out. At this time, the ring will detect and judge and automatically stop charging. Please remove the ring and unplug the power adapter.
- Keep the charging interface dry and clean before charging to prevent short circuit or other risks.
- Due to the magnet inside the charger, it is inevitable to adsorb metal foreign matter and other foreign matter. Please pay attention to clean up when using. Forbid the charger to contact with high temperature for a long time, resulting in demagnetization of the charger and other abnormalities.

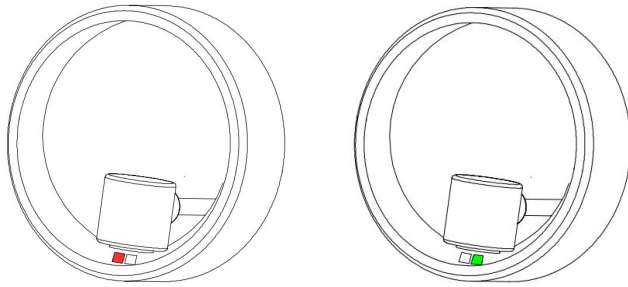2. Wear a ring

- Wear the purchased smart ring on the size of your finger and choose a comfortable size; It is recommended to wear on the index finger, middle finger and ring finger and heart rate IC and heart rate blood oxygen lamp on the inside of the finger.
- Ensure that the sensor part of the smart ring is tightly fitted to the finger to ensure the accuracy of the data;
- Avoid wearing too tight or too loose, so as not to affect the comfort and use effect;
- Ensure that the monitoring unit area is attached to the skin, and there is no foreign body occlusion;

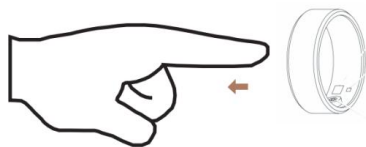Before use:

1. Power on the ring

The ring automatically powers on when it is charged and powers off when it runs out of power. Connect the ring device to the magnetic charging cable. When the red light is on and the green light is flashing, the device starts.

Note: When the ring is out of power, connect it to the original charging device for charging, and the ring will automatically turn on. If you wish to turn off the ring manually, you can go to "Device" > "Other Settings" > "Ring Off" to complete the operation.

2. Download APP

- Please turn on the Bluetooth of the mobile phone, and make sure the mobile phone is

connected to the WLAN or open the mobile network.

- Please download the latest version of APP first. The APP shall be based on Android 6.0 or above, and Apple 9.0 or above

Way one: Search and download

Android: You can download it by searching for "Da Rings" in the App store and App Store

IOS: Search the APP Store for "Da Rings" to download

Way two: Scan and download

Scan the QR code below and use your mobile phone to scan and download.

3. Bluetooth connection

1 > Pair the ring with the APP

- Please turn on the Bluetooth of the mobile phone, and make sure the mobile phone is connected to the WLAN or open the mobile network

Manual connection: open the APP, select "Device" --> Select "Add Device"

Note: (1) The device will automatically synchronize data each time it is rebonded after disconnection or unbinding.

(2) Manually refresh the synchronization data by pulling down from the main page of the APP

2 > The ring is unpaired with the APP

Android: Open Da Rings APP --> Go to "Device" --> Click "Remove Device"

IOS: Open Da Rings APP --> Go to "Device" --> Click "Remove Device"

The main page of the APP

1. Activities

Wear the ring correctly, and the ring will automatically record your daily activity consumption of calories, distance, time and other related data.

- The connection between the ring and the APP is synchronized. On the main page of the APP, the distance of the day's activity, the calories consumed and the duration of the activity can be viewed in the activity card. Click on the activity card to enter the step count, you can switch the daily, weekly and monthly exercise data details, you can view the physical activity, the past 7 days goal achievement and exercise record. Click the icon in the top right corner of the activity page to view the historical activity goal achievement data for the corresponding date.

2. Physical state

Wearing the ring properly, the ring will automatically record your daily activity and goal achievement percentage, sleep quality score, resting heart rate, average heart rate variability, and average stress, and AI will comprehensively analyze the body state of the previous day and give recommendations. The body status card provides a comprehensive view of activity, sleep, sleep ratio, and stress ratio. Click on the icon  in the top right corner of the body status page to view the historical body status data for the corresponding date.

3. Sleep

When worn properly, the ring can automatically track sleep duration, identify deep sleep, light sleep, rapid eye movement (REM) and wakefulness, and score sleep quality based on sleep status.

The ring is connected to the APP. In the APP sleep page, you can see the sleep details of the previous day, sleep legend, sleep score, sleep onset time, wake up time, total sleep duration, sporadic nap duration, rapid eye movement duration, deep sleep duration, light sleep duration, awake duration, and sleep advice. Toggle to view daily, weekly and monthly sleep data details. Click the icon  in the top right corner of the sleep page to view the historical sleep data for the corresponding date.

4. Exercise Records

If you wear the ring properly, you can click on the exercise record card to initiate the exercise and select the corresponding item based on the exercise you intend to perform. During the exercise, the duration, calories burned, pace, number of steps, and heart rate will be recorded.

5. Heart rate

Heart rate was measured at a single time

The ring is connected to the APP. On the main page of the APP --> Heart rate card, the top right part shows the latest heart rate value, and the heart rate fluctuation diagram of sleep and daytime time is displayed in the middle --> Click the bottom to measure the heart rate --> Start the measurement, the heart rate light at the ring end lights up, and wait for the heart rate measurement value --> Click the heart rate card to switch between day, week, month heart rate data details; In the top right corner of the heart rate page  the icon  and  , you can view the historical heart rate data and the historical sleep data for the corresponding date.

Continuous heart rate measurement

The ring is connected to the APP, and the ring automatically measures the heart rate every half hour.

6. Oxygen saturation

A single measurement of blood oxygen

The ring is connected to the APP. On the main page of the APP, the last blood oxygen value is displayed on the top right of the blood oxygen card. Click below to measure blood oxygen --> Start measurement, the blood oxygen light at the ring end lights up, and wait for the blood oxygen value (click below the page to show the last seven blood oxygen measurements). Click to

enter the blood oxygen page and click  on the top right corner of the page to view historical blood oxygen data.

7. Heart rate variability

The ring is connected to the APP. On the main page of the APP, the top right of the heart rate variability card shows the latest heart rate variability value. Click the bottom to measure the heart rate --> Start the measurement, and the heart rate light at the ring end lights up. Click to enter the heart rate variability page, click  and  on the top right corner of the page and to view the historical heart rate variability data and the historical heart rate variability data of the corresponding date.

8. Pressure

The ring is connected to the APP. On the main page of the APP, the pressure value of the last test is displayed on the top right of the pressure card. Click the bottom to measure the pressure --> Start the measurement. Click to enter the stress page, click  and  on the top right corner of the page and to view the historical heart rate variability data and the historical heart rate variability data of the corresponding date.

9. Your physiological cycle

The ring is connected with the APP, and on the main page of the APP, click on the physiological cycle card to enter the physiological cycle page to set the physiological cycle, and predict the physiological period, ovulation date and ovulation period. Notifications can be set.

APP meditation page

Click on the APP side meditation page, there are options to reduce stress and help sleep. Click "Stress Reduction" --> Click "Start" to follow the instructions for stress reduction exercises; --> Click the default bonfire icon below to select multiple sound scenes; Click "Sleep Aid" --> Click "Start" to follow the instructions for stress reduction exercises; --> Click the default bonfire icon below to select multiple sound scenes;

APP device page

1, add device/remove device, check the battery percentage

- Pair the ring with the APP, open the APP, select "Device" --> Select "Add Device"
- When the ring is normally connected to the APP, the device can be removed and the battery percentage can be checked in the device page of the APP.

2. Goal setting

- When the ring is normally connected to the APP, click on the device page of the APP to enter the goal setting. You can set the target step number, calories, exercise time and exercise day goal.

3. Information Settings

- When the ring is normally connected to the APP, click on the device page of the APP to enter the information Settings, and you can set the gender, height, weight, birthday and step length (the APP will estimate the recommended step length based on the height).

4. Other Settings

- When the ring is normally connected to the APP, click Other Settings in the device page of the APP to set the unit format, temperature unit, and daytime heart rate /HRV measurement interval; You can choose to turn on and off low battery warning, all-day pressure, exercise heart rate warning, blood oxygen during sleep, ring off, and period display.

5. Gesture control to take photos

- When the ring is normally connected to the APP, click the gesture control in the device page of the APP to take photos. According to the gesture instructions, you can use the gesture control to take photos on the APP side.

6. Firmware upgrade

- When the ring is normally connected to the APP, click Firmware Upgrade in the device page of the APP to view the firmware information and restore the factory and firmware records.

7. Apple Health

- When the ring is normally connected to the APP, click Apple Health (the "Health" APP in iPhone) on the device page of the APP to access and update your health data.

8. Strava

- Click on Strava in the APP device detail page to connect to Strava.

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FCC Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two

conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

RF Exposure Information

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.