

VIGGLAB AVILON



QUICK START GUIDE

IMPORTANT: THIS IS NOT A MEDICAL DEVICE.

THE VIGGLAB™ AVILON MONITOR IS DESIGNED TO SUPPORT BABY CARE. IT IS **NOT** INTENDED FOR DIAGNOSIS, CURE, TREATMENT, OR PREVENTION OF DISEASE. PLEASE CONSULT A DOCTOR FOR MEDICAL CONCERNS.

THIS DEVICE IS DESIGNED TO HELP MONITOR A BABY'S MOVEMENTS BUT DOES NOT REPLACE PARENTAL SUPERVISION, NOR IS IT A REPLACEMENT FOR PROFESSIONAL MEDICAL ADVICE OR JUDGMENT. ALWAYS CONSULT A HEALTHCARE PROFESSIONAL IF YOU HAVE CONCERNS ABOUT A MEDICAL CONDITION.

THIS MONITOR IS INTENDED FOR USE ONLY WHEN YOUR BABY IS LYING IN A MOTIONLESS CRIB OR BED. DO NOT USE WHEN CO-SLEEPING, OR IN CAR SEATS, STROLLERS, SWINGS, OR ANY MOVING ENVIRONMENT, AS VIBRATIONS OR MOVEMENT CAN AFFECT ITS ACCURACY IN DETECTING ABDOMINAL MOVEMENTS.

THIS DEVICE COMPLIES WITH PART 15 OF THE FCC RULES. OPERATION IS SUBJECT TO THE FOLLOWING TWO CONDITIONS: (1) THIS DEVICE MAY NOT CAUSE HARMFUL INTERFERENCE, AND (2) THIS DEVICE MUST ACCEPT ANY INTERFERENCE RECEIVED, INCLUDING INTERFERENCE THAT MAY CAUSE UNDESIRE OPERATION. CHANGES OR MODIFICATIONS NOT EXPRESSLY APPROVED BY THE PARTY RESPONSIBLE FOR COMPLIANCE COULD VOID THE USER'S AUTHORITY TO OPERATE THE EQUIPMENT.

THIS EQUIPMENT HAS BEEN TESTED AND FOUND TO COMPLY WITH THE LIMITS FOR A CLASS B DIGITAL DEVICE, PURSUANT TO PART 15 OF THE FCC RULES. THESE LIMITS ARE DESIGNED TO PROVIDE REASONABLE PROTECTION AGAINST HARMFUL INTERFERENCE IN A RESIDENTIAL INSTALLATION. THIS EQUIPMENT GENERATES, USES, AND CAN RADIATE RADIO FREQUENCY ENERGY AND, IF NOT INSTALLED AND USED IN ACCORDANCE WITH THE INSTRUCTIONS, MAY CAUSE HARMFUL INTERFERENCE TO RADIO COMMUNICATIONS. HOWEVER, THERE IS NO GUARANTEE THAT INTERFERENCE WILL NOT OCCUR IN A PARTICULAR INSTALLATION. IF THIS EQUIPMENT DOES CAUSE HARMFUL INTERFERENCE TO RADIO OR TELEVISION RECEPTION, WHICH CAN BE DETERMINED BY TURNING THE EQUIPMENT OFF AND ON, THE USER IS ENCOURAGED TO TRY TO CORRECT THE INTERFERENCE BY ONE OR MORE OF THE FOLLOWING MEASURES:

- REORIENT OR RELOCATE THE RECEIVING ANTENNA.
- INCREASE THE SEPARATION BETWEEN THE EQUIPMENT AND RECEIVER.
- CONNECT THE EQUIPMENT INTO AN OUTLET ON A CIRCUIT DIFFERENT FROM THAT TO WHICH THE RECEIVER IS CONNECTED.
- CONSULT THE DEALER OR AN EXPERIENCED RADIO/TV TECHNICIAN FOR HELP.

RF EXPOSURE INFORMATION: THIS DEVICE HAS BEEN EVALUATED TO MEET THE FCC'S GENERAL RF-EXPOSURE REQUIREMENTS AND MAY BE USED IN PORTABLE EXPOSURE CONDITIONS WITHOUT RESTRICTION.




Vigglab Avilon
Model: **AVL-002**
FCC ID: **2BNFI-AVL002**

PLACEMENT



SIMPLE

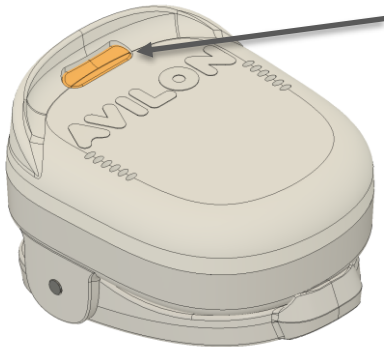

HOLD BUTTON FOR 2 SECONDS TO POWER ON/OFF


MOVEMENT IS SHOWN WITH MOVING LIGHTS.

LIGHTS SPREAD OUT IF MOVEMENT IS NOT DETECTED.
LIGHTS GET NARROWER WHEN MOVEMENT IS
DETECTED.

AVILON WILL ALERT IF MOVEMENTS ARE TOO SMALL OR
STOP FOR TOO LONG.

**THIS PRODUCT IS NOT A REPLACEMENT FOR ATTENTIVE
PARENTING.**

STAY CONNECTED

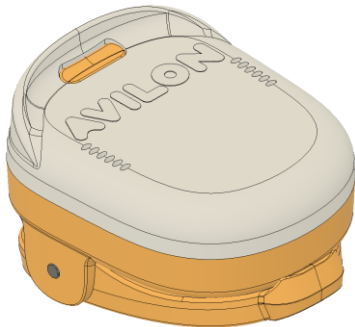

ORANGE SHOWS BABY MOVEMENT, AND WARNS IF MOVEMENT STOPS.

RELIABLE CONNECTION THAT WARNS YOU IF THE SIGNAL IS LOST, AND RECONNECTS AUTOMATICALLY.

POWER



RECHARGE WITH A USB-C CABLE.

RED WHEN CHARGING.

GREEN WHEN FULL.



BATTERY WARNING LIGHT:

- FLASHES SLOWLY IF LESS THAN 12H BATTERY REMAINING.
- FLASHES QUICKLY IF LESS THAN 6H REMAINING.

FULL RECHARGE IN 2 HOURS.

- OVER 40 HOURS OF BATTERY LIFE.

FAQ

Can I clip Aviron to the outside of a sleeping bag or other covers?

- It is better to clip to the diaper. Thick fabric will reduce movement and you may get false alarms.

Can I use a swaddle?

- The swaddle will reduce movement, so you may get false alarms.

Can I use Aviron in the car or in a swing?

- No, it will detect the motion of the car or the swing, so you will fail to get alerts when needed.

What about co-sleeping? Can I sleep in the same bed with my baby.

- No, it will detect your movement, so you will not to get alerts when needed.

What is Ailon made of?

- Medical grade materials. Biocompatible, hypoallergenic, designed for medical devices and pharmaceutical applications.
- ISO10993, USP Class VI. Makrolon 2458 and Cycloy HC1204HF

©2025 Vigglab ApS Denmark. All rights reserved. Vigglab™, the Vigglab™ logo, Vigglab Ailon™, and other Vigglab™ marks may be registered. Trademarks are owned by Vigglab ApS Denmark. All other trademarks are the property of their respective owners.