


116mm

Intelligent Steam Eye Protector

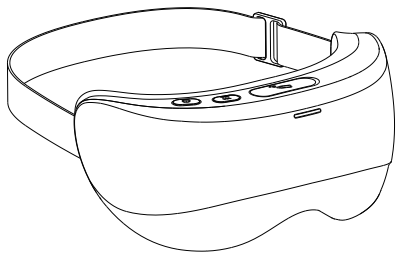


Intelligent Voice Control


		directive commands	Broadcast/Function
wake word	01	hello aimee hi aimee	hi,I am here.
	02	cold mist	cold mist
	03	light mist	light mist
	04	warm mist	warm mist
	05	cleaning mode	cleaning mode
	06	power off	power off
	07	play the music	play the music
	08	stop the music	stop the music
	09	previous song	previous song
	10	next song	next song
	11	turn up the volume	turn up the volue
	12	turn dow the volue	turn dow the volume

130mm

PRODUCT ACCESSORIES

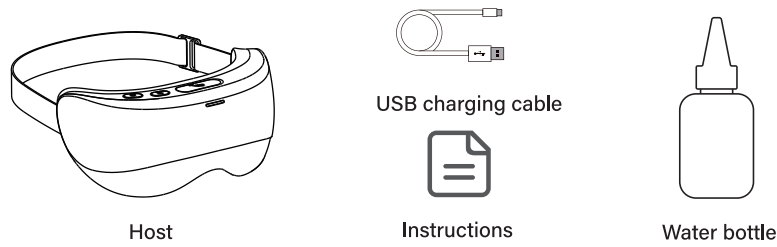

PRODUCT INTRODUCTION

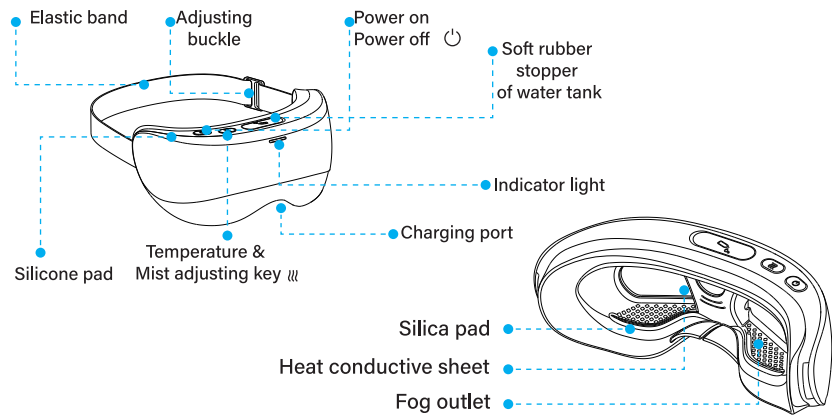

* The illustration is for reference only, and the specific appearance accessories should be subject to the objects

METHOD OF USE

1 USE

- Adjust the elastic band to the appropriate position, wear the eye mask and start using it.
- Please keep your eyes closed during use:**
After a single operation of 10 minutes, the device will automatically shut down. If you need to use it again, you can repeat the 2nd and 3rd actions.



After use, please remove the residual water in the water tank and wipe the device, and store it in a cool and dry place.


Warm Tips

- If the intelligent steam eye protector needs to be stopped halfway during work, it can be turned off by long pressing the power key.
- Please do not turn the equipment upside down when the device is working.
*** When the water volume in the water tank is too small(about less than 2.5 ml), do not make it lie flat for use.**
- If the intelligent steam eye protector is not used for along time, you shal remove all the residual water in the water tank to avoid breeding bacteria.
- The intelligent steam eye protector has a perfect safety protection mechanism, which will give an alarm in time and shut down when the working temperature is abnormal.

2 WATER INJECTION

- Add water to the water bottle.**
Do not add hot water, essential oils and other liquids.
- Open the soft rubber stopper of the water tank upwards, inject proper amount of water into the water tank with a water bottle, and then fasten the soft rubber stopper.
The water tank should not be too full every time, and it is recommended not to exceed 5ml.


! Note: Do not add hot water to avoid scalding!

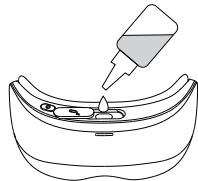

3 POWER ON/POWER OFF

Power on:Long press the power key * for 3s, and the device will be started.
Power off:Long press the power key * for 3s, and the device will be shut down.


Long press the key for starting/shutting down.

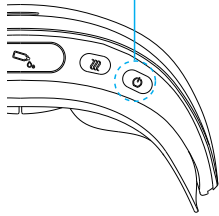

4 TEMPERATURE & MIST MODE ADJUSTMENT


- Mode 1.**
Cool Mist: ambient temperature, only cool mist.
- Mode 2.**
Low warm mist: about 38 degrees, medium warm mist.
- Mode 3.**
High warm mist: about 42 degrees, high warm mist.

Mode adjustment

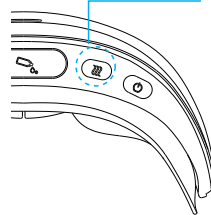

5 MOBILE PHONE MUSIC PLAYING

After the device is power on, turn on the Bluetooth of the mobile phone and search for the connected device"Steam eye mask ELS-02". The device can play mobile phone music, and the mobile phone can turn off/turn on of the music. The mobile phone can adjust the music volume. Note: When switching modes, the music on the phone will pause, and you need to click on the phone to play music again.

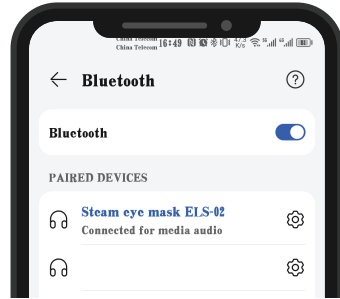

INDICATOR DESCRIPTION


Startup :	The blue & red light is on for 1 minute, and it will automatically go out. If you touch it again, you canwake up the operation.
Charging tips:	The red light is always on during charging; and the red light goes out when charging is completed.

CHARGING INSTRUCTIONS

This product uses a rechargeable lithium battery and should usean adapter passing 3C, with the rated input of the device is **5V 1A**
The charging time is about 3 hours, and after full charge, it can be used about 9 times (affected by ambient temperature, the times are slightly different)

Charging state

When charging, the indicator's "red light" is always on; and after full charging, the "red light" goes out.

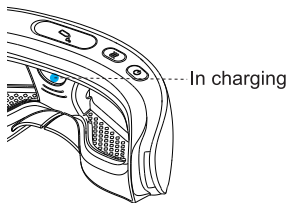

Warm Tips

- Please choose the power adapter that meets the rated voltage of this product.
- Please do not start the product for use during charging.
- This product contains built-in rechargeable battery, so you shall not put it into fire or heat it. (Otherwise, there may be risks of liquid leakage in the battery, overheating or explosion of the battery.)

CLEANING METHOD

- Silicone pad can be disassembled and washed separately. In addition, the whole device is not waterproof, and it cannot be immersed in water or washed with water.
- If there is dust or stain, wipe it with clean and dry cloth or wet cloth that is wrung to dry.
- When cleaning, Please make sure that the device is turned off and the charging cable is unplugged.
- Do not use any detergents or cleansers for cleaning.

Clean the device body

Soak a piece of clean and dry cloth in clean water, and after wringing it dry, wipe the device body with it.

! Do not use alcohol, paint thinner, gasoline or nail polish remover etc, for cleaning.




Clean the components

Remove the silicone pad from the device body, clean it with water, and dry it with a dry cloth. After aired, it will be reinstalled on the device body.

REMOVE THE SILICONE PAD

If cosmetics stick to the body, siicone pad or elastic band, please wipe off immediately, otherwise it may cause discoloration. After cleaning and airing, it shall be stored (Please refer to"Storage").

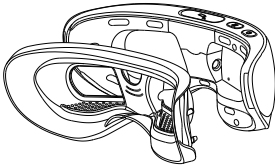

Installation method of the silicone pad

Clamp the silicone pad buckles fully into the 9 connection points (holes)on the device body.



NOTICE FOR USE

- Please keep your eyes closed during use.
- If there is any discomfort in the skin or the eyes during use,please stop using this product immediately and consult the doctor.
- This product is not suitable to diagnosis, treatment and other purposes.
- Please adjust the use time and frequency according to your skin conditions.
- Do not wear contact lenses when using the product.
- Do not use this product with eye mask patches.
- Please do not use this product with makeup on your face (your makeup may be destroyed).
- If you have used eye drops, please wait 30 minutes before using this product.
- If the operation isn't carried out according to the instructions and any problem arises, the Company will not bear any legal liability and related expenses.
- Do not use the product in the following situations:
 - When you feel unwell or have a fever;
 - When there is a wound,eczema or swelling in the skin;
 - When there are abnormalities like dermatitis or excessive sunburn, etc. in the skin.
- This product should not be used by the following personnel:
 - Those with allergic constitution;
 - Those with skin diseases or atopic dermatitis;
 - Those with sensitive skin;(Otherwise, it may cause problems in the skin or the body.)
 - Pregnant women,(It may cause skin discomfort due to hormonal unbalance.)
- Do not use this product for the following persons:
 - Those with facial nerve disorder or skin sensory disorder;
 - Those injecting metal, silicon or plastic materials into the face due to fracture or plastic surgery, and their skin is not capable of sensing heat;
 - Those having head or eye sequelae due to traffic accidents,etc.;
 - Those with high myopia (diopter exceeding-6D);
 - Those who have or had had eye diseases (cataract, glaucoma, retinal detachment, fundus hemorrhage floaters, etc.);
 - Those who have had a history of treatment for eye diseases (those who have recently undergone laser eye surgery should consult medical personnel before using this product).

PRECAUTIONS

- Please store the device beyond the contact range of children and infants. Do not allow children and infants to operate this product.
- Do not use this device in a bathroom or a shower room.
- Please do not refit or repair the device without permission. If you need repair, please contact the after-sales service personnel.
- In case of abnormality and failure, please stop using it immediately, disconnect the power supply and contact the after-sales personnel.

[Abnormal or malfunction conditions]

- Deformation or abnormal heating of the device, adapter or charging cable.
- The device, adapter or charging cable emits a burnt smell.
- The device, adapter or charging cable makes abnormal sound during use or charging.

PARAMETERS

Product name: Intelligent Steam Eye Protector

Product model: ELS-02

Appearance dimensions: 160 X 66 x 30 mm

Product weight: 175g

Rated input: 5V~1A

Operating temperature: It is recommended to use at 18C -40C

Storage/transportation conditions: Temperature: 5°C-45C

Relative humidity: <= 80% RH -without condensate

Atmospheric pressure:(50.0 106.0) kpa

Working hours: 10 min/time

Water tank capacity: 6ml

Charging duration: 3 hours

PRECAUTIONS

Unable to charge	Low power Others	Please charge Please contact after-sales customer service personnel
Unable to charge	Wrong USB cable	Please check the connection of USB cable and adapter
	Wrong adapter	Check whether the adapter parameters are consistent with the product requirements
	Others	Please contact after-sales customer service personnel
No fog	There are bubbles in the water	Shake up and down or restart the device
	Fog spraying	Turn on self-cleaning mode
	Calcium substances in water clog the meshes of the atomizing film	1. Drop 4-5 drops of household white vinegar in the water tank 2. Make the device operate for about 2 minutes and the fog will be recovered
	Others	3. Turn off the device, pour out the water in the water tank, add mineral water in and clean it, and then turn it on Please contact after-sales customer service personnel

Made in China



FCC Warning:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 0cm between the radiator and your body.