

User Manual
MG03

Please read this user manual carefully before use and keep it for future reference.

3.2 Wear

Personal Sound Amplifier L/R Recognition



3.3 Power on/off

Power on: Open the charging box and take out the Personal Sound Amplifier. It will automatically turn on, or long press the button for 5s to turn on.

Power off: Put the Personal Sound Amplifier into the charging box, close the charging box cover, the Personal Sound Amplifier will automatically shut down, or long press the button for 5s to shutdown the Personal Sound Amplifier.


3.4 Volume adjustment

Short press the button each time, the volume increases one level, if the volume is at maximum volume, press the button again, the volume would be turned to minimum volume.



6.2 Viva Hear APP operating guide

(1) Connect Bluetooth:

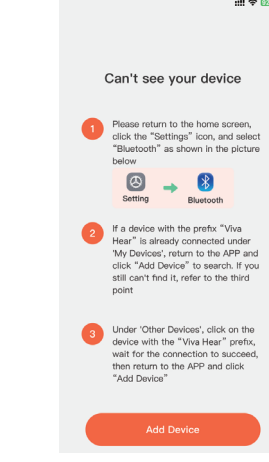

(2) Bluetooth settings:



(3) Device attachment:



(7) Acoustic environment checking: (if too noisy, the test cannot be carried out)



(8) Hearing loss testing:



(9) Test report:



(10) Home page:



(11) Listening fine-tuning: (To be done)



(12) Pre-test preparation:



(13) Start Hearing Test:



(14) Next Step:

(15) Start Hearing Test:

(16) Next Step:

(17) Start Hearing Test:

(18) Next Step:

(19) Start Hearing Test:

(20) Next Step:

(21) Start Hearing Test:

(22) Next Step:

(23) Start Hearing Test:

(24) Next Step:

(25) Start Hearing Test:

(26) Next Step:

(27) Start Hearing Test:

(28) Next Step:

(29) Start Hearing Test:

(30) Next Step:

(31) Start Hearing Test:

(32) Next Step:

(33) Start Hearing Test:

(34) Next Step:

(35) Start Hearing Test:

(36) Next Step:

(37) Start Hearing Test:

(38) Next Step:

(39) Start Hearing Test:

(40) Next Step:

(41) Start Hearing Test:

(42) Next Step:

(43) Start Hearing Test:

(44) Next Step:

(45) Start Hearing Test:

(46) Next Step:

(47) Start Hearing Test:

(48) Next Step:

(49) Start Hearing Test:

(50) Next Step:

(51) Start Hearing Test:

(52) Next Step:

(53) Start Hearing Test:

(54) Next Step:

(55) Start Hearing Test:

(56) Next Step:

(57) Start Hearing Test:

(58) Next Step:

(59) Start Hearing Test:

(60) Next Step:

(61) Start Hearing Test:

(62) Next Step:

(63) Start Hearing Test:

(64) Next Step:

(65) Start Hearing Test:

(66) Next Step:

(67) Start Hearing Test:

(68) Next Step:

(69) Start Hearing Test:

(70) Next Step:

(71) Start Hearing Test:

(72) Next Step:

(73) Start Hearing Test:

(74) Next Step:

(75) Start Hearing Test:

(76) Next Step:

(77) Start Hearing Test:

(78) Next Step:

(79) Start Hearing Test:

(80) Next Step:

(81) Start Hearing Test:

(82) Next Step:

(83) Start Hearing Test:

(84) Next Step:

(85) Start Hearing Test:

(86) Next Step:

(87) Start Hearing Test:

(88) Next Step:

(89) Start Hearing Test:

(90) Next Step:

(91) Start Hearing Test:

(92) Next Step:

(93) Start Hearing Test:

(94) Next Step:

(95) Start Hearing Test:

(96) Next Step:

(97) Start Hearing Test:

(98) Next Step:

(99) Start Hearing Test:

(100) Next Step:

(101) Start Hearing Test:

(102) Next Step:

(103) Start Hearing Test:

(104) Next Step:

(105) Start Hearing Test:

(106) Next Step:

(107) Start Hearing Test:

(108) Next Step:

(109) Start Hearing Test:

(110) Next Step:

(111) Start Hearing Test:

(112) Next Step:

(113) Start Hearing Test:

(114) Next Step:

(115) Start Hearing Test:

(116) Next Step:

(117) Start Hearing Test:

(118) Next Step:

(119) Start Hearing Test:

(120) Next Step:

(121) Start Hearing Test:

(122) Next Step:

(123) Start Hearing Test:

(124) Next Step:

(125) Start Hearing Test:

(126) Next Step:

(127) Start Hearing Test:

(128) Next Step:

(129) Start Hearing Test:

(130) Next Step:

(131) Start Hearing Test:

(132) Next Step:

(133) Start Hearing Test:

(134) Next Step:

(135) Start Hearing Test:

(136) Next Step:

(137) Start Hearing Test:

(138) Next Step:

(139) Start Hearing Test:

(140) Next Step:

(141) Start Hearing Test:

(142) Next Step:

(143) Start Hearing Test:

(144) Next Step:

(145) Start Hearing Test:

(146) Next Step:

(147) Start Hearing Test:

(148) Next Step:

(149) Start Hearing Test:

(150) Next Step:

(151) Start Hearing Test:

(152) Next Step:

(153) Start Hearing Test:

(154) Next Step:

</

FCC Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

RF Exposure Information

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.