

amicoipet Product Manual

Smart Collar/Harness for monitoring Pet's vital signs

The monitoring device utilizes Bluetooth connectivity and is paired with a dedicated collar/harness equipped with a sensing antenna. Through the collar's/Harness's harmless low-power radar waves, pets can be monitored for their vital data without the need to shave their fur. Simply by wearing the collar, pets can have their health data monitored.




Accessories




- ① Package
- ② amicoipet Device
- ③ Sensor Collar/Harness
- ④ Type C cable
- ⑤ Product Manual

Product Design

Device

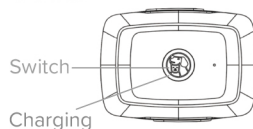

Sensor Collar/Harness

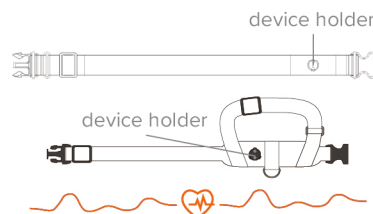

Precautions

- ① A new battery requires 3-5 full charge and discharge cycles to reach optimal performance. Please turn off the device during charging.
- ② It is recommended to use the provided USB Type-C charging cable to charge the device. The charging current does not exceed 1A.
- ③ Each device corresponds to one pet, and there is a possibility of data confusion or failure to recognize if monitoring pets is changed arbitrarily.
- ④ It is recommended to monitor the pet during rest, with the device worn for at least 30 minutes to 1 hour each time.
- ⑤ Due to the self-discharge characteristics of the battery, to avoid abnormal charging caused by excessive discharge, it is recommended to recharge immediately after each complete battery depletion. If the device is not used for an extended period, it is advisable to recharge at least every 3 months to maintain battery performance.
- ⑥ The product has an Ingress Protection Grade 54. It is not suitable for immersion in water or seawater.
- ⑦ This product is not a medical device. The monitoring of pet vital data is primarily intended to facilitate the convenience of pet care for owners and cannot replace medical diagnosis.

Indicator Lights Instruction

Device Power On

Press the power button, successful power-on is indicated by the illumination of the blue light, which continues to flash.

Device Power Off

Long-press the power button, successful power-off is indicated by the blue light turning off.

Indicator Lights	Status
Constant Red Light in the Bottom Left Corner	Charging
Bottom Left Red Light Off	Charging Complete
Blue Light Flashing	Device Powering On
No Indicator Light	Device Powering Off

User Instructions

① Download amicoipet App

Search for [amicoipet] on Google Play or Apple Store, and you can download the app by clicking on it.

Amicoipet



② Registration and Bluetooth Pairing

After the completion of downloading the amicoipet app

- a. Please ensure that your internet and Bluetooth functions are enabled, then power on the device.
- b. Open the amicoipet app, and choose "Register Member."
- c. Follow the app instructions to scan the QR code on the device, complete Bluetooth pairing, member registration, and add a pet.

③ Assembly and Wearing

Ensure that the device is powered on and successfully connected via Bluetooth. In the app, select the pet to be monitored and the desired mode. After putting the assembled collar/harness on your pet, monitoring will begin, and you can view your pet's physiological data through the Amicoipet app.

Note: Each time the device is powered on for monitoring, it needs to collect data for approximately 1-2 minutes before the app starts displaying the data and waveform.

Customer Service:
Mon. to Fri. 09:00-12:30 / 13:30-17:00
Customer Service Email: service@amicoipet.com
amicoipet Website:
<https://www.amicoipet.com.tw/>
amicoipet Line: @amicoipet
amicoipet Facebook:
<https://www.facebook.com/amicoipet.official>
amicoipet Instagram: @amicoipet



amicoipet Facebook



amicoipet Instagram



amicoipet LINE



amicoipet Website