

aaboRing⁴®

Your circle of health



Setup Guide

Charge the Ring



Charge for 60 - 180 minutes

Connect charging case




Plug USB to Power Adapter, Green LED light will blink while charging & it will remain solid green when fully charged

Keep the Ring in Charging Case
till it is Bound with aaboRing® App

Download aaboRing® App




Login aaboRing® App

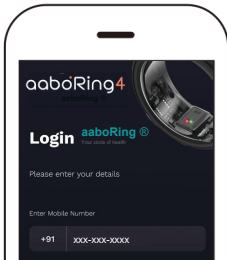


Bind Ring with app

Ring must be in the charging case

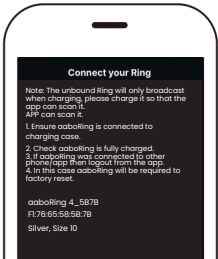

Wear aaboRing⁴

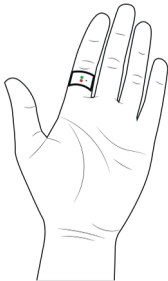

Ensure sensors
touch palm side
of finger

All Set!



Check your mindfulness now



Experience aaboRing⁴

Anytime, Anywhere!

Data capturing starts
within 30 mins

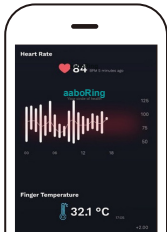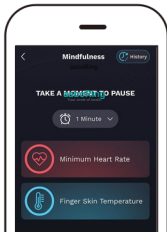

FCC Caution.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

***RF warning for Portable device:**

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.



218-250738

FCC ID : 2BRUZ-AR4

The product is too small to be engraved.

**Enter the circle of
health & wellness**



aabo

aarrooggaa bohilabhaam
Good Health for All

aabo® & aaboRing® are registered trademarks
of 1008 Digital Health Pvt. Ltd. India