


vitahunt

Smart Ring

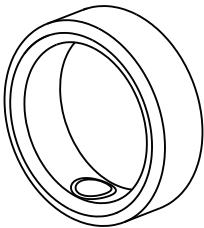

Quick Guide

Important safety instructions

This product contains electronic devices. To ensure your safety, please carefully read the following safety instructions before you use it:

1. Do not disassemble the housing of the product.
2. Do not replace its built-in battery.
3. Do not use corrosive cleaning agents to clean it.
4. Do not charge when it is damp or not completely dry.
5. Do not expose it to direct sunlight for a long time.
6. Do not use it in a sauna or steam bath.
7. Do not place it in an environment with extreme temperature or humidity.
8. Do not use it if it is deformed or internally damaged.
9. Do not put it in a fire as the battery may explode.
10. Please use Vitahunt's standard charging compartment to charge the ring, and use a

certificated power adapter to charge the charging compartment.


11.If the product temperature rises too high or becomes seriously hot during charging, please suspend charging immediately.

12. This product has a built-in lithium battery. Being out of use for a long time could cause capacity loss, expansion, aging, circuit board corrosion, or other safety hazards to the battery or product damage due to the properties of the battery.It is recommended that you charge the battery to about 70% before stopping use, store it in a charging compartment to maintain power, and charge your ring and the charging compartment regularly.

Package contains



Vitahunt Smart ring

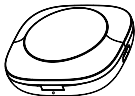

Smart Charging Case




Type-C Cable



Quick Guide

This information is for reference only and does not constitute any form of commitment. Please view the actual product for exact details (including but not limited to color, size, etc). This information is subject to change without prior notice.

Function indication and charging

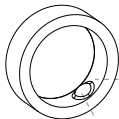


Ring indicator light

Charging: Red light flashes slowly


Fully charged: Solid green light

Ring indicator light



PPG Sensor area

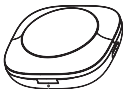

Charging indicator light

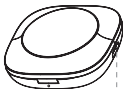

Type-C charging port

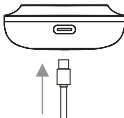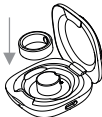

Bluetooth Connection

Ring Bluetooth Name: W596

Download and Install the Vitahunt App

Scan the "Vitahunt App" QR code below to download and install the Vitahunt App. Both Android and iOS devices are supported.

Compatible Systems:

iOS 12.0 and above,

Harmony 2.0 and above, 5.0 and below,

Android 5.0 and above.



Connecting the Ring

1. Binding the Ring

Make sure your phone's Bluetooth is turned on. Open the Vitahunt app, log in, select Bind Device, and select Bind Vitahunt. Follow the prompts to bind the ring.

Note:

- If you are connecting a new ring for the first time, place the ring in the charging dock and connect the charging cable to activate it.
- When the red light on the ring is solid, activation is successful. You can then start connecting.

2. Connecting the Ring

After successfully binding, the ring will automatically connect to the ring upon subsequent power-ups. There is no need to bind again. Note: If the device connected to the ring goes out of range, the Bluetooth connection will automatically disconnect. It will automatically reconnect when it is brought closer.

3. If you are unable to connect, please troubleshoot the following:

- Make sure the ring is activated.
- Check that the Bluetooth on your phone is turned on.
- Check that the Bluetooth permission in the app is enabled.
- Keep your phone as close to the ring as possible.


Correct Ring Wearing

To ensure accurate heart rate and blood oxygen measurements, wear the ring on your index finger as is customary. You can also wear it on other fingers.

Choose a ring size that fits your finger and avoid wearing it too tight or too loose. If redness or swelling develops around the ring, discontinue use immediately and consult a doctor. Remain still while measuring heart rate and blood oxygen for more accurate data.

If the ring becomes stuck when you remove it, try wetting your finger with cold water and mild soap and slowly twisting it to remove it.

Wear the ring on your index finger as you normally would. Wear it with the touch area and indicator light facing your thumb.

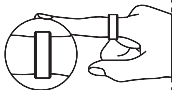

Sensor area

Wear it with the ring sensor bump facing your fingertips.

Product information

Battery Capacity Of The Ring:

3.7V/16mAh(6#-8#)3.7V/18mAh(9#-13#)

Charging Method/Time: Wireless

charging/approx. 90 minutes for full charge

Transmission Mode: Bluetooth 5.0

Compatibility: iOS12.2,Android7.0,Bluetooth4.0 and above

Offline Storage: 7 days

Material: Ceramic

Ring Weight (Net): 3-4grams

Working Temperatur: -10°C~50°C/14°F~122°F

Product Name: Smart Charging Case

Charging Method: Type-C Cable & rated 5V/1A

Charging Case Weight(Net): 40 grams

Instructions for Use

Legal Notice

Copyright © Vitahunt. All rights reserved. This notice may be updated from time to time due to product upgrades or other reasons. Unless otherwise agreed, this notice is for reference only and does not constitute any form of commitment. Please refer to the actual product (including but not limited to color, weight, size, etc.) for details.

Disclaimer

Thank you for purchasing this product. The information in this notice relates to your safety and your legal rights and responsibilities.

By using this product, you and your guest users have carefully read this notice, understand, and agree to abide by its terms and conditions. You and your guest users agree to bear full responsibility for your actions and the consequences of using this product. You and your guest users agree to use this product only for

appropriate purposes and in compliance with all applicable local laws, regulations, rules, and all terms, notices, etc. established or to be established by this product.

This product assumes no obligation or liability for your or your guest users' access to or use of third-party content, or any content or functionality contained therein. You or your guest users will be held responsible for any personal or property damage resulting from improper use or use in accordance with these notices.

You or your guest users must comply with all recommendations and guidelines for any third-party software you use to avoid potential risks. Your or your guest users' obligations and rights regarding the use and access of third-party content are governed solely by the agreement between you and the provider of such content.

Medical Disclaimer

Warning: This product is not a medical device. Vitahunt and related applications are not

intended to diagnose, treat, cure, or prevent any disease or medical condition. Please consult a medical professional before adjusting your exercise, sleep, or diet program.

Important Safety Instructions:

1. If you have any pre-existing health conditions that may be affected by the use of this product, please consult your doctor.
2. This product is designed to record your heart rate, blood oxygen, sleep, and other related data. Prolonged use of the product may cause skin irritation. If you have any redness, swelling, itching, or other allergic symptoms, please consult your doctor and seek medical advice promptly.
3. Please keep the ring at the right level of looseness so that it is not too tight to prevent poor blood flow.
4. Do not use or view the product while driving or in situations that require concentration.
5. This product is not a toy and is not intended for

use by minors. It contains small parts that may pose a choking hazard if swallowed. Please keep it out of the reach of minors.

6. The performance of this product may be affected by external factors. Please do not place the ring in an environment with excessive heat, cold, or outdoor bright light. Therefore, the performance data in this information is for reference only.

7. Do not pierce the ring or the battery in the center as this could risk the product short circuiting.

8. This product may pose a risk to the users with certain symptoms or certain medical conditions.

If you have any of the following conditions, please stop using this product and consult your doctor for advice:

- Epilepsy or sensitivity to light;
- Poor blood circulation or abrasions;
- Heart disease or other medical conditions
- If you are taking any photosensitizing medication.

Usage and Maintenance

1. The size of your finger changes at different times of the day, and sometimes it may be difficult to remove the rings. If you have difficulty removing the ring, please try the following methods:

- Wet your finger with mild soapy water or cold water and slowly turn it to remove the ring;
- Raise your hand above your chest until your blood pressure drops and then try to remove it;
- In case of emergency and if you feel uncomfortable and are unable to remove the ring on your own, please seek medical advice immediately.

2. Remove the ring to avoid wear and tear when performing strength training, holding a shovel or other heavy tool, or holding objects made of metal, ceramic, or stone.

3. Please do not put the ring into an ultrasonic cleaner to avoid damage to some of the

components.

4. The product is waterproof in daily life. You can wear it in swimming pools and other shallow water activities. Please avoid using it in hot springs, saunas (steam baths), diving, and any other environment with extreme temperature or humidity.

5. While using this product, do not apply too many skincare products to your hands to avoid corrosion of its surface.

6. Clean the ring regularly with water and mild soap, and wipe the ring with a soft cloth, paying special attention to the parts in contact with the skin.

7. Remove your ring from time to time to give your skin a chance to relax, and you may use this time to recharge your ring.

8. The ring may be scratched and may scratch the skin that is in close contact with the ring. Jewelry made of softer materials such as gold and silver, and electronic products with soft coatings should be used with care and

protection to avoid scratching by the ring.

Regulatory USA



FCC compliance statement for the VITAHUNT RING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B

digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a

circuit different from that to which the receiver is connected.

- Consult the dealer or an experienced radio/TV technician for help.

To maintain compliance with FCC's RF Exposure guidelines, This equipment should be installed and operated with minimum distance between 20cm the radiator your body: Use only the supplied antenna.

CHARGER

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a

circuit different from that to which the receiver is connected.

-Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.



RoHS

Made in China

W691A-A02-说明书

尺寸: 58x68mm

材质: 80g铜版

成型: 折页

印刷: 彩色双面印刷

此框内容只备注不印刷

vitahunt

Email: service@vitahunt.com

EU	REP
----	-----

UK	REP
----	-----