

AIVELA RING PRO

Quick Start Guide

WHAT'S IN THE BOX

Your AIVELA package includes the following items:



Charging Case



Quick Start Guide




Ring



Type-C Charging Cable


CHARGE THE RING

To activate the ring, **please charge it first**. The ring charges magnetically using the charging post inside the case. **Align the magnetic charging pin** on the inner side of the ring with the charging post. **A blue LED** indicator on the charging case **will light up** to confirm that charging has started

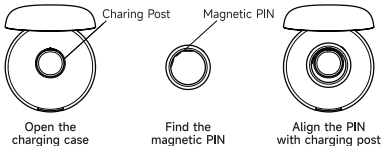


KNOW YOUR RING

AIVELA Ring Pro is equipped with **PPG**, **IMU**, **Temp** and **OFN sensors**. Find the OFN sensor (a black diamond-shaped window) on the outer surface, and the health sensors (little square windows) on the inner surface.

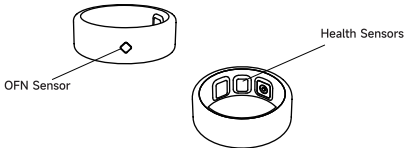

WEAR THE RING

For convenient use of the touch features, we recommend wearing **the ring on your index finger**. Make sure the PPG sensors face the pad of your finger and are directly underneath it, and that **the OFN Sensor faces your thumb** and **is aligned at 90° to your finger**, on its exact side.

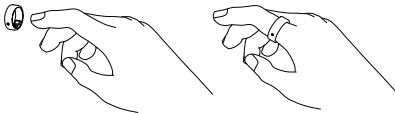

DOWNLOAD THE APP

Please search "AIVELA" in **App Store/Google Play** or use your phone to scan the **QC code** below.



aivela.com/pages/app-download




Download on the
App Store



GET IT ON
Google Play

PAIR THE RING

1. Place your ring inside the charging case to make it discoverable.
2. Ensure your phone is near the ring and **that Bluetooth is enabled.**
3. Open the **AIVELA App** and **follow the on-screen** instructions to complete the pairing process.

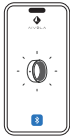

Congratulations! Your ring is ready to use. For more information, please go to User Guide & FAQ in the App or visit aivela.com/pages/user-guide.



AIVOLA

© 2025 AI Plus Appliances Inc. All rights reserved.