

尿液报警器使用说明



感谢您购买我们的产品使用前，
请务必仔细阅读此说明书并妥善保管，以备将来参考。

产品组成



尿液报警器由三个部分构成:
发射器(尿液检测): 检测尿液通过无线方式发送报警信息给接收器
接收器: 接收器用于接收发射器报警信息，通过震动或者声音方式通知用户。
充电盒(收纳盒): 用于发射器，接收器的充电。不充电时盖上盖子后方便用户收纳，携带。

开关机



接收器报警模式设置
震动模式设置:
按开机后默认无定时 小图标亮起，双击开关键解除无定时，然后单击开关键（设置模式）至震动 小图标不停闪烁，按确认键，震动模式设置完成。

定时模式设置

双击开关键解除无定时，继续单击至，30、60、90分钟的图标亮起，然后按确认键，定时设置完成。



报警提示

当发出报警提示时，用户单击开关键或确认键就可以解除报警提示



尿液检测模块佩戴方式1




尿液检测模块佩戴方式2



产品充电

将接收器，发射器放入充电盒，插入USB充电。当充满电时，充电盒的两个红色指示灯会熄灭。

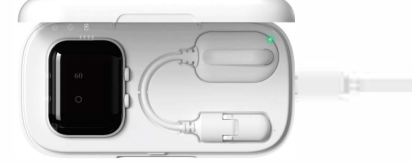

注意事项

1. 使用接收器温度功能时接收器不得接触或放置于如口袋等特殊环境，不将将温度传感器孔位遮挡，接收器在充电过程中温度功能无效。
2. 产品充电采用标准DC5V供电
3. 充电完成后，产品不使用时可将接收器关机以节省电能消耗。

产品介绍

重要

在孩子入睡大约1至2个小时，家长想要完全弄醒孩子是不容易的，而报警将会帮助孩子阻止睡得太熟的孩子不会被叫醒。通过白天让孩子小便一会，确保孩子在晚上睡觉不会筋疲力尽这样将会有助于报警器将孩子叫醒。
每天晚上提醒您的孩子，当听到报警声音时，须停止排尿并上厕所。
在晚上上床睡觉前，没必要限制液体的摄入量。

重要

在最初的几周里，当报警器响起时孩子自己可能不会被报警器叫醒，所以孩子需要家长的帮助；叫醒孩子并上厕所。
孩子参与的过程是孩子取得成功不可缺少的。不要为孩子按下报警器，让孩子自己感觉报警器响起并自己取下或断开报警器为更好。家长需要帮助孩子上厕所并完全。在最初的一到两周里，孩子可能会在未上完全尿湿但家长仍需要继续提醒上厕所的重要性并且鼓励他们。问问孩子第二天是否记得报警器响起和他们做了些什么。

重要

家长的参与对于这个计划的成功也是必不可少的。在最初的几周里，孩子需要家长的帮助和鼓励。这个计划可能缓慢和沮丧的，但家长坚持每天晚上使用报警器，每次训练过程报警器都随正常响起，家长要帮助叫醒孩子。
你需要继续使用报警器3至4周，亦确保孩子在最后一次报警后不会复发。复发是指孩子在一个月或更多月之后又连续尿湿。如果这种情况发生，请立刻开始再次使用报警器，这样复发的间隔将会变得越来越长。

计划

在报警器发生响声后，要给孩子重新连接报警器，同时要给孩子换上一条新的内裤。在孩子的房间使用一个通宵灯，这样有助于孩子醒来时更容易醒来，并且能够找到去卫生间的道路。
每次当孩子通过一个小时的睡眠或感到醒来，一定要表扬他们。陪伴重要的是不要因为孩子没有感到报警器而责骂他们。因为改变孩子的行为要花费一些时间。此外，使用进度图是很有帮助的，以便孩子能参与到计划的进展中。

保养/清洁

在每次有尿湿部分后，记住要清理传感器夹子。如果有需要的话，可一并取下传感器。不要把报警器本体浸在水里，因为这样会损坏内部零件。传感器夹子可在水里清洗。
一定要经常检查电池，以确保报警器仍在工作。这可以通过在夹子中间插入一块湿布来检查。当很长时间不使用报警器时，请取出电池。
本公司报警器是为单个用户设计的，如果你想另一个孩子使用报警器，请联系购买另外一个报警器。

警告：使用中报警后，请将发射端夹子打开并用干的纸巾擦干净，方可使用。

BED WETTING ALARM CY-090



Thank you for purchasing our products.
Please read the manual carefully before use and keep it properly for future reference.

Product Composition



The urine alarm consists of three parts:
Transmitter (urine detection): Detect urine and send alarm information to receiver by wireless way.
Receiver: The receiver is used to receive the alarm information from the transmitter and notify the vibration or sound.
Charging box (storage box): used for charging transmitter and receiver. When not charging, it is convenient for users to store and carry after being covered.

Switch on / off



Receiver alarm mode setting

Vibration mode setting:
After pressing the power button, the default No Timer icon lights up. Double-click the power button to disable the No Timer mode, then single-click the power button (to enter setting mode) until the "vibration" icon flashes continuously. Press the confirmation button, and the vibration mode setup is complete.



Timer Mode Setup

Double-click the power button to disable the No Timer mode. Continue single-clicking until the 30-minute, 60-minute, or 90-minute icon lights up, then press the confirmation button to complete the timer setup.




Alarm prompt



When an alarm notification is triggered, the user can press either the power button or the confirmation button to deactivate it.



Wearing Method of the Urine Detection Module-1

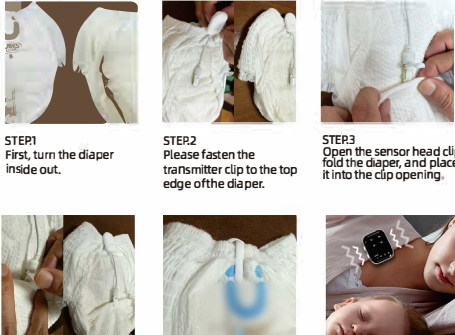


Wearing Method of the Urine Detection Module-2

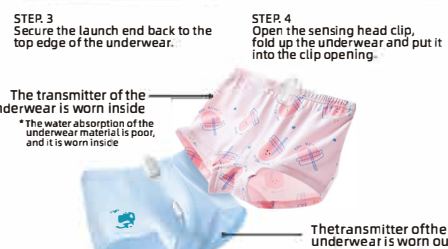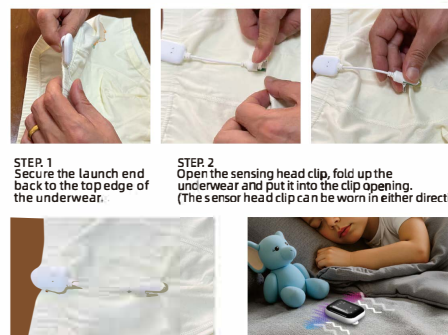

Product Charging

Please place the receiver and transmitter into the charging case and plug in the USB cable for charging. When fully charged, the two red indicator lights on the charging case will turn off.

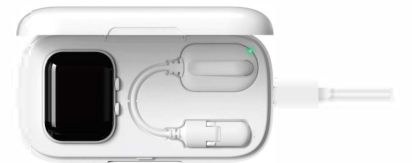

Matters needing attention

1. When using the receiver's temperature and humidity function, the receiver shall not contact or be placed in a special environment such as a pocket, and the hole position of the temperature and humidity sensor shall not be blocked. The receiver's temperature and humidity function is invalid during charging.
2. The product is charged by standard DC5V
3. After charging, the receiver can be turned off when the product is not in use to save power consumption.

Product introduction

Importance

It is not easy for parents to fully wake a child who has been asleep for about 1 to 2 hours, and the alarm will help prevent deep-sleeping children from falling to be awakened.
Allowing the child to take a short nap during the day to ensure they are not exhausted at night will help the alarm wake them up effectively.
Remind your child every night that when they hear the alarm, they must stop urinating and go to the toilet immediately.
There is no need to restrict fluid intake before going to bed at night.

Importance

In the first few weeks, when the alarm sounds, the child may not be woken up by the alarm, so the child needs the help of parents: wake up the child and let him go to the toilet.
The process of children's participation is indispensable for their success. Don't unplug the alarm for your child. It's better you're your child to notice that the alarm sounds and remove or disconnect the alarm by himself.
Parents need to help the child to go to the bathroom to discharge urine. In the first one to two weeks, children may be completely wet in bed, but parents still need to continue to emphasize the importance of going to the toilet and encourage them. As the children if they remember the alarm and what they did the next day.

Importance

Parents' participation is also essential for the success of this program. In the first few weeks, children need parents' help and encouragement. This plan may be slow and frustrating, but parents insist on using the alarm every night, and the alarm can sound normally during every training course, so parents should help to wake up their children.
You need to continue using the alarm for 3 to 4 weeks, and make sure that your child will not relapse after the last alarm. Recurrence refers to a child's continuous wet urine after one or more months. If this happens, please start using the alarm again immediately, so that the recurrence interval will become longer and longer.

Plan

After the alarm sounds, reconnect the alarm to the child and put on a new pair of underwear. Use an all-night light in the child's room, which will help the child wake up more easily when the alarm sounds and find the way to the bathroom.
Be sure to praise your children every time they spend a dry night or notice the alarm. It is also important not to scold children because they don't notice the alarm sound, because it takes some time to change their behavior. In addition, it is helpful to use the progress chart so that children can participate in the progress of the plan.

Maintenance/cleaning

Remember to clean the sensor clip after every wet part. If necessary, clean the alarm once and for all. Do not immerse the alarm body in water, because it will damage the internal parts, and the sensor clip can be washed in water.
Always check the battery to ensure that the alarm is still working. This can be checked by inserting a damp cloth in the middle of the clip. When the alarm is not used for a long time, please take out the battery.
The alarm of our company is designed for a single user. If you want another child to use the alarm, please contact to buy another alarm.



Warning: After use, please open the transmitter clip and wipe it with dry paper towel before use.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Radiation Exposure statement

The device has been evaluatec to meel general RF exposure requirement. The device can be used in portable exposure condition without restriction.