



User Manual



For your safety and convenience, please read user manual before using product.
Please keep the included warranty information for your records.

About BeCozU

What's Included




01 Charging Cradle



02 Earbuds



03 Ear tip

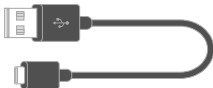

04 USB Type-C
Charging Cable



05 User Manual &
Product Warranty


Caution

User only authentic BeCozU components. If product is damaged due to use of inauthentic components, warranty is void. Please use provided USB-Type C cable only. If the ear tips are damaged, you can purchase replacements by contacting customer support email.


Parts of the BeCozU

Charging Cradle

Front


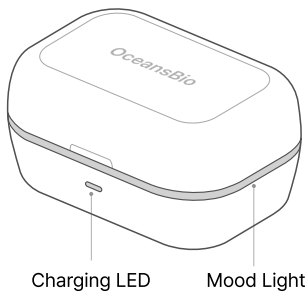


Back

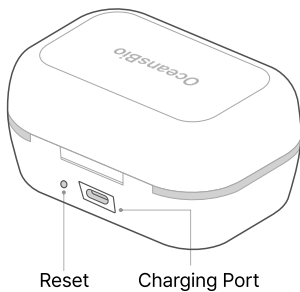

Parts of the BeCozU

Earbuds

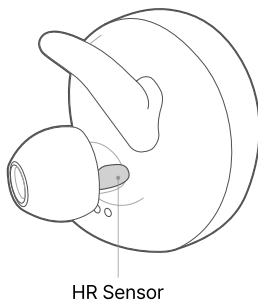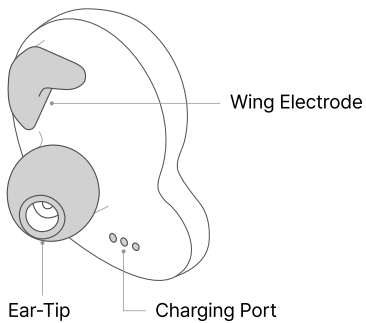

How to wear earbuds



- 01 Take the earbud out of the charging cradle and fit it properly in your ears.



- 02 Adjust the earbud to fit your ears. Effectiveness of the device may be limited if the electrodes do not make proper contact due to variations in ear shape.

App installation

Android users can download BeCozU from Google Play, and iPhone users can download it from the App Store. Search for 'BeCozU' or scan the QR code below.



Android



iOS

Caution

Only Android 13 / iOS 16 or higher is supported, and Bluetooth 5.3 and location information setup is required on mobile devices.

Note

To connect with another device, follow the steps above again.


How to connect



- 01 Turn on Bluetooth on the smartphone or device you want to connect.



- 02 Open the charging case lid while the earbuds are inside.

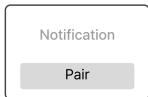


- 03 When the fast notification appears on your smartphone, tap to pair connect.




- 04 If the connection pop-up does not appear, select and connect to [BeCozU] on your smartphone's Bluetooth menu.


How to connect - iOS

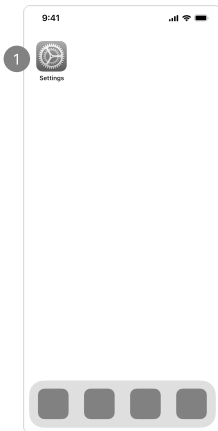

1. Open settings

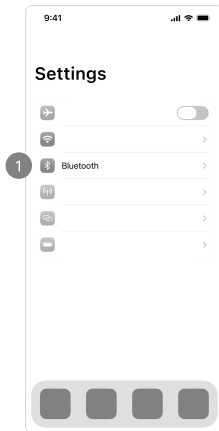

1. Select Bluetooth from the menu

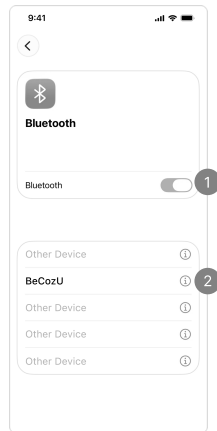

1. Toggle Bluetooth on
2. Select BeCozU from the list to pair

How to connect - Android



1. Select Bluetooth from the menu




1. Tap pair with a new device




1. Select BeCozU from the list to pair

How to Connect

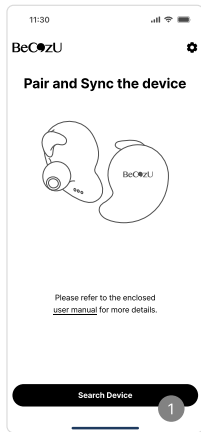


1. Search device

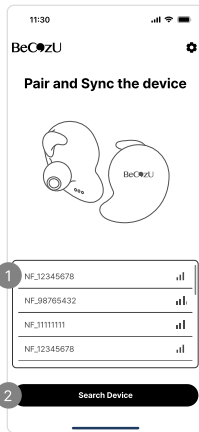

1. Select device and pair
2. Select 'BeCozU(8-digit device serial number)' from the search list to complete the connection.
 - If the device is already registered to another account, it cannot be registered again.
(One device per account)


How to Connect

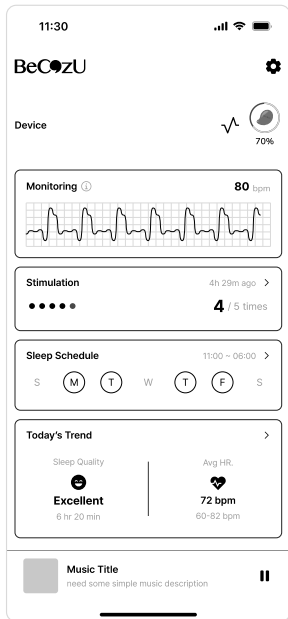


Note

- The device connects automatically.
- Battery level is displayed when the device is connected.
- Pulse indicators and HR monitoring activate when the device is worn.

How to use the app

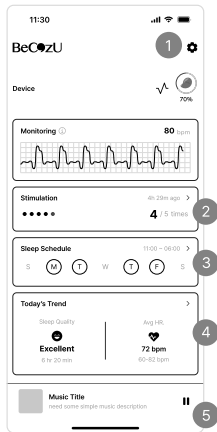


1. Settings
2. Stimulation
3. Sleep & Stimulation alarm schedule
4. Daily & Weekly stimulation and sleep records
5. Music playback

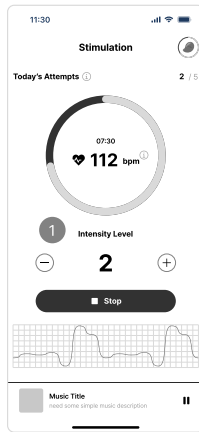


1. Stimulation intensity selection
- This screen is accessible only when the device is connected.
 - When the device is disconnected, you will be returned to the previous screen.

How to use the app

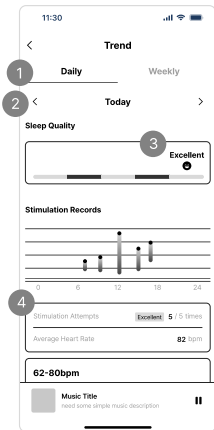

1. Daily & Weekly view
2. Date navigation
3. Sleep quality
4. Stimulation log & Status

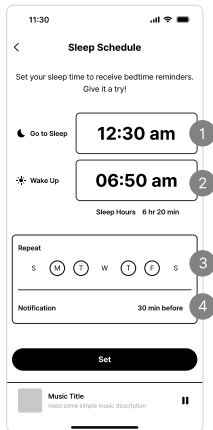

1. Set bedtime
2. Set wake-up time
3. Alarm repeat settings
4. Set alarm time

Charging

The charging cradle shows its charging status when connected to power.

Note

- The earbud and charging cradle charge at the same time.
- The charging case connector is not designed to provide power to other devices, such as smartphones.
- When the charging cradle battery is depleted, the earbud cannot be turned on or off. Please charge the charging cradle.
- Ensure the earbuds are inserted correctly into the corresponding left and right slots in the charging case, then connect the charger.
- Improper connection of the charger may cause damage to the charging case. Damage resulting from incorrect usage is not covered under the warranty.
- Unplug the charger when not in use. As the charger does not have a power switch, it must be disconnected from the power outlet to prevent unnecessary power consumption.
- Moisture or sweat on the charging terminals may lead to corrosion. If there is any moisture or sweat on the charging terminals or earbuds, wipe them dry before placing the earbuds in the charging case.
- The earbuds and charging case may become warm during charging. This is a normal condition and does not affect the product's performance or service life. If the temperature exceeds a certain threshold, charging may stop automatically as a safety precaution.

Safety information

WARNING

- Do not place heavy objects on product.
- Do not drop product from high places.
- Do not allow the product to come into contact with water, alcohol, or benzene.
- Do not store the product in damp or dusty environments.
- Do not put the product in water or use it in very humid places.

CAUTION

- Do not disassemble, modify, or repair the product yourself.
- Keep the product away from excessive heat and flammable materials.
- Use only a 5V, 1A certified charger. Do not use fast chargers, as they may damage the product.
- Small parts attached to magnets and products are dangerous if swallowed.
Keep them out of reach of children.
- Do not allow foreign objects to enter the charging contacts or power plug.
This may cause a fire or electric shock.
- Users with medical devices should consult a doctor before using the product.
- The built-in battery cannot be removed or replaced by the user. For your safety, do not attempt to remove it. Improper or careless battery replacement may cause an explosion.
- Lithium batteries are hazardous and can cause injury if mishandled.
- Replacing the battery by an unqualified person may damage the product.
- Dispose of batteries properly, following local regulations; do not discard carelessly.
- Use only certified authorized charging adapters.
- Disposing of batteries in fire or hot ovens, or mechanically crushing or cutting them, may cause them to explode. Extremely low air pressure may cause the battery to explode or flammable liquids or gases to leak.

Battery indicator

Status	Indication	Description
Cradle Charging	Mood Light ON	USB-C
	Red LED ON	
Cradle Fully Charged	Green LED ON	Battery Full
Cradle Low Battery	Red LED 5 Blinks	Low Battery
Earbud Charging	Green LED 5 Blinks	Earbuds Charging
Cradle Standby	Red Blinks 5 Sec	MCU Normal

Note

Earbud charging status via App

Product specification

Item	Explanation
Model Name	OC-NF01
Bluetooth Specification	Bluetooth 5.3
Charging Input	DC 5V 1A
Playing Time	30min
Charging Time	Charging Cradle : 750mAh Earbud : 100mAh Fully Charged : 4hours
Size/Weight	Earbud : 28×37×25 / 13g Charging Cradle : 62×90×45 / 93g Total Weight : 120g
Frequency	25Hz
Verification No.	FCC ID : 2BTFO-OC-NF01
Customer Support	service@oceansbio.com



Safety Manual

For your safety and convenience, please read user manual before using product.
Please keep the included warranty information for your records.

Precautions for Use

- The device provides electrical stimulation with a maximum output of 2mA, which is within the safe range for human use.
- During use, you may feel a mild electrical sensation. While some individuals may find this slightly uncomfortable, the device uses a lower current than typical stimulators and does not pose a significant safety risk.
- Depending on the individual, symptoms such as tinnitus, pain, abnormal sensations, skin itching, or erythema may occur.
- If you experience any discomfort while using the product, stop using it and consult a doctor or visit a medical facility.
- This product is a microcurrent stimulation device. Because it uses microcurrent, some users may not feel any stimulation. If you wish to feel the stimulation, clean your ear with a wet tissue before use, and start the stimulation while the ear is slightly moist. This can help enhance the sensation of stimulation.

Important Safety Instruction

- Do not disassemble the product.
- Keep out of reach of children. The product contains small parts that pose a choking hazard.
- Do not expose the battery to extreme temperatures above 50°C (122°F) or below -10°C (14°F). Use the product at room temperature to maintain maximum battery capacity.
- Battery capacity decreases when used in low temperatures.
The battery can be charged at temperatures between 0°C and 45°C.

- Do not use the product in excessively high or low temperature environments, and do not expose it to strong sunlight or extremely humid conditions.
- Ensure no liquid, oil, or water gets into the charging port during charging.
- Do not submerge the earbuds in water or allow water to splash onto them.
- This product generates a strong permanent magnetic field that may interfere with pacemakers, implantable cardioverter-defibrillators (ICDs), and other implanted devices. Always maintain a minimum distance of 6 inches/15 cm between product components containing magnets (earbuds and charging case) and pacemakers, ICDs, or other implanted devices.
- The Bluetooth word mark, and logos are registered trademarks owned by Bluetooth SIG, Inc. All applicable trademarks are subject to the license agreement. Android, Google Play, the Google Play logo, and other related marks and logos are trademarks of Google LLC.
- The Bluetooth BeCoZU stereo device is a wireless transceiver. It communicates with Bluetooth-equipped mobile devices by transmitting and receiving radio frequency(RF) electromagnetic fields(ultra-short waves) in a frequency band of 2400 MHz to 2483 MHz when operating.
The following use environments can affect the range and sensitivity of Bluetooth communication.
- When an object like a wall, metal, or person is between the Bluetooth device and connected device.
- When using the device near other devices that use 2.4GHz frequencies (e.g. wireless LAN devices, wireless phones, microwaves)
- When using devices other than smartphones that support Bluetooth (like laptops), playback quality can be affected depending on device support standards/performance.

Hearing Safety

- To prevent hearing damage, avoid listening at high volumes for extended periods.
- Do not turn the volume up so high that you cannot hear ambient sounds.
- Do not use the product while driving.
- Exercise caution or temporarily stop use in potentially hazardous situations.
- If you experience ringing in your ears or difficulty hearing speech, get your hearing checked. Higher volumes can cause hearing damage more quickly.
- This product is not intended for listening to music. Audio quality is guaranteed only for the 25 sleep-supporting audio tracks included in the application. The sound quality of other audio files cannot be ensured.



Note

This product has been tested and found to comply with the limits for Class B digital device, pursuant Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product. Off and on, the user is encouraged to try to correct the interference by on more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

WARNING

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note : This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generate, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Warranty

Oceansbio. Co., Ltd will warrant the equipment according to the consumer reward policy as follows :

- Contact nearest dealer or sales representative from whom the equipment was purchased in case that the equipment found defective. (Contact Oceansbio directly if no dealer or representative was found near you)
- Warranty service will be performed by Oceansbio or Oceansbio's authorized service.

Product Name : BeCozU

Serial No.

Model No : OC-NF01

Dealer Name

Date of Purchase : Year Month Date

Amount You Paid

- Service of product in-warranty: The defective equipment will be serviced at free in case the equipment is in-warranty period
(the equipment is defective within one year from the date of original purchase)
- Service of the product out of warranty or service to be charged

1) Warranty of product with no defect found

- The warranty of the product with no defect found will be charged and you are recommended to read the user manual carefully before the request of warranty service.
- Service for the defect caused by the transportation or re-location of the equipment will be charged.
- The defect of the accessories in the manufacturer's materials and workmanship will be excluded.
- The equipment is requested to be reinstalled after having been installed by the user based on the purchase condition at the time of purchase.

2) Warranty of the product caused by misuse or fault of the user shall be void

- The defect resulted from misuse, alteration, unauthorized disassembly, repair and installation.
- The defect resulted from the improper or inadequate repair by user or service performed by anyone other than Oceansbio or a dealer authorized by Oceansbio.
- The defect resulted from the connection of improper electric utilities.
Ex) When the 220V device is connected to the 110V electric utility.
- The defect or damage resulted from the careless movement or transportation after installation.
- The defect resulted from the use of accessories, consumables or options not manufactured nor approved by Oceansbio for use with this equipment.

3) Others

- Defect and damage due to the accident or force majeure conditions like fire, damage from seawater, flood, freezing and earthquake consumable over-used its life time.
- In case of finishing life cycle for consumable parts.

BeCozU