



Follow the official
account



Scan to get the
user manual.

instruction manual



Manufacturer: Zhejiang Junkong Intelligent Technology Co., LTD

Address: Room 401, Building 13, No.36, Changsheng South Road, Economic and Technological
Development Zone, Jiaxing City, Zhejiang Province

Thank you for using this product

Please read this manual carefully before use and keep it properly for future reference

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement

Safety precautions

Before use, please read this "Safety Precautions" carefully, which will help you to use the product safely and correctly.

The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction;

children being supervised not to play with the appliance ;

The appliance incorporates an earth connection for functional purposes only;

During use, make sure the power plug of the household appliance is fully inserted into the power socket and maintains good contact;

When using an extension cord power socket, use one that meets national standard requirements;

Disconnect the power when leaving the site;

Stop using immediately if the product makes unusual noise, emits unusual odors, becomes excessively hot, or malfunctions; turn off the switch and disconnect the power;

Do not supply power to household appliances from a lamp head with a socket;

The total rated current value of household appliances connected to an extension cord power socket should not exceed the rated current value of the original fixed socket or circuit;

Do not use this product near flammable or explosive items or similar environments;

Users should not disassemble household appliances themselves, nor change internal wiring, components, and protective devices;

When unplugging the plug from the socket, do not pull on the power cord;

Users should not replace the fuse of household appliances and power supply with ordinary conductors or fuses exceeding the specified current capacity.



warn



**Do not
disassemble**

Except for the repair technicians designated by the company, the product shall not be disassembled and repaired or modified. If there is any problem, please contact the dealer in time.



pay attention to



prohibit

1. Do not use this product near the instrument or equipment that generates magnetic field or microwave, do not share the same power supply with the instrument or equipment that can generate magnetic field or microwave;
2. This instrument is a precision equipment, in the process of transportation, storage and use should not be inverted, vibration, collision;
3. When replacing the attachment, it is necessary to use the special attachment of Junkong. The specified size range of the moxa stick (diameter 4.5/5cm; length ≤ 10 cm);
4. Do not use sharp objects to stab the main engine and accessories, do not splash liquid on the machine;

The above situation prohibits the continued use of this product!



force

1. If the product is moved from a low temperature place to a warm place for use, please wait for an hour before use to avoid the possibility of water vapor condensing on the machine and causing problems;
2. The optimal duration of each use is 45 minutes;
3. Do not plug or unplug the power cord with wet hands. Use a separate power outlet and ensure reliable grounding;
4. Do not place the product outdoors, should be kept ventilated, avoid direct sunlight, avoid using the product in damp, dusty places;
5. People who are not sensitive to heat must be accompanied by normal people when using it;
6. Please use this product in an environment where the room temperature is higher than 30 degrees and lower than 16 degrees;
7. To avoid the risk of electric shock, this equipment must be connected to a power supply network with protective grounding.
8. Do not use in the humid environment of high tide, so as to avoid electrical shock caused by damp electrical appliances;


	9.Do not use other instruments while using this machine; 10.In the process of use, if the body appears or feels abnormal, please stop using immediately and consult a doctor; 11.Do not use for therapeutic purposes; 12.If the power cable is damaged, please contact the dealer to consult and replace the qualified power cable to supply power to this product;
 Unplug	1.When performing maintenance, unplug the power cord from the device when it is off. In addition, do not plug or unplug the plug with wet hands to avoid electric shock and injury. 2.In case of power failure, please immediately cut off the power switch and unplug the power plug from the socket.

product configuration

Your robot contains the following:

order number	name	quantity
1	main engine	1 set
2	power line	1 root
3	Heat resistant silicone sleeve	One
4	telecontroller	One
5	Pre-filter element	Twenty
6	Product specification (electronic version)	1 part
7	work certificate	1 sheet
8	Product warranty card	1 sheet
Note: The product will be adjusted when it leaves the factory. Please refer to the attached configuration list for details.		


Product appearance and structure description

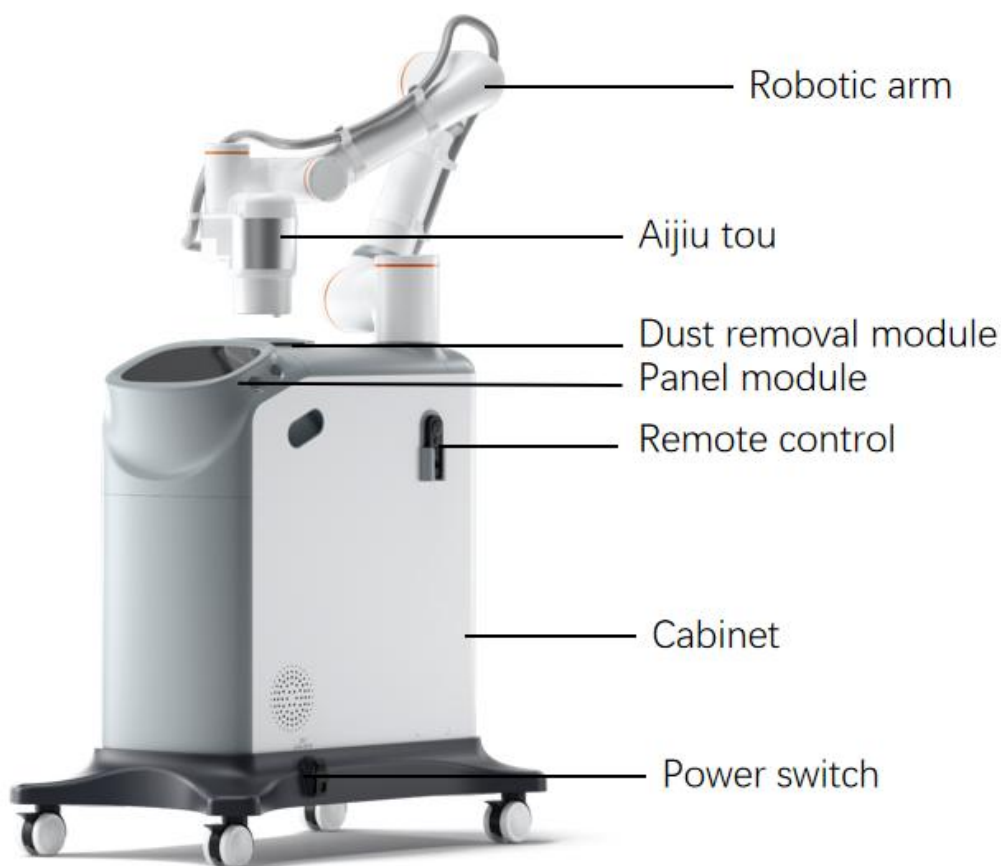

parameter declaration

product name :	JMR robot	product model :	JFR6
Rated input power:	220V~/50Hz	Weight of the device:	106KG
power rating :	300VA	operative norm	
Timing function	5~60min, step: 5min	The minimum distance between the acupuncture point and the human body	30mm
Interval range of powder ash:	1~10min adjustable	Function of remote control button:	Temperature +, temperature-, volume +, volume-, pause, call, switch music

method of operation

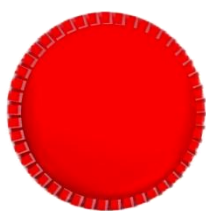


scram button



pause button



Open (close) button

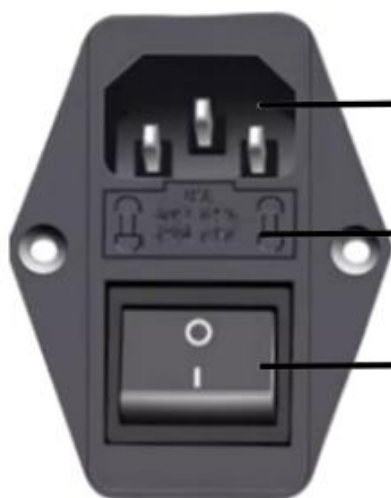
 Click the Confirm
button


three-prong plug end

fuse box

toggle switch

1. First, insert one end of the power cord into the device —— (pin-shaped plug end);
2. Insert the other end of the power cord into the nearest socket on the device;
3. Move the boat-shaped switch above the power filter to "I" for a normally open state;
4. Click the power button on the device to start the device;
5. Press the OK button according to the voice prompt of the device;
6. The mobile terminal enters the "smart screen" to place an order;
7. The specific operation process is as follows: activate the robot → the robot fixes the hole → load the moxa stick → start automatic service → end service.

The meaning of the lights

Panel indicator light color	meaning
-----------------------------	---------

white	Power on
blue	Self-inspection of movement, service completed
	Status normal, idle
green	unlimber
red	equipment failure

Mobile terminal information display

character display		define
task status	waiting	New service tasks are received and waiting to be received and processed
	In service	The automatic service is in continuous service
	It has been suspended	The automated service has been suspended and needs to be addressed
	to be confirmed	The automated service has ended and the robot needs to be released
Service duration	Total duration	Automatically total duration
	Remaining duration	The automatic remaining service duration is changing continuously
Acupuncture and acupoints	moxibustion therapy	It includes hovering moxibustion, sparrow pecking moxibustion, rotary moxibustion, reciprocating moxibustion and meridian moxibustion
device name		Users can customize the device name to distinguish the robot

About battery scrapping

The product comes with a hand controller with a battery inside. When replacing the battery, you must remove the battery from the product. The specific steps are as follows



Step 1: Use a screwdriver to rotate the back cover of the remote control counterclockwise

Step 2: Remove the rear cover and remove the battery



Step 3: Put in the new battery, paying attention to the positive and negative terminals of the battery (face outward)

Step 4: Cover the rear cover and rotate it clockwise to secure it in place

Robot maintenance item

The robot needs to be maintained on time during its use. Maintenance is divided into daily cleaning and wear parts processing.

Ash storage box: a large amount of ash is stored, long-term accumulation may lead to odor generation, ash overflow and other situations;

Tool head: used to measure the distance and temperature of the human body. In case of ash, the moxibustion head may be abnormal and the robot cannot start;

Smoke pipe: the internal accumulation of ash oil, long time without replacement will lead to smoke duct blockage, affecting or causing smoking function failure;

E-cigarette filter: contains a large amount of filtering material, long time without replacement may lead to smoking and odor filtering failure.



position	name	act on	life span	process mode
filter	Pre-filter element	Filter out fine ash	45min	renewal

	Smoke pipe	The Aisian passage	Eighty hours	renewal
	Physical adsorption components	Filter out the smoke	Eighty hours	renewal
Tool head	The head is the main body	The head bearing area	Eight hours	maintenance
	cranidium	The head holds the area	2 months / 8 hours	Replacement/maintenance
	Shield	Where the moxa stick is placed	2 months / 8 hours	Replacement/maintenance
	Smoking cavity		2 months / 8 hours	Replacement/maintenance
	Filter channel	flue	3 days / 20 hours	Change the oil and water
	The visor lens	discriminate	Seven days	Clean once
telecontroller	Remote control battery	Power supply for remote control	Six months	renewal

The names and contents of harmful substances in this product are listed in the following table

Component name	hazardous material					
	lead (Pb)	mercury (Hg)	cadmium (Cd)	hexavalent chromium (Cr)	Polychlorinated biphenyls (PBB)	Polychlorinated dibenzene (PBDE)
hull	○	○	○	○	○	○
Circuit board components	×	○	○	○	○	○
metalwork	○	○	○	○	○	○
Non-metallic parts	○	○	○	○	○	○
Internal cable	○	○	○	○	○	○

This table is prepared in accordance with SJ/T11364

Note: Circuit board components include printed circuit board machine components, such as resistors, capacitors, integrated circuits, connectors, etc.

○: It indicates that the content of the harmful substance in all homogeneous materials of the component is below the limit requirements specified in GB/T26572.

X: indicates that the content of the harmful substance in at least one homogeneous material of the component exceeds the limit requirements specified in GB/T26572.



pay attention to

- ▶ **Do not pull the plug of the power adapter with wet hands.**
- ▶ **Do not rinse the product with water or let water or other liquids into the product.**
- ▶ **Do not use benzene or diluent and other corrosive cleaning agents to clean the product, otherwise it may cause the product to fade, product structure corrosion, cracking.**
- ▶ **The product needs to be cleaned and maintained regularly**

Environmental protection explanation

The equipment smoke purification adsorption consumables are made of environmentally friendly materials, and the adsorption consumables after use need to be buried after solid waste splitting.

Common problems and solutions

Start the fault

Fault information (voice prompt)	fault treatment
Fault occurs, fault code: XXXX. Please press the confirm button to continue	Method 1: Clear the fault according to the prompt and press the confirm key; Method 2: Restart the robot, the fault disappears and continue to work; Method 3: Replace other idle robots, cut off power to the robot and contact

after troubleshooting	technical support;
The robot arm cannot move to the initial position	Method 1: Press the confirm key and try to move back to the initial position; Method 2: Hold down the pause button and drag the robot arm to the initial position; Method 3: Replace other idle robots, cut off power to the robot, and contact technical support;

Fault in service

fault message	fault treatment
The gray box is missing or installed abnormally	Method 1: Adjust the position of the grey box to ensure normal installation, click to clear the fault; Method 2: If method 1 cannot solve the problem, please replace other idle robots, disconnect power to the robot, and contact technical support for abnormal processing;
Robot arm collision detection	Method 1: The robot is impacted or collides with other objects. Please ensure the safety of the working environment, click to clear the fault. If the abnormality disappears, continue to work; Method 2: Restart the robot, the exception disappears and continue to work; Method 3: If method 1 and method 2 cannot be solved, please replace other idle robots, disconnect the power of the robot, and contact technical support for processing;
The robot arm collided with the cabinet	Method 1: Click to clear the fault. If the exception disappears, click to start locating and restoring the robot posture;
Cobot arm collision detection	Method 2: After adjusting the distance between the cabinet and the customer, adjust the fixed position or adjust the height of the part, and continue to work; Method 3: Restart the robot and check whether it is normal. If it is normal, continue to work; Method 4: If the above method cannot solve the problem, please replace other idle robots, disconnect the power of the robot, and contact technical support for processing;
The panel flashes	When WIFI is not connected or the signal is less than -65, the flashing frequency is 1hz; when WIFI is not connected to the store service, the flashing frequency is 3hz;

If there is still an abnormal phenomenon after the above treatment, please contact your local dealer for maintenance service

conformally connected space

Dear customers, thank you for your support. The products you purchased are manufactured under strict quality control and standard inspection process. This warranty card refers to the provision of maintenance and service to customers in accordance with the provisions of the warranty card when the products fail under normal use.

Warranty statement

1. From the date of purchase, if the product fails under normal use, you can enjoy one year's free warranty service, except for damage caused by human destruction (such as improper disassembly of parts, damage treatment, etc.); In addition, damage caused by natural damage (such as fire, earthquake, flood, rodent damage, etc.) is not covered by the free warranty;
 2. The warranty card will not be reissued once it is lost, so please take good care of it;
 3. Free warranty service will not be provided for the faults caused by consumers' improper use;
 4. If the product has failed after the warranty period, if there are spare parts of our company, we will provide warranty service and charge a certain part fee and handling fee;
 5. This warranty card is valid only within the territory of the People's Republic of China;
-

The following situations are not covered by the free warranty

1. Damage caused by prohibited operation not in accordance with the instructions;
2. Counterfeit and inferior products without purchase invoice or receipt, without valid warranty certificate;
3. The outer packaging of the returned products is not complete, the accessories are not complete, or the user disassembled and repaired by himself;
4. Damage caused by human factors: such as falling, water immersion, violent impact;
5. The warranty card is altered or the information on the warranty card does not match the actual product.



保修条例	
产品名称: NAME	尊敬的客户: 感谢您购买本公司的产品。感谢您! 为了保护您的合法权益, 免除您的后顾之忧, 同时为改善对客户的服务, 特制定本保修条例。请您认真阅读并尽快告知您的意见和建议。 一、本公司产品实行自出厂之日起12个月的免费保修服务, 如产品自出厂之日起超过12个月的保修期限, 本公司将视产品故障情况酌情收取零件成本费用。 二、送修产品请妥善保管, 送修过程有损毁或遗失, 本公司恕不负责。 三、在免费保修期内, 发生下列事项者, 本公司有权拒绝服务或酌收材料、服务费。 1、用户操作不当造成产品故障或损害; 2、电压或安装不当导致烧毁之事故; 3、非经服务人员修理造成损坏的; 4、产品标签被破坏者; 四、安装使用本公司产品前请认真阅读本产品说明书。 五、保修卡务必随机携带以便盖章及注明日期, 以确保您的权益。 Dear customer Thank you for purchasing our company's products. Thank you very much! To protect your legitimate rights and interests. To ensure your service throughout the after-sales service for customers, we have specially formulated this insurance policy. Please read the relevant regulations carefully and provide valuable opinions and suggestions. 1. Our company's products enjoy a 12-month free warranty service from the date of sale. Our company will consider if exceed the cost of parts will be charged as appropriate sale, our company will consider if exceed the cost of parts will be charged as appropriate sale, our company will consider if exceed the cost of parts will be charged as appropriate sale. 2. Please pack and transport the returned products properly. If there is any damage or loss during transportation, please note that the company is not responsible. 3. During the free maintenance period, if any of the following events occur, our company has the right to refuse service. Alternatively, material and service fees may be charged as appropriate. 1. Accidents caused by electric shock or improper installation leading to burning. 2. Damage caused by disassembly and repair by non-maintenance service personnel. 4. Damaged product label. 5. Please carefully read this product manual before installing and using our company's products. 6. The warranty card must be stamped and dated by the selling unit to ensure your rights.
产品型号: MODEL	
生产日期: PRODUCT DATE	
售出日期: SELLING DATE	
售出单位: SELLING UNITS (请加盖公章) SEAL	
售后电话: TEL	400-016-8356