

01 Safety precautions and system requirements for use

■ Precautions for safe use:

- To achieve an optimal performance and ensure your pet's safety, please read this manual carefully.
- Please make sure you keep the product out of kid's reach, and be trained before using the otoscope for clinical application.
- It is recommended to use alcohol wipe for disinfection of the outer case of the main unit and accessories. Please do not use oxidizing gases, steam, or any liquid disinfectants. Make sure to avoid presence of flammable anesthetics.
- First, install the software completely onto the phone/tablet, and then connected to the host via wifi.
- Be careful when inserting the probe into the ear canal.
- Please contact Vericare Technology Service support for troubleshooting issues, official accessories purchase, or any after-sales service.

■ System Requirements:

- Android/iOS system on Phone/Pad

02 Product Description and Technical Specifications

- Vericare otoscope system provides a image resolution up to 3840*2160 (according to different size of screen display) with 50X magnification, to deliver an enhanced clarity and efficiency compared to traditional otoscope images.



Power on/Off
Press for 3 Seconds to turn on
Brightness Adjustment



4 High-Intensity LED Lamp Beads



Zoom for magnification



Type-C port for charging (5V-A)

Optional ear tips for replacement :

1#	2#	3#	4#
Size:3 * 22 mm	Size:4 * 22 mm	Size:5.5 * 35 mm	Size:7 * 40 mm

03 Software installation and operation steps

① Software download.

- Please kindly confirm your system version and consult your supplier for a latest QR code for downloading the Smart Vericare.


④ Find the hotspot (Wi-Fi) starting with JM, connect it and display a- /



② Insert the corresponding ear tips



⑤ Register with the basic information, and log in with your registered account.

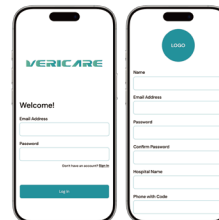


③ Press and hold the power button on the host for 3 seconds, and the light will turn on




⑥ open the app and Select-Wireless otoscope



- 1 Select the correspondent device from Vericare, to start the scanning for images or videos.

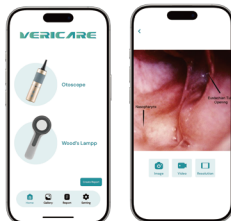


- 2 Create a report for your patient or re-edit on the previous patient record.

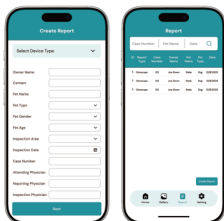


- 3 Select the images from your mobile album to complete the report.



- 4 Diagnostic options allows for template preset and enhancing your diagnostic efficiency

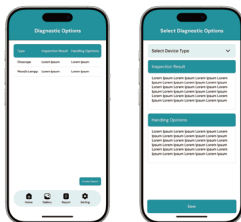

- 5 Draft and preview your report before sharing with your vet or pet-owners

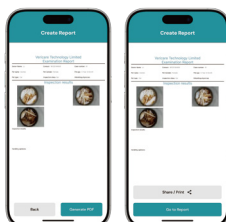

Vericare Veterinary Wireless Video Otoscopy & Smart Vericare App User Manual



Safety precautions and system requirements for use	01
Product Description and Technical Specifications	02
Software installation and operation steps	03

FCC Warning:

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:-- Reorient or relocate the receiving antenna.

-- Increase the separation between the equipment and receiver.

-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-- Consult the dealer or an experienced radio/TV technician for help.

Warning: changes or modifications not expressly approved by the party responsible for compliance could

void the user's authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement, the device can be used in

portable exposure condition without restriction.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference

received, including interference that may cause undesired operation.