



3.2 Wear

Personal Sound Amplifier L/R Recognition



3.3 Power on/off

Power on: After the charging case is closed and the Personal Sound Amplifier is fully charged, turn on, or long press the button for 3 seconds to turn on.

Power off: Put the Personal Sound Amplifier into the charging case, close the charging case, and the Personal Sound Amplifier will be powered off automatically.

3.4 Volume adjustment


Press the volume button each time, the volume increases one level. If it is volume at maximum value, press the button again, the volume will be reset to minimum value.

4.2 Via Hear APP operating guide

(1) Connect Bluetooth



(7) Acoustic environment checking (If no noisy, the test can not be carried out)

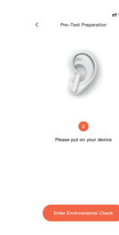

1. Attachment

Personal Sound Amplifier x 2

Charging case x 1

Charging cable x 1

Ear cap x 4

User manual x 1

3.5 Pairing

(1) When you first open the charging case of the Personal Sound Amplifier, the device will automatically enter the Bluetooth pairing mode.

(2) Press the Bluetooth pairing button on the back of the Personal Sound Amplifier for 3 seconds to enter the pairing mode.

(3) Press the Bluetooth pairing button on the back of the Personal Sound Amplifier for 3 seconds to enter the pairing mode.

(4) Press the Bluetooth pairing button on the back of the Personal Sound Amplifier for 3 seconds to enter the pairing mode.

3.6 Replace the ear cap

(1) Hold your fingers on the hearing aid body and pull out the ear cap.

(2) Push the indicated ear cap slightly to ensure that it is fully inserted.

(3) Repeat the steps above to replace the other ear cap.

(4) Repeat the steps above to replace the other ear cap.

(5) Repeat the steps above to replace the other ear cap.

(6) Repeat the steps above to replace the other ear cap.

(7) Repeat the steps above to replace the other ear cap.

(8) Repeat the steps above to replace the other ear cap.

(9) Repeat the steps above to replace the other ear cap.

(10) Repeat the steps above to replace the other ear cap.

(11) Repeat the steps above to replace the other ear cap.

(12) Repeat the steps above to replace the other ear cap.

(13) Repeat the steps above to replace the other ear cap.

(14) Repeat the steps above to replace the other ear cap.

(15) Repeat the steps above to replace the other ear cap.

(16) Repeat the steps above to replace the other ear cap.

(17) Repeat the steps above to replace the other ear cap.

(18) Repeat the steps above to replace the other ear cap.

(19) Repeat the steps above to replace the other ear cap.

(20) Repeat the steps above to replace the other ear cap.

(21) Repeat the steps above to replace the other ear cap.

(22) Repeat the steps above to replace the other ear cap.

(23) Repeat the steps above to replace the other ear cap.

(24) Repeat the steps above to replace the other ear cap.

2. Product structure

2.1 Personal Sound Amplifier

2.2 Charging case

2.3 Type-C charging port

2.4 Troubleshooting guide

2.5 Technical parameters

2.6 APP guide

2.7 Product safety information

2.8 Product safety information

2.9 Product safety information

2.10 Product safety information

2.11 Product safety information

2.12 Product safety information

2.13 Product safety information

2.14 Product safety information

2.15 Product safety information

2.16 Product safety information

2.17 Product safety information

2.18 Product safety information

2.19 Product safety information

2.20 Product safety information

2.21 Product safety information

2.22 Product safety information

2.23 Product safety information

2.24 Product safety information

2.25 Product safety information

2.26 Product safety information

2.27 Product safety information

2.28 Product safety information

2.29 Product safety information

2.30 Product safety information

2.31 Product safety information

2.32 Product safety information

2.33 Product safety information

2.34 Product safety information

2.35 Product safety information

2.36 Product safety information

2.37 Product safety information

2.38 Product safety information

2.39 Product safety information

2.40 Product safety information

3. Guide to operation

3.1 How to charge the Personal Sound Amplifier and charging case

3.2 How to charge the Personal Sound Amplifier

3.3 How to charge the Personal Sound Amplifier

3.4 How to charge the Personal Sound Amplifier

3.5 How to charge the Personal Sound Amplifier

3.6 How to charge the Personal Sound Amplifier

3.7 How to charge the Personal Sound Amplifier

3.8 How to charge the Personal Sound Amplifier

3.9 How to charge the Personal Sound Amplifier

3.10 How to charge the Personal Sound Amplifier

3.11 How to charge the Personal Sound Amplifier

3.12 How to charge the Personal Sound Amplifier

3.13 How to charge the Personal Sound Amplifier

3.14 How to charge the Personal Sound Amplifier

3.15 How to charge the Personal Sound Amplifier

3.16 How to charge the Personal Sound Amplifier

3.17 How to charge the Personal Sound Amplifier

3.18 How to charge the Personal Sound Amplifier

3.19 How to charge the Personal Sound Amplifier

3.20 How to charge the Personal Sound Amplifier

3.21 How to charge the Personal Sound Amplifier

3.22 How to charge the Personal Sound Amplifier

3.23 How to charge the Personal Sound Amplifier

3.24 How to charge the Personal Sound Amplifier

3.25 How to charge the Personal Sound Amplifier

3.26 How to charge the Personal Sound Amplifier

3.27 How to charge the Personal Sound Amplifier

3.28 How to charge the Personal Sound Amplifier

3.29 How to charge the Personal Sound Amplifier

3.30 How to charge the Personal Sound Amplifier

3.31 How to charge the Personal Sound Amplifier

3.32 How to charge the Personal Sound Amplifier

3.33 How to charge the Personal Sound Amplifier

3.34 How to charge the Personal Sound Amplifier

3.35 How to charge the Personal Sound Amplifier

3.36 How to charge the Personal Sound Amplifier

3.37 How to charge the Personal Sound Amplifier

3.38 How to charge the Personal Sound Amplifier

3.39 How to charge the Personal Sound Amplifier

3.40 How to charge the Personal Sound Amplifier

4. Troubleshooting guide

4.1 Problem

4.2 Solution

4.3 Solution

4.4 Solution

4.5 Solution

4.6 Solution

4.7 Solution

4.8 Solution

4.9 Solution

4.10 Solution

4.11 Solution

4.12 Solution

4.13 Solution

4.14 Solution

4.15 Solution

4.16 Solution

4.17 Solution

4.18 Solution

4.19 Solution

4.20 Solution

4.21 Solution

4.22 Solution

4.23 Solution

4.24 Solution

4.25 Solution

4.26 Solution

4.27 Solution

4.28 Solution

4.29 Solution

4.30 Solution

4.31 Solution

4.32 Solution

4.33 Solution

4.34 Solution

4.35 Solution

4.36 Solution

4.37 Solution

4.38 Solution

4.39 Solution

4.40 Solution

5. Technical parameters

5.1 Acoustic parameters

5.2 Electrical parameters

5.3 Mechanical parameters

5.4 Environmental parameters

5.5 Safety parameters

5.6 Safety parameters

5.7 Safety parameters

5.8 Safety parameters

5.9 Safety parameters

5.10 Safety parameters

5.11 Safety parameters

5.12 Safety parameters

5.13 Safety parameters

5.14 Safety parameters

5.15 Safety parameters

5.16 Safety parameters

5.17 Safety parameters

5.18 Safety parameters

5.19 Safety parameters

5.20 Safety parameters

5.21 Safety parameters

5.22 Safety parameters

5.23 Safety parameters

5.24 Safety parameters

5.25 Safety parameters

5.26 Safety parameters

5.27 Safety parameters

5.28 Safety parameters

5.29 Safety parameters

5.30 Safety parameters

5.31 Safety parameters

5.32 Safety parameters

5.33 Safety parameters

5.34 Safety parameters

5.35 Safety parameters

5.36 Safety parameters

5.37 Safety parameters

5.38 Safety parameters

5.39 Safety parameters

5.40 Safety parameters

6. APP guide

6.1 Download the Via Hear APP

6.2 Open the Via Hear APP

6.3 Connect the Via Hear APP

6.4 Adjust the volume

6.5 Adjust the volume

6.6 Adjust the volume

6.7 Adjust the volume

6.8 Adjust the volume

6.9 Adjust the volume

6.10 Adjust the volume

6.11 Adjust the volume

6.12 Adjust the volume

6.13 Adjust the volume

6.14 Adjust the volume

6.15 Adjust the volume

6.16 Adjust the volume

6.17 Adjust the volume

6.18 Adjust the volume

6.19 Adjust the volume

6.20 Adjust the volume

6.21 Adjust the volume

6.22 Adjust the volume

6.23 Adjust the volume

6.24 Adjust the volume

6.25 Adjust the volume

6.26 Adjust the volume

6.27 Adjust the volume

6.28 Adjust the volume

6.29 Adjust the volume

6.30 Adjust the volume

6.31 Adjust the volume

6.32 Adjust the volume

6.33 Adjust the volume

6.34 Adjust the volume

6.35 Adjust the volume

6.36 Adjust the volume

6.37 Adjust the volume

6.38 Adjust the volume

6.39 Adjust the volume

6.40 Adjust the volume

7. Product safety information

7.1 Read the instructions carefully before using the product

7.2 Do not use the product in high temperature, strong light and strong magnetic field environments

7.3 Do not use the product in high temperature, strong light and strong magnetic field environments

7.4 Do not use the product in high temperature, strong light and strong magnetic field environments

7.5 Do not use the product in high temperature, strong light and strong magnetic field environments

7.6 Do not use the product in high temperature, strong light and strong magnetic field environments

7.7 Do not use the product in high temperature, strong light and strong magnetic field environments

7.8 Do not use the product in high temperature, strong light and strong magnetic field environments

7.9 Do not use the product in high temperature, strong light and strong magnetic field environments

7.10 Do not use the product in high temperature, strong light and strong magnetic field environments

7.11 Do not use the product in high temperature, strong light and strong magnetic field environments

7.12 Do not use the product in high temperature, strong light and strong magnetic field environments

7.13 Do not use the product in high temperature, strong light and strong magnetic field environments

7.14 Do not use the product in high temperature, strong light and strong magnetic field environments

7.15 Do not use the product in high temperature, strong light and strong magnetic field environments

7.16 Do not use the product in high temperature, strong light and strong magnetic field environments

7.17 Do not use the product in high temperature, strong light and strong magnetic field environments

7.18 Do not use the product in high temperature, strong light and strong magnetic field environments

7.19 Do not use the product in high temperature, strong light and strong magnetic field environments

7.20 Do not use the product in high temperature, strong light and strong magnetic field environments

7.21 Do not use the product in high temperature, strong light and strong magnetic field environments

7.22 Do not

FCC Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

RF Exposure Information

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.