


尺寸: 70*100mm 157g 单铜 折页

adomii

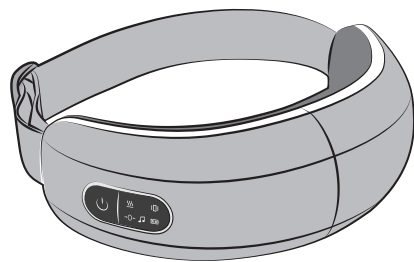

Wireless Eye Massager

User Manual

AD-MASG-EYE-WH

Table of Contents

1. Package Contents	02
2. Product Overview	03
3. Operation Guide	04
4. How to Charge	04
5. Bluetooth Pairing	04
6. Precautions	05
7. Product Specification	06
8. Troubleshooting	07


Thank you for purchasing your new adomii Wireless Eye Massager. This user manual will help provide you with all the information needed to get started. Please read this manual carefully before using the product and reference it for later use.

Enjoy your new product.

1. Package Contents



Massager

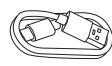


USB Charge Cable



User Manual

-2-

2. Product Overview

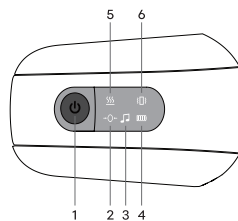

- 1. Power Button**
Long Press ON/OFF
Short Press Switches Modes
- 2. Air Compression Indicator**
- 3. Audio Mode Indicator**
Music (White Light)
Bluetooth (Blue Light)
- 4. Battery Indicator**
Flashing: Low Power
Breathing: Charging
Solid: Fully Charged
- 5. Heating Indicator**
- 6. Vibration Indicator**

Music Audio Mode

Mode	Air Compression	Heating	Vibration	Music	Timer (Minutes)
Vitality	✓	✓	✓	✓	10
Relaxation	✓	✓		✓	15
Surfing	✓	✓			15
Sleep	✓			✓	20

Bluetooth Audio Mode

Mode	Air Compression	Heating	Vibration	Bluetooth	Timer (Minutes)
Vitality	✓	✓	✓	✓	10
Relaxation	✓	✓		✓	15
Surfing	✓	✓		✓	15
Sleep	✓			✓	20


-3-

3. Operation Guide

- Long press the power button to power on the Eye Massager.
- You will be greeted by a voice prompt "Power On".
- Vitality Mode begins by default.
- Short press the power button to cycle through the 4 modes:
Vitality → Relaxation → Surfing → Sleep
- When the battery power is low, the Battery LED will flash, and the beep sound will chime. When the battery power is exhausted, the massager will auto shut down, requiring a charge before usage.

4. How to Charge

- Connect the charge cable to a power adapter (not included).
- Connect to the USB-C port on the bottom of the Eye Massager

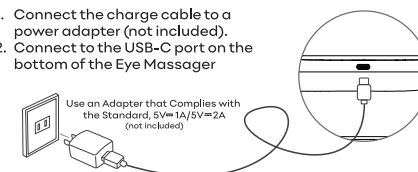

5. Bluetooth Pairing

- With the Eye Massager powered off, long press the power button to power on.
- Enable Bluetooth on your device.
- Search for "adomii Eye Massager" and tap to connect the Eye Massager with your device.
- After connecting via Bluetooth successfully, the Audio Mode Indicator on the Eye Massager will change to blue.



-4-

adomii Limited Warranty

Warranty Coverage:

We stand behind our products with confidence, offering warranties designed to give you peace of mind and a worry-free experience. Our commitment covers most adomii products with specific warranty periods based on product type such as cables, cases, wall chargers, and more. Either way, you're protected from defects in materials and workmanship.

• Consumer Electronics: 24 Months

30-Day Money-Back Guarantee:

You may return undamaged products within 30 days for a full refund. Products must include all accessories and original packaging. Shipping costs for non-quality-related returns are the responsibility of the buyer.

Replacement Warranty:

If a defect arises within the warranty period, adomii will repair, replace, or issue a credit for the product. Repaired or replacement products are covered for 90 days or the remainder of the original warranty, whichever is longer.

Need Assistance?

For more information or to file a warranty claim, please visit adomii's Warranty Support. (<https://www.adomii.com/pages/warranty-information>)
(866) 686-1800
support@adomii.com
Mon - Fri: 9:00 AM - 5:30 PM (PST)
www.adomii.com

FCC Compliance Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference.

(2) This device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 10cm between the radiator and your body.



WARNING: This product can expose you to chemicals known to the State of California to cause cancer.
For more information, go to www.P65Warnings.ca.gov.

-5-

7. Product Specification

Charging Voltage	5V=1A
Power	5W
Battery Capacity	1200mAh

-6-

6. Precautions

Attention: This massager is designed for private use, not intended for medical or commercial purposes. Children should be supervised to ensure that they do not use this massager as a toy.

- Adjust elastic band to a comfortable setting
- Use while seated, leaning against a sofa or chair for comfort and stability
- Suggested use of this massager is once or twice per day, with each time being 10-20 minutes
- Do not open your eyes when using this massager
- Do not use this massager while charging
- It is recommended to use the included charging cable and standard adapter (not included) for charging
- Do not clean the massager with water
- Please clean surface of this massager with a damp towel or cotton ball with 75% medical alcohol.
- Do not put this device close fire, corrosive liquid, or other adverse environment.
- When using it for the first time, or the battery power is depleted after the product has been shelved for a long time, please follow the instruction manual to charge it to ensure that it can work properly
- Do not use with wet face/eyes
- Not recommended for children under 7 years old

Please consult with a doctor before using the Eye Massager if you:

- Are pregnant
- Are shortsighted by over 800 degrees
- Are experiencing glaucoma or cataract
- Are sequela of apoplexy
- Are experiencing an eye injury, or in the period of rehabilitation of finished eye surgery (double eyelid)
- Had glaucoma, cataract, retinal detachment, and other serious ocular region surgery
- Are suffering from serious eye inflammation

8. Troubleshooting

If you encounter problems when using the product, please refer to the following tips and solutions. If the product is still in a malfunctioning state or cannot be used, please turn off the power. Users should not disassemble or repair the product by themselves. Damage to the product caused by disassembly or repair will invalidate the product warranty, and the company will not be responsible for any damage caused by the user.

Problem	Possible Solutions
Unable to Start Massager	Please confirm that the battery is charged; please confirm that the ON/OFF button has been pressed.
Unable to Charge or No Charge Lamp Display	Please check whether the USB-C charging cable is connected to the massager's charging port.
The Massager Suddenly Stops During Operation	The battery is low; Please charge to 100%. The time set by the automatic timing device has passed.

Danger	The battery used in this product is a special rechargeable battery and cannot be used in other products. Do not charge the rechargeable battery that has been removed; Do not throw it into fire to cause excessive heat; Do not pierce it with nails, impact it with gravity, or disassemble it for modification; Do not directly connect the positive and negative poles of the battery with metal; Do not store it with necklaces, hairpins, etc. Do not charge, use, or place it near a stove or in the hot sun. The above situations may cause heating, fire, rupture, etc.
Warning	Please place the removed rechargeable battery out of reach of infants and young children to prevent accidental swallowing; If accidental swallowing occurs, please seek medical attention immediately; When the battery fluid leaks, do not touch the liquid with bare hands, and handle it as follows: When the battery fluid enters the eyes, there is a risk of blindness. Do not rub the eyes directly with your hands. Immediately wash with clean water and seek medical attention; When the battery fluid gets on the body or clothing, it may cause skin inflammation or injury. Please rinse with clean water. If you still feel uncomfortable, seek medical attention immediately. The product can only be disassembled when it is discarded. Failure to do so may result in fire, electric shock, injury, etc.

-7-