

NoctIQ™
Calmer Nights, Brighter Days



User Manual
Smart Anti-Snore Pillow
Powered by nitetronic 

WELCOME

Thank you for choosing the
NoctIQ Smart Anti-Snore Pillow.

Need Help?

The **NoctIQ Care Team** is ready to
support your sleep journey.

Support



Live Chat: noctiq.com

Phone: +1-888-598-7231 | Hours: 12PM-8PM EST

Email: care@noctiq.com

NoctIQTM
Calmer Nights, Brighter Days

You're about to begin a simple journey toward quieter nights, deeper rest, and more connected mornings.

The **NoctIQ Smart Anti-Snore Pillow** is designed to work with your natural sleep patterns over time. As it learns, adapts, and responds, you'll begin to experience a more personalized path to better sleep.

Because sleep is personal, your NoctIQ experience begins with a **21-day sleep calibration** that helps the AI learn your natural snoring patterns, gives your body time to adapt, and refines **Snore Intervention** for your comfort and results.

Your 21-Day Sleep Calibration

Week 1: Establish Your Baseline - App Only

For the first 7 nights, continue using your current pillow. Use the **NoctIQ app** to track your snoring patterns. This allows our AI to learn your natural baseline so **Snore Intervention** can be better personalized once it begins.

Weeks 2-3: Adaptation Period - Pillow + App

Switch to your **NoctIQ Smart Anti-Snore Pillow** and activate Snore Intervention. During this period, the pillow gently encourages side head movement while the AI personalizes its response. Most users naturally adapt to the pillow's gentle movement within about **14 days**.

Be patient with the process. By allowing this three-week calibration, you are building a foundation for years of better rest.

If you need help along the way, the **NoctIQ Care Team** is ready to support your sleep journey.

IMPORTANT SAFETY INSTRUCTIONS

READ BEFORE USE

Read all instructions before using this product. Failure to follow these instructions may result in damage to the product or risk of injury. The **NoctIQ Smart Anti-Snore Pillow, powered by Nitetronic**, contains precision components and an integrated Wi-Fi module. Handle with care.

General Safety

- **Handling:** The NoctIQ Smart Anti-Snore Pillow contains precision sensors and integrated electronics. Do **NOT** excessively squeeze, fold, bend, or strike the pillow.
- **Intended Use:** Use this product only for its intended use as described in this manual. Do **NOT** use attachments not recommended by the manufacturer.
- **Environment:** Keep away from heat sources such as electric blankets, heating pads, space heaters, or any appliance that may cause excessive heat.
- **Environment & Moisture:** Keep the pillow in a dry, indoor environment. Do **NOT** expose the unit to moisture, liquids, corrosive gases, or excessive dust. (See **Electrical Safety** for water-related hazards).
- **Damage & Malfunction:** If the pillow is not functioning properly, emits an unusual odor, or appears damaged (including tears in the inner casing or exposed wires), stop use immediately. Disconnect the power and contact the NoctIQ Care Team. Using a damaged electronic product may result in fire or electric shock.
- **Modification:** Do **NOT** attempt to disassemble, repair, or modify the product. There are no user-serviceable parts inside. Unauthorized tampering voids the warranty and may create a safety hazard.

Electrical Safety

- **Power Supply:** Use only the original power cable and the **Shenzhen Sunshine Technological Co., Ltd. / XSC-0502000AUSU** power adapter included with the product.
- **Certified Power Sources:** Use only with a USB port or power supply that is certified for Class 2 operation, such as the included power adapter.
- **Foreign Objects:** Do **NOT** insert any objects into the air tube, air ports, or any openings in the control unit. Doing so may cause mechanical failure, puncture the internal air bladder, or result in electric shock.
- **Cord Care:** Do **NOT** pull the power cord forcefully. This may damage the cord and create an electrical hazard. Position cords to avoid tripping hazards.
- **Damage Protocol:** If the power cord or adapter is damaged, disconnect power immediately and stop use. Contact the NoctIQ Care Team for a manufacturer approved replacement.
- **Accessibility:** Plug the power adapter into an easily accessible outlet so it can be quickly disconnected in an emergency.
- **Prohibitions:** Do **NOT** use the product near water (e.g., bathtubs or sinks) or in environments where it may become wet. Do **NOT** obstruct internal ventilation by placing the control unit under heavy bedding or in tightly enclosed spaces.

Children & Pets

- This product is not a toy. Keep away from infants, young children, and pets.
- Do **NOT** allow children to play with or sleep on the product unsupervised.

User Capability Requirements

This product should be used only by individuals who can safely operate household electronic devices. Users with reduced physical, sensory, or cognitive capabilities should use the product only with supervision or after receiving instruction from a person responsible for their safety.

Unsuitable Users

Do **NOT** use this product if you:

- Have intracranial electrodes
- Have head or neck injuries
- Have an implanted pacemaker or other implanted electronic medical device

Users with **severe sleep apnea** should consult a physician before using this product.

Non-Medical Device Disclaimer

The **NoctIQ Smart Anti-Snore Pillow, powered by Nitetronic**, is a General Wellness device designed to help reduce snoring disturbances by gently encouraging side head movement during sleep.

- This product is **NOT** a medical device.
- It does **NOT** diagnose, treat, cure, or prevent sleep apnea or any medical condition.
- It is **NOT** intended for medical monitoring and is not a substitute for professional medical evaluation.

Care & Maintenance

Proper care ensures the longevity of your smart components.

- **Outer Pillowcase:** The outer pillowcase is removable and washable. Please follow the washing instructions on the care label. The outer pillowcase must be installed during use.
- **Inner Pillowcase:** The inner pillowcase is not removable. Do **NOT** attempt to open or remove it.
- **Cleaning the Unit:** The entire pillow unit and internal components **CANNOT** be washed or dry-cleaned. To clean, wipe with a dry or slightly damp cloth only.
- **Storage:** When not in use, store the pillow in a cool, dry place. Do **NOT** place heavy objects on top of the pillow during storage to avoid damaging the air-inflation system.

Specifications

Type	Electrical Rating (Product)	Power Adapter Rating
• ASP04	• Input: 5V $\equiv$ 2A • Power: 10W Max	• Input: 100-240V ~ 50/60 Hz 0.4 A Max • Output: 5V $\equiv$ 2A

INSIDE THE BOX



Power Adapter
Provides power to the control unit.



USB Power Cable
Connects the power adapter to the control unit.



Smart Anti-Snore Pillow
Adaptive Inflation Technology that helps reduce snoring.



Cooling Pillowcase
Soft, breathable fabric that helps keep you cool and comfortable.



Height Adjustment Pad
Add or remove to adjust pillow height.



User Manual
Detailed instructions for setup, use, care, and safety.

Before First Use

- ☒ Remove all packaging
- ☒ Air out components if desired
- ☒ Optional: Wash pillowcase for a fresher feel

Getting Started Guide

App Setup & Account Registration

Before first use, you must set up the **NoctIQ: Anti-Snore Pillow App**. The app works together with your Smart Anti-Snore Pillow to detect snoring, track your sleep, and deliver personalized snore reduction.

Step 1: Download the NoctIQ: Anti-Snore Pillow App

Scan the QR code with your phone to download the app. You'll be directed to the App Store or Google Play based on your device.



Scan to download

Scan QR code or search **NoctIQ: Anti-Snore Pillow** in the App Store or Google Play

Step 2: Create Your Account

- 1 Open the NoctIQ app and tap **Signup**.



- 2 Enter your **Personal Details**—helps personalize your recommended starting settings.



- 3 Tap **Continue**.



- 4 Check your email for the **6-digit verification code**.



5

Enter the verification code in the app and tap **Verify**.



6

Create and confirm your password.



7

Accept the **Terms & Conditions** and **Privacy Policy**.



8

Tap **Create Account**.



Step 3: Phone Settings for Accurate Tracking

To ensure uninterrupted tracking:

iPhone:

- **Low Power Mode:** OFF
- **Background App Refresh:** ON
- **Focus / Sleep:** Allow NoctIQ app notifications

Android:

- **Battery Saver:** OFF
- **Battery Usage:** Settings → Apps → **NoctIQ app** → Battery → Unrestricted

Optional for better accuracy:

- **Sound Recognition / Sound Detection:** OFF
- **Bluetooth Audio Devices:** Disconnect headphones or speakers before tracking, as they may interfere with snore detection.

⚠️WARNING

Turning off Sound Recognition / Sound Detection may prevent your phone from alerting you to safety alarms, including smoke detectors, fire alarms, or other emergency notifications. If you rely on these alerts, keep Sound Recognition / Sound Detection enabled. The NoctIQ Smart Anti-Snore Pillow is not a medical, safety, or emergency alert device.


Week 1: Baseline Tracking – Without Snore Intervention

For the first 7 nights, use the **NoctIQ app** with your **current pillow** before switching to the **Smart Anti-Snore Pillow**. This helps create a more accurate personal baseline and gives the AI time to learn your natural snoring patterns before Snore Intervention begins.

If you do not snore on some nights, you may continue tracking until you have 7 nights with snoring data.

Phone Placement & Power Guidelines

For accurate detection, place your phone as shown below.

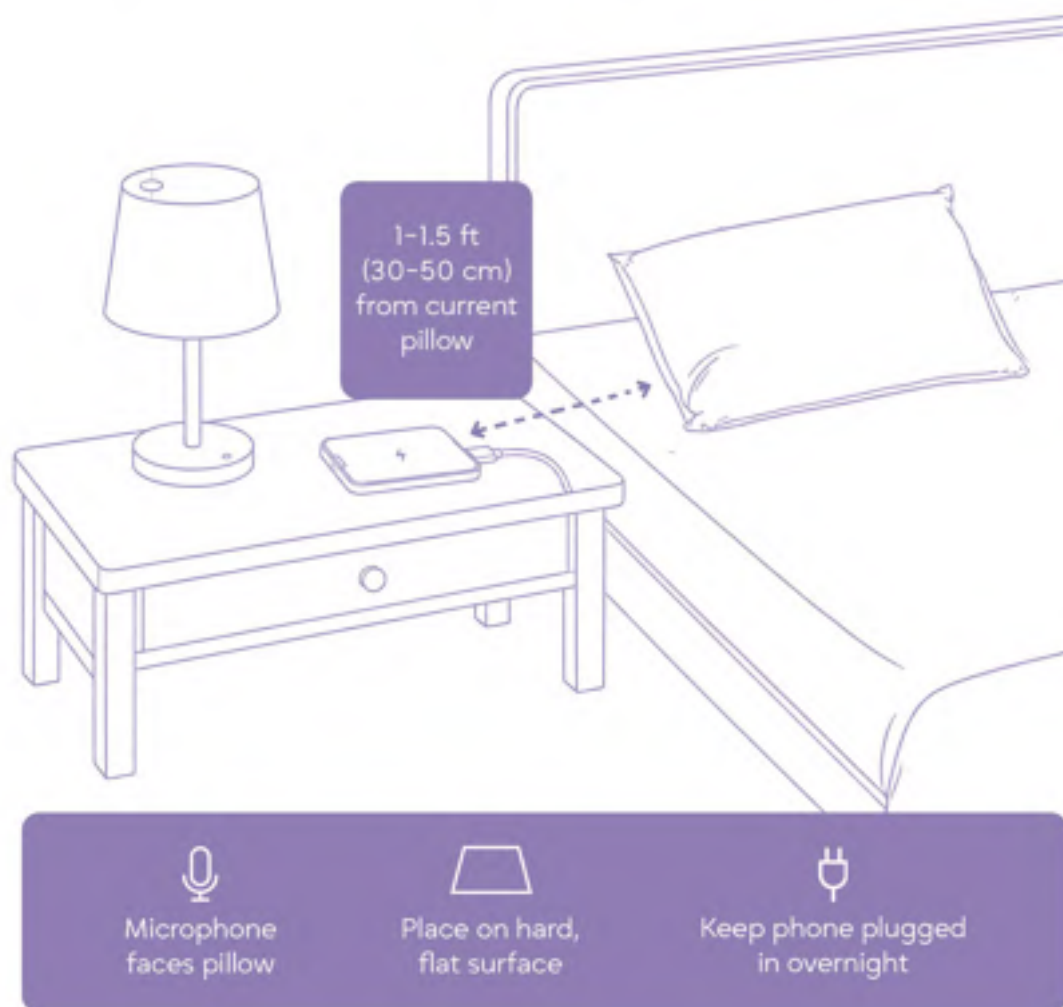

⚠️WARNING

Placing your phone on the bed, under a blanket, or under a pillow may block the microphone or cause the device to overheat. Place the phone on a hard, flat surface during use.

Nightly Routine – App Only

- 1 Open the **NoctIQ** app.




- 2 Tap **Start Tracking**.



- 3 Place phone on **nightstand** near pillow.



- 4 Allow your **phone screen to turn off** automatically, or press the **power button** to turn off the screen / sleep normally.

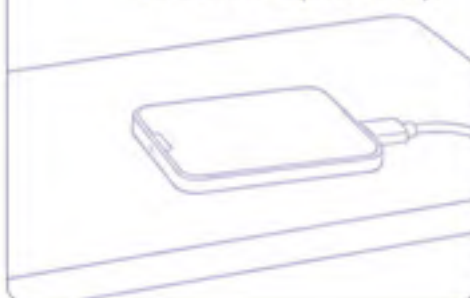

- 5 In the Morning, tap **Stop Tracking**.



- 6 Review **Snore Report**. See **Understanding Your Snore Report**, page 22.



Privacy Note

Microphone access is used only for snoring detection during active tracking sessions. Snoring recordings stay on your device and are not shared.

Important

After tapping Start Tracking, do **NOT** use other audio-playing apps or sleep-detection software during the night, as this may affect snoring detection accuracy.


Week 2+: Snore Reduction & Adaptation – Snore Intervention ON


After completing baseline tracking, set up your Smart Anti-Snore Pillow and activate Snore Intervention.

This begins a 2-week adaptation period as your body adapts to gentle pillow movement and the AI personalizes Snore Intervention over time.

Step 1: Physical Setup & Cable Routing

Follow these steps to prepare your Smart Anti-Snore Pillow for setup.

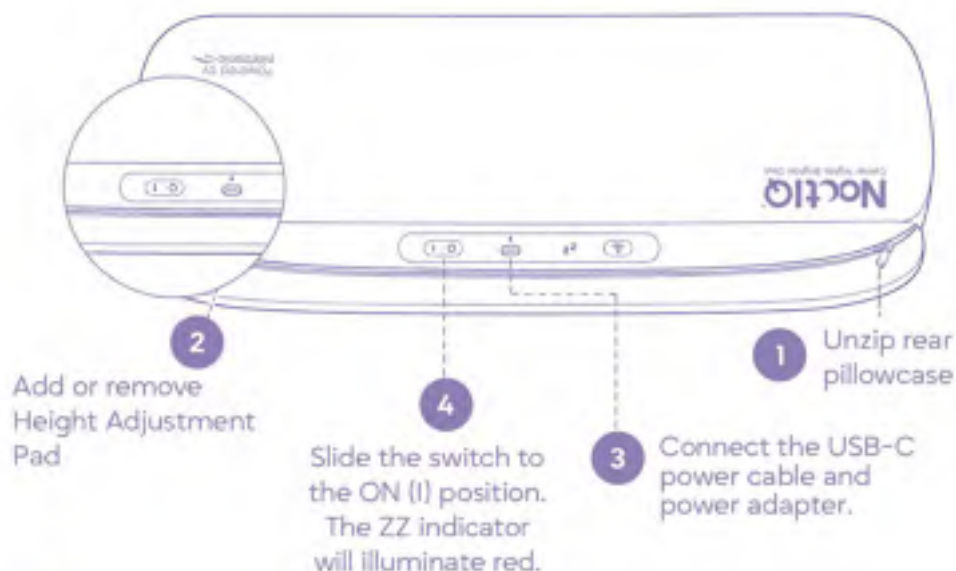

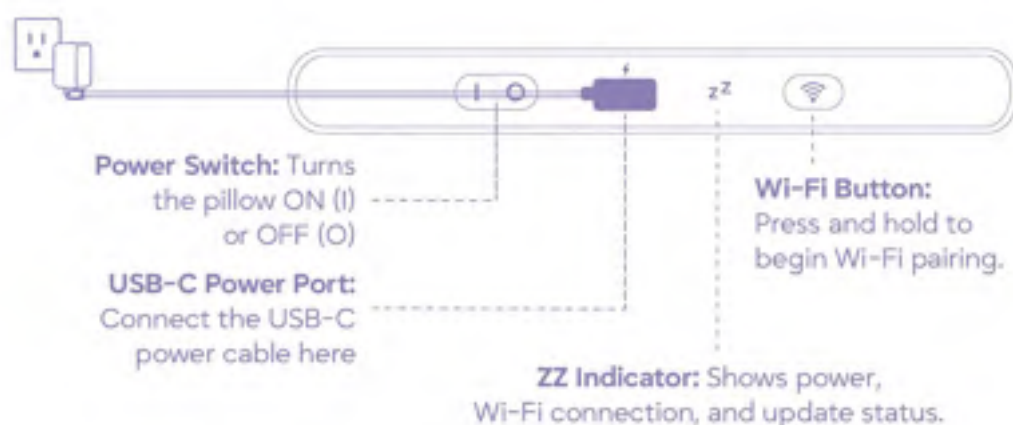

ZZ Indicator Light Guide

ZZ Indicator	Status	What It Means
Red	Solid	Pillow is powered on but not connected to Wi-Fi.
Blue	Fast blinking	Connecting to Wi-Fi.
Blue	Slow blinking	Firmware update in progress.
Blue	Solid	Connected to Wi-Fi.

Notice

The Smart Anti-Snore Pillow **does not contain a battery**. It must remain **connected to power during overnight use**.

Cable Management

Route the cable properly to protect the power connection.

1. Ensure the L-shaped plug is fully inserted into the control panel.
2. Rotate the plug so the cable lies flat and points toward the outlet.
3. Guide the cable through the nearest exit notch.
4. Zip the pillowcase closed with the cable exiting cleanly.

L-shaped plug

Connects securely to the control panel and helps the cable lie flat.




Exit notch

Allows the cable to exit cleanly from the pillowcase corner.



Cable routing direction

Route the cable toward the nearest wall outlet and away from the sleeping area.

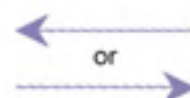

Step 2: Connect Your Pillow to Wi-Fi


Wi-Fi Pairing

To use your Smart Anti-Snore Pillow, you must pair it with the **NoctIQ App**. This connection enables AI features, sleep tracking, and smart, personalized Snore Intervention.

- 1 Open the app → Profile → My Pillow → Add Pillow.



- 2 Press and hold the Wi-Fi button until zZ flashes blue.

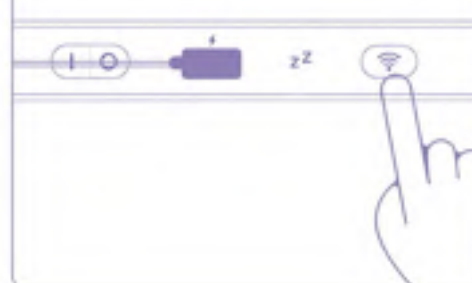

- 3 Check the box "I see the blue light flashing."



- 4 Tap Continue.



- 5 Enter Wi-Fi password.



- 6 Allow setup to reach 100%.



- 7 Confirm your pillow has successfully paired.



Permission Note: If prompted, allow required permissions so the app can find and connect to your pillow.

Connection Tip

If the app shows **We Couldn't Connect Your Pillow** or the ZZ indicator does **NOT** turn solid blue, follow the Wi-Fi Pairing Troubleshooting steps below.

Wi-Fi Pairing Troubleshooting

If your pillow does not connect:

- **Check your password**
Confirm the Wi-Fi password is correct, then retry pairing.
- **Use a 2.4 GHz network**
The pillow requires **2.4 GHz Wi-Fi**. 5 GHz networks are not supported.
- **Use your primary router during setup**
If you use a mesh system or Wi-Fi extender, **connect your phone to the primary 2.4 GHz network** before pairing. **Avoid extender or 5 GHz network names** such as **_EXT** or **_5G**.


Note


After setup is complete and the ZZ indicator is solid blue, you may reconnect your phone to your preferred network.

Step 3: Confirm Proper Setup

Pillow Self-Test

Confirm the pillow is working properly before your first night.

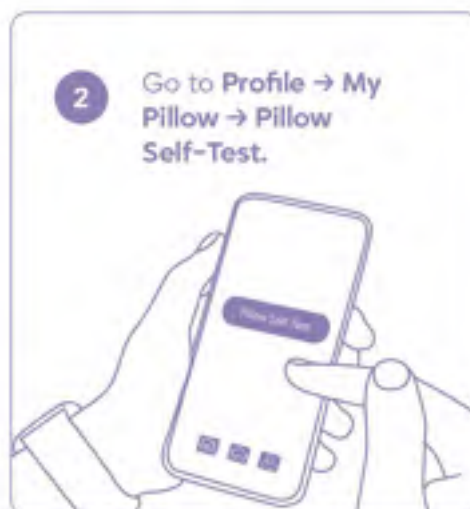

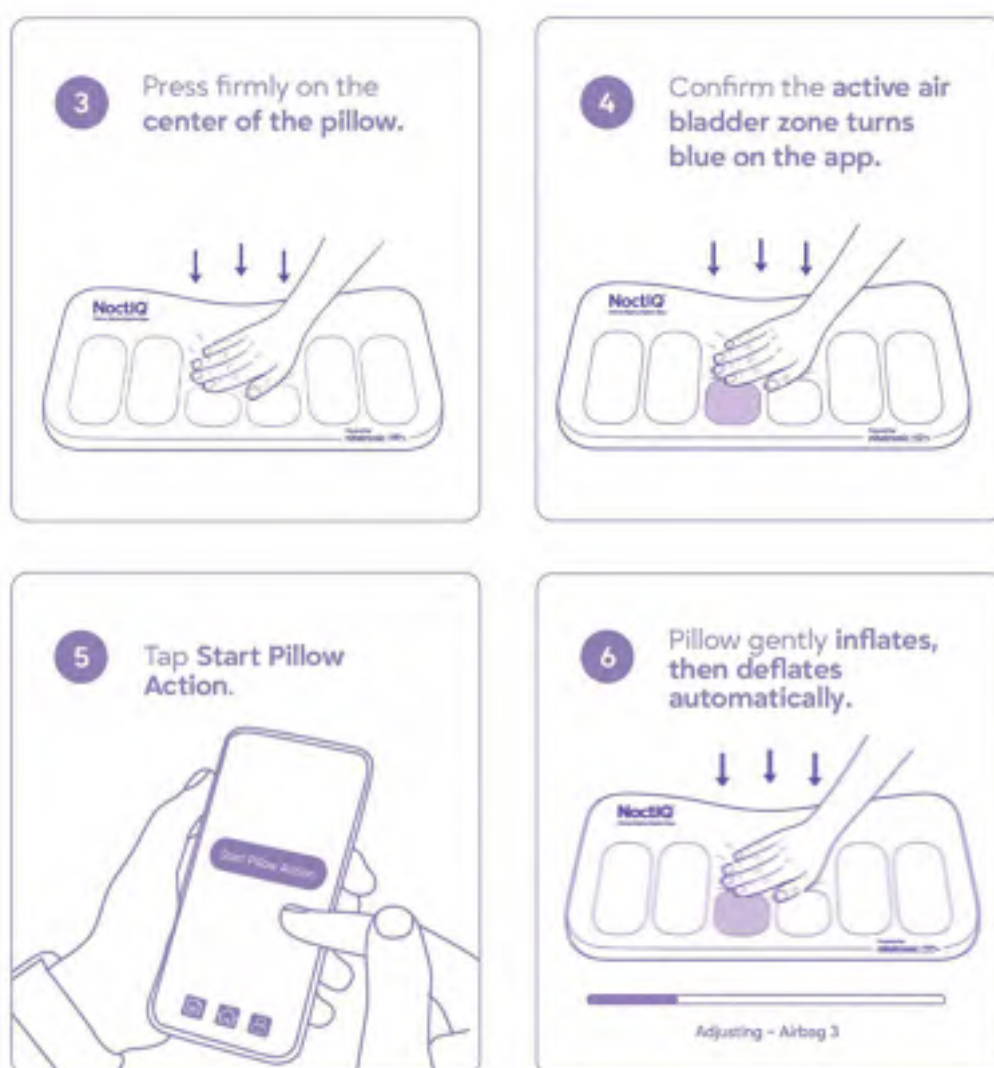

Reset Head Position Sensor (If Needed)

If Snore Intervention is ON but no intervention reports are being generated, the head position sensor may need to be reset.


1. Place the **Smart Anti-Snore Pillow** in its intended position on the bed.
2. Turn the **pillow ON (I)**.
3. In the app, go to **Profile → My Pillow → Reset Sensor**.
4. Wait approximately **1 minute** for the sensor to reset automatically.

Important

Do **NOT** lie on the pillow or place any objects on it during sensor reset.


Step 4: Activate Snore Intervention & Confirm Settings

Follow the steps below to turn **Snore Intervention ON** and confirm your starting settings.

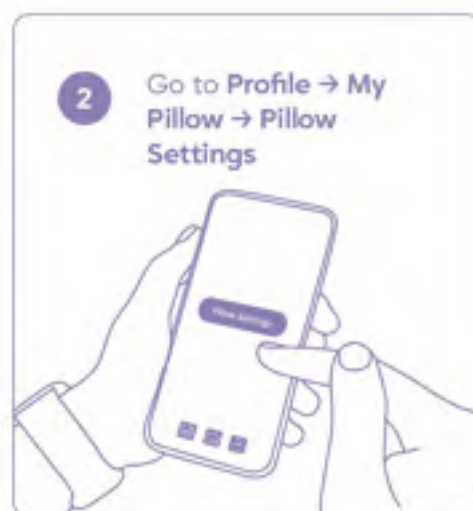

Max Inflation Height controls how high the pillow inflates during **Snore Intervention**; the AI may use less height when less response is needed.

Inflation Speed controls how fast the pillow inflates during **Snore Intervention**; the AI may adjust the speed and timing of each inflation based on snoring activity.

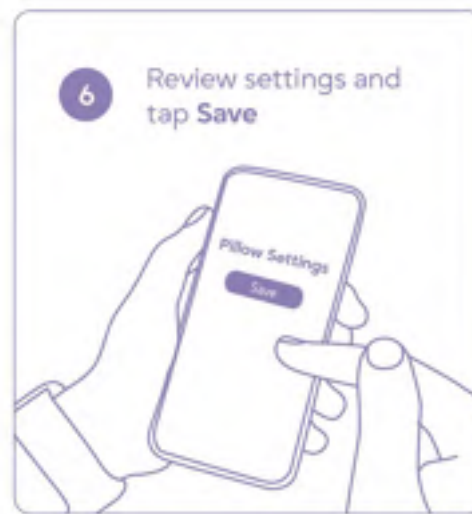

Profile Tip

Recommended settings are available only after your profile information has been completed. For best results, we recommend adding your profile information, such as **age, gender, height, and weight**. If you prefer not to provide this information, you can start with the default settings (**Max Inflation Height: Level 3, Inflation Speed: Level 2**). You can add or update your profile at any time in Profile.

Important

If **Snore Intervention feels too strong** or causes a significant sleep disturbance, reduce **Max Inflation Height** immediately.

Your Smart Anti-Snore Pillow will now begin active **Snore Intervention**.

Step 5: Your Next 2 Weeks – Adaptation & Fine-Tuning

Most users need about 2 weeks to adapt and experience the best snore-reduction results. During this time, your body adapts to gentle pillow movement while the AI learns your snoring patterns and personalizes Snore Intervention.

What Happens During This Time


Week 2: Getting Started

You may notice gentle pillow movement as the system begins responding to snoring. Keep your recommended or default settings unchanged during the first few nights.


Week 3: Adaptation & Optimization

Most users begin noticing the movement less and less each night as the AI continues personalizing the Snore Intervention. If needed, small setting adjustments can be made during this stage to improve comfort and performance. See **Understanding Your Snore Report (page 22)** for guidance on reviewing recommendations and adjusting settings.


Step 6: Nightly Routine – Pillow + App

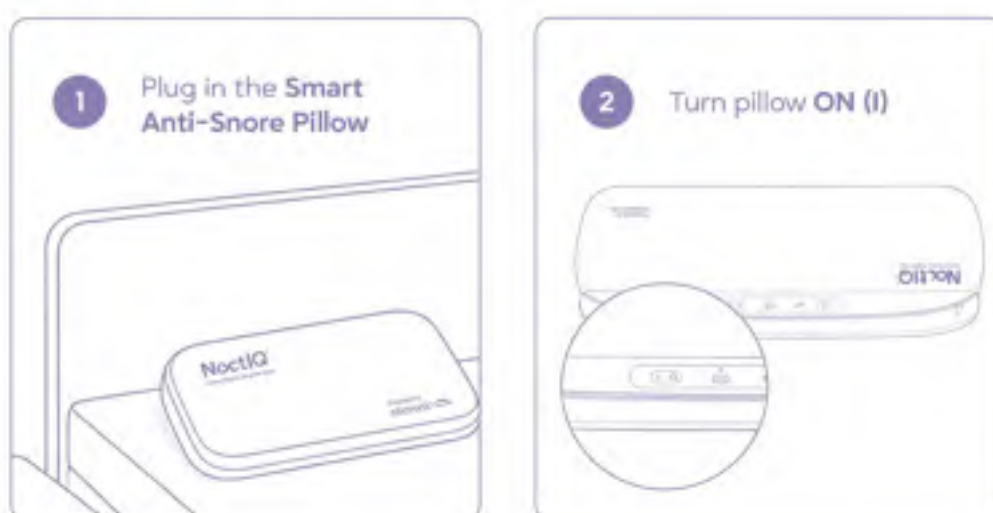

- 3 Confirm the **ZZ indicator** is solid blue.

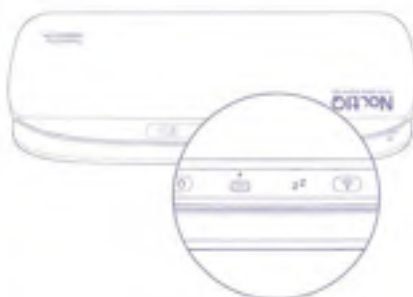

- 4 Keep your phone plugged in overnight.

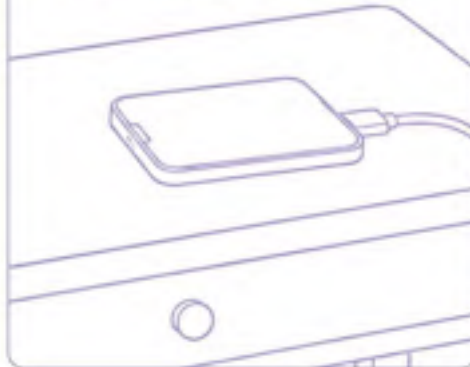


- 5 Place your phone on your **nightstand**, 1-1.5 ft away, **microphone facing pillow**.



- 6 Open the app → Tap **Start Tracking**.



- 7 Allow your **phone screen to turn off** automatically, or press the **power button to turn off** the screen / sleep normally.

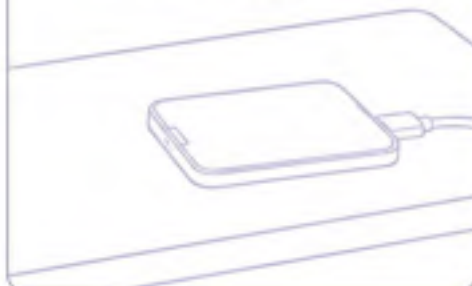

- 8 In the Morning, tap **Stop Tracking**.






Important

After tapping Start Tracking, do **NOT** use other audio-playing apps or sleep-detection software during the night, as this may affect snoring detection accuracy.

Understanding Your Snore Report

Your **Snore Report** helps you understand your snoring trends, compare **Baseline** nights with **Snore Intervention** nights, and track changes over time.

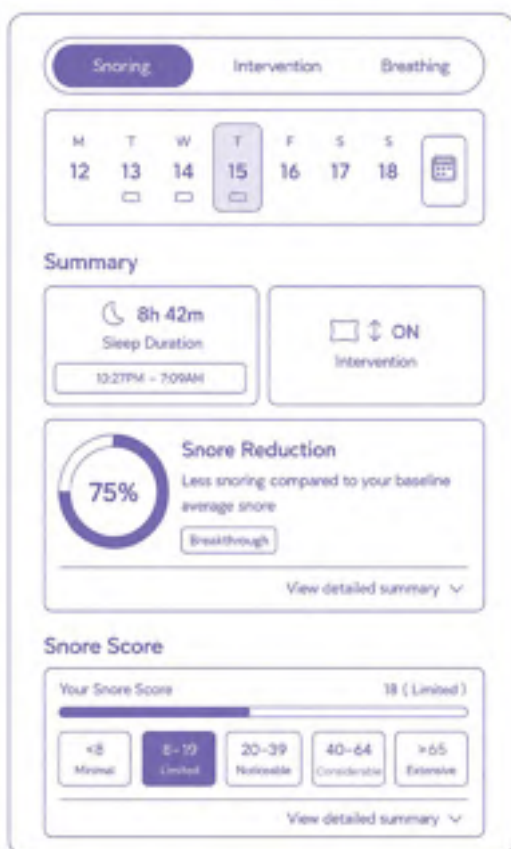

How to Use This Report

1. Review **Snore Report Metrics** to understand your snoring patterns and track changes over time.
2. Use the **Settings Recommendation Guide** to determine whether your current settings should be adjusted based on your results.

Snore Report Metrics

Your report compares two key periods:

- **Baseline – Snore Intervention OFF:** Your natural snoring patterns before active intervention.
- **Snore Intervention ON:** Results with active intervention.

This comparison helps show how the Smart Anti-Snore Pillow supports quieter, more restful sleep. Nights without detected snoring are excluded to help keep trends accurate.

The app highlights your **Snore Score** and **Snore Reduction**, with additional details available throughout the report. Tap a metric in the app to learn more.

- **Snore Score:** overall snoring level for the night.
- **Snore Reduction:** improvement compared to baseline nights.
- **Snore Duration:** total time spent snoring.
- **Snore Events:** number of individual snoring sounds detected.
- **Snore Intensity:** how often snoring was detected within each minute.
- **Snore Decibels:** average and range of snoring volume.
- **Snore Type Metrics:** additional insight into the types of snoring detected throughout the night.
 - **Soft Snoring:** lighter snoring sounds.
 - **Mid-Range Snoring:** moderate snoring sounds.
 - **Deep Snoring:** louder or more intense snoring sounds.
- **Snore Recordings:** optional playback stored only on your phone.

Note

Snore Reduction reporting begins after at least **5 snore nights with Snore Intervention OFF** and **3 snore nights with Snore Intervention ON**. Until then, use the other report metrics to understand your snoring patterns and progress.

Before You Adjust

Before making any adjustment, review your **Snore Report** and **Snore Reduction** result.

To change your settings: open the **NoctIQ** app and go to **Profile** → **My Pillow** → **Pillow Settings**.

For most setting changes, use the new setting for **3 nights** before making another adjustment. Do **NOT** adjust settings based on one unusual night affected by alcohol, congestion, illness, extreme fatigue, or an unusual sleep position.

Adjust immediately if:

- **Head feels pushed too high / neck pressure:** reduce **Max Inflation Height** by 1 level.
- **Pillow movement wakes you repeatedly / feels too sudden:** reduce **Inflation Speed** by 1 level.

Settings Recommendation Guide

Use this guide together with your **Snore Reduction** result.

If **Snore Reduction** is not yet available, continue collecting sleep data. The app may continue gathering information before providing a **Snore Reduction** result.

Some recommendations in this guide may vary based on body profile.

- **Lower BMI:** below 25
- **Higher BMI:** 25 or higher

If you are unsure, contact the **NoctIQ Care Team** for assistance.

Remember

Higher settings are not always better. The best setting is the one that reduces snoring while remaining comfortable.

Snore Reduction Result	Recommended Setting Adjustment
Green/ Light Green	Keep your current settings. If the pillow feels uncomfortable, use the Comfort Adjustment Guide below.
Yellow	Make a slight adjustment for 3 nights: • Higher BMI: Increase Max Inflation Height by 1 level. • Lower BMI: Reduce Max Inflation Height by 1 level. • If Max Inflation Height changes do not help after 2 levels, adjust Inflation Speed by 1 level.
Orange	Make a slight adjustment for 3 nights: • Higher BMI: Increase Max Inflation Height by 1 level. • Lower BMI: Reduce Max Inflation Height by 1 level. If results do not improve, return to your original setting, then increase Max Inflation Height by 1 level. Keep the setting that works better.
Red	Try a stronger adjustment for 3 nights: • Increase Max Inflation Height by 2 levels. OR • Increase Max Inflation Height by 1 level and Inflation Speed by 1 level. • If your result remains red, try the other option for 3 nights. If stronger settings do not improve results or feel too disruptive, return to your last comfortable setting and contact the NoctIQ Care Team. Some users may also choose to consult a qualified healthcare professional.

Comfort Adjustment Guide

If the pillow feels uncomfortable, adjust based on what you feel:

What You Feel	What to Adjust
Head pushed too high / neck pressure / neck strain	Reduce Max Inflation Height by 1 level.
Pillow movement wakes you / feels too sudden	Reduce Inflation Speed by 1 level.
General discomfort, not sure why	Reduce Max Inflation Height first. If discomfort continues, reduce Inflation Speed .

Pro Tip

Change only one setting at a time, tap **Save**, and allow **3 nights** before making another change, unless the app recommends a combined adjustment.

Advanced Sleep Wellness & Optional Membership Benefits

The **NoctIQ app** membership will provide a deeper view of your nighttime environment and long-term sleep wellness patterns.

Subscription features may include:

Respiratory Pattern Insights: Wellness-oriented analysis of breathing pattern consistency throughout the night.

Expanded Sleep History: Access to longer-term trends and past reports.

Enhanced Audio Review: Additional playback options for snoring and nighttime sound events.

Advanced Wellness Metrics: AI-supported insights to help you better understand your sleep environment.

Learn More or Manage Your Membership:

- Visit **noctiq.com** to learn more about optional membership benefits
- Open the app → **Profile** → **Manage Subscription**
- Contact the **NoctIQ Care Team**: Live Chat: **noctiq.com** |
Phone: **+1-888-598-7231** | Hours: **12PM-8PM EST** | Email: **care@noctiq.com**

Important Notice

Advanced sleep wellness features are for general wellness and informational purposes only. They are not medical features, are not intended for medical monitoring, and do not diagnose, treat, cure, mitigate, or prevent sleep apnea or any other medical condition. If you have concerns about your breathing, sleep quality, or symptoms, consult a qualified healthcare professional.

Quick Tips & FAQs

Pillow movement wakes you?

Lower **Inflation Speed** by 1 level.

Neck discomfort?

Lower **Max Inflation Height** by 1 level.

Prefer a higher or lower pillow?

Add or remove the **Height Adjustment Pad** to increase comfort.

Snoring not improving?

Adjust gradually. Change only **one setting at a time**, tap **Save**, and allow **3 nights** for your body to adapt before making another change. Follow the recommended Settings Recommendation Guide.

Why is there a 2-week adaptation period?

Your body adapts to the pillow's support structure and gentle movement while the AI learns your snoring patterns and personalizes **Snore Intervention**.

Why don't I have a Snore Reduction result yet?

Snore Reduction reporting becomes available after at least 5 snore nights with Snore Intervention OFF and 3 snore nights with Snore Intervention ON. Continue tracking your sleep and reviewing your Snore Report while sufficient data is collected.

Will the Smart Anti-Snore Pillow wake me?

Pillow movement is designed to be gentle. Most users adapt over time.

Should I start at the highest Max Inflation Height and Inflation Speed?

No. Start with the app's recommended settings. If no profile information was provided, use the default settings: **Max Inflation Height Level 3** and **Inflation Speed Level 2**.

Does the pillow need to stay plugged in?

Yes. The Smart Anti-Snore Pillow does not contain a battery and must remain connected to power during overnight use.

Can I use the pillow while traveling?

Yes. If 2.4GHz Wi-Fi is available.

Still need help?

Contact the **NoctIQ Care Team** for personalized assistance.

Returns, Warranty & Regulatory Information

Return Policy

30-Night Sleep Trial

We want you to love your sleep-and the peace it brings to your bedroom. If the product isn't right for you, you can return it within **30 nights of delivery**. You can feel confident that we will honor your return.

Contact Us for Help

You made an investment in better sleep and a quieter, more connected night with your partner.

Before deciding to return your product, we strongly encourage you to reach out to the **NoctIQ Care Team**. Most customers see meaningful improvement - and real snore-reduction results - after we help fine-tune their settings.

We're here to help you get the benefits - and the peaceful nights - you deserve.

Live Chat: noctiq.com

Phone: +1-888-598-7231 | **Hours:** 12PM-8PM EST

Email: care@noctiq.com

How to Start a Return

If you ultimately decide to return the product, contact the NoctIQ Care Team. We'll guide you through the simple steps.

Return Shipping

For non-defective returns, a return shipping fee will be deducted from your refund. Defective items are returned at no cost.

Refund Timing

Refunds are issued within **7 business days** after the product is received.

Condition of Returned Items

Returned products must include all components and be in good condition.

Limited Warranty (2 Years)

This product is covered by a two-year limited warranty from the date of original purchase. During the warranty period, NoctIQ Inc. will repair or replace, at its discretion, any defects in materials or workmanship related to the internal electronic components, sensors, air system, and control unit.

What Is Covered

- Internal electronics
- Sensors
- Air pump and air chambers (manufacturing defects only)
- Control unit and internal wiring
- Functional failures not caused by misuse or improper care

What Is Not Covered

This warranty does not cover:

- Stains, discoloration, or odor caused by normal use
- Fabric wear, pilling, or cosmetic changes over time
- Foam compression or changes in firmness
- Damage caused by liquids, improper cleaning, or misuse
- Normal aging of materials
- Unauthorized repairs or modifications

How to Obtain Service

Proof of purchase is required. Please contact **NoctIQ Care Team** for troubleshooting and warranty service instructions.

Limitations

This limited warranty is non-transferable and applies only to products purchased from authorized retailers. NoctIQ is not responsible for incidental or consequential damages. Some jurisdictions do not allow limitations on implied warranties or exclusion of incidental damages; these limitations may not apply to you.

Optional Extended Warranty (+1 Year)

For added peace of mind, **you can extend your coverage to a full 3 years with our optional +1-year extended warranty.** This upgrade protects your investment longer and ensures uninterrupted support and coverage.

Extend Your Coverage

Visit noctiq.com to learn more or add the optional **+1-year extended warranty.**

Regulatory Information

ETL Listing (U.S. & Canada)

This product is ETL Listed for the United States and Canada under the following standards:

- **UL 60335-1**
- **CSA C22.2 No. 60335-1**

The ETL mark indicates that the product has been evaluated to applicable safety standards for household electrical appliances.

FCC Compliance (United States)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the manufacturer could void the user's authority to operate this equipment.

Note

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Compliance (Canada)

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

CAN ICES-003(B) / NMB-003(B)

This Class B digital apparatus complies with Canadian ICES-003.

Textile & Flammability Compliance

This product meets the following textile and flammability requirements:

- **16 CFR 1610** - Flammability of Clothing Textiles
- **TB 117-2013** - Smoldering Resistance of Upholstered Furniture Components
- **Flame Retardant Claim Compliance Verification**

Chemical Safety & California Proposition 65

This product has been tested and meets applicable chemical safety requirements, including those established under California Proposition 65.

WEEE / Recycling Information

Do not dispose of this product in household trash.

This product contains electronic components and must be recycled or disposed of according to local regulations for waste electrical and electronic equipment (WEEE).

- Follow local guidelines for electronic waste recycling
- Contact your municipality or recycling center for proper disposal
- Do not incinerate or expose electronic components to extreme heat

Wishing you **calmer nights**
and brighter days.

Thank you for choosing **NoctIQ.**





Customer Support

Live Chat: noctiq.com

Phone: +1-888-598-7231

Email: care@noctiq.com

Distributed by NoctIQ Inc.,
Glen Allen, VA 23059 USA

©2026 NoctIQ Inc. All rights reserved

