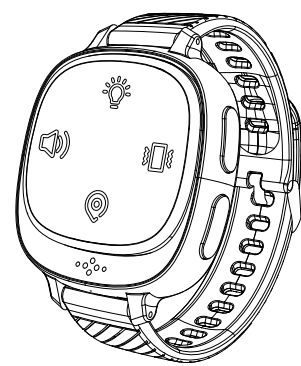
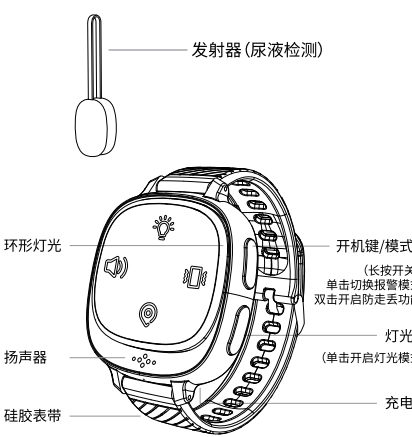
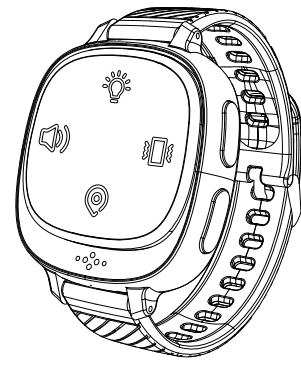
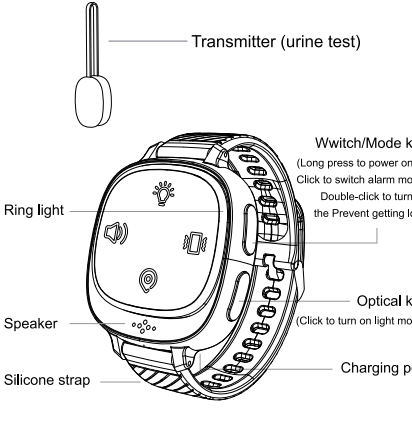
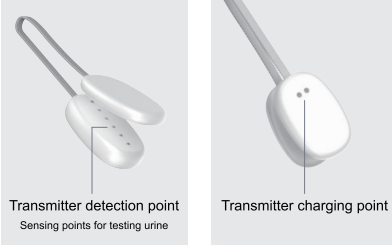
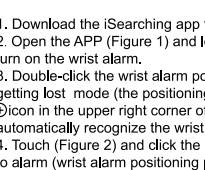
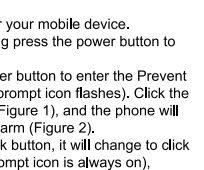
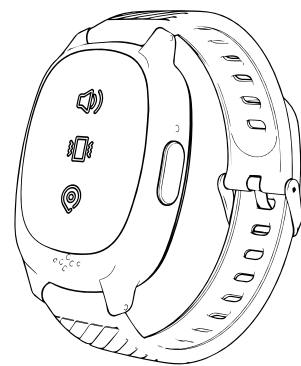
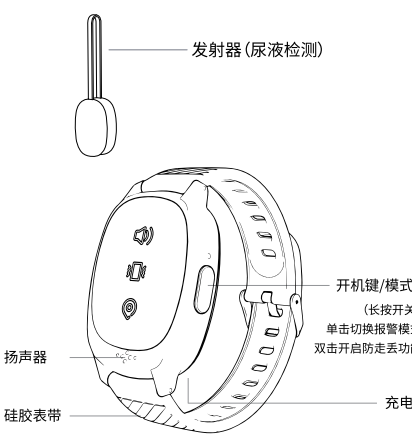
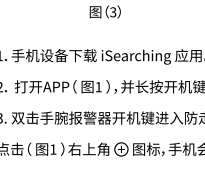
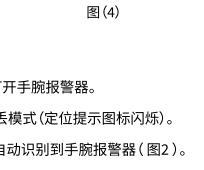
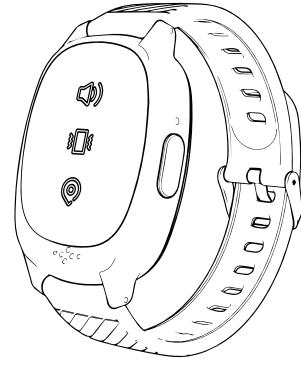
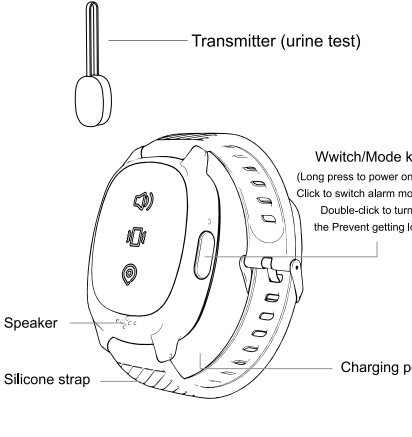
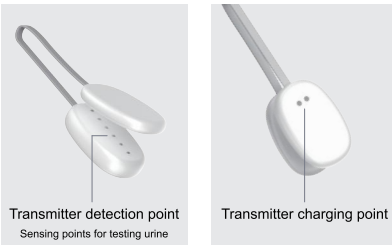
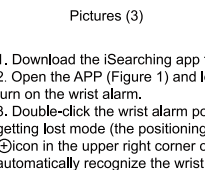
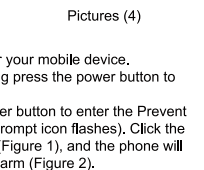


说明书 (折页)
尺寸:63.3*90mm
材质:100克铜版纸
工艺:双面4色印刷

尿湿报警器 使用说明  ● 感谢您购买我们的产品 使用前请务必仔细阅读此说明书并妥善保管,以备将来使用参考。	产品组成 尿液报警器由两个部分构成: 发射器(尿液检测):检测尿液通过无线方式发送报警信息给接收器。 接收器(手腕报警器):接收器用于接收发射器报警信息,通过振动或者声音方式通知用户。  发射器(尿液检测) 发射器检测点 检测尿液的位置点 发射器充电点 开机键/模式键 (长按开机关机 单击切换报警模式, 双击开启防走失功能) 灯光键 (单击开启灯光模式) 扬声器 硅胶表带 充电口	发射器介绍  发射器检测点 检测尿液的位置点 发射器充电点 特别提醒:将发射器扣在能检测尿液的地方,发射器检测点的金属片要接触到尿液才能报警。 显示屏图标功能介绍 <table><tr><th>图标</th><th></th><th></th><th></th><th></th></tr><tr><th>功能</th><td>报警声音提示</td><td>振动提示</td><td>灯光提示</td><td>防走失提示</td></tr></table>	图标					功能	报警声音提示	振动提示	灯光提示	防走失提示	报警器正确使用说明 主功能: 声音提醒+振动提醒+灯光模式+防走失模式。 功能使用说明: - 长按右侧上键(开机键)1.5 秒开机。 - 开机后声音和灯光图标同时亮起,幻彩灯光环绕一周之后和灯光图标同时熄灭,声音图标常亮。 - 单击右侧上键(开机键)切换振动模式,对应振动灯光图标常亮。 - 双击右侧上键(开机键)切换到防走失模式,防走失模式灯光图标常亮,进入配对状态,打开手机浏览器下载ISEARCHING防走失器进行配对,配对名为ITAG,连接上后防走失灯光图标常亮。 - 双击右侧上键(开机键),可退出防走失模式。 - 单击右侧下键(灯光键),开启幻彩灯光模式,共5种灯光模式。单击灯光键切换灯光模式,每组灯光模式时间20秒自动关闭,中途可长按右侧下键(灯光键)1.5 秒关闭灯光。 - 低电时自动关机,充电时声音图标闪烁。 - 额定输入:5V/2A。 发射器: - 报警器开机后,发射器自动配对。 - 发射器遇水报警(蓝色灯光闪烁),报警器声音/振动 报警。 - 充电时红灯微闪烁,充满电后灯光熄灭。	手机APP操作说明 尿湿报警器APP图标  ● 苹果设备请扫描二维码下载APP,或在APP Store中搜索ISEARCHING关键词进行下载。 ● 安卓设备请扫描二维码下载APP,或在设备浏览器中搜索ISEARCHING关键词进行下载。 操作步骤:  图 (1)  图 (2)	手机APP操作说明  图 (3)  图 (4) - 手机设备下载 ISeaching 应用。 - 打开APP (图1),并长按开机键打开手腕报警器。 - 双击手腕报警器开机键进入防走失模式(定位提示图标闪烁)。点击(图1)右上角图标,手机会自动识别到手腕报警器(图2)。 - 触碰(图2)点击链接按钮,会变为点击报警(手腕报警器定位提示图标常亮),表示手机与设备已经链接成功(图3)。 - 点击设置图标(图3),打开蓝牙断开报警按钮(图4),手腕报警器距离手机20米外,手机自动响铃报警。 - 手腕报警器关机后设备与手机断开连接,下次连接需要再次进入防走失模式即可从新链接。	产品介绍 重要作用 尿湿报警器是一种用于监测尿液泄漏的电子设备,主要服务于婴幼儿、老年人或失禁患者: 1.及时提醒护理人员 当使用者(如婴儿、卧床患者)排尿时,报警器通过温度传感器检测到尿液,立即触发声音/振动提醒看护人员,及时更换尿布或床单,避免长时间潮湿刺激皮肤。 2.预防尿布疹和感染 减少尿液或粪便与皮肤长时间接触,降低红疹、湿疹、尿路感染等风险,尤其对婴幼儿和行动不便的老人尤为重要。 3.改善睡眠质量 对于夜间护理,报警器可避免频繁手动检查尿布,让看护者和使用者(如父母或老人)获得更好的睡眠。 4.辅助如厕训练 在训练婴幼儿自主如厕时,报警器可帮助孩子建立排尿与如厕的关联,通过提醒逐步养成习惯。 5.减轻护理负担 对长期照护失禁患者的家庭或机构,报警器能减少人工监控压力,提升护理效率。	产品信息 型号:bj03 电池容量:报警器电池容量300mAh,发射器电池容量50mAh。 电源要求:产品充电采用标准5V/2A供电。 ■ 需定期清洁传感器,避免误报。 保养/清洁 在每次使用本产品后,记住要清理发射器。如有需要的话可和报警器一并清洗干净,不要把报警器本体浸在水里,容易损坏内部零件。在清理发射器时,请先关机。发射器可在水里冲洗,清洗干净擦干即可。本公司报警器是为单个用户设计的,如果你想另一个孩子使用报警器,请联系购买另外一个报警器。  警告:在使用中途发生报警时,请将发射器磁吸扣打开并用干质的纸巾擦拭干净后,方可使用。
图标																	
功能	报警声音提示	振动提示	灯光提示	防走失提示													
URINE ALARM  ● Thank you for purchasing our product Please read this manual carefully before use and keep it properly for future reference.	Product composition The urine alarm consists of two parts: Transmitter (urine test): Detect urine and send alarm information to the receiver wirelessly. Receiver (wrist alarm): The receiver is used to receive alarm information from the transmitter and notify the user through vibration or sound.  Transmitter (urine test) Transmitter detection point Sensing points for testing urine Transmitter charging point Switch/Mode key (Long press to power off Click to switch alarm mode, Double-click to turn on the Prevent getting lost) Ring light Speaker Silicone strap Optical key (click to turn on light mode) Charging port	Transmitter introduction  Transmitter detection point Sensing points for testing urine Transmitter charging point Special reminder: Place the transmitter somewhere that can detect urine. The metal piece at the transmitter detection point Contact with urine is required to call the police.d. Icon introduction <table><tr><th>Icon</th><th></th><th></th><th></th><th></th></tr><tr><th>Function</th><td>Alarm sound prompt</td><td>Vibration reminder</td><td>Light indication</td><td>Prevent getting lost</td></tr></table>	Icon					Function	Alarm sound prompt	Vibration reminder	Light indication	Prevent getting lost	Proper instructions for using the alarm Main functions: Sound reminder + vibration reminder + light mode + Prevent getting lost mode. Function instructions: - Press and hold the upper button on the right side (power on button) for 1.5 seconds to turn on the phone. - After long pressing to power on, the sound and light icons will light up at the same time. The magic light circles around in a circle and turns off at the same time as the light icon. The sound icon will always be on. - Click the upper button on the right (power on button) to switch to vibration mode, and the corresponding vibration light icon will be on. - Double-click the upper button on the right side (power on button) to switch to Prevent getting lost mode. The light icon of Prevent getting lost mode flashes slowly. Enter the pairing state, open the mobile browser to download ISEARCHING Prevent getting lost device for pairing. The pairing name is ITAG. After connecting, the Prevent getting lost light icon is always on. Double-click the upper button on the right (power button) to exit the Prevent getting lost mode. - Click the lower button on the right (light button) to turn on the magic lighting mode. There are 5 lighting modes in total. Click the light key switches the light mode. Each group of light mode will automatically turn off after 20 seconds. You can turn off the light by pressing and holding the right down button (light button) for 1.5 seconds. - It automatically shuts down when the battery is low, and the sound icon flashes when charging. - Rated input: 5V/2A. Transmitter: - After the alarm is turned on, the transmitter pairs automatically. - When the transmitter encounters water, it will alarm (blue light flashes), and the alarm will sound/vibrate. - The red light flashes slowly when charging, and turns off when fully charged.	Mobile APP operation instructions Humidity alarm APP icon  ● For Apple devices, please scan the QR code to download the APP, or search in the APP Store Searching keywords to download. ● For Android devices, please scan the QR code to download the APP, or search in the device browser Searching keywords to download. Operating Instructions:  Step ①  Step ② Pictures (1) Pictures (2)	Mobile APP operation instructions  Pictures (3)  Pictures (4) - Download the ISeaching app for your mobile device. - Open the APP (Figure 1) and long press the power button to turn on the wrist alarm. - Double-click the wrist alarm power button to enter the Prevent getting lost mode (the positioning prompt icon flashes). Click the in the upper right corner of (Figure 1), and the phone will automatically recognize the wrist alarm (Figure 2). - Touch (Figure 2) and click the link button. It will change to click to alarm (wrist alarm positioning prompt icon is always on), indicating that the mobile phone and the device have been successfully linked (Figure 3). - Click the settings icon (Picture 3), turn on the Bluetooth disconnection alarm button (Picture 4), and the wrist alarm When the phone is 20 meters away, the phone will automatically ring the alarm. - After the wrist alarm is turned off, the device is disconnected from the mobile phone, and you need to re-enter the connection next time. Prevent getting lost mode can be re-linked.	Product introduction Important Role Wetness alarm is an electronic device used to monitor urine leakage. It mainly serves infants, young children, Elderly or incontinent patients. 1. Remind the nursing staff promptly When the user (such as an infant or a bedridden patient) urinates, the alarm detects the urine through the humidity sensor. Liquid, immediately triggers sound/vibration to remind caregivers to change diapers or sheets in time to avoid long Moisture in the room irritates the skin. 2. Prevent diaper rash and infections Reduce the prolonged contact of urine or feces with the skin and reduce the risk of rash, eczema, urinary tract infection, etc. It is especially important for infants and young children and the elderly with limited mobility. 3. Improve sleep quality For nighttime care, the alarm eliminates the need for frequent manual diaper checks, allowing caregivers and users (like parents or the elderly) get better sleep. 4. Assist toilet training When training infants and young children to go to the toilet independently, the alarm can help the child establish the relationship between urination and going to the toilet. Build habits gradually through reminders. 5. Reduce the burden of care For families or institutions that care for incontinence patients for a long time, alarms can reduce the pressure of manual monitoring and improve Nursing efficiency.	Product information Model:bj03 Battery capacity: Alarm battery capacity is 300mAh, transmitter battery capacity is 50mAh. Power requirements: Product charging uses standard 5V/2A power supply. ■ Sensors need to be cleaned regularly to avoid false alarms. Cleaning and maintenance Please remember to clean the transmitter after each use of this product. Call the police if necessary Wipe it off together. Do not immerse the alarm body in water, otherwise the internal parts may be damaged, clean up When using the transmitter, please turn off the power first. The transmitter can be rinsed, cleaned and dried with water. Our alarms are designed for a single user, if you would like another child to use the alarm, Please contact us to purchase another alarm.  Warning: When an alarm occurs during use, please open the magnetic buckle of the transmitter.Wipe clean with dry paper towel before use.
Icon																	
Function	Alarm sound prompt	Vibration reminder	Light indication	Prevent getting lost													

尿湿报警器 使用说明  ● 感谢您购买我们的产品 使用前请务必仔细阅读此说明书并妥善保管,以备将来使用参考。	产品组成 尿液报警器由两个部分构成: 发射器(尿液检测):检测尿液通过无线方式发送报警信息给接收器。 接收器(手腕报警器):接收器用于接收发射器报警信息,通过振动或者声音方式通知用户。  发射器(尿液检测) 发射器检测点 检测尿液的位置点 发射器充电点 开机键/模式键 (长按开机关机 单击切换报警模式, 双击开启防走失功能) 扬声器 硅胶表带 充电口	发射器介绍  发射器检测点 检测尿液的位置点 发射器充电点 特别提醒:将发射器扣在能检测尿液的地方,发射器检测点的金属片要接触到尿液才能报警。 显示屏图标功能介绍 <table><tr><th>图标</th><th></th><th></th><th></th><th></th></tr><tr><th>功能</th><td>报警声音提示</td><td>振动提示</td><td>灯光提示</td><td>防走失提示</td></tr></table>	图标					功能	报警声音提示	振动提示	灯光提示	防走失提示	报警器正确使用说明 主功能: 声音提醒+振动提醒+防走失模式。 功能使用说明: - 长按右侧开机键1.5 秒开机。 - 开机后默认声音报警功能开启,声音灯光图标常亮。 - 单击右侧开机键切换振动模式,对应振动灯光图标常亮。 - 双击右侧开机键切换到防走失模式,防走失模式灯光图标常亮,进入配对状态,打开手机浏览器下载ISEARCHING防走失器进行配对,配对名为ITAG,连接上后防走失灯光图标常亮。 - 双击右侧开机键,可退出防走失模式。 - 连接上后防走失灯光图标常亮,双击右侧开机键可退出防走失模式。 - 低电时自动关机,充电时声音/振动/防走失图标依次闪烁。 - 额定输入:5V/2A。 发射器: - 报警器开机后,发射器自动配对。 - 发射器遇水报警(蓝色灯光闪烁),报警器声音/振动报警。 - 充电时红灯微闪烁,充满电后灯光熄灭。	手机APP操作说明 尿湿报警器APP图标  ● 苹果设备请扫描二维码下载APP,或在APP Store中搜索ISEARCHING关键词进行下载。 ● 安卓设备请扫描二维码下载APP,或在设备浏览器中搜索ISEARCHING关键词进行下载。 操作步骤:  图 (1)  图 (2)	手机APP操作说明  图 (3)  图 (4) - 手机设备下载 ISeaching 应用。 - 打开APP (图1),并长按开机键打开手腕报警器。 - 双击手腕报警器开机键进入防走失模式(定位提示图标闪烁)。点击(图1)右上角图标,手机会自动识别到手腕报警器(图2)。 - 触碰(图2)点击链接按钮,会变为点击报警(手腕报警器定位提示图标常亮),表示手机与设备已经链接成功(图3)。 - 点击设置图标(图3),打开蓝牙断开报警按钮(图4),手腕报警器距离手机20米外,手机自动响铃报警。 - 手腕报警器关机后设备与手机断开连接,下次连接需要再次进入防走失模式即可从新链接。	产品介绍 重要作用 尿湿报警器是一种用于监测尿液泄漏的电子设备,主要服务于婴幼儿、老年人或失禁患者: 1.及时提醒护理人员 当使用者(如婴儿、卧床患者)排尿时,报警器通过温度传感器检测到尿液,立即触发声音/振动提醒看护人员,及时更换尿布或床单,避免长时间潮湿刺激皮肤。 2.预防尿布疹和感染 减少尿液或粪便与皮肤长时间接触,降低红疹、湿疹、尿路感染等风险,尤其对婴幼儿和行动不便的老人尤为重要。 3.改善睡眠质量 对于夜间护理,报警器可避免频繁手动检查尿布,让看护者和使用者(如父母或老人)获得更好的睡眠。 4.辅助如厕训练 在训练婴幼儿自主如厕时,报警器可帮助孩子建立排尿与如厕的关联,通过提醒逐步养成习惯。 5.减轻护理负担 对长期照护失禁患者的家庭或机构,报警器能减少人工监控压力,提升护理效率。	产品信息 型号:bj02 电池容量:报警器电池容量300mAh,发射器电池容量50mAh。 电源要求:产品充电采用标准5V/2A供电。 ■ 需定期清洁传感器,避免误报。 保养/清洁 在每次使用本产品后,记住要清理发射器。如有需要的话可和报警器一并清洗干净,不要把报警器本体浸在水里,容易损坏内部零件。在清理发射器时,请先关机。发射器可在水里冲洗,清洗干净擦干即可。本公司报警器是为单个用户设计的,如果你想另一个孩子使用报警器,请联系购买另外一个报警器。  警告:在使用中途发生报警时,请将发射器磁吸扣打开并用干质的纸巾擦拭干净后,方可使用。
图标																	
功能	报警声音提示	振动提示	灯光提示	防走失提示													
URINE ALARM  ● Thank you for purchasing our product Please read this manual carefully before use and keep it properly for future reference.	Product composition The urine alarm consists of two parts: Transmitter (urine test): Detect urine and send alarm information to the receiver wirelessly. Receiver (wrist alarm): The receiver is used to receive alarm information from the transmitter and notify the user through vibration or sound.  Transmitter (urine test) Transmitter detection point Sensing points for testing urine Transmitter charging point Switch/Mode key (Long press to power off Click to switch alarm mode, Double-click to turn on the Prevent getting lost) Ring light Speaker Silicone strap Optical key (click to turn on light mode) Charging port	Transmitter introduction  Transmitter detection point Sensing points for testing urine Transmitter charging point Special reminder: Place the transmitter somewhere that can detect urine. The metal piece at the transmitter detection point Contact with urine is required to call the police.d. Icon introduction <table><tr><th>Icon</th><th></th><th></th><th></th><th></th></tr><tr><th>Function</th><td>Alarm sound prompt</td><td>Vibration reminder</td><td>Light indication</td><td>Prevent getting lost</td></tr></table>	Icon					Function	Alarm sound prompt	Vibration reminder	Light indication	Prevent getting lost	Proper instructions for using the alarm Main functions: Sound reminder + vibration reminder + Prevent getting lost mode. Function instructions: - Press and hold the power button on the right side for 1.5 seconds to turn on the phone. - After power on, the sound alarm function is turned on By default, the sound icon is always on. - Click the power button on the right to switch to vibration mode, and the corresponding vibration light icon will stay on. - Double-click the power button on the right side to switch to Prevent getting lost mode. The light icon of Prevent getting lost mode flashes slowly. Enter the pairing state, open the mobile browser to download ISEARCHING Prevent getting lost device for pairing. The pairing name is ITAG. After connecting, the Prevent getting lost light icon is always on. Double-click the power button on the right to exit the Prevent getting lost mode. - It automatically shuts down when the battery is low, and the sound/vibration/Prevent getting lost icons flash in sequence during charging. - Rated input: 5V/2A. Transmitter: - After the alarm is turned on, the transmitter pairs automatically. - When the transmitter encounters water, it will alarm (blue light flashes), and the alarm will sound/vibrate. - The red light flashes slowly when charging, and turns off when fully charged.	Mobile APP operation instructions Humidity alarm APP icon  ● For Apple devices, please scan the QR code to download the APP, or search in the APP Store Searching keywords to download. ● For Android devices, please scan the QR code to download the APP, or search in the device browser Searching keywords to download. Operating Instructions:  Step ①  Step ② Pictures (1) Pictures (2)	Mobile APP operation instructions  Pictures (3)  Pictures (4) - Download the ISeaching app for your mobile device. - Open the APP (Figure 1) and long press the power button to turn on the wrist alarm. - Double-click the wrist alarm power button to enter the Prevent getting lost mode (the positioning prompt icon flashes). Click the in the upper right corner of (Figure 1), and the phone will automatically recognize the wrist alarm (Figure 2). - Touch (Figure 2) and click the link button. It will change to click to alarm (wrist alarm positioning prompt icon is always on), indicating that the mobile phone and the device have been successfully linked (Figure 3). - Click the settings icon (Picture 3), turn on the Bluetooth disconnection alarm button (Picture 4), and the wrist alarm When the phone is 20 meters away, the phone will automatically ring the alarm. - After the wrist alarm is turned off, the device is disconnected from the mobile phone, and you need to re-enter the connection next time. Prevent getting lost mode can be re-linked.	Product introduction Important Role Wetness alarm is an electronic device used to monitor urine leakage. It mainly serves infants, young children, Elderly or incontinent patients. 1. Remind the nursing staff promptly When the user (such as an infant or a bedridden patient) urinates, the alarm detects the urine through the humidity sensor. liquid, immediately triggers sound/vibration to remind caregivers to change diapers or sheets in time to avoid long Moisture in the room irritates the skin. 2. Prevent diaper rash and infections Reduce the prolonged contact of urine or feces with the skin and reduce the risk of rash, eczema, urinary tract infection, etc. It is especially important for infants and young children and the elderly with limited mobility. 3. Improve sleep quality For nighttime care, the alarm eliminates the need for frequent manual diaper checks, allowing caregivers and users (like parents or the elderly) get better sleep. 4. Assist toilet training When training infants and young children to go to the toilet independently, the alarm can help the child establish the relationship between urination and going to the toilet. Build habits gradually through reminders. 5. Reduce the burden of care For families or institutions that care for incontinence patients for a long time, alarms can reduce the pressure of manual monitoring and improve Nursing efficiency.	Product information Model:bj02 Battery capacity: Alarm battery capacity is 300mAh, transmitter battery capacity is 50mAh. Power requirements: Product charging uses standard 5V/2A power supply. ■ Sensors need to be cleaned regularly to avoid false alarms. Cleaning and maintenance Please remember to clean the transmitter after each use of this product. Call the police if necessary Wipe it off together. Do not immerse the alarm body in water, otherwise the internal parts may be damaged, clean up When using the transmitter, please turn off the power first. The transmitter can be rinsed, cleaned and dried with water. Our alarms are designed for a single user, if you would like another child to use the alarm, Please contact us to purchase another alarm.  Warning: When an alarm occurs during use, please open the magnetic buckle of the transmitter.Wipe clean with dry paper towel before use.
Icon																	
Function	Alarm sound prompt	Vibration reminder	Light indication	Prevent getting lost													

FCC Warning:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 0cm between the radiator and your body.