




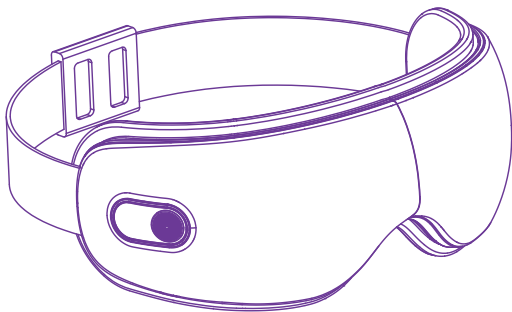

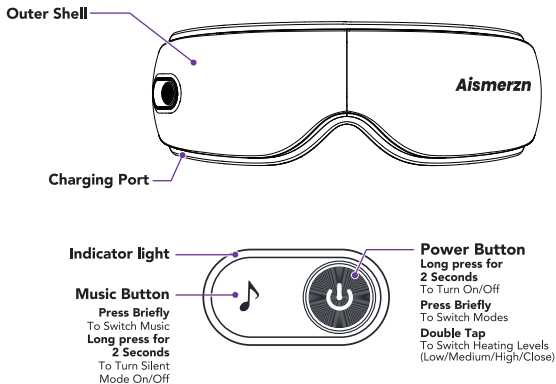

Internal Components
(Air Compression + Hot Compress)

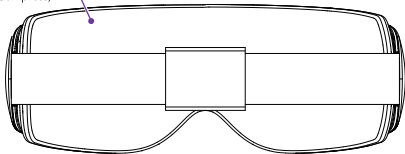

* Images and diagrams referenced in this manual are for illustrative purposes only,
The actual appearance may differ.

Basic Product Parameters

Product Name	Eye Massager	Model	ASME02
Rated Input Power	5W	Rated Input Voltage	5V $\overline{\text{---}}$
Battery Capacity	1500mAh/5.55Wh		

*Before using this product, please first adjust the strap to an appropriate length, then wear it on the eye area, and adjust the position and tightness as needed to experience the massage effect.

Operating Steps:

- 1) Long press the power button for two seconds to turn the machine on/off. Press it briefly to switch among three working modes: Vitality Mode, Soothing Mode and Sleep Mode. The default mode is Vitality Mode. Double tap the power button to switch hot compress levels(Low/Medium/High/Close). Long press the music button to turn off music, and long press it again to turn it on.
- 2) Manual pairing and connection is required for the first time with Bluetooth devices that have never been connected to the product. The product's Bluetooth name is ASME02 . Devices that have previously been connected before will automatically reconnect to the Bluetooth when the product is turned on. A voice announcement "Bluetooth Connected" will be made after a successful connection.
- 3) When the product's battery is low, a voice prompt will indicate, "Low Battery, Please Charge," and the indicator will glow red. The product will shut down automatically if it is used further.
- 4) The product's charging icon will remain red while charging, and it will turn green when fully charged.
- 5) The product will shut down automatically while charging if it is on and cannot be turned on while charging.

Working Mode

Vitality Mode: Air Compression 100% + Heat + Music

Soothing Mode: Air Compression 80%+Heat + Music

Sleep Mode: Heat+Music

- When the surface of the massager is dirty, wipe it with a clean, soft cloth dampened with a little water. Do not wipe with volatile liquids like alcohol. The product cannot be disassembled or washed in water.
- Please cover the product when not in use for an extended period of time to prevent dust accumulation.
- Prevent water or other liquids from entering the machine during cleaning to avoid short circuits during use.

Troubleshooting Instructions

If the product does not work properly, please check if the following situations have occurred. If our suggestions do not solve the problem, please return the product to the dealer for repair.

Condition	Possible Cause	Solution
Unable to Power On	Internal Battery Low	Charge the Product
	Product is Charging	Use After Fully Charged
Unable to Charge /No Charging Indicator	Charging Cable Not Properly Connected to Main Unit and Power Adapter	Reconnect Charging After Unplugging
There is a mechanical or air sound coming from the device.	Normal sound of air inflation or air compression.	No action required.
The device turns off automatically. / Buttons are not responsive. / Air compression is weak.	Low battery.	Charge the device.
	Built-in 15mins auto-off timer is activated.	Turn on the device again.

Please support the recycling of rechargeable batteries (lithium-ion batteries)

 Danger	This product's battery is a dedicated rechargeable battery that should not be used in other products. Do not charge the battery once it has been removed. • Do not throw into the fire or expose to excessive heat. • Do not puncture with nails, subject to heavy impact, or disassemble • Do not directly connect the battery's positive and negative terminals with metal • Do not store it with necklaces, hairpins, etc. • Do not charge, use, or store in high-temperature locations such as near a stove or under direct sunlight • All of the abovementioned scenarios can result in overheating, fire, or bursting.
 Warning	Place the removed rechargeable battery out of the reach of infants to prevent accidental ingestion • If swallowed, seek medical attention immediately. • If battery leakage occurs, do not touch the liquid with bare hands. • Follow the safety instructions below: • There is the risk of blindness if battery fluid contacts the eyes. Do not rub the eye with your hands. Rinse with clean water immediately and seek medical attention. If battery liquid gets on the body or clothing, it may cause skin inflammation or injury. Rinse with clean water, and if discomfort persists, seek medical attention mediatly. • Disassembly is only permitted when the product is being discarded. Failure to comply may result in fire, electric shock, injury, or other accidents. Please contact the retailer where you purchased the product for repairs (except battery replacement).

The lithium-ion battery used in this product is a valuable recyclable resource. Please remove the lithium-ion battery prior to disposal and hand over it to a professional organization capable of recycling and processing rechargeable batteries.

- The battery must be removed from the device before discarding the appliance.
- The device must be turned off before removing the battery.
- The battery should be disposed of safely.

Removal Steps:

1. Pry open the main unit's inner and outer shells using professional tools.
2. Unscrew the fixed screws of the internal bracket using a screwdriver.
3. Unscrew the screws at the circuit board.
4. Cut the battery leads using scissors.
5. Remove the battery.



Caution

- Read and follow the operation guidelines in the User Manual carefully before using this product. The company is not responsible for any consequences resulting from incorrect use or failure to comply with the User Manual.
- At the first use, charge the battery in accordance with the User Manual to ensure its normal operation, as the battery may have run out after not being used for an extended period.
- As we are constantly improving our products, technical specifications are subject to change without prior notice.
- Please use this product with caution after applying facial makeup, as it may ruin your makeup. Cosmetics can also contaminate the product's contact surface. If such cases occur, wipe the contaminated surface with clean water.
- Do not use the eye massager while driving.
- Do not use the eye massager in excessively high temperatures or excessively humid conditions.
- Do not use the product if you have a history of stroke, hypertension, or hypotension.
- Do not use the product if you have eye conditions such as glaucoma, cataracts, retinal detachment, or pinkeye.
- Avoid using the eye massager after consuming alcohol.
- The product is prohibited for individuals who have undergone eye surgery or have heart conditions.
- Remove contact lenses, glasses, false eyelashes, hairpins, and accessories before use.
- This product is a personal care massage product and is not intended to replace pharmaceuticals and medical devices.
- Children, pregnant women, and women during menstruation should use the eye massager with caution.
- The product contains a lithium battery. Do not place it in high temperatures.
- For sensitive skin, facial allergic dermatitis, and those with plastic surgery materials on the face, use under the guidance of a doctor.
- In case of any discomfort during use, please stop using the product immediately.
- Do not massage anywhere on your body other than the eyes.
- Do not disassemble the product yourself.
- The product is not waterproof. Do not wash it or use it in high humidity environments.
- If the lining is stained, gently wipe it with a damp cloth.
- Please use this product in a normal operating environment
Temperature: 0°C–40°C;
Relative humidity: <80%RH, non-condensing
Atmospheric pressure: 50Kpa–106Kpa
- Please transport and store the product under normal conditions
Ambient temperature range: -20°C–+55°C;
Relative humidity range: ≤95%RH, non-condensing;
Atmospheric pressure: 50Kpa–106Kpa.

In the following situations, use the product with extreme caution. It is recommended to consult a doctor before use; otherwise, accidents or physical discomfort may occur:

- People under the age of 18 (children, adolescents, etc.)
- People of 70 years and older
- Patients who are being treated by a doctor and feel abnormal
- Patients with infectious diseases, heart or cranial nerve abnormalities
- Patients with impaired skin perception or abnormal skin
- Patients with a temperature of more than 38 degrees (fever period)
- Patients who are taking medication
- Patients with abnormal blood pressure
- Patients who are unconscious, have severe cognitive disorders, or are not free to express themselves
- Patients with a high level of peripheral circulatory dysfunction caused by (e.g.) diabetes
- Patients with acute suppurative inflammation
- Patients with metal implant
- Patients with spastic paralysis
- Patients with coagulopathy
- Patients with spasticity after a cerebrovascular accident
- Patients with nerve irritability
- People with osteoarticular contracture
- The product has a heating surface, and users who are insensitive to heat shall take care when using this product

Disclaimer: For the above circumstances, please ensure that this product is suitable for your health condition. The company is not responsible for any discomfort or injury caused by failure to comply with this statement and the instructions and recommendations in the User Manual

This device complies with part 15 of the FCC Rules. Operation is subject to the condition that this device does not cause harmful interference (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance

could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure condition without restriction

Manufacturer: Ningbo Ishimei Electronic Technology Co.,

Ltd.Room 1303-5, No. 158 Tai'an Middle Road, Yinzhou District,
Ningbo City, Zhejiang Province, People's Republic of China,Postal
Code: 315100



Made in China