⚠ 警告 WARNING

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The following action must be taken to properly receive the transmitter signal of the correct patient on the receiving instrument. Otherwise, there may be signal loss or signals may mix causing a serious accident, such as monitoring a different patient.

- Assign a channel administrator in the hospital and only he or she should manage channel assignment on his or her responsibility.
- The channel administrator must manage the channels in the facility so that there is no signal interference.
- When the transmitter channel is changed, the channel administrator must check that the channel on the receiving monitor is also changed and the signal is properly received.
- The channel administrator must replace the channel number label on the transmitter with the new one after changing the channel.

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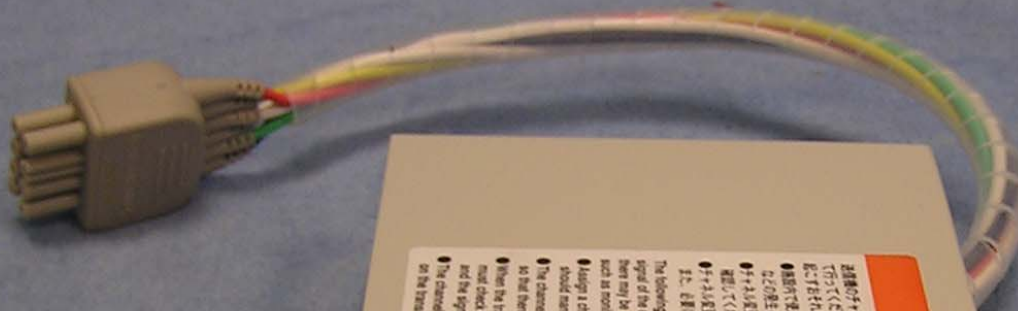
The following action must be taken to properly receive the transmitter signal at the correct patient on the receiving instrument. Otherwise, there may be signal loss or signals may mix causing a serious accident, such as monitoring a different patient.

- Assign a channel administrator in the hospital and only he or she should manage channel assignment on his or her responsibility.
- The channel administrator must manage the channels in the facility as he or she is to signal interference.
- When the transmitter channel is changed, the channel administrator must check that the channel on the receiving monitor is also changed and the signal is properly received.
- The channel administrator must replace the channel number label on the transmitter with the new one after changing the channel.

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