

User Manual

Body Fat Analyzer GBF-2008-B



Guangdong Transtek Medical Electronics Co., Ltd.
Zone A, No.105, Dongli Road, Torch Development District,
528437 Zhongshan, Guangdong, China
TEL:86-760-88282982 www.transtek.cn

- Thank you very much for selecting the Transtek Body Fat Analyzer GBF-2008-B.
- Please do read the user manual carefully and thoroughly so as to ensure the safe usage of this product, and keep the manual well for further reference in case you have problems.

Table of Contents

Safety Information	
Safety and Usage Information	2
Tips	3
Overview	
Device Components	4
List	4
LED Display.....	5
Initial Start-Up	
General Instructions	6
Insert the Batteries	6
Install App and Pair-up	7
Select Measurement Unit	8
Start Measuring	
Daily Measurement	8
Manage Your Health	
Data Transmission	9
Troubleshooting	
Error Prompt	10
When Measuring	10
When Data Transmitting... ..	12
Specifications	13
Maintenance	14
Warranty	14
FCC Regulations	15
Appendix	
Warning	16
EMC Guidance	16

Safety Information

♥ Safety and Usage Information

The warning signs and symbols are essential to ensure your correct and safe use of this product and protect you and others from injury. Please kindly find the meanings of the warning signs and symbols, which you may encounter in the label and user manual, as follows:

	Symbol for "THE OPERATION GUIDE MUST BE READ"		Symbol for "MANUFACTURER"
	Symbol for "RECYCLE"		Caution: These notes must be observed to prevent any damage to the device.
	Symbol for "TYPE BF APPLIED PARTS"		Symbol for "ENVIRONMENT PROTECTION – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice"
	Symbol for "SERIAL NUMBER"		
	Symbol for "DIRECT CURRENT"		
	Symbol for "MANUFACTURE DATE"		

CAUTION

The Body Fat Analyzer GBF-2008-B offers you a seamless way to manage your health. Please be aware that this device is designed for adults' self-measuring and self-monitoring body fat level. Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening. This device should not be used by anyone who is acutely or chronically ill, suffering from a disease or taking medications that affect your water levels. The accuracy of readings for these patients has not been verified. Specific medical advice should be obtained from a physician.

Before diet and exercise, you must first find a professional medical guide.

The body fat analyzer is equipped with data transmission function. It may emit electromagnetic energy so as to perform its intended function. Nearby portable and mobile RF communications equipment can affect the performance of the device.

Portable and mobile RF communications equipment can affect the measuring accuracy of the body fat analyzer.

Kindly note that the use of accessories, transducers or cables other than those specified, with the exception of transducers and cables sold by the manufacturer as replacement parts for internal components, may result in increased EMISSIONS or decreased IMMUNITY of the device.

Be aware that misuse of electrical equipments can cause electric shock, burns, fire and other hazards. Please do not use the device in flammable gas environment.

Warning that the body fat analyzer should not be used adjacent to or stacked with other equipment.

Manufacturer will make available on request circuit diagrams, component parts list, etc.

WARNING: No modifications of this equipment is allowed.

During using the analyzer, the patient will contact with the upper surface and the electrode of the Analyzer. The upper surface and the electrode of the analyzer is type BF Applied part.

The device doesn't need to be calibrated in five years lifetime.

Do not touch the output of batteries when using the analyzer.

The patient is an intended operator. The patient can perform all the operations in the manual, such as measurement, data transmitting, changing batteries.

Safety Information

INDICATIONS FOR USE

- The Transtek Body Fat Analyzer measures weight and uses bio-electrical impedance analysis(BIA) technology to estimate body fat in generally healthy children 10-17 years old and healthy adults.
- It is intended for use in the domestic setting only.

CONTRAINDICATIONS

1. This device is contraindicated for any female subject who may be suspected of, or is pregnant. Besides provided inaccurate readings, the affects of this device on the fetus are unknown.
2. This device is contraindicated for any person who is connected to a wearable or implantable electronic device or instrument such as a pacemaker or defibrillator.

♥ Your Scale and Its Environment

To ensure your safety as well as the service life of your scale, please avoid using the scale under the following circumstances:

- Slippery floor such as tile floor
- Jumping onto the platform immediately after bath or with wet feet
- Near a cell phone or microwave oven

Avoid storage in the following locations:


- Where there is water
- Where the device may be exposed to extreme temperatures, humidity, moisture, direct sunlight, dust, or salt air
- Where there is risk of shock or drop
- Where you store chemicals or full of corrosive gases
- Where in reach of the infants or children, otherwise there is a risk of swallowing batteries.

♥ Efficient Use of Your Scale

To ensure the accuracy of measurement, please follow below instructions when you start measurement.

- Place the scale on a flat, hard surface. Soft surface such as carpet will affect the performance of the scale.
- Step onto the platform with bare feet. Stand still and keep full contact with the scale until the measurement is complete.
- Start measurement at least two hours after Getting up or Dinning.
- Avoid measurement immediately after strenuous exercise, sauna or bath, drinking, and dining.
- Always start measurement in the same time slot and on the same scale located on the same flat, hard surface.
- This device is not intended for use in an oxygen rich environment.


♥ Device Components

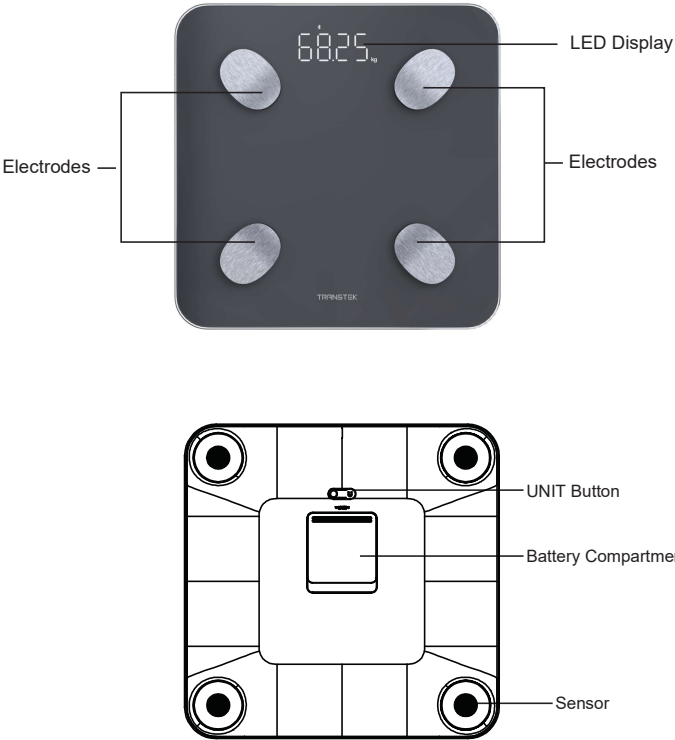

♥ List

- 1. Body Fat Analyzer GBF-2008-B
- 2. Four AAA-size Batteries (1.5V per each)
- 3. User Manual

♥ LED Display

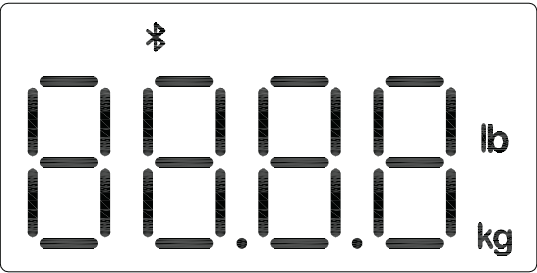

lb	Pound	Bluetooth symbol	Bluetooth Connection Indication
kg	Kilogram		

Initial Start-Up

♥ General Instructions

Transtek Body Fat Analyzer GBF-2008-B applies BIA (Bio-impedance Analysis) technology. A small amount of weak current flows through the human body so as to detect the bio-impedance and estimate body fat. The electrical current is small and may not be felt.

This BIA technology is cheap, safe, non-invasive, toxic-free and harmless. It also possesses the characteristics of simple operation and abundant information.

The current mentioned above is less than 0.5mA. However, please be aware that anyone with an wearable or implantable medical electronic instrument, such as a pacemaker, must avoid using this device.

The intended use of this device is for healthy children 10-17 years old and healthy adults.


♥ Insert the Batteries

- Open the battery door in the back of the scale.
- Insert the batteries (4 x 1.5V AAA) into the battery compartment according to the polarity indications marked inside the compartment.

* The digits "8888" will be shown on the LED.



- Close the battery door and wait until the digits "00.00" are shown on the LED.

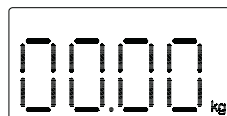

CAUTION

- When the symbol "Lo" appears, the device will power off in about 4s. Then you shall replace with a new set of batteries. Please replace all four batteries at the same time. Do NOT mix the old batteries with the new one.
- Worn batteries are hazardous waste. Do NOT dispose of them together with the household garbage.
- Please refer to the local ordinances and recycling instructions regarding disposal of the worn batteries and scrapped device.
- If you do not intend to use this unit for a prolonged period of time, it is advisable to remove the batteries before storing.
- Do not mix alkaline, standard(carbon-zinc), or rechargeable(Ni-Cd, Ni-MH, etc) batteries.
- If you do not intend to use this unit for a prolonged period of time, it is advisable to remove the batteries before storing.
- Use the scale three times a day, and the lifetime of battery is one year.

Initial Start-Up

♥ Install APP and Pair-Up


With the advanced Bluetooth 5.0 technology applied, the mobile or portable equipments, which are equipped with Bluetooth function in line with BLE Technical Specifications as well as BLP Protocol established by global organization Bluetooth SIG, are capable to receive your personal health data.

Just simply install the Tuya Smart App and pair up your scale with your mobile or portable equipments. Then you may enjoy the comprehensive health solution provided by Tuya.


Note: It is only available for the device which supports Bluetooth 5.0 technology.

- For the first time, the scale will automatically enter bluetooth pairing mode.
- Turn on bluetooth and the Tuya Smart App. Make sure both are ON when pair-up is proceeding.

If SUCCEED, symbol "✱" will be shown on the LED.

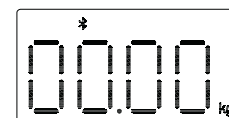

If FAIL, symbol "✱" will not be shown on the LED.

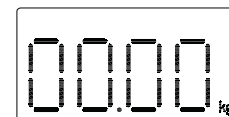


If symbol "✱" blinks, it means the device is being searched.

Start Measuring


♥ Select Measurement Unit

With batteries correctly installed, press "UNIT" button in the back of the scale to select measurement unit. The default measurement unit is "kg". You may press "UNIT" button to choose between kilogram and pound.

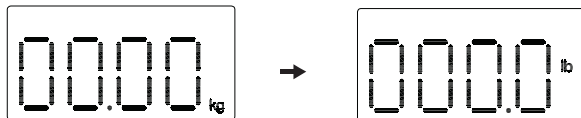

♥ Daily Measurement

- With original SENSE ON patent technology, GBF-2008-B will switch on as you step on the platform barefooted.
- Stand still and keep full contact with the electrodes until the LED stops blinking digits.

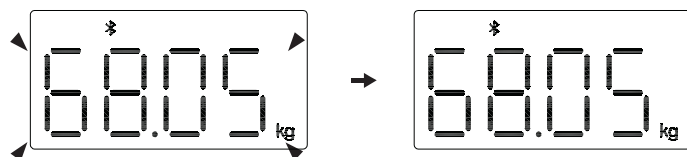

Manage Your Health

♥ Data Transmission

- With GBF-2008-B successfully pair-up with your smartphone, the measurement data will be transmitted to your mobile via Bluetooth.



- ONLY when the measuring results are attached to a specific User ID will it be transmitted to your mobile after measurement.
- The symbol  will display after successful data transmission, and you may check your personal health data stored in your smartphone.
- If the data transmission fails, the symbol  doesn't display. The pending measurement data will be temporarily kept in the scale and transmitted to your smartphone when next measurement is complete.

CAUTION


- Interference may occur in the vicinity of equipment marked with the following symbol . And the Analyzer may interfere vicinity electrical equipment.
- To enable the data transmission function, this product should be paired to a Bluetooth 5.0 or above end.

How to mitigate possible interference?

1. The range between the Analyzer and the Bluetooth end should be reasonably close, from 1 meter to 10 meters. Please ensure no obstacles between the Analyzer and the Bluetooth end so as to obtain quality connection.
2. To avoid interference, other electronic devices (particularly those with Bluetooth transmission / Transmitter) should be kept at least 1 meter away from the Analyzer.

Troubleshooting

♥ Error Prompt

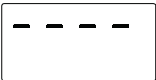Error	Description	Solution
	Overload.	Stop using this scale for measurement.
	Low Battery.	Replace all four batteries in the same time. Please purchase the authorized batteries for replacement.

♥ When Measuring ...

Problem	Root Cause	Solution
Abnormal measuring results: - Too high; - Too low; - Huge difference between two recent measurement.	Incorrect posture	Please step on the platform barefooted and stand still.
	The device is located on the soft ground such as a carpet or on a rugged surface.	Please place the device on a flat, hard surface.

Troubleshooting

Problem	Root Cause	Solution
Abnormal measuring results: - Too high; - Too low; - Huge difference between two recent measurement.	Your feet are too dry.	Wipe your feet with a damp cloth, keeping them slightly damp when starting measurement.
No display on LED when the device powers on.	Batteries not yet installed.	Install the batteries. (Please refer to Insert the Batteries)
	Worn batteries.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.
CANNOT proceed to analyze body fat.	Step onto the platform wearing socks or shoes.	Please keep barefooted during the measurement, and keep full contact with the electrodes as well.
	The system cannot identify the possible User ID with most similar data.	Please assign a User ID again on the APP.
The device powers off.	Low battery.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.

Troubleshooting

♥ When Data Transmitting ...

Problem	Root Cause	Solution
Data transmission failed.	Bluetooth is OFF.	Turn ON the Bluetooth via Settings.
	The APP is OFF.	Press the icon to turn ON your APP.
	Out of range of Bluetooth transmission.	Place your smartphone closer to the scale.
	None of the user ID is assigned.	Please assign a User ID on the APP.

About the Accuracy of This Product

- This product passes strict inspection before delivery and therefore its accuracy is guaranteed by the manufacturer. Please refer to the above table for the descriptions on accuracy.
- This product is specially designed for body fat analysis as well as weight measurement. It should NOT be used by anyone during the process of transaction for verification of goods' weight.

Specifications

♥ Specifications

Product Name	Body Fat Analyzer (GBF-2008-B)
Dimension	Scale: 300x300x25.5mm (Approximately)
Net Weight	Approximately 1.7kg (Including the dry cells)
Display	Digital LED
Measurement Unit	Pound / Kilogram
Measurement Range	11lb to 397lb / 5kg to 180kg
Division	5-100kg: 0.05kg / 0.1lb; 100-180kg: 0.1kg / 0.2lb
Accuracy	5-50: ±0.2kg; 50.05-100: ±0.3kg; 100.1-150: ±0.4kg; 150.1-180: ±0.5kg;
Working Environment	Temperature: 5℃ to 40℃ Relative Humidity: ≤90% RH Atmospheric pressure: 86kPa to 106kPa
Storage Environment	Temperature: -20℃ to 60℃ Relative Humidity: ≤90% RH Atmospheric pressure: 50kPa to 106kPa
Power Source	6V (Four AA-size Batteries)
Auto-ON	SENSE ON technology
Auto-OFF	The scale will turn off after about 30s if there is no operation.
Accessories	1. Four AA-size Batteries 2. User Manual
Mode of Operation	Continuous Operation
Protection Against Ingress of Water	IP21(The first number 2: Protected against solid foreign objects of 12.5 mm Φ and greater. The second number 1: Protection against vertically falling water drops.)
Software Version	V1.4
Degree of protection	Type BF applied part

♥ Maintenance

When carrying out usual maintenance, please ensure practice of the following Do's and Don'ts:

- DO use a dry soft cloth to wipe the dust.
- DO use a wet soft cloth, dipped into water and wrung out, to wipe the dirt. Then use a dry soft cloth to dry up the device.
- DON'T wash the device with water or immerse it in water.
- DON'T use propellant, abrasive or other chemicals to wipe the dirt in avoidance of discolor or malfunction.
- DON'T disassemble this device. If you have any problems, please contact Transtek.
(Please refer to Warranty for contact information)

♥ Warranty

- Transtek warrants its products free of defects in materials and workmanship in normal use for a period of FIVE years from the date of retail purchase.
- This warranty does NOT cover damages caused by misuse or abuse, including but not limited to:
 - Failure caused by unauthorized repairs or modifications;
 - Damage caused by shock or drop during transportation;
 - Failure caused by improper operation inconsistent with the instructions stated in this user manual;
 - Malfunction or damage from failure to provide the recommended maintenance;
 - Damage caused by improper use of power supply.
- Should this device require maintenance (or replacement at our option) under warranty, please deliver the original package to GUANGDONG TRANSTEK MEDICAL ELECTRONICS CO., LTD prepaid. Please return the store receipt (with the retail purchase date) and a note with reasons to return on it as well.

Guangdong Transtek Medical Electronics Co., Ltd.
Zone A, No.105, Dongli Road, Torch Development District, 528437 Zhongshan, Guangdong, China
Tel: 86-760-88282982
Website: <http://www.transtek.cn>

♥ FCC Regulations

FCC User Guide Information

Radio Frequency Interface Requirements - FCC

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna;
- Increase the separation between the equipment and receiver;
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected;
- Consult the dealer or an experienced radio / TV technician for help.

Radio Transmitters (Part 15)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC RF Exposure Guidelines

Safety Information

Reducing RF Exposure - Use Properly

Only operate the device in accordance with the instructions supplied.
This device complies with FCC radiation exposure limits set forth for an uncontrolled environment.

♥ Warning

Warning: Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation."

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment GBF-2008-B, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

♥ EMC Guidance

Guidance and manufacturer's declaration -electromagnetic emissions and Immunity

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Not applicable
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable

Table 2

Guidance and manufacturer's declaration - electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	Not applicable	Not applicable
Surge IEC 61000-4-5	Not applicable	Not applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Not applicable	Not applicable
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	Not applicable	Not applicable
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz
NOTE U _T is the a.c. mains voltage prior to application of the test level.		

Appendix

Table 3

Guidance and manufacturer's declaration - electromagnetic Immunity							
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
	385	380 – 390	TETRA 400	Pulse modulation 18 Hz	1,8	0.3	27
	450	430 – 470	GMRS 460, FRS 460	FM ± 5 kHz deviation 1 kHz sine	2	0.3	28
	710	704 – 787	LTE Band 13, 17	Pulse modulation 217 Hz	0,2	0.3	9
	745						
	780						
	810	800 – 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18 Hz	2	0.3	28
	870						
	930						
	1720	1 700 – 1 990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	2	0.3	28
	1845						
	1970						
	2450	2 400 – 2 570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28



Guangdong Transtek Medical Electronics Co., Ltd.


广东乐心医疗电子股份有限公司

Document Name: 文件名称:	USER MANUAL
Document No 文件编号:	LS-USER MANUAL-GBS-2012-B
Rev.: 版本号:	A/1
Page: 页码:	Total 17 pages

Prepared 编制		Reviewed 审核		Approved 批准	
Coco Chen		Aimee Li		Rachel	
◆ Revision History 修订履历					
NO.	Date 日期	Update Description 变更描述	Prepared 编制	Reviewed 审核	Approved 批准
A/0	2023/12/22	Initial release	林丽丽 Lily Lin	黎萍芳 Aimee Li	徐托 Alex Xu
A/1	2026/01/30	Update the manufacturer's name.	陈伟璇 Coco Chen	黎萍芳 Aimee Li	胡丽琼 Rachel Hu

TeleRPM

Body Scale
GBS-2012-B

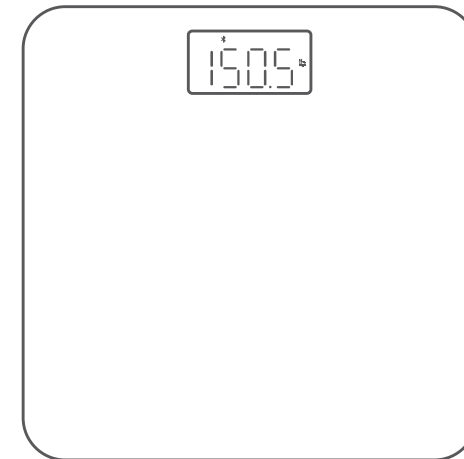

User Manual



Guangdong Transtek Medical Electronics Co., Ltd.
Zone A, No.105, Dongli Road, Torch Development District,
528437 Zhongshan, Guangdong, China

version: A/1

- Thank you for choosing Transtek Body Scale!
- Please read this manual thoroughly to know how to use your Body Scale safely and correctly. Keep the manual well for future reference.

Table of Contents

Symbol Meaning	2
Safety Information	
Your Scale and Its Environment	3
Efficient Use of Your Scale	3
Overview	
Device Components	4
List	4
LED Display	5
Initial Start-Up	
General Instructions	6
Insert the Batteries	6
Install APP and Pair-Up	7
Select Measurement Unit	7
Start Measuring	
Daily Measurement	8
Troubleshooting	
Error Prompt	9
When Measuring	9
When Data Transmitting	10
Specifications	11
Maintenance	12
Warranty	12
FCC Regulations	13
Appendix	
Warning	14
EMC Guidance	14

Symbol Meanings

The symbols below may be in user manual, labeling or other component. They are the standard requirement of using.

	Refer to instruction manual/booklet To signify that the instruction manual/ booklet must be read.		Manufacturer
	Serial Number		Date of Manufacture
	Direct Current		General symbol for recovery/recyclable
	Type BF applied part		Temperature limit
	Humidity limitation		Atmospheric pressure limitation
	Caution Indicates that caution is necessary when operating the device or control close to where the symbol is placed, or that the current situation needs operator awareness or operator action in order to avoid undesir- able consequences.		The symbol indicates that the product should not be discarded as unsorted waste but must be sent to separate collection facilities for recovery and recycling.

CAUTION

The Body Scale offers you a way to manage your health. Please be aware that this device is designed for children (10–17 years old) and adults' self-measuring. Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening.

The Body Scale is equipped with data transmission function. It may emit electromagnetic energy so as to perform its intended function. Nearby portable and mobile RF communications equipment can affect the performance of the device.

Kindly note that the use of accessories, transducers or cables other than those specified, with the exception of transducers and cables sold by the manufacturer as replacement parts for internal components, may result in increased EMISSIONS or decreased IMMUNITY of the device.

Be aware that misuse of electrical equipments can cause electric shock, burns, fire and other hazards. Please do not use the device in flammable gas environment.

Warning that the Body Scale should not be used adjacent to or stacked with other equipment.

WARNING: No modifications of this equipment is allowed.

The top surface of product is considered as Type BF applied part.

The device do not need to be calibrated in five years lifetime.

Do not touch the output of batteries when using the scale.

The patient is an intended operator. The patient can perform all the operations in the manual, such as measurement, changing batteries, no third-party operator is required.

Please use the device under the environment which is provided in the user manual. Otherwise, the performance and lifetime of the device will be impacted and reduced.

Store your device in a clean and dry place. Never place any heavy objects on it.

Not intended for use in oxygen rich environment, Not intended for use with flammable anaesthetics, not intended for use in conjunction with flammable agents

Dispose of accessories, detachable parts, and the device according to the local guidelines.

Degeneration or looseness of sensors and electrodes might reduce accuracy.

Stop using the Scale if any discomfort or allergic reaction occurs.

Safety Information

INDICATIONS FOR USE

- The Transtek Body Scale is intended for adults and children (10–17 years old) to measure weight at home.

Your Scale and Its Environment

To ensure your safety as well as the service life of your scale, DO NOT use the scale:

- on the slippery floor such as tile floor,
- immediately after bath or with wet feet,
- near a cell phone or microwave oven.

Avoid storage in the following locations:

- Where there is water.
- Where the device may be exposed to extreme temperatures, humidity, moisture, direct sunlight, dust, or salt air.
- Where there is risk of shock or drop.
- Where you store chemicals or full of corrosive gases.
- Where in reach of the infants or children.


Efficient Use of Your Scale

To ensure the accuracy of measurement, please follow below instructions when you start measurement.

- Place the scale on a flat, hard surface. Soft surface such as carpet will affect the performance of the scale.
- Start measurement at least two hours after getting up or dinning.
- Avoid measurement immediately after strenuous exercise, sauna or bath, drinking, and dining.
- Always start measurement in the same time slot and on the same scale located on the same flat, hard surface.
- This device is not intended for use in an oxygen rich environment.
- Any information provided by this scale is in no way meant to treat, cure or prevent any disease or illness from happening. If in doubt, contact your physician.
- The device is designed for use in the home healthcare environment.

Overview

Device Components

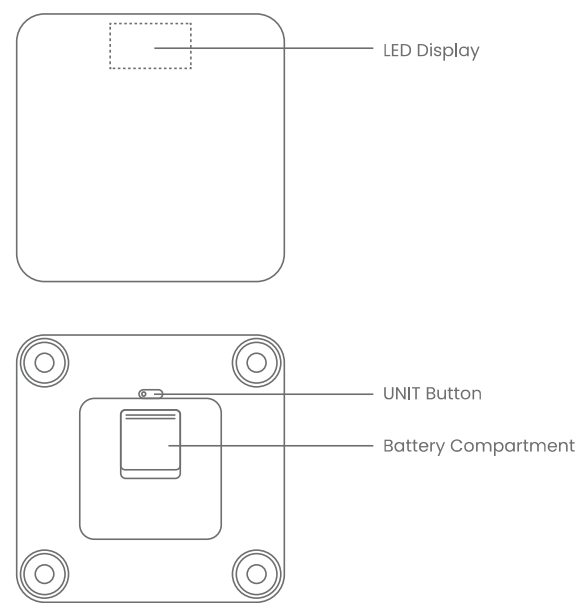


List

- 1. Body Scale GBS-2012-B
- 2. Four AAA-size Batteries (installed in the device)
- 3. User Manual

Overview

LED Display

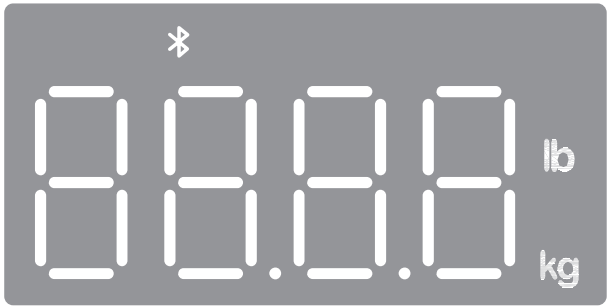

lb	Pound	Bluetooth symbol	Bluetooth Connection Indication
kg	Kilogram		

Initial Start-Up

General Instructions

Transtek Body Scale GBS-2012-B applies SENSE ON technology. It's not necessary to turn on or turn off the scale, just stand on the scale for measurement. It is delicate and trendy, convenience for using.

Insert the Batteries

- Open the battery door in the back of the scale.
- Remove insulating strip from beneath the battery (if fitted) or insert batteries (4 x 1.5V AAA) observing the polarity signs (+ and -) inside the battery compartment. The all digits will be shown on the LED display.



- Close the battery door and wait until the digits "0.0lb" are shown on the LED display.



CAUTION

- When symbol "⚡" shows, the device will power off. Then you shall replace with a new set of batteries. Please replace all four batteries at the same time. Do NOT mix the old batteries with the new ones.
- Worn batteries are hazardous waste. Do NOT dispose of them together with the household garbage.
- Please refer to the local ordinances and recycling instructions regarding disposal of the worn batteries and scrapped device.
- Do not mix alkaline, standard (carbon-zinc), or rechargeable (Ni-Cd, Ni-MH, etc) batteries.
- If you do not intend to use this unit for a prolonged period of time, it is advisable to remove the batteries before storing.
- Use the scale three times a day, and the lifetime of battery is one year.

Install APP and Pair-Up

You can measure your weight and send measurement data to a smartphone with Bluetooth wireless connectivity and MedM Health application.


1. Download the MedM Health app from APP Store or Google Play, install the APP, and register an account. Then set your personal information (Gender, Birthday, Height, Name and so on).
2. Turn on Bluetooth and the MedM Health application. Make sure both are ON when pair-up is proceeding.
3. When the scale is powered on for the first time, the symbol "⚡" on the LED display blinks, indicating pair-up is proceeding.
4. If pair-up is successful, the Bluetooth symbol will stop blinking.

Note: Hold the unit button for 5 seconds and the pair-up process will be restarted.



Select Measurement Unit

With batteries correctly installed, press "UNIT" button in the back of the scale to select measurement unit. The default measurement unit is "lb". You may press "UNIT" button to choose between pound and kilogram.

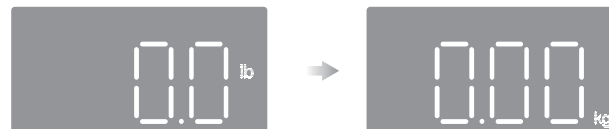

Start Measuring

Daily Measurement

- With original SENSE ON patent technology, the scale will automatically switch on as you step on the platform barefooted.



- Stand still and keep full contact with the scale, the weight data on the LED display will blink twice, then the weight is locked.
- When your scale is successfully connected to your smartphone, it will process data transmission automatically. After the data transmission, the LED display goes off automatically.
- If the data transmission fails, the pending measurement data will be temporarily kept in the scale and transmitted to your smartphone when next measurement is complete.



Troubleshooting

Error Prompt

Error	Description	Solution
	Overload. The device will power off.	Stop using this scale for measurement.
	Low Battery. The device will power off soon.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.

When Measuring ...

Problem	Root Cause	Solution
No display on LED when the device powers on.	Batteries not yet installed.	Install the batteries. (Please refer to Insert the Batteries)
	Worn batteries.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.
The device powers off automatically.	Low battery.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.

Troubleshooting

When Data Transmitting ...

Problem	Root Cause	Solution
Data transmission failed.	Bluetooth is OFF.	Turn ON the Bluetooth via Settings.
	App is OFF.	Press the icon to turn ON your app.
	Out of range of Bluetooth transmission.	Place your smartphone closer to the scale.

CAUTION

- Interference may occur in the vicinity of equipment marked with the following symbol . And the Scale may interfere vicinity electrical equipment.
- To enable the data transmission function, this product should be paired to a Bluetooth 5.0 or above end.
How to mitigate possible interference?
 - The range between the Scale and the Bluetooth end should be reasonably close, from 1 meter to 10 meters. Please ensure no obstacles between the Scale and the Bluetooth end so as to obtain quality connection.
 - To avoid interference, other electronic devices (particularly those with Bluetooth transmission / Transmitter) should be kept at least 1 meter away from the Scale.

About the Accuracy of This Product

- This product passes strict inspection before delivery and therefore its accuracy is guaranteed by the manufacturer. Please refer to "Specifications" chapter for the descriptions on accuracy.
- This product is specially designed for human weight measurement. It should NOT be used by anyone during the process of transaction for verification of goods' weight.
- At least 2 hours required for the device to warm / cool from the minimum / maximum storage temperature between uses until it is ready for intended use when the ambient temperature is 20°C.

Specifications

Product Name	Body Scale
Dimension	300x300x25.1mm
Net Weight	Approximately 1.67kg (Excluding the dry cells)
Display	Digital LED
Measurement Unit	Pound/Kilogram
Measurement Range	11 lb to 440 lb / 5 kg to 200 kg
Division	5~100kg, 0.05kg/0.2lb; 100.1~200kg, 0.1kg/0.2lb
Accuracy	Temperature requirement: 25±5°C 50 ± 0.2 kg; 100 ± 0.3 kg; 150 ± 0.4 kg; 200 ± 0.5 kg;
Working Environment	Temperature: 5°C to 40°C Relative Humidity: ≤90% RH Atmospheric pressure: 86kPa to 106kPa
Storage Environment	Temperature: -20°C to 60°C Relative Humidity: ≤90% RH Atmospheric pressure: 50kPa to 106kPa
Power Source	6V (Four AAA-size Batteries)
Auto-ON	SENSE ON technology
Auto-OFF	About 15 seconds while showing 0.0; About 15 seconds after the weight data is locked;
Accessories	1. Four AAA-size Batteries 2. User Manual
Mode of Operation	Continuous Operation
Degree of protection	Type BF applied part
Software Version	Version 1.4
Protection Against Ingress of Water	IP21

Maintenance

When carrying out usual maintenance, please ensure practice of the following Do's and Don'ts:

- Clean the device after use.
- DO use a dry soft cloth to wipe the dust.
- DO use a wet soft cloth, dipped into water and wrung out, to wipe the dirt. Then use a dry soft cloth to dry up the device.
- DON'T wash the device with water or immerse it in water.
- DON'T use propellant, abrasive or other chemicals to wipe the dirt in avoidance of discolor or malfunction.
- Do not attempt to repair the device or replace components (e.g., transducers, LED, etc.) by yourself. If you suspect incorrect operation, contact customer support for further assistance. At the request of authorized service personnel, circuit diagrams, component part lists, descriptions, and calibration procedures will be made available by the manufacturer or distributor.
- Do Not disassemble this device. If you have problems, contact customer support (Please refer to the section of "Warranty" for customer support information)
- No service and maintenance while in use.

Warranty

- Transtek warrants its products free of defects in materials and workmanship in normal use for a period of TWO years from the date of retail purchase.
- This warranty does NOT cover damages caused by misuse or abuse, including but not limited to:
Failure caused by unauthorized repairs or modifications;
Damage caused by shock or drop during transportation;
Failure caused by improper operation inconsistent with the instructions stated in this user manual;
Malfunction or damage from failure to provide the recommended maintenance;
Damage caused by improper use of power supply.
- Should this device require maintenance (or replacement at our option) under warranty, please deliver the original package to GUANGDONG TRANSTEK MEDICAL ELECTRONICS CO., LTD prepaid. Please return the store receipt (with the retail purchase date) and a note with reasons to return on it as well.

Guangdong Transtek Medical Electronics Co., Ltd.
Zone A, No.105, Dongli Road, Torch Development District, 528437 Zhongshan, Guangdong, China
Website: www.transtekcorp.com

FCC Regulations

FCC User Guide Information

Radio Frequency Interface Requirements – FCC

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna;
- Increase the separation between the equipment and receiver;
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected;
- Consult the dealer or an experienced radio / TV technician for help.

Radio Transmitters (Part 15)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC Radiation Exposure Statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Appendix

Warning

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments.

- Warning: Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.
- Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
- Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.
- Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

EMC Guidance

Guidance and manufacturer's declaration–electromagnetic emissions and Immunity.

Table 1

Guidance and manufacturer's declaration – electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Not applicable
Voltage fluctuations/flicker emissions IEC 61000-3-3	Not applicable

Appendix

Table 2

Guidance and manufacturer's declaration–electromagnetic immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	Not applicable	Not applicable
Surge IEC 61000-4-5	Not applicable	Not applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Not applicable	Not applicable
Power frequency magnetic field	30 A/m 50 Hz/60 Hz	30 A/m 50 Hz/60 Hz
IEC 61000-4-8		
Conducted RF IEC 61000-4-8	Not applicable	Not applicable
Radiated RF IEC 61000-4-3	10 V/m 80 MHz-2.7 GHz 80% AM at 1 kHz	10 V/m 80 MHz-2.7 GHz 80% AM at 1 kHz
NOTE U T is the a.c. mains voltage prior to application of test level		

Appendix

Table 3

Guidance and manufacturer's declaration – electromagnetic immunity								
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Maximum Power (W)	Distance (m)	IEC 60601-1-2 Test Level (V/m)	Compliance level (V/m)
	385	380–390	TETRA 400	Pulse modulation 18 Hz	1.8	0.3	27	27
	450	430–470	GMRS 460, FRS 460	FM ±5k Hz deviation 1 kHz sine	2	0.3	28	28
	710	704–787 800–960	LTE Band 13, 17 GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 217 Hz Pulse modulation 18 Hz	0.2 2	0.3 0.3	9 28	9 28
	745							
	780							
	810							
	870							
	930							
	1720	1700–1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	2	0.3	28	28
	1845							
	1970							
	2450	2400–2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28	28
	5240	5100–5800	WLAN, 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9	9
	5500							
	5785							