MD10W2

10.1" Medical Android Tablet

ARM Genio 700, Octa-core, up to 2.2 GHz

User Manual

Version 1.0

Document Part No. 91521110118I

Contents

Chapter 1: Read Me First	3
1.1 Advisory Conventions	3
1.2 Intended Use	4
1.3 Safety Information	5
1.4 Handling and Maintenance	8
1.5 Declaration of Conformity (DoC)	9
1.6 Copyright Notice	10
1.7 Trademark Acknowledgement	10
1.8 Disclaimer	10
Chapter 2: Getting Started	11
2.1 Introduction	11
2.2 Unpacking	11
2.3 Description of Parts	12
2.4 First-Time Operation	13
2.4.1 Inserting Battery Pack	14
2.4.2 Inserting MicroSD Card	15
2.4.3 Connecting to Power Supply	15
2.5 Turning On and Off the Device	16
2.5.1 Turning On	16
2.5.2 Turning Off	16
2.6 Resetting to Factory Defaults	16
Chapter 3: Using the Tablet	20
3.1 Using the Touch Screen	20
3.2 Waking Up the Tablet	21
3.3 Using the Home Screen	22
3.3.1 Personalizing Home Screen	23
3.3.2 Viewing other parts of the Home Screen	26
3.4 Status and Notification Icons	28
3.4.1 Notification Icons	28
3.4.2 Status Icons	28
3.5 Using the On-Screen Keyboard	29
3.5.1 Entering Text	29
3.5.2 Editing Text	31
3.6 Display Settings	34
3.7 Audio Profile Settings	35
3.7.1 Sound Settings	35
3.7.2 Adjusting Volume	37
3.8 Language and Keyboard Settings	38

3.9 Changing Touchscreen Modes	42
3.10 Extreme power saving setting	44
3.11 About Phone	46
Chapter 4: Winmate Kiosk Manager	47
4.1 Lock Applications Launcher	47
4.2 Unlock Applications Launcher	53
4.3 How to Change Password	56
4.4 White List App	57
Chapter 5: Networks and Wireless Devices	60
5.1 Connecting to WLAN Networks	60
5.1.1 Turning on WLAN and Connecting to WLAN Networks	60
5.1.2 Adding WLAN Network	63
5.2 Connecting to BT Devices	65
5.2.1 Turning BT On or Off	65
5.2.2 Connecting the BT Device	66
5.2.2 Disconnecting BT Device	67
5.3 Connecting to Computer via USB	68
5.3.1 USB Computer Connection	68
5.3.2 USB Debugging Settings	69
Chapter 6: Using Camera	73
6.1 Capturing Photos or Video	73
6.1.1 Opening Camera	73
6.1.2 Taking Pictures	73
6.1.3 Viewing pictures	74
6.1.4 Shooting Videos	74
6.2 The Gallery	74
Appendix	75
Appendix A: Basic Troubleshooting	75

Chapter 1: Read Me First

1.1 Advisory Conventions

Four types of advisories are used throughout the user manual to provide helpful information or to alert you to the potential for hardware damage or personal injury. These are Notes, Important, Cautions, and Warnings. The following is an example of each type of advisory.

**Note:**

A note is used to emphasize helpful information

**Important:**

An important note indicates information that is important for you to know.

Caution A Caution alert indicates potential damage to hardware and explains how to avoid the potential problem.

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Warning! An Electrical Shock Warning indicates the potential harm from electrical hazards and how to avoid the potential problem.

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1.2 Intended Use

1. Intended Use Statement

The Winmate MD10W2 is a rugged medical-grade tablet computer intended to serve as an accessory to Medical Electrical (ME) Equipment or an ME System.

Its primary function is to act as a **Human-Machine Interface (HMI)** and **Data Processing Unit** for:

1. **Visualizing Clinical Data:** Displaying Electronic Medical Records (EMR), medical imaging (PACS), vital signs, and nursing workflows.
2. **System Control:** Interfacing with and controlling connected active medical devices (e.g., mobile medical carts, diagnostic stations, or patient monitoring systems).

Note: The device itself does not perform direct automated diagnosis or therapeutic functions. It is intended to support healthcare professionals in clinical decision-making by providing a reliable platform for medical software applications.

2. Detailed Medical Indication & Specification

A. Intended Medical Indication

- **Device Classification:** Active Medical Device Accessory / Medical Data Display.
- **Indication for Use:** The device is indicated for use in general clinical workflows, data entry, and the reviewing of medical information. It supports the management of patient data for diseases ranging from chronic conditions (e.g., diabetes management software) to acute care monitoring.
- **IEC 62304 (Software) Consideration:** The device provides a verified **Operating Environment** (Android OS) for **Software as a Medical Device (SaMD)**. The medical software running on the device determines the specific clinical indication (e.g., reviewing an X-ray for bone fractures).

B. Patient Population

- **Target Population:** The device is used to manage data for patients of all ages (neonatal, pediatric, adult, geriatric), genders, and weight categories.
- **Direct Interaction:** The device is **not** intended for direct use on the patient. Patient contact is considered incidental.

C. Part of the Body or Type of Tissue Interacted With)

- **Applied Parts: None.** The device is defined as having **No Applied Parts**.
- **Patient Environment:** The device is designed to be safe for use within the **Patient Vicinity** (within 1.5m of the patient). It meets **MOPP (Means of Patient Protection)** standards for electrical isolation to prevent leakage current hazards if accidentally touched by a patient or by a nurse touching the patient.
- **Accessible Parts:** The touchscreen and enclosure are "Accessible Parts" intended for the operator (Healthcare Provider) and meet **MOOP (Means of Operator Protection)** standards.

D. User Profile

- **Intended Users:**
 - **Clinical Professionals:** Physicians, Nurses, Surgeons.

- **Allied Health Professionals:** Radiologists, Lab Technicians.
- **Biomedical Engineers:** IT staff responsible for device configuration and maintenance.
- **Competence:** Users are expected to have professional medical training and basic computer literacy.
- **Lay Users:** The device is **not** intended for use by patients or laypersons in a home healthcare setting unless specifically configured and validated for a "Home Healthcare Environment" under IEC 60601-1-11.

E. Use Environment

- **Professional Healthcare Facility Environment:** Hospitals, Intensive Care Units (ICU), Operating Rooms (OR), Emergency Rooms (ER), and Clinics.
- **Electromagnetic Compatibility:** The device is designed to function in environments with typical hospital RF interference (e.g., near electrosurgical units or MRI control rooms), compliant with **CISPR 11 Class B** (or A) and **IEC 60601-1-2 Edition 4.0** immunity levels.

3. Risk Management & Safety Limitations

- **Contraindication (Life Support):** This device must **not** be used as the *sole* monitoring display for life-critical alarms (e.g., ventilator disconnection) where a delay in software processing or Wi-Fi transmission could result in patient injury or death. It is a *secondary* or *adjunct* viewer.
- **Latency Warning:** Real-time waveform transmission via WLAN/Bluetooth may be subject to environmental latency. Clinical judgment should not rely exclusively on wireless real-time data without verifying the primary monitor.
- **Hygiene & Sterilization:** The device supports cleaning with hospital-grade disinfectants (e.g., 70% Isopropyl Alcohol, Sodium Hypochlorite). It is **not autoclavable**. Do not immerse the device in liquid unless all ports are sealed and the IP rating (IP65) is verified intact.
- **Cybersecurity (IEC 81001-5-1):** As a network-connected device, the Responsible Organization (Hospital IT) must enforce security policies (e.g., "Kiosk Mode", strong passwords, disabled USB debugging) to prevent unauthorized access to Protected Health Information (PHI).

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Do Not Disassemble, Modify, or Insert Sharp Objects Into the Battery Pack. Electrolyte leakage, generation of heat, ignition or rupture may result.

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Do Not Short the Positive (+) and Negative (-) Terminals of the Battery. Generation of heat, ignition or rupture may result. When carrying or storing the device, do not place the battery pack together with articles that may contact electricity.

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Do Not Apply Shocks to the Battery by Dropping It or Applying Strong Pressure to It. Electrolyte leakage, generation of heat, ignition or rupture of battery may occurs.

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Do Not Charge the Battery Other Than Methods Those Specified In This Manual. If the battery is not charged using one of the specified methods, electrolyte leakage, generation of heat, ignition or rupture may result.

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When the Battery Pack Has Deteriorated, Replace It with a New Battery. Continued use of a damaged battery pack may result in heat generation, ignition or battery rupture.

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Do Not Expose the Battery Pack to Excessive Heat, or Extreme Heat (Near Fire, in Direct Sunlight for example) Generation of heat, ignition or rupture may result

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Do Not Use the Battery Pack With Any Other Products. The battery pack is rechargeable and is only designed for this specific product. If it is used with a product other than the specified product, electrolyte leakage, generation of heat, ignition or rupture may result.

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Do Not Use This Product With Battery Other Than the One Specified.

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Discontinue using device immediately if you notice the battery is swelling, or looks puffy or misshapen. Do not attempt to charge or use the device. Do not handle a damaged or leaking battery. Do not let leaking battery fluid come into contact with your eyes, skin or clothing.

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Do not block any ventilation openings. Do not place the product against a wall or in an enclosed space as doing so might impede the flow of cooling air.

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Do not place the product near any heat sources such as radiators, warm air intakes, stoves or other heat-generating devices (including amplifiers or other heat producing equipment).

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Protect the power cord from being crushed or pinched, particularly in the area of the plugs, any receptacles used, and the point of exit from the apparatus. Do not break off the ground pin of the power cord.

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Do not disregard the safety that a polarized or grounded connection can provide. A polarized connection consists of two plugs, one of which is larger than the other. A grounded connection has two blades and a third one that is grounded. If the main plug provided cannot be inserted into your outlet, consult an electrician to replace your obsolete outlet.

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Use only accessories specified by the manufacturer.

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This product should not be used near water such as bathtubs, sinks, swimming pools, wet basement, etc.

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Use only with its transport carriage, stand, tripod, console or table specified by the manufacturer or sold with the product. When using a carriage, move the carriage assembly with the product carefully to avoid damage from falling over.

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Unplug this apparatus during lightning storms or when it is not being used for long periods.

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Refer all servicing to qualified service personnel. A technical intervention is required when the apparatus has been damaged in any way, for example if the power cord or plug is damaged, liquid has been spilled or objects have fallen inside the unit, if the product has been exposed to rain or moisture, if it does not operate normally or if it has been dropped.

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General Guideline **Directive générale**

It is recommended to reboot the device when some functions are defect or inactive. If it still can't solve the problems, please contact your dealer or agent.

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1.4 Handling and Maintenance

General Guideline *Guide Général*

Keep the product at least 13 cm (5 inches) away from electrical appliances that can generate a strong magnetic field such as TV, refrigerator, motor, or a large audio speaker.

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Do not place heavy objects on top of the tablet computer as this may damage the display.

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The screen surface can be easily scratched. Use with care. Do not use paper towels to clean the display but use the included cleaning cloth. Avoid touching it with a pen or pencil.

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To maximize the life of the backlight in the display, allow backlight to automatically turn off as a result of power management. Avoid using a screen saver or other software that prevents the power management from working.

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Cleaning Guideline *Guide de nettoyage*

Use a soft cloth moistened with water or a non-alkaline detergent to wipe the exterior of the product.

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Gently wipe the display with a soft, lint-free cloth. Do not use alcohol or detergent on the product.

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Use your finger or the stylus on the display. Using sharp or metallic object may cause scratches and damage the display, thereby causing errors.

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Use a soft cloth to remove dirt on the display. The touchscreen surface has a special protective coating that prevents dirt from sticking to it. Not using a soft cloth may cause damage to the special protective coating on the touchscreen surface. To clean off hard to remove dirt, blow some steam on the particular area and gently wipe away using a soft cloth.

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Do not use excessive force on the display. Avoid placing objects on top of the display as this may cause the glass to break.

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When there is a noticeable discrepancy in the operation of the touchscreen function (your touch does not correlate to the proper location on the screen), recalibrate the touchscreen display.

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1.5 Declaration of Conformity (DoC)

FCC Statement

This device complies with part 15 FCC rules.

Operation is subject to the following two conditions:

- This device may not cause harmful interference.
- This device must accept any interference received including interference that may cause undesired operation

This equipment has been tested and found to comply with the limits for a class "B" digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at him own expense.

CE Notice (European Union)

Electromagnetic Compatibility Directive (2014/30/EU)

- EN55035:2017/A11:2020
 - IEC61000-4-2: 2008
 - IEC61000-4-3: 2020
 - IEC61000-4-4: 2012
 - IEC61000-4-5: 2014/A1:2017
 - IEC61000-4-6: 2013+COR1:2015
 - IEC61000-4-8: 2009
 - IEC61000-4-11: 2020
 - EN55032: 2015+A11:2020
 - EN61000-3-2:2019+A1:2021
 - EN61000-3-3:2013+A2:2021
- Low Voltage Directive (2014/35/EU)
- EN 62368-1:2014+A11:2017

This equipment is in conformity with the requirement of the following EU legislations and harmonized standards. Product also complies with the Council directions.

1.6 Copyright Notice

No part of this document may be reproduced, copied, translated, or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the prior written permission of the original manufacturer.

1.7 Trademark Acknowledgement

Brand and product names are trademarks or registered trademarks of their respective owners.

1.8 Disclaimer

Winmate Inc. reserve the right to make changes, without notice, to any product, including circuits and/or software described or contained in this manual in order to improve design and/or performance. We assume no responsibility or liability for the use of the described product(s) conveys no license or title under any patent, copyright, or masks work rights to these products, and make no representations or warranties that these products are free from patent, copyright, or mask work right infringement, unless otherwise specified. Applications that are described in this manual are for illustration purposes only. We make no representation or guarantee that such application will be suitable for the specified use without further testing or modification.

Chapter 2: Getting Started

This chapter tells you the step-by-step to how to get the tablet computer up and running.

2.1 Introduction

Congratulations on purchasing the Winmate® 10.1-inch Rugged Android Tablet MD10W2.

The MD10W2 is a rugged tablet with robust features designed to withstand medical use while providing high-tech solutions that increase productivity, improve safety, and reduce operational costs. The processing power comes from the ARM Genio 700, Octa-core (up to 2.2 GHz) processor for high performance. The tablet features a brilliant, in-plane switching with anti-glare screen protection film, projective capacitive touch screen, which is outdoor viewable and offers 1920 x 1200 pixel resolution, and weighing at just 1.0 kg, the MD10W2 delivers lightweight mobility in a rugged tablet.


Highlights

- ARM Genio 700, Octa-core, up to 2.2 GHz
- 10.1” TFT LCD Display
- 8.0 MP Camera in Front Side
- 13.0 MP Rear Camera (with LED auxiliary light, Auto Focus)
- 4 GB RAM + 64 GB eMMC
- Android™ 13.0

2.2 Unpacking

Check and identify the standard supplied accessories.

Standard factory shipment list:

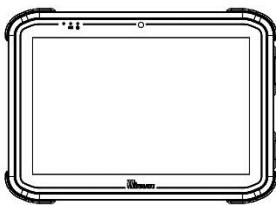Tablet Computer

Standard Battery

User Manual

Power Adapter


Power Cord and Plugs

2.3 Description of Parts

Dimensions: 268 x 187 x 22.2 mm


Front:

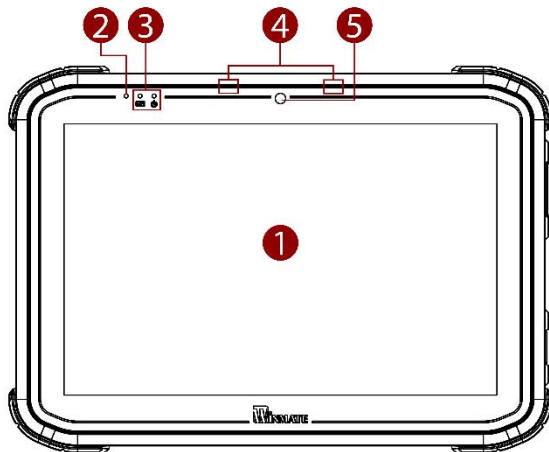Left:

Right:


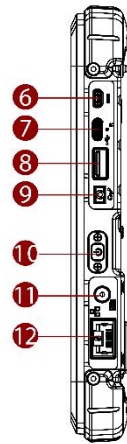
No	Item	Description
1	10.1" Touch Screen	Display and act as one of the inputs for the Tablet Computer
2	Light Sensor	To automatically adjust the screen brightness based on the surrounding light.
3	LED Indicators	Show the current status of the Tablet Computer Orange: Charging Green: Full Charge, Blue: Turning on.
4	Microphone	Enable audio recording for communications, music vocals, speech and sound recording.
5	Front Camera	8 Megapixel Webcam

No	Item	Description
6	Micro HDMI	HDMI output to external display
7	USB Type-C*	USB port connects to USB devices, such as USB flash drive, mouse, keyboard, printer, etc.
8	USB Type-A	USB port connects to USB devices, such as USB flash drive, mouse, keyboard, printer, etc.
9	Audio Jack	Headset or headphone can be connected (Mic is only supported with 3.5mm 4-pin headset)
10	Power Button	Power on/off tablet
11	Power Jack	DC power input jack.
12	RJ45 Connector	Connects the LAN cable.

Front LED Indicator Details:

* **Note:** In accessibility setting, enable the "Connect to PC" to allow USB type-C host function and disabled USB type-A connection. In newer systems might support this feature automatically without additional enable the "Connect to PC" steps.

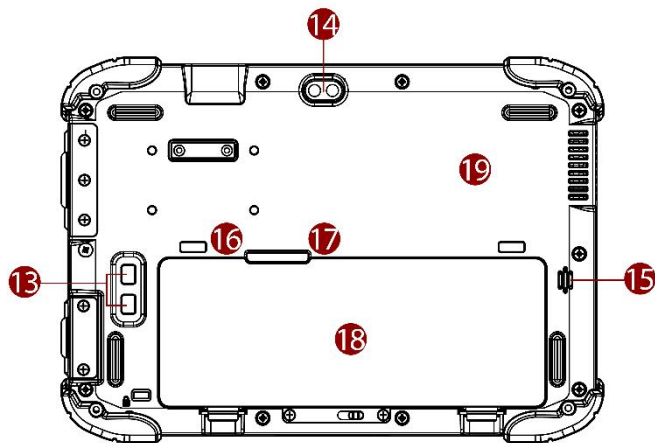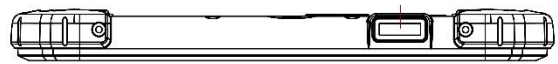Rear:**Top:****Bottom:**

No	Item	Description
13	Function Buttons	Programmable function button. Default: Menu
14	Rear Camera	13 MP Autofocus Camera. Captures photos and videos.
15	Speaker	Send out sound and voice from the computer.
16	SIM card Slot	Service window to insert SIM card.
17	MicroSD Card Slot	Memory card reader port for SD Card
18	Battery	11.4 V, typ. 3500 mAh. Li-Polymer Battery (3S1P)

2.4 First-Time Operation

To get started:

1. Install the battery.
2. Install memory card.
3. Charge your tablet.
4. Start the tablet.

2.4.1 Inserting Battery Pack

Warning! Discontinue using device immediately if you notice the battery is swelling, or looks puffy or misshapen. Do not attempt to charge or use the device. Do not handle a damaged or leaking battery. Do not let leaking battery fluid come into contact with your eyes, skin or clothing.

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
Caution To keep battery health and prevent possible issues from battery over-discharged, please charge the tablet built-in back-up battery once every two months.

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

To install the standard battery pack:

1. Disconnect any cords connected to the tablet computer and turn the tablet over, back facing up.
2. Locate the Battery Lock and slide it to the "Unlocked" position.
3. Insert the battery with the connector-pins side first, and then lower the battery into the battery compartment. The battery will be flat against the housing once fully inserted.
4. Slide the Battery Lock to the "Locked" position.

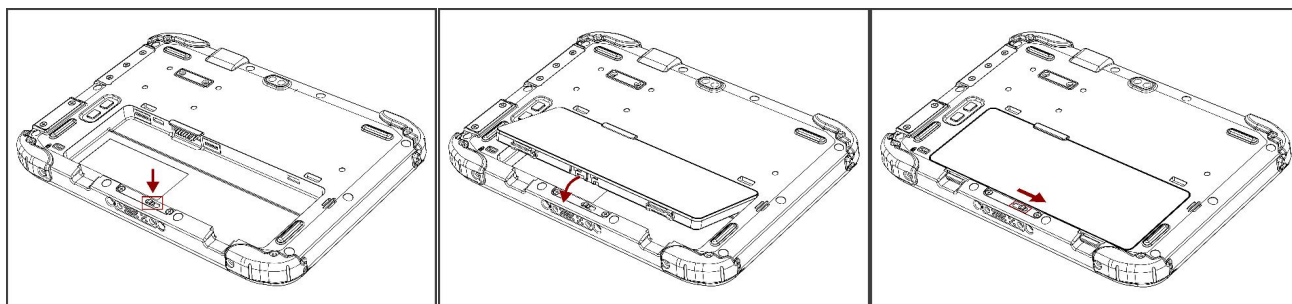To remove the battery pack:



1. Shut down the tablet (or put it into sleep mode). Disconnect any cords connected to the computer.
2. Turn the tablet computer over, back facing up.
3. Locate the Battery Lock and slide it to the "Unlocked" position.
4. Push and hold the **Battery Release Latch** while lifting the battery out of the compartment.

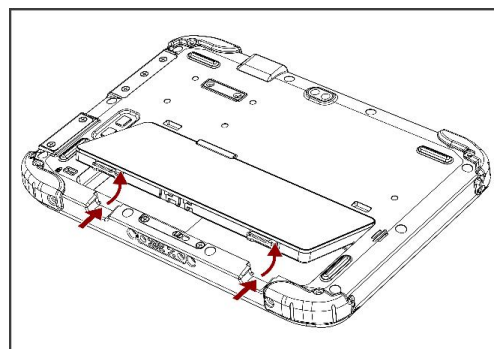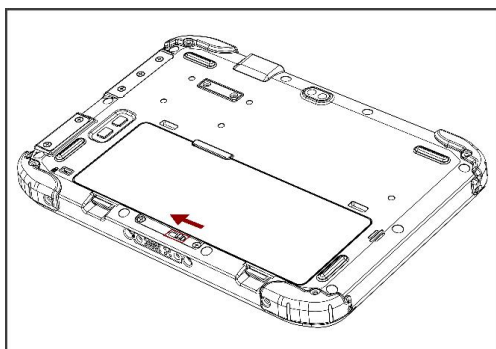2.4.2 Inserting MicroSD Card

To install MicroSD Card:

1. Remove the battery, as described in the section 2.4.1 *Inserting Battery Pack*.
2. Insert microSD card into slot with contacts facing down.
3. Press the card in to ensure that it is seated properly.
4. Put back the battery, as described in the section 2.4.1 *Inserting Battery Pack*.

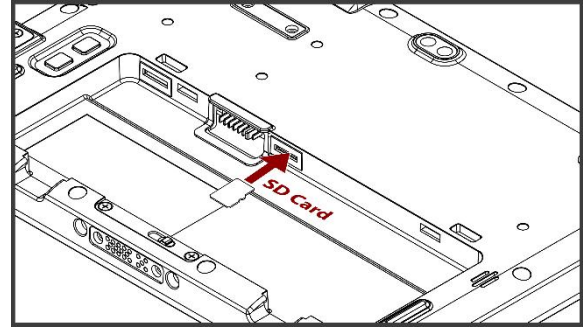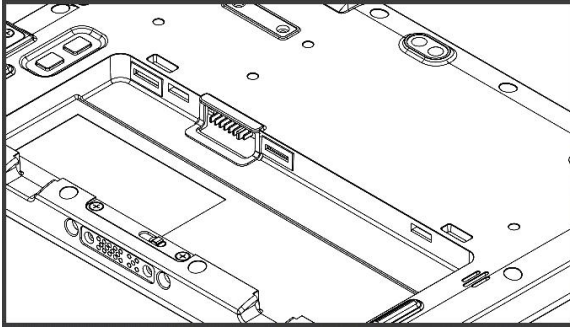To remove Micro SD Card:

1. Remove the battery, as described in the section 2.4.1 *Inserting Battery Pack*.
2. Push and release the microSD card.
3. Pull out the microSD card.
4. Put back the battery, as described in the section 2.4.1 *Inserting Battery Pack*.

2.4.3 Connecting to Power Supply

The Tablet Computer operates either on external AC power or internal battery power.

It is recommended to use AC power when you start up the Tablet Computer for the first time.

Caution Use only the AC adapter included with your Tablet Computer. Using other AC adapters may damage the Tablet Computer and the Battery.

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Caution For users start up the Tablet Computer for the first time with internal battery power, please make sure the battery is fully charged before using.

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Perform the following to connect the Tablet computer to AC Power:

1. Open I/O protective cover by pulling out the latch.
2. Plug the DC cord of the AC adapter into the power jack of the tablet computer.
3. Plug the AC power cord into the AC adapter.
4. Plug the AC power cord into the electrical outlet.

When the AC Adapter is connected, the battery pack starts charging. The battery LED indicator on the tablet computer would show different states. Please refer to description of LED indicator on Description of Parts.

Standard Battery	
Battery Operating Time	10 Hours
Battery Charging Time	2.5 Hours

2.5 Turning On and Off the Device

This section describes how to turn on/ off the device.

2.5.1 Turning On


1. Press and hold the power button for **3 seconds**, located on the right side of the unit.
2. When the device vibrates, release the power button. The desktop screen is displayed.

Note: For power saving purpose, please enable Connect to PC in Accessibility settings. In newer systems might support this feature automatically without additional enable the Connect to PC steps.

2.5.2 Turning Off

1. To turn off the device, press and hold the power button until the phone options menu appears, tap the **Power Off** .

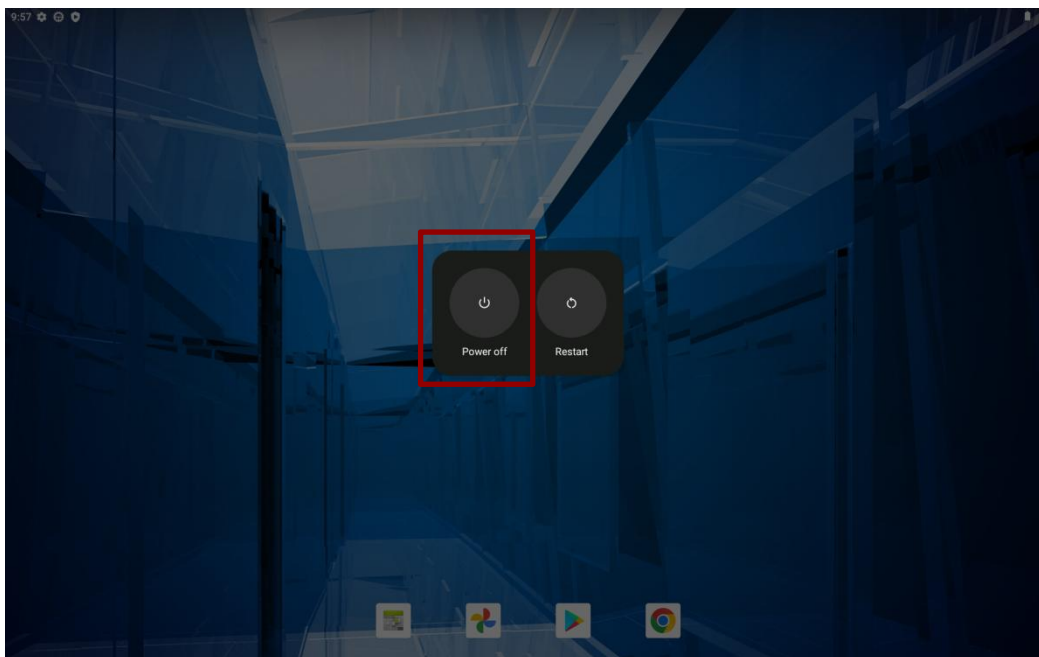
2. User can also use the Phone options menu to put the device into silent or airplane mode, see features table.


2.6 Resetting to Factory Defaults

To reset the device to factory defaults:


1. Slide from the bottom side, then open **Settings**  app.
2. Tap **System** , user can select to back up application data, WLAN passwords and other settings to **Google** servers by tapping **Backup My Data**.

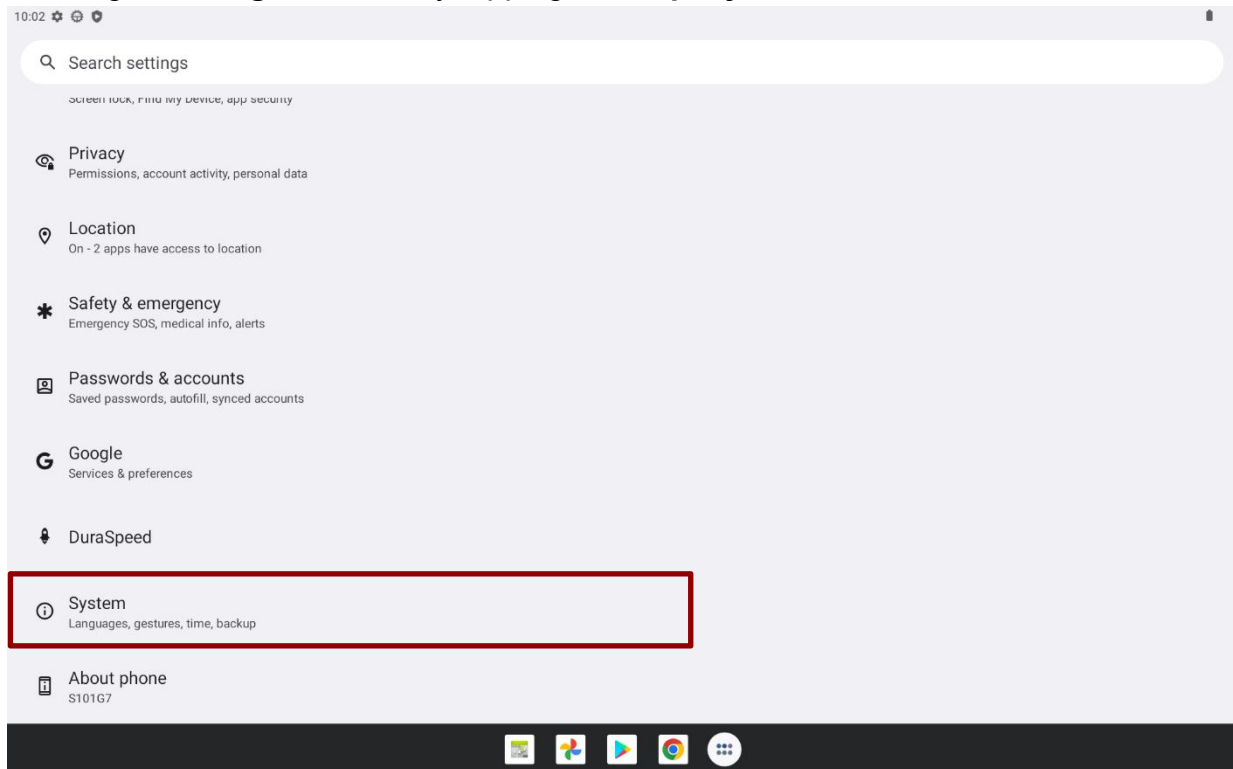3. Tap **Reset options**.

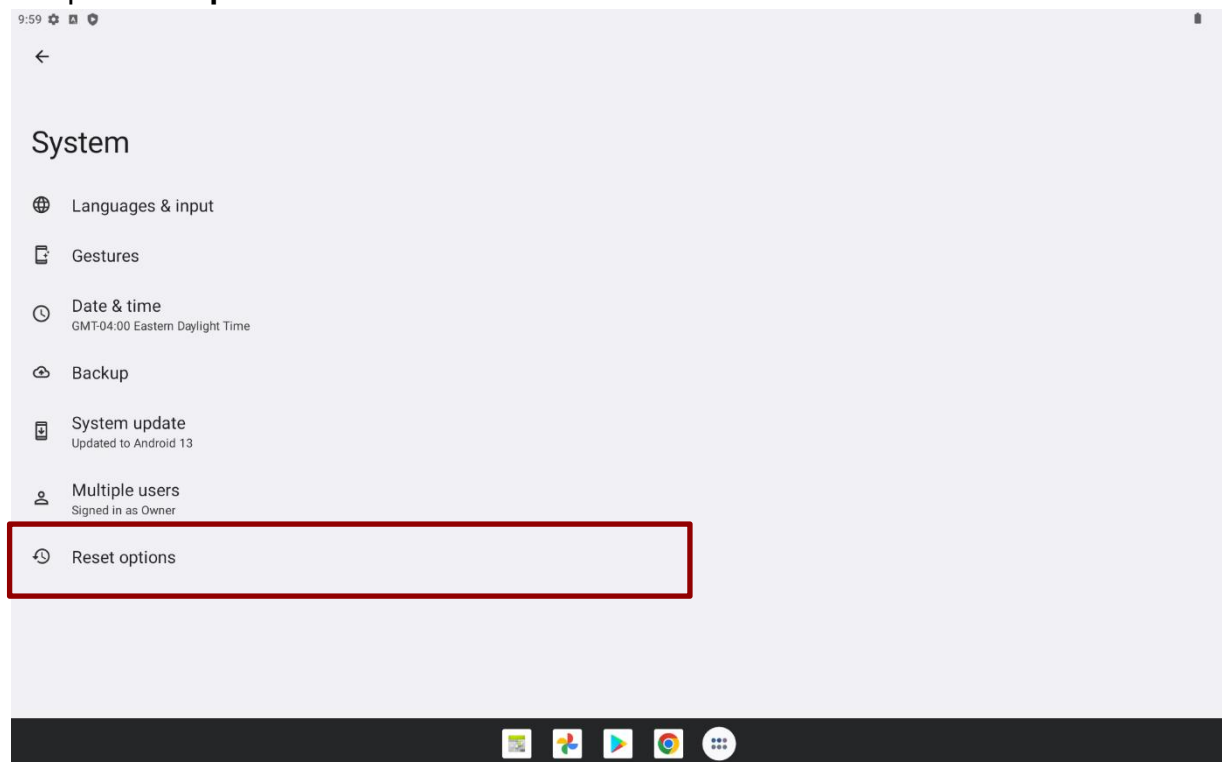4. Tap **Erase all data (factory reset)** to reset the tablet.

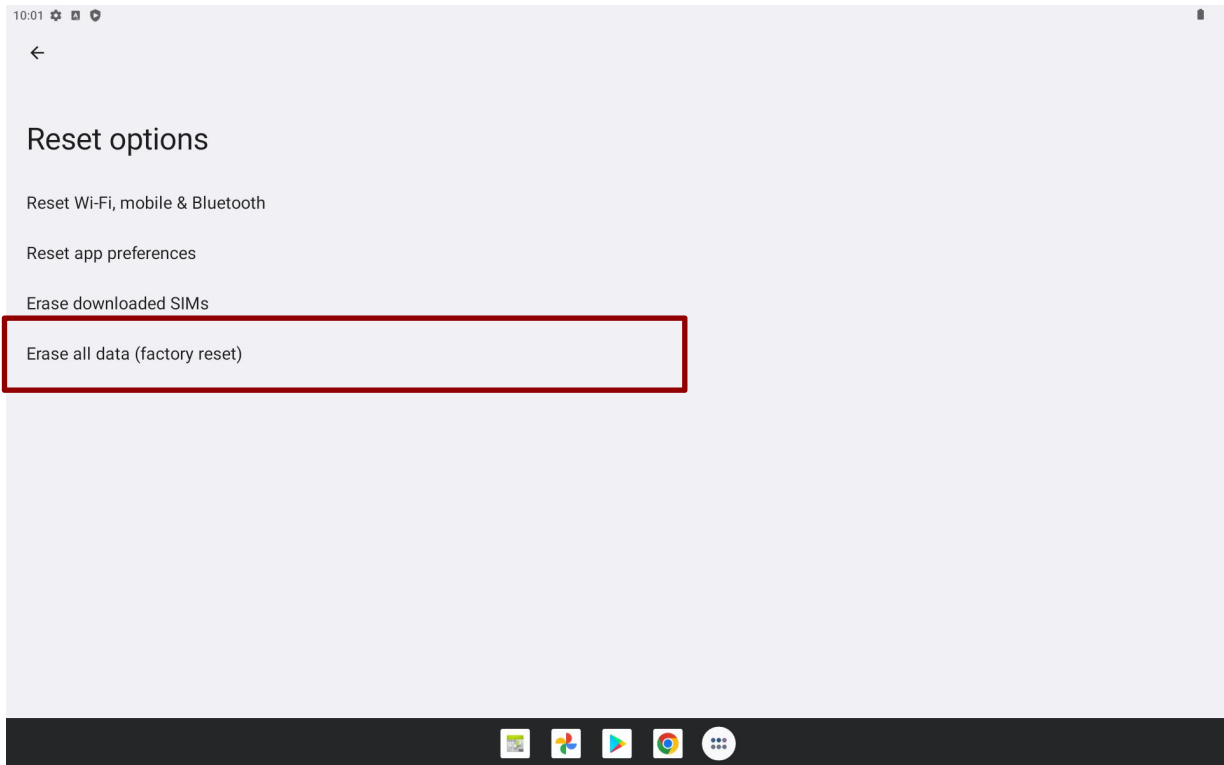5. Read all the available information and then press **Erase all data** to delete all the data in storage such as music or photo.

Erase all data (factory reset)

This will erase all data from your phone's **internal storage**, including:

- Your Google Account
- System and app data and settings
- Downloaded apps
- Music
- Photos
- Other user data

6. Press **Erase all data** to remove all your personal information and downloaded apps.

Erase all data?

All of your personal information and downloaded apps will be deleted.
You can't undo this action.

ERASE ALL DATA

Chapter 3: Using the Tablet

This chapter explains how to use the buttons, status icons, and controls of the MD10W2 Rugged Tablet running on Android 13.0 operating system and provides instructions on entering and capturing data.

3.1 Using the Touch Screen


The touchscreen supports the core gestures for Android.

Gesture	Action	Description
Tap		To select or activate something
Tap & Hold		Touch and hold an item on the screen by touching it and not lifting your finger until it responds
Touch, Hold & Drag		Touch and hold an item for a moment and then without lifting your finger, move your finger on the screen until you reach the target position. i.e. moving the apps around the home screen
Swipe or slide		Quickly move the finger across the surface of the screen. i.e. slide a lock screen to open the lock
Double-tap		Tap quickly twice on a webpage, screen, or App to zoom
Scale		Place two fingers on the screen and pinching them together

3.2 Waking Up the Tablet


If you do not use the tablet for a while, the screen dims and then darkens to conserve the Battery. To wake up the tablet:

1. Press the Power button on the right side of the tablet.
2. Drag the screen upwards to unlock.

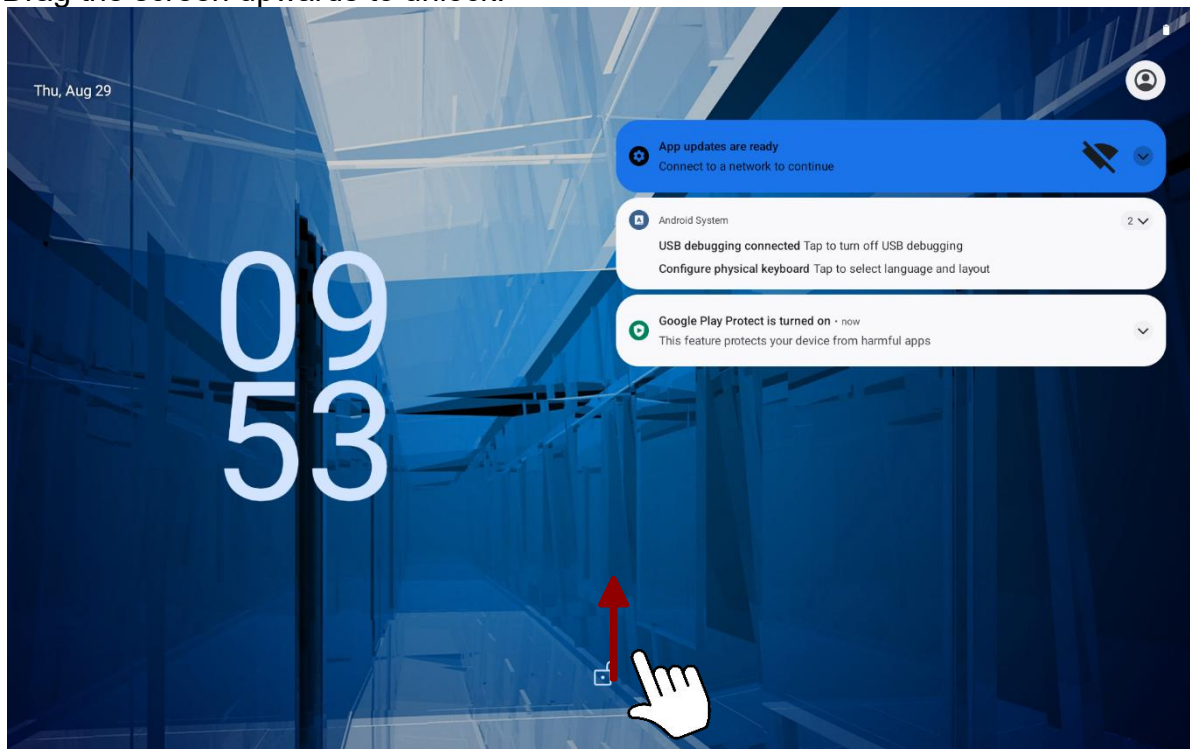3.3 Using the Home Screen

Home screen is the starting point to access all the features on your tablet Device. It displays application icons, widgets, shortcuts, and other features. User can personalize the home screen. See “*Personalizing home screen*”.

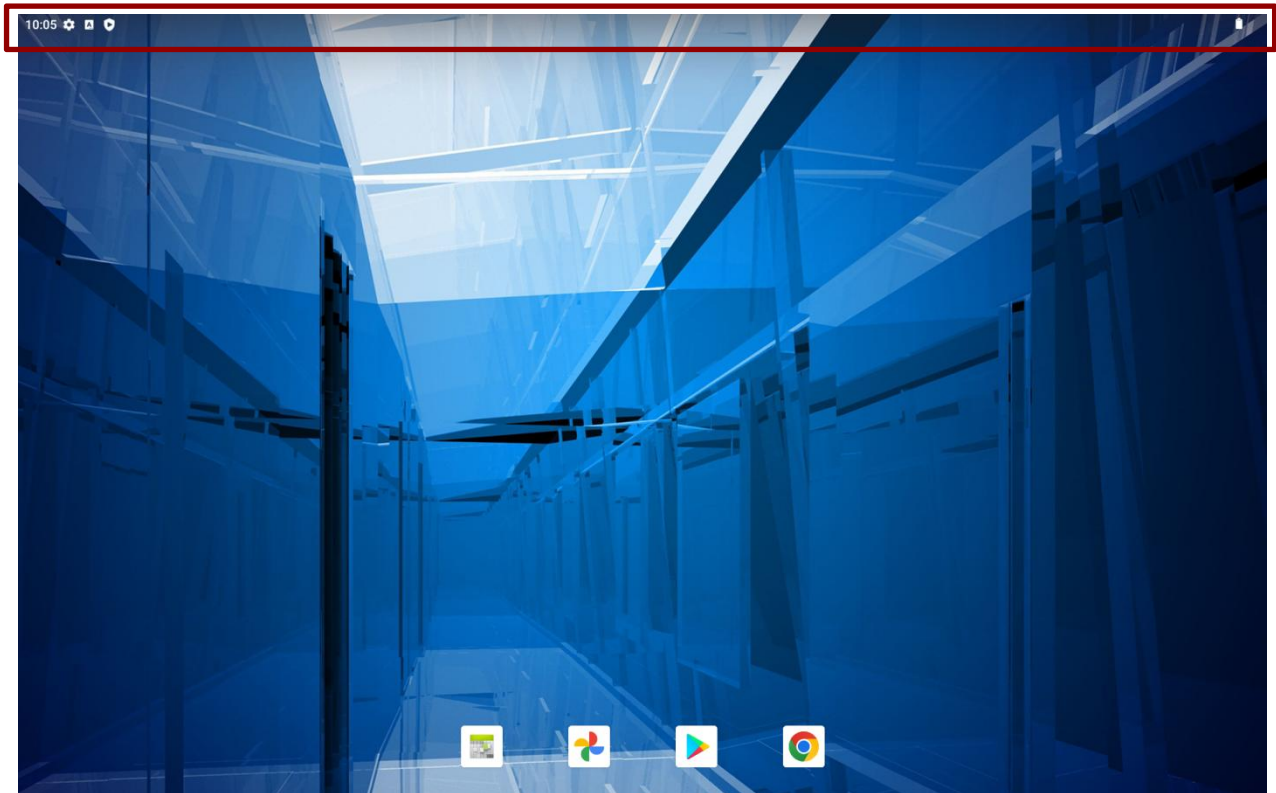
No.	Description
1	Status Bar – Shows the time, signal strength, battery status, and other notifications. It also displays notification icons.
2	Application – Shows applications that user can open directly on the home screen.

3.3.1 Personalizing Home Screen


User can add application icons, shortcuts, widgets, and other items to any part of any Home screen where there is free space. User can also switch the wallpaper.

To change wallpaper on the HOME SCREEN:

1. Swipe from the bottom of the main desktop, tap **Settings** ⚙️, and then tap **wallpaper** 🖼️.

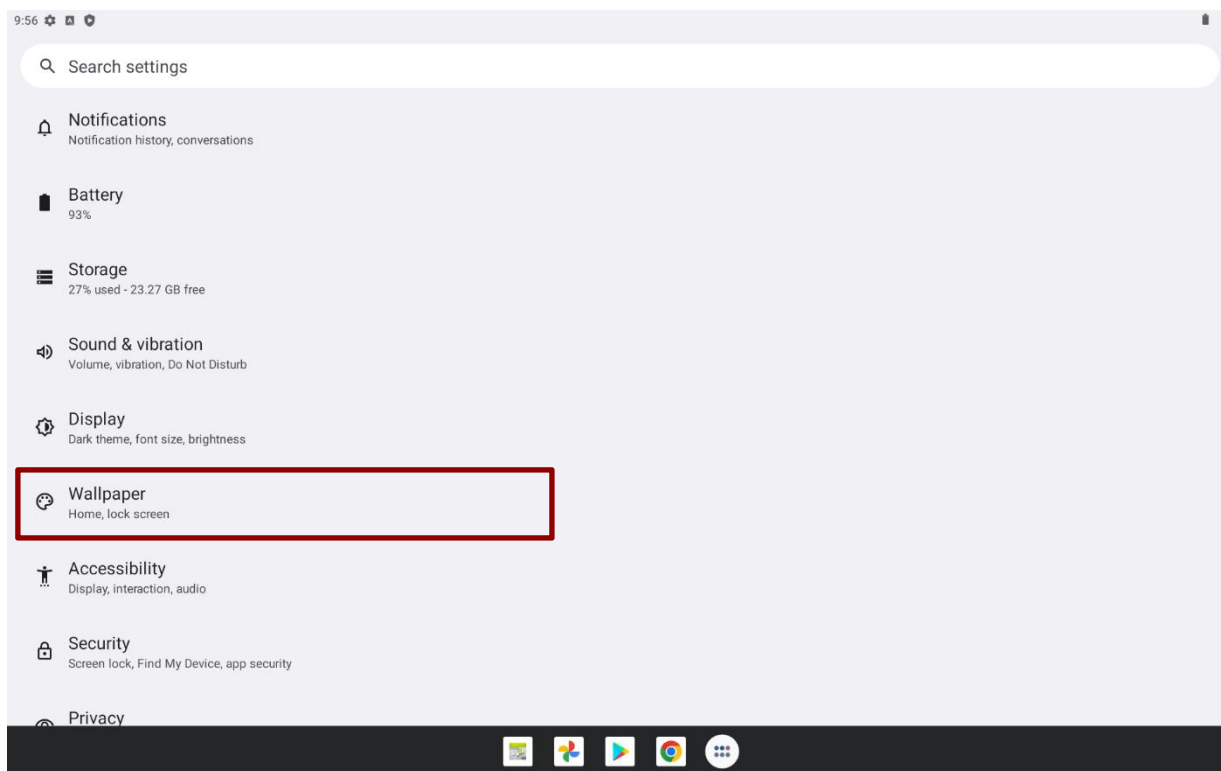
2. In the Wallpaper menu select the location of wallpaper, and then tap the type of item you want to add.

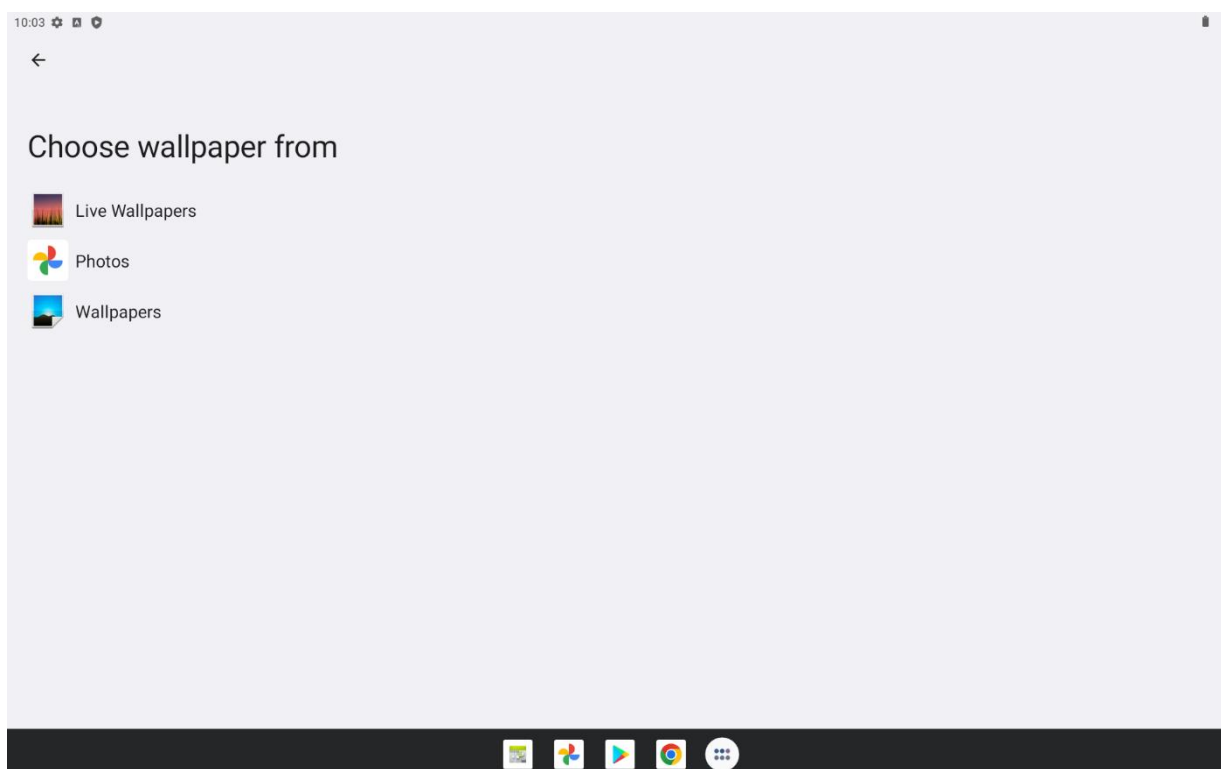
To move an item on the HOME SCREEN:

1. Tap and hold the item you want to move, until it can be moved. Drag it to a new location on the screen.

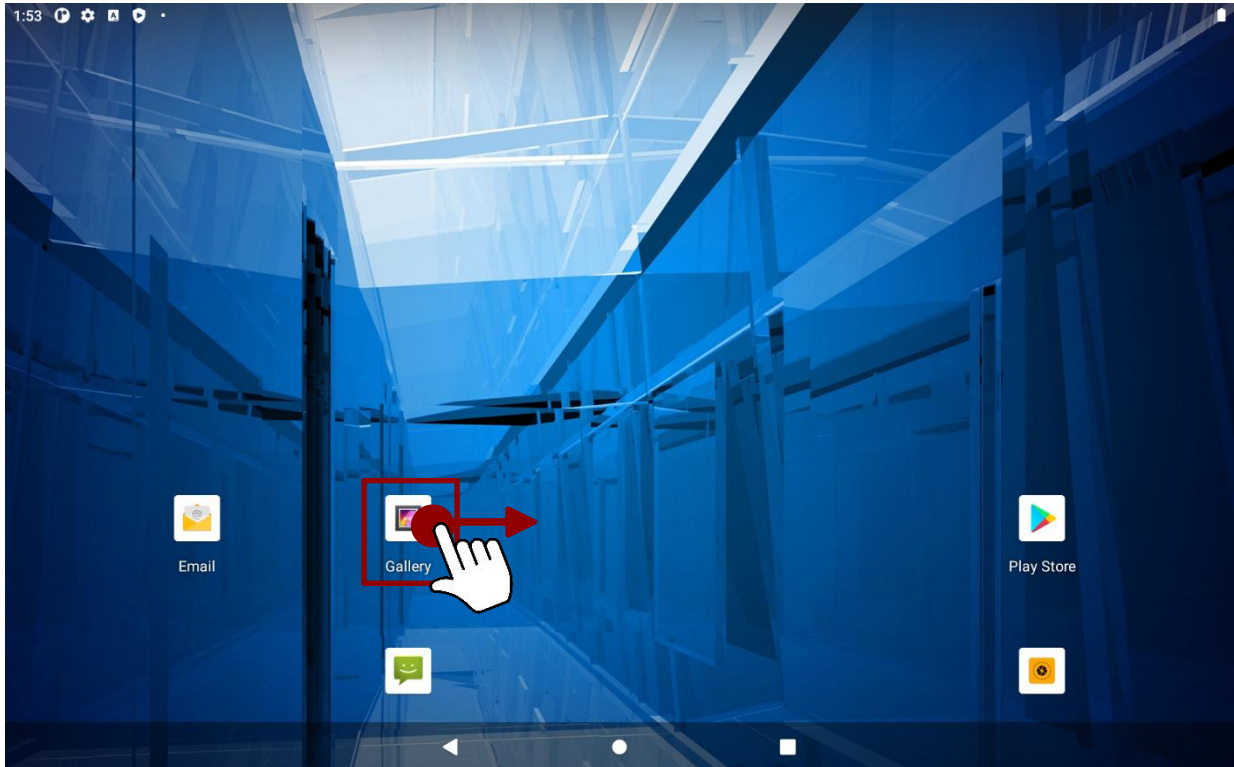
2. Place the item where you want to put it and lift your finger.


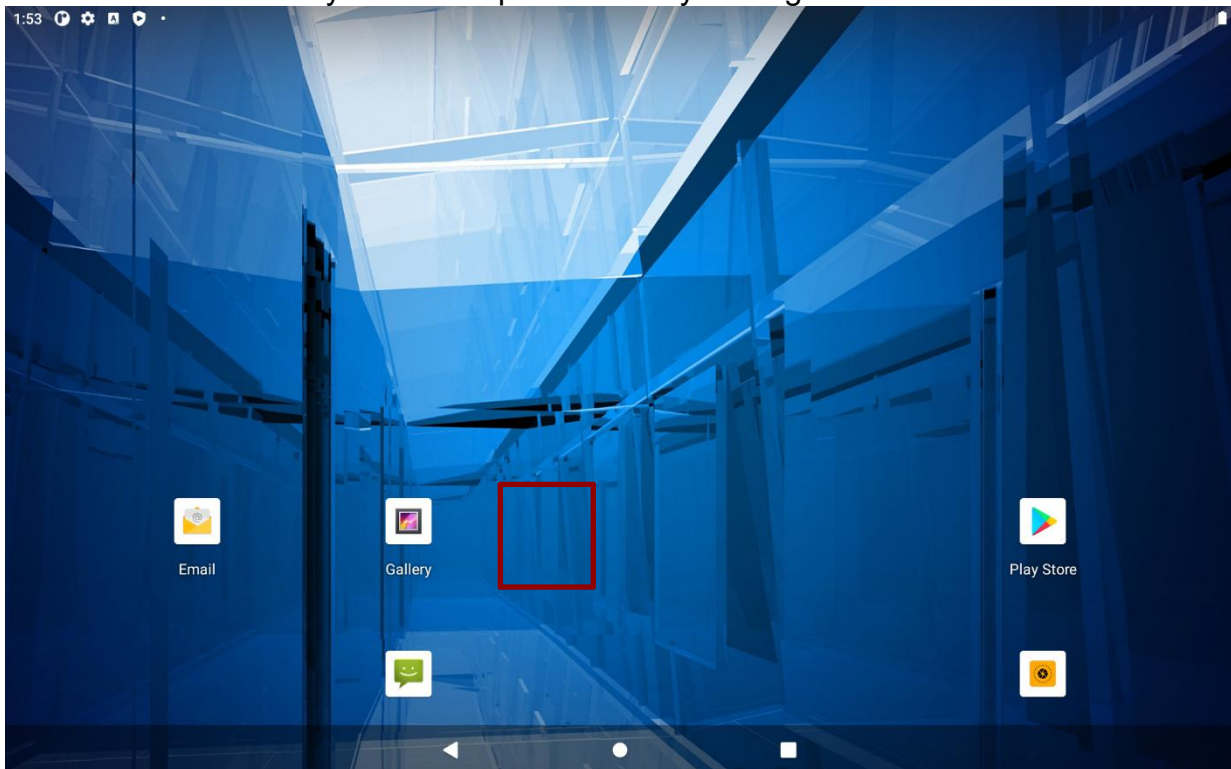
To remove items from the HOME SCREEN:


1. Tap and hold the item you want to remove. The launcher icon will show **Remove** on the top of the screen.
2. Drag the item you want to remove to the Remove icon.
3. When the icon turns red, release your finger.

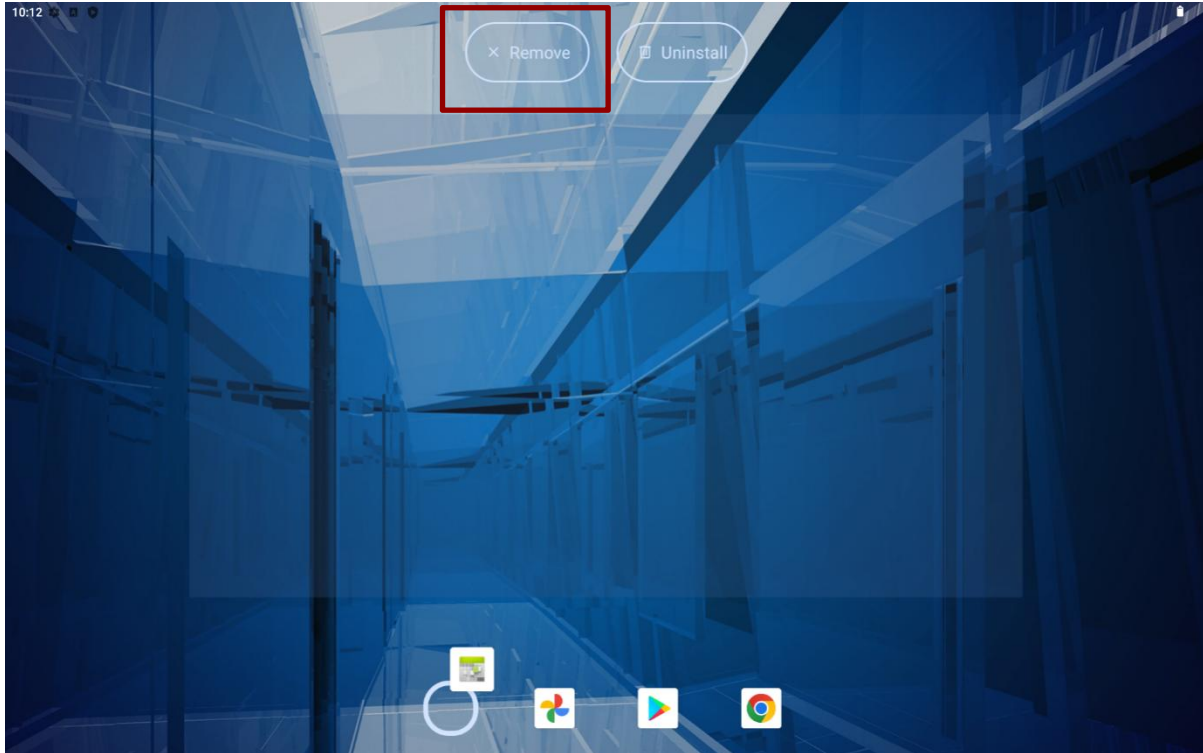
3.3.2 Viewing other parts of the Home Screen

1. Slide your finger left or right across the Home screen. Extensions to the Home screen provide extra space for Apps, shortcuts, and other items.

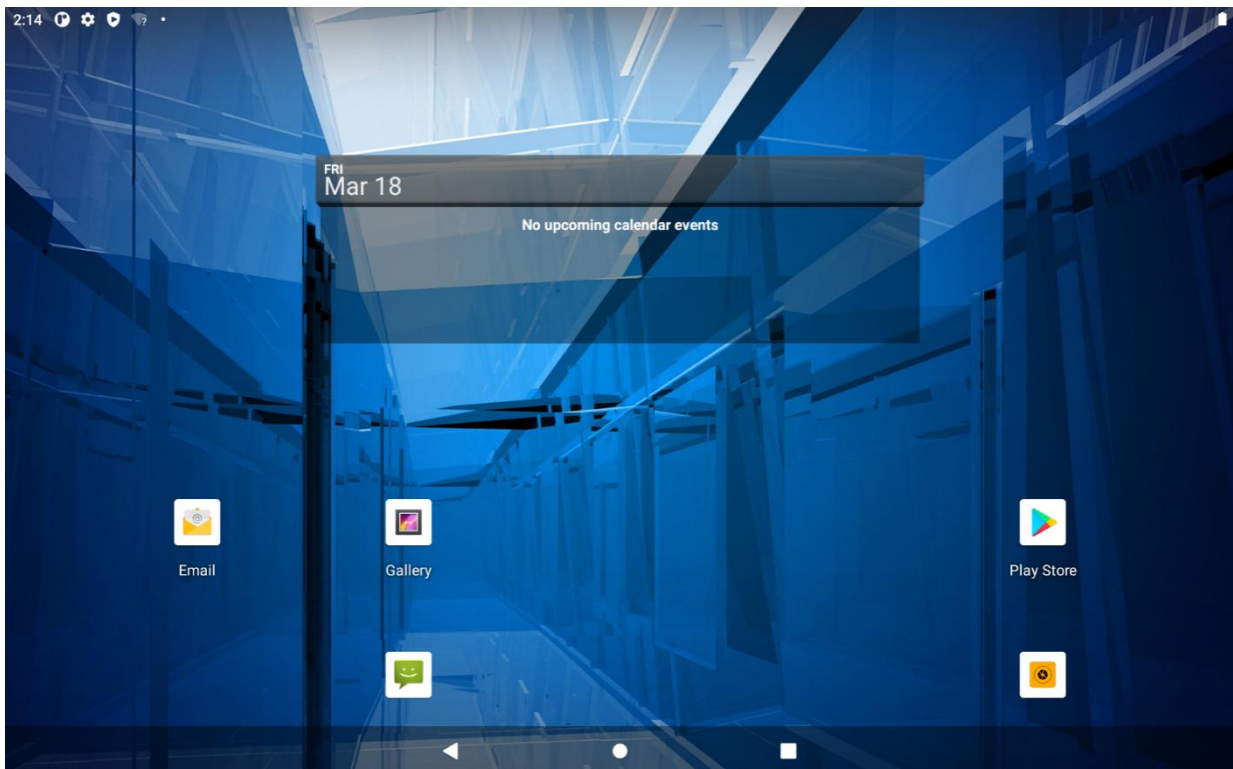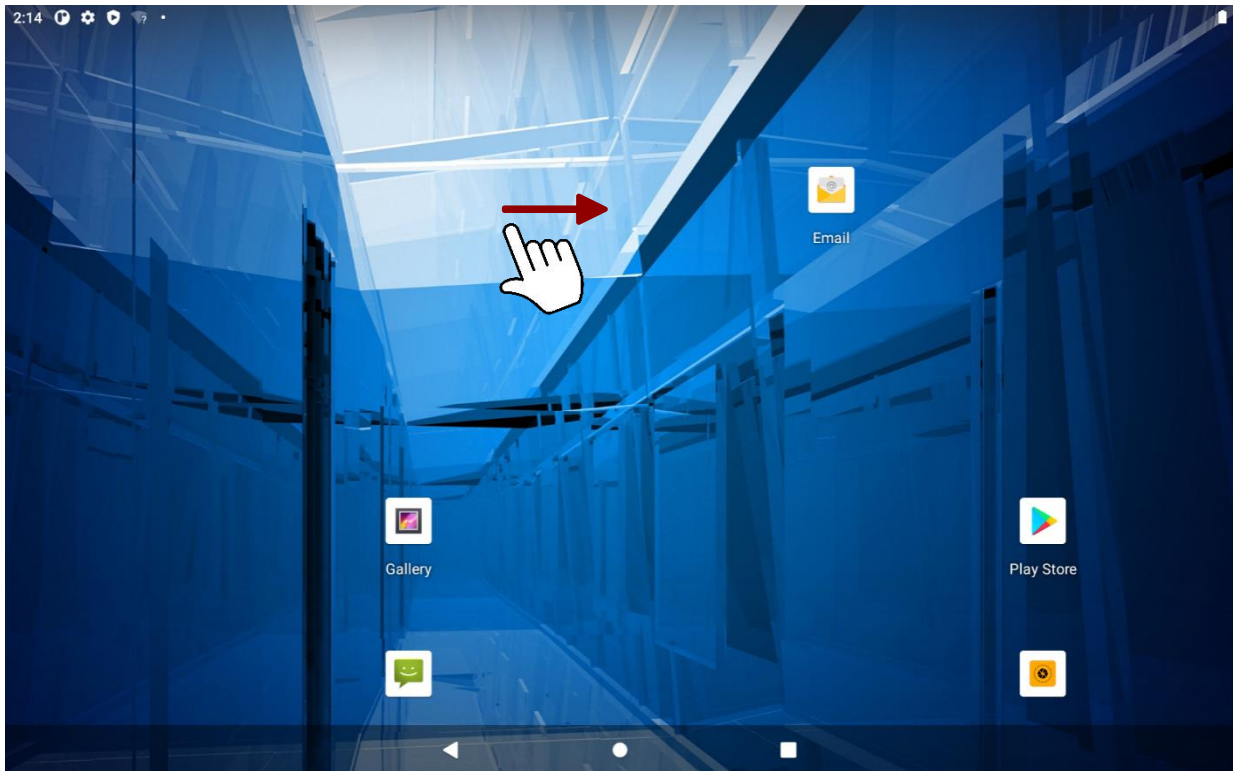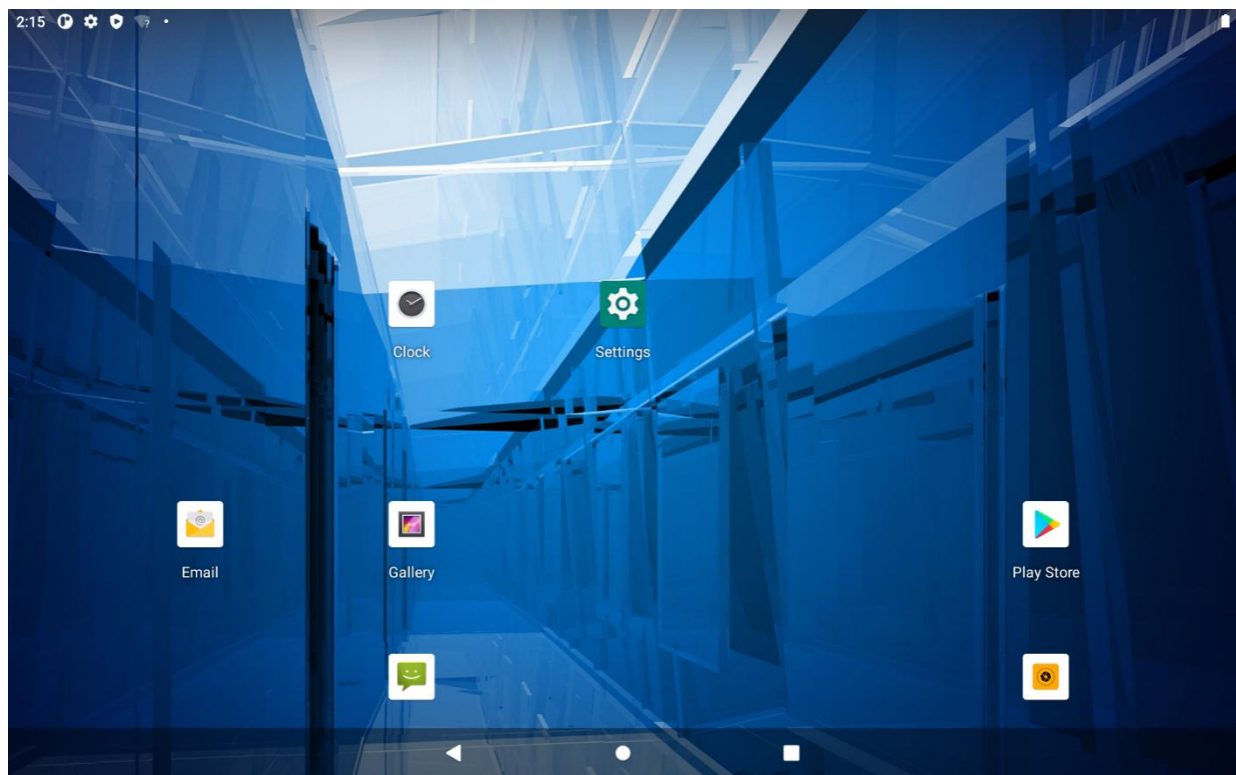
2. Swipe from the bottom of the main desktop to view the entire APPS on the device.

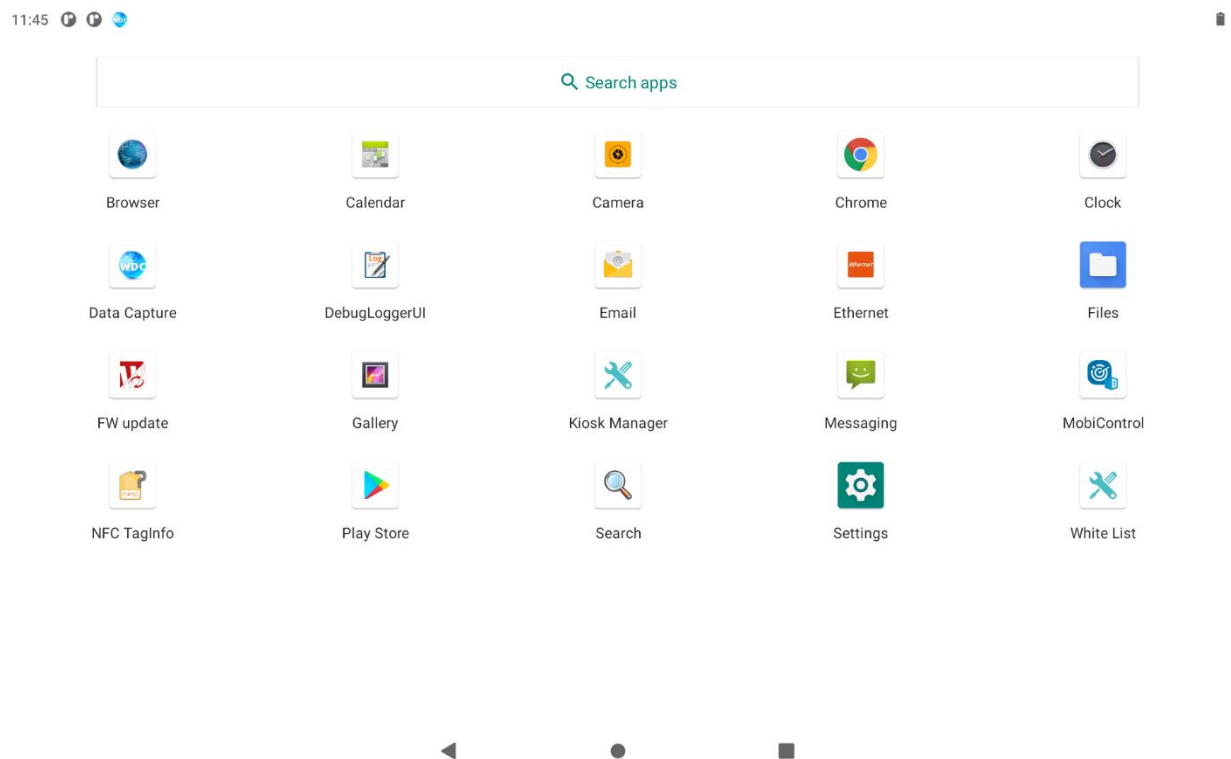3.4 Status and Notification Icons

The status bar appears at the top of every screen. On the left are notification icons along with the current time, and on the right are status icons,.

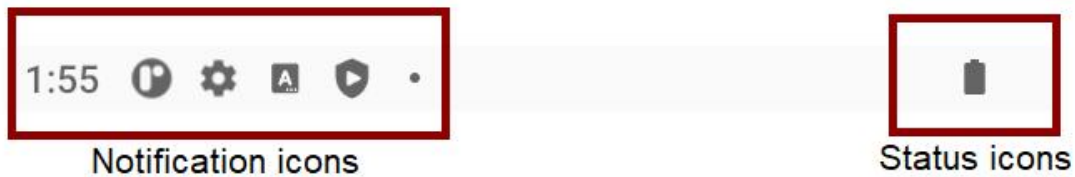3.4.1 Notification Icons

Pull down from the left top of the screen to get the notifications area. The following icons indicate that you have received notification.

	New email message		Incoming call
	New text or multimedia message		Missed call
	Alarm is set		Speakerphone is activated
	Upcoming event		Phone microphone is muted
	Data is syncing		Downloading data
	Problem with sign-in or sync		Phone is connected via USB cable
	A screenshot is captured		

3.4.2 Status Icons

Pull down from the right top of the screen to get the notifications area. The following icons indicate the status of your device, for example if users are connected to a WLAN network that is not connected to the internet, the network icons are appears in white color.

Status icons that appear in the status bar are as follows:


	Connected to GPRS mobile network		Bluetooth is on
	Connected to 1xRTT mobile network		Connected to a Bluetooth device
	Connected to EDGE mobile network		No SIM card installed
	Connected to 3G mobile network		Vibrate mode
	Connected to 4G mobile network		Ringer is silenced
	Mobile network signal strength		Battery is very low
	Searching for a network signal		Battery is low
	No signal		Battery is partially drained
	Connected to a Wi-Fi network		Battery is full
	Roaming		Battery is charging
	Airplane mode		Receiving location data from GPS

3.5 Using the On-Screen Keyboard

User will be able to enter text using the On-Screen Keyboard. Some applications will launch the keyboard automatically. In others, user needs to touch a text field where user wants to enter text to launch the keyboard. User will be able to enter text using the On-Screen Keyboard. Some applications will launch the keyboard automatically. In others, user needs to touch a text field where user wants to enter text to launch the keyboard.

3.5.1 Entering Text

1. Open a text field, in this example we used Messaging  to compose a message. The On-Screen Keyboard opens.

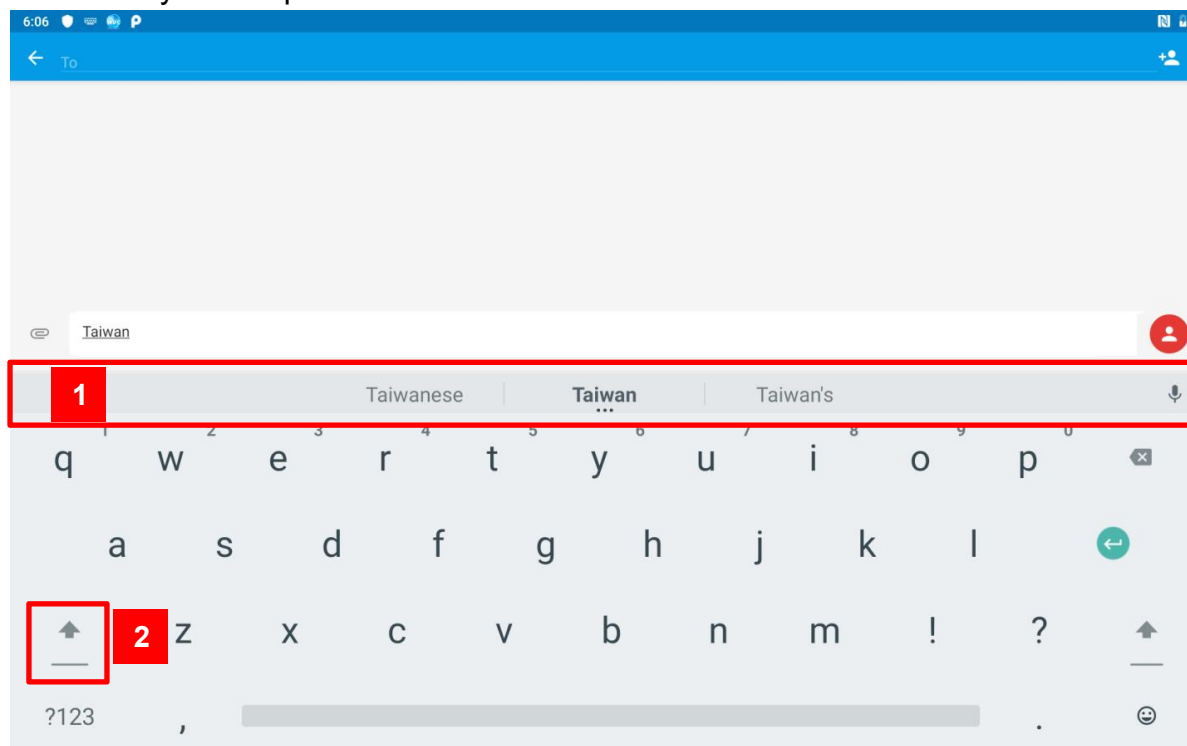
No.	Description
1	Suggestions – Tap a suggestion to enter it in the place of the underlined.
2	Capitalize – Tap once to capitalize the alphabet.


**If user taps in a text field which already contains text, the insertion point is set where user tap and the cursor appears. User can drag the cursor to pick exactly where you want to add more text. User can also select a word to copy, cut and paste it. See 3.5.2 Editing Text.*

2. Tap the delete key  to erase characters to the left of the cursor.


3. Tap **?123** to enter numbers. Tap **ABC** to return to letters.

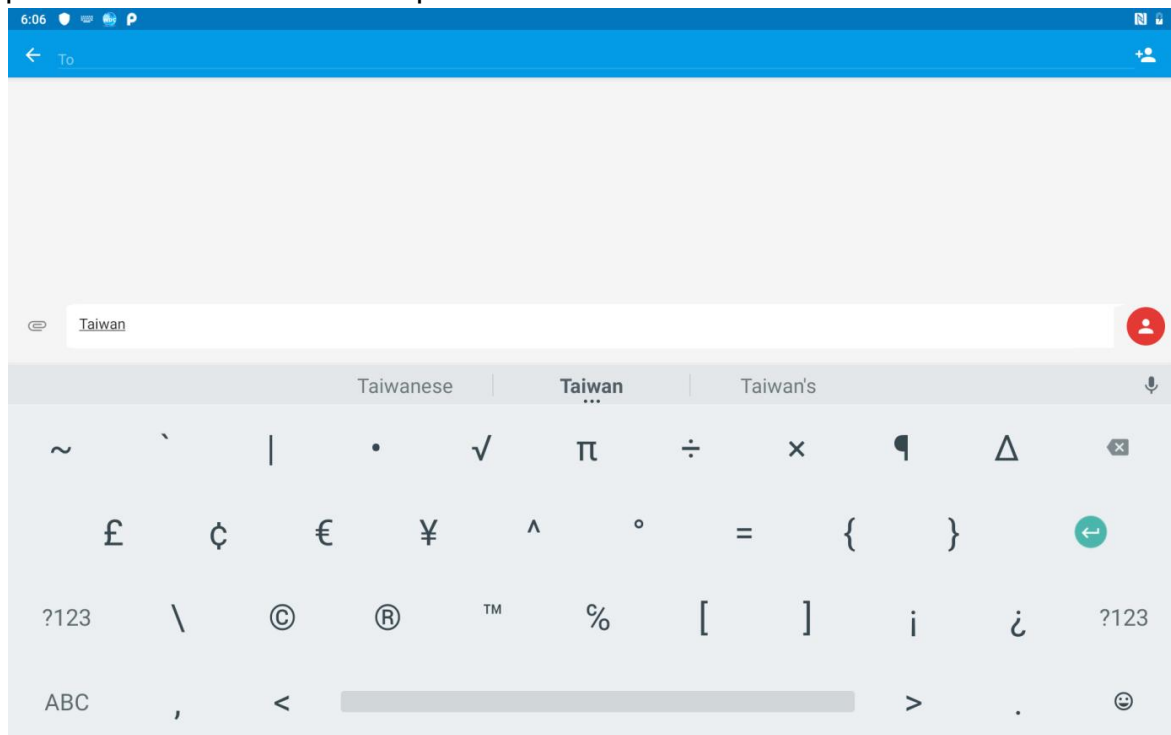4. To enter symbols and other characters, tap and hold one of the top-row keys and slide your finger up to the list of the options. Lift your finger when you find your selection.

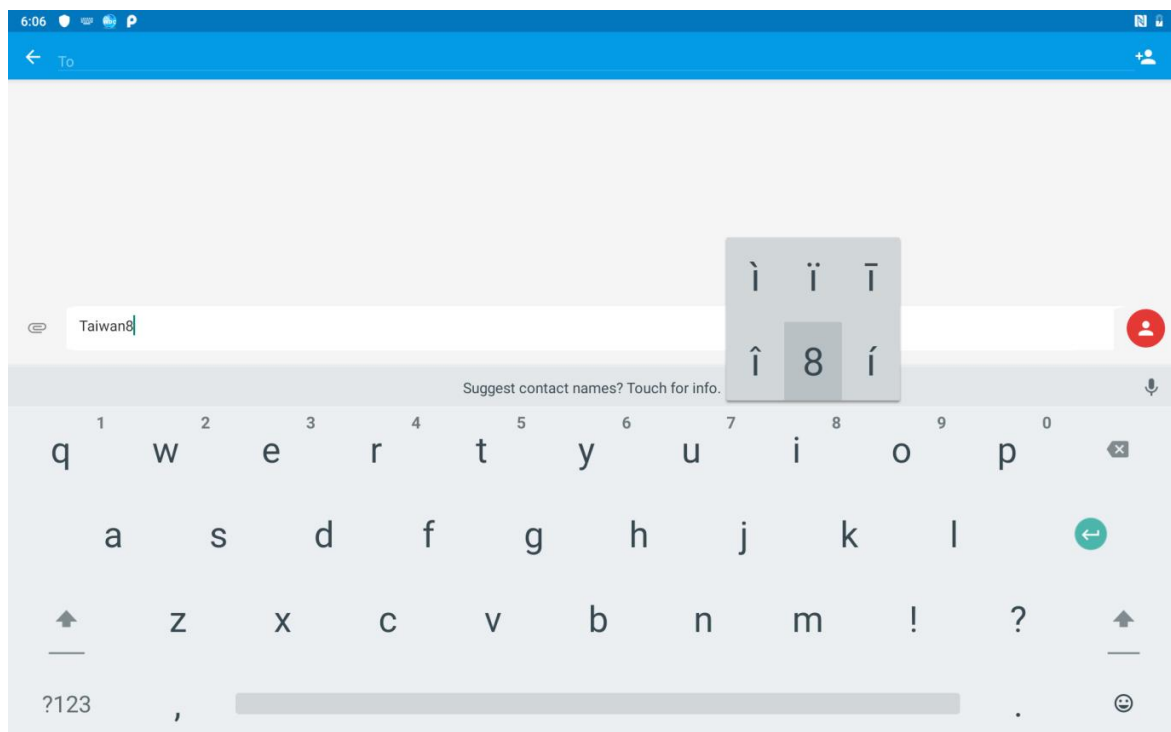5. To view additional symbols, tap  on the numeric keyboard.


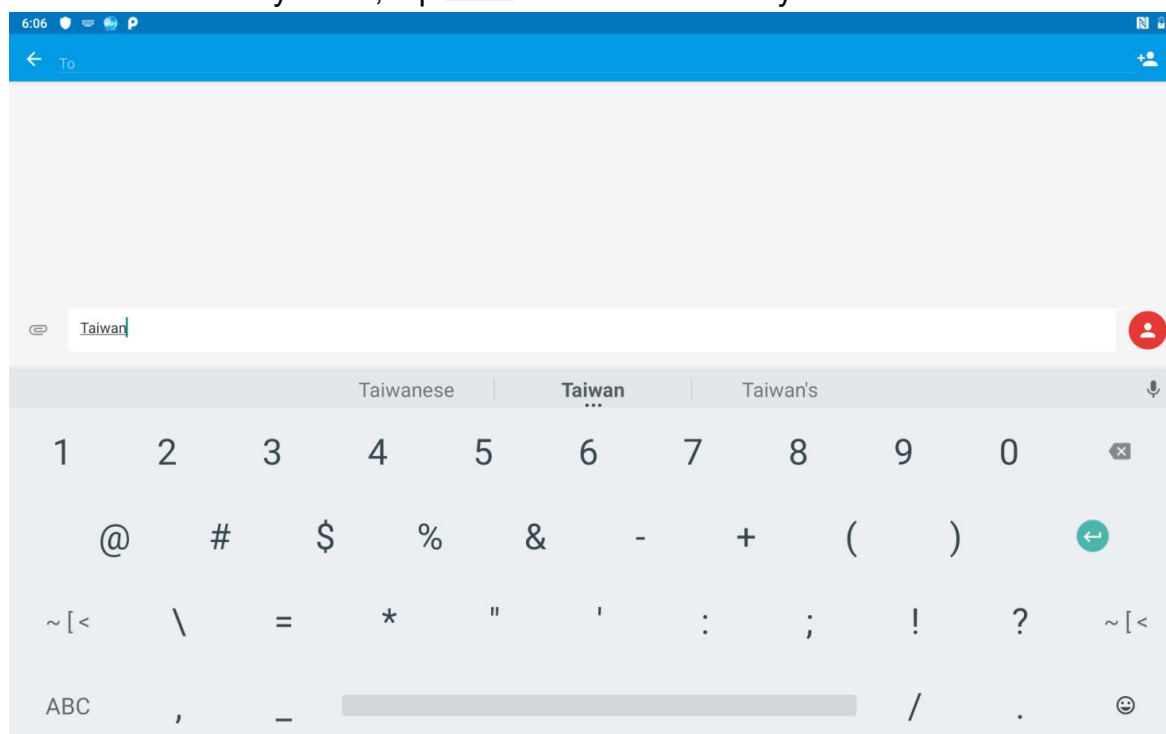
6. To close the On-Screen keyboard, press the back button  on the top or **Fn** button.

3.5.2 Editing Text

To select text:

1. Double-tap the text you want to select.
2. Tap  to select all the sentences in the menu that appears, the selected text is highlighted in green, with a tab at each end of the selected word.
3. Drag selection tab to expand or reduce the range of the selected text.

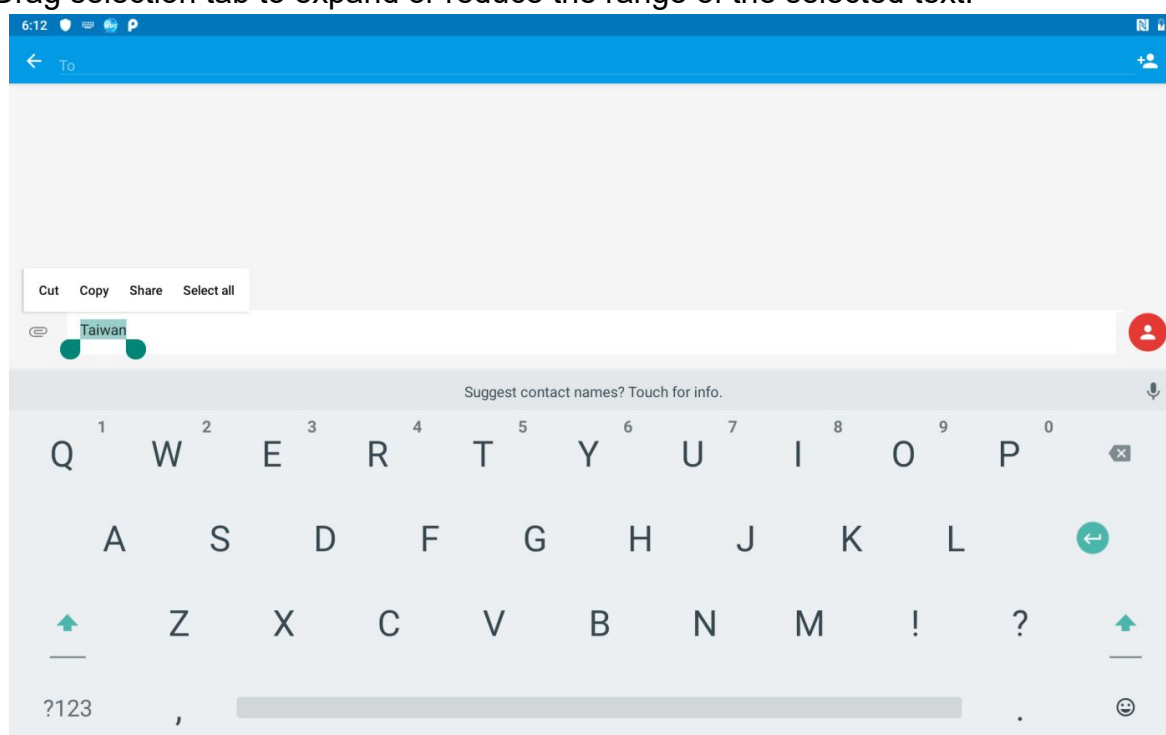
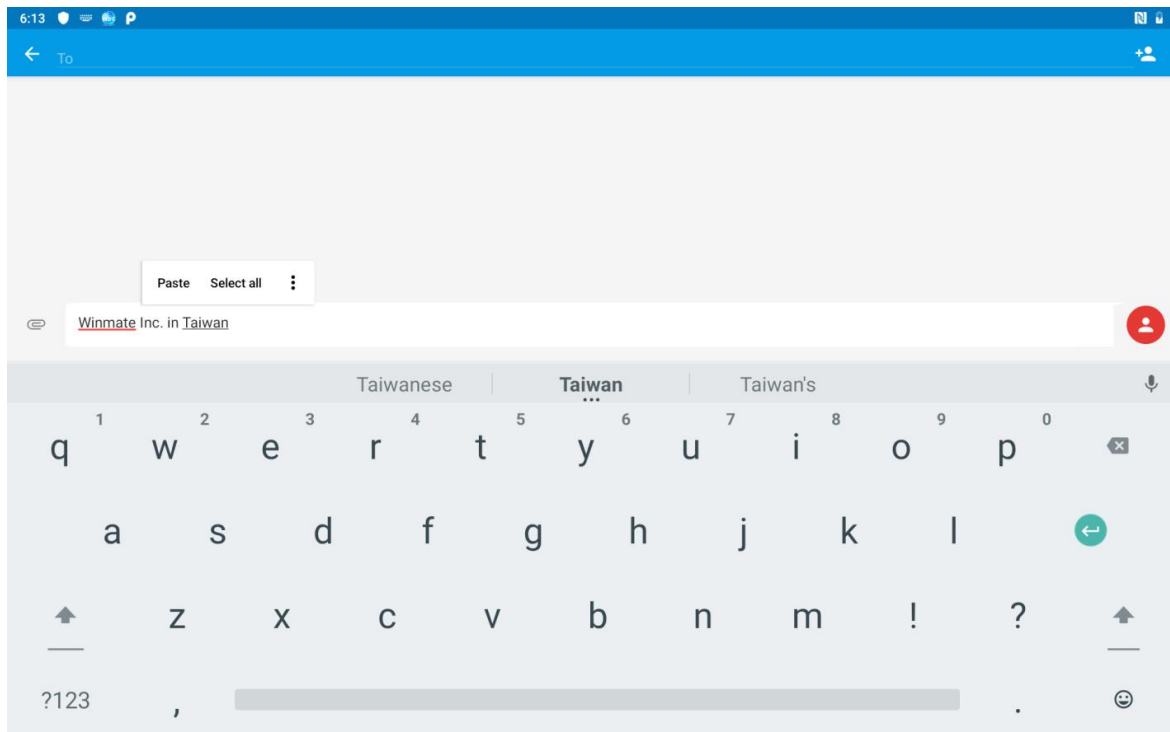
To Cut or Copy Text:

1. Select the text to cut or copy, see the section above [“To select text”](#).
2. Tap **Cut** or **Copy** in the Edit text menu.
3. If you select **Cut**, the selected text is removed from the text field; it is stored in temporary area so that you can paste it into different text field.

To Paste Text:

1. After Cut or Copy the text, see the section above [“To Cut or Copy Text”](#).
2. Place the cursor in the text field where you want to paste the text, user can paste text that you copied from one application into a text field in any application.
3. Tap and hold into the location that you want to paste. Tap **Paste** in the menu that opens to replace the existing text with copied text.

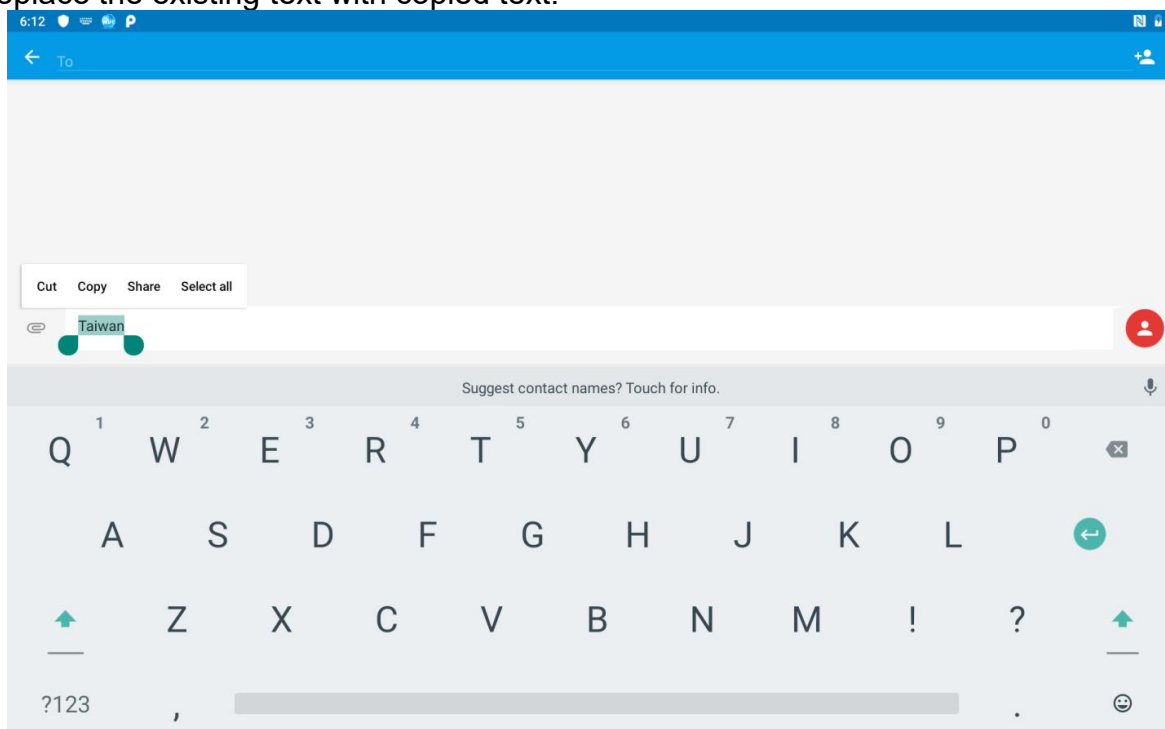
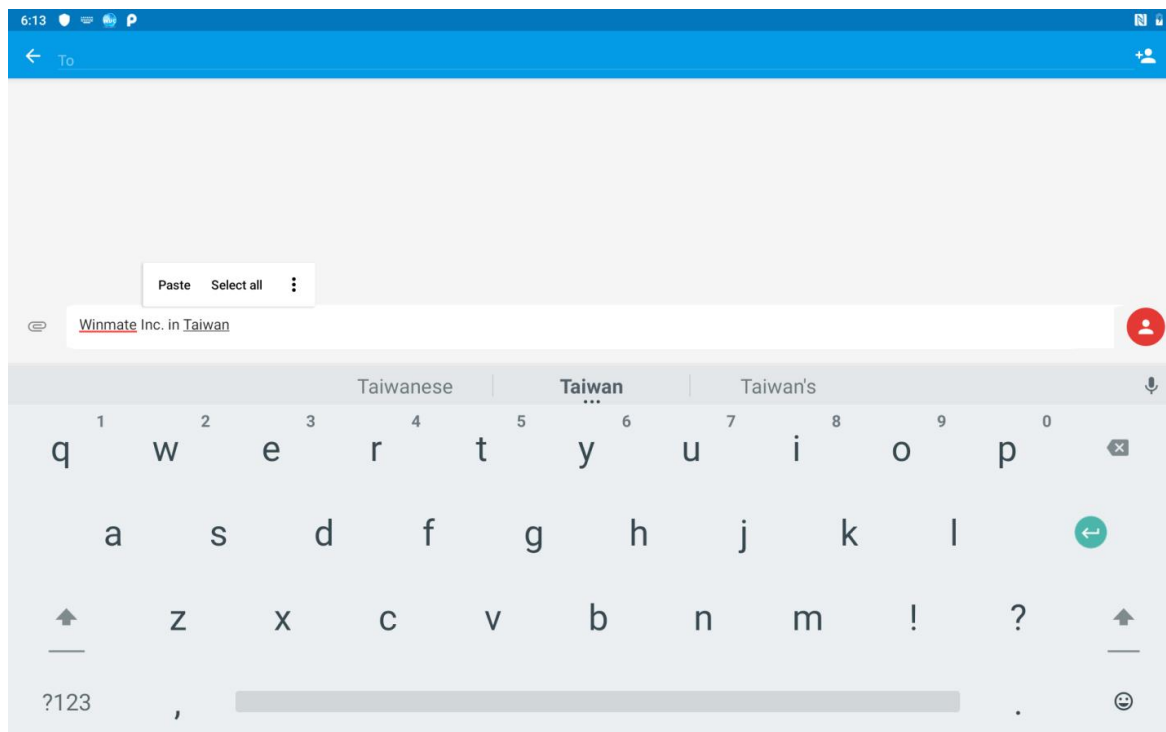
3.6 Display Settings

Use the display settings to configure brightness and other screen settings.

1. Swipe from the bottom of the main desktop, tap **Settings** , and then tap **Display** .

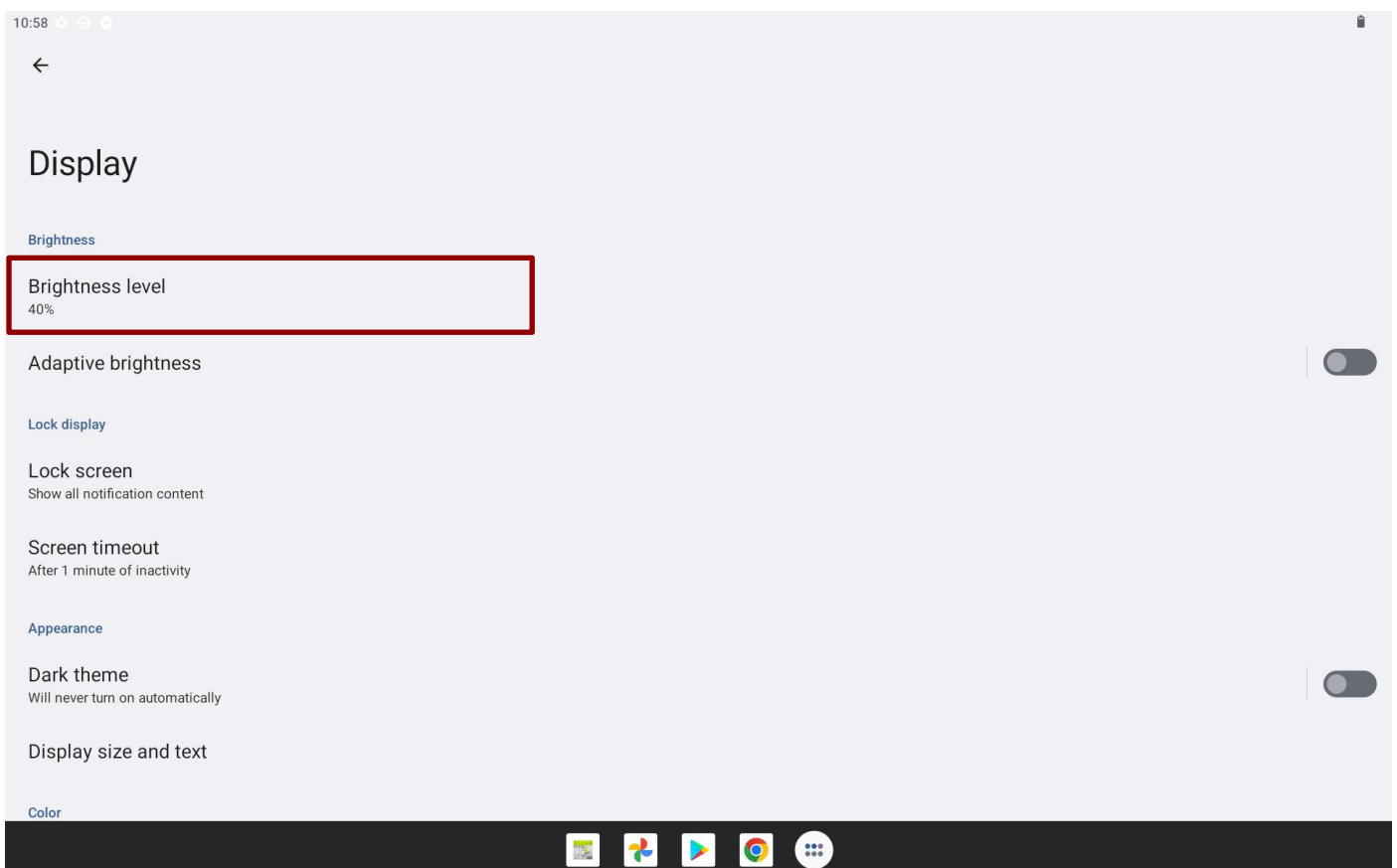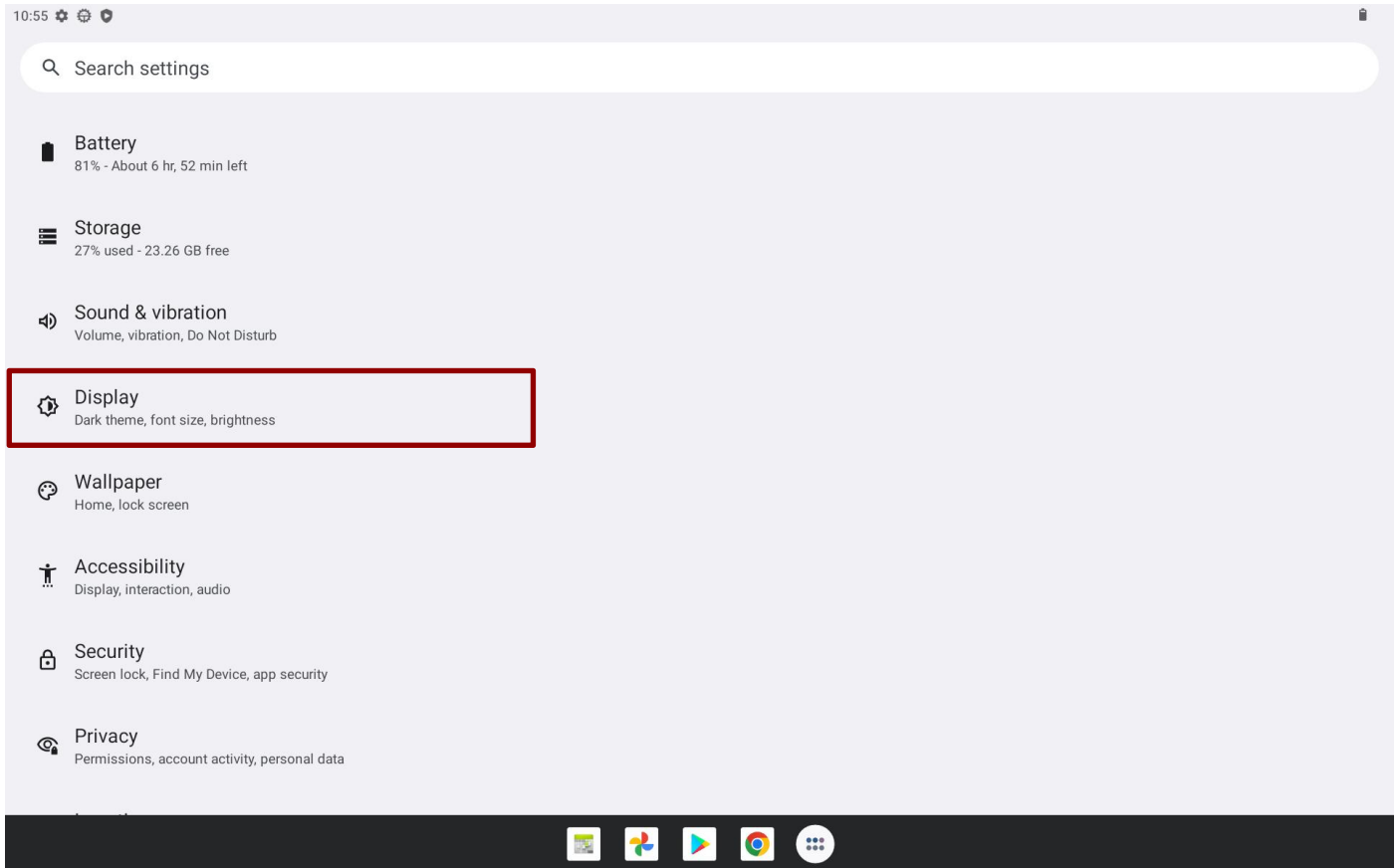Display

- **MiraVision™** – to display picture quality optimization options.
- **Brightness level** – for adjusting the brightness of the screen.
- **Dark theme** – uses black background to help keep battery alive longer on some screens.
- **Night light** – tints your screen amber. This makes it easier to look at your screen or read in dim light.
- **Adaptive brightness** – to optimize brightness level for available light
- **Wallpaper** – to select wallpaper
- **Screen timeout** – to set up the device’s sleep mode after the specified period of inactivity
- **Auto-rotate screen** – for switching the orientation of the screen when user turns the device sideways or upright.
- **Font size** – Make the text on screen smaller or larger
- **Display size** – Make the items on your screen smaller or larger. Some apps on your screen may change position.
- **Screen saver** – to control what happens when the tablet is docked and/or sleeping, turn screen saver on
- **Lock screen** – to control all notification content
- **HDMI setting** – to configure HDMI


3.7 Audio Profile Settings

Use the sound settings to configure how and at what volume the tablet rings, vibrates or alerts user in other ways when user receive a call, receive a notification or an alarm sound.

3.7.1 Sound Settings

Use the sound settings to configure how and at what volume the tablet rings, vibrates or alerts user in other ways when user receive a call, receive a notification or an alarm sound.

1. Swipe from the bottom of the main desktop, tap **Settings** , and then tap **Sound**  .

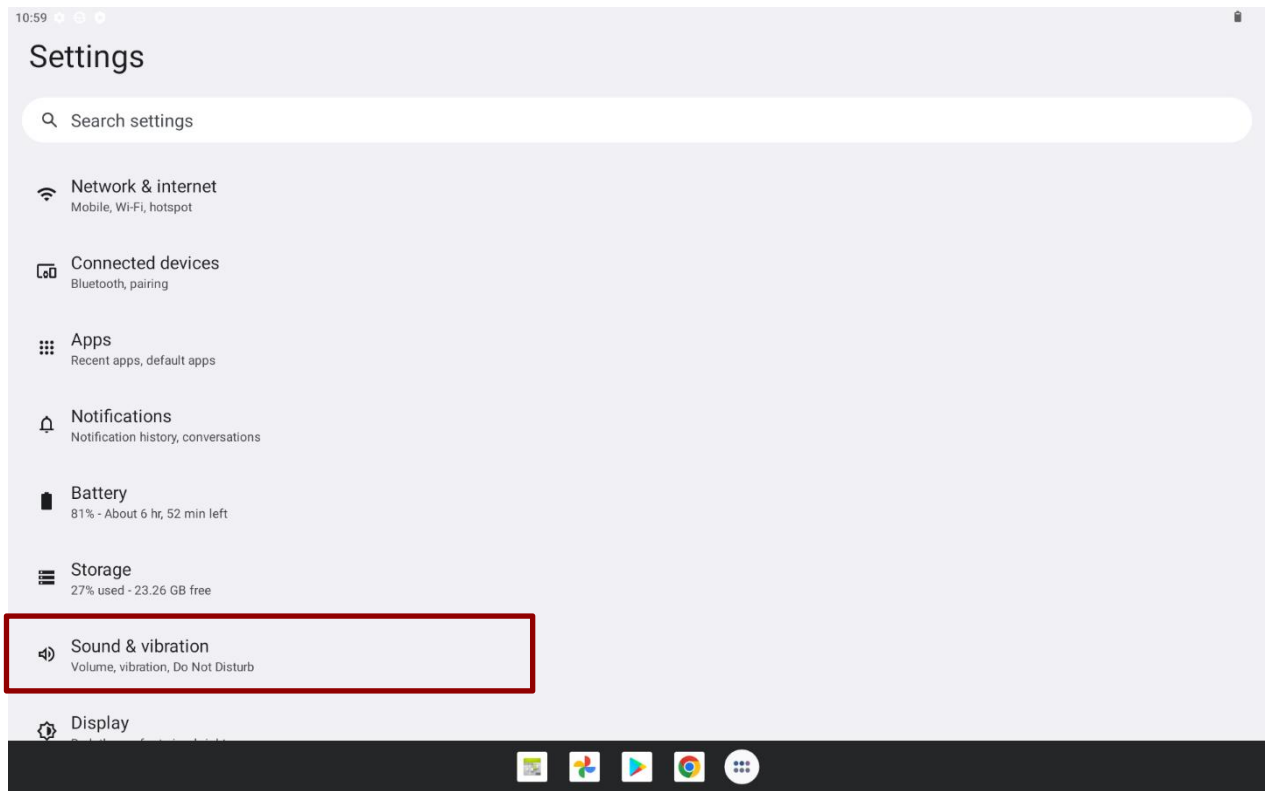
2. In the **Sound** user can select Media volume, Call volume, ring volume, alarm volume and other sound settings.


3.7.2 Adjusting Volume

1. Pull down the notification area and select the Menu Mode **Menu Mode** **Scan Mode** first.

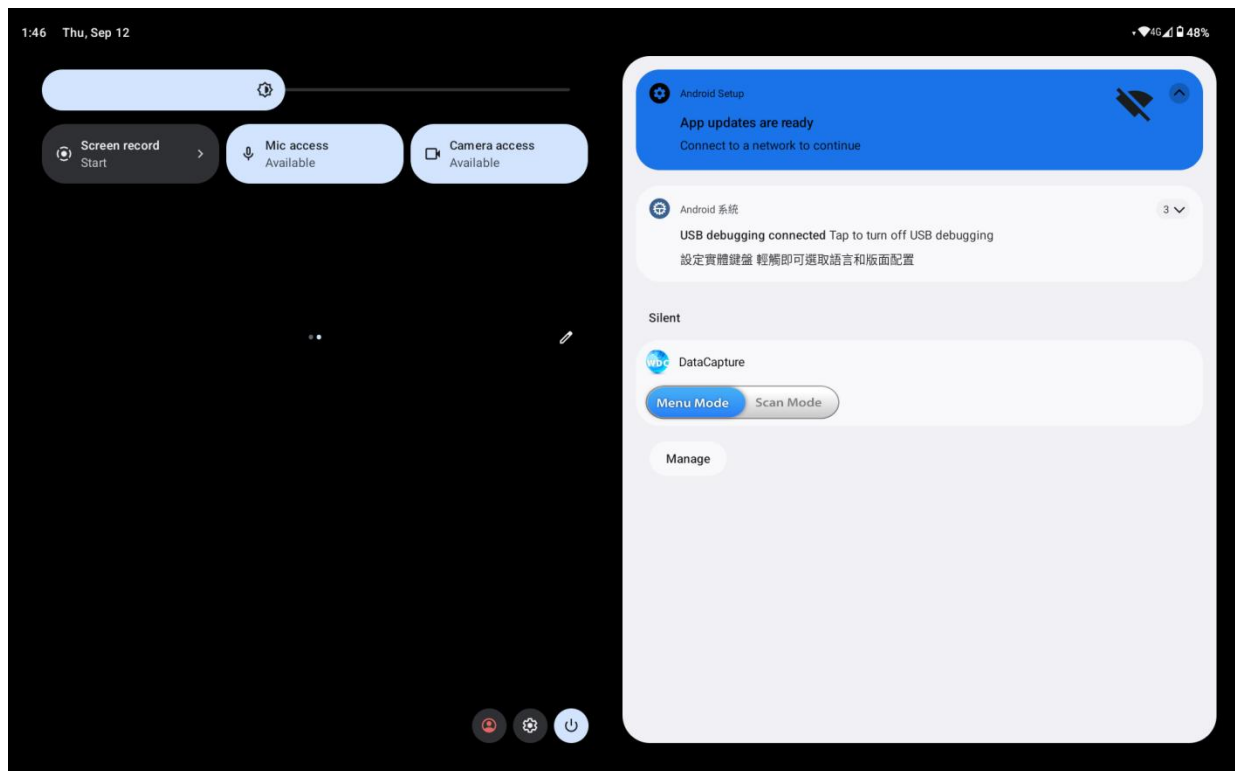
2. Press the volume keys to change the media volume (On Home screen) or earpiece volume (during a call). Press the volume key again to set up Vibration mode.

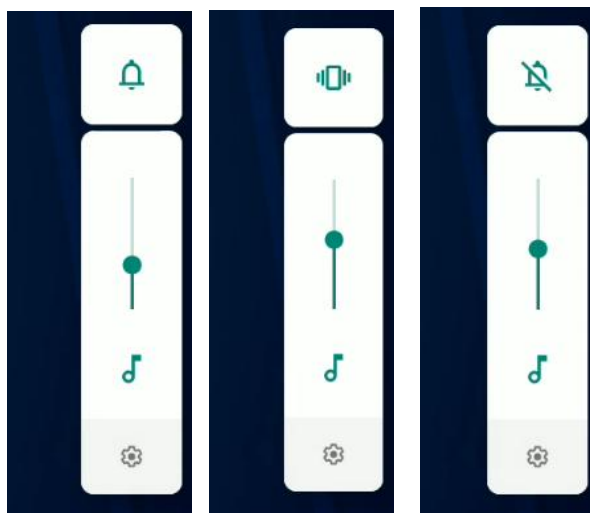
3.8 Language and Keyboard Settings

Use the Language & keyboard settings to select the language for the text on your tablet and for configuring On-Screen Keyboard, including words that you have added to its dictionary.

1. Swipe from the bottom of the main desktop, and then tap **Settings** , and then tap **System** .


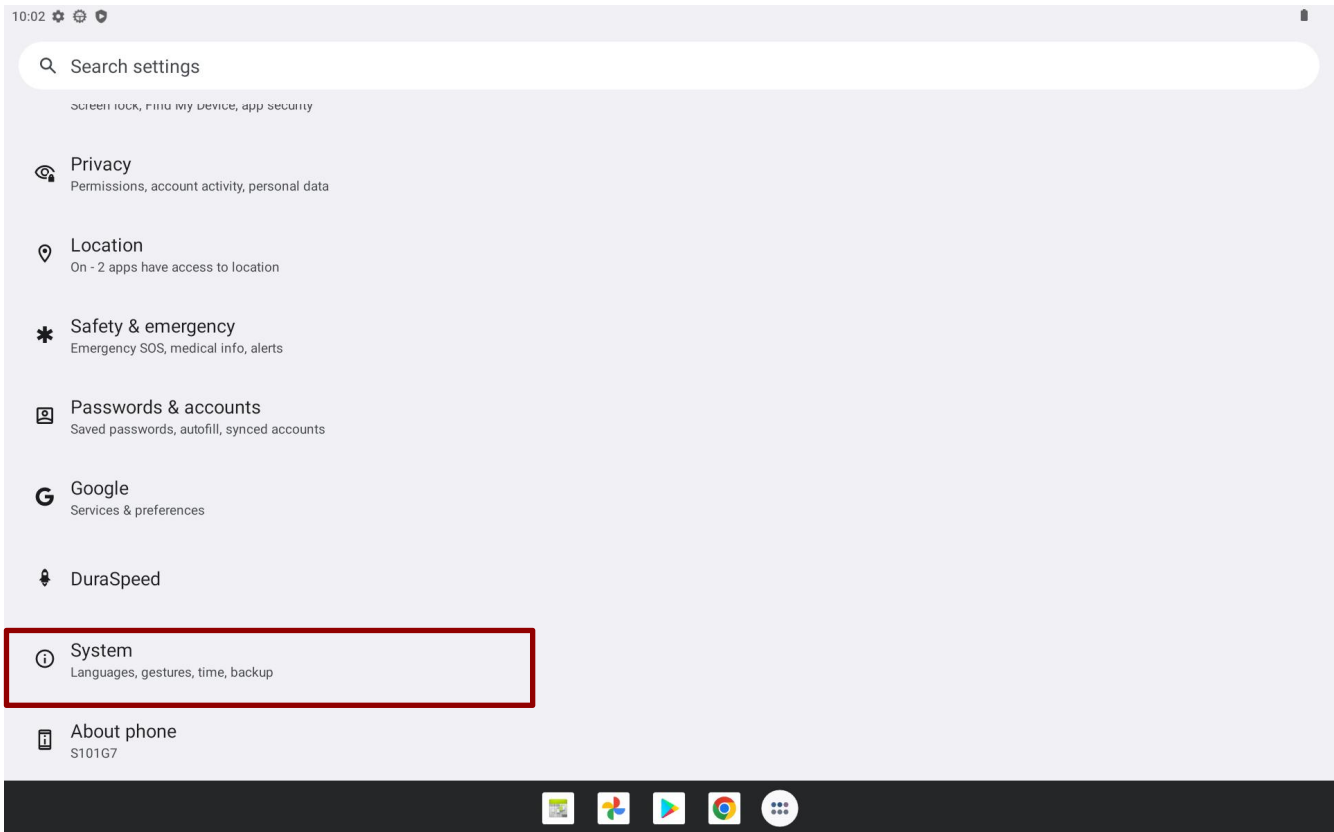
2. In the **Personal** tab, select **Language & Input** .

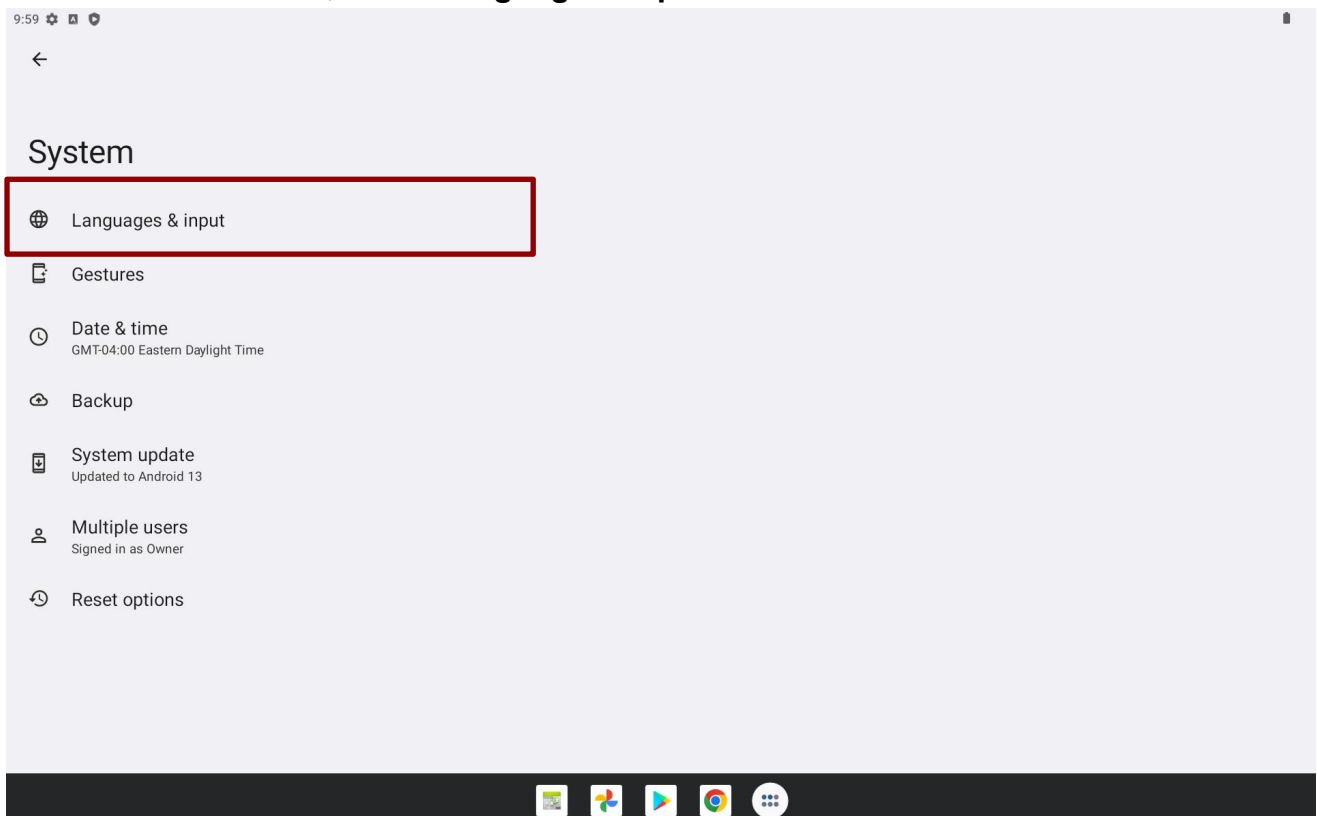- **Keyboard & input methods**

- **On-screen Keyboard** – shows default Android Keyboard (AOSP).
- **Physical Keyboard** – allows the user to apply physical keyboard settings if connected.

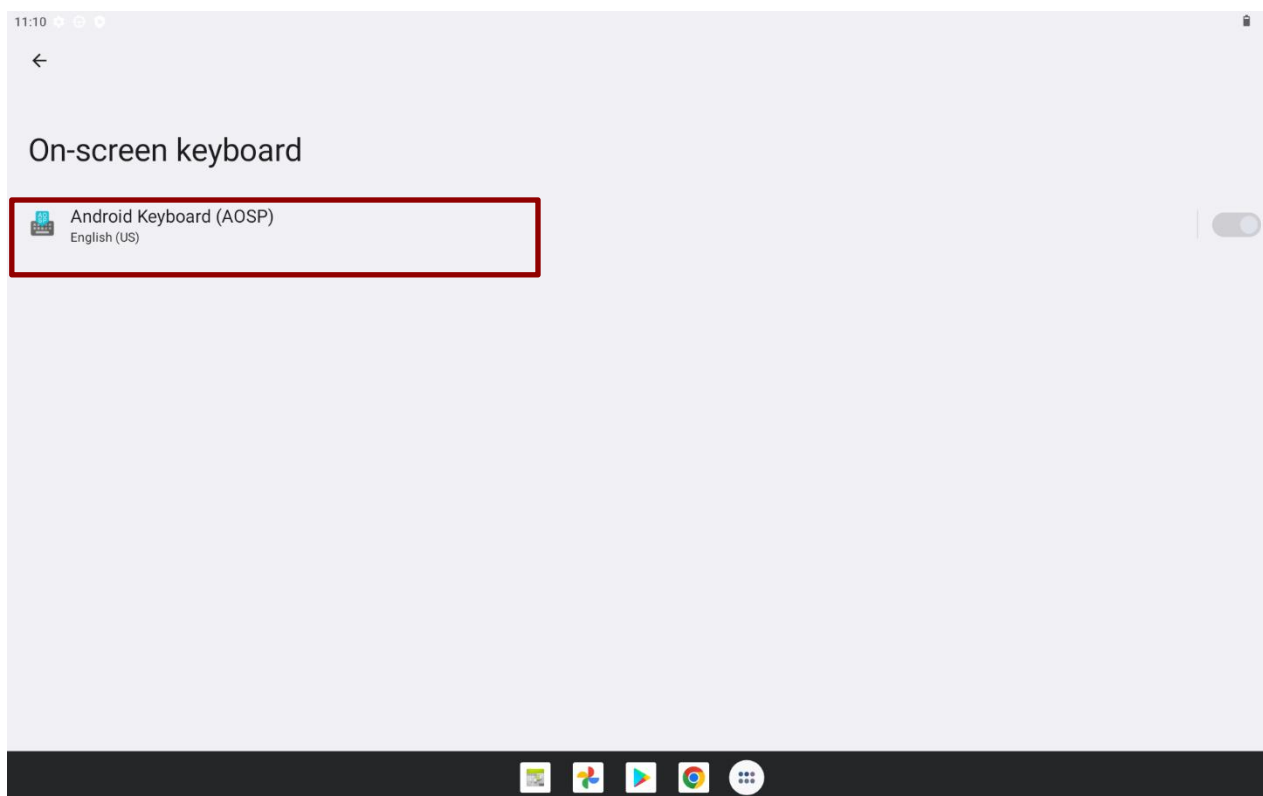
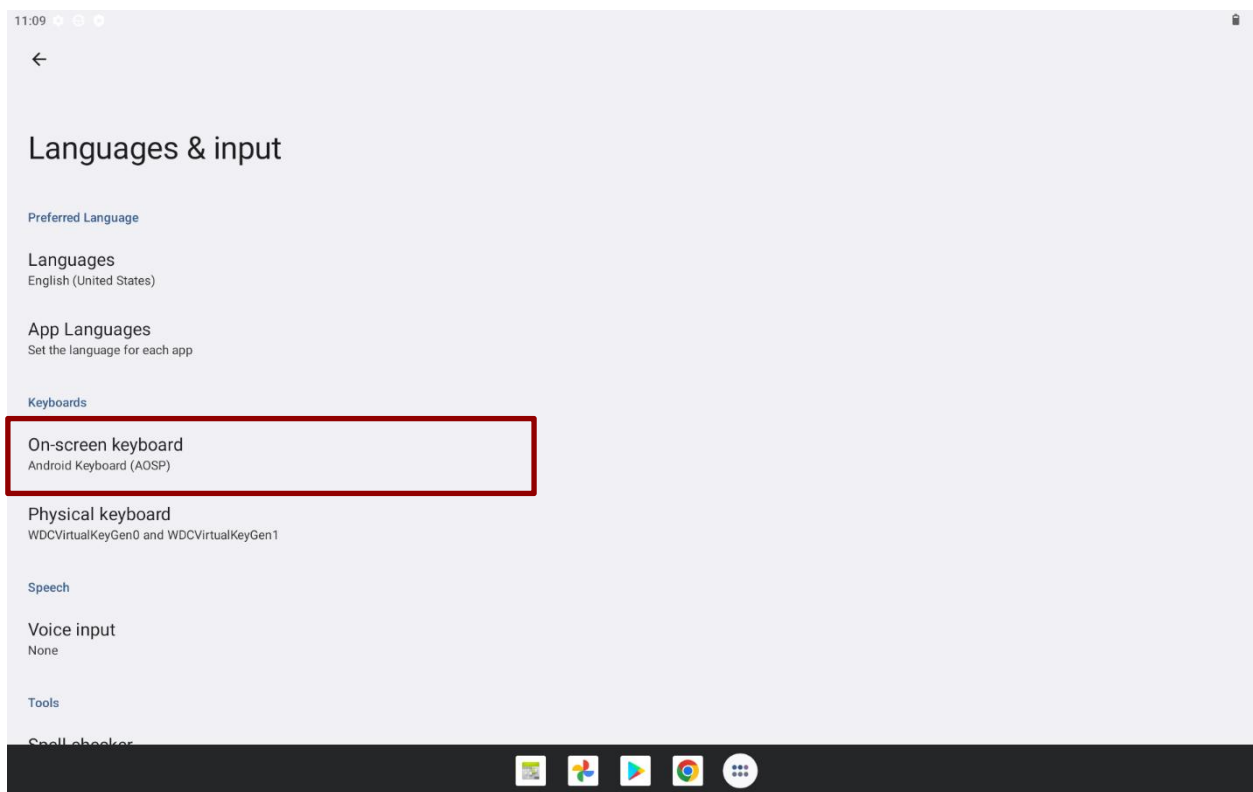
- **Appearance & layouts**

- **Theme** – to select the theme of layout.
- **Custom input styles** – user can select custom input style, such as QWERTY keyboard.
- **Enable split keyboard** – toggle to activate split keyboard layout.

- **Text Correction**

- **Personal dictionaries**
- **Add-on dictionaries**
- **Block offensive words** – disable suggesting potentially offensive words
- **Auto-correction** – spacebar and punctuation automatically correct mistyped words.
- **Show correction suggestions** – display suggested words while typing.
- **Personalized suggestions** – learn from your communications and typed data to improve suggestions.
- **Suggest contact names** – use names from Contacts for suggestions and corrections.
- **Next-word suggestion** – use previous word in making suggestion.

- **Advanced**

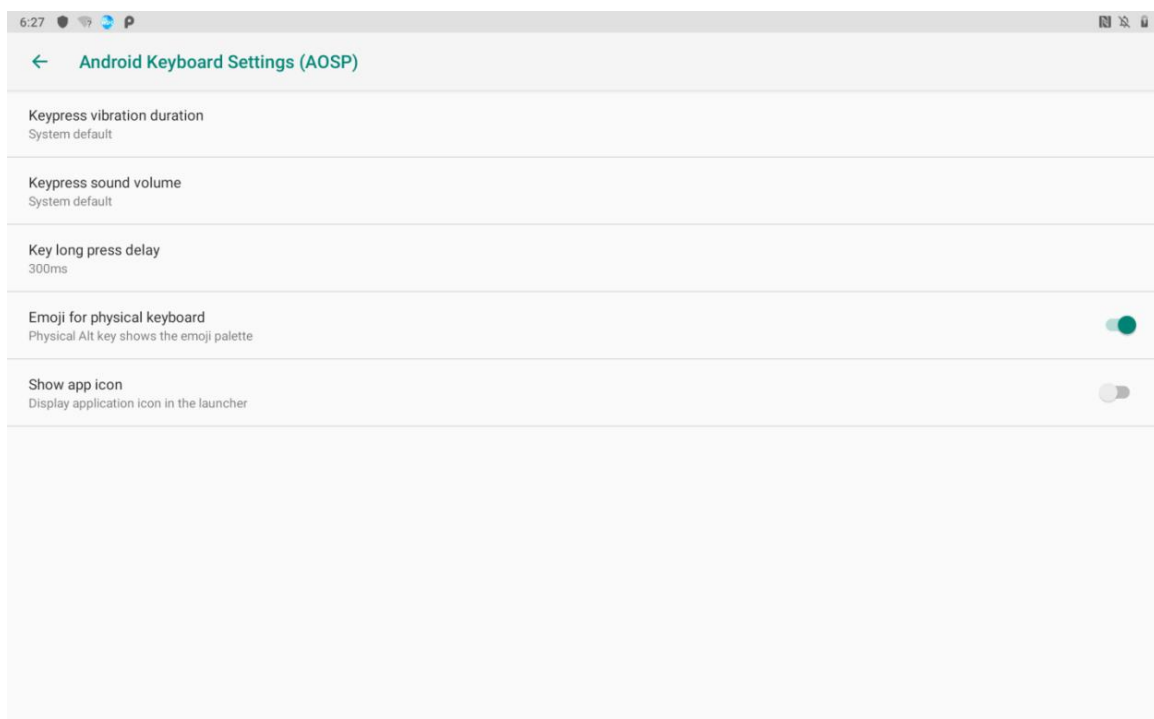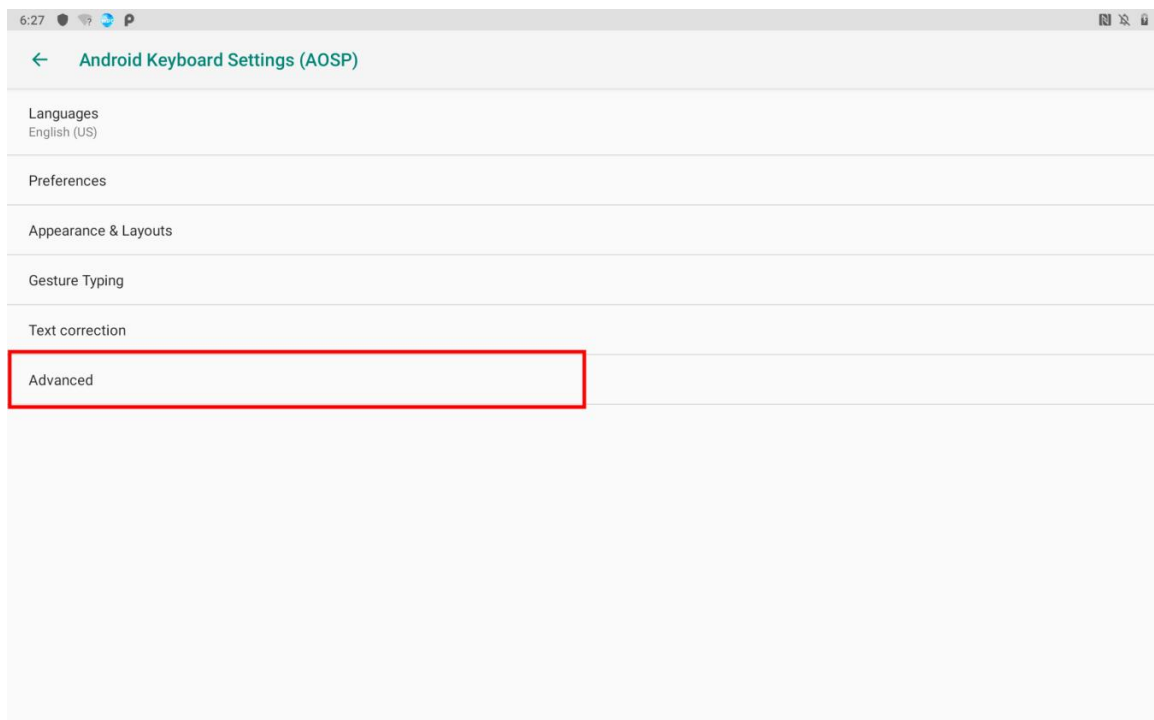- ***Keypress vibration duration***
- ***Keypress sound volume***
- ***Key long press delay***
- ***Emoji for physical keyboard***
- ***Show app icon***

Changing the Keyboard Language

The current keyboard language is displayed on the Space key and user can switch languages.

1. Tap and hold the space key.
2. Without lifting the finger, slide it to the left or right until the language you want to use appears in the center of the small window above the Space Key.
3. Lift the finger to change the keyboard language.


3.9 Changing Touchscreen Modes

The touch screen of the MD10W2 Series Rugged Tablet supports hand/ Water Resistant, Stylus and glove modes. You can change the touch mode in Accessibility settings.


Mode	Description
Hand/ Water Resistant Mode (Default)	In <i>Hand/ Water Resistant</i> mode the touchscreen rejects false inputs from liquid drops. When the liquid is removed, the touchscreen accepts touch inputs from bare fingertips. This is a default mode setting. <i>Supports input by:</i> • <i>Rubber tip passive stylus</i> • <i>Bare fingertips</i>
Stylus Mode	In <i>Stylus</i> mode, the touchscreen supports inputs from the fine-tip passive capacitive stylus (shipped with tablet). <i>Support input by:</i> • <i>Fine-tip passive stylus with palm rejection</i> • <i>Rubber tip passive stylus</i> • <i>Bare fingertips</i>
Glove Mode	In <i>Glove</i> mode, the touchscreen supports inputs while the user is operating the touchscreen with gloves on. <i>Supports input by:</i> • <i>Glove</i> • <i>Rubber tip passive stylus</i> • <i>Bare fingertips</i>

To change touch mode:


1. Open **Settings > Accessibility > Touch Mode**.

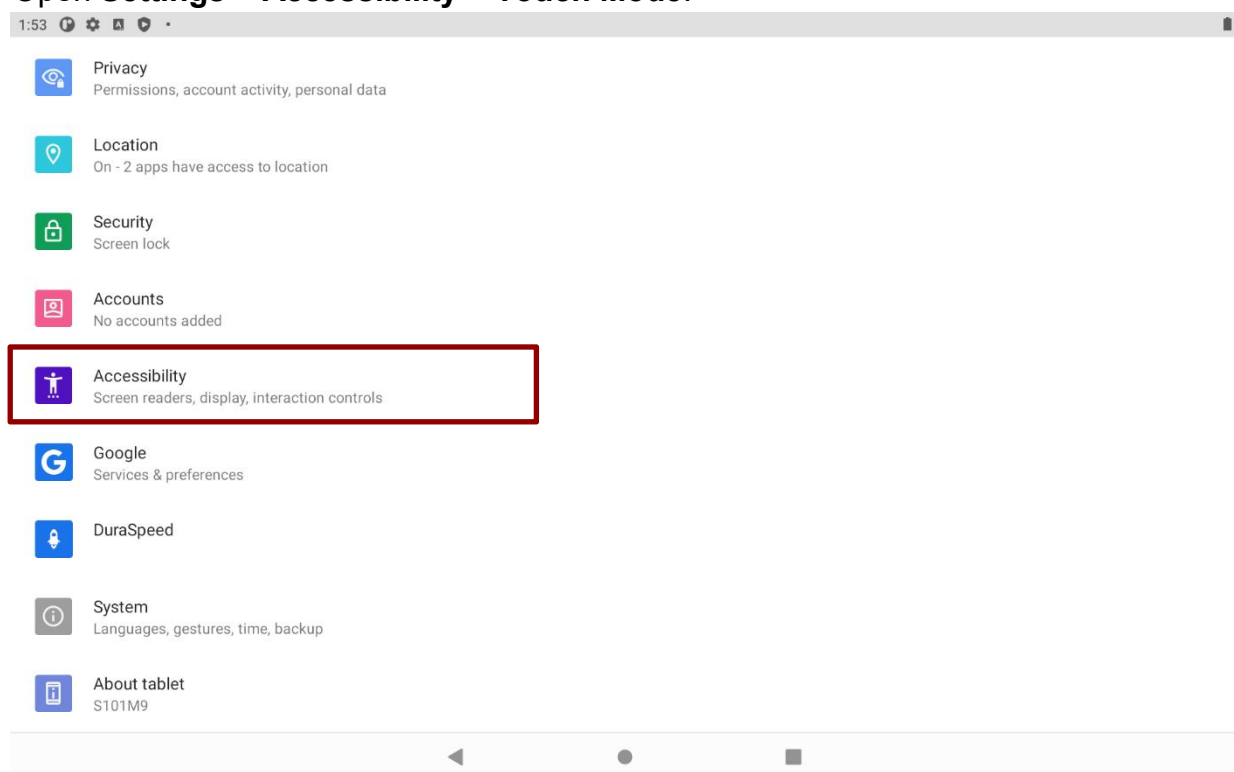2. Select touch mode. Follow the table below for more information about each touch mode.

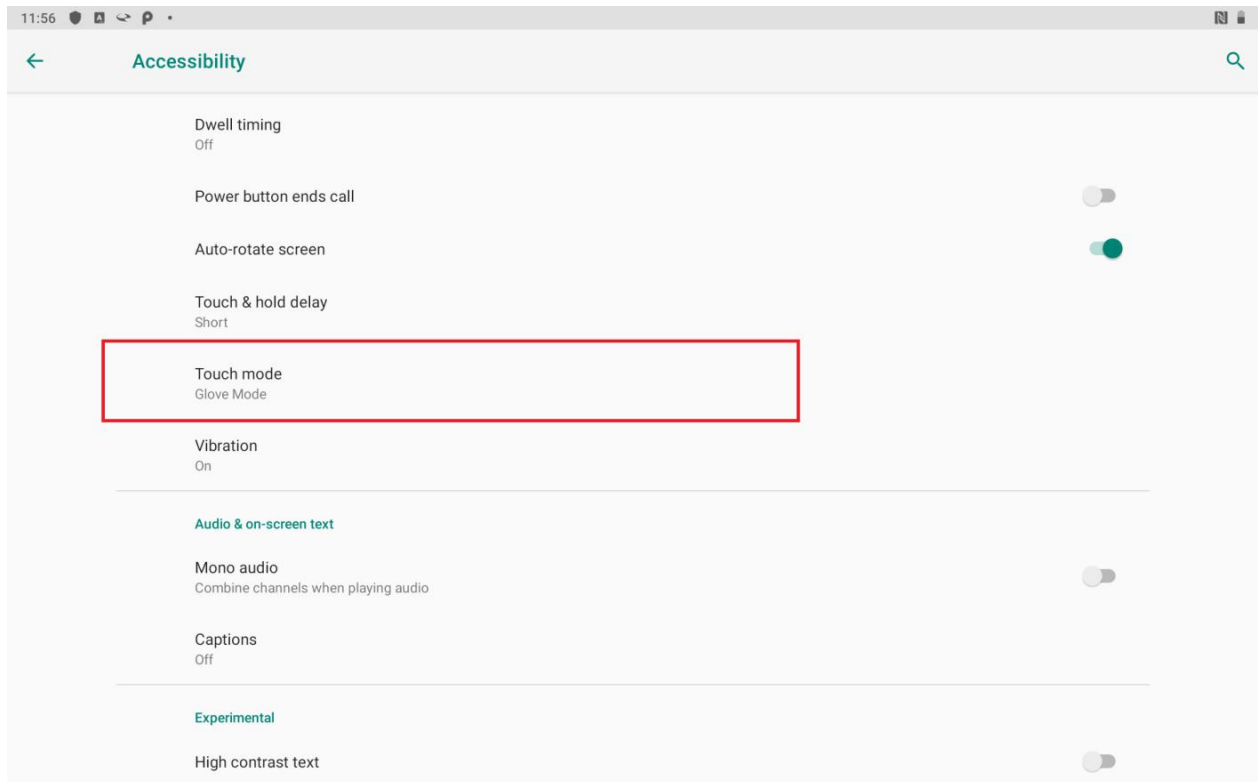
3. Select the touch mode and leave the Touch Mode settings.

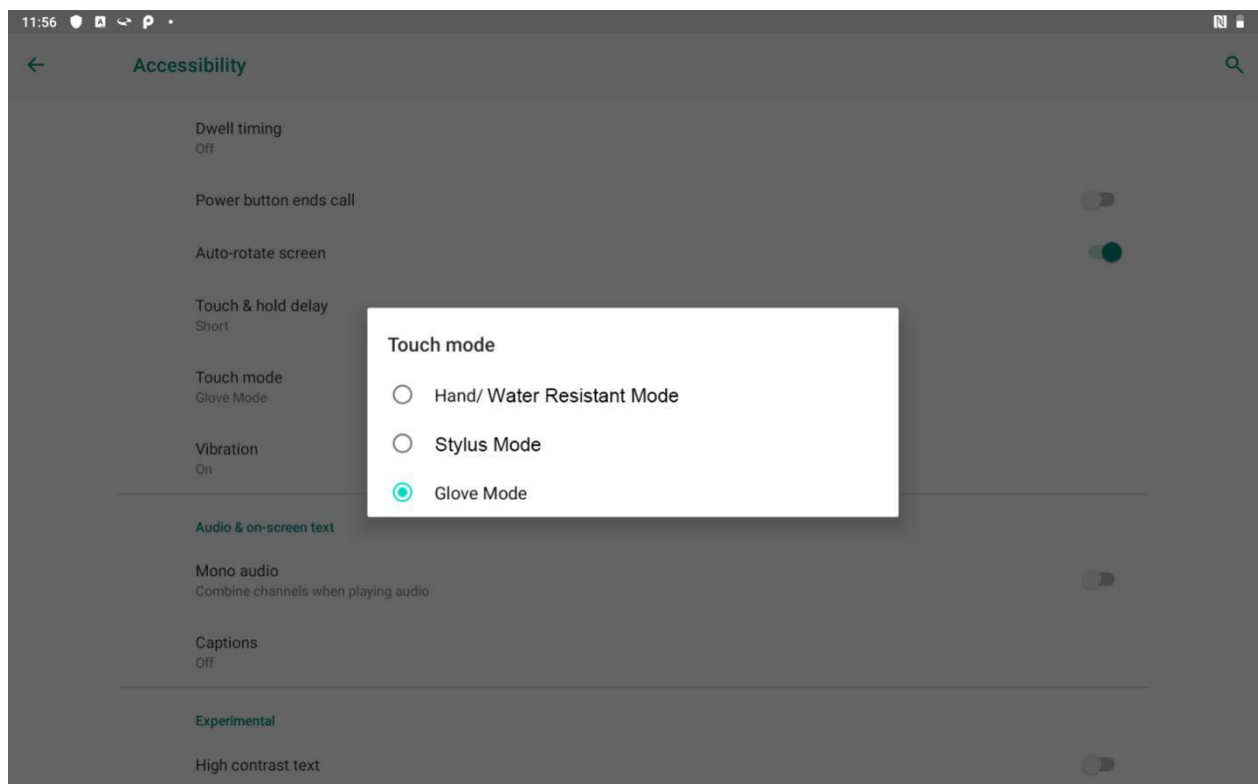
3.10 Extreme power saving setting

After the MD10W2 Tablet enters sleep mode, LTE functions will disconnect unless the "Extreme power saving" is set to "never"


1. Swipe from the bottom of the main desktop, tap "Settings".
2. Select battery mode

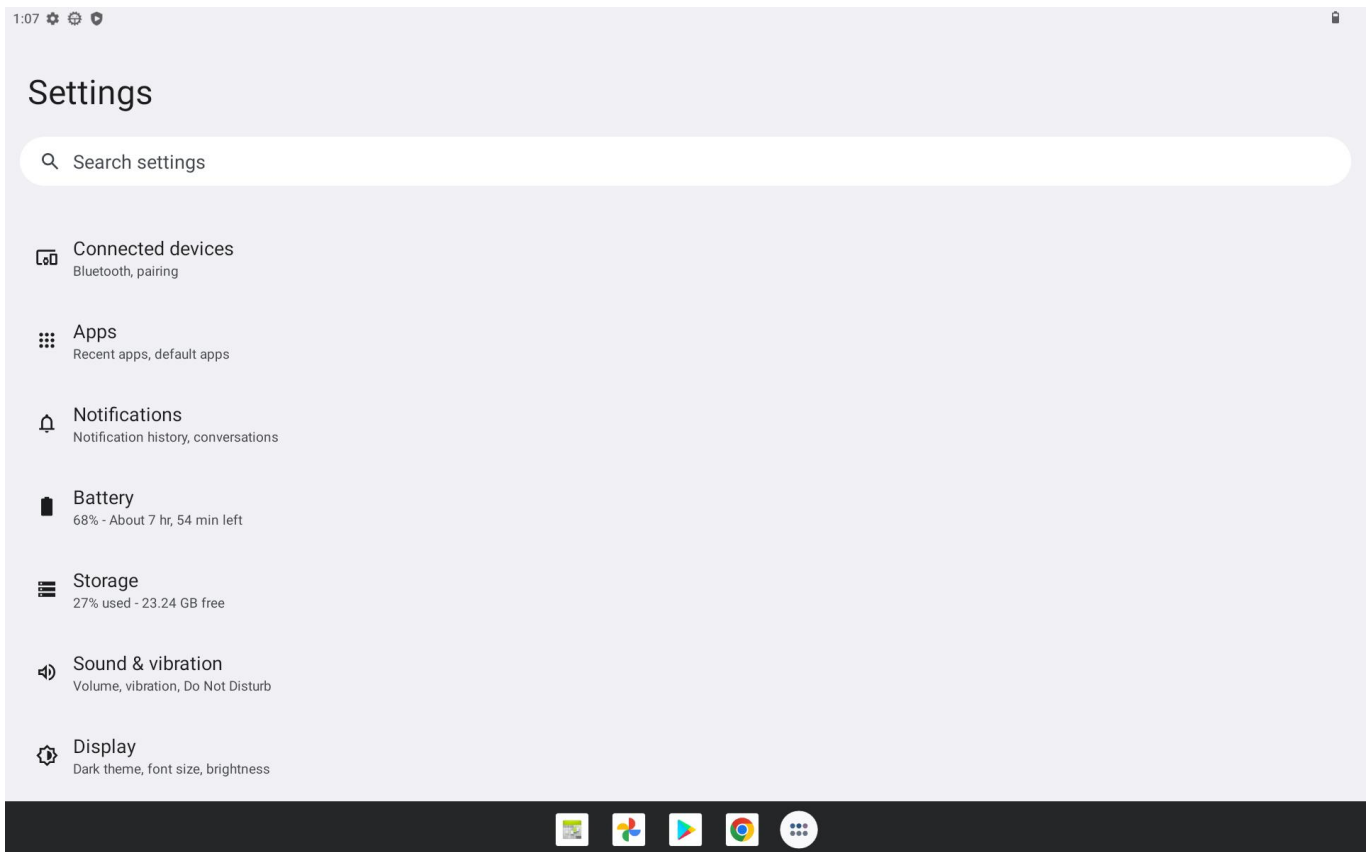3. In the battery mode, select Extreme Power Saving

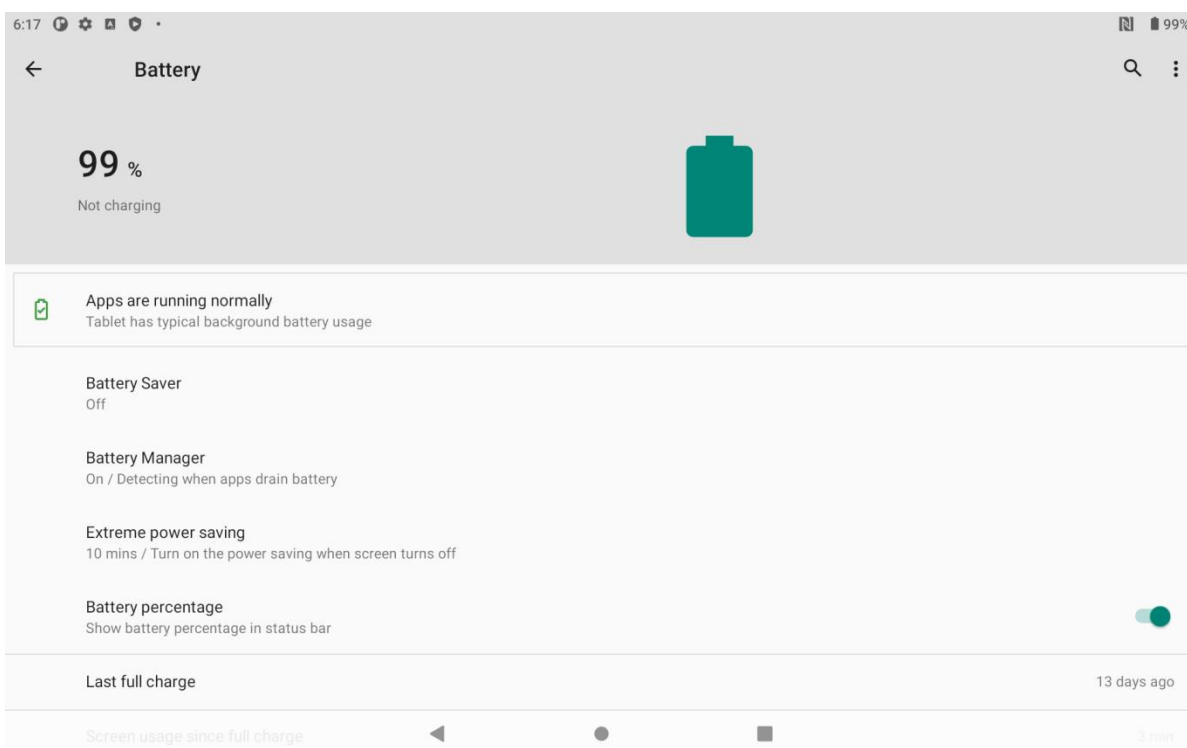4. Set up for extreme power saving to never

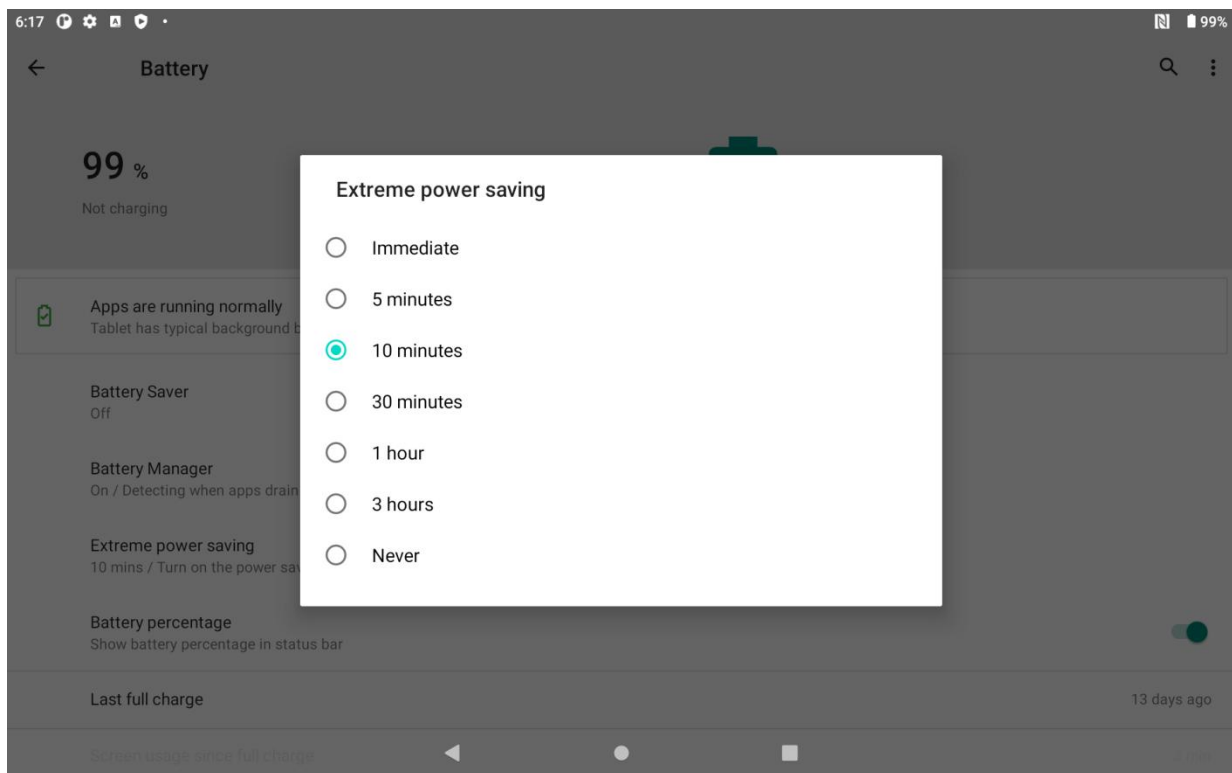
3.11 About Phone

About Phone includes information about user's device.

- **System Update** — provides information regarding the local or online update of your devices.
- **Status** — opens the status screen with a list of information about battery, network, signal strength and other details.
- **Legal Information** — lists details about the model number, Android Version, Baseband version etc. this information is useful if user need to work with support team.
- **Model number** – displays the model name of the current device.
- **Android version** – displays the current Android operating system version the device running

To access the information about the phone:

1. Swipe from the bottom of the main desktop, tap **Settings** ⚙️, tap **About Phone** ⓘ.


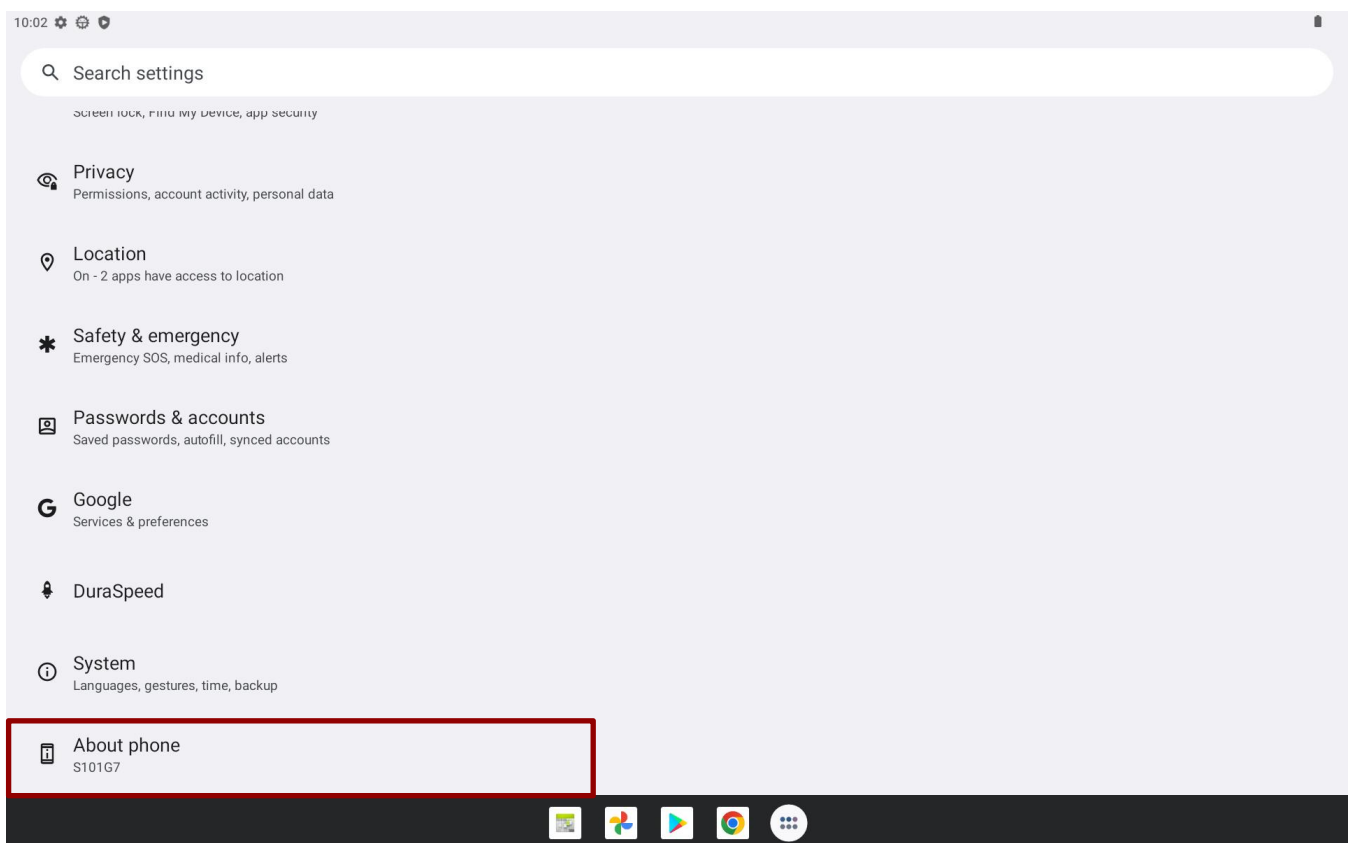
Chapter 4: Winmate Kiosk Manager


The MD10W2 Series Rugged Tablet comes with preinstalled APP named Winmate® Kiosk Manager (WKM). This application lets the device manager to lock the launcher. When the launcher is locked the user can only use one application assigned by the device manager.

4.1 Lock Applications Launcher


1. In the main APPS window, select the **Kiosk Manager** . Select **Kiosk Manager** to start your settings. Go to **Password Manager** and set your own password for this APP. The default password is “password”.

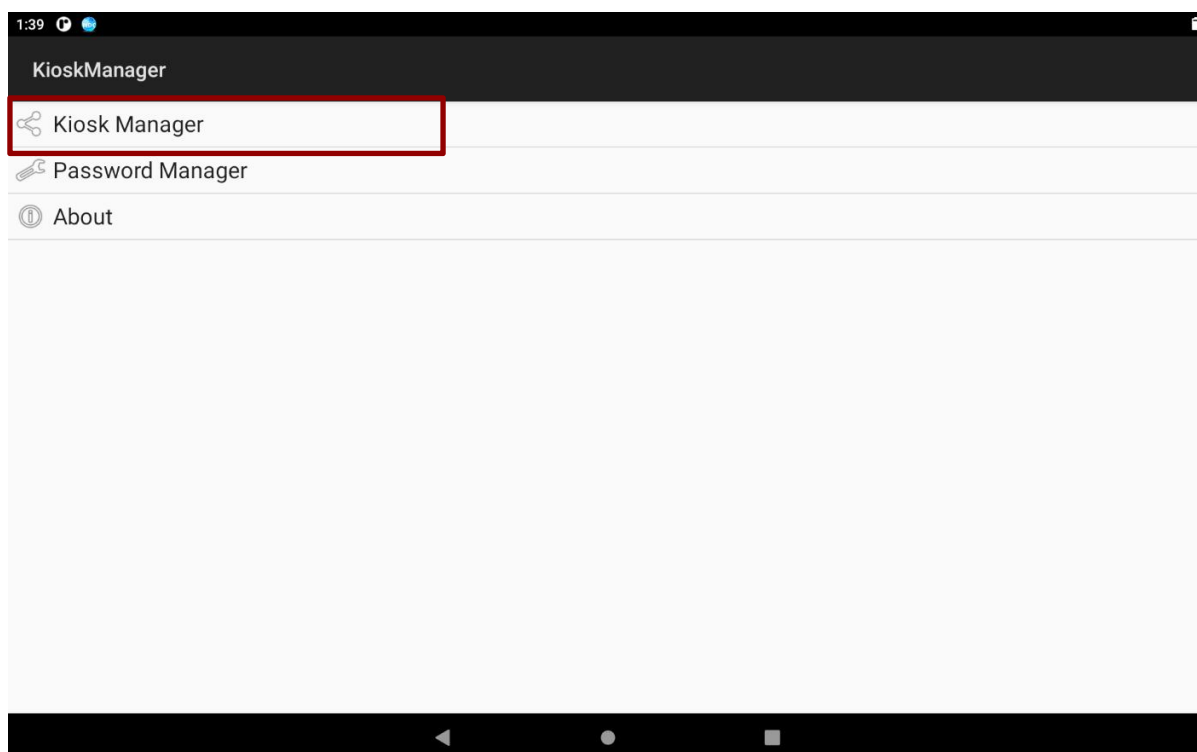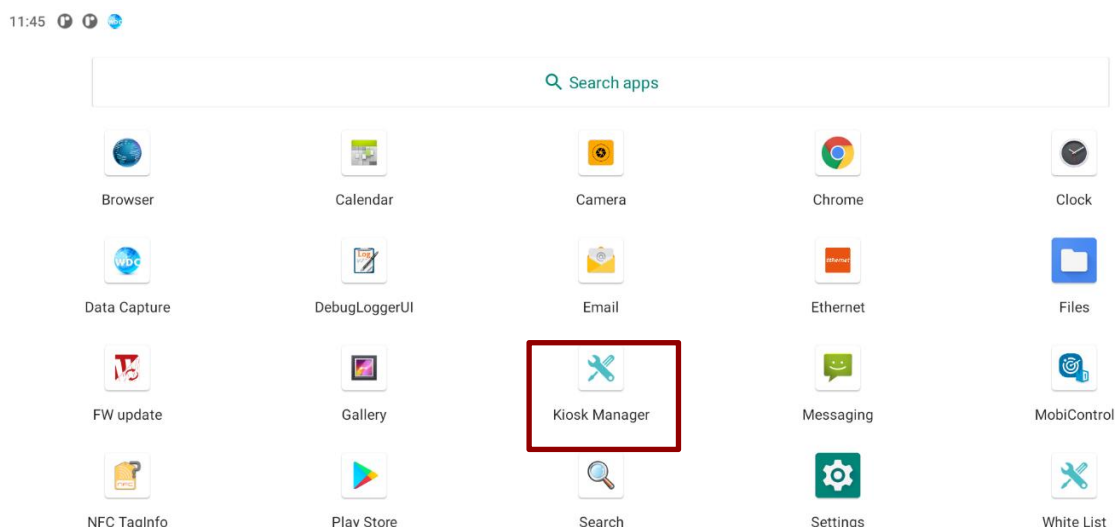
2. Press **Clear and Reset Default Home**.

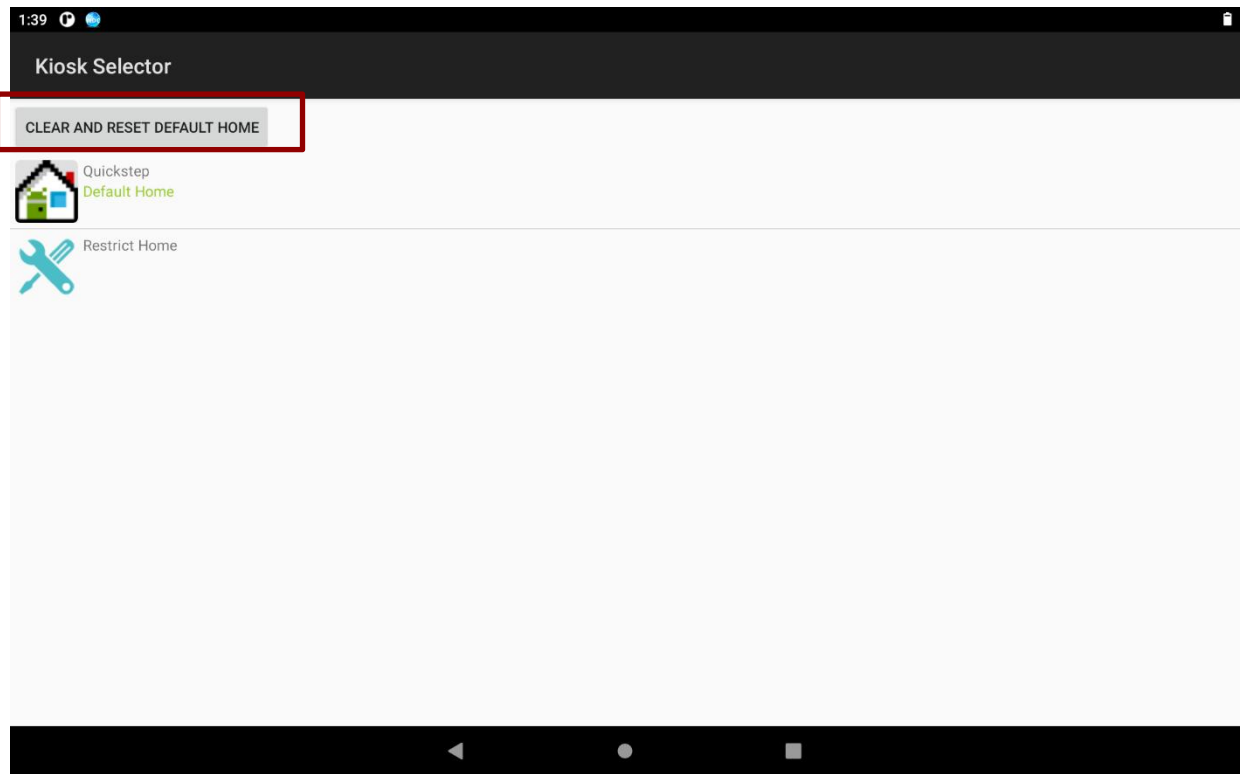
3. Select **Restrict Home** then **Always**. The system will automatically be locked on selected Apps.

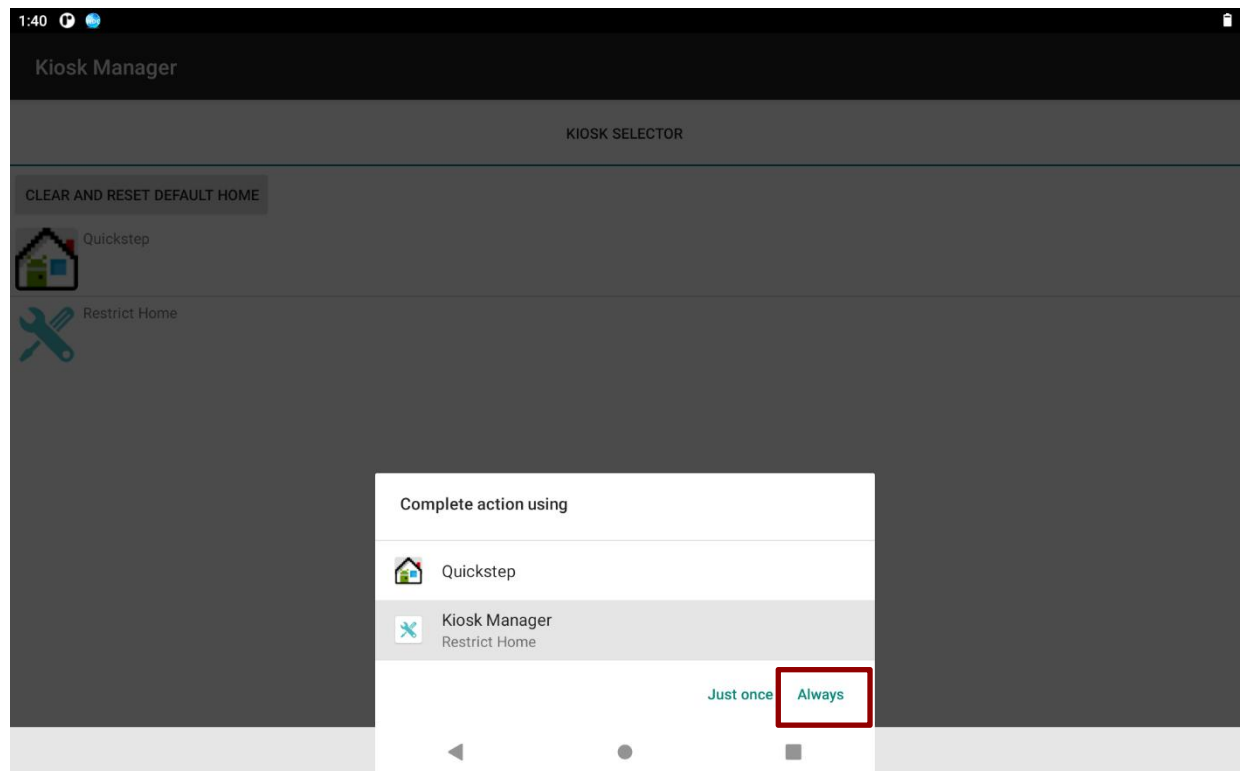
4. Click the gear icon on the upper right corner. Then choose **Change UI**.

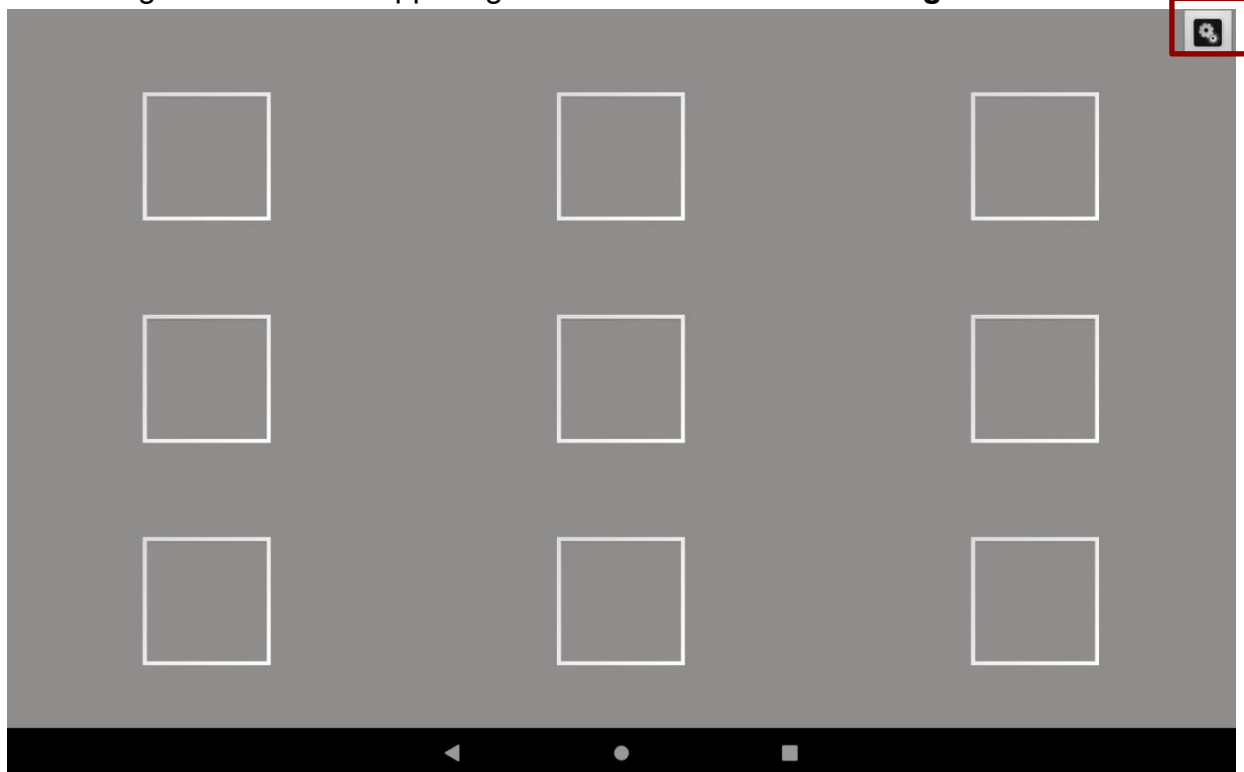
5. Input the default password: **password**

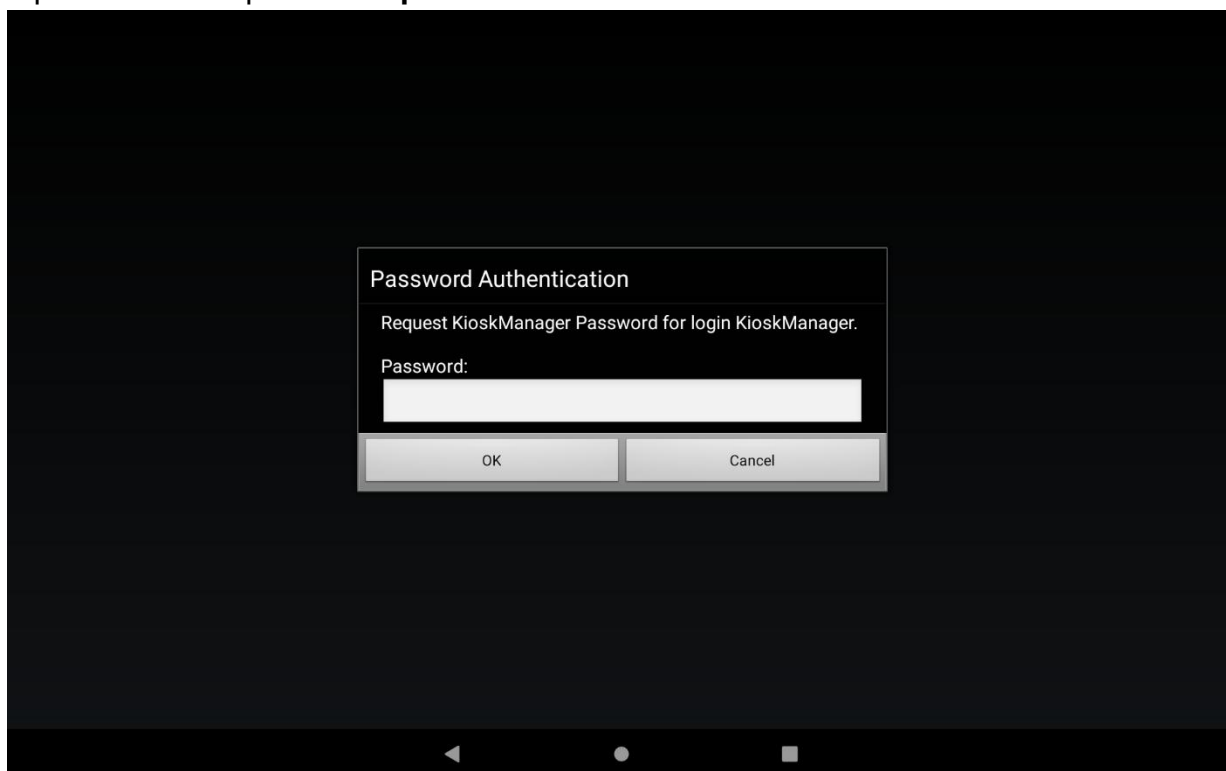
6. Choose the layout UI, then press the empty slot to add/remove apps.

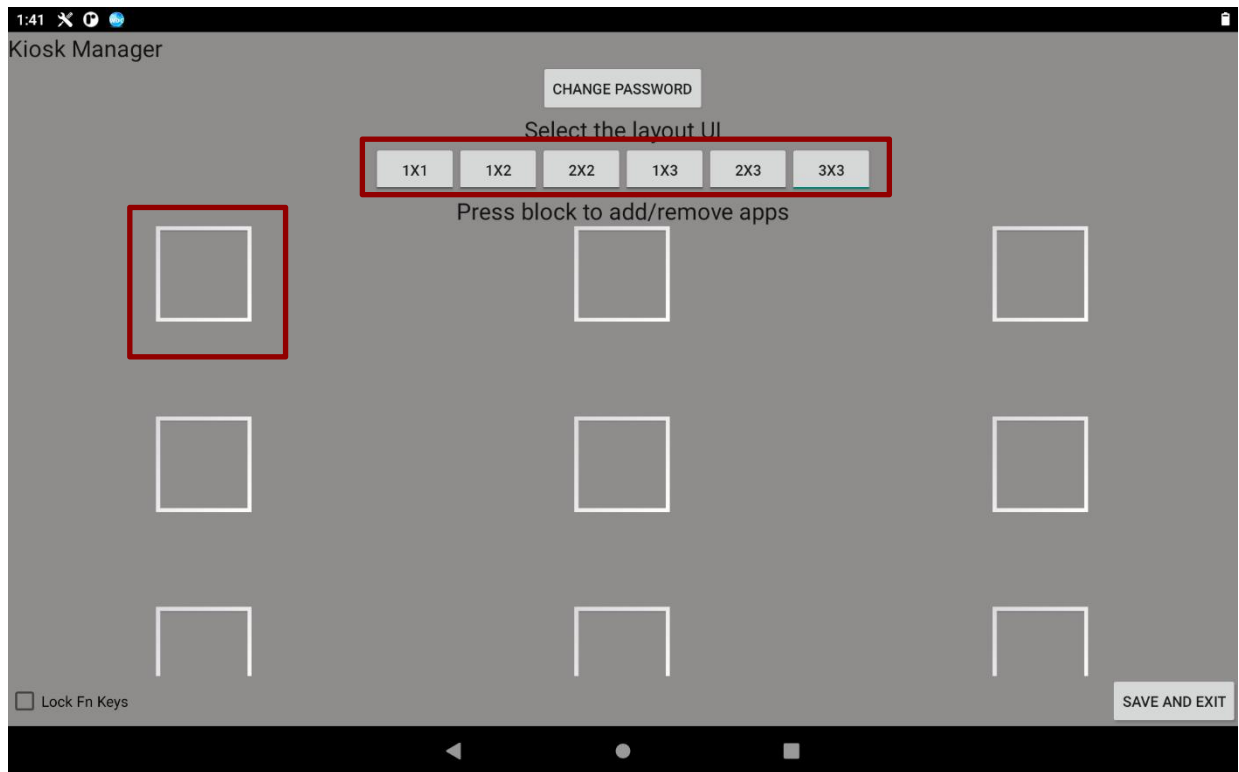
7. Click **Add App**.


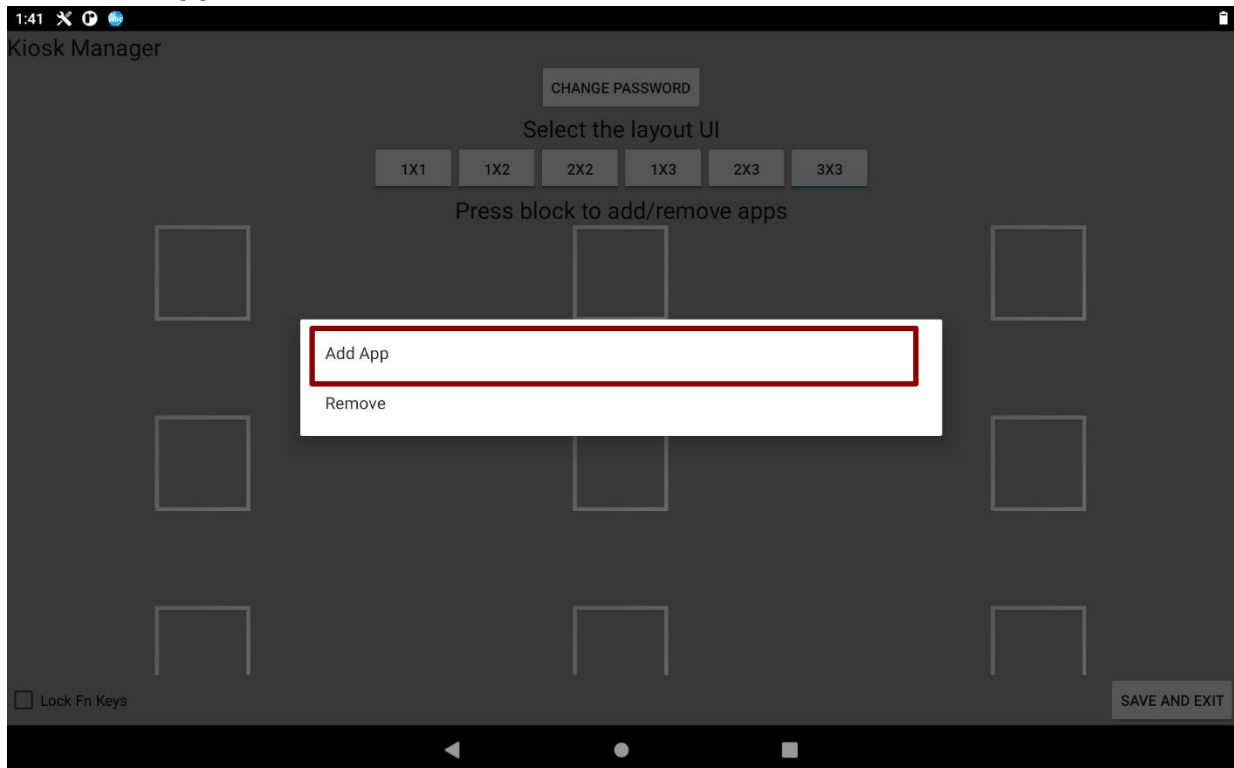
8. Select applications you want to lock. For example, **Ethernet**.

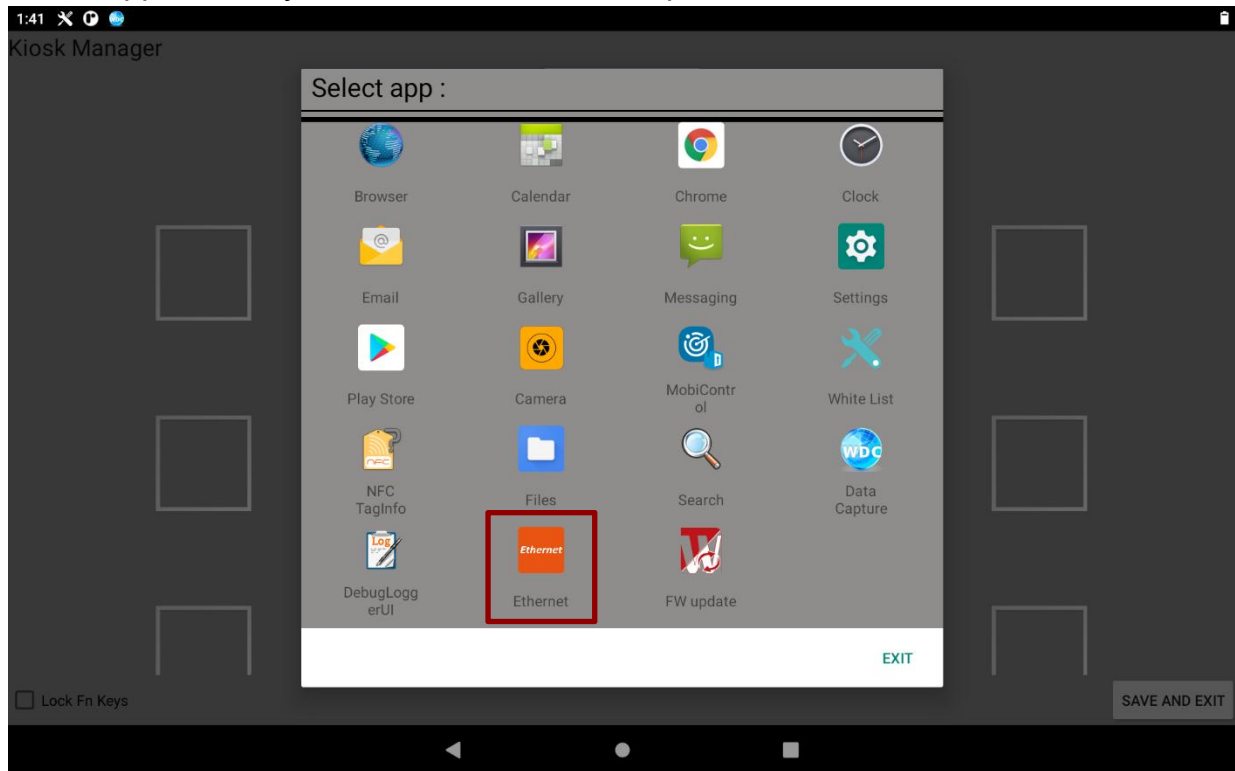9. Click **SAVE AND EXIT**.


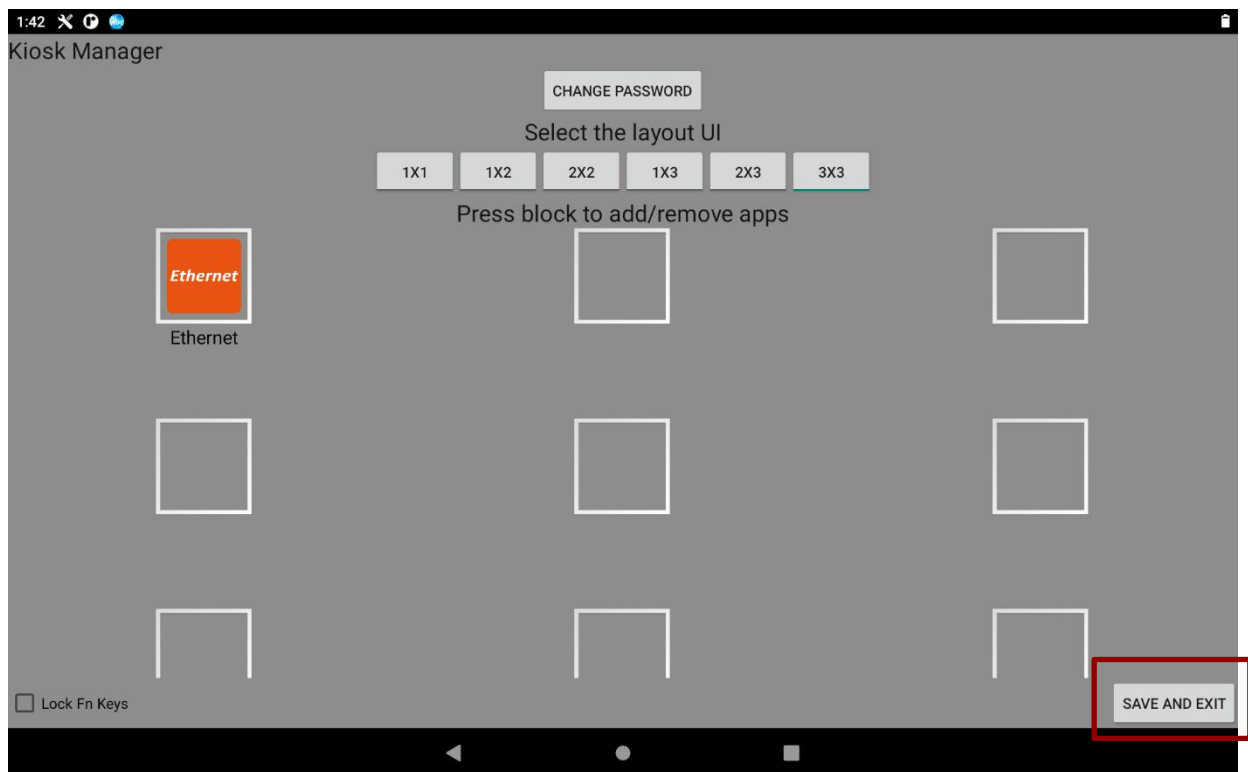
10. Click the selected Apps.

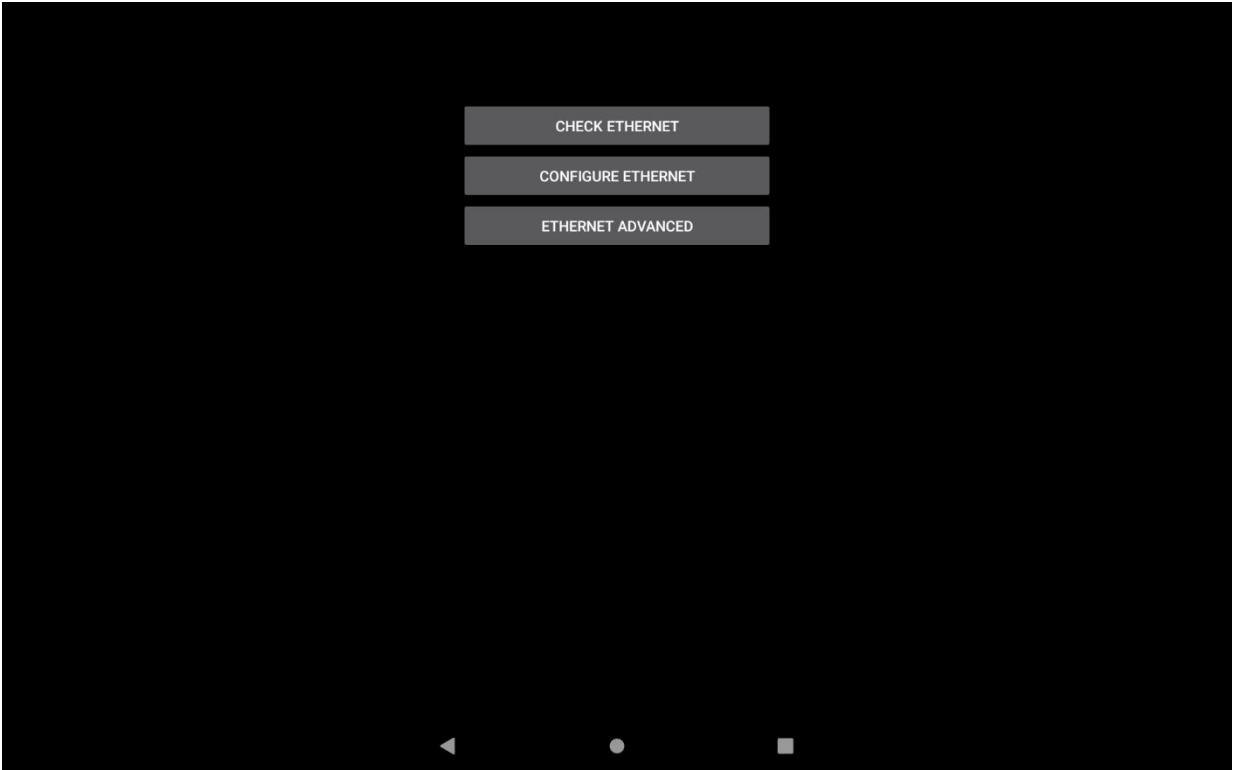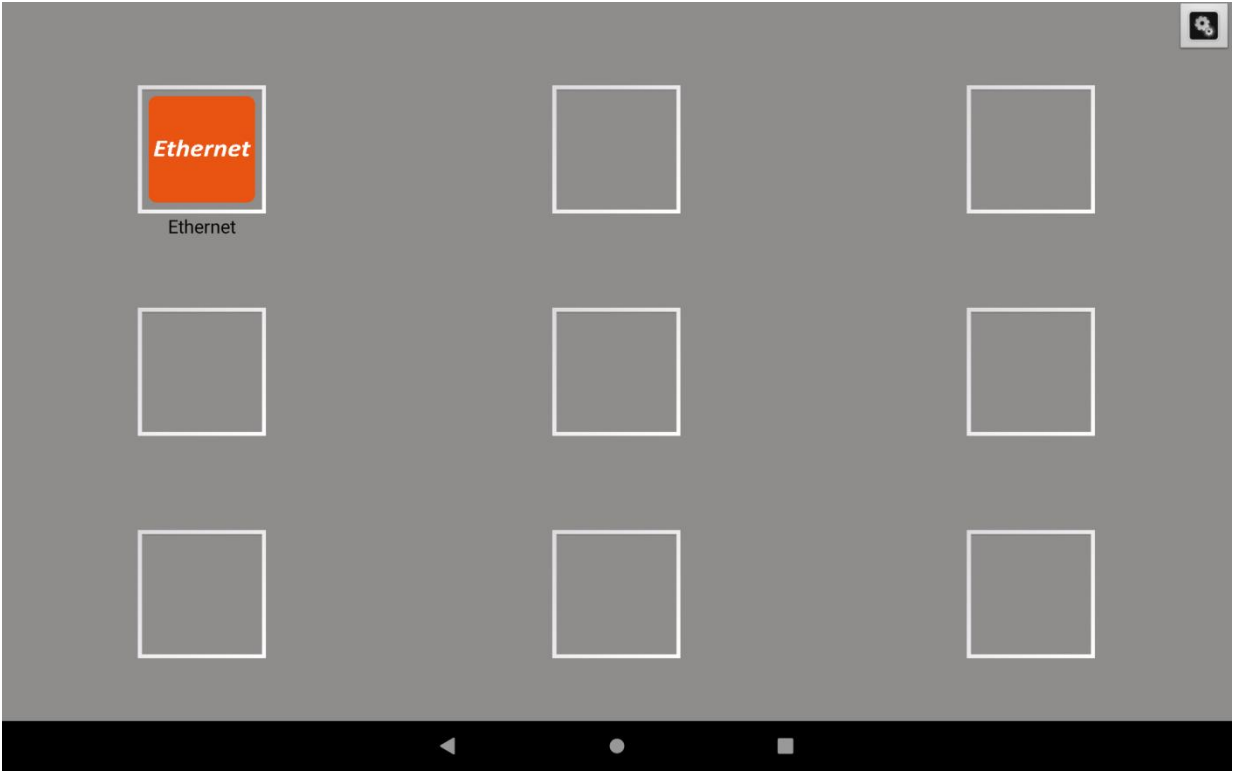
4.2 Unlock Applications Launcher


1. Slide down the top bar and press **Kiosk Manager** icon, click **Exit** and please enter the password to login (the default password is **password**). Press **Ok**.

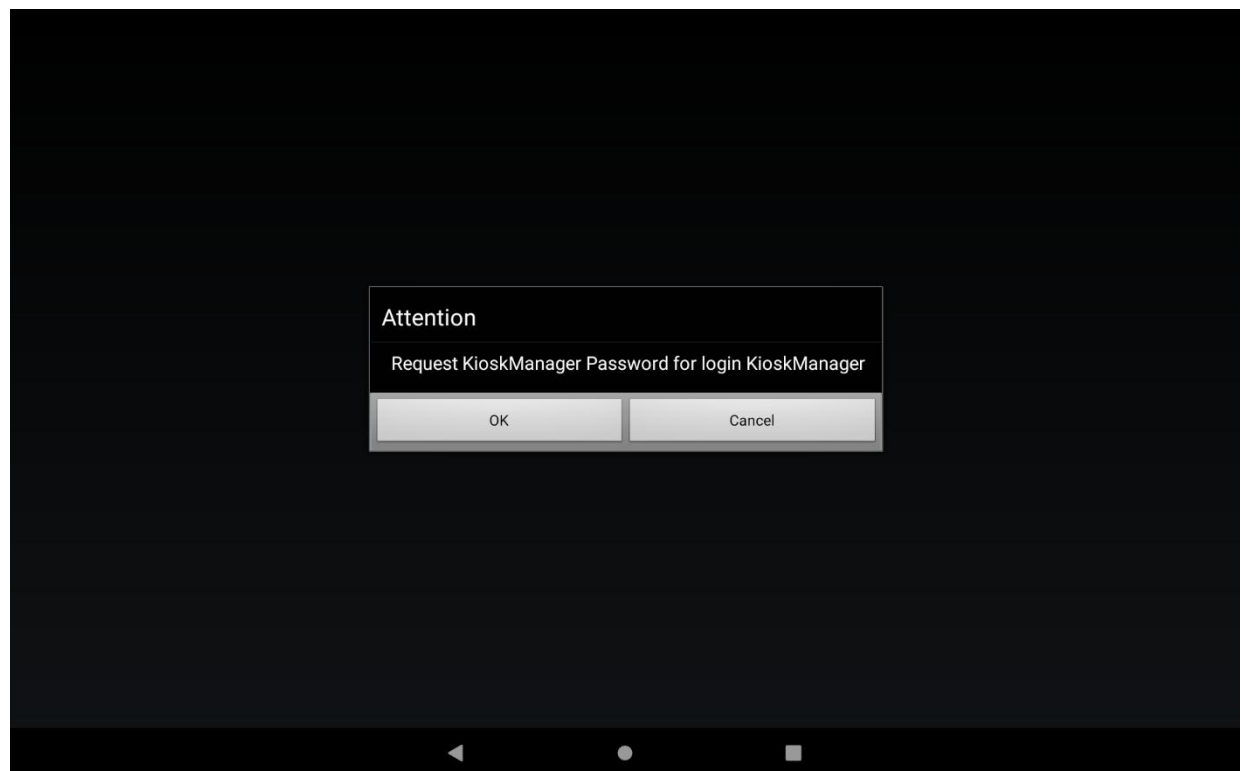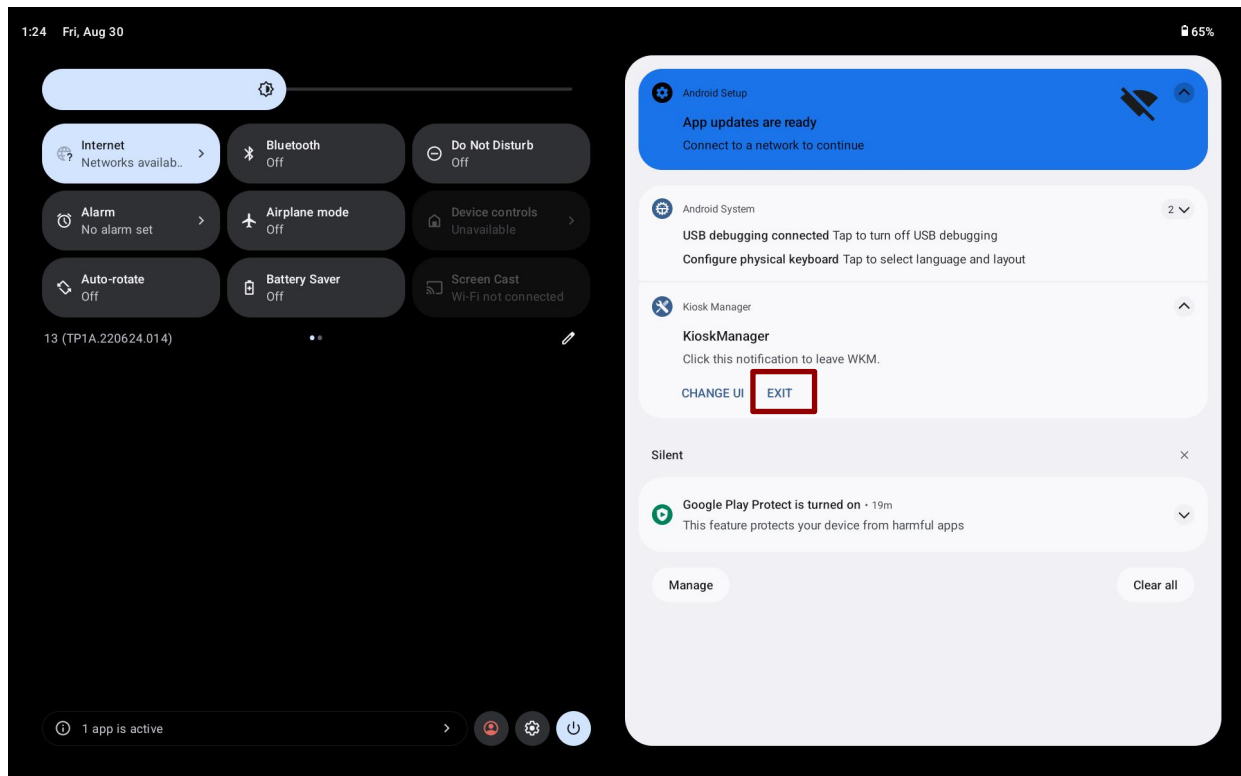
2. Click **Clear and Reset Default Home**.

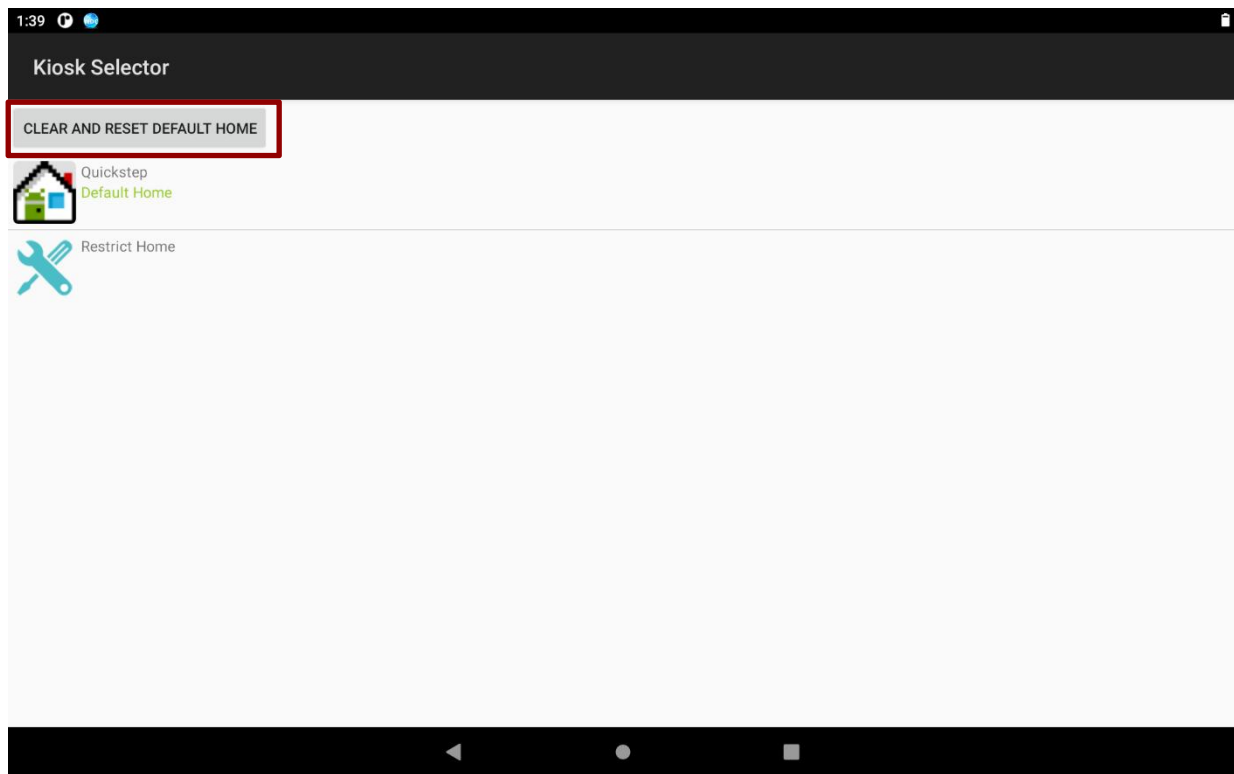3. Select **Quickstep**, and press **ALWAYS**.

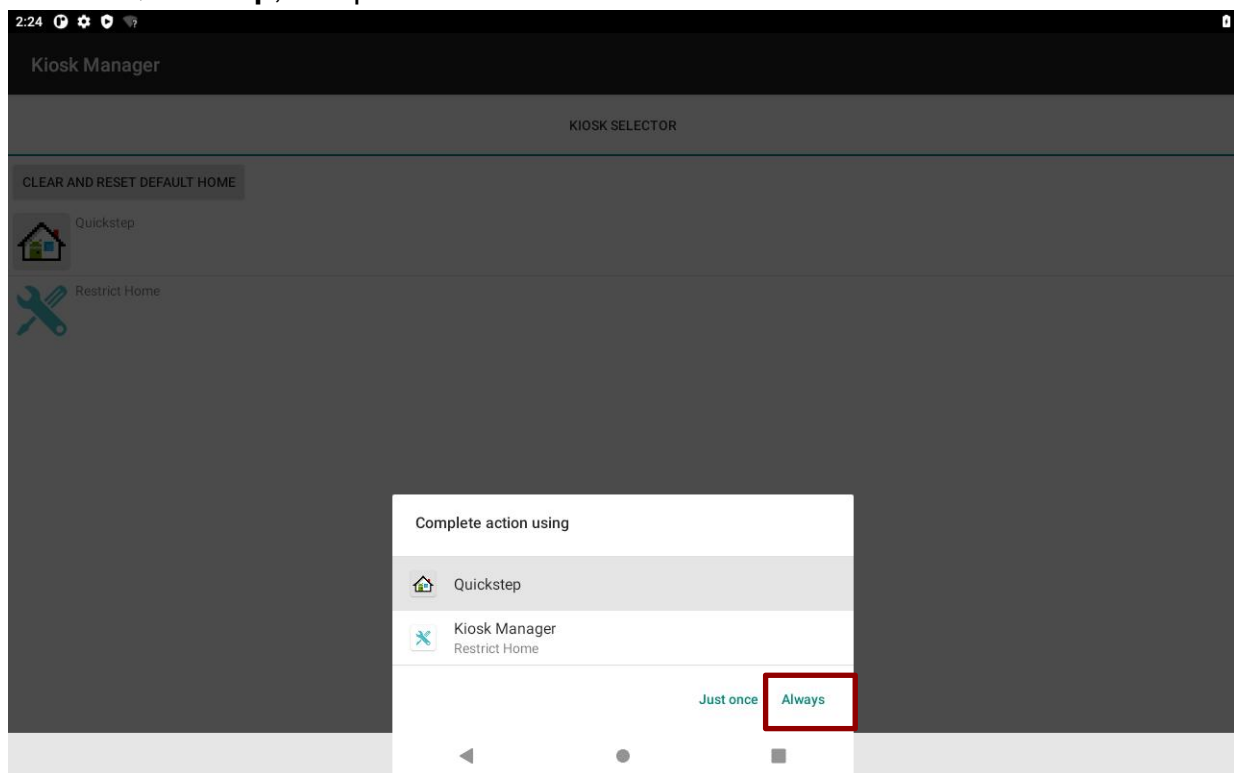
4. Restart the system and you will go back to the normal Launcher.

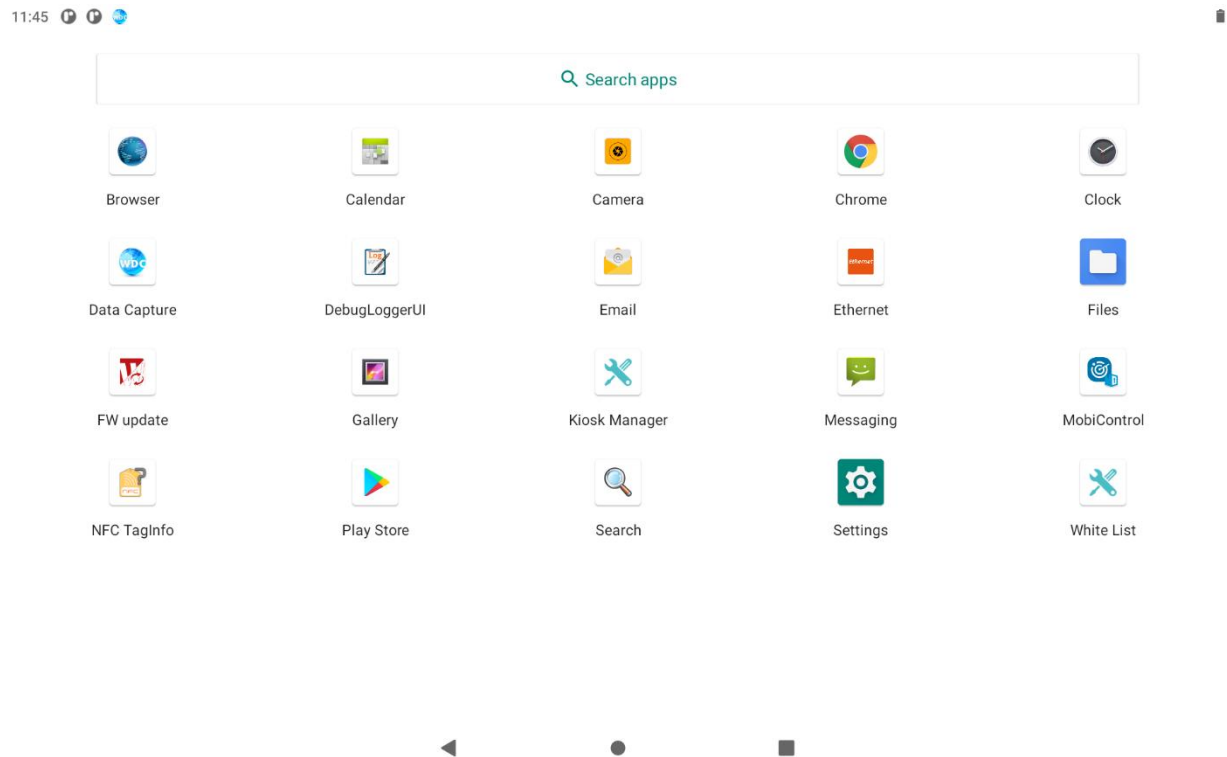
4.3 How to Change Password

1. In the main Apps window, select the **Kiosk Manager** . Select **Password Manager** to start your settings.

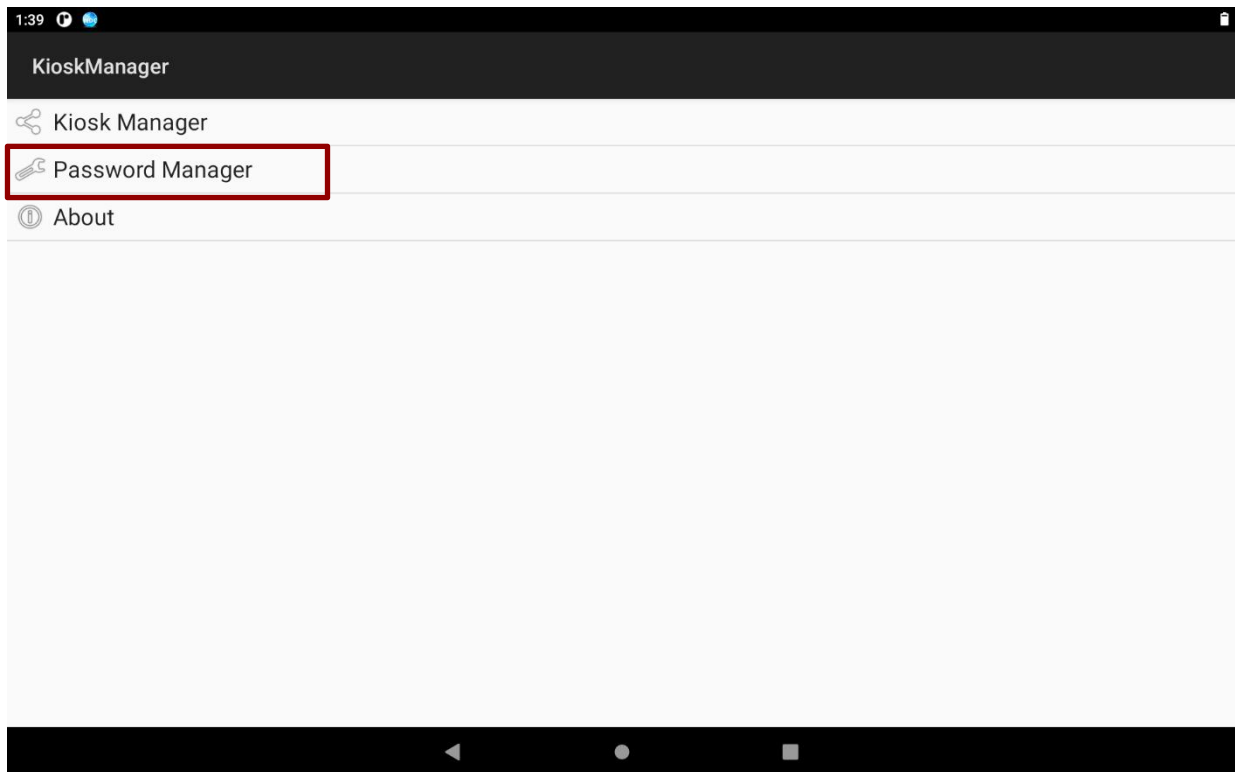
2. The default password is “**password**”. Enter a new password in **New Password** field, and then tap **Change Password**.

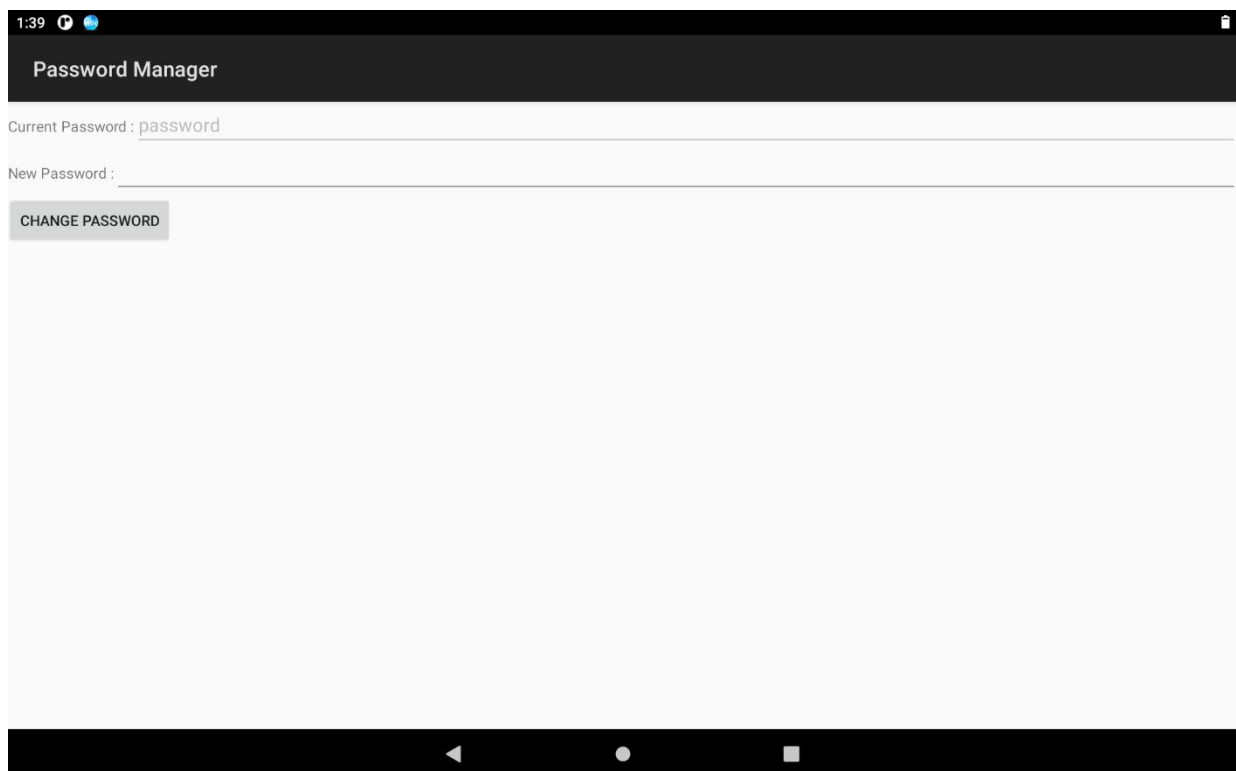
4.4 White List App


1. Open the **White List** .

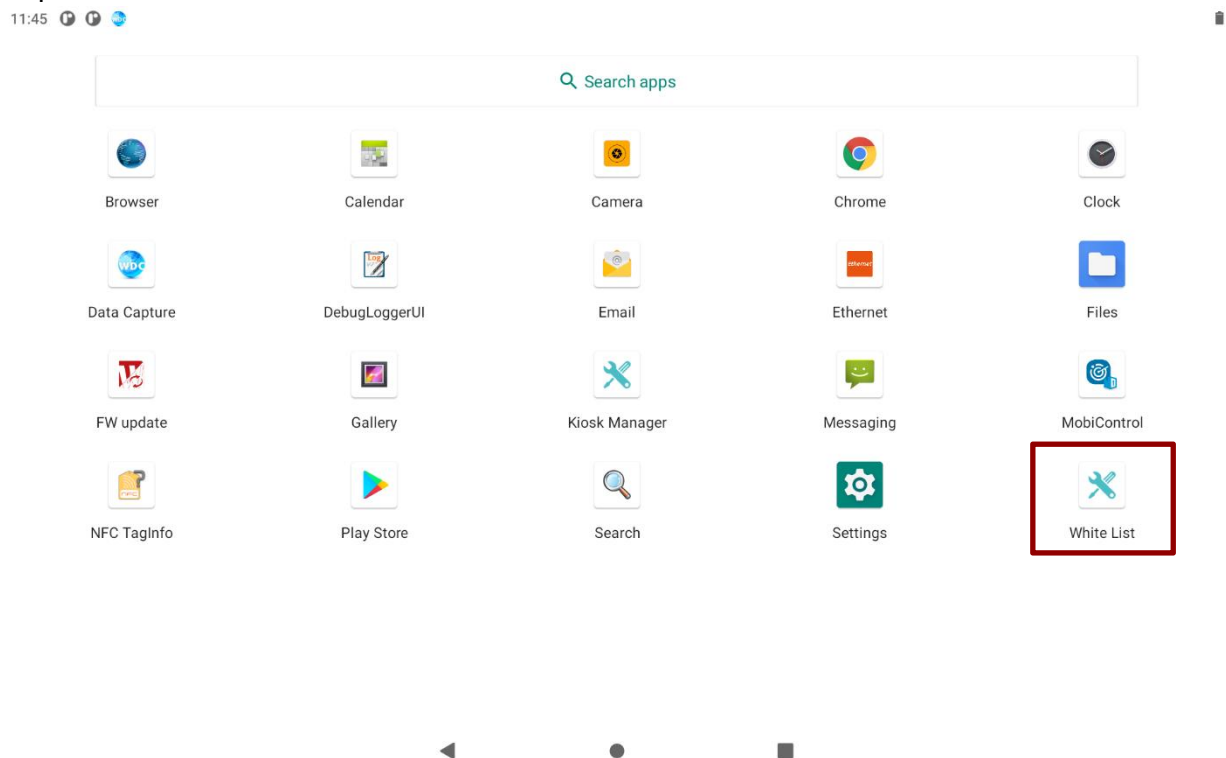
2. Select an application of **White List**.

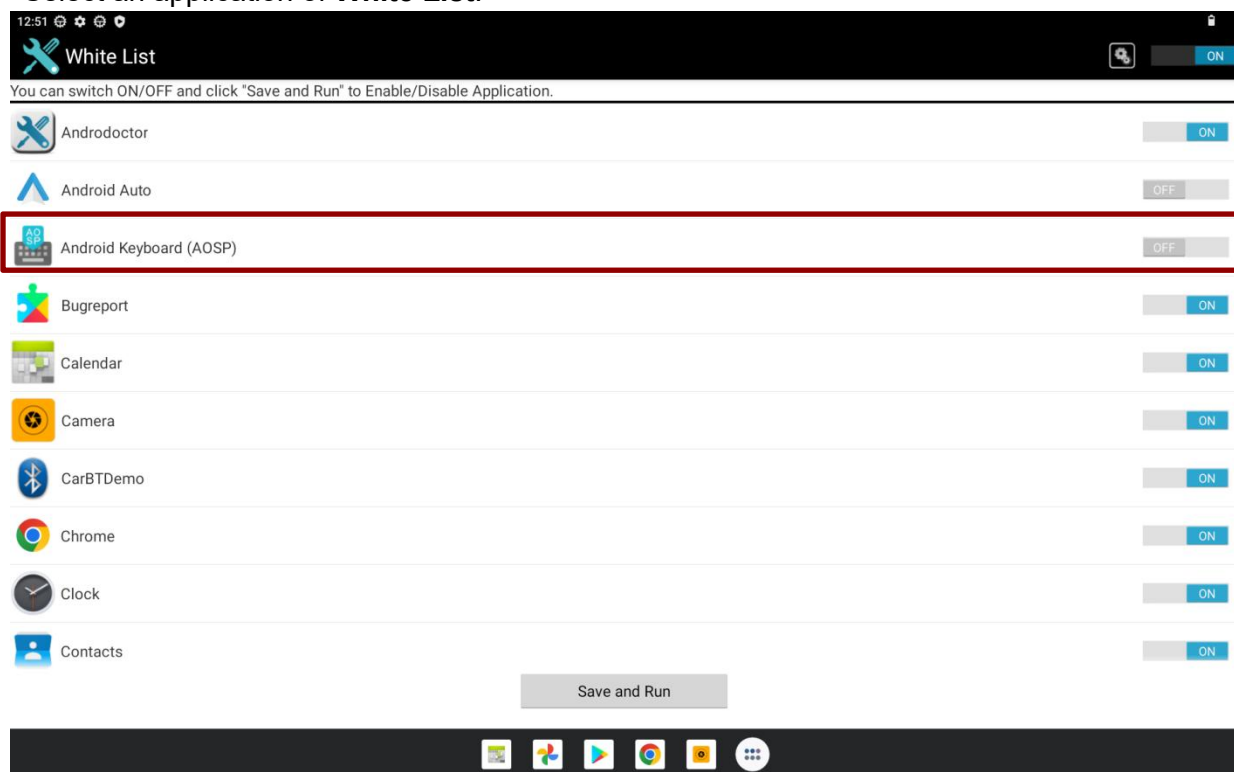3. Enter the password. The password default value is “**password**”.

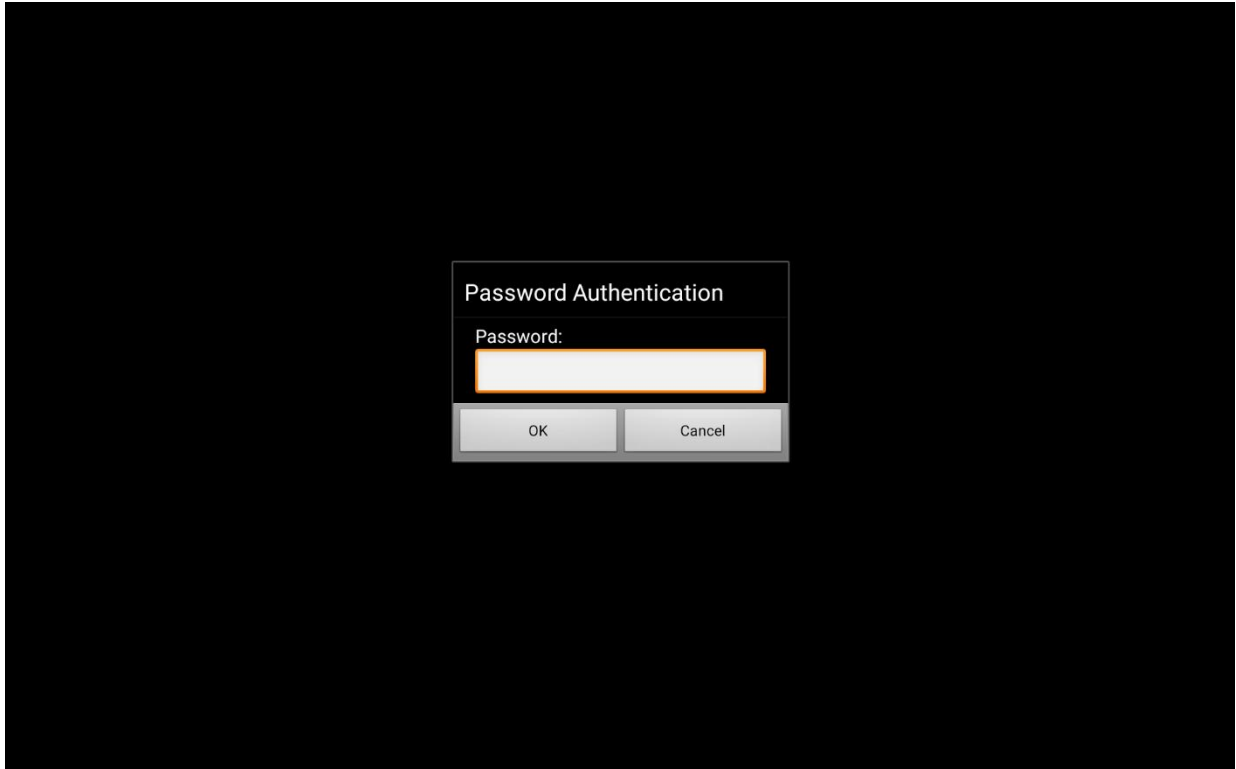
4. To change the password, click the setting icon  on the White List, choose **Password Manager**. Input the current and new password, then click **Change Password**.

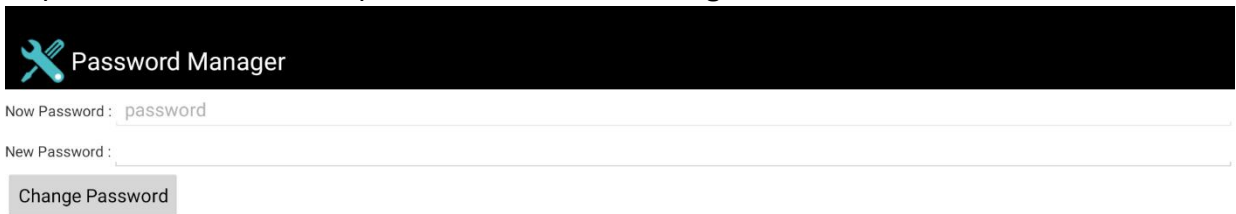
A screenshot of the 'Password Manager' interface. It features a dark header bar with a blue wrench icon and the text 'Password Manager'. Below the header, there are two text input fields. The first field is labeled 'Now Password :' and contains the text 'password'. The second field is labeled 'New Password :'. Below these fields is a button labeled 'Change Password'.

Disable or Enable Application:

1. Click list item to switch App Status

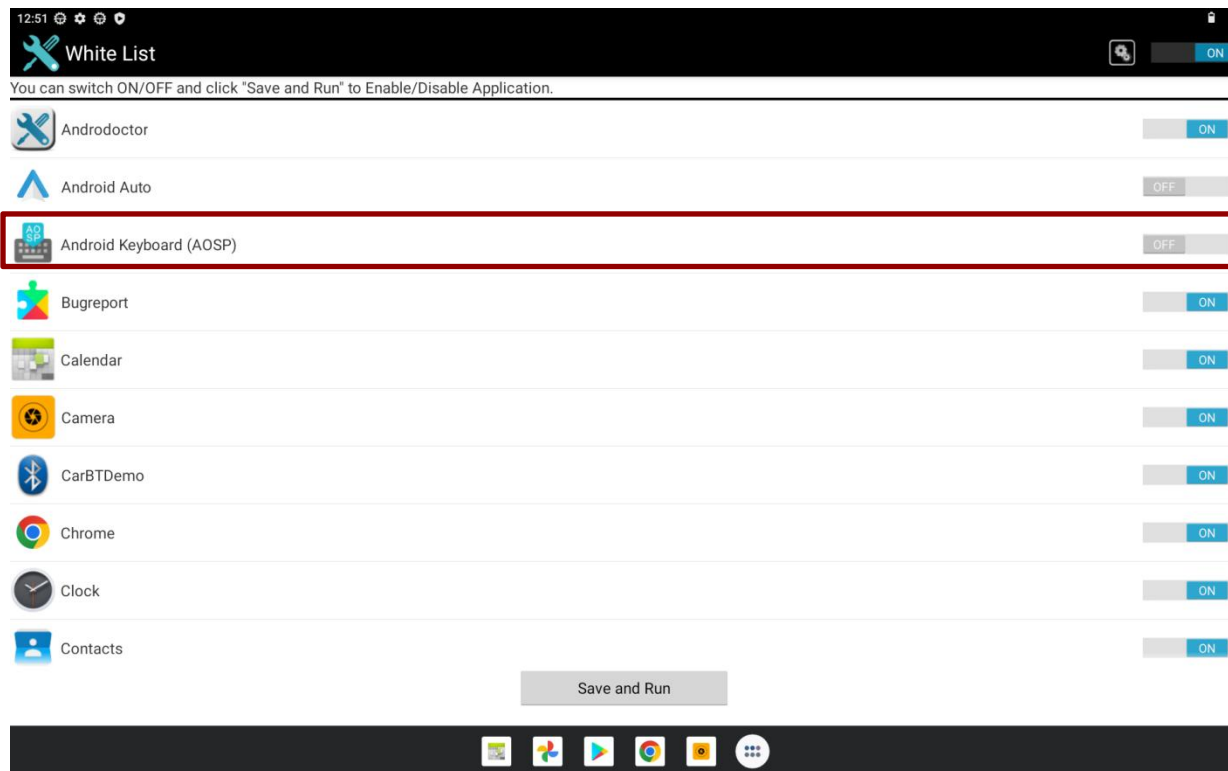2. Click **Save and Run**.

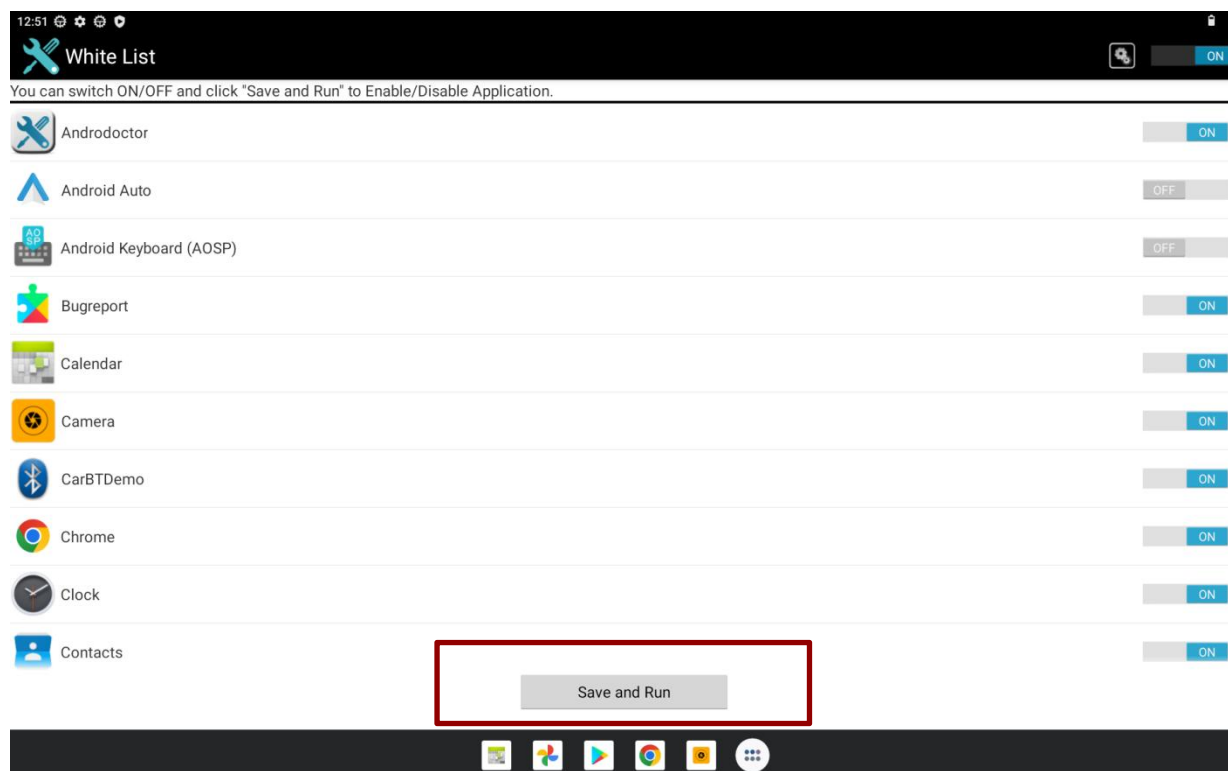
Chapter 5: Networks and Wireless Devices

This chapter explains how to connect the MD10W2 Series Rugged Tablet to network and peripheral devices.


5.1 Connecting to WLAN Networks

5.1.1 Turning on WLAN and Connecting to WLAN Networks

1. Swipe from the bottom of the main desktop, tap **Settings** .
2. Tap **Network & Internet**  and then toggle the on/off switch on the top-right hand side to turn on the **WLAN**.

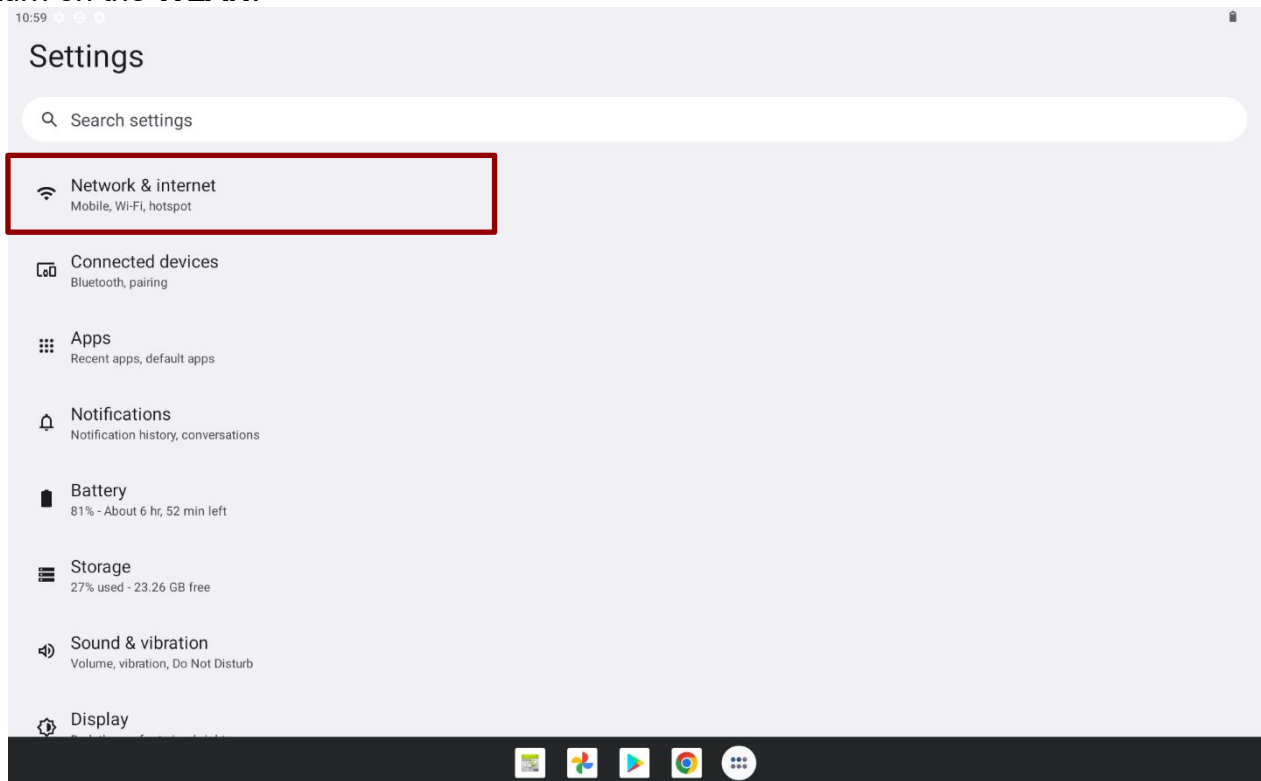
3. Tap **Internet** . The device scans for available WLAN networks and displays the list of it. Secured networks are indicated with a lock icon. If the tablet finds a network that you have connected previously, it connects to it automatically.

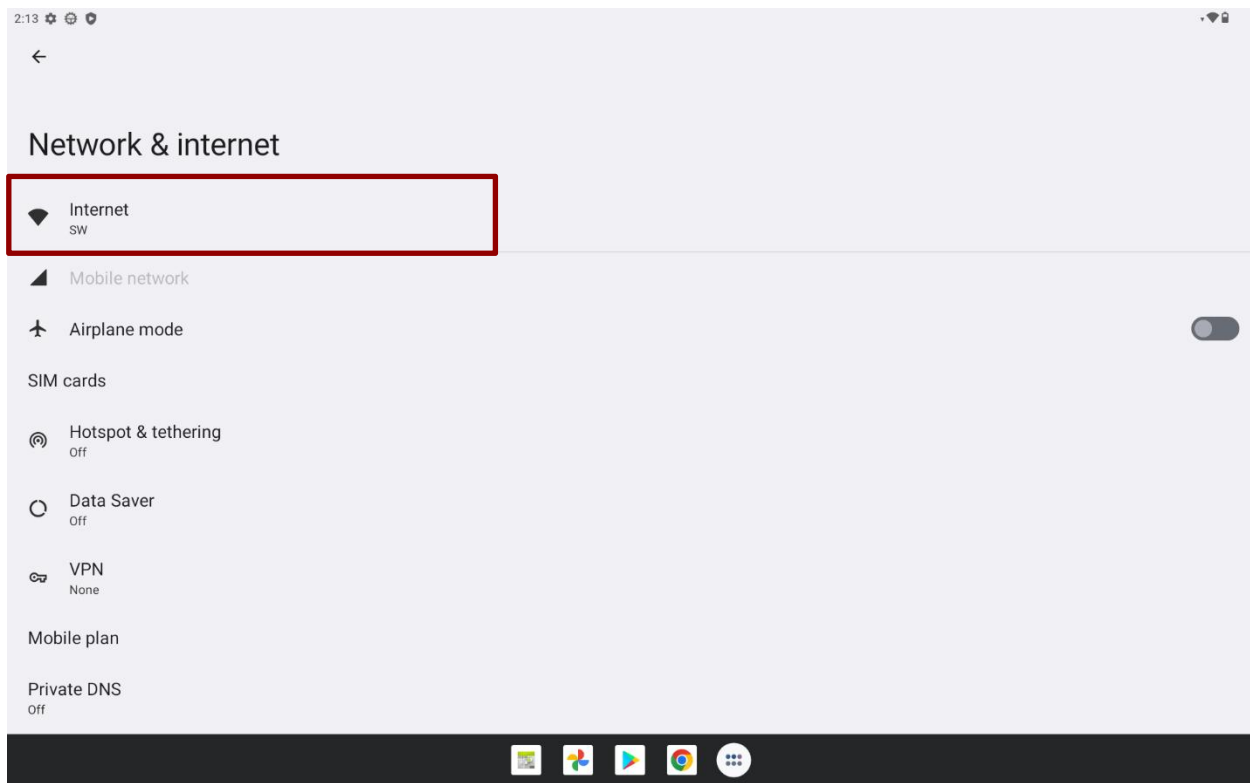
4. Tap a network to connect to it:

- If the network is open, you are prompted to confirm that you want to connect it by tapping **Connect**.
- If the network is secured, users are prompted to enter a password or other credentials.



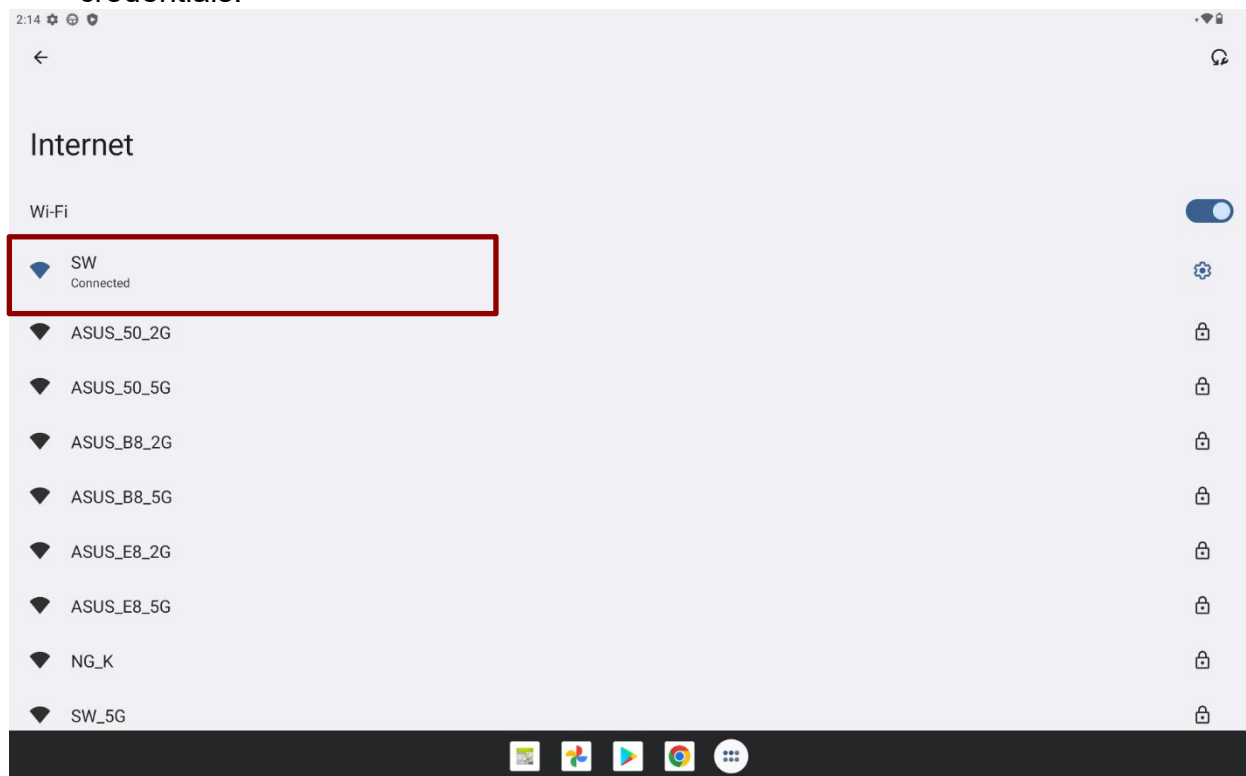
5. To check the details about the speed, security, address and related settings, tap its name in the WLAN settings screen.

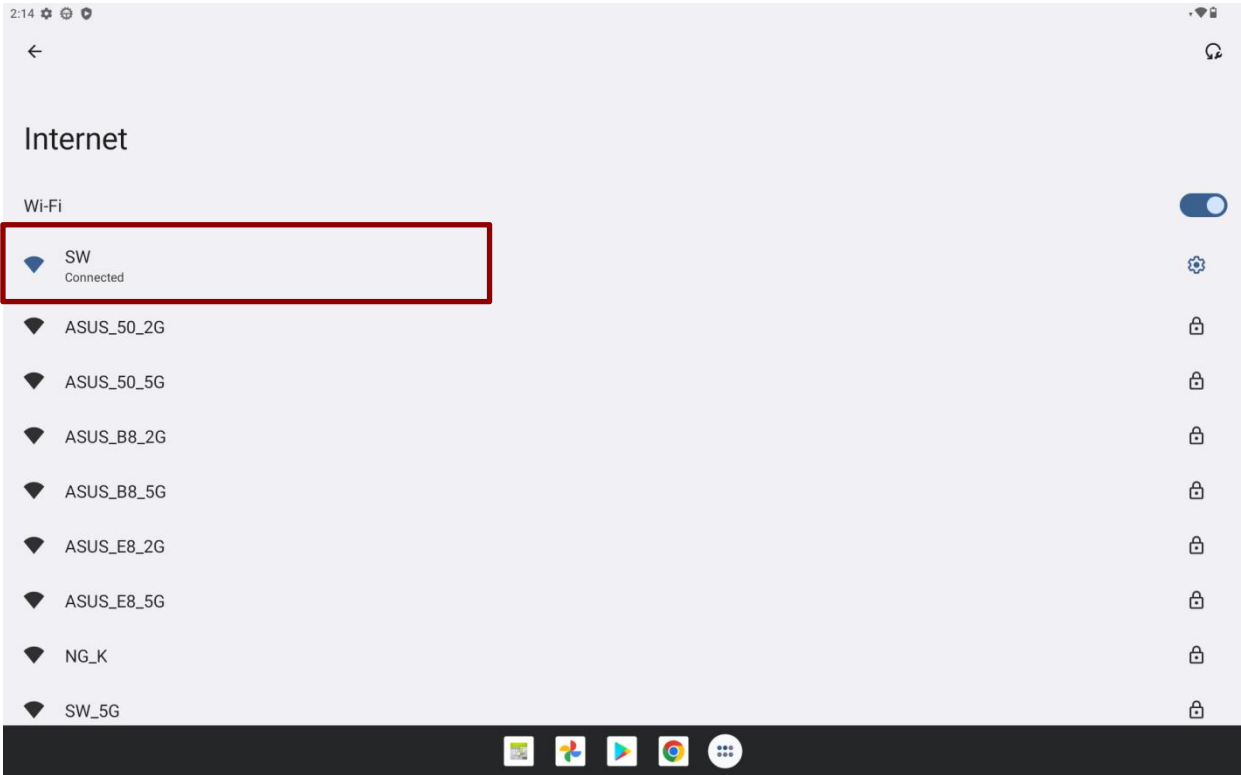
5.1.2 Adding WLAN Network

1. Turn on WLAN, in the WLAN settings screen, tap **+ Add network** (at the bottom of the list of discovered networks).

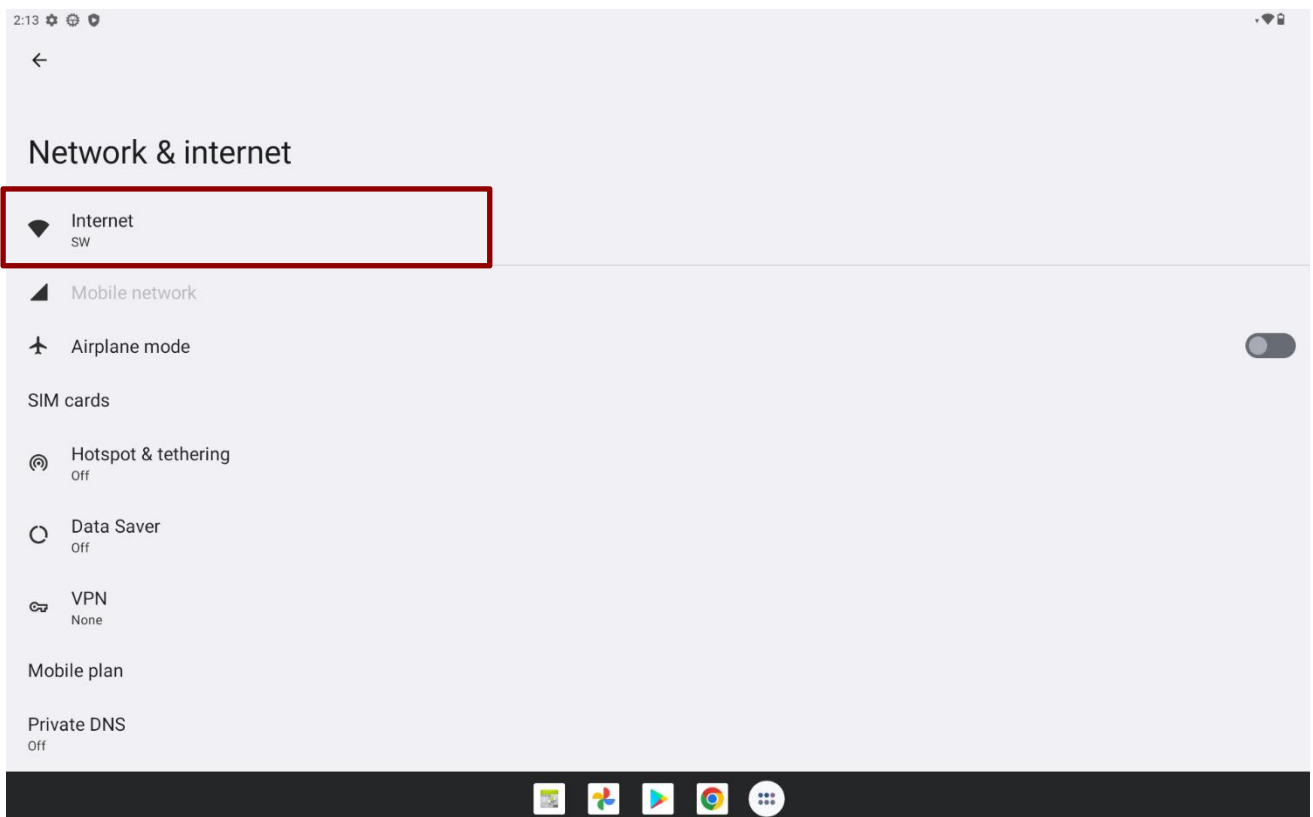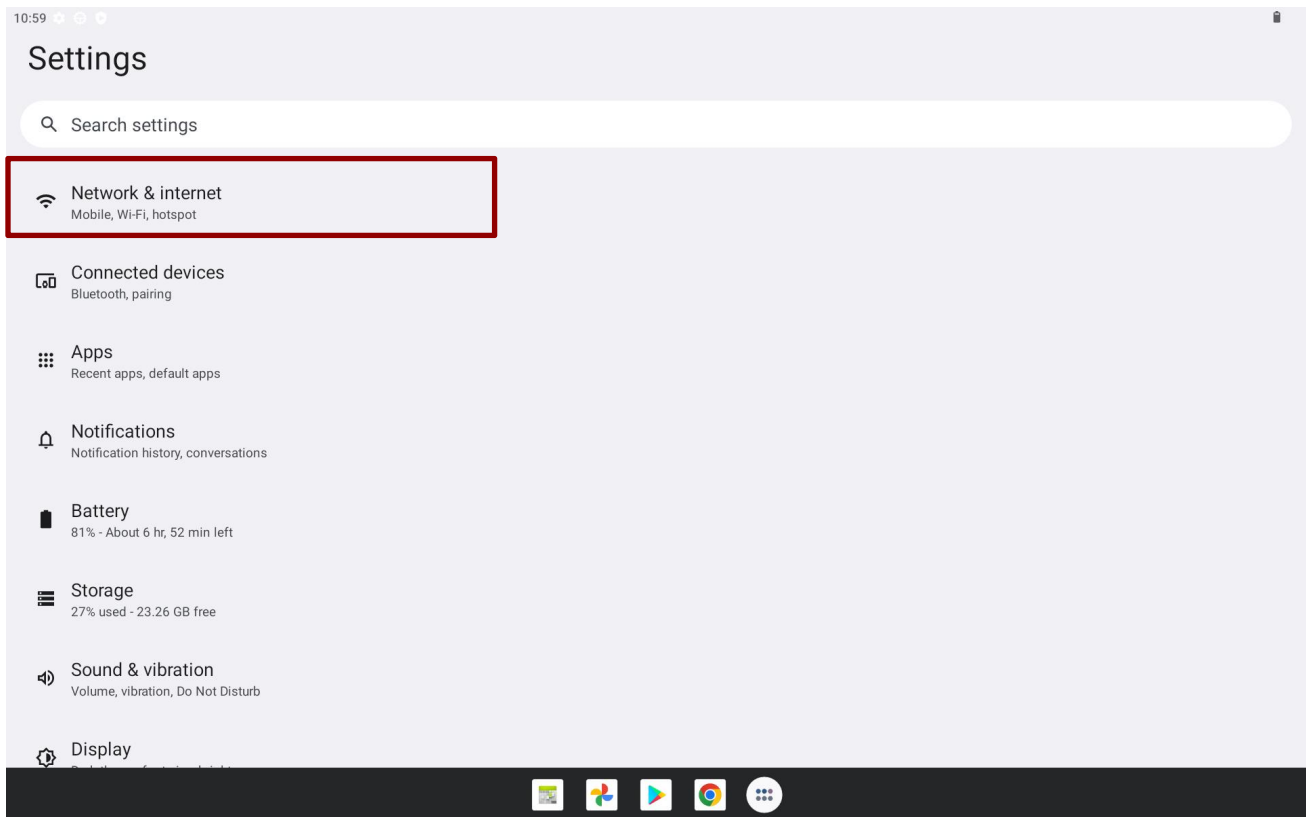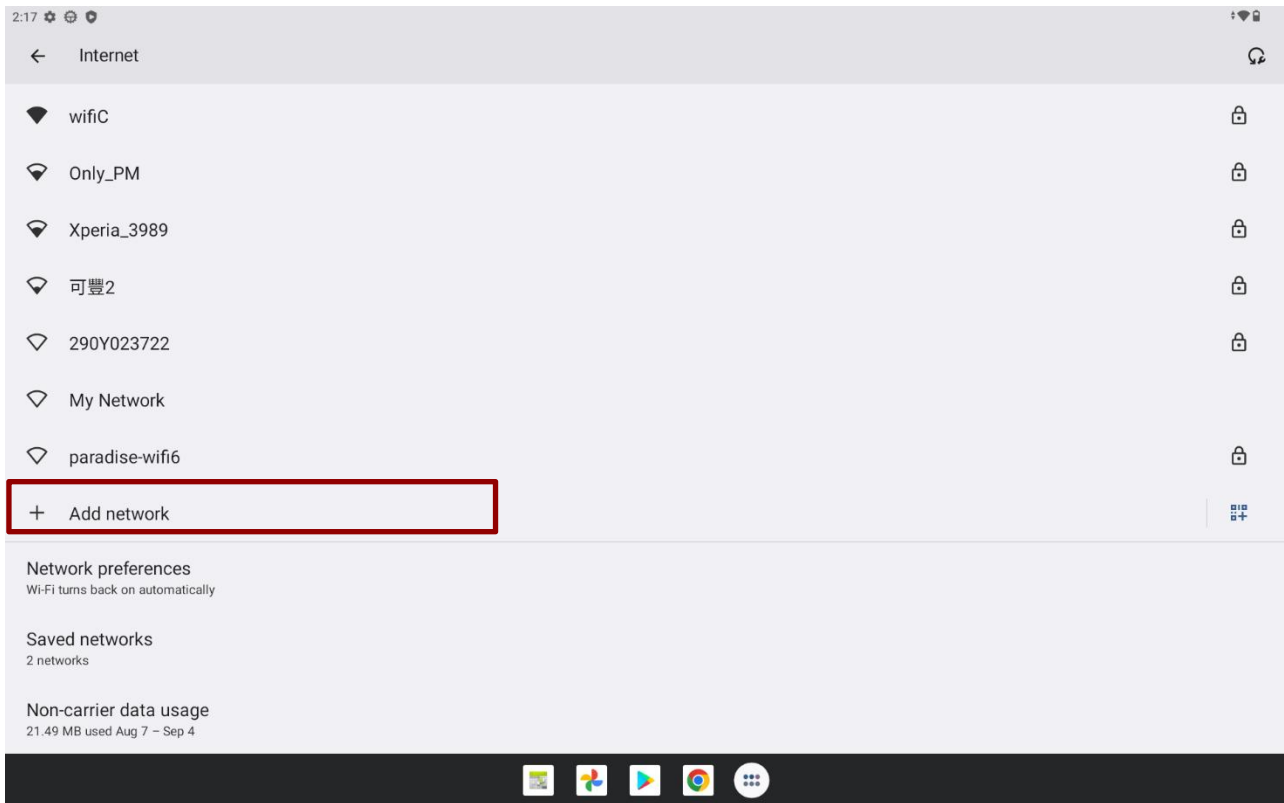2. Input the SSID (name) of the network, if the network is secured, tap the security menu and then tap the type of security deployed on the network.

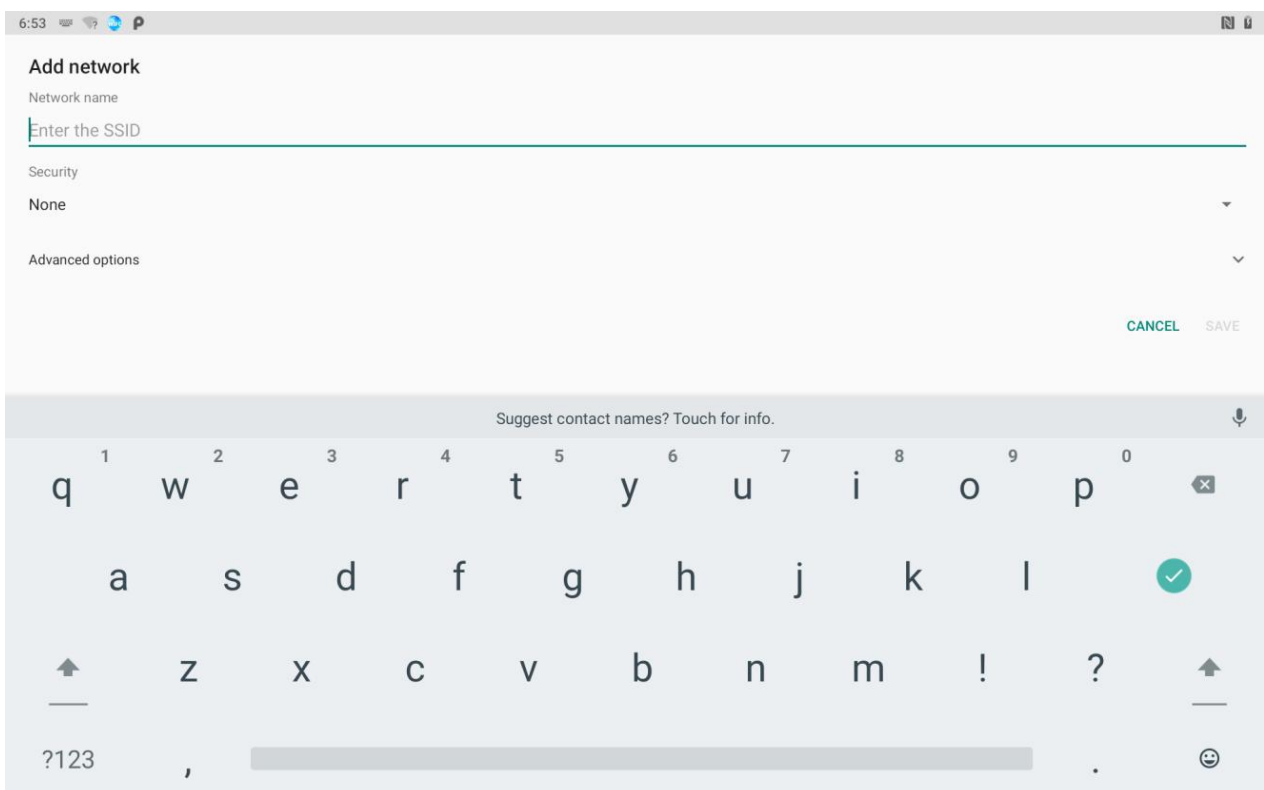
3. Enter any required security credentials.
4. Tap **Save**.

5.2 Connecting to BT Devices

5.2.1 Turning BT On or Off


1. Swipe from the bottom of the main desktop, tap **Settings** .
2. Tap **Connected Devices** .

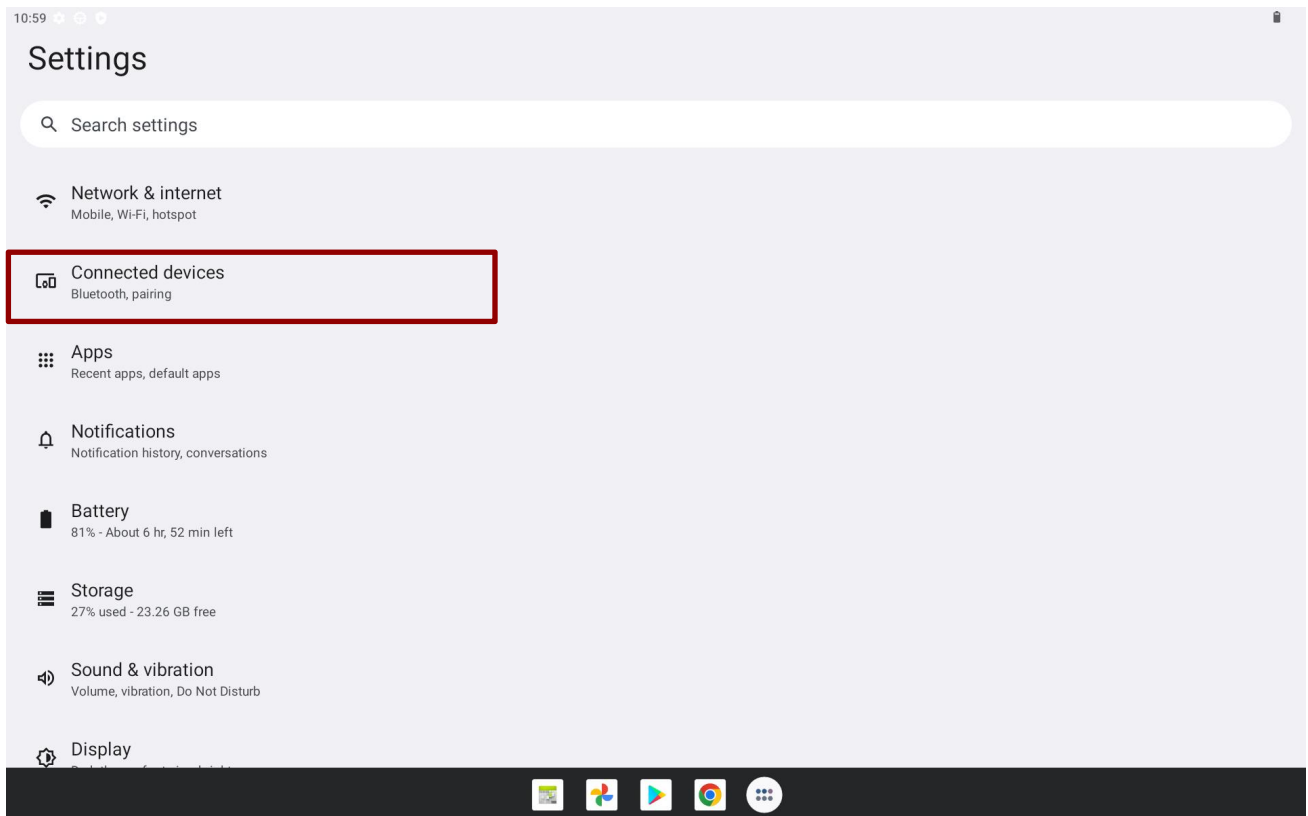
3. Toggle the on/off switch on the top-right hand side to turn on the BT.

5.2.2 Connecting the BT Device

1. Swipe from the bottom of the main desktop, tap **Settings** ⚙️.
2. Tap **Connected Devices** 📶.
3. Tap **Pair new device**. The system will automatically scan and display the IDs of all available BT device in range.

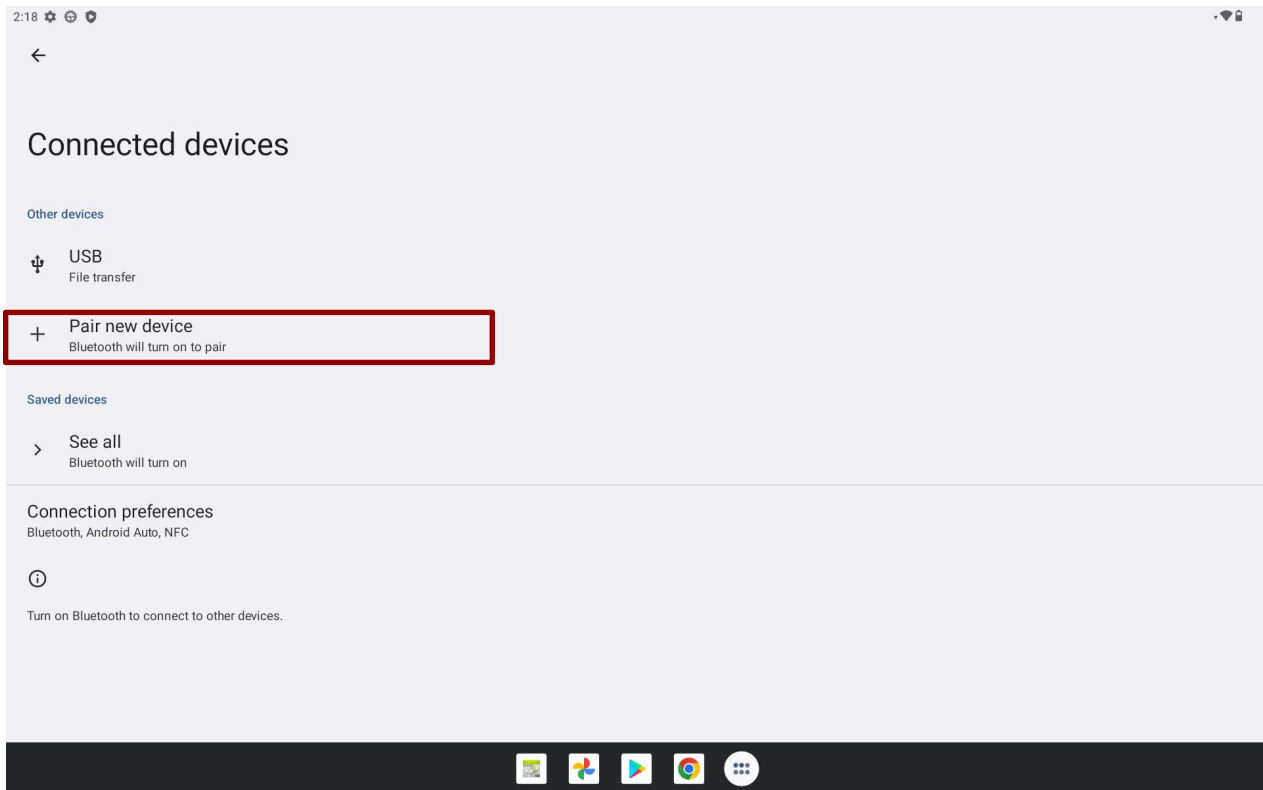
4. Select the device to be paired.

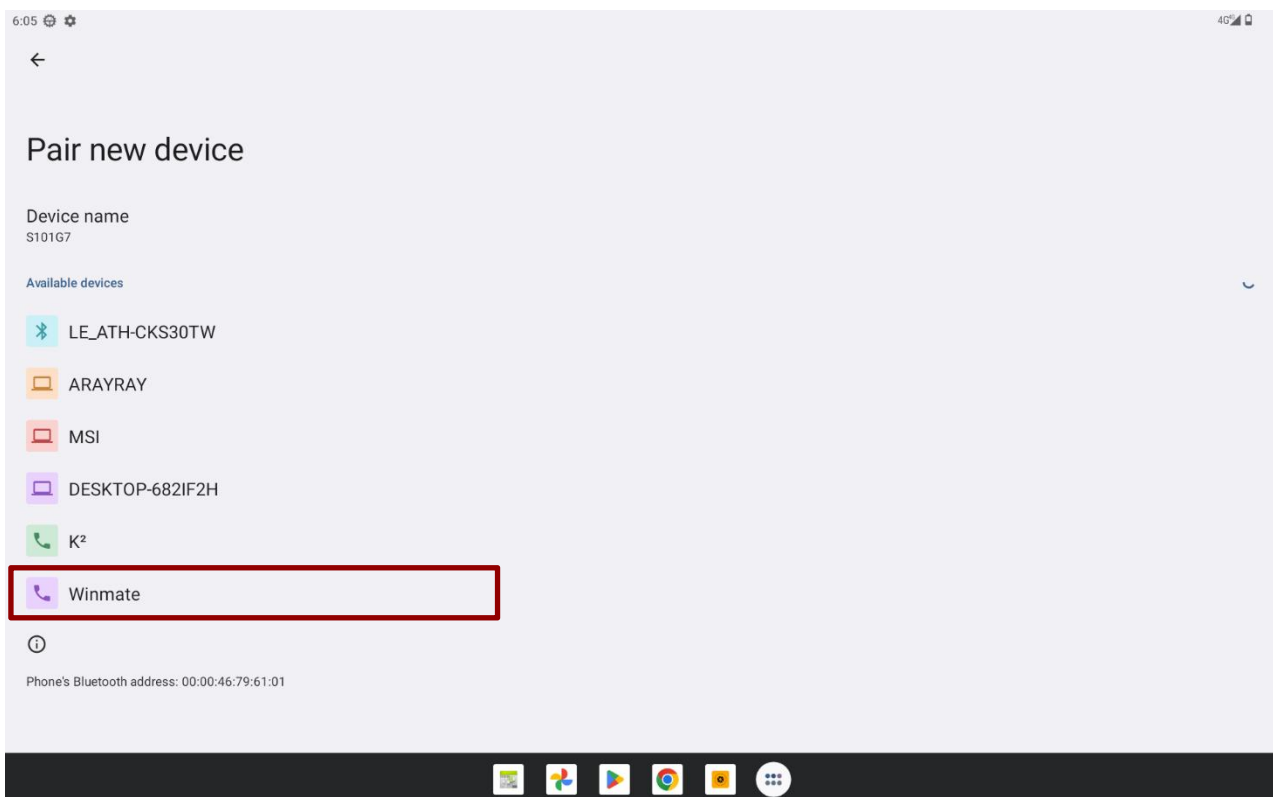
5.2.2 Disconnecting BT Device

1. Tap the device in the **BT** settings screen.

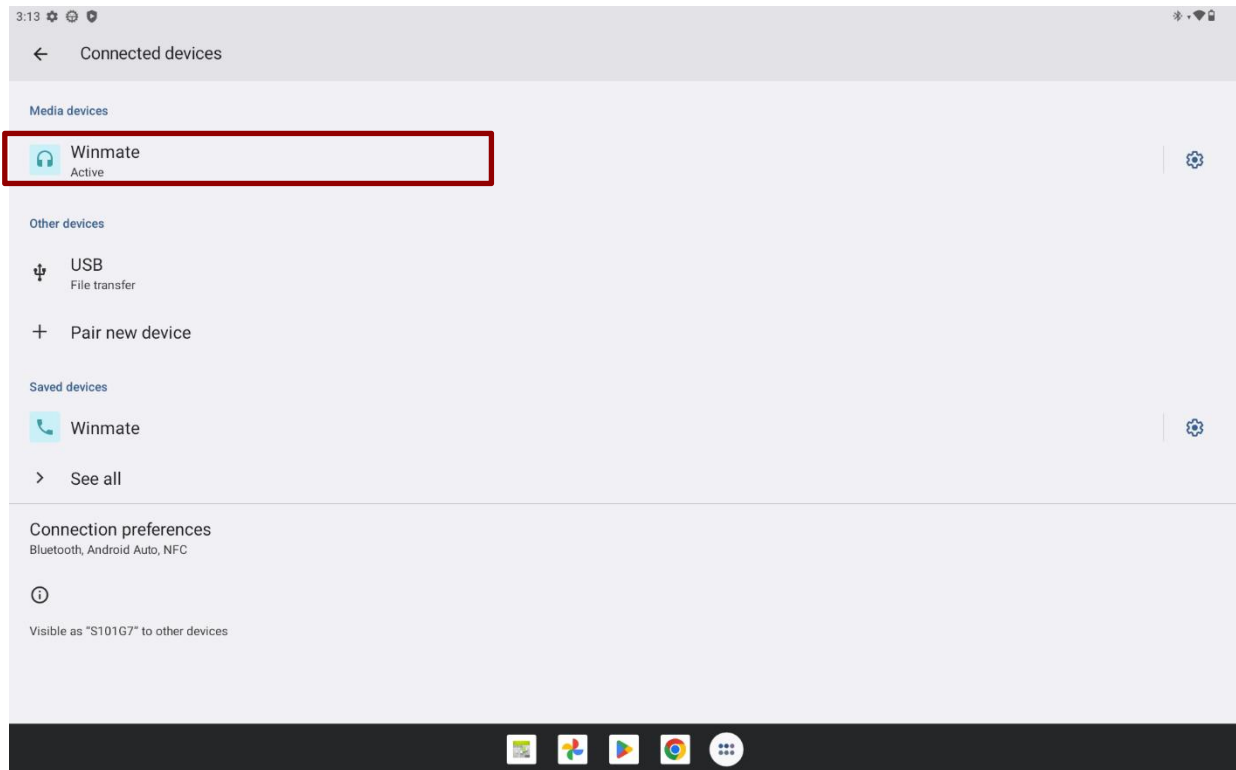2. Tap **FORGET** to disconnect.

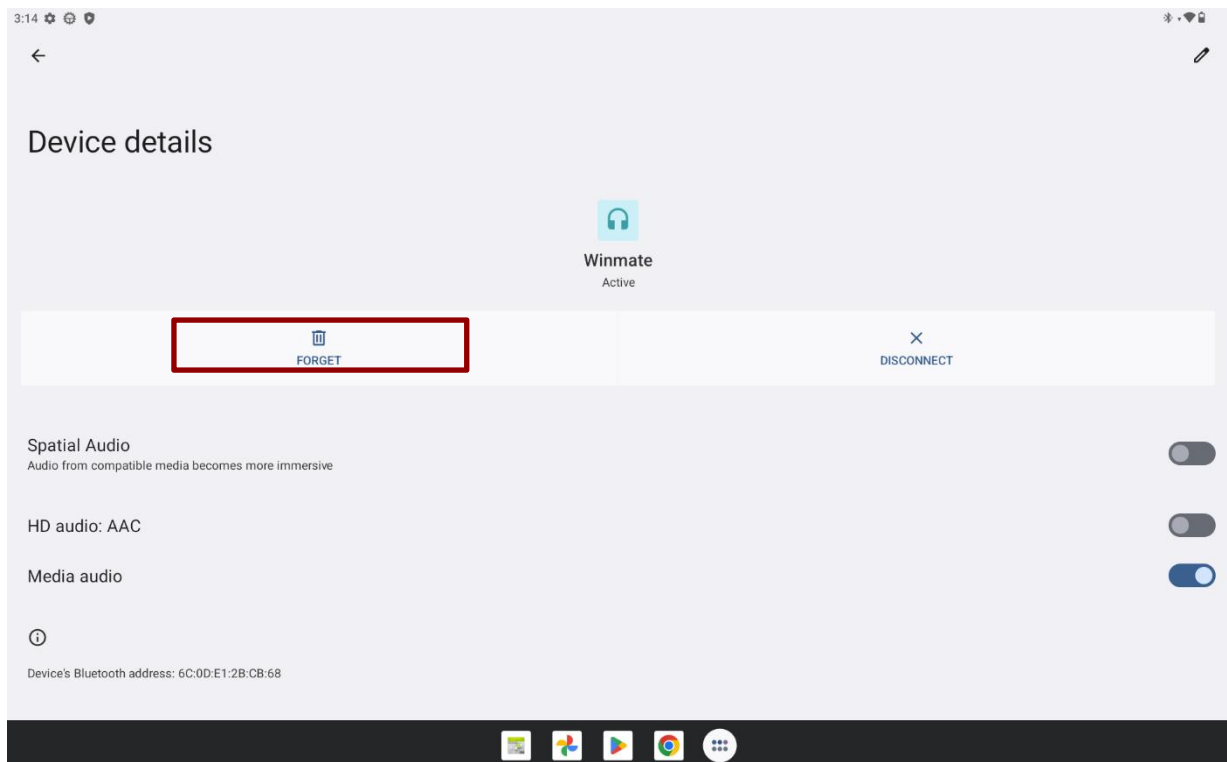5.3 Connecting to Computer via USB

User can connect their tablet to a computer with USB cable to transfer files between tablet and computer.


To allow device connection to PC: Click **Setting** > **Accessibility** > enable **Connect to PC**.


Important: When connecting your tablet to a computer and inserting its USB storage or SD card, user must follow the computer's instructions for connecting and disconnecting USB device, to avoid damaging or corrupting files. Notice that the MD10W2 Series Tablet PC cannot be charged via USB.

5.3.1 USB Computer Connection

1. Connect tablet to the computer via USB connector.
2. Pull down the notification area to open notification panel, tap **Android System**.

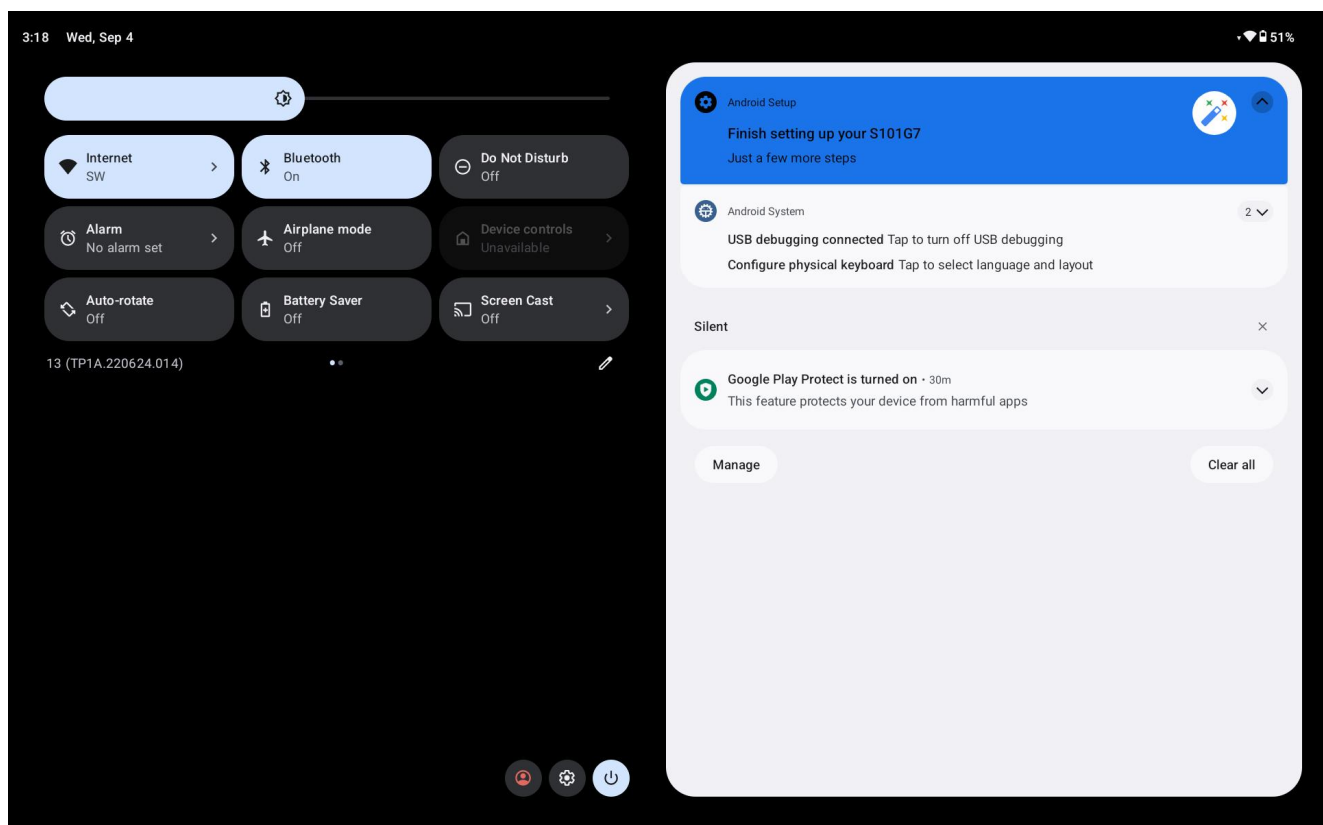
3. In the **USB Preferences** choose the type of USB connection.
 - **File Transfer** – let user transfers files between computer and tablet.
 - **PTP** – let users transfer photos or media files if MTP is not supported.
 - **MIDI** – let users use device as MIDI.
 - **No data transfer**– disable data transfer.


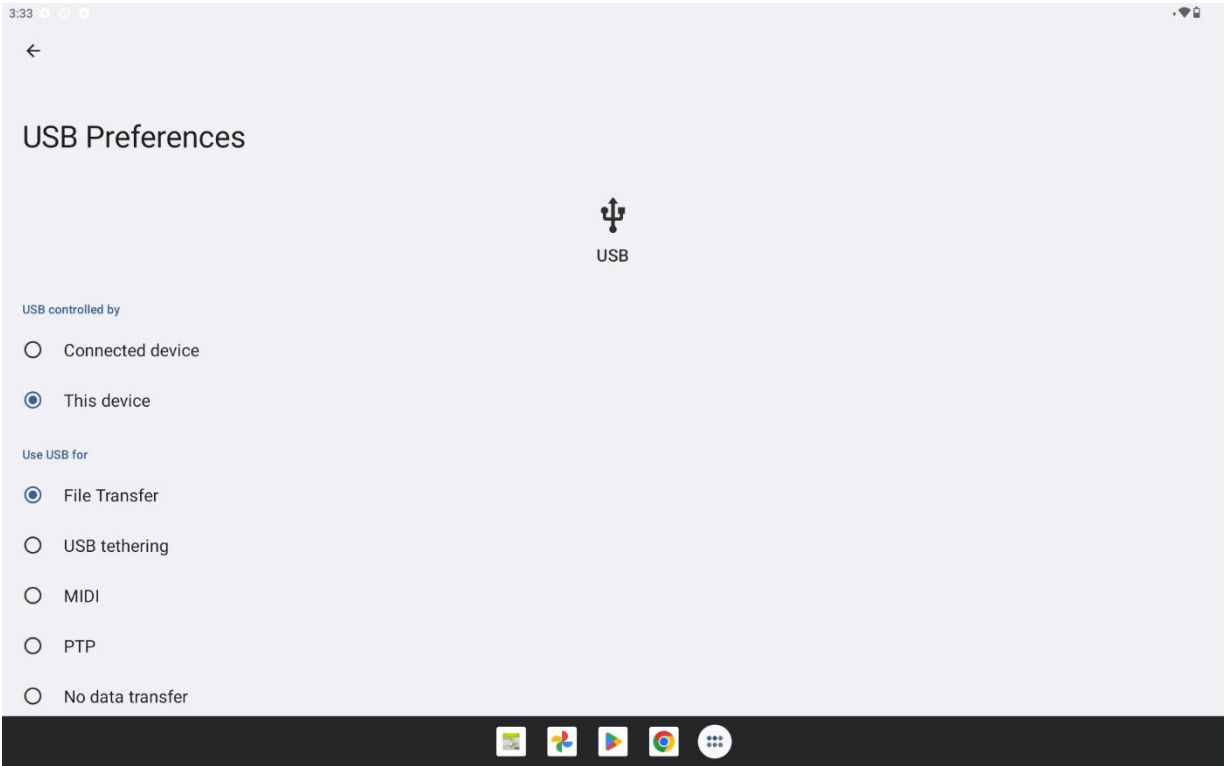
5.3.2 USB Debugging Settings


To permit debugging tools on a computer to communicate with User’s tablet via USB Connection.

1. Connect tablet to the computer via USB cable.
2. Swipe from the bottom of the main desktop, tap **Settings**, tap **About phone**.

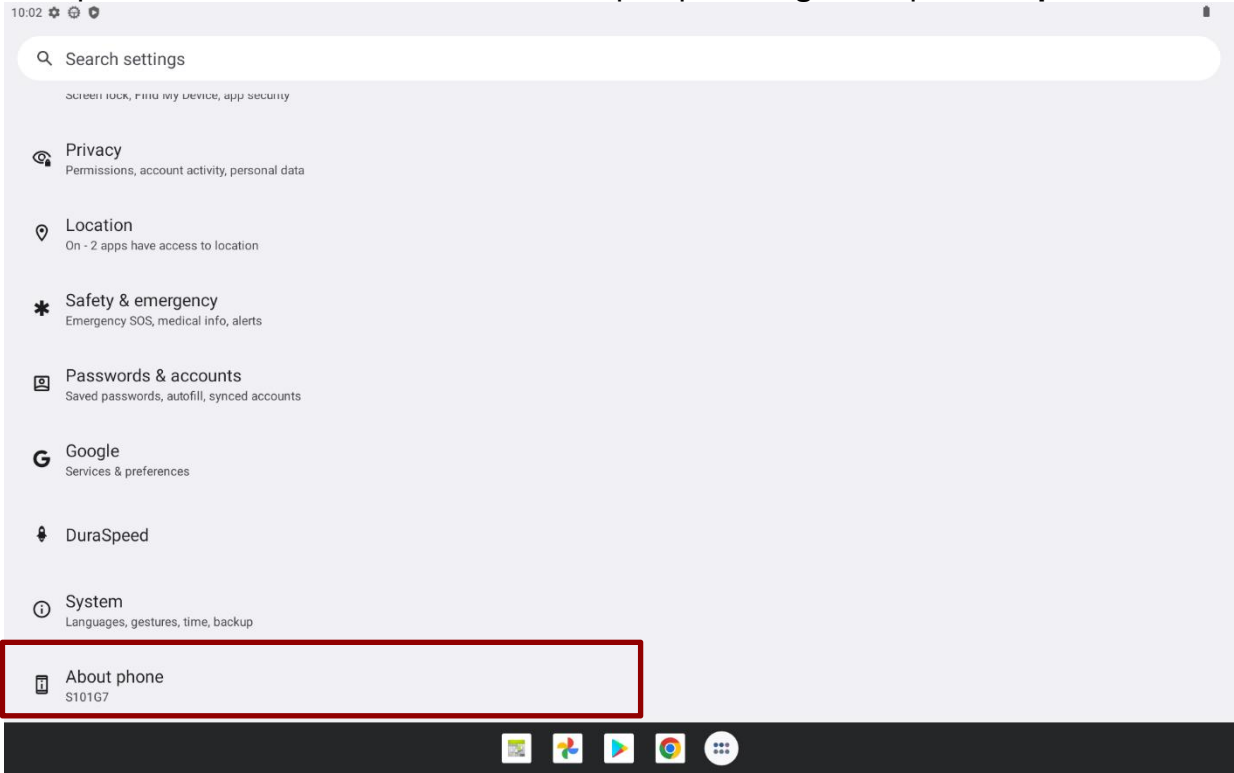
3. Tap **Build number**.

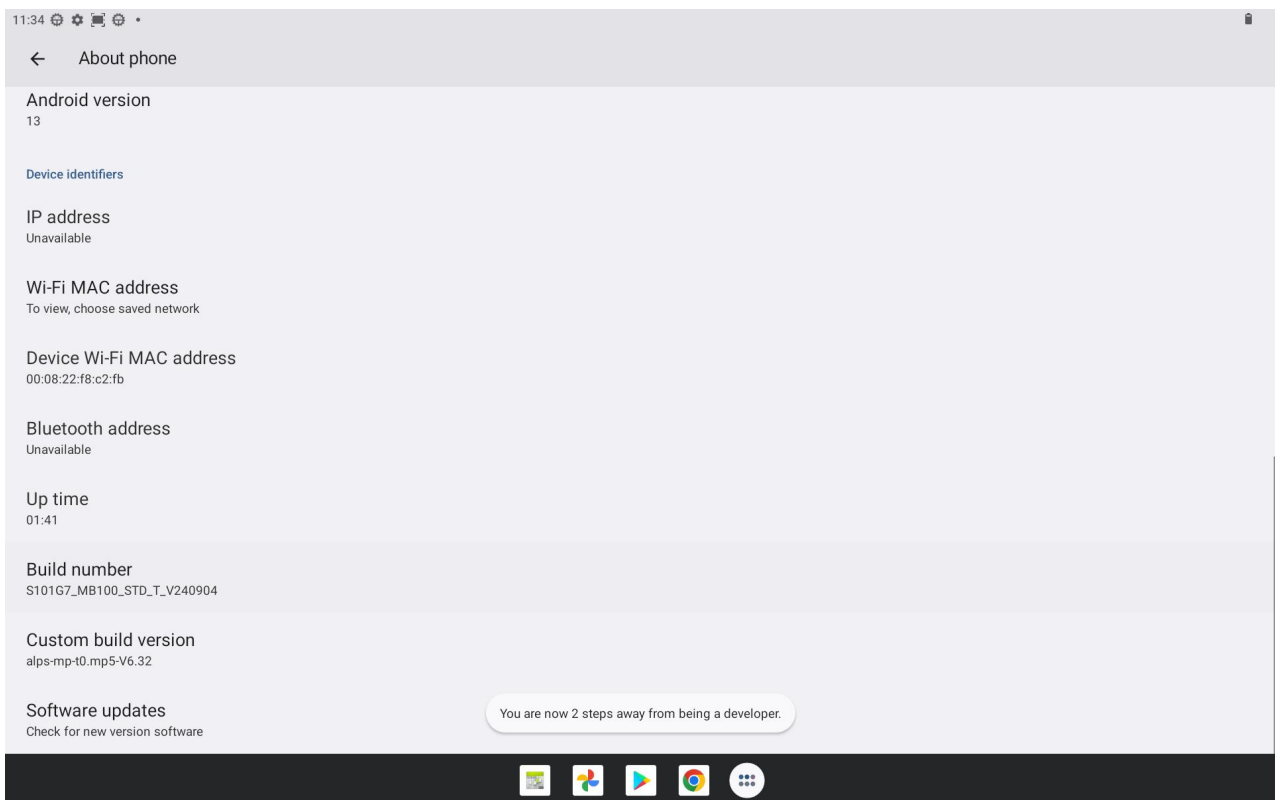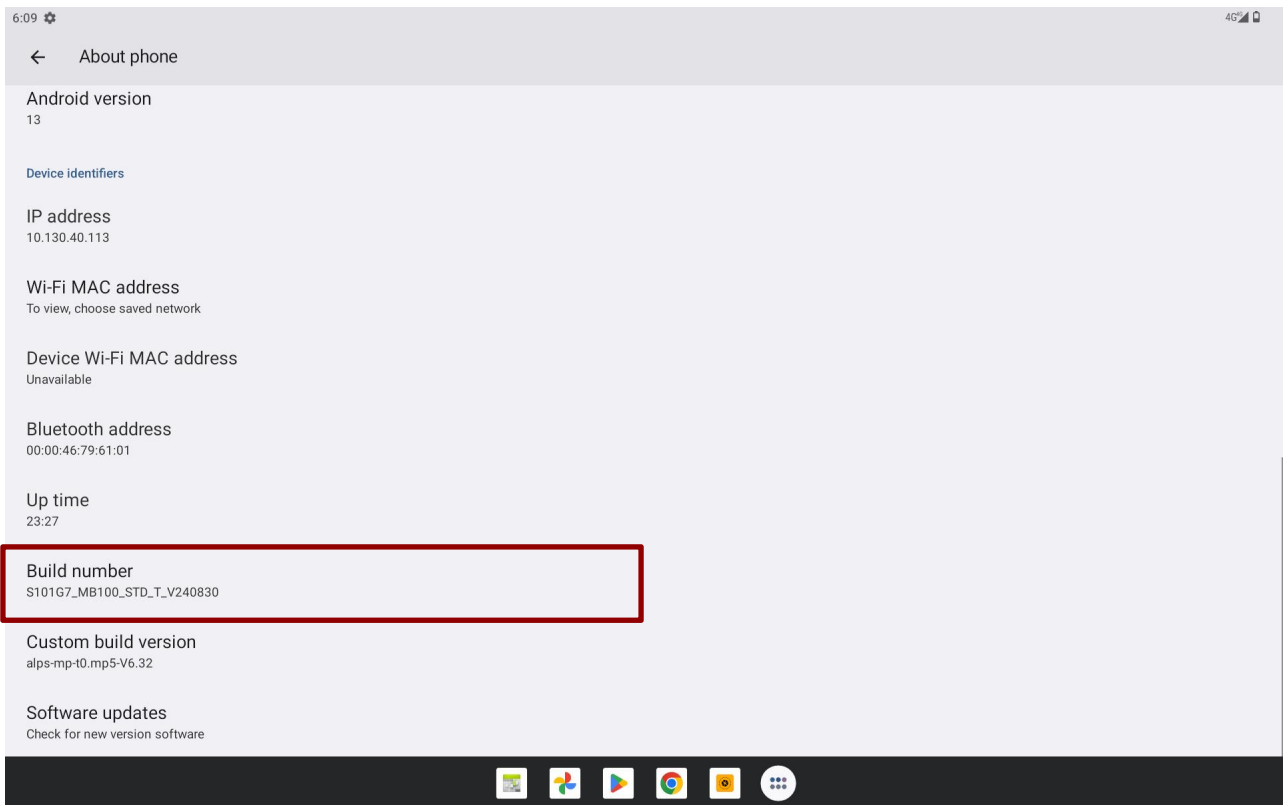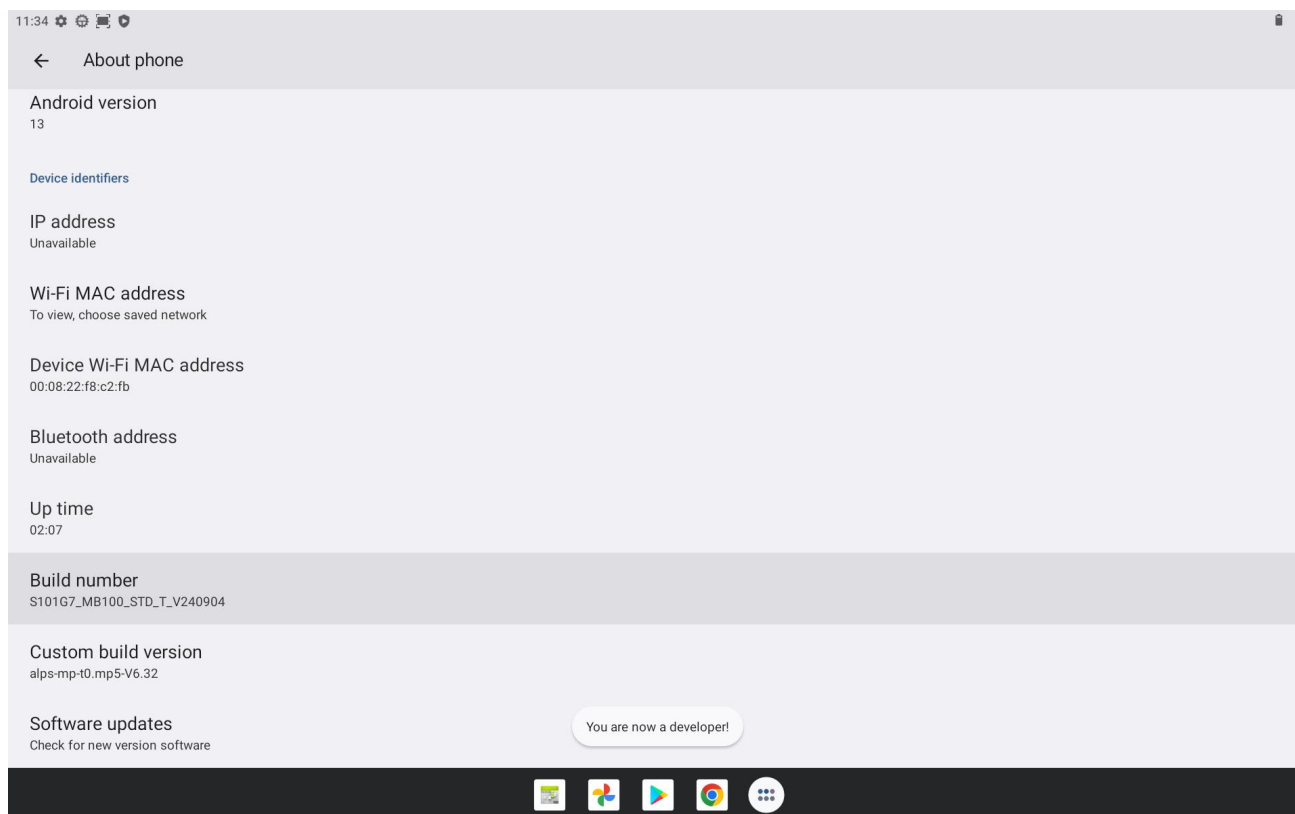
4. Pull down the screen to open notification panel. Tap **USB debugging connected**.

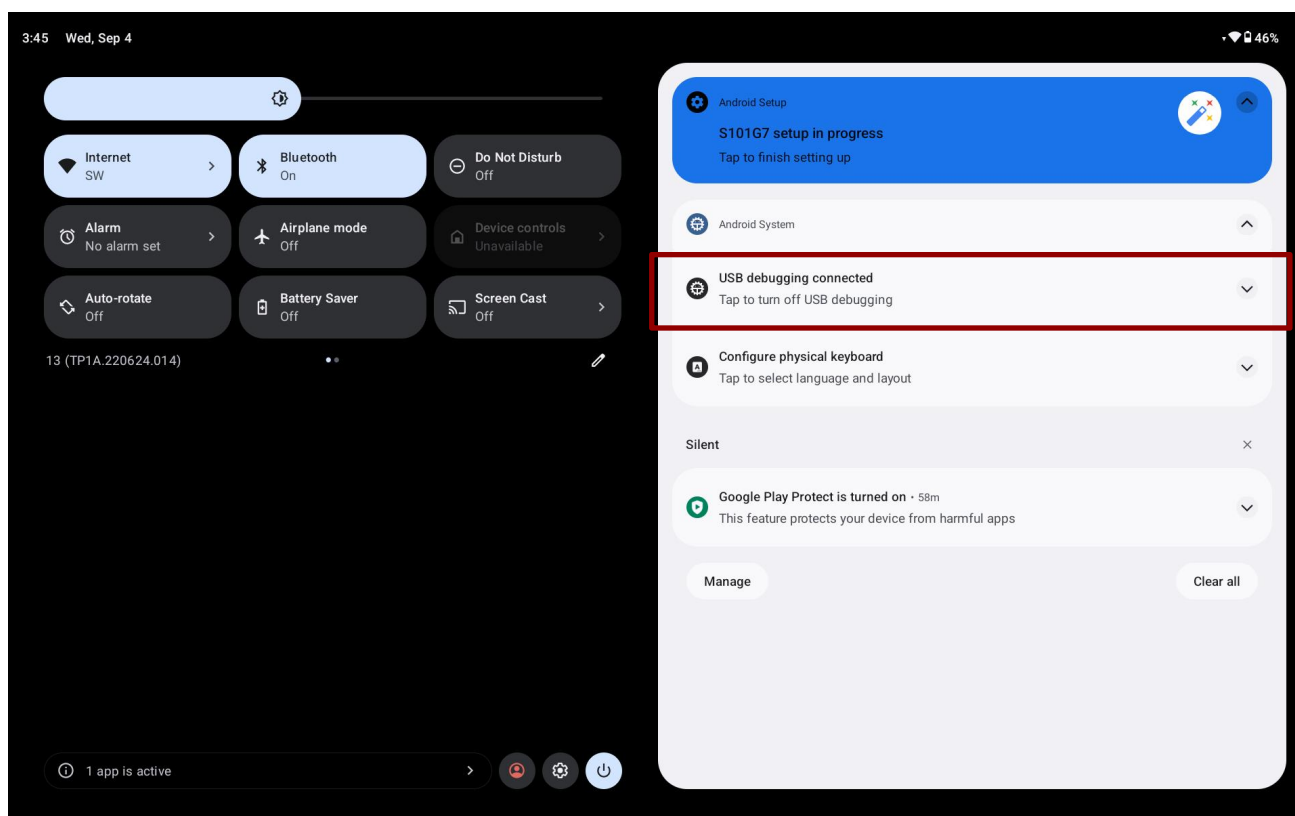5. In the **Developer options** user can choose the action to perform.

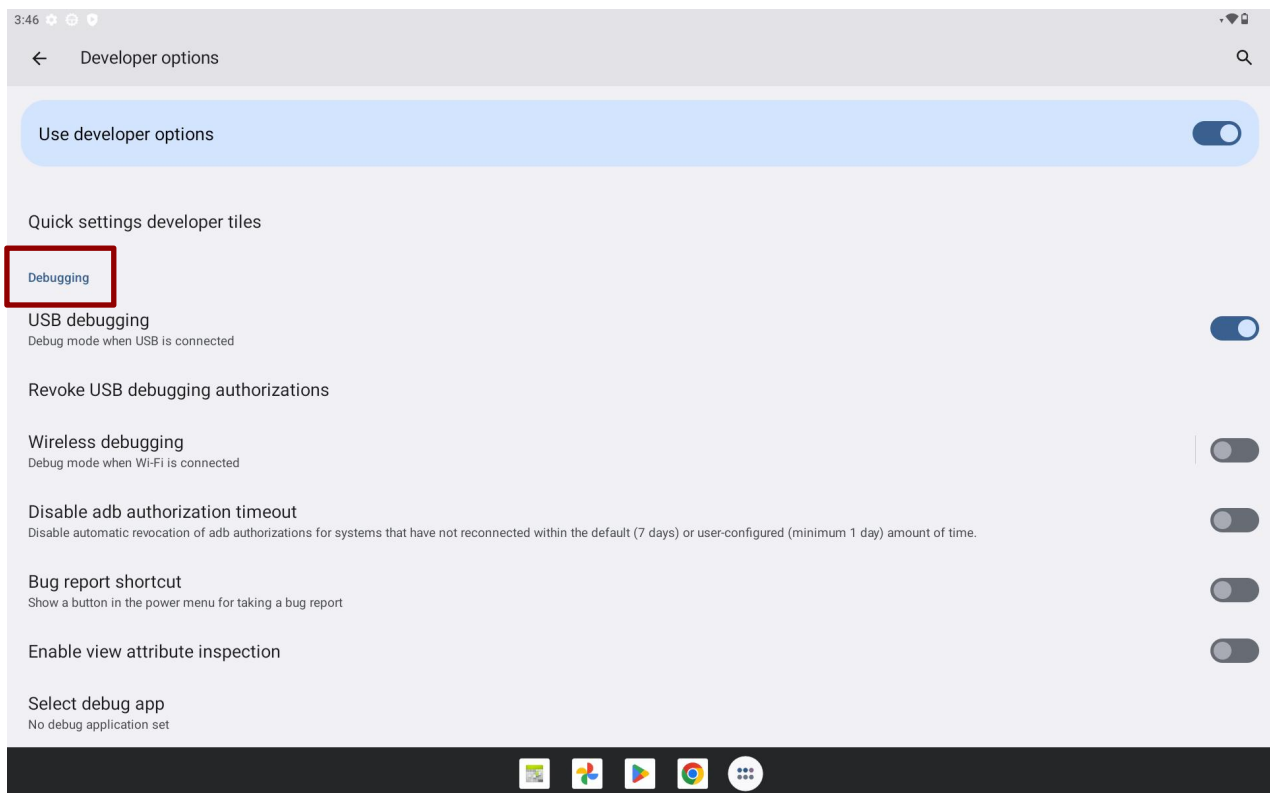Important: Carefully follow your computer's instructions to un-mount the USB storage or SD card. Disconnect USB device not properly can cause a loss of important data or information.

Chapter 6: Using Camera

This chapter explains how to use camera to capture photos and videos.

6.1 Capturing Photos or Video

Winmate® MD10W2 is equipped with an 13MP Camera on the rear side and an 8 MP Camera on the front side, which user can use to capture photos or video.

6.1.1 Opening Camera


1. Tap the launcher and then tap **Camera** . Camera opens in Portrait mode, ready to take a picture or shoot a video.

Important: Make sure the protective lens cover is clean before taking pictures. A lens cover with smudges from fingers can cause blurry pictures.

6.1.2 Taking Pictures

1. Adjust the exposure, flash, and other necessary settings or leave it automatic.
2. Tap the Shutter icon .
3. The camera brings the image into focus, when the image is in focus, the focus indicators in each corner will turn green and ready to capture.
4. Tap , then  to change camera settings.

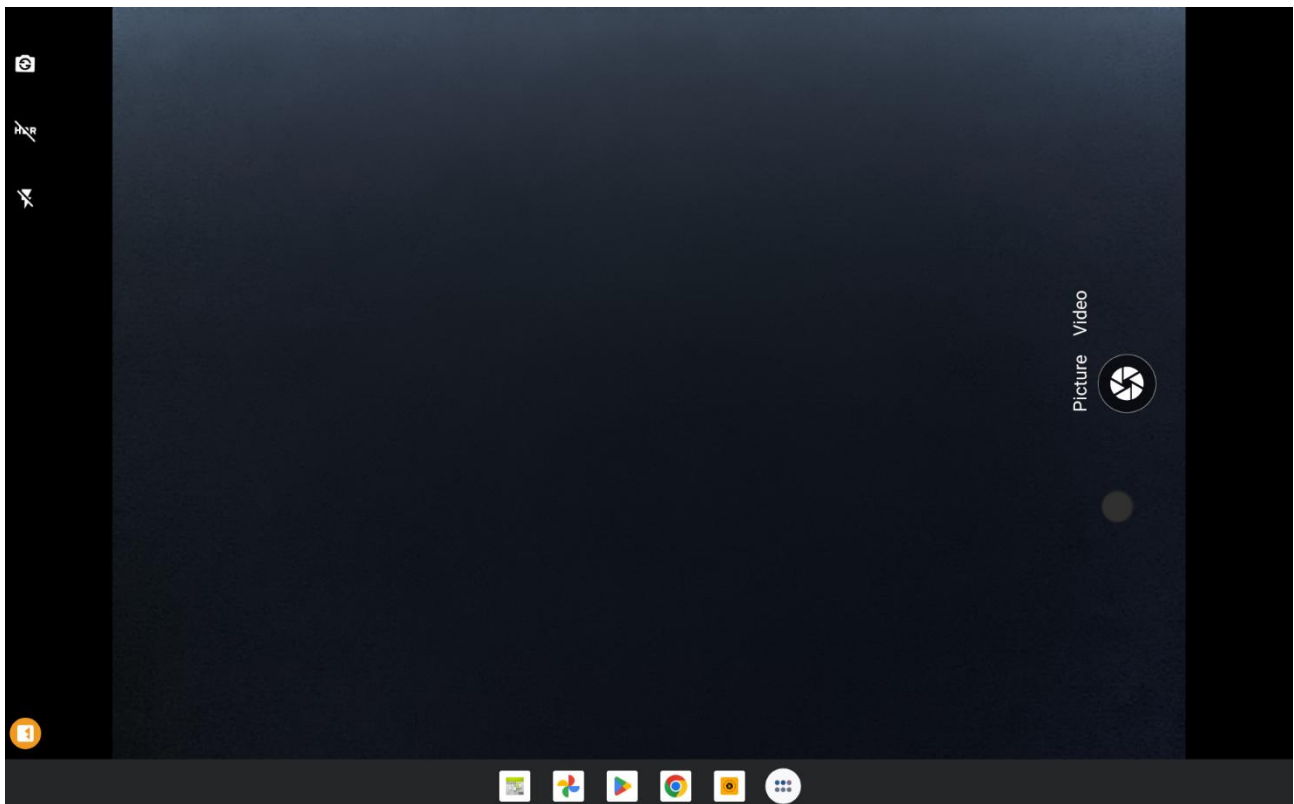
6.1.3 Viewing pictures

1. After taking a photo, a thumbnail of the picture is shown on the corner of the capture screen. Tap the thumbnail to view it.
2. Tap the onscreen button:
 -  to share a photo.
 -  to delete a photo.
 -  to edit the photo
 -  For more Camera options such as Details, Print.
3. Swipe right or left to view more photos.

6.1.4 Shooting Videos

1. Point the lens to record the scene where you want to start.
2. Tap  to start recording video.
3. The duration shoot is located on the top.
4. Tap the  to stop shooting.

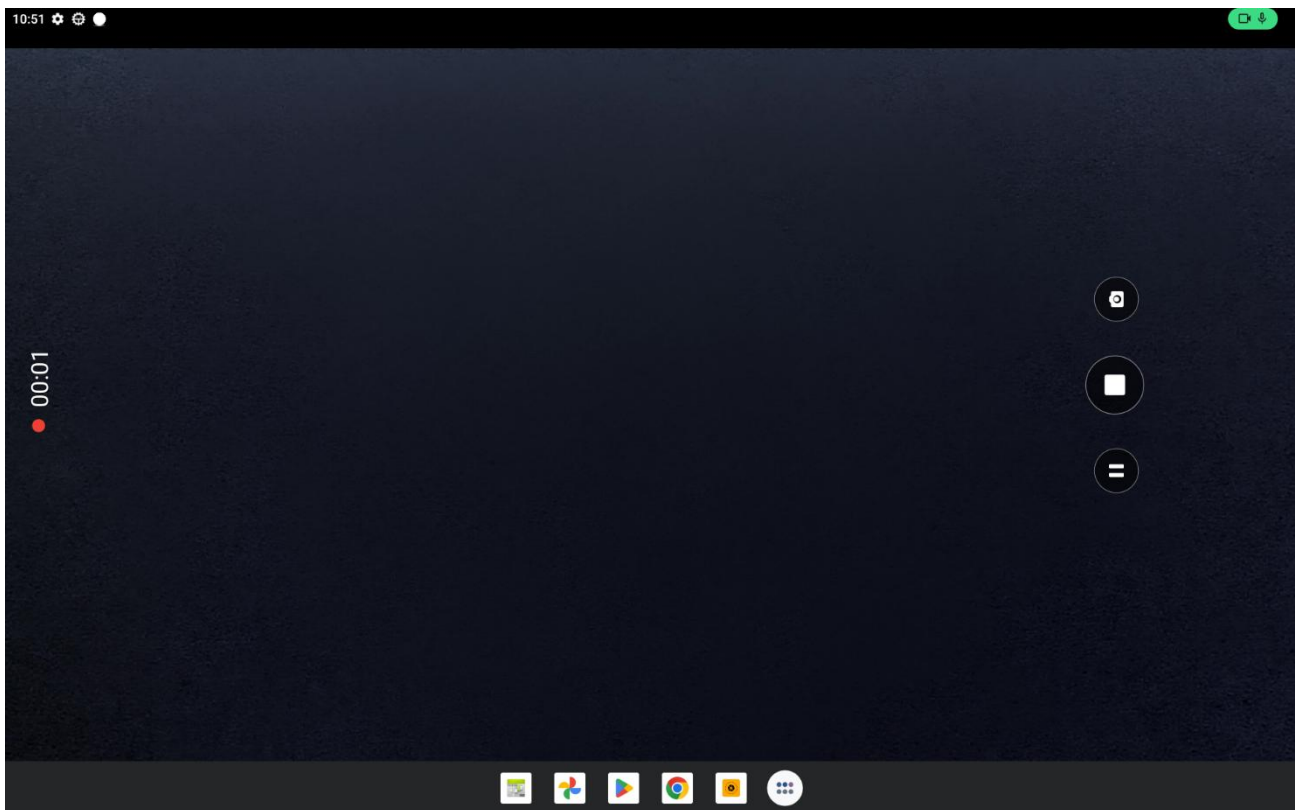The time remaining to shoot counts down at the left; when less than a minute is left, the countdown turns red.

6.2 The Gallery

1. Go to the launcher and then tap **Gallery** .
2. Tap a folder to view the pictures or video inside in full screen mode.
3. Flick the pictures left or right to view the previous or next picture.

Appendix

Appendix A: Basic Troubleshooting

Problem	Solution
The removable battery is not charging.	Make sure the AC adapter is properly connected, and that the removable battery is properly inserted.
Nothing appears on the screen.	- The tablet computer may be in sleep mode. Press the power button to wake up the display. - The brightness level may be too low. Increase brightness.
The screen suddenly turns blank and the power indicator goes off.	- Make sure the AC adapter is properly connected. - If running on battery power, the battery may be depleted, charge the battery. - Make sure that LED Indicator shows orang.
Wireless LAN signal quality is poor.	- The tablet computer may be out of range. Move the tablet computer closer to the access point or another wireless LAN device associated with it. - Check if your surroundings have interferences, such as microwave ovens, cordless phones. Move the tablet computer away from objects causing interference.
I cannot connect to another wireless LAN device.	- Make sure that wireless LAN feature is powered on. - Make sure that the SSID setting is the same for every wireless LAN device in the network. - Restart the tablet computer. - Make sure the IP address or subnet mask setting is correct.
I cannot connect to another device with BT.	- Make sure that both devices have BT powered on (For the Tablet Computer, check Settings to check if BT is powered on). - Make sure that the distance between the two devices is within 10 meters and there are no walls or large obstructions between the devices. - Make sure the both devices are powered on and configured to be discoverable. - Make sure that both devices are compatible.
The GPS application cannot find my location.	- Make sure that GPS is turned on in Settings. - Make sure that GPS receiver and the location sensor are enabled in Control Panel. - The tablet computer may be out of coverage area. Move your tablet computer outdoors or closer to a window. - Check if your surroundings have interferences, such as microwave ovens, cordless phones. Move the tablet computer away from objects causing interference.

Problem	Solution
No sound is produced.	• Adjust volume of the tablet computer by pressing on the volume rocker. • If using an external audio device, make sure that the cable is properly connected.
I cannot project the tablet computer display to another monitor.	• Make sure the VGA driver is installed. • Make sure the output setting selected is correct.
Battery leaking or swollen	• The lithium battery used in this tablet may swell after being charged/ discharged many times, but there is no problem regarding safety. • Continued use in such a state may result in the battery getting jammed, so please stop using the battery, and purchase a new one.

NOTES

FCC Statement

a、 § 15.19 Labelling requirements.

This device complies with part 15 of the FCC Rules. Operation is subject to the condition that this device does not cause harmful interference.(10.1 "Medical Android Tablet License")

a、 § 15.19 Labeling requirements.

This device complies with part 15 of the FCC Rules. Operation is subject to the following

two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.(10.1 "Medical Android Tablet License")

b、 § 15.21 Changes or modification warning.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

c、 § 15.105 Information to the user.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

-Reorient or relocate the receiving antenna.

-Increase the separation between the equipment and receiver.

-Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-Consult the dealer or an experienced radio/TV technician for help.

d、 Specific Absorption Rate (SAR) information:

This 10.1" Medical Android Tablet meets the government's requirements for exposure to radio waves. The guidelines are based on standards that were developed by independent

scientific organizations through periodic and thorough evaluation of scientific studies. The standards include a substantial safety margin designed to assure the safety of all persons regardless of age or health.

e、 FCC RF Exposure Information and Statement

The SAR limit of USA (FCC) is 1.6 W/kg averaged over one gram of tissue. Device types:

MD10W2 (FCC ID: RZ5-MD10W2) has also been tested against this SAR limit. when properly worn on the body is 0.527W/kg. This device was tested for typical body-worn operations with the back of the handset kept 0mm from the body. To maintain

compliance with FCC RF exposure requirements, use accessories that maintain a 0mm separation distance between the user's body and the back of the handset. The use of belt clips, holsters and similar accessories should not contain metallic components in its assembly. The use of accessories that do not satisfy these requirements may not comply with FCC RF exposure requirements, and should be avoided.

Body-worn Operation

This device was tested for typical body-worn operations. To comply with RF exposure requirements, a minimum separation distance of 0mm must be maintained between the user's body and the handset, including the antenna. Third-party belt-clips, holsters, and similar accessories used by this device should not contain any metallic components.

Body-worn accessories that do not meet these requirements may not comply with RF exposure requirements and should be avoided. Use only the supplied or an approved antenna.

For 5150-5250 frequency band, add below statement for indoor use warning statement on manual, Operations in the 5.15-5.25GHz band are restricted to indoor usage only.

*** RF warning for Portable device:**

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Contact Us

Winmate Inc.

No.18, Zhongxing South Rd.,

San-Chung District, New Taipei City 241017, Taiwan

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Website: www.winnate.com.cn

Winmate JP Office/ HPC System Inc.

LOOP-X 8F,3-9-15 Kaigan,
Minato-ku, Tokyo 108-0022, Japan

Tel: +81-3-5446-5535

Fax: +81-3-5446-5550

Website: www.hpc.co.jp

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Winmate USA Inc.
TL Electronic GmbH