

Pure 13 BT primax

User Guide | Master

preliminary



Hearing
Systems

SIEMENS

Content

Welcome	4
Your hearing instruments	5
Instrument type	5
Getting to know your hearing instruments	5
Components and names	6
Controls	8
Settings	9
Batteries	10
Battery size and handling tips	10
Replacing batteries	11
Daily use	12
Turning on and off	12
Inserting and removing the hearing instruments	14
Adjusting the volume	18
Changing the hearing program	19
Further adjustments (optional)	19
Special listening situations	20
On the phone	20
Audio streaming and control with iPhone	21
Bluetooth off	21

Maintenance and care	22
Hearing instruments	22
Ear pieces	23
Professional maintenance	25
Further information	26
Safety information	26
Accessories	26
Symbols used in this document	26
Troubleshooting	27
Country-specific information	28
Service and warranty	31

Welcome

Thank you for choosing our hearing instruments to accompany you through everyday life. Like anything new, it may take you a little time to become familiar with them.

This guide, along with support from your Hearing Care Professional, will help you understand the advantages and greater quality of life your hearing instruments offer.

To receive the most benefit from your hearing instruments, it is recommended that you wear them every day, all day. This will help you get adjusted to them.



CAUTION

It is important to read this user guide and the safety manual thoroughly and completely. Follow the safety information to avoid damage or injury.

Your hearing instruments



This user guide describes optional features that your hearing instruments may or may not have. Ask your Hearing Care Professional to indicate the features that are valid for your hearing instruments.

Instrument type

Your hearing instruments are a RIC (Receiver-in-Canal) model. The receiver is placed within the ear canal and connected to the instrument via a receiver cable. The instruments are not intended for children under the age of 3 years or persons with a developmental age of under 3 years.

Your hearing instruments support direct streaming technology allowing for easy and seamless audio streaming and data exchange with your iPhone*.

Getting to know your hearing instruments


It is recommended to familiarize yourself with your new hearing instruments. With the instruments in your hand try using the controls and note their location on the instrument. This will make it easier for you to feel and press the controls while wearing the hearing instruments.

*iPhone is a trademark of Apple Inc., registered in the U.S. and other countries.

i

If you have problems pressing the controls of your hearing instruments while wearing them, you can use a remote control to control your instruments. With a dedicated smartphone app you have even more control options.

Components and names

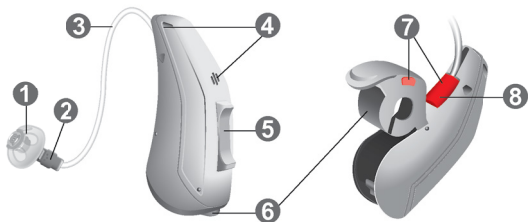

-
- | | |
|---------------------------|---|
| ① Ear piece | ⑥ Battery compartment (on/off switch) |
| ② Receiver | ⑦ Side indicator (red = right ear, blue = left ear) |
| ③ Receiver cable | ⑧ Receiver connection |
| ④ Microphone openings | |
| ⑤ Rocker switch (control) | |
-

You can use the following standard ear pieces:


Standard ear pieces	Size
A small, white, cylindrical ear piece with a flared base.	Click Sleeve (vented or closed)
A small, white, dome-shaped ear piece with a central opening.	Click Dome™ single (open or closed)
A small, white, dome-shaped ear piece with a semi-circular opening.	Click Dome semi-open
A small, white, dome-shaped ear piece with two circular openings.	Click Dome double

You can easily exchange the standard ear pieces.
Read more in section "Maintenance and care".

Custom-made ear pieces		
☐		Custom shell
☐		Click Mold™

Controls

With the rocker switch you can, for example, adjust the volume or switch hearing programs. Your Hearing Care Professional has programmed your desired functions to the rocker switch.

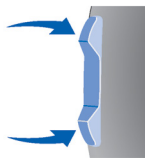

Rocker switch function	L	R
Press briefly:		
Program up/down	☐	☐
Volume up/down	☐	☐
Tinnitus therapy signal level up/down	☐	☐
Transmitter mode on/off	☐	☐
Press for about 2 seconds:		
Program up/down	☐	☐
Volume up/down	☐	☐
Tinnitus therapy signal level up/down	☐	☐
Transmitter mode on/off	☐	☐
Press for more than 3 seconds:		
Standby/turn on	☐	☐
Bluetooth on/off	☐	☐

L = Left, R = Right



You can also use a remote control to change hearing programs and adjust the volume or your hearing instruments. With our smartphone app you have even more control options.

Settings

Hearing programs

1

2

3

4

5

6

Read more in section "Changing the hearing program".

Features

- ☐ **Power-on delay** enables whistle-free insertion of hearing instruments.
Read more in section "Turning on and off".
- ☐ The **tinnitus therapy feature** generates a sound to distract you from your tinnitus.
- ☐ The **e2e wireless** functionality enables advanced audiological features and remote control options.

Batteries

When the battery is low the sound becomes weaker or you will hear an alert signal. The battery type will determine how long you have to replace the battery.

Battery size and handling tips

Ask your Hearing Care Professional for recommended batteries.

Battery size:	13
---------------	----

- Always use the correct battery size for your hearing instruments.
- Remove the batteries if you intend not to use the hearing instruments for several days.
- Always carry spare batteries.
- Remove empty batteries immediately and follow your local recycling rules for battery disposal.

Replacing batteries

Removing the battery:

- Open the battery compartment.



- Use the magnet stick to pull out the battery. The magnet stick is available as an accessory.



Inserting the battery:

- If the battery has a protective film, remove it only when you are ready to use the battery.



- Insert the battery with the "+" symbol facing up (refer to the picture).



- Carefully close the battery compartment. If you feel resistance, the battery is not inserted correctly.

Do not attempt to close the battery compartment by force. It can be damaged.

Daily use

Turning on and off

You have the following options to turn your hearing instruments on or off.

Via battery compartment:

- ▶ Turning on: Close the battery compartment.
The **default** volume and hearing program are set.
- ▶ Turning off: Open the battery compartment to the first catch.

Via rocker switch:

- ▶ Turning on or off: Press and hold the rocker switch for several seconds. Refer to section "Controls" for the settings of the rocker switch.
After turning on, the **previously** used volume and hearing program remain set.

When wearing the hearing instruments, an alert tone can indicate when an instrument is being turned on or off.

When the **power-on delay** is activated, the hearing instruments turn on after a delay of several seconds.

During this time you can insert the hearing instruments into your ears without experiencing unpleasant feedback whistling.

The "power-on delay" can be activated by your Hearing Care Professional.

Inserting and removing the hearing instruments

Your hearing instruments have been fine-tuned for your right and left ear. Colored markers indicate the side:


- red marker = right ear
- blue marker = left ear



Inserting a hearing instrument:

- For Click Sleeves, take care that the bend of the Click Sleeve is in line with the bend of the receiver cable.


Correct:

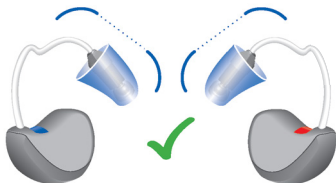

False:



- ▶ Hold the receiver cable at the bend closer to the ear piece.
- ▶ Carefully push the ear piece in the ear canal **1**.
- ▶ Twist it slightly until it sits well.
Open and close your mouth to avoid accumulation of air in the ear canal.
- ▶ Lift the hearing instrument and slide it over the top of your ear **2**.

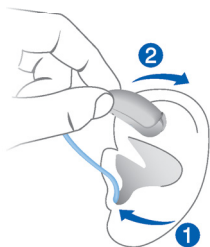

CAUTION

Risk of injury!

- ▶ Insert the ear piece carefully and not too deeply into the ear.



- It may be helpful to insert the **right** hearing instrument with the **right** hand and the **left** hearing instrument with the **left** hand.
- If you have problems inserting the ear piece, use the other hand to gently pull your earlobe downwards. This opens the ear canal and eases insertion of the ear piece.


The optional retention strand helps to securely retain the ear piece in your ear. To position the retention strand:

- Bend the retention strand and position it carefully into the bottom of the bowl of your ear (refer to the picture).



Removing a hearing instrument:

- ▶ Lift the hearing instrument and slide it over the top of your ear ①.
 - ▶ If your hearing instrument is equipped with a custom shell or a Click Mold, remove it by pulling the small removal cord toward the back of your head.
 - ▶ For all other ear pieces: Grip the receiver in the ear canal with two fingers and pull it carefully out ②.
- Do not pull the receiver cable.

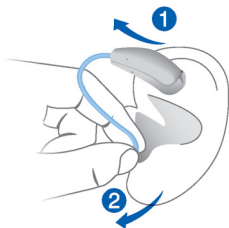

CAUTION

Risk of injury!

- ▶ In very rare cases the ear piece could remain in your ear when removing the hearing instrument. If this happens, have the ear piece removed by a medical professional.

Clean and dry your hearing instruments after usage.

Read more in section "Maintenance and care".

Adjusting the volume

Your hearing instruments automatically adjust the volume to the listening situation.

- If you prefer manual volume adjustment, press the rocker switch.

Refer to section "Controls" for the setting of the rocker switch.

An optional signal tone can indicate the volume change. Once the maximum or minimum volume has been reached, you may hear an optional signal tone.

Changing the hearing program

Depending on the listening situation, your hearing instruments automatically adjust their sound.

Your hearing instruments may also have several hearing programs which allow you to change the sound, if needed. An optional signal tone can indicate the program change.

- To change the hearing program, press the rocker switch.

Refer to section "Controls" for the setting of the rocker switch. Refer to section "Settings" for a list of your hearing programs.

Further adjustments (optional)

You can also use the hearing instrument controls to change other features, for example, the level of the tinnitus therapy signal.

Refer to section "Controls" for the setting of your controls.

Special listening situations

On the phone

When you are on the phone, hold the telephone receiver slightly above your ear. The hearing instrument and telephone receiver have to be aligned. Turn the receiver slightly so that the ear is not completely covered.



Telephone program

You may prefer a certain volume when using the phone. Ask your Hearing Care Professional to configure a telephone program.

- Switch to the telephone program whenever you are on the phone.

If a telephone program is configured for your hearing instruments, it is listed in section "Settings".

Audio streaming and control with iPhone

Your hearing instruments are **Made for iPhone** hearing instruments. This means you can get phone calls and listen to music from your iPhone directly into your hearing instruments. And you can control your hearing instruments with the myControl App.

For more information about compatible iOS devices, pairing, audio streaming and other useful functions, contact your Hearing Care Professional.

Bluetooth off

Your hearing instruments feature Bluetooth®* low energy technology. In areas where the use of Bluetooth technology is not allowed (e.g. on some airplanes), you can turn Bluetooth off temporarily. The hearing instruments still work without Bluetooth, however, direct audio streaming will not be possible.

- To turn Bluetooth off or on, use the smartphone app or use the controls of your hearing instruments.

Refer to section "Controls" for more information.

An optional alert tone can indicate when Bluetooth is being turned on or off.

*The Bluetooth word mark and logos are owned by the Bluetooth SIG, Inc., and any use of such marks the legal manufacturer of this product is under licenses. Other trademarks and trade names are those of their respective owners.

Maintenance and care

To prevent damage it is important that you take care of your hearing instruments and follow a few basic rules, which will soon become a part of your daily routine.

Hearing instruments

Drying and storage

- ▶ Dry your hearing instruments overnight.
- ▶ Ask your Hearing Care Professional for recommended drying products.
- ▶ For longer periods of non-use, store your hearing instruments with open battery compartment and batteries removed in a drying system to avoid the adverse effects of moisture.

Cleaning

Your hearing instruments have a protective coating. However, if not cleaned regularly you may damage the hearing instruments or cause personal injury.

- ▶ Clean your hearing instruments daily with a soft, dry tissue.
- ▶ Never use running water or immerse the devices in water.
- ▶ Never apply pressure while cleaning.



- Ask your Hearing Care Professional for recommended cleaning products, special care sets, or more information on how to keep your hearing instruments in good condition.

Ear pieces


Cleaning

Cerumen (ear wax) may accumulate on the ear pieces. This may affect sound quality. Clean the ear pieces daily.

For all types of ear pieces:


- Clean the ear piece right after removal with a soft and dry tissue.

This prevents that cerumen becomes dry and hard.

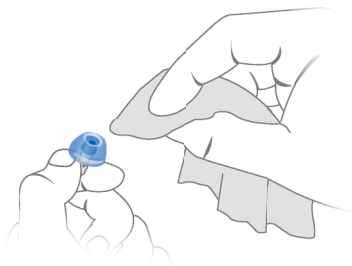

For Click Domes or Click Sleeves:

- Squeeze the tip of the ear piece.

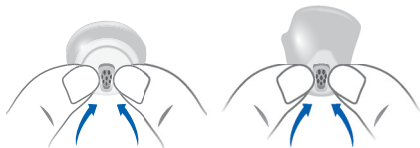


Exchanging

Exchange standard ear pieces approximately every three months. Replace them sooner if you notice cracks or other changes. The procedure for exchanging standard ear pieces depends on the type of ear piece. In section "Components and names", your Hearing Care Professional has marked your type of ear piece.

Exchanging Click Domes

- ▶ To remove the old Click Dome, use the removal tool or follow the instructions on the packaging of the Click Domes. Be careful not to pull the receiver cable.

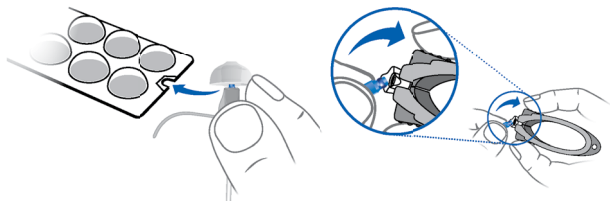

The removal tool is available as an accessory.

- ▶ Take particular care that the new Click Dome clicks noticeable in place.
- ▶ Check the firm fit of the new Click Dome.



Exchanging Click Sleeves

- ▶ Grip the Click Sleeve with your fingers and pull it off.
- ▶ Take particular care that the new Click Sleeve clicks noticeably in place.



Professional maintenance

Your Hearing Care Professional can perform a thorough professional cleaning and maintenance.

Custom-made ear pieces and wax filters should be exchanged as needed by your Hearing Care Professional.

Ask your Hearing Care Professional for your individual recommendation of maintenance intervals and for support.

Further information

Safety information

For more safety information, refer to the safety manual supplied with the device.

Accessories

You can use a smartphone app for convenient remote control of your hearing instruments. Additionally, you can create personal settings, for example, deactivate signal tones.

You can also use a remote control. Consult your Hearing Care Professional for further information.

Symbols used in this document



Points out a situation that could lead to serious, moderate, or minor injuries.



Indicates possible property damage.



Advice and tips on how to handle your device better.



"Made for iPhone" means that an electronic accessory has been designed to connect specifically to iPhone and has been certified by the developer to meet

Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone may affect wireless performance.

Troubleshooting

Problem and possible solutions

Sound is weak.

- Increase the volume.
 - Replace the empty battery.
 - Clean or exchange the ear piece.
-

Hearing instrument emits whistling sounds.

- Reinsert the ear piece until it fits securely.
 - Decrease the volume.
 - Clean or exchange the ear piece.
-

Sound is distorted.

- Decrease the volume.
 - Replace the empty battery.
 - Clean or exchange the ear piece.
-

Hearing instrument emits signal tones.

- Replace the empty battery.
-

Problem and possible solutions

Hearing instrument does not work.

- Turn the hearing instrument on.
 - Gently close the battery compartment completely.
 - Replace the empty battery.
 - Make sure the battery is inserted correctly.
 - Power-on delay is active. Wait several seconds and check again.
-

Consult your Hearing Care Professional if you encounter further problems.

Country-specific information

Risk mitigations (for internal use only)

M18_1/2 | replace battery11

M31_2/3, M48_2/3, M95-2/3 | not too deeply 15

M31_3/3, M48_3/3, M89_2/2, M95_3/3 | remain in ear | HI 17

M44_1/2 | sound is weak27

M44_2/2 | sound is distorted.....27

M49 | ask HCP28

M103 | cleaning22

Service and warranty

Serial numbers	Instrument	Receiver unit
----------------	------------	---------------

Left:		
-------	--	--

Right:		
--------	--	--

Service dates

1:	4:
----	----

2:	5:
----	----

3:	6:
----	----

Warranty

Date of purchase:

Warranty period in months:

Your Hearing Care Professional



Legal Manufacturer

Signia GmbH

Henri-Dunant-Strasse 100

91058 Erlangen

Germany

Phone +49 9131 308 0

Manufactured by Signia GmbH

under Trademark License of Siemens AG.

Document No. **preliminary**-99T##-#### ##

Order/Item No. ### ### ##

Master Rev01, 02.2017

© Signia GmbH, ##.2017

www.signia-hearing.com