Bluetooth® Blood Pressure Monitor

Instruction
Booklet

Clinically
Proven
Accurate

Model #BP3KV1-5WPC

BP3KV1-5WPC-0325

Important Product and Safety Information

Follow Instructions for Use. This document provides important product operation and safety information regarding this Blood Pressure Monitor. Please read this document thoroughly before using the device and keep for future reference.

This blood pressure monitor is an automatic digital blood pressure measuring device for use by adults, on the upper arm at home or in the doctor's office. It enables a very fast and reliable measurement of the systolic and diastolic blood pressure, as well as the pulse rate, by way of the oscillometric method during cuff inflation. This device detects the appearance of irregular heartbeats during measurement and provides a warning signal when the irregular heartbeat is detected. The device can be connected via Bluetooth® to a smart mobile device running the App.

This blood pressure monitor IS intended to be used:

- For self-measurement/monitoring of blood pressure and pulse in adults; whereas the person being measured may be the user/operator of the device.
- Within a home healthcare environment.
- With a cuff located upon the user's upper arm; ½ inch above the elbow over the artery as indicated on the cuff.

This blood pressure monitor IS NOT intended for use with:

- Pregnancy
- Preeclampsia
- Children under 12
- Neonatal patients

WARNING – Potentially hazardous situation that if not avoided may result in serious injury or death.

- Self-measuring means monitoring, not diagnosis or treatment. Unusual values must always be discussed with your doctor. Under no circumstances should you independently alter the dosages of any drugs prescribed by your physician.
- Consult your physician before using this device if any of the following or similar conditions are present: arrhythmias such as atrial or ventricular premature beats or atrial fibrillation, arterial sclerosis, poor perfusion, diabetes, age, pregnancy, preeclampsia, renal diseases. Motion during measurement, including trembling or shivering may affect the measurement.

- Though not for use with children under 12, ensure that any children around this device are supervised; some parts are small enough to be swallowed and any present tubes or cables may provide a risk of strangulation.
- Ensure the cuff tubing is not kinked during use as harmful injury may occur due to the effect of blood flow interference caused by high pressure in the cuff not releasing.

DO NOT:

- ✗ Use this device if you think it is damaged or anything appears unusual (e.g. sporadic operation, open packaging upon purchase).
- ✗ Use the displayed pulse for checking the frequency of heart pacemakers as this device is not suitable for this action.
- ✗ Open/modify this device; inaccuracy and/or harmful injury may result.
- ✗ Conduct many frequent measurements as this may result in harmful injury due to blood flow interference.
- ✗ Place the Cuff over a wound as this may cause additional injury.
- ✗ Place and pressurize the Cuff over/near any present intravascular access or therapy, or arteriovenous shunt, as this may cause blood flow interference and result in harmful injury.
- ✗ Place and pressurize the Cuff over a limb near the side of a mastectomy as this may cause harmful injury.
- ✗ Use the Blood Pressure Monitor on a limb simultaneously with other medical equipment on the same limb due to possible interference with such medical equipment.
- ✗ Maintain pressure in the Cuff applied to the limb for a prolonged amount of time. Ensure that circulation in the limb is not impaired by checking circulation if prolonged/repeated exposure to pressure occurs.
- ✗ Service the device when being used or when power is supplied. When the device is not in use and power is removed, cuff assembly, batteries and AC adapter (if applicable) may be replaced by the user with Microlife supplied replacements. No other parts/components are accessible.

⚠ CAUTION – Potentially hazardous situation that if not avoided may result in minor/moderate injury, property damage, and/or damage to the device.

- This device contains sensitive electronics components. Avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g., mobile telephones, microwave ovens). These can lead to temporary impairment of the measuring accuracy. Move the device to another location if interference is determined.

- When not using the Blood Pressure Monitor for extended lengths of time, remove the batteries to avoid potential battery leakage and damage to the monitor.
- When replacing the batteries, ensure all the batteries are replaced simultaneously to avoid battery damage and potential damage to the monitor. Microlife USA does not recommend using rechargeable batteries.
- To avoid inaccurate measurements and to lessen any discomfort from Cuff pressure, ensure the Cuff is placed correctly on the limb and fits correctly when snug (not tight), as indicated by markings with the Cuff.
- Consult your physician in cases of frequent irregular heartbeat detections.
- This Risk Indicator feature is provided in order to help you understand your potential blood pressure risk. However, this feature is neither a diagnosis nor a substitute for a medical examination. It is important to consult with your physician to determine your risk.

DO NOT:

- ✗ Drop this device or expose it to strong vibrations; sensitive components may be affected resulting in inaccuracies and/or operational issues.
- ✗ Use the Blood Pressure Monitor outside of its specified operation temperature and humidity rating, or if stored outside of its specified storage temperature and humidity rating. Avoid storage in direct sunlight.
- ✗ Use this device in a moving vehicle; inaccurate measurements may result.
- ✗ Use third party accessories. Only use Microlife authorized accessories, such as cuffs or AC adapters, as those not approved for use with the device may provide inaccurate measurements, injury, and/or damage the device.

STANDARDS:

In addition to the standards stated in the Instruction Manual:

- This medical device is compliant with medical device and non-invasive blood pressure monitor standards IEC 60601-1, IEC 60601-1-2, IEC 60601-1-11, and AAMI/ANSI/IEC 80601-2-30, and
- Electromagnetic standards IEC 60601-1-2 along with FCC Part 15, and
- Clinical Testing per standard ISO 81060-2 was conducted on blood pressure device using the same measurement technology.

TYPE:

IP20: Protected against solid foreign particles with a diameter of more than 12.5 mm, no protection against water. Keep Dry. 

Type BF Applied Part

Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

Expected Life:

Monitor: 5 Years

Cuff: 2 Years

FCC

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

* RF warning for Portable device:

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Bluetooth® Blood Pressure Monitor**Table of contents****1. Introduction**

- 1.1. Your automatic blood pressure monitor

2. How is blood pressure measured?

- 2.1. What is blood pressure, and why does it fluctuate?
- 2.2. Measurement guidelines
- 2.3. What can I do to change my blood pressure?

3. Components of your blood pressure monitor**4. How do I get started?**

- 4.1. Inserting the batteries
- 4.2. Protective screen film
- 4.3. Setting the date and time
- 4.4. Tube connection
- 4.5. Select the user
- 4.6. Select the measuring mode: Standard or Measurement Averaging Mode (MAM)
- 4.7. Measurement Averaging Mode technology (MAM)

5. How do I take a measurement?

- 5.1. Preparing to take a measurement
- 5.2. Common errors
- 5.3. Using the cuff
- 5.4. Taking a measurement
- 5.5. Irregular heartbeat detector
- 5.6. Hypertension risk indicator
- 5.7. MyTrend analysis
- 5.8. Memory: reviewing readings
- 5.9. Stopping a measurement
- 5.10. Battery change indicator
- 5.11. Using the AC power adapter

6. Bluetooth® functions and application (App)

- 6.1. Activating your Bluetooth®
- 6.2. Pairing your smartphone

7. Error messages/troubleshooting

8. Care and maintenance

9. Limited warranty

10. Technical specifications

11. How to contact us

1. Introduction

1.1. Your automatic blood pressure monitor

Thank you for purchasing a fully automatic blood pressure monitor. Your monitor is designed to provide fast and reliable digital readings of your pulse, and systolic / diastolic blood pressure using the oscillometric method on your upper arm. It offers clinically proven accuracy and has been designed to be user friendly.

Before using your blood pressure monitor, please read this instruction manual carefully to ensure correct use. If you have additional questions regarding blood pressure measurements please contact your doctor.

2. How is blood pressure measured?

2.1. What is blood pressure, and why does it fluctuate?

Your level of blood pressure is determined in the circulatory center of the brain and adjusts to a variety of situations through feedback from the nervous system. To adjust blood pressure, the strength and speed of the heart (Pulse), as well as the width of circulatory blood vessels is altered. Blood vessel width is controlled by fine muscles in the blood vessel walls.

Your level of arterial blood pressure changes periodically during heart activity: During the "blood ejection" (Systole) the value is highest (systolic blood pressure value). At the end of the heart's "rest period" (Diastole) pressure is lowest (diastolic blood pressure value).

Blood pressure values must lie within certain normal ranges in order to prevent particular diseases.

2.2. Measurement guidelines

Blood pressure is very high if your systolic blood pressure is over 140 mmHg and/or your diastolic pressure is above 90 mmHg, **while at rest**. In this case, please consult your doctor immediately. Long-term values at this level endanger your health due to continual damage to the blood vessels in your body.

Should the systolic blood pressure values lie between 130 mmHg and 139 mmHg and/or the diastolic blood pressure values lie between 80 mmHg and 89 mmHg, consult your doctor. Regular self-checks will be necessary.

If you have blood pressure values that are too low, (i.e., systolic values under 105 mmHg and/or diastolic values under 60 mmHg), consult your physician.

Even with normal blood pressure values, a regular self-check with your blood pressure monitor is recommended. You can detect possible changes in your values early and react appropriately.

If you are undergoing medical treatment to control your blood pressure, keep a record of values along with time of day and date. Show these values to your physician. **Never use the results of your measurements to independently alter the drug doses prescribed by your physician.**

Which values are normal?

The following standards for assessing high blood pressure in adults have been established by the American Heart Association (AHA) and American College of Cardiology (ACC) in 2017.

Category	Systolic (mmHg)	Diastolic (mmHg)
Normal	<120	and <80
Elevated	120-129	and <80
Hypertension		
Stage 1 Hypertension	130-139	or 80-89
Stage 2 Hypertension	140-179	or 90-119
Hypertensive Crisis	≥180	or ≥120

Additional information

- This chart reflects a 2017 update to blood pressure standards.
- If your values are mostly normal under resting conditions but exceptionally high under conditions of physical or psychological stress, it is possible that you are suffering from so-called “labile hypertension”. Consult your doctor.

2.3. What can I do to change my blood pressure?

- Consult your doctor.
- Increased blood pressure values (various forms of hypertension) are associated with considerable health risks over time. Arterial blood vessels in your body are endangered due to constriction caused by deposits in the vessel walls (Arteriosclerosis). A deficient supply of blood to important organs (heart, brain, muscles) can result from arteriosclerosis. Furthermore, the heart will become structurally damaged with increased blood pressure values.
- There are many different causes of high blood pressure. We differentiate between the common primary (essential) hypertension, and secondary hypertension. The latter group can be ascribed to specific organ malfunctions. Please consult your doctor for information about the possible origins of your own increased blood pressure values.
- There are measures which you can take to reduce and even prevent high blood pressure. These measures must be permanent lifestyle changes.

1) Eating habits

- Strive for a normal weight corresponding to your age. See your doctor for your ideal weight.
- Avoid excessive consumption of common salt.
- Avoid fatty foods.

2) Previous illnesses

- Consistently follow all medical instructions for treating illness such as:
 - Diabetes (Diabetes mellitus or sugar diabetes)
 - Fat metabolism disorder
 - Gout

3) Habits




- Give up smoking completely.
- Drink only moderate amounts of alcohol.
- Restrict your caffeine consumption (coffee, tea, chocolate).

4) Your physical condition


- After a medical examination, and with your doctor's approval and direction; exercise.
- Choose sports which require stamina and avoid those which require strength.
- Avoid reaching the limit of your performance.
- With previous illnesses and/or an age of over 40 years, please consult your doctor before beginning your exercise routine. Speak with your doctor regarding the type and extent of exercise that is appropriate for you.

3. Components of your blood pressure monitor

a) Measuring unit

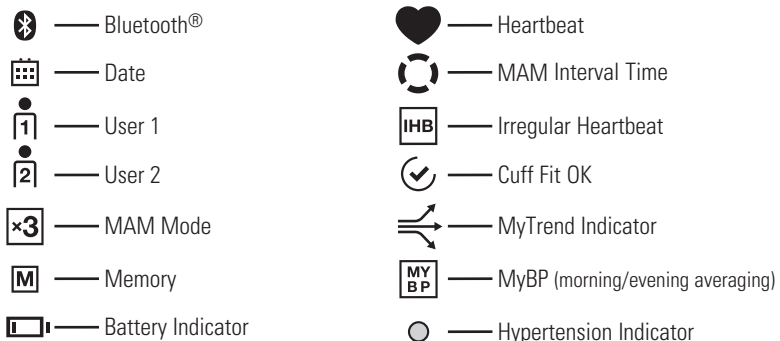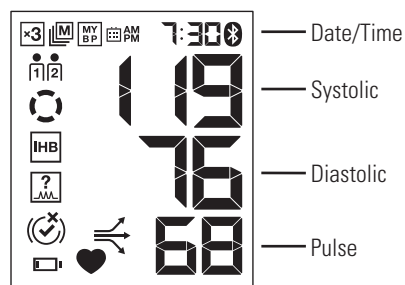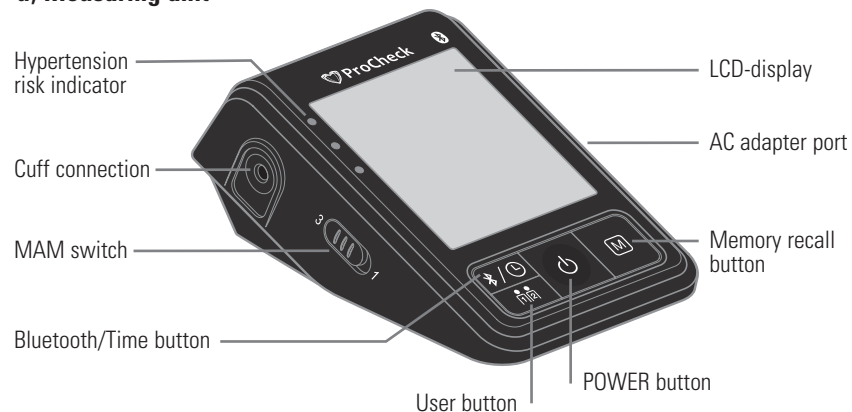
b) Cuff holder

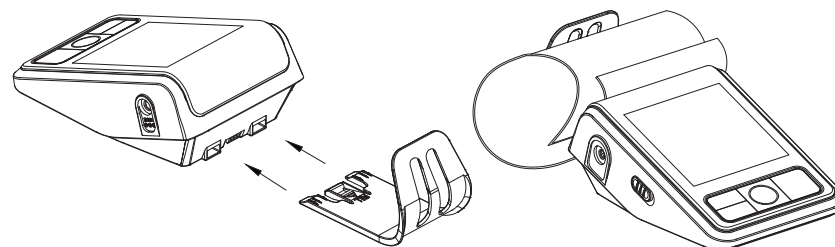c) Wide Range Easy-Fit Contoured Cuff:

For arm circumference 8.7" - 16.5" (22 - 42 cm)

If you ever need to buy a replacement cuff call customer service at 1-866-289-6611 for information.

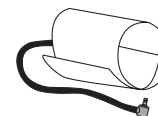Please Note:

Arm circumference should be measured with a measuring tape in the middle of the relaxed upper arm. Do not force cuff connection into the opening. Make sure the cuff connection is not pushed into the AC adapter port. If the cuff is too small, call 1-866-289-6611 for additional information.

4. How do I get started?

4.1. Inserting the batteries

Insert batteries immediately after unpacking the unit. The battery compartment is located on the bottom of the device.

- Remove cover.
- Insert the batteries (4 x size AA 1.5V), observing the indicated polarity. When replacing the batteries, replace all at the same time.
- If a battery warning appears in the display, the batteries are discharged and must be replaced.

Low Battery Indicator

Notes:


- We do not recommend using rechargeable batteries.
- If the blood pressure monitor is not used for long periods, remove the batteries from the device.

4.2. Protective screen film

There is a protective piece of static film covering the monitor's screen. Please remove it before use.

4.3. Setting the date and time


1. After the new batteries are activated, the year number flashes in the display. You can advance the year by pressing the "M" button. To confirm and then set the month, press the "BT/TIME" button.
2. You can now set the month using the "M" button. To confirm and then set the day, press the "BT/TIME" button.
3. Please follow the instructions above to set the day, hour and minutes.
4. Once you have set the last minute and pressed the "BT/TIME" button, the date and time are set and the time is displayed.
5. If you want to change the date and time, press and hold the "BT/TIME" button down for approx. 8 seconds or until the year number starts to flash. Now you can enter the new values as described above.

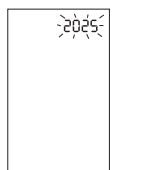Memory Button

BT/TIME Button

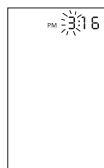
4.4. Tube connection

Make sure cuff connector is securely inserted into the side of your blood pressure monitor (a soft "click" may be heard when fully inserted).


Please Note:


- A loose connection will result in inaccurate readings, and an error message (Err 3).

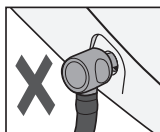Not engaged

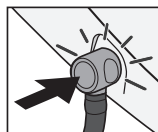No stem visible

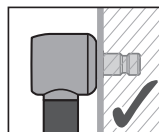4.5. Select the user

This blood pressure monitor is designed to store 99 measurements for each of two users. Before taking a measurement, be certain that the correct user has been selected.

- a) Press the "user button" to toggle between users, user 1 or user 2.

4.6. Select the measuring mode:

Standard or Measurement Averaging Mode (MAM)

This instrument enables you to select either Standard (single measurement) or Measurement Averaging Mode (automatic triple measurement).

- a) To select Standard mode, push the switch on the side of the monitor to Position "1".
- b) To select Averaging mode, push the switch to Position "3". Note the "MAM" icon will illuminate on the screen.

If you select 1, then only one measurement will be taken. If you select 3, the unit will inflate and deflate three times resulting in one final averaged measurement.

4.7. Measurement Averaging Mode technology (MAM)

- In Averaging mode, three measurements are automatically taken in succession and the result then automatically analyzed and displayed. Because your blood pressure constantly fluctuates, a result determined in this way is more reliable than one produced by a single measurement.
- After pressing the POWER button the selected Averaging mode appears in the display.
- The bottom, right hand section of the display shows a 1, 2 or 3 to indicate which of the 3 measurements is currently being taken.
- There is a break of 15 seconds between the measurements (15 seconds are adequate according to "Blood Pressure Monitoring, 2001, 6:145-147" for oscillometric instruments). A count down indicates the remaining time before the 2nd and 3rd readings will begin.
- The individual results from measurements 1, 2 and 3 are not initially displayed.
- Do not remove the cuff between measurements.

5. How do I take a measurement?

Please note: You should always be seated before and during measurement.

5.1. Preparing to take a measurement:

- Avoid eating and smoking as well as all forms of exertion directly before measurement. These factors influence the measurement result. Find time to relax by sitting in an armchair in a quiet atmosphere for about 5 minutes before taking a measurement.
- Remove any garment that fits closely to your upper arm.
- Always measure on the same arm (normally left).
- Always compare measurements taken at the same time of day, since blood pressure changes during the course of the day, as much as 20-40 mmHg.

5.2. Common errors:

Note: Comparable blood pressure measurements always require the same conditions! Conditions should always be quiet.

- All efforts by the user to support the arm can increase blood pressure. Make sure you are in a comfortable, relaxed position and do not flex any of the muscles in the measurement arm during the measurement. Use a cushion for support if necessary.
- If the arm artery lies considerably lower or higher than the heart, an erroneously high or low blood pressure will be measured! Each 15 cm (6 inch) difference in height between your heart and the cuff results in a measurement error of 10 mmHg!
- Cuffs that are too narrow or too short result in false measurement values. Selecting the correct cuff is extremely important. Cuff size is dependent upon the circumference of the arm (measured in the center). The permissible range is printed on the cuff. If this is not suitable for your use, please call 1-866-289-6611.
- A loose cuff or a sideways protruding air pocket causes false measurement values.
- With repeated measurements, blood accumulates in the arm, which can lead to false results. Consecutive blood pressure measurements should be repeated after a 1 minute pause or after your arm has been held up in order to allow the accumulated blood to flow.

If you decide to take your Averaging Mode measurement again, be sure to wait at least one minute beforehand.

5.3. Using the cuff

- Pass the end of the cuff through the flat metal ring so that a loop is formed. The hook and loop material must be facing outward. (Ignore this step if the cuff has already been prepared.)
- Slide the cuff onto your upper left arm. The rubber tube should be on the inside of your arm extending downward to your hand. Make certain the cuff lies approximately $\frac{1}{2}$ " (1 to 2 cm) above the elbow.

IMPORTANT: The red strip on the edge of the cuff (Artery Mark) must lie over the artery which runs down the inner side of the arm.




TIP: Align red artery mark to pinky finger.

- To secure the cuff, wrap it around your arm and press the hook and loop material together. Check the size range indicator on the cuff. The arrow should point within the size range on the inside of the cuff.
- There should be little free space between the arm and the cuff. You should be able to fit 2 fingers between your arm and the cuff. Clothing must not restrict the arm. Any piece of clothing which does must be removed. Cuffs that don't fit properly result in false measurement values. Measure your arm circumference if you are not sure of proper fit.
- Lay your arm on a table so the cuff is at the same height as your heart. Make sure the tube is not kinked.
- Remain seated quietly for five minutes before you begin the measurement.

Comment:

If it is not possible to fit the cuff to your left arm, it can also be placed on your right arm. However, all measurements should be made using the same arm.

Comparable blood pressure measurements always require the same conditions (relax for several minutes before a reading).

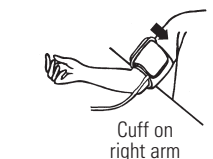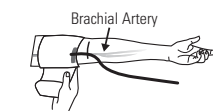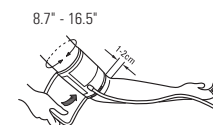5.4. Taking a measurement

After the cuff has been appropriately positioned on the arm and connected to the blood pressure monitor, the measurement can begin. Remain still during your measurement, do not flex muscles.

- Press the POWER button. The pump begins to inflate the cuff. In the display, the increasing cuff pressure is continually displayed.
- When the device has detected your pulse, the heart symbol in the display begins to blink.
- After automatically reaching an individual pressure, the pump stops and the pressure slowly falls. The cuff pressure is displayed during the measurement.
- When the cuff is positioned on the arm correctly, and is sufficiently tight, the "optimal" (cuff fit OK) cuff symbol will appear. When the cuff position is less than ideal, the "sub-optimal" symbol will appear.
- When the measurement has been concluded, the measured systolic and diastolic blood pressure values, as well as the pulse are now displayed.
- The measurement results are displayed until you switch the device off. If no button is pressed for 1 minute, the device switches off automatically.
- When the unit is set to the Averaging mode setting, 3 separate measurements will take place in succession after which your result is calculated and displayed as a single, averaged measurement. There is a 15 second resting time in-between each measurement. A count down indicates the remaining time before the 2nd and 3rd readings will begin.

If one of the measurements causes an error message, it will be repeated one more time. If any additional error occurs, the measurement will be discontinued and error code displayed.

Pumping
Pressure

Measuring

Systolic

Diastolic


Pulse

Measurement
complete

5.5. Irregular heartbeat detector

This symbol  indicates that certain pulse irregularities were detected during the measurement. In this case, the result may deviate from your normal basal blood pressure – repeat the measurement. In most cases, this is no cause for concern. However, if the symbol appears on a regular basis (e.g. several times a week with measurements taken daily), we advise you to tell your doctor. Please show your doctor the following explanation:

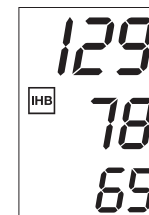Information for the doctor on frequent appearance of the Irregular Heartbeat Symbol.

This instrument is an oscillometric blood pressure monitor device that also analyzes pulse frequency during measurement. The instrument is clinically tested.

If pulse irregularities occur during measurement, the irregular heartbeat symbol is displayed after the measurement. If the symbol appears more frequently (e.g. several times per week on measurements performed daily) or if it suddenly appears more often than usual, we recommend the patient to seek medical advice.

The instrument does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

5.6. Hypertension risk indicator

The colored LED indicator lights on the left-hand edge of the display has been designed to provide a quick visual representation of your blood pressure. Once a measurement has been completed, one of the LED lights will illuminate and display onscreen next to your measurement. The color light that illuminates will show if the measurement is within the normal (green), borderline (orange) or danger (red) range.

This classification is based on standards established by the American Heart Association (AHA) and American College of Cardiology (ACC) in 2017.

The LED indicator lights:

- **Green** light, your measurement is “Normal”.
- **Orange** light, it is “Elevated” or “Stage 1 Hypertensive”.
- **Red** light, it is “Stage 2 Hypertensive” or “Hypertensive Crisis”.

5.7. MyTrend analysis

After each measurement, the MyTrend icon will illuminate on the screen to show how the current measurement is trending compared with the stored average measurement. Depending on the most recent measurement a series of arrows will appear on screen to show if the current measurement is trending up, down or flat compared to the stored averaged measurement. See MyTrend Guidelines below:

- If the Systolic or Diastolic measurement is more than 5 mmHg above the stored average value, then an upward trend arrow will illuminate.
- If the Systolic or Diastolic measurements are less than 5 mmHg below the stored average value, then a downward trend arrow will illuminate.
- If the Systolic or Diastolic measurements do not differ by more than 5 mmHg from the stored average value, then a flat/straight arrow will illuminate.
- If the Systolic and Diastolic measurements trend in different directions from the stored average value, then the Systolic reading will flash along with an up or down arrow for a few seconds followed by the Diastolic reading will flash along with an up or down arrow for a few seconds.

5.8. Memory: reviewing readings

At the end of a measurement, this monitor automatically stores each result with date and time. This unit stores 99 memories for each of 2 users.

Viewing the MyBP blood pressure average

Press the “M” button when device is off to turn on the memory mode. MyBP average includes only readings taken in the morning or evening of the most recent 3 to 7 days for computation, per clinical guideline suggestion. MyBP average is displayed only when sufficient number of readings meeting the criteria (below) have been obtained, to provide a more clinically relevant average. If data in the memory are insufficient, MyBP average is not displayed. Note: Measurements with suboptimal cuff fit are excluded from the average computation.

- Only measurement that were performed in the morning between (5:00 am-10:59 am) or in the evening (5:00 pm-10:59 pm) in the last 28 days are evaluated.
- Measurements performed in MAM mode or single standard mode are both treated as single measurement for MyBP average.
- A maximum of 4 readings per day are evaluated for the average (2 from the morning and 2 from the evening). MyBP average requires at least 12 readings (3 days), and evaluates up to 7 days of readings.

Viewing the stored values

With the unit off, press the “M” button. The display first shows “28 A”, which is an average of all measurements in the last 28 days. Pressing the “M” button again displays the previous value. To view a particular stored memory, press and hold the “M” button to scroll to that stored reading.

Please note: Measurements for each user are averaged and stored separately. Be certain that you are viewing the measurements for the correct user.


Memory full

When the memory has stored 99 results a new measured value is stored by overwriting the oldest value.

Clear all values

If you are sure that you want to permanently remove all stored values, hold down the “M” button (the instrument must have been switched off beforehand) until “CL” appears and then release the button. If you do not want to clear the values, press the POWER button. To permanently clear the memory, press the “BT/TIME” button while “CL” is flashing.

Individual values cannot be cleared.

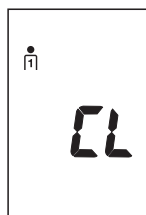
5.9. Stopping a measurement

If it is necessary to interrupt a blood pressure measurement for any reason (e.g. the patient feels unwell), the POWER button can be pressed at any time. The device then immediately lowers the cuff pressure automatically.

5.10. Battery change indicator

Batteries almost discharged

When the batteries are approximately 75% used the battery symbol will flash as soon as the instrument is switched on (if at least one of the batteries still has some charge). Although the instrument will continue to measure reliably, you should obtain replacement batteries.

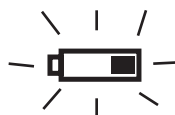Batteries discharged – replacements required

When the batteries are discharged, the battery symbol will appear unblinking as soon as the instrument is switched on. You cannot take any additional measurements and must replace the batteries.

1. Open the battery compartment at the bottom of the instrument.
2. Replace all of the batteries at the same time – ensure they are correctly connected, as shown on the symbols in the compartment.
3. The memory retains all values although date and time (and possibly also set alarm times) must be reset - the year number therefore flashes automatically after the batteries are replaced.
4. To set date and time, follow the procedure described in Section 4.3.

Notes:

Use four new, longlife 1.5V AA batteries. Do not use batteries beyond their expiration date. If the monitor is not going to be used for a prolonged period the batteries should be removed.

5.11. Using the AC power adapter (included)

You may also operate this monitor using the included AC adapter. Use only the included AC adapter to avoid damaging the unit.

- a) Ensure that the AC adapter and cable are not damaged.
- b) Plug the adapter cable into the AC adapter port on the right side of monitor.
- c) Plug the AC adapter into a 110 V power socket (U.S. or Canada).
- d) Test that power is available by pressing the POWER button.

Notes:

- No power is taken from the batteries while the AC adapter is connected to the instrument.
- If the power is interrupted during a measurement (e.g., by removal of the adapter from the wall socket), the instrument must be reset by removing the plug from the instrument.

If you have any questions regarding an AC adapter, call us at 1-866-289-6611.

6. Bluetooth® functions and application (App)

Please download the Microlife Connected Health+ US App

(referred to as App in this document) from Apple's® App Store® or Google Play® before pairing your devices:

Questions?

- To access help, go to the App home screen (Microlife Connected Health+ US) and press the "Menu" button (in the top left corner of the screen).
- Press "FAQ/Help" and select your topic.
- If you have any issues with your setup, please contact customer service toll free: 1-866-289-6611 for additional support.


6.1. Activating your Bluetooth®

The Bluetooth® icon indicator on your blood pressure monitor, located in the upper right portion of the screen, is designed to provide information about the connection between your device and the monitor. Please note the following:

- The Bluetooth® icon on your blood pressure monitor will flash initially when ready to connect with a device, and continue flashing for 2 minutes until a connection occurs.
- To manually activate the Bluetooth® on your monitor, press the BT/TIME button or until the Bluetooth® icon starts blinking.
- If a connection does not occur within 2 minutes, the Bluetooth® function on the monitor will automatically turn off.
- Set the date and time on your blood pressure monitor before pairing to your phone.

Note: The App will not be able to receive readings unless the date and time have been properly set on your blood pressure monitor.

Bluetooth® icon

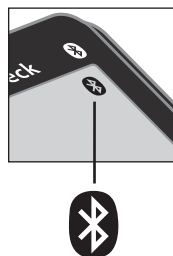6.2. Pairing your smartphone

- Download and install the App on your phone.
- Confirm phone's Bluetooth® function is turned on.
- Open the App on your phone.
 - Select your User Profile.
 - Select the blood pressure measurement page.
- Press Start on your monitor.

Once the measurement has been completed, the monitor will automatically transfer your reading to the App.
- If you have any issues please contact customer service toll free: 1-866-289-6611.

Setting up your blood pressure monitor

Confirm the correct user is selected on your blood pressure monitor:

- Confirm the user icon on the top left hand side of the monitor screen is set to user 1 or 2.
- To change users on your monitor, see section 4.5.

7. Error messages/troubleshooting

If an error occurs during a measurement, the measurement is discontinued and a corresponding error code is displayed. We recommend you switch the device off and take a new measurement making sure you allow time for relaxation (example: Error no. 2).

Err
2

Error No.	Possible cause(s)/Solutions
ERR 1 	The tube may have loosened, or no pulse was detected.* Ensure cuff connections are tight with proper cuff placement. See section 5.3.
ERR 2 	Unnatural pressure impulses influenced the measurement result. Reason: The arm was moved during the measurement (artefact). Repeat measurement, keeping still and quiet.
ERR 3 	Inflation of the cuff takes too long. The cuff is not correctly seated or the hose connection is not tight. Re-position cuff and repeat the measurement. Please follow the procedure described in "4.4 Tube connection" to ensure the cuff is properly connected.
ERR 5	The measured readings indicated an unacceptable difference between systolic and diastolic pressures. Take another reading following directions carefully. Contact your doctor if you continue to get unusual readings.
ERR 6	There were too many errors during the measurement in MAM mode, making it impossible to obtain a final result. Read through the checklist for performing reliable measurements and then repeat the measurement.*
HI	The cuff pressure is too high. Relax for 5 minutes and repeat the measurement.*
LO	The pulse is too low (less than 40). Repeat the measurement.*
BLUETOOTH BLINKS RAPIDLY	Bluetooth® connection error. Turn off the device Bluetooth® and close the app on the smartphone. Wait for 1 minute, open the app on the smartphone and manually activate Bluetooth® on the device, to re-try Bluetooth® connection and data transfer.
ERR BT	Bluetooth® self check error. Bluetooth® is malfunctioning. To solve the problem, please contact us at 1-866-289-6611.

***If this or any other problem occurs repeatedly, please consult your doctor.**

Other possible errors and their solutions

If problems occur when using the device, the following points should be checked:

Malfunction	Remedy
The display remains blank when the instrument is switched on although the batteries are in place.	• Check battery installation or AC adapter connection is correct and wall outlet is working. • If the display is unusual, remove the batteries and then exchange them for new ones. Check polarity.
The pressure does not rise although the pump is running.	• Check the connection of the cuff tube and connect properly.
The device frequently fails to measure, or the values measured are too low or high.	1. Ensure the cuff fits properly. (see 5.3) 2. Before starting measurement make sure that the cuff is not too tight and that clothing is not exerting pressure on the arm. Take articles of clothing off if necessary. 3. Measure blood pressure again in complete peace and quiet.
Every measurement results in different values, although the device functions normally and normal values are displayed.	• Please read the following information and the points listed under "Common sources of error" (5.2). Repeat the measurement. • Blood pressure changes constantly. The observed readings may accurately reflect your pressure.
Blood pressure values differ from those measured by my doctor.	• Record daily values and consult your doctor. • Pressure readings in your doctor's office may be higher due to anxiety.
After the instrument has inflated the cuff the pressure falls very slowly, or not at all. (No reasonable measurement possible).	• Check cuff connections. Ensure the unit has not been tampered with.

8. Care and maintenance

- a) The cuff contains a sensitive air-tight bubble. Handle this carefully and avoid all types of straining through twisting or buckling.
- b) Clean the device with a soft, dry cloth. Do not use gas, thinners or similar solvents. Spots on the cuff can be removed carefully with a damp cloth and soapsuds. **Do not wash the cuff in a dishwasher, clothes washer, or submerge in water.**
- c) Handle the tube carefully. Do not pull on it. Do not allow the tubing to kink and keep it away from sharp edges.
- d) **Never open the monitor.** This voids the manufacturer's warranty.

9. Limited warranty

Your Automatic Blood Pressure Monitor is warrantied for the estimated life of the monitor by ProCheck, against manufacturer defects for the original purchaser only, from date of purchase.

The lifetime warranty applies to the monitor only. The following accessories are warranted for 1 year: cuff, AC adapter.

Batteries are not covered by this warranty. There are no user serviceable parts inside.

The warranty does not apply to consequential and incidental damages, or damage caused by batteries, improper handling, and accidents. Professional use, not following the operating instructions, and alterations made to the monitor or accessory by third parties, are also not included in this warranty. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

ProCheck will investigate your concern. A monitor or accessory as defined by this warranty, determined to be out of specification, will be replaced and shipped to you at no cost. A monitor or accessory as defined by this warranty, determined to be within specification, will be returned to you with a report of findings, at no cost.

Please use the below customer service contact information to reach ProCheck regarding any warranty concerns. We ask that you please contact us before sending any product back in order to better identify, and more quickly process, your concern.

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

10. Technical specifications

Weight:	415 g (with batteries)
Size:	157.5 (W) x 105 (L) x 61.5 (H) mm
Storage temperature:	-20° C to +55° C (-4° F to +131° F)
Humidity:	15 to 90% relative humidity maximum
Operation temperature:	10° C to 40° C (50° F - 104° F)
Display:	LCD-Display (Liquid Crystal Display)
Measuring method:	Oscillometric
Pressure sensor:	Semiconductor
Measuring range:	
SYS:	60 to 255 mmHg
DIA:	40 to 200 mmHg
Pulse:	40 to 199 per minute
Cuff pressure display range:	0–299 mmHg
Memory:	Automatically stores the last 99 measurements for 2 users
Measuring resolution:	1 mmHg
Accuracy:	pressure within ± 3 mmHg or 2% of reading > 200 mmHg pulse ± 5 % of the reading
Power source:	a) 4 AA batteries, 1.5 V b) AC adapter 6 V DC 600 mA (voltage 4.5 V DC to 6 V DC)
Accessories:	Wide-range cuff for arm circumference 8.7" - 16.5" (22 - 42 cm)

Technical alterations reserved.

Monitoring Device Assembled In Vietnam; Cuff Made In China

11. How to contact us

Distributed by: ProCheck
1617 Gulf to Bay Blvd.
Clearwater, FL 33755

Toll Free Customer Support: 1-866-289-6611

Trademark Usage:

Apple, the Apple logo, iPad, and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

Android and Google Play are both trademarks of Google Inc.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks in this Blood Pressure Monitor is under license. Other trademarks and trade names are those of their respective owners.

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Android y Google Play son marcas comerciales de Google Inc.