


TANITA[®]

Instruction Manual

Body Composition Monitor

BC-401F

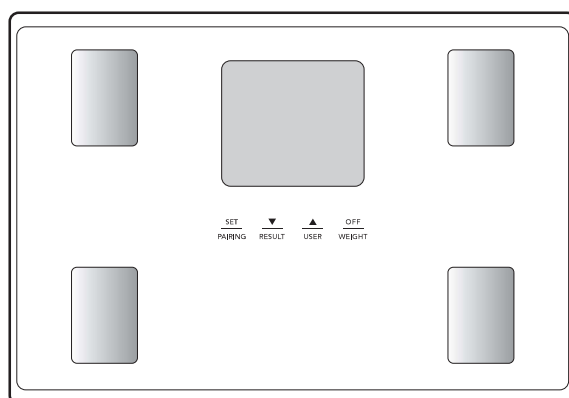

Images in this manual are for illustrative purposes only.

Read this Instruction Manual carefully and keep it for future reference.

Table of Contents

Introduction	3
For Your Safety	4
Features and Functions	6
Preparations before Use	6
Changing the Weight Mode	7
Automatic Shut-down Function	7
Getting Accurate Readings	7
Initial Setup (Type 1) via Smartphone	8
Initial Setup (Type 2) from Body Composition Monitor	10
Taking Measurements	11
Setting and Storing Personal Data in Memory	12
Taking Weight and Body Fat Readings	14
Programming the Guest Mode	16
Taking Weight Measurement Only	16
Changing Personal Data	16
Erasing Personal Data	16
Troubleshooting	17

Introduction

Thank you for purchasing this TANITA body composition monitor.
 This monitor is one of a wide range of home health care products produced by TANITA.
 This Instruction Manual will guide you through the setup procedures and outline the key features of this monitor. Additional information on healthy living can be found on our website www.tanita.com
 TANITA products incorporate the latest clinical research and technological innovations.
 All data is monitored by the TANITA Medical Advisory Board to ensure accuracy.



The TANITA body Composition Monitor is intended for use in the measurement of weight and impedance, and the estimation of body mass index (BMI), body fat percentage, total body water percentage, muscle mass, physique rating, bone mass, visceral fat rating with healthy range, basal metabolic rate (BMR), and metabolic age.

Safety Precautions

Do not use the body composition monitor if you have an electronic medical implant such as a pacemaker, as it passes a low-level electrical signal through the body which may interfere with the operation of the implant.
 Pregnant women should only use the weight function. Other functions are not intended for use when pregnant.

Important Notes for Users

This body composition monitor is intended for adults and children (age 7-17) who have inactive to moderately active lifestyles, and adults with athletic lifestyles.

Thanks to major investment in the latest BIA Technology and sports science research, TANITA has upgraded the Athlete Mode function to make measurements more accurate and suitable for a wider range of users.

Who should use Athlete Mode?

Adults aged 18 years and over who either:

- Train or exercise for 12 hours or more a week and have been doing so for at least six months.
- Are body builders.
- Are professional athletes who want to monitor their progress at home.
- Have a lifetime history of fitness and who used to do more than 12 hours a week but do less now.

Recorded data may be lost if the unit is used incorrectly or is exposed to electrical power surges.

TANITA cannot accept any responsibility for losses incurred due to the loss of recorded data.

TANITA cannot accept any responsibility for damages or losses incurred through the use of this product or any third-party claims.

Note: Body fat percentage estimates will vary with the amount of water in the body, and can be affected by dehydration or over hydration due to such factors as alcohol consumption, menstruation, illness, intense exercise, etc.

For Your Safety

This section explains precautionary measures to be taken to avoid injury to the users of this device and others, and to prevent damage to property. Please familiarize yourself with this information to ensure safe operation of this equipment.

	Warning	Failure to follow instructions highlighted with this mark could result in death or severe injury.
	Caution	Failure to follow instructions highlighted with this mark could result in injury or damage to property.
 Prohibited	This mark indicates actions that are prohibited.	
 Required	This mark indicates instructions that must always be followed.	

Warning

 Prohibited	This monitor must not be used by people with pacemakers or other medical implants. This monitor passes a weak electrical current through the body which could interfere with and cause the malfunction of electrical medical implants, resulting in serious harm to the user.
	Interpretation of measurements (e.g. evaluation of measurements and formulation of exercise programs based on the measurements) must be performed by a professional. Implementing weight loss measures and exercise programs based on self-analysis could be harmful to your health. Always follow the advice of a qualified professional.
	Do not step on the edge of the platform. Do not jump onto the monitor. Doing so may cause you to fall or slip, resulting in serious injury.
	Do not place this monitor on a slippery surface such as a wet floor. Doing so may cause you to fall or slip, resulting in serious injury. It may also cause water to get inside the monitor resulting in damage to internal components.
	Do not drop any objects onto the monitor. Avoid excessive impact or vibration to the monitor. Doing so may damage the glass, resulting in injury.

Caution

 Prohibited	Do not stand on the display screen.
	Do not insert your fingers into gaps or holes. Doing so may result in injury.
	Do not use rechargeable batteries. Do not use old batteries together with new batteries, or a mix of different types of batteries at the same time. Doing so may cause battery fluid to leak or the batteries to become excessively hot and rupture, resulting in damage to the monitor or injury.
	Never allow an infant or small child to carry the monitor. The child may drop the monitor, resulting in injury.
	This monitor is intended for home use only. This monitor is not intended for professional use including hospitals, medical or fitness facilities. It is not equipped with the components required for heavy usage. Using the monitor in a professional location will invalidate the warranty.



Required

Place the monitor on a hard, flat surface where there is minimal vibration to ensure safe and accurate measurement.

Avoid excessive impact or vibration to the monitor.

Doing so may damage the monitor or cause it to malfunction, or may cause settings to be erased.

Do not use or store this monitor in the following locations.

- Locations that are exposed to direct sunlight, near to heaters, or exposed directly to air from air conditioning units, etc.
- Locations where there are severe temperature changes or a large amount of vibration.
- Locations with high levels of humidity or moisture.

Doing so may damage the monitor or cause it to malfunction.

Do not step on the monitor when your body is wet (immediately after a shower or bath, etc.).

Be sure to thoroughly dry your body and feet before using the monitor. Failure to do so may cause you to slip or fall, resulting in injury. Also, the monitor cannot take accurate measurements if your body or feet are wet.

Keep away from water.

Do not use boiling water, benzene, thinner, etc. to clean the monitor.

If the monitor becomes dirty, soak a soft cloth in water or neutral household detergent, wring it out thoroughly and wipe the monitor clean, then use a dry cloth to wipe the monitor dry.

Do not attempt to disassemble the monitor. Do not modify this equipment in any way.

Do not wash the monitor in water.

Do not submerge this monitor in water.

Doing so may damage the monitor or cause it to malfunction.

Be sure to insert the batteries with the polarities (+, -) in the correct position.



If the batteries are incorrectly inserted, battery fluid may leak and damage the floor. (If you do not intend to use this monitor for a long period of time, it is advisable to remove the batteries before storage.)

Do not use the monitor near other products such as microwave ovens that emit electromagnetic waves.

Doing so may cause the monitor to malfunction, preventing accurate measurements.



Another person should assist persons with disabilities who may not be able to take a measurement alone.

Features and Functions

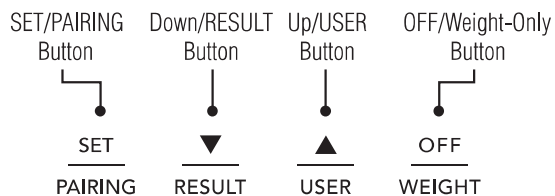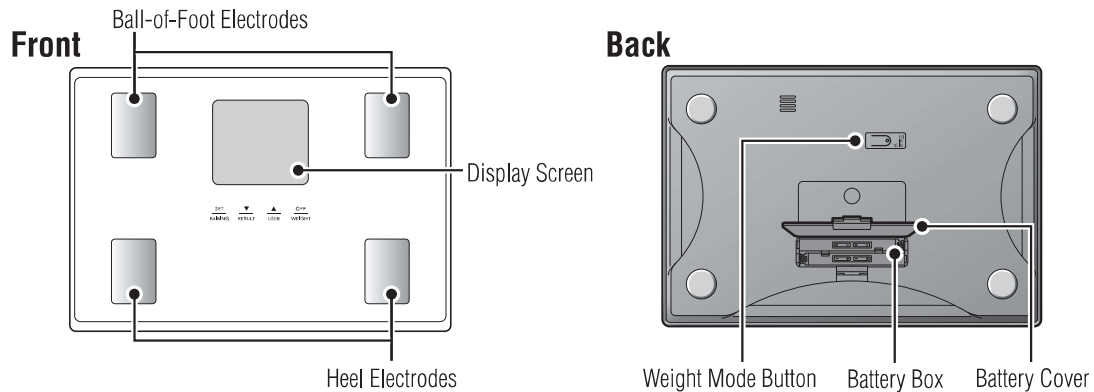

Accessories

Location Code Instruction Manual AAA-Size Batteries (4)

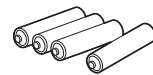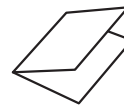


Preparations before Use

Inserting the Batteries

Open the battery cover on the back of the measuring platform. Insert the supplied AAA-batteries as indicated.

Note: Be sure that the polarity of the batteries is set properly. If the batteries are incorrectly positioned the fluid may leak and damage floors. If you do not intend to use this unit for a long period of time, it is advisable to remove the batteries before storage.

Please note that the included batteries from factory may have decreased energy levels over time.

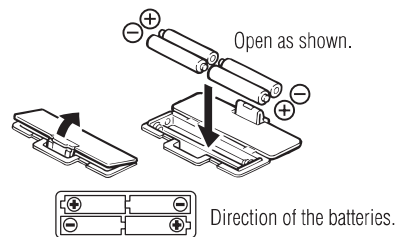

Positioning the Monitor

Place the measuring platform on a hard, flat surface where there is minimal vibration to ensure safe and accurate measurement.

Note: To avoid possible injury, do not step on the edge of the platform.



Handling Tips

This monitor is a precision instrument utilizing state-of-the-art technology. To keep the unit in the best condition, follow these instructions carefully:

- Do not attempt to disassemble the measuring platform.
- Store the unit horizontally, and place it so that the buttons will not be pressed accidentally.
- Avoid excessive impact or vibration to the unit.
- Place the unit in an area free from direct sunlight, heating equipment, high humidity, or extreme temperature change.
- Never submerge in water. Use alcohol to clean the electrodes and glass cleaner (applied to a cloth first) to keep them shiny; avoid soaps.
- Do not step on the platform when wet.
- Do not drop any objects onto the platform.
- Do not jump onto or off the unit.

Changing the Weight Mode

You can change the weight unit by using the button on the back of the scale.

Press the OFF button when the display is blank. When  appears in the screen, press the Weight Mode button on the underside of the unit to switch the unit in which weight is displayed.


Note: If weight mode is set to pounds or stone-pounds, the height programming mode will be automatically set to feet and inches. Similarly, if kilograms is selected, height will be automatically set to "cm".

Automatic Shut-Down Function


The automatic shut-down function shuts off the power automatically in the following cases:

- If you interrupt the measurement process. The power will shut down automatically within 30 seconds, depending upon the type of operation.
- During programming, if you do not touch any of the keys or buttons within 60 seconds.
- After you have completed the measuring process.

Getting Accurate Readings

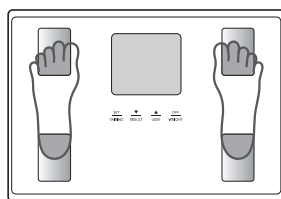

Heels centered on electrodes

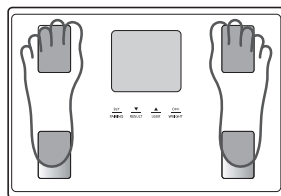

Toes may overhang measuring platform

To ensure accuracy, readings should be taken without clothing and under consistent conditions of hydration. If you do not undress, always remove your socks or stockings, and be sure the soles of your feet are clean before stepping on the measuring platform.

Be sure that your heels are correctly aligned with the electrodes on the measuring platform. Don't worry if your feet appear too large for the unit – accurate readings can still be obtained if your toes overhang the platform.

It is best to take readings at the same time of day. Try to wait about three hours after rising, eating, or hard exercise before taking measurements.

While readings taken under other conditions may not have the same absolute values, they are accurate for determining the percentage of change as long as the readings are taken in a consistent manner. To monitor progress, compare weight and body fat percentage taken under the same conditions over a period of time.

Note: An accurate reading will not be possible if the soles of your feet are not clean, or if your knees are bent or you are in a sitting position.

Initial Setup (Type 1) via Smartphone

1. Download the the app and turn Bluetooth® on.


Visit the App Store/Google Play Store to download the "Healthy Edge Mobile" app.
Turn the Bluetooth setting on.

- Note:**
- If you are not using the app, enter settings manually in the body composition monitor (pg. 10) before using.
 - If you purchase a new smartphone, be sure to adjust the registration of your body composition monitor accordingly.
 - Approximate signal range: 16 feet in areas without obstructions




2. Set your "Profile"

Open the app.
Open "Settings" and then register your profile details.
Touch "<Settings" to return to the previous screen.

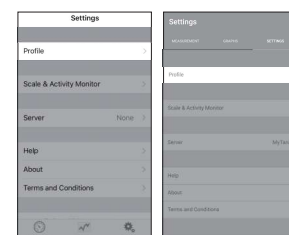


3. "Server"

Fill out your info registered at my.tanita.com in "Server" under "Settings".
Measurement results can be uploaded to my.tanita.com after measurement.
Touch "Save" to return to the previous screen.

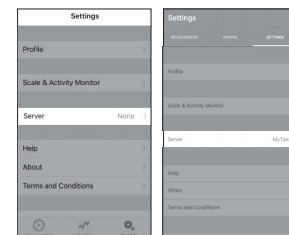

4. "Scale" settings

Touch "+" and follow the instructions in the app to register your platform.
Return to the previous screen.

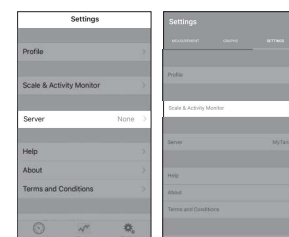


- Note:**
- Compatible apps may be changed or added, in which case screen displays and operation methods will vary accordingly. Contact www.tanita.com with any questions or concerns.
 - The screen displays and instructions may differ depending on the app specifications.


5. Operate the body composition monitor as indicated by the app.
If the body composition monitor is not powered on, press and hold SET(enter button) for 3 seconds or more, then enter your information as indicated by the app.

Note: • Depending on signal strength, it may take a few seconds to connect to the body composition monitor.

- The  icon is lit while data is being transmitted.
- To weigh yourself while holding a smartphone, first enter the amount of weight you want to subtract into "Weight of Clothing, etc." The input amount will be subtracted from your weight display.
- If 0.4 lb is entered into "Weight of Clothing, etc.," the following will be displayed before a measurement is taken.

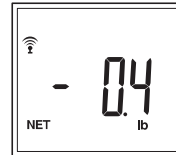


Display before measurement

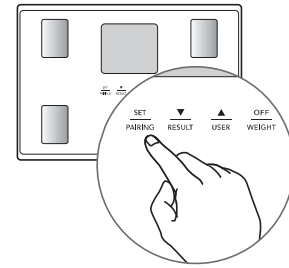

6. When the body composition monitor is activated,  is displayed.

7. After  is displayed, carefully step barefoot onto the measuring platform.

Note: Please refrain from moving during measurement.

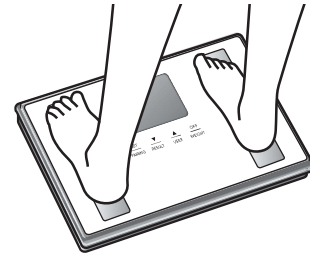

Registration of your body composition monitor is completed once measurement is completed.


Initial Setup (Type 2) from Body Composition Monitor

If the BC-401F is already paired with a smartphone, the location, date and time are already set so please skip this procedure. You may begin taking measurements right away.

Different locations across the world have different gravitational pulls. These differences affect the accuracy of the weight measurement. You can obtain the most accurate weight readings by selecting the correct gravity setting on the scale according to your geographical location.

1. Set your location code

When the power is off, press the SET, ▼ or ▲ button. Refer to the map enclosed, and press the ▼ or ▲ button to select your location code. Press the SET button to enter and save the setting.

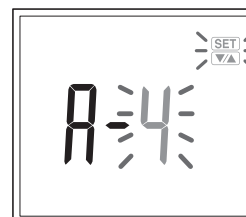


2. Set Date

Press the ▼ or ▲ button to set year and then press the SET button. Set the month and date in the same way.




3. Set Hour

Press the ▼ or ▲ button to set hour and then press the SET button.

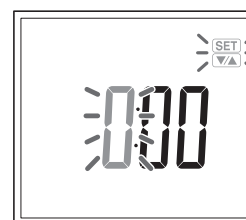

4. Set Minute

Press the ▼ or ▲ button to set minute and then press the SET button.

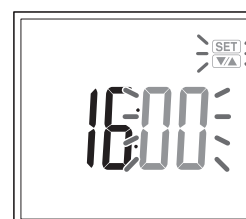

Note: 60 seconds or more of inactivity will cause automatic shut off during this process. No data will be saved as a result.

If you make a mistake or wish to turn the unit off before you complete the programming, press the OFF button to turn the unit off. Settings will not be saved in this case.

To change your location code, press and hold the OFF button for 3 seconds when the power is off, and then follow the procedure above to change your location code.

To change the date or time settings, press and hold the ▼ button for 3 seconds when the power is off, and then follow the procedure above to change the date or time settings.

You can fast-forward numbers if you keep holding down the ▼ or ▲ button.

When you replace batteries, press the SET, ▼ or ▲ button and then follow the above steps 2 to 4.

Taking Measurements

- When the monitor power is off, touch  and then follow the instructions in the app. When  (0.0 kg, 0.0 stlb) is displayed, step onto the monitor.




- Measurement data is automatically sent to the smartphone.

Note: When taking measurements while holding your smartphone, enter the weight that you want to subtract into the "Smartphone's weight" section in the dedicated app before taking measurements.

The weight that you entered to be subtracted from your weight is displayed.

The "Smartphone's weight" can be entered up to -4.4 lb (-2.0 kg, -4 lb).

The maximum weight (weighing capacity) is reduced by the amount entered for the "Smartphone's weight".

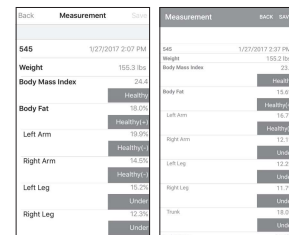

If you do not have your smartphone when taking measurements.


Simply step onto the monitor. The auto-recognition feature will automatically recall your personal data. Your monitor will store up to 10 sets of measurements.

Saved data is sent to the smartphone together with new measurements the next time the smartphone is connected to the monitor when taking measurements.

Note: Up to 10 measurements can be stored in the device per ID number.

Each new data entered after this amount causes the oldest data to be deleted.

 will notify you if pairing is complete and there is unsent data.

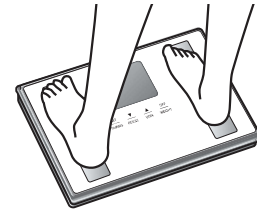

Setting and Storing Personal Data in Memory

The unit can be operated only if data has been programmed into one of the personal data memories.

1. Turn on the Power

Press the SET button to turn on the unit.

The unit will beep to confirm activation, the Personal data number will be displayed, and the display will flash.

Note: If you don't operate the unit for 60 seconds after the unit has been turned on, the unit will turn off automatically.
If you make a mistake or want to turn the unit off before you have finished programming it, press the OFF button to force quit.

2. Select a Personal data Number


Press the ▼ or ▲ button to select a Personal data number (1, 2, 3, 4 and 5). Once you reach the Personal data number you wish to use, press the SET button. The unit will beep once to confirm.



3. Set your Birthday

Set the date of your birthday (Year/Month/Day) using the ▼ or ▲ button. Press the SET button to confirm at each programming. The unit will beep once to confirm.

Note: Year settings can be configured from 1900 onwards.

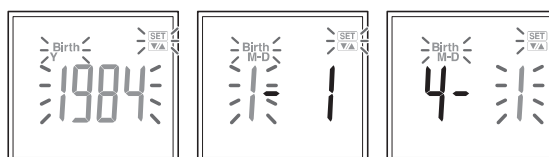

4. Select Female or Male

Use the ▼ or ▲ button to scroll through Female (♀) or Male (♂), Female/Athlete (♀🏃) and Male/Athlete (♂🏃) settings, then press the SET button. The unit will beep once to confirm.

Who should use Athlete Mode?

Adults aged 18 years and over who either:

- Train or exercise for 12 hours or more a week and have been doing so for at least six months.
- Are body builders.
- Are professional athletes who want to monitor their progress at home.
- Have a lifetime history of fitness and who used to do more than 12 hours a week but do less now.



5. Specify the Height

The display defaults to 5' 7.0" (170 cm) (range of height is from 3'4.0" –7'3.0" or 100 cm – 220 cm). Use the ▼ or ▲ button to specify Height and then press the SET button. The unit will beep once to confirm, and then the display shows  and is ready for measurement.



6. Measuring your weight and body fat

Step onto the scale after  is shown on the display within 60 seconds.

Do not step off until "0000" disappears.

The readings will be shown for 40 seconds.

The unit will then shut off and the readings will be stored in the memory.



Note: If you step onto the scale before  appears the display will show *Err* and you will not obtain a reading. Furthermore, if you do not step onto the measuring scale within about 60 seconds after  appears, the power is shut-off automatically. Please try the setting and storing personal data process again.

Warning

Do not use the body fat reading feature of this product if you have a pacemaker or other implanted electronic medical device.


A diagram of a hand with a sensor on the palm and a central display unit. The sensor is a rectangular block on the palm, and the display unit is a rectangular block in the center. Below the display unit are four small icons: a person, a person with a plus sign, a person with a minus sign, and a person with a plus/minus sign.

C
M
Y
CM
MY
CY
CMY
K


1. Weight




2. BMI

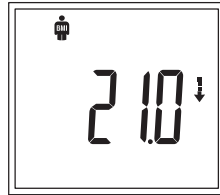


3. Body Fat %

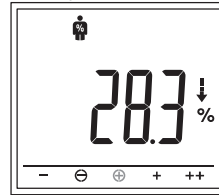


4. Muscle mass

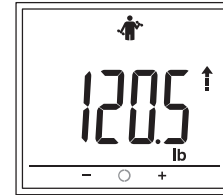


5. Physique Rating

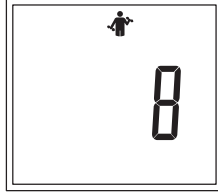


6. Bone mass

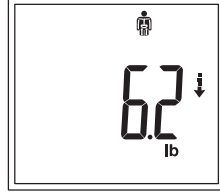


7. Visceral fat

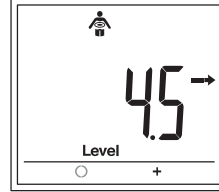


8. BMR

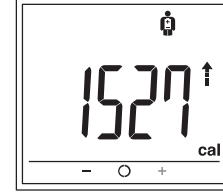

9. Metabolic Age

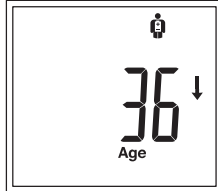

10. Total Body Water %

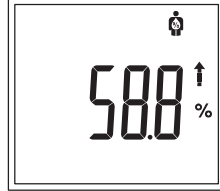

	-	⊖	○	⊕	+	++
Body Fat Percentage Assessment	Under	Standard -		Standard +	Over	Obese
Muscle Mass Assessment	Under		Average		High	
Visceral Fat Level Assessment			Average		High	
BMR Assessment	Low		Average		High	

↑ → ↓ Arrows indicate the trend compared to the previous time that a measurement was taken.

Note: For children (age 7–17), the monitor only displays weight, BMI and Body Fat %.

For children (age 6 or less), the monitor only displays the weight.

"Err" is displayed if the body fat percentage exceeds the measurable range.

You can use the ▼ button to scroll through the results.

Press and hold the ▼ button for fast scrolling. After reading the results that you want to see, release the button to see the rest of the results.

Press the OFF button to cancel measurement or turn the power off.


Do not move the monitor for approximately 3 seconds after taking measurements and stepping off.

To prevent accidents such as falling off the monitor, always step off the monitor before operating the buttons.

Programming the Guest Mode

Guest mode allows you to program the unit for a one-time use without reading a personal data number.

Press the $\blacktriangle$ button to activate the unit. Press the $\blacktriangledown$ or $\blacktriangle$ button to select $-0-$. Press the SET button. Then follow the steps: Set the age, gender (male/female/male athlete/female athlete) and height. Step onto the scale after 00_{16} is shown on the display. Your weight will be shown first, continue standing on the scale until the body fat % reading is shown on the display. The reading will be shown for 40 seconds.

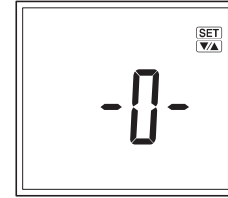

Note: Do not step off until “0000” disappears.

Note: If you step onto the scale before 00_{16} appears the display will show E_{rr} and you will not obtain readings.

Taking Weight Measurement Only

Press the OFF button. Step onto the scale after 00_{16} is shown on the display. The body weight value will appear. The display will light for about 20 seconds, and then the power will automatically turn off.


Changing Personal Data

1. When the power is off, press the SET button.
2. Press the $\blacktriangledown$ or $\blacktriangle$ button to select a personal data number (1, 2, 3, 4 and 5). Once you reach the personal data number you wish to use, press the SET button. The unit will beep once to confirm.
3. Follow the steps as listed on P.12 to change your personal data.

Note: If you change the birthday and/or gender, the unit will show E_{Lr} .


If you select YES using the $\blacktriangledown$ or $\blacktriangle$ button and pressing the SET button to confirm, the existing personal data will be deleted.

If you select no using the $\blacktriangledown$ or $\blacktriangle$ button and pressing the SET button to confirm, no new data will be saved (existing data will be unchanged) then the unit will shut off.

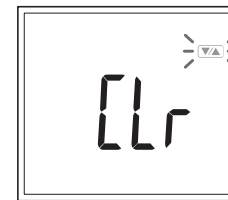

Erasing Personal Data

1. When the power is off, press and hold the SET button.
2. Press the $\blacktriangledown$ or $\blacktriangle$ button to select a personal data number (1, 2, 3, 4 and 5). Once you reach the personal data number you wish to use, press the SET button for more than 3 seconds. The unit will beep once to confirm. E_{Lr} is displayed.
3. Select YES using the $\blacktriangledown$ or $\blacktriangle$ button.
4. Press the SET button.

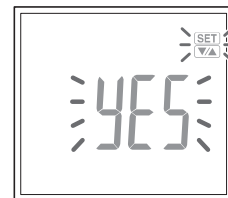







Troubleshooting







If the following problems occur ...then...

- **A wrong weight format appears followed by kg, lb, or st-lb.**
Refer to "Changing the Weight Mode" on P.7.
- **Lo appears on the display, or all the data appears and immediately disappears.**
Batteries are low. When this message appears, be sure to replace the batteries immediately, since weak batteries will affect the accuracy of your measurements. Change all the batteries at the same time with new AAA-size batteries.
Note: Your settings will not be erased from the memory when you remove the batteries.
- **Err appears while measuring.**
Please stand on the platform keeping movement to a minimum. The unit cannot accurately measure your weight if it detects movement.
When the Body Fat % is more than 75%, readings cannot be obtained from the unit.
When the BMI is less than 5, readings cannot be obtained from the unit.
- **The Body Fat Percentage measurement does not appear.**
Your personal data has not been programmed. Follow the steps on P.12.
Make sure socks or stockings are removed, and the soles of your feet are clean and properly aligned with the guides on the measuring platform.
- **OL appears while measuring.**
Reading cannot be obtained if the weight capacity is exceeded.
- **-OL appears on the display.**
Did you step onto the body composition monitor before  was displayed?
- **Errb is displayed.**
Please set the correct date and time.
When not in use, the body composition monitor will periodically calibrate in preparation for the next weighing. If the unit has been stored horizontally, with items placed on top or moved from one location to another, select your personal data number as explained on p.12 prior to taking your measurement.

Cannot connect my smartphone to the body composition monitor

- **ErrP or Errd appears on the display**
Are batteries inserted correctly in the body composition monitor?
Are the batteries running low?
Are the smartphone and body composition monitor within communication range?
The communication range for the body composition monitor is 5 m in an unobstructed location.
Is the "Bluetooth" setting in the smartphone "Settings" turned on?
- **ErrP appears repeatedly.**
Delete the pairing information from "Settings" > "Bluetooth" on the smartphone, and pair the smartphone with the body composition monitor again.
- **ErrU appears repeatedly.**
Delete the pairing information from "Settings" > "Bluetooth" on the smartphone, and pair the smartphone with the body composition monitor again.
- **Errg appears on the display**
Remove and replace the batteries from the body composition monitor, and set the date and time again.
- **The date and time settings have been changed**
When using the body composition monitor linked with a smartphone, the date and time settings on the body composition monitor are overwritten with those from the smartphone.
- **The "Input Data" button in the app does not work**
Is the "Bluetooth" setting in the smartphone "Settings" turned on?
Linking between the app and the body composition monitor is only enabled when the body composition monitor is turned off.
Check that the power of the body composition monitor is turned off and try again.

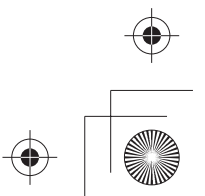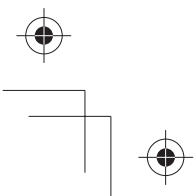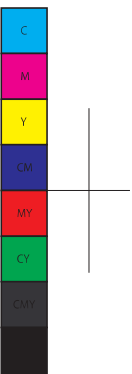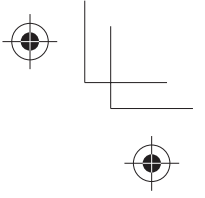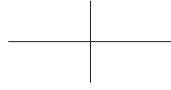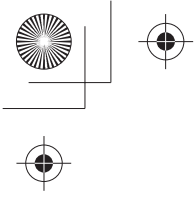

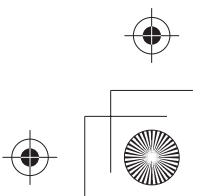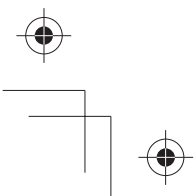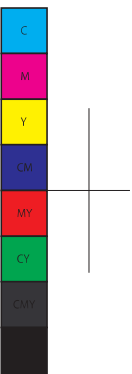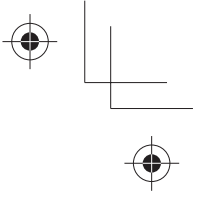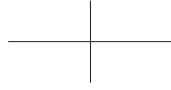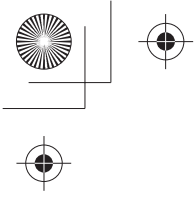

Specifications

BC-401F

Weight Capacity

330 lb (150 kg) (23 st 8 lb)

Weight Increments

0 – 330 lb/0.2 lb (0 – 150 kg/0.1 kg) (0 – 23 st 8 lb/0.2 lb)

Body Fat %

7 – 99 years old*

Body Mass Index

7 – 99 years old

Muscle Mass

18 – 99 years old*

Physique Rating

18 – 99 years old

Bone Mass

18 – 99 years old

Visceral Fat

18 – 99 years old*

Basal Metabolic Rate

18 – 99 years old*

Metabolic Age

18 – 99 years old

Total Body Water %

18 – 99 years old

Personal Data

5 memories

Power Supply

4AAA alkaline batteries (included)

Power Consumption

50 mA maximum

Measuring Current50 kHz, 5 kHz, 100 μ A**Communication Method**

Bluetooth Low Energy technology

Frequency Band

2402 – 2480 MHz

* When age exceeds 99 years old, it is just for reference.

USA**Federal Communications Commission**

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio or television technician for help.

Modifications

The FCC requires the user to be notified that any changes or modifications made to this device that are not expressly approved by Tanita Corporation may void the user's authority to operate the equipment.

- The product design and specifications may be changed at any time without prior notice.
- Apple, the Apple logo, and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.
- Google Play and the Google Play logo are trademarks of Google LLC.
- Android is trademark of Google LLC.
- The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by TANITA Corporation is under license.

Manufacturer**TANITA Corporation**

1-14-2 Maeno-cho, Itabashi-ku,
Tokyo 174-8630 Japan
TEL: +81-(0)3-3968-7048
www.tanita.co.jp

EU Representative**TANITA Europe B.V.**

Hoogoorddreef 56-E, 1101 BE
Amsterdam, the Netherlands
TEL: +31-(0)20-560-2970
FAX: +31-(0)20-560-2988
www.tanita.eu

TANITA Corporation of America, Inc.

2625 South Clearbrook Drive,
Arlington Heights, Illinois 60005
U.S.A.
TEL: +1-847-640-9241
FAX: +1-847-640-9261
www.tanita.com

TANITA Health Equipment H.K. Ltd.

Unit 301-303, Wing On Plaza, 3/F., 62 Mody
Road,
Tsimshatsui East, Kowloon, Hong Kong
TEL: +852-2834-3917
FAX: +852-2838-8667
www.tanita.asia

TANITA (Shanghai) Trading Co., Ltd.

Room 8005, 877 Huai Hai Zhong Lu,
Shanghai,
The People's Republic of China
TEL: +86-21-6474-6803
FAX: +86-21-6474-7901
www.tanita.com.cn