



Ecolab Hand Hygiene Program Compliance Monitoring System

HHCM 915 Bed Beacon 92053072
Installation and Users' Guide

Table of Contents

1. Introduction to the Bed Beacon and Mat Antenna.....	3
2. Required Materials.....	5
2.1 Installation Materials	5
2.2 Bed Beacon Components.....	6
3. Bed Beacon Installation Functions.....	6
3.1 Bed Mode (Default) (Internal antenna only)	8
3.2 Bed Mode, Using External Antenna Only.....	8
3.3 Bed Beacon, Using Internal and External Antenna Mode (toggling back and forth between the internal and external antenna).....	9
3.4 Stretcher Mode, using Internal Antenna Only	9
3.5 Stretcher Mode, using External Antenna Only	9
3.6 Stretcher Mode Using Internal and External Antenna Mode (toggling back and forth between the internal and external antenna)	10
4. Installation.....	10
4.1 Mounting Bracket 9205-2337 Installation - Patient Bed	10
4.2 Bed Beacon Installation	11
4.3 AC Line Power Sensor.....	12
5. Cable Management – Hill Rom VersaCare P3200	12
6. Adjusting the Patient Zone	14
6.1 Range Adjustment on Beds.....	14
6.2 Test Badge Range Adjustment	14
6.3 Ecolab Dashboard Test with the Network Badge and Bed Beacon	15
7. Metadata Entry	17
7.1 New Bed Beacon Metadata Entry.....	17
7.2 Replacing an Existing Bed Beacon – Swapping Metadata	18
8. Troubleshooting.....	19
8.1 No Bed Beacon Heartbeats on the Dashboard.....	19
8.2 The Bed Installation Does Not Communicate with a Badge	19

9. System Component Care and Maintenance	20
9.1 Cleaning the Components	20
9.2 Handling the Bed Beacon	20
9.3 Replacement Batteries	20
Appendix A - Certification and Safety Approvals	20
FCC Statement	20
Industry Canada	21
Appendix C – Bed Label	21
Appendix D- Z Bracket (9205-2236)	23
Appendix E - Patient Bed Hill-Rom Advanta Installation	24
1. Installation	25
1.1 Step 1: Bed Mode Configuration	25
1.2 Step 2: Installing the Bed Beacon (92053072)	25
1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord	25
1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range	26
2. Cable Management	26
2.1 Step 5: AC Line Power Sensor Cable Routing	26
3. Step 6: Final Test	27
3.1 Disconnect and Reconnect the Bed's AC Power Cord	27
4. Bed Installation Checklist	27
Appendix F - Patient Bed Hill-Rom 1000 & Advanta 2 Installation	28
1. Installation	29
1.1 Step 1: Bed Mode Configuration	29
1.2 Step 2: Installing the Bed Beacon (92053072)	30
1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord	30
1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range	31
2. Cable Management	31
2.1 Step 5: AC Line Power Sensor Cable Routing	31
3. Step 6: Final Test	32
3.1 Disconnect and Reconnect the Bed's AC Power Cord	32
4. Bed Installation Checklist	32
Appendix G - Patient Bed Stryker Installation	33
1. Installation	34
1.1 Step 1: Bed Mode Configuration	34
1.2 Step 2: Installing the Bed Beacon (92053072)	34
1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord	35
1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range	35
2. Cable Management	36

2.1 Step 5: AC Line Power Sensor Cable Routing.....	36
3. Step 6: Final Test	37
3.1 Disconnect and Reconnect the Bed's AC Power Cord	37
4. Bed Installation Checklist.....	37
Appendix H - Patient Bed Hill-Rom TotalCare Installation	38
1. Installation	39
1.1 Step 1: Bed Mode Configuration	39
1.2 Step 2: Installing the Bed Beacon (92053072).....	39
1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord.....	40
1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range	40
2. Cable Management.....	40
2.1 Step 5: AC Line Power Sensor Cable Routing.....	40
3. Step 6: Final Test	41
3.1 Disconnect and Reconnect the Bed's AC Power Cord	41
4. Bed Installation Checklist.....	42
Appendix I - Patient Bed Hill-Rom Affinity 2 Installation	43
1. Installation	44
1.1 Step 1: Bed Mode Configuration	44
1.2 Step 2: Installing the Bed Beacon (92053072).....	44
1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord.....	44
1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range	45
2. Cable Management.....	45
2.1 Step 5: AC Line Power Sensor Cable Routing.....	45
3. Step 6: Final Test	46
3.1 Disconnect and Reconnect the Bed's AC Power Cord	46
4. Bed Installation Checklist.....	46
Appendix J – Stryker Renaissance Stretcher.....	47
1. Installation	48
1.1 Step 1: Bed Mode Configuration.....	48
1.2 Step 2: Installing the Bed Beacon (92053072).....	49
1.4 Step 3: Adjusting the Bed Beacon to Proper Field Range	49
2. Stretcher Installation Checklist.....	49


1. Introduction to the Bed Beacon and Mat Antenna

In full hospital implementations of the Ecolab™ Hand Hygiene Program Compliance Monitoring System (the **System**), the **Bed Beacon**, which contains an internal low frequency antenna, is mounted to a patient bed. Its function is to: communicate with the system **Badge** worn by a HCW, when the HCW comes into close proximity with the patient bed (“bed event”); collect bed event and hand hygiene status data from that **Badge**; and then

transmit the collected data to the **Ecolab Proprietary Wireless Network**, where it will ultimately be sent to an offsite server for processing and archiving. Similarly, the “bed beacon” may be attached to stretchers, cribs, infusion chairs, etc.: anywhere you wish to establish the requirement for hand hygiene before and after patient contact. For simplicity, the examples of beds and stretchers will be used in this document, but the principals may be applied to installation around other patient-centered areas.




Ecolab Hand Hygiene Compliance Monitoring System: Bed Beacon, Mounting Bracket and AC Line Power Sensor

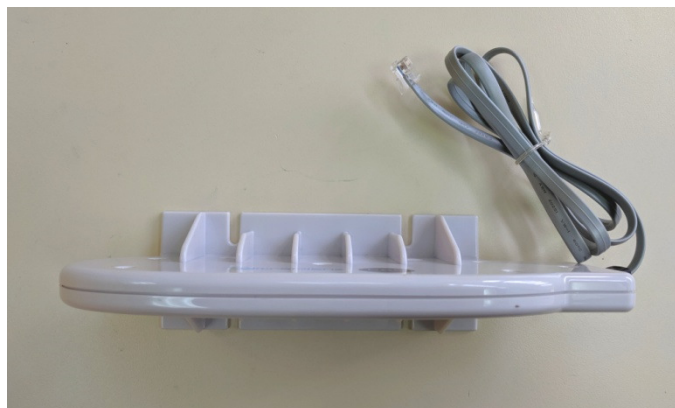

Ecolab Bed Beacon External Antenna (P/N: 92052232)

2. Required Materials

2.1 Installation Materials

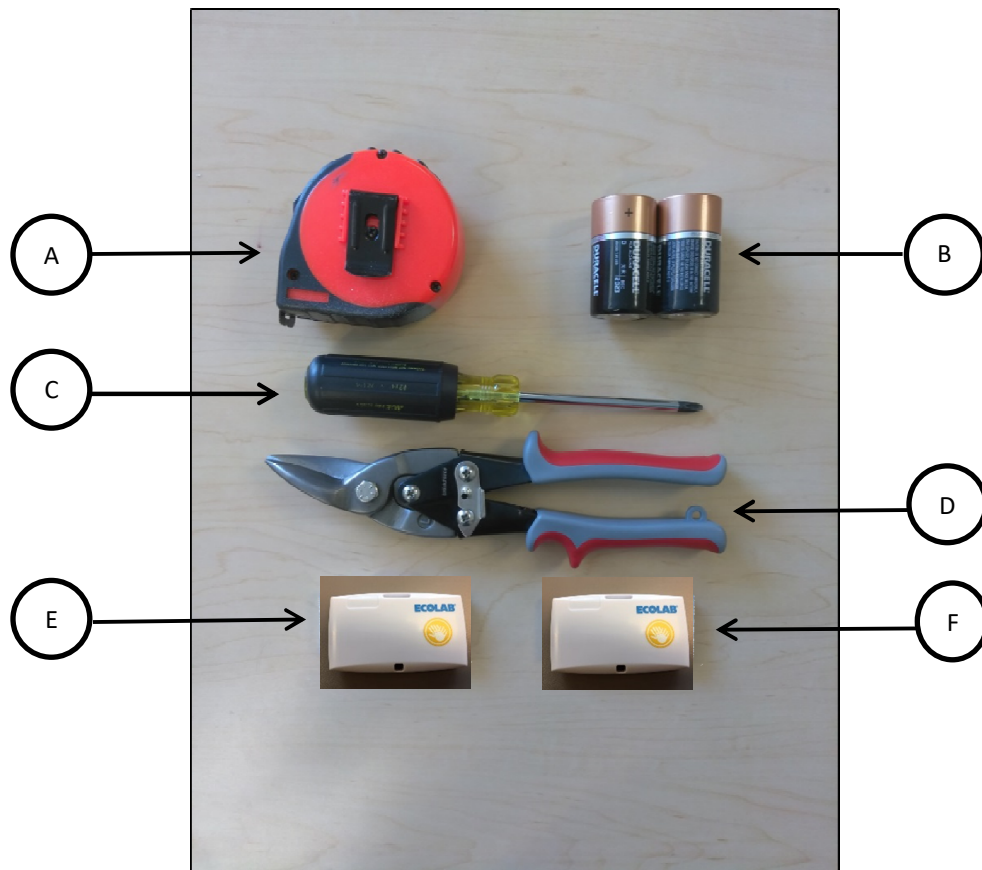

(A) Tape Measure	(B) Duracell Coppertop Alkaline D cell batteries (2 ea.)
(C) Phillips screwdriver	(D) Wire cutters
(E) Range Test Badge with Network	(F) Range Test Badge Without Network

Items not pictured:

(G) Cable Tie 4in long (10ea.)
(H) Cable Tie Holder (10ea.)
(I) Cable Tie 18in long (4ea.)
(B) Alcohol Wipe
(C) Scotch-Brite™ (or similar) scouring pad
(D) Spray bottle with a mixture of liquid dish soap and water (1 oz. soap per gallon of water)

2.2 Bed Beacon Components

The patient bed installation is comprised of four main components. The first component is the **Bed Beacon**, which is used to create and control the **Patient Zone** around the patient bed. The second component is the **AC Line Power Sensor** that senses when the bed is connected (turns the **Patient Zone** on) or disconnected (turns the **Patient Zone** off) from AC power. The third component is an **External Antenna** that will be mounted on the Patient Bed and create another Patient Zone around the antenna. This antenna provided by Ecolab is specific to this equipment and may not be replaced by any other antenna models. The last component is the quick release **Mounting Bracket** clip for attaching the **Bed Beacon** to a bed or stretcher.



Figure 1. *Bed Beacon Installation Components*

3. Bed Beacon Installation Functions

The **Bed Beacon** creates a modulated magnetic field around a patient bed, referred to as the **Patient Zone**. When a healthcare worker (HCW) wearing a System **Badge** enters the **Patient Zone**, the **Bed Beacon** wirelessly communicates with the **Badge** and signals to the **Badge** that it has entered the **Patient Zone** around a specific patient bed. The **Bed Beacon** will emit an audible beep and flash its status LED (Light Emitting Diode) green when it has successfully communicated with a **Badge**. If communication with the **Badge** is successful, the **Bed Beacon** sends information about its interaction with the **Badge** to the server via the System's proprietary wireless network. For more information about how the **Badge** will behave when interacting with the **Patient Zone**, please read the document entitled "*Ecolab™ Hand Hygiene Program Compliance Monitoring System Healthcare Worker Badge User's Guide*".

To verify which mode the **Bed Beacon** is in, press and hold both range buttons down. Immediately, the **Beacon** will flash three (3) times the green LED. With each flash, there will be an audible alert (beep). If the range buttons are released during or immediately after the three green flashes, the bed beacon is set to the default range (10 = maximum). One (1) second later, the **Bed Beacon** will flash red and beep once. The single red LED flash and audible alert corresponds to the **Bed Beacon** being set for mode number 1: bed mode. If both range buttons continue to be held down, the bed beacon will continue to cycle through the operation modes, at a rate of 1 every 500 ms. The number of red LED flashes and audible beeps corresponds to the operation mode (see list below.) **In order to select the mode, release the range buttons after the desired mode is reached.**

An alternative method of setting the bed mode is with the Ecolab **Beacon Installation Tool** (9250-3074). For more information about how the **Beacon Installation Tool** can change the **Beacon** mode, please read the document entitled, “*Ecolab™ Hand Hygiene Program Compliance Monitoring System Beacon Installation Tool User’s Guide*”.

Operation Modes:

- 1 = Bed Mode using Internal Antenna Only
- 2 = Bed Mode using External Antenna Only
- 3 = Bed Mode using Internal and External Antenna (toggle between internal and external antenna)
- 4 = Stretcher Mode using Internal Antenna Only
- 5 = Stretcher Mode using External Antenna Only
- 6 = Stretcher Mode using Internal and External Antenna (toggle between internal and external antenna)

Note:

Bed Modes (1-3) - AC Power Cord Sensor must detect AC power to broadcast the 125kHz signal.

Stretcher Modes (4-6) – the 125kHz signal is always broadcast. No AC Power Cord Sensor required.

The default operation mode is 1. To select a different mode, press and hold both range buttons simultaneously. The LEDs will flash and the piezo will beep to indicate the mode.

***	3 green flashes and beeps – default NFC range, which is 10 (maximum)
*	1 red flash and beep – operation mode 1
**	2 red flashes and beeps – operation mode 2
***	3 red flashes and beeps – operation mode 3
****	4 red flashes and beeps – operation mode 4
*****	5 red flashes and beeps – operation mode 5
*****	6 red flashes and beeps – operation mode 6

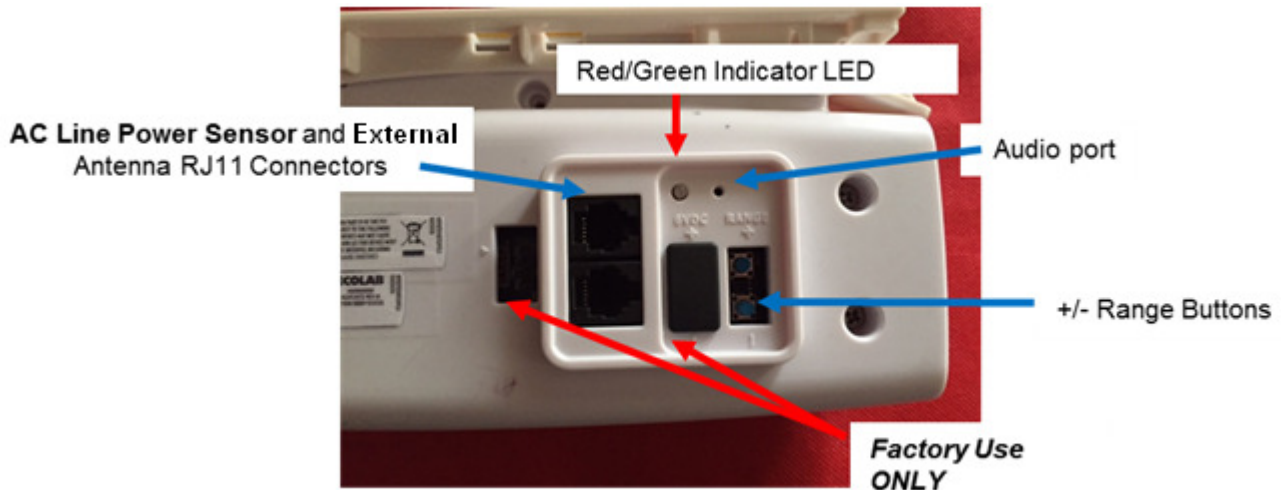


Figure 2. Bed Beacon Details

3.1 Bed Mode (Default) (Internal antenna only)

* 1 red flash and beep – operation mode 1

In “Bed Mode”, the **Bed Beacon** is mounted on a patient bed and connected to an **AC Line Power Sensor** through either RJ11 connector. The **Bed Beacon** is then able to detect when the bed's AC power cord is connected to a wall outlet, causing the **Bed Beacon** to emit a beep and the LED to flash green. The state of the AC Line Power Sensor is included in the information sent when the **Bed Beacon** sends out a heartbeat message over the System's proprietary wireless network. If the bed's power cord is disconnected from AC wall power, the **Patient Zone** will be turned off and will remain off while the bed is in transit. When the **AC Line Power (ALP) Sensor** senses that the power is disconnected, the **Bed Beacon** will emit a beep and flash the red LED. When the bed's power cord is plugged in again, the **Bed Beacon** will emit a beep and flash the green LED (may take up to 3S for the Bed Beacon to acknowledge the change in state).

3.2 Bed Mode, Using External Antenna Only

** 2 red flashes and beeps – operation mode 2

To set the mode, press and hold both range buttons down. Immediately, the **Beacon** will flash the green LED three (3) times. With each flash, there will be an audible alert (beep). Continue to hold the range buttons down while the bed beacon continues to cycle through the operation modes, until you reach two (2) red LED flashes / 2 audible beeps. Release the range buttons to set.

In this mode, the patient zone is created by the External Antenna only. The Bed Beacon itself is used to communicate with badges and the System network hubs. In “Bed Mode, External Antenna Only”, the **Bed Beacon** is mounted on a patient bed and connected to the **AC Line Power Sensor** through either RJ11 connector. The **External Antenna** is connected to the Bed Beacon using the remaining RJ11 connector and is then also attached to the bed, in a central location, as it will generate the patient zone. The **External Antenna** (Ecolab P/N: 92052232) provided by the Ecolab is specific to this equipment and may not be replaced by any other antenna model. The **Bed Beacon** is then able to detect when the bed's AC power cord is connected to a wall outlet, causing the **Bed Beacon** to emit a beep and the LED to flash green. If the

bed's power cord is disconnected from AC wall power, the **Patient Zone** will be turned off and will remain off while the bed is in transit. When the **AC Line Power (ALP) Sensor** senses that the power is disconnected, the **Bed Beacon** will emit a beep and flash the red LED. When the bed's power cord is plugged in again, the **Bed Beacon** will emit a beep and flash the green LED.

3.3 Bed Beacon, Using Internal and External Antenna Mode (toggling back and forth between the internal and external antenna)

*** 3 red flashes and beeps – operation mode 3

To set the mode, press and hold both range buttons down. Immediately, the **Beacon** will flash three (3) times the green LED. With each flash, there will be an audible alert (beep). Continue to hold the range buttons down while the bed beacon continues to cycle through the operation modes, until you reach three (3) red LED flashes / 3 audible beeps. Release the range buttons to set.

In this mode, **External Antenna** (Ecolab P/N: 92052232) is plugged into either RJ11 connector on the Bed Beacon, and both the Bed Beacon and the External Antenna are attached to the bed. The patient zone is created by alternating between the Internal and the External Antenna, allowing for extension of the patient zone area. The Bed Beacon communicates with badges and the System network hubs. The **AC Line Power Sensor** is connected to the Bed Beacon through the remaining RJ11 connector. The **External Antenna** is connected to the Bed Beacon using the remaining RJ11 connector and is then also attached to the bed, in a central location, as it will generate the patient zone. The **Bed Beacon** is then able to detect when the bed's AC power cord is connected to a wall outlet, causing the **Bed Beacon** to emit a beep and the LED to flash green. If the bed's power cord is disconnected from AC wall power, the **Patient Zone** will be turned off and will remain off while the bed is in transit. When the **AC Line Power (ALP) Sensor** senses that the power is disconnected, the **Bed Beacon** will emit a beep and flash the red LED. When the bed's power cord is plugged in again, the **Bed Beacon** will emit a beep and flash the green LED.

3.4 Stretcher Mode, using Internal Antenna Only

**** 4 red flashes and beeps – operation mode 4

“Stretcher Mode” is used when the **Bed Beacon** is mounted directly under a patient stretcher and there is no **AC Line Power Sensor** connected to it. In stretcher mode, the **Bed Beacon** will always have the **Patient Zone** “on”. To set the **Bed Beacon** mode, press and hold both range buttons down. Immediately, the **Beacon** will flash three (3) times the green LED. With each flash, there will be an audible alert (beep). Continue to hold the range buttons down while the bed beacon continues to cycle through the operation modes, until you reach four (4) red LED flashes / 4 audible beeps. Release the range buttons to set.

3.5 Stretcher Mode, using External Antenna Only

***** 5 red flashes and beeps – operation mode 5

“Stretcher Mode using External Antenna” is used when the **Bed Beacon** and the **External Antenna** (Ecolab P/N: 92052232) are both mounted directly under a patient stretcher and there is no **AC Line Power Sensor** connected to the Bed Beacon. The **External Antenna** is connected to the Bed Beacon using the remaining RJ11 connector and is then also attached to the bed, in a central location, as it will generate the patient zone. In stretcher mode, the **Bed Beacon** will always have the **Patient Zone** “on”, but in this mode the patient zone is created by the External Antenna only. To set the mode, press and hold both range buttons down. Immediately, the **Beacon** will flash three (3) times the green LED. With each flash, there will be an audible alert (beep). Continue to hold the range buttons down while the bed beacon continues to cycle

through the operation modes, until you reach five (5) red LED flashes / 5 audible beeps. Release the range buttons to set.

3.6 Stretcher Mode Using Internal and External Antenna Mode (toggling back and forth between the internal and external antenna)

***** 6 red flashes and beeps – operation mode 6

“Stretcher Mode using Internal and External Antenna” is used when the **Bed Beacon** and the **External Antenna** (Ecolab P/N: 92052232) are mounted directly under a patient stretcher and there is no **AC Line Power Sensor** connected to it. In this stretcher mode, the **Bed Beacon** will always have the **Patient Zone** “on”. In this mode the patient zone is created alternately by the External Antenna and the Bed Beacon’s internal antenna. This is used to extend the effective area of the patient zone. To set the mode, press and hold both range buttons down. Immediately, the **Beacon** will flash three (3) times the green LED. With each flash, there will be an audible alert (beep). Continue to hold the range buttons down while the bed beacon continues to cycle through the operation modes, until you reach six (6) red LED flashes / 6 audible beeps. Release the range buttons to set.

4. Installation

4.1 Mounting Bracket 9205-2337 Installation - Patient Bed

The **Mounting Bracket** can be mounted in the center of the lower bed frame (see Figure 3).

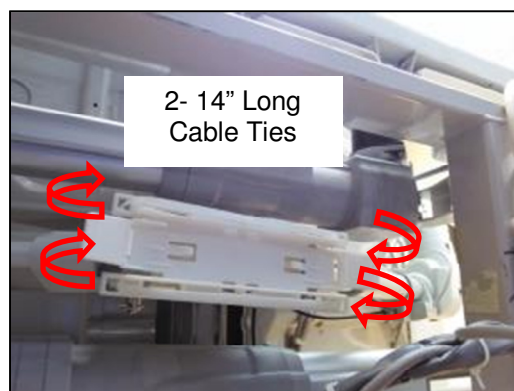


Figure 3. **Mounting Bracket on Patient Bed Model Hill-Rom Versa Care P3200**

Make sure the **Mounting Bracket** location surface is clean and free of dirt and oil. Spray a mixture of liquid dish soap (1% soap per gallon of water) and scrub the mounting surface with a plain Scotch-Brite scouring pad. Wipe the surface with a lint-free cloth to clean off any contaminants. Next, use a soft cloth dampened with isopropyl alcohol to further clean the mounting surface. Repeat the process as necessary. Ensure that the mounting surface is dry before attaching the **Mounting Bracket**. Peel off the double-sided tape liner on the **Mounting Bracket**. Firmly press and hold the **Mounting Bracket** in place for several minutes. Secure the **Mounting Bracket** with the two (2) 14 inch long cable ties, one (1) on each side of the **Bracket**. Finally, flush cut off the excess cable tie ends.

4.2 Bed Beacon Installation

Remove the front battery cover, using a Phillips head screw driver. Correctly position battery polarity and insert two new Duracell® Coppertop alkaline D cell batteries (see Figure 4). After the **Bed Beacon** batteries are installed, the **Bed Beacon** will perform a POST (Power On Self-Test).

(a)*** - Three green LED flashes (Range Switch test). If only one green flash, + range switch is stuck closed (pressed). If only two green flashes, - range switch is stuck closed (pressed).

(b) Pause

(c)*** - Alternate between green and red LED 2 times (LED test)

(d) BEEP - Piezo test

(Note: If the POST LED flashes red three times, then the **Bed Beacon** is defective.)

Replace the **Bed Beacon's** battery door and re-tighten screw. Set the **Bed Beacon** in the desired Bed/Stretcher Mode.



Figure 4. Installing Batteries

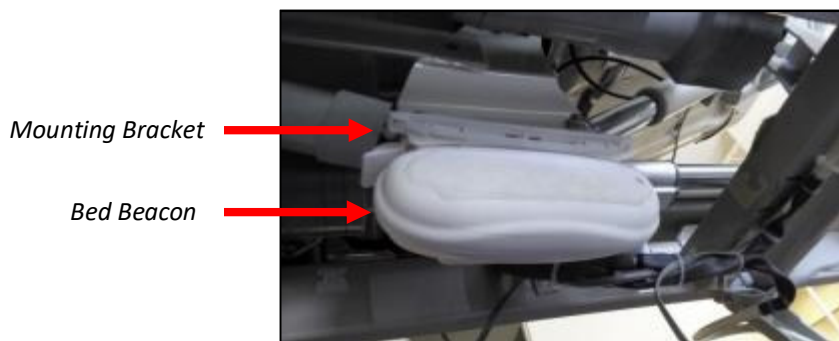


Figure 5. Version 1 Bed Beacon (9205-2185) Mounted on Patient Bed Model Hill-Rom Versa Care P3200

(Note: Photographs showing the mounting of the Bed Beacon on Hospital Beds will depict an earlier version of the Bed Beacon, 9205-2185.)

Before attaching the **Bed Beacon** to the **Mounting Bracket**, position the **Bed Beacon** based on **AC Line Power Sensor** cable management. Press the **Bed Beacon** (flat side) into the **Bracket** until a click is heard. Pull downward on the **Bed Beacon** to verify that the attachment is secure.

To unclip the **Beacon**, push up on either **Mounting Bracket** tab to release (See Figure 6).



Figure 6. Mounting Bracket Release Tabs

4.3 AC Line Power Sensor

The **AC Line Power Sensor** (see Figure 7) is strapped to the bed's AC power cord approximately three (3) inches from the bed cord's power connector. The **AC Line Power Sensor** senses if the bed is connected to AC power, and will activate or deactivate the **Bed Beacon** accordingly. The **AC Line Power Sensor** is connected to the **Bed Beacon**. It mounts to bed's power cord with two (2) cable ties. Make sure a service loop is left in the **AC Line Power Sensor's** cable run, allowing the bed to be raised and lowered. **Note:** Stretchers are installed without the **AC Line Power Sensor** because they don't have an AC power cord. In stretcher mode; the stretcher field is continuously on.



Figure 7. AC Line Power Sensor Attached to Power Cord

Cable Tie mounting

Attach the **AC Line Power Sensor** to the patient bed's AC power cord using 2 cable ties (see Figure 7). Use cable tie holders for routing the wire to the **Bed Beacon**, and make sure the mounting surface is clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean both surfaces.

5. Cable Management – Hill Rom VersaCare P3200

1. Place cable tie holders and cable ties, then route cable to the **Bed Beacon** as shown (see Figure 7, 8 & 9). Make sure the mounting surfaces are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean the mounting surfaces. Avoid routing cables near the bed's patient fall detect sensor IR beam path.

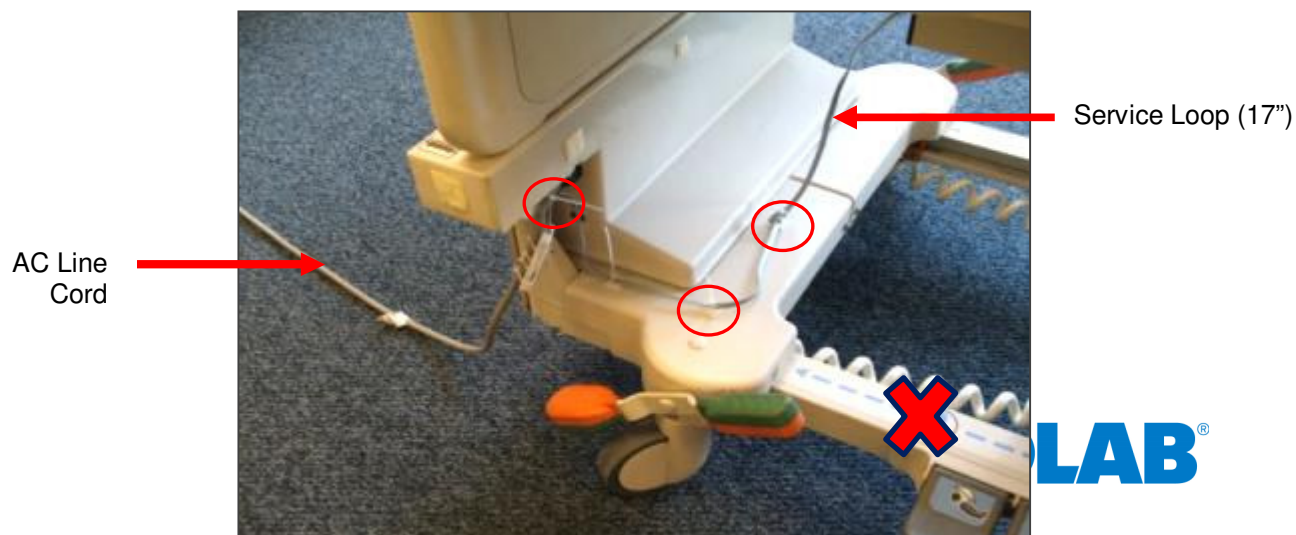


Figure 8. AC Line Power Sensor Cable Routing and Cable Tie Holder Placement

2. Place cable tie holders and route the **AC Line Power Sensor** cable along the inside rail and tie down any excess cable (service loop) as shown (see Figure 8).



Figure 9. AC Line Power Sensor Cable Routing and Cable Tie Holder Placement

3. Place a cable tie holder and route the **AC Line Power Sensor** cable as shown (see Figure 10).



Figure 10. AC Line Power Sensor Cable Routing to Bed Beacon and Cable Tie Holder Placement

6. Adjusting the Patient Zone

6.1 Range Adjustment on Beds

The **Bed Beacon** communication range is adjusted during installation to optimize communication with the healthcare worker **Badge**. To adjust the communication range, remove the rubber range cover (see *Figure 11*) located on the back of the **Bed Beacon**. After the range cover has been removed, the **Bed Beacon** communication range is adjusted by repeatedly pushing either the “+” increment or the “-” decrement range buttons (see *Figure 11*) to achieve the desired range. The total adjustable range of the **Bed Beacon** is from 1 to 10 increments between the lowest and the highest setting. Each push of the “+” button will produce an audible beep and the green LED will flash once indicating one (1) increment up. When the highest setting is reached, two (2) audible beeps will be heard and two (2) green LED flashes will be visible. Each push of the “-” button will produce an audible beep and the red LED will flash once indicating one (1) increment down. When the lowest setting is reached, two (2) audible beeps will be heard and two (2) red LED flashes will be visible. To reset a **Bed Beacon’s** range back to the factory default setting, push and hold both the “+” and “-” buttons at the same time. Three (3) audible beeps will be heard and the LEDs will simultaneously flash green three (3) times. The factory default range setting for the **Bed Beacon** is position 10. Replace the range button cover when the adjustment is complete.

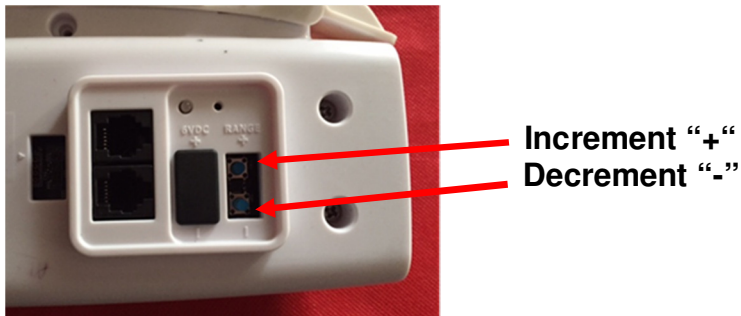


Figure 11. *Bed Beacon Range Adjustments*

6.2 Test Badge Range Adjustment

1. Bed Installs - Adjust the bed and mattress height to the average height for beds in the facility.
Stretcher Installs - Adjust the stretcher and mattress height to the average height for stretchers in the facility.
2. Bed Installs - While holding the Test **Badge** (without **Network**) 50 inches from the floor and 18 inches from the bed's handrail center location (side of bed) (see *Figures 12 & 13*), the Test **Badge's** green and yellow LEDs should flash once a second, indicating that the proper activation range and communication have been successfully established. If range is not 18" adjust the **Bed Beacon's** range buttons until the **Badge** activation range is 18 inches (see *Figure 12*). Once the range has been set, use the Test **Badge** (with **Network**) to verify the **Bed Beacon** communication to the **Proprietary wireless network** by placing the **Badge** in the field.

Stretcher Installs - While holding the Test **Badge** (Network) 50 inches from the floor (*Figure 12*) and at the edge of the stretcher handrail center location (side of stretcher), the Test **Badge's** green and yellow LEDs should flash once a second, indicating that the proper activation range and communication have been successfully established. If range is not at the edge of the handrail adjust the **Bed Beacon's** range buttons until the **Badge** activation range is at the edge. Once the range has been set, use the Test **Badge** (with **Network**) to verify the **Bed Beacon** communication to **Network** by placing the **Badge** in the field.



Figure 12. Test Badge Height Is 50 Inches From the floor



Figure 13. Test Badge Activation Distance from Bed's Handrail Center is 18" (at edge for stretcher installations)

6.3 Ecolab Dashboard Test with the Network Badge and Bed Beacon

If needed, please refer to "Ecolab™ Hand Hygiene Program Compliance Monitoring System Dashboard Directions for Use" for further information on how to use the **Dashboard**.

1. Using a HCW **Badge**, verify that the **Bed Beacon** and **Badge** activate at 18 x 50 inches. Also verify that the **Bed Beacon's** device address appears on the **Dashboard**.
2. Using a web browser, type the facility's **Dashboard** URL and you should see login screen shown in *Figure 14*.
3. Type in the correct email address and password to login to the **Dashboard**.

A screenshot of a web login form. It has two input fields: "Email Address:" with a placeholder "Email address" and "Password:" with a placeholder "Password". Below the fields is a blue button labeled "Sign in".

Figure 14. Dashboard Entry Screen

4. When access is granted you should see the screen shown in *Figure 15*.

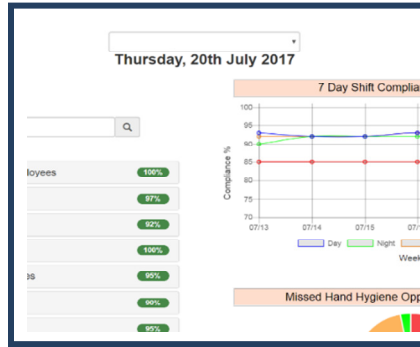


Figure 15. Main Dashboard Page

- Click on “System Menu” under the System Tab as shown in Figure 16.

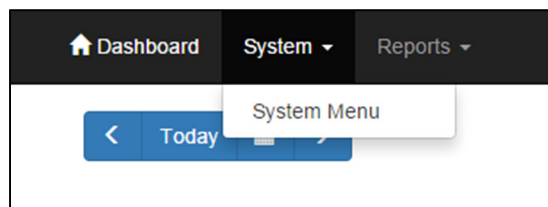


Figure 16. Choose System Menu Screen

- Click on **Bed Beacons** and a list of all current active **Bed Beacons** will be displayed as shown in Figure 17.

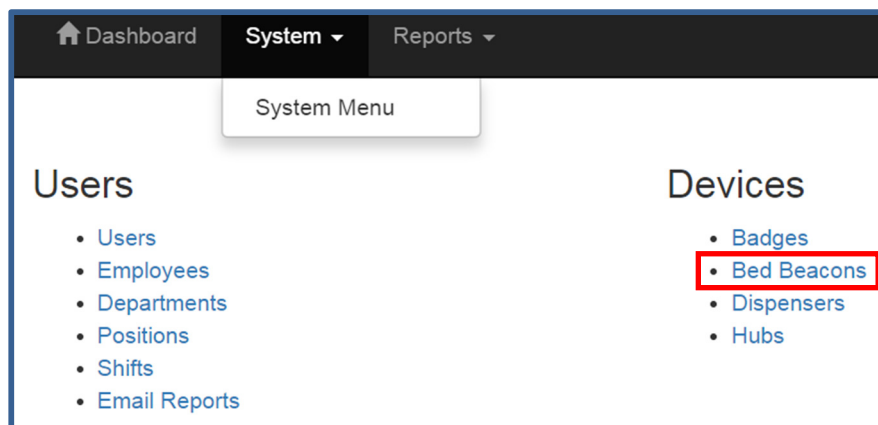


Figure 17. System Menu Screen

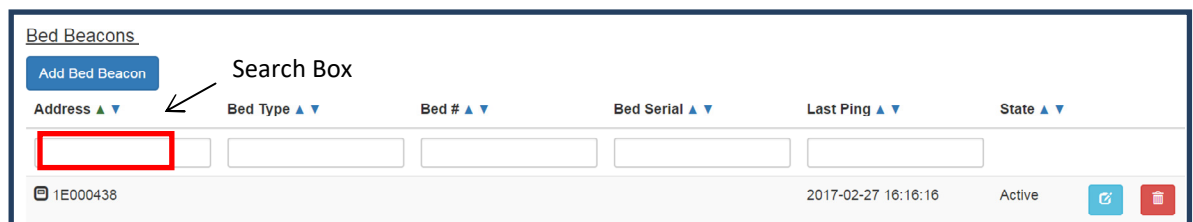


Figure 18. Search for Bed Beacon address

7. Enter the **Bed Beacon's** device address number as shown in *Figure 18* and click the enter key. The list will automatically display the results.

Address ▲ ▼	Bed Type ▲ ▼	Bed # ▲ ▼	Bed Serial ▲ ▼	Last Ping ▲ ▼	State ▲ ▼
1E000438				2017-02-27 16:16:16	Active

Figure 19. Device Screen Results

8. If the **Bed Beacon** address appears in the results and has a recent “Ping Time”, this confirms that the **Bed Beacon** successfully communicated with the **Network** (see *Figure 19*). If the **Bed Beacon** address does not appear in the results, retest with a **Badge**.

There may be circumstances which require a **Bed Beacon's** range to be adjusted, after the initial installation, to properly communicate with a **Badge**. Examples of why a **Bed Beacon's** range may need to be adjusted are explained below.

Cases for **Bed Beacon** Range Adjustment

1. Two beds with **Bed Beacons** in one room.
2. Two beds with **Bed Beacons** in separate but adjacent rooms separated by a wall.

In both cases, the range of communication of a given **Bed Beacon** may overlap or interfere with a nearby **Bed Beacon**. If this occurs, the HCW's **Badge** may communicate with a nearby **Bed Beacon** instead of the intended **Bed Beacon**. This miscommunication can result in the **System** reporting contact with the wrong patient. To eliminate the possibility of miscommunication, the range of **Bed Beacons** in close proximity may be decreased, at any time, to ensure that a HCW's **Badge** will communicate with the intended **Bed Beacon**.

7. Metadata Entry

7.1 New Bed Beacon Metadata Entry

Activate the **Bed Beacon** with a **Badge**, twice, at the hospital. The events will be transmitted to the server to be viewed on the **Dashboard**. To locate the active **Bed Beacon**, login to the **Dashboard** and select the “**Bed Beacon**” option under the **System** menu.

The metadata can only be accessed by restricted logins. Read the **Bed Beacon** address on the product label. The **Bed Beacon's** hexadecimal address will start with a 1E. The 1E at the beginning of the address is the device type and tells the **Dashboard** that this is a **Bed Beacon**. Enter the full address into the “Address” column search window. In the example below, the address is 1E000767 (see *Figure 20*). Click the “Click to Edit Device” icon under the “State” column to open the device's metadata page.

Address ▲ ▼	Bed Type ▲ ▼	Bed # ▲ ▼	Bed Serial ▲ ▼	Last Ping ▲ ▼	State ▲ ▼
1E000438				2017-02-27 16:16:16	Active

Figure 20. Active Bed Beacon List



- Device Information (see Figure 21)
The **Bed Beacon** device information is optional. The device information can be found on the product label. An example from the product label can be located on the same page (see Figure 21)
- Bed Type: example – “Hill Rom/VersaCare” (optional metadata). The bed manufacturer and model information can assist staff in identifying a bed
- Bed Number: example – “001” (required metadata). The installer can choose to place a numbered label at the bed’s footboard frame to easily identify beds since the **Bed Beacon** label may not be in a convenient location to easily view.
- Bed Serial Number: example – “123456” (optional metadata). The bed serial number information can assist staff in identifying a bed.
- State:
 - Active- Device is currently active within the **System**.
 - Inactive- Device will be removed from any reports.
 - RMA- Repair

The figure shows a web form for 'Bed Beacon Information' and a product label for 'UltraClenz'.

Bed Beacon Information Form:

- Address: 1E000767
- Part No: PartNo
- Revision: Revision
- Serial: Serial
- Bed Type: Hill Rom/VersaCare/Bed
- Bed #: BED 199
- Bed Serial: Bed Serial
- State: Active
- Save Beacon button

UltraClenz Product Label:

- QR Code
- UltraClenz
- 0A003CBA
- Model: FA51527-02
- Rev. 2.0
- 1616-001B-00027E
- FCC ID: Z10-FA51527, IC: 10060A-FA51527
- THIS DEVICE COMPLIES WITH PART 15 OF THE FCC RULES. OPERATION IS SUBJECT TO THE FOLLOWING TWO CONDITIONS: (1) THIS DEVICE MAY NOT CAUSE HARMFUL INTERFERENCE, AND (2) THIS DEVICE MUST ACCEPT ANY INTERFERENCE RECEIVED, INCLUDING INTERFERENCE THAT MAY CAUSE UNDESIRABLE OPERATION.

Red lines connect the label information to the form fields: '0A003CBA' to 'Address' and 'Rev. 2.0' to 'Revision'.

Figure 21. Bed Beacon Metadata

7.2 Replacing an Existing Bed Beacon – Swapping Metadata

Remove the **Bed Beacon** and read the address from the product label. Select the “Replace Device” function under system menu (see Figure 22). The replacement tab will open as shown in Figure 23. Enter or scan via bar code reader the old **Bed Beacon’s** address into the first box labeled “old device address”. Next, enter or scan via bar code reader the new **Bed Beacon’s** address into the second box labeled “new device address”. Finally, click the “submit button” to complete replacing the existing **Bed Beacon**.

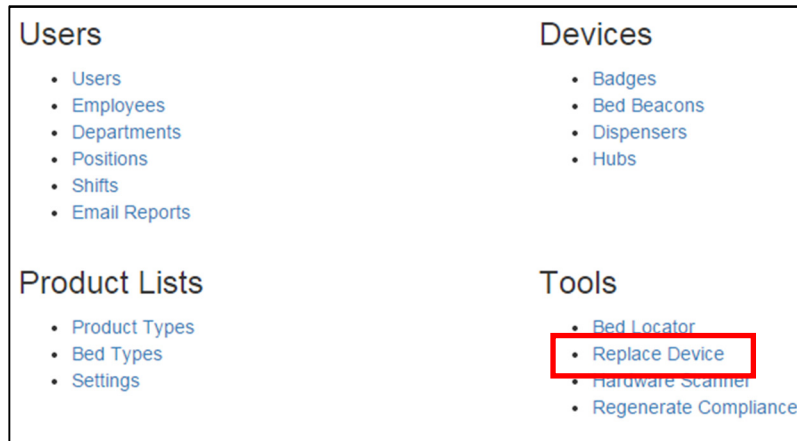


Figure 22. Replace Device

The screenshot shows the 'Replace Device' form. It has two input fields: 'Old Device Address' and 'New Device Address', each with a 'Scan' button next to it. Below these fields is a 'Submit' button. To the right of the form is a QR code for the device. The QR code contains the following information: 'LiltronCensor', '0A003CBA', 'Model: FA1627-02', 'Rev: 2.0', '1616-001B-00027E', and 'PCC ID: 290-FA51027-IC19080A-FA51027'. A red line points from the 'Address' label to the QR code.

Figure 23. Replace Device Tab

8. Troubleshooting

8.1 No Bed Beacon Heartbeats on the Dashboard

If the **Bed Beacon** has stopped transmitting heartbeats to the **Dashboard**, there are three possible causes:

1. The batteries have died. Place a new Duracell Coppertop alkaline D set of batteries in the **Bed Beacon** and test the **Bed Beacon's** other function, communications with a **Badge** as described in section 6.3. If the unit tests well, then give the **Bed Beacon** a couple of hours to transmit a few heartbeats to the **Dashboard**.
2. If new batteries do not fix the problem, there may be a problem with the **Bed Beacon** hardware. Remove and replace the **Bed Beacon**. Refer to section 7.2 – *Replacing an Existing Bed Beacon – Swapping Metadata* -- for instructions on how to replace/update device information on the **Dashboard**.
3. If replacing the **Bed Beacon** does not remedy the problem, there may be a problem with **the Network**. Contact your Ecolab representative for further assistance.

8.2 The Bed Installation Does Not Communicate with a Badge

Verify that the bed's AC power cord is plugged into an active power source and that the **AC Line Power Sensor** is connected to the bed's power cord.

Note: If any of the components of the Bed Installation are replaced, the bed's range will have to be re-checked and readjusted as needed.

9. System Component Care and Maintenance

9.1 Cleaning the Components

The components may be cleaned by wiping with a soft cloth. The cloth should be damp but not wet. A soft cloth dampened with isopropyl alcohol works well to clean the external surfaces. Only the exterior of the components may be cleaned. Do not attempt to clean any interior surface of these components as this will damage the component's circuitry. Do not use abrasive cleaners or cleaning products in aerosol cans as they will damage the component's finish.

9.2 Handling the Bed Beacon

The **Bed Beacon** is an electronic device and should be handled with care. Like other electronic devices such as cell phones, the **Bed Beacon** must be protected from extreme heat, cold and moisture. Avoid handling the **Bed Beacon** with wet hands or exposing it to rain. Avoid dropping or tossing the **Bed Beacon**. The shock can damage the **Beacon's** sensitive internal electronics.

9.3 Replacement Batteries

Replacement batteries for **Bed Beacons** must be Duracell® Coppertop™ alkaline D batteries. Use of any other brand batteries will result in suboptimal performance.

Appendix A - Certification and Safety Approvals

FCC Statement

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment causes harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try and correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

WARNING: Changes or modifications not expressly approved by Ecolab could void the user's authority to operate the equipment.

The antenna provided by the manufacturer is specific to this equipment and may not be replaced by any other antenna models.

RF EXPOSURE:

"This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator and your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter."

Industry Canada

Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that necessary for successful communication.

Conformément à la réglementation d'Industrie Canada, le présent émetteur radio peut fonctionner avec une antenne d'un type et d'un gain maximal (ou inférieur) approuvé pour l'émetteur par Industrie Canada. Dans le but de réduire les risques de brouillage radioélectrique à l'intention des autres utilisateurs, il faut choisir le type d'antenne et son gain de sorte que la puissance isotrope rayonnée équivalente (p.i.r.e.) ne dépasse pas l'intensité nécessaire à l'établissement d'une communication satisfaisante.

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

This radio transmitter (IC:10060A-92053072) has been approved by Industry Canada to operate with the antenna types listed below with the maximum permissible gain and required antenna impedance for each antenna type indicated. Antenna types not included in this list, having a gain greater than the maximum gain indicated for that type, are strictly prohibited for use with this device.

Le présent émetteur radio (IC:10060A-92053072) a été approuvé par Industrie Canada pour fonctionner avec les types d'antenne énumérés ci-dessous et ayant un gain admissible maximal et l'impédance requise pour chaque type d'antenne. Les types d'antenne non inclus dans cette liste, ou dont le gain est supérieur au gain maximal indiqué, sont strictement interdits pour l'exploitation de l'émetteur.

Magnetic Loop Antenna

Ceramic Chip Antenna: -2.5 dBi

RF Exposure: This equipment complies with Industry Canada radiation exposure limits set forth for an Uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator and your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Exposition aux radiofréquences : Cet équipement est conforme aux limites d'exposition par rayonnements définies par l'industrie du Canada pour une utilisation dans un environnement non clos. Cet équipement doit être utilisé à une distance minimale de 20 cm entre l'émetteur de radiation et votre corps. Cet émetteur ne doit pas être situé au même endroit qu'un autre émetteur et ne doit pas être connecté à une antenne différente.

The antenna provided by the manufacturer is specific to this equipment and may not be replaced by any other antenna models.

L'antenne fournie par le fabricant est spécifique à cet équipement et ne peut pas être remplacée par aucun autre modèle d'antenne..

Appendix C – Bed Label

The bed label (see *Figure 1*) is a unique bed number that will be linked to the **Ecolab Dashboard Bed Beacon** metadata. During installation, the bed labels are located on or near the foot-board frame of the patient bed (see *Figure 2*). The bed number helps to quickly identify the **Bed Beacon** address that is on the bed.



Figure 1. Bed Installation Label



Figure 2. Bed Label Location

Appendix D- Z Bracket (9205-2236)

The **Z Bracket** (see *Figure 1*) is an optional mounting brace for the **Bed Beacon** to be mounted onto the VersaCare model bed. The **Z Bracket** helps to elevate the **Bed Beacon** from the ground when the bed is at its lowest point.

Mounting

Mount the **Z Bracket** under the bed on the metal surface shown in *Figure 2* with double-sided tape and cable tie. Make sure the mounting surface is clean and free of dirt and oil. First, spray a mixture of liquid dish soap (1% of soap per gallon of water) and scrub the mounting surface with a plain Scotch-Brite scouring pad. Wipe the surface with a lint free cloth to clean off any contaminants. Second, use a soft cloth dampened with isopropyl alcohol to further clean the mounting surface. Repeat the process as necessary. Ensure the mounting surface is dry before attaching the **Z Bracket** and firmly push and hold in place for several minutes (see *Figure 2*). Once the **Z Bracket** is firmly mounted, secure the **Bracket** by cable tying it to the bed frame (see *Figure 2*), use small wire cutters to carefully flush cut off and dispose of all cable tie ends. Finally, place the **Bed Beacon** onto the **Z Bracket** as shown in *Figure 2* with double-sided tape.



Figure 1. Z Bracket

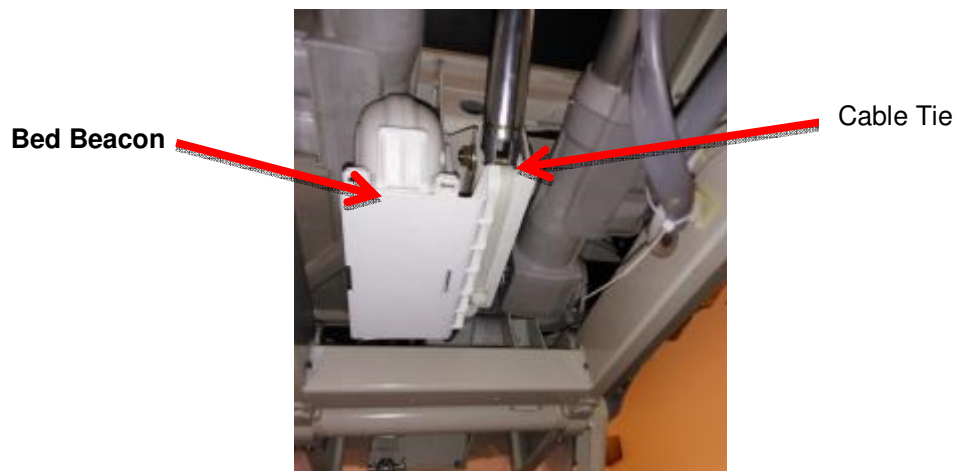


Figure 2. Mounting the Z

Appendix E - Patient Bed Hill-Rom Advanta Installation



1. Installation

1.1 Step 1: Bed Mode Configuration

Install 2 Duracell Coppertop alkaline D cell batteries into the **Bed Beacon**. The **Bed Beacon** needs to be configured to bed mode before installing the **Bed Beacon** onto the bed. To verify the **Beacon's** mode, hold both range buttons down. Immediately, the **Beacon** will alternate four (4) times between green and red flashes. With each flash, there will be an audio alert. Two (2) seconds later, the **Bed Beacon** should flash green and beep once. The green flash and audio alert confirms that the **Beacon** is set for bed mode. If there is one (1) red LED flash and one (1) beep, the **Bed Beacon** was already in bed mode and is now in stretcher mode. Repeat the process to get the **Bed Beacon** back into bed mode.

1.2 Step 2: Installing the Bed Beacon (92053072)

The patient bed must be raised to a height of 32 inches from the floor to gain access to all the mounting locations.

Bed Beacon mounting

Mount the **Bed Beacon** underneath the bed frame with double-sided tape as shown in *Figure 1*. Make sure the **Bed Beacon** and the mounting surface are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean both surfaces. Firmly push and hold in place for several minutes (see *Figure 1*).

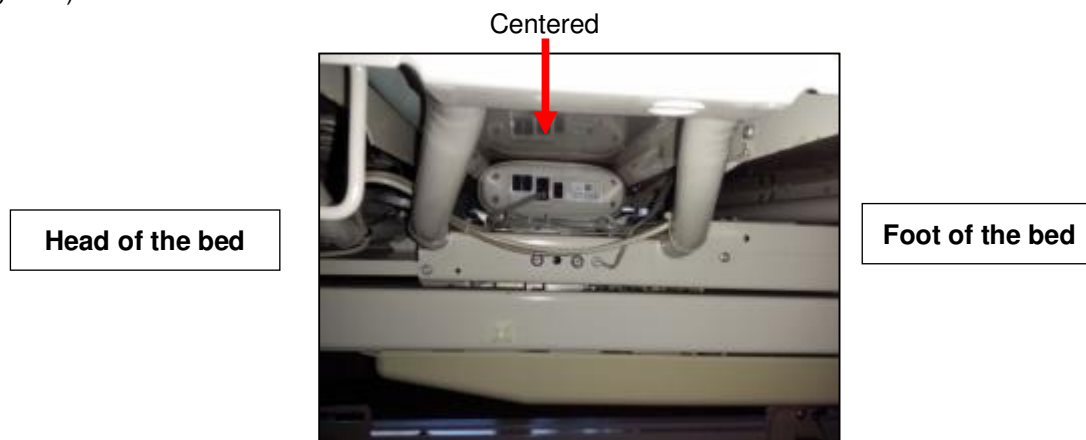


Figure 1. **Bed Beacon** Installation Location on the Bed

1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord

The **AC Line Power Sensor** (see *Figure 2*) is mounted onto the bed's gray AC power cord located underneath the bed. Use the supplied cable ties to secure the **AC Line Power Sensor** to the bed's power cord as shown. The **Bed Beacon** Patient Zone is only enabled when the **Sensor** detects AC power on the bed's line cord after it is plugged into an AC outlet. The **Bed Beacon** Patient Zone will be disabled while the bed is unplugged. The **Bed Beacon** will send heartbeats in both the enabled and disabled modes.

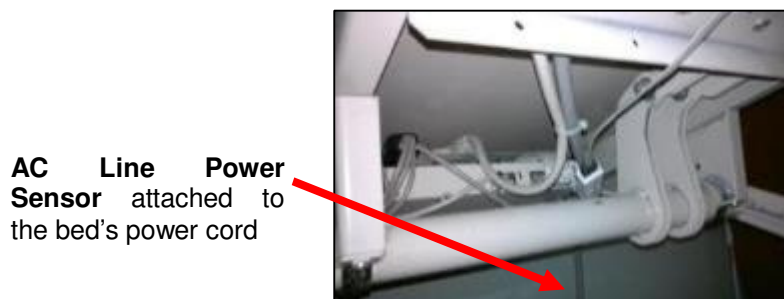


Figure 2. The **AC Line Power Sensor** Installation Location on the Bed

1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range

1. Connect the **AC Line Power Sensor** to the **Bed Beacon** via any RJ11 jacks as shown in *Figure 3*.

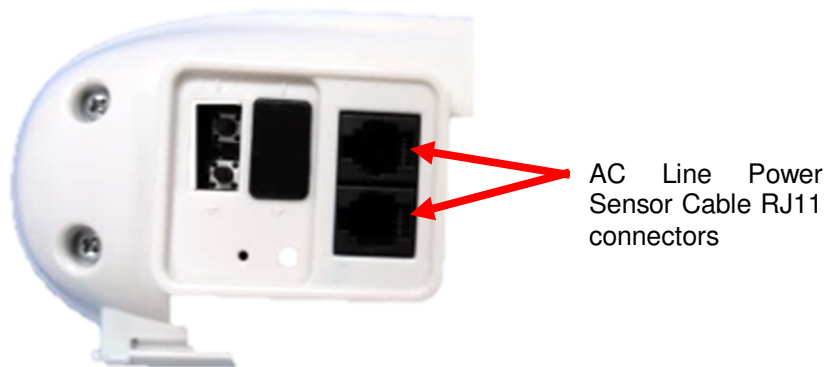


Figure 3. Connecting the **AC Line Power Sensor** to the Bed Beacon via the RJ11 connector

2. Plug the patient bed's AC power cord into a 120vAC power source. This will activate the **Bed Beacon** so that the range setting and all testing can be performed. Refer to section 6 - *Adjusting the Patient Zone* -- in the main section of this document for instructions on adjusting the communication range.

2. Cable Management

2.1 Step 5: AC Line Power Sensor Cable Routing

1. Place cable tie holders, cable ties and route **AC Line Power Sensor** cable to the **Bed Beacon** as shown (see *Figure 4 & 5*). Make sure the mounting surfaces are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean mounting surfaces.

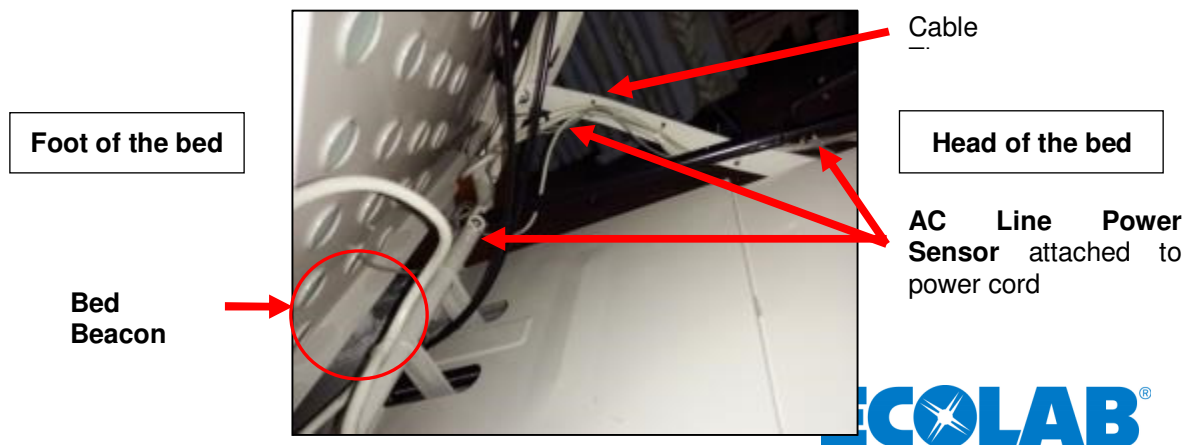


Figure 4. **AC Line Power Sensor** Cable Routing to the **Bed Beacon**

2. Route the **AC Line Power Sensor** cable underneath the headboard along the metal rail frame. Secure the cable with cable ties shown in *Figure 4 & 5*.



Figure 5. **AC Line Power Sensor** Cable Routing

3. Place cable tie holders, route and tie down any excess cable from the **AC Line Power Sensor** as shown (see *Figure 5*). After all cable management is complete, use small wire cutters to carefully flush cut off and dispose of all cable tie ends.

3. Step 6: Final Test

3.1 Disconnect and Reconnect the Bed's AC Power Cord

Confirm that the bed's AC power cord is connected into a 120vAC power source. Disconnect the bed's AC power cord and within 5 seconds you will hear one beep and the **Bed Beacon** LED will flash red once. Reconnect the bed's AC power cord and within 5 seconds you will hear one beep and the **Bed Beacon** LED will flash green once. This will confirm that everything is working properly. Complete the installation by following the instructions on metadata entry in section 7 – *Metadata Entry* -- in the main section of this document.

4. Bed Installation Checklist

Organize all required equipment & tools

1. For installation adjust the bed mattress height to 32 inches
2. For safety reasons disconnect the bed AC power cord from 120vAC source
3. Step 1 - Bed Beacon Mode Configuration
4. Step 2 - Installing the Bed Beacon (92053072)
5. Step 3 - Install the AC Line Power Sensor (9205-2186)

- | | |
|----|--|
| 6. | Step 4 - Adjusting the Bed Beacon to Proper Field Range Refer to section 6 |
| 7. | Step 5 - AC Line Power Sensor Cable Management |
| 8. | Step 6 - Final Test: Disconnect and Reconnect the Bed's AC Power Cord |
| 9. | Step 7 - Metadata Entry |

Appendix F - Patient Bed Hill-Rom 1000 & Advanta 2 Installation



1. Installation

1.1 Step 1: Bed Mode Configuration

Install 2 Duracell Coppertop alkaline D cell batteries into the **Bed Beacon**. The **Bed Beacon** needs to be configured to bed mode before installing the **Bed Beacon** onto the bed. To verify the **Beacon's** mode, hold both range buttons down. Immediately, the **Beacon's** LED will alternate four (4) times between green and red flashes. With each flash, there will be an audio alert. Two (2) seconds later, the **Bed Beacon** should flash green and beep once. The green flash and audio alert confirms that the **Beacon** is set for bed mode. If

there is one (1) red LED flash and one (1) beep, the **Bed Beacon** was already in bed mode and is now in stretcher mode. Repeat the process to get the **Bed Beacon** back into bed mode.

1.2 Step 2: Installing the Bed Beacon (92053072)

The patient bed must be raised to a height of 32 inches from the floor to gain access to all of the mounting locations.

Bed Beacon mounting

Mount the **Bed Beacon** on the center rail and secure the **Bed Beacon** with two 18" natural cable ties as shown in *Figure 1*. Make sure the **Bed Beacon** and the mounting surface are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean both surfaces.

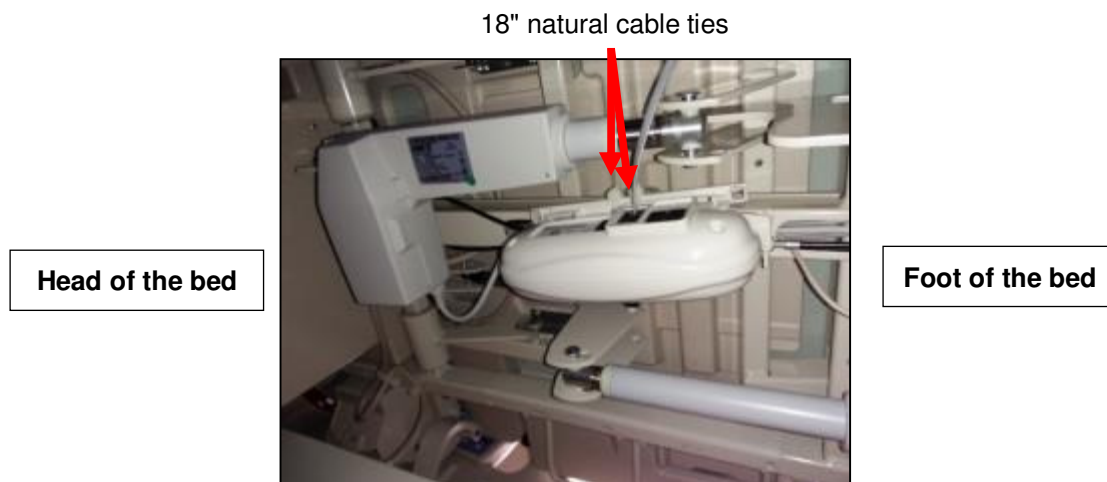


Figure 1. **Bed Beacon** Installation Location on the Bed

1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord

The **AC Line Power Sensor** (see *Figure 2*) is mounted onto the bed's black AC power cord located underneath the bed. Use the supplied cable ties to secure the **AC Line Power Sensor** to the bed's power cord as shown. The **Bed Beacon** Patient Zone is only enabled when the **Sensor** detects AC power on the bed's line cord after it is plugged into an AC outlet. The **Bed Beacon** Patient Zone will be disabled while the bed is unplugged. The **Bed Beacon** will send heartbeats in both the enabled and disabled modes.

AC Line Power Sensor attached to the bed's power cord



Figure 2. The **AC Line Power Sensor** Installation Location on the Bed

1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range

1. Connect the **AC Line Power Sensor** to the **Bed Beacon** via any RJ11 jack as shown in Figure 3.

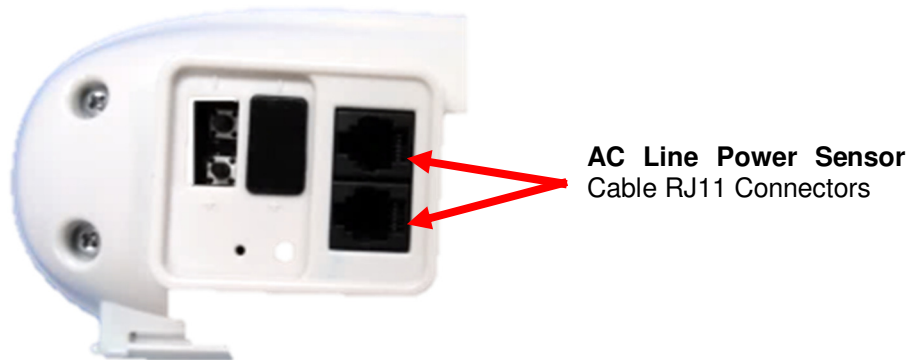


Figure 3. Connecting the **AC Line Power Sensor** to the **Bed Beacon** via the RJ11 connector

2. Plug the patient bed's AC power cord into a 120vAC power source. This will activate the **Bed Beacon** so that the range setting and all testing can be performed. Refer to section 6 - Adjusting the **Patient Zone** in the main section of this document for instructions on adjusting the communication range.

2. Cable Management

2.1 Step 5: AC Line Power Sensor Cable Routing

1. Place cable tie holders, cable ties and route **AC Line Power Sensor** cable to the **Bed Beacon** as shown (see Figure 4 & 5). Make sure the mounting surfaces are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean mounting surfaces.



Figure 4. **AC Line Power Sensor** Cable Routing to the **Bed Beacon**

2. Route the **AC Line Power Sensor** cable underneath the bed along the metal rail frame. Secure the cable with cable ties shown in *Figure 4 & 5*.

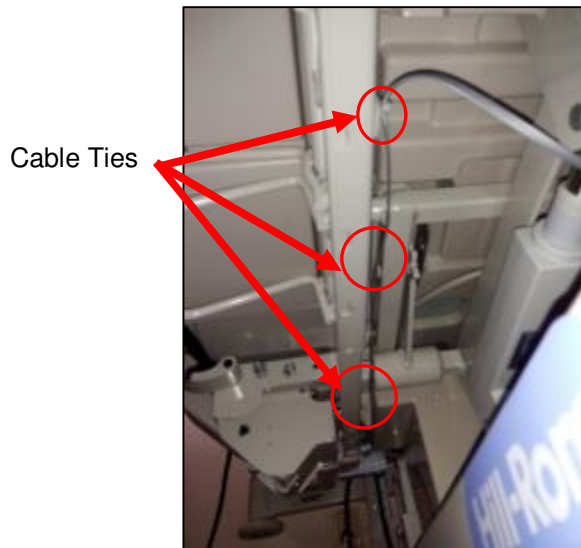


Figure 5. AC Line Power Sensor Cable Routing

3. Place cable tie holders, route and tie down any excess cable from the **AC Line Power Sensor** as shown (see *Figure 5*). After all cable management is complete, use small wire cutters to carefully flush cut off and dispose of all cable tie ends.

3. Step 6: Final Test

3.1 Disconnect and Reconnect the Bed's AC Power Cord

Confirm that the bed's AC power cord is connected into a 120vAC power source. Disconnect the bed's AC power cord and within 5 seconds you will hear one beep and the **Bed Beacon** LED will flash red once. Reconnect the bed's AC power cord and within 5 seconds you will hear one beep and the **Bed Beacon** LED will flash green once. This will confirm that everything is working properly. Complete the installation by following the instructions on metadata entry in section 7 – Metadata Entry in the main section of this document.

4. Bed Installation Checklist

Organize all required equipment & tools

1. For installation adjust the bed mattress height to 32 inches
2. For safety reasons disconnect the bed AC power cord from 120vAC source
3. Step 1 - Bed Beacon Mode Configuration
4. Step 2 - Installing the Bed Beacon (92053072)
5. Step 3 - Install the AC Line Power Sensor (9205-2186)
6. Step 4 - Adjusting the Bed Beacon to Proper Field Range Refer to section 6
7. Step 5 - AC Line Power Sensor Cable Management
8. Step 6 - Final Test: Disconnect and Reconnect the Bed's AC Power Cord
9. Step 7 - Metadata Entry



Appendix G - Patient Bed Stryker Installation

1. Installation

1.1 Step 1: Bed Mode Configuration

Install 2 Duracell Coppertop alkaline D cell batteries into the **Bed Beacon**. The **Bed Beacon** needs to be configured to bed mode before installing the **Bed Beacon** onto the bed. To verify the **Beacon's** mode, hold both range buttons down. Immediately, the **Beacon** will alternate four (4) times between green and red flashes. With each flash there will be an audio alert. Two (2) seconds later, the **Bed Beacon** LED should flash green and beep once. The green flash and beep confirms that the **Beacon** is set for bed mode. If there is one (1) red LED flash and one (1) beep, the **Bed Beacon** was already in bed mode and is now in stretcher mode. Repeat the process to get the **Bed Beacon** back into bed mode.

1.2 Step 2: Installing the Bed Beacon (92053072)

The patient bed must be raised to a height of 32 inches from the floor to gain access to all the mounting locations.

Bed Beacon mounting

Mount the **Bed Beacon** directly double-sided tape as shown in and the mounting surface are cloth dampened with isopropyl surfaces.



underneath the bed's headboard with *Figure 1*. Make sure the **Bed Beacon** clean and free of dirt and oil. A soft alcohol works well to clean both



Figure 1. **Bed Beacon** Installation Location on the Bed

1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord

The **AC Line Power Sensor** (see *Figure 2*) is mounted onto the bed's black AC power cord located underneath the bed. Use the supplied cable ties to secure the **AC Line Power Sensor** to the bed's power cord as shown. The **Bed Beacon** field is only enabled when the **Sensor** detects AC power on the bed's line cord after it is plugged into an AC outlet. The **Bed Beacon** field will be disabled while the bed is unplugged. The **Bed Beacon** will send heartbeats in both the enabled and disabled modes.

AC Line Power Sensor attached to the bed's power cord



Figure 2. The **AC Line Power Sensor** Installation Location on the Bed

1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range

1. Connect the **AC Line Power Sensor** to the **Bed Beacon** via any RJ11 jacks as shown in *Figure 3*.

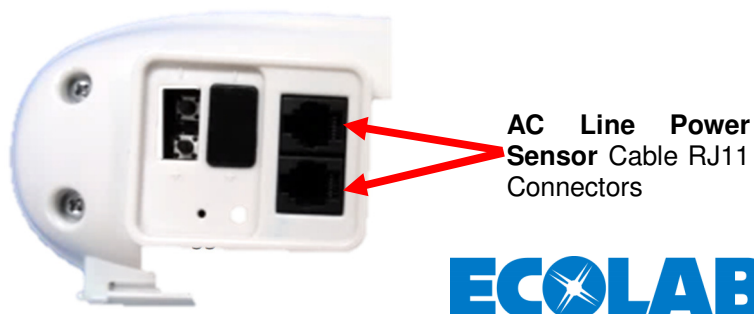


Figure 3. Connecting the **AC Line Power Sensor** to the **Bed Beacon** via the RJ11 connector

2. Plug the patient bed's AC power cord into a 120vAC power source. This will activate the **Bed Beacon** so that the range setting and all testing can be performed. Refer to section 6 - Adjusting the **Patient Zone** in the main section of this document for instructions on adjusting the communication range.

2. Cable Management

2.1 Step 5: AC Line Power Sensor Cable Routing

1. Place cable tie holders, cable ties and route **AC Line Power Sensor** cable to the **Bed Beacon** as shown (see *Figure 4 & 5*). Make sure the mounting surfaces are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean mounting surfaces.

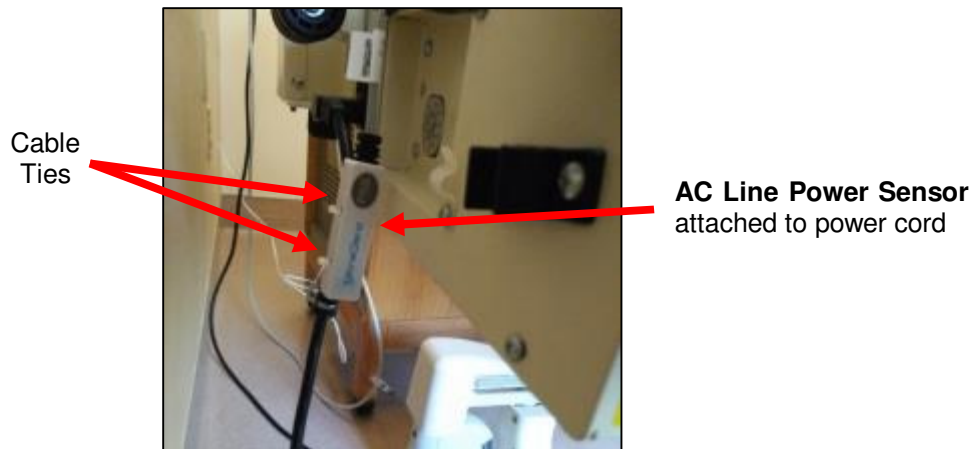


Figure 4. **AC Line Power Sensor** Cable Routing to the **Bed Beacon**

2. Route the **AC Line Power Sensor** cable underneath the headboard along the metal rail frame. Secure the cable with cable ties shown in *Figure 4 & 5*.



Figure 5. **AC Line Power Sensor** Cable Routing

3. Place cable tie holders, route and tie down any excess cable from the **AC Line Power Sensor** as shown (see Figure 5). After all cable management is complete, use small wire cutters to carefully flush cut off and dispose of all cable tie ends.

3. Step 6: Final Test

3.1 Disconnect and Reconnect the Bed's AC Power Cord

Confirm that the bed's AC power cord is connected into a 120vAC power source. Disconnect the bed's AC power cord and within 5 seconds you will hear one audio alert and the **Bed Beacon** LED will flash red once. Reconnect the bed's AC power cord and within 5 seconds you will hear one audio alert and the **Bed Beacon** LED will flash green once. This will confirm that everything is working properly. Complete the installation by following the instructions on metadata entry in section 7 – Metadata Entry in the main section of this document.

4. Bed Installation Checklist

Organize all required equipment & tools

- | | |
|----|--|
| 1. | For installation adjust the bed mattress height to 32 inches |
| 2. | For safety reasons disconnect the bed AC power cord from 120vAC source |
| 3. | Step 1 - Bed Beacon Mode Configuration |
| 4. | Step 2 - Installing the Bed Beacon (92053072) |
| 5. | Step 3 - Install the AC Line Power Sensor (9205-2186) |
| 6. | Step 4 - Adjusting the Bed Beacon to Proper Field Range Refer to section 6 |
| 7. | Step 5 - AC Line Power Sensor Cable Management |
| 8. | Step 6 - Final Test: Disconnect and Reconnect the Bed's AC Power Cord |
| 9. | Step 7 - Metadata Entry |

Appendix H - Patient Bed Hill-Rom TotalCare Installation



1. Installation

1.1 Step 1: Bed Mode Configuration

Install 2 Duracell Coppertop alkaline D cell batteries into the **Bed Beacon**. The **Bed Beacon** needs to be configured to bed mode before installing the **Bed Beacon** onto the bed. To verify the **Beacon's** mode, hold both range buttons down. Immediately, the **Beacon** will alternate four (4) times between green and red flashes. With each flash there will be an audio alert. Two (2) seconds later, the **Bed Beacon** should flash green and beep once. The green flash and audio alert confirms that the **Beacon** is set for bed mode. If there is one (1) red LED flash and one (1) beep, the **Bed Beacon** was already in bed mode and is now in stretcher mode. Repeat the process to get the **Bed Beacon** back into bed mode.

1.2 Step 2: Installing the Bed Beacon (92053072)

The patient bed must be raised to a height of 32 inches from the floor to gain access to all of the mounting locations.

Bed Beacon mounting

Mount the **Bed Beacon**, with double-sided tape, centered (10in x 9in) under the bed on the plastic surface as shown in *Figure 1*. Make sure the **Bed Beacon** and the mounting surface are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean both surfaces. Firmly push and hold in place for several minutes (see *Figure 2*).

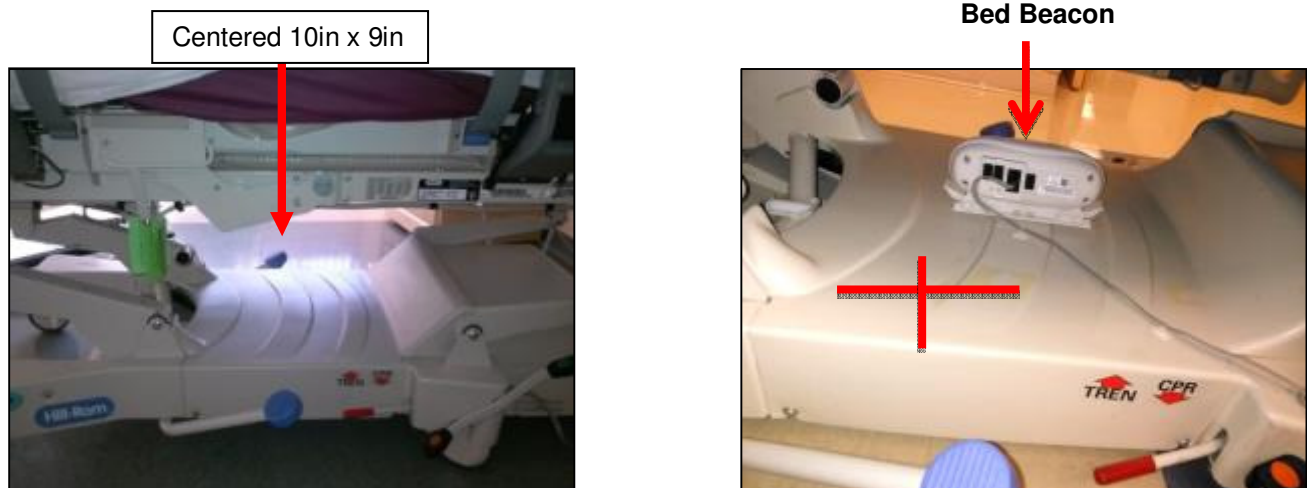


Figure 1 & 2. **Bed Beacon** Installation Location on the Bed

1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord

The AC Line Power Sensor (see *Figure 3*) is mounted onto the bed's gray AC power cord located underneath the bed. Use the supplied cable ties to secure the AC Line Power Sensor to the bed's power cord as shown. The Bed Beacon field is only enabled when the Sensor detects AC power on the bed's line cord after it is plugged into an AC outlet. The Bed Beacon field will be disabled while the bed is unplugged. The Bed Beacon will send heartbeats in both the enabled and disabled modes.



AC Line Power Sensor
attached to the bed's
power cord

Figure 3. The **AC Line Power Sensor** Installation Location on the Bed

1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range

1. Connect the **AC Line Power Sensor** to the **Bed Beacon** via any RJ11 jacks as shown in *Figure 4*.



AC Line Power Sensor
Cable RJ11 Connectors

Figure 4. Connecting the **AC Line Power Sensor** to the **Bed Beacon** via the RJ11 connector

2. Plug the patient bed's AC power cord into a 120vAC power source. This will activate the **Bed Beacon** so that the range setting and all testing can be performed. Refer to section 6 - Adjusting the **Patient Zone** in the main section of this document for instructions on adjusting the communication range.

2. Cable Management

2.1 Step 5: AC Line Power Sensor Cable Routing

1. Place cable tie holders, cable ties and route **AC Line Power Sensor** cable to the **Bed Beacon** as shown (see *Figure 5 & 6*). Make sure the mounting surfaces are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean mounting surfaces.



Figure 5. **AC Line Power Sensor** Cable Routing to the **Bed Beacon**

2. Place cable tie holders and route the **AC Line Power Sensor** cable on the bottom frame of the bed shown in *Figure 5 & 6*.



Figure 6. **AC Line Power Sensor** Cable Routing

3. Place cable tie holders, route and tie down any excess cable from the **AC Line Power Sensor** as shown (see *Figure 6*). After all cable management is complete, use small wire cutters to carefully flush cut off and dispose of all cable tie ends.

3. Step 6: Final Test

3.1 Disconnect and Reconnect the Bed's AC Power Cord

Confirm that the bed's AC power cord is connected into a 120vAC power source. Disconnect the bed's AC power cord and within 5 seconds you will hear one audio alert and the **Bed Beacon** LED will flash red once. Reconnect the bed's AC power cord and within 5 seconds you will hear one audio alert and the **Bed Beacon**

LED will flash green once. This will confirm that everything is working properly. Complete the installation by following the instructions on metadata entry in section 7 – Metadata Entry in the main section of this document.

4. Bed Installation Checklist

Organize all required equipment & tools

1. For installation adjust the bed mattress height to 32 inches
2. For safety reasons disconnect the bed AC power cord from 120vAC source
3. Step 1 - Bed Beacon Mode Configuration
4. Step 2 - Installing the Bed Beacon (92053072)
5. Step 3 - Install the AC Line Power Sensor (9205-2186)
6. Step 4 - Adjusting the Bed Beacon to Proper Field Range Refer to section 6
7. Step 5 - AC Line Power Sensor Cable Management
8. Step 6 - Final Test: Disconnect and Reconnect the Bed's AC Power Cord
9. Step 7 - Metadata Entry

Appendix I - Patient Bed Hill-Rom Affinity 2 Installation



1. Installation

1.1 Step 1: Bed Mode Configuration

Install 2 Duracell Coppertop alkaline D cell batteries into the **Bed Beacon**. The **Bed Beacon** needs to be configured to bed mode before installing the **Bed Beacon** onto the bed. To verify the **Beacon's** mode, hold both range buttons down. Immediately, the **Beacon** will alternate four (4) times between green and red flashes. With each flash there will be an audio alert. Two (2) seconds later, the **Bed Beacon** should flash green and beep once. The green flash and audio alert confirms that the **Beacon** is set for bed mode. If there is one (1) red LED flash and one (1) beep, the **Bed Beacon** was already in bed mode and is now in stretcher mode. Repeat the process to get the **Bed Beacon** back into bed mode.

1.2 Step 2: Installing the Bed Beacon (92053072)

The patient bed must be raised to a height of 32 inches from the floor to gain access to all of the mounting locations.

Bed Beacon mounting

Mount the **Bed Beacon** directly centered underneath the bed with double-sided tape as shown in *Figure 1*. Make sure the **Bed Beacon** and the mounting surface are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean both surfaces.



Figure 1. **Bed Beacon** Installation Location on the Bed

1.3 Step 3: Installing the AC Line Power Sensor on the Bed's AC Power Cord

The **AC Line Power Sensor** (see *Figure 2*) is mounted onto the bed's gray AC power cord located underneath the bed. Use the supplied cable ties to secure the **AC Line Power Sensor** to the bed's power cord as shown. The **Bed Beacon** field is only enabled when the **Sensor** detects AC power on the bed's line cord after it is plugged into an AC outlet. The **Bed Beacon** field will be disabled while the bed is unplugged. The **Bed Beacon** will send heartbeats in both the enabled and disabled modes.



AC Line Power Sensor attached to the bed's power cord



Figure 2. The **AC Line Power Sensor** Installation Location on the Bed

1.4 Step 4: Adjusting the Bed Beacon to Proper Field Range

1. Connect the **AC Line Power Sensor** to the **Bed Beacon** via any RJ11 jacks as shown in *Figure 3*.



AC Line Power Sensor
Cable RJ11 Connectors

Figure 3. Connecting the **AC Line Power Sensor** to the **Bed Beacon** via the RJ11 connector

2. Plug the patient bed's AC power cord into a 120vAC power source. This will activate the **Bed Beacon** so that the range setting and all testing can be performed. Refer to section 6 - Adjusting the **Patient Zone** in the main section of this document for instructions on adjusting the communication range.

2. Cable Management

2.1 Step 5: AC Line Power Sensor Cable Routing

1. Place cable tie holders, cable ties and route **AC Line Power Sensor** cable to the **Bed Beacon** as shown (see *Figure 4 & 5*). Make sure the mounting surfaces are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean mounting surfaces.

AC Line Power Sensor
attached to power cord



Cable
Ties

ECOLAB[®]



Figure 4. **AC Line Power Sensor** Cable Routing to the **Bed Beacon**

2. Route the **AC Line Power Sensor** cable underneath the headboard along the metal rail frame. Secure the cable with cable ties shown in *Figure 4 & 5*.



Figure 5. **AC Line Power Sensor** Cable Routing

3. Place cable tie holders, route and tie down any excess cable from the **AC Line Power Sensor** as shown (see *Figure 5*). After all cable management is complete, use small wire cutters to carefully flush cut off and dispose of all cable tie ends.

3. Step 6: Final Test

3.1 Disconnect and Reconnect the Bed's AC Power Cord

Confirm that the bed's AC power cord is connected into a 120vAC power source. Disconnect the bed's AC power cord and within 5 seconds you will hear one audio alert and the **Bed Beacon** LED will flash red once. Reconnect the bed's AC power cord and within 5 seconds you will hear one audio alert and the **Bed Beacon** LED will flash green once. This will confirm that everything is working properly. Complete the installation by following the instructions on metadata entry in section 7 – Metadata Entry in the main section of this document.

4. Bed Installation Checklist

Organize all required equipment & tools

1. For installation adjust the bed mattress height to 32 inches
2. For safety reasons disconnect the bed AC power cord from 120vAC source
3. Step 1 - Bed Beacon Mode Configuration
4. Step 2 - Installing the Bed Beacon (92053072)
5. Step 3 - Install the AC Line Power Sensor (9205-2186)

- | | |
|----|--|
| 6. | Step 4 - Adjusting the Bed Beacon to Proper Field Range Refer to section 6 |
| 7. | Step 5 - AC Line Power Sensor Cable Management |
| 8. | Step 6 - Final Test: Disconnect and Reconnect the Bed's AC Power Cord |
| 9. | Step 7 - Metadata Entry |

Appendix J – Stryker Renaissance Stretcher



1. Installation

1.1 Step 1: Bed Mode Configuration

Install 2 Duracell Coppertop alkaline D cell batteries into the **Bed Beacon**. The **Bed Beacon** needs to be configured to stretcher mode before installing the **Bed Beacon** onto the stretcher. To get the **Bed Beacon** into stretcher mode, hold both range buttons down. Immediately, the **Bed Beacon** will alternate 4 times between green and red flashes with each flash there will be an audio alert. Two (2) seconds later, the **Bed Beacon** should flash red and beep twice. The two red flashes and audio alert confirms that the **Bed Beacon** is set for stretcher mode. If there is one (1) green LED flash and one (1) beep, the **Bed Beacon** was already in stretcher mode and is now in bed mode. Repeat the process to get the **Bed Beacon** back into stretcher mode.

1.2 Step 2: Installing the Bed Beacon (92053072)

The patient stretcher must be raised to a height of 32 inches from the floor to gain access to all of the mounting locations.

Bed Beacon mounting

Mount the **Bed Beacon** directly centered underneath the stretcher with double-sided tape as shown in *Figure 1*. Make sure the **Bed Beacon** and the mounting surface are clean and free of dirt and oil. A soft cloth dampened with isopropyl alcohol works well to clean both surfaces.



Figure 1. Bed Beacon Installation Location on the Bed

1.4 Step 3: Adjusting the Bed Beacon to Proper Field Range

1. Refer to section 6 - Adjusting the **Patient Zone** in the main section of this document for instructions on adjusting the communication range.

2. Stretcher Installation Checklist

Organize all required equipment & tools

1.	For installation adjust the stretcher mattress height to 32 inches
2.	Step 1 - Bed Beacon Mode Configuration
3.	Step 2 - Install the Bed Beacon (9205-3072)
4.	Step 3 - Adjusting the Bed Beacon to Proper Field Range Refer to section 6
5.	Step 4 - Metadata Entry