

RENPHO



Embrace Cup

Electric Cupping Massager

User Manual

Model: R-WMV01

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
2. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
3. Children should be supervised to ensure that they do not play with the device.
4. WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons.
 - a) This device should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - c) Close supervision is necessary when this device is used by, on, or near children, invalids, or disabled persons.
 - d) Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - e) Never operate this device if it is not working properly, if it has been dropped or damaged, or dropped into water.
 - f) Do not carry this device by supply cord or use cord as a handle.
 - g) Keep the cord away from heated surfaces.
 - h) Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - i) Never drop or insert any object into any opening.
 - j) Do not use outdoors.
 - k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
5. DANGER: Do not use while bathing or in a shower.
6. DANGER: Do not place or store device where it can fall or be

pulled into a tub or sink.

7. DANGER: Do not place in or drop into water or other liquid.

8. WARNING : Never operate on a soft surface such as a bed or couch where the air openings may be blocked.

Precautions:

1. DO NOT USE if you have an implantable medical electronic device (such as a cardiac pacemaker), a life-sustaining medical device (such as an artificial heart or lung), or other internal medical implants. The device may interfere with or damage these instruments.
2. DO NOT USE if you are connected to medical electronic monitoring or support equipment (such as an electrocardiograph).
- 3. Please consult your doctor before using if you are one of the following persons:**
 - 1) Do not use on swollen, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombophlebitis, varicose veins).
 - 2) Do not use on areas of the skin that are insensitive to heat or have poor circulation.
 - 3) Do not use on the face, neck, or throat area.
 - 4) Consult your doctor before use if you are pregnant, have a medical condition (especially heart disease, blood clotting disorders, or cancer), or have a pacemaker or other electronic medical device implanted.
4. The device uses suction. Do not use for more than 10 minutes on a single area to prevent tissue damage or severe bruising.
5. Immediately discontinue use and consult your doctor if you experience severe pain, dizziness, or skin irritation.
6. The infrared function produces slight heat. Do not use on sensitive skin areas or on areas where you cannot sense temperature. Do not fall asleep while the device is in use.

About the Battery:

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short- circuited by metal objects.

3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

ABOUT THE DEVICE

1. What's in the Box

- 1 x Embrace Cup
- 3 x Cup (S, L, Curved)
- 1 x Storage Pouch
- 1 x USB-C Charging Cable

2. Technical Specifications

Product Name: Electric Cupping
Massager

Model: R-WMV01

Battery Voltage: DC 3.7V

Battery Capacity: 1200mAh 4.44Wh

Infrared Wavelength: 850nm

Auto-off Timer: 3/5/8 min


Charging Time: ≤2.5hrs

Operating Time: ≥70mins

Dimension: 3.35 x 3.35 x 2.52 in /
85 x 85 x 64 mm

Weight: 0.46 lb / 208.5 g

3. Overview

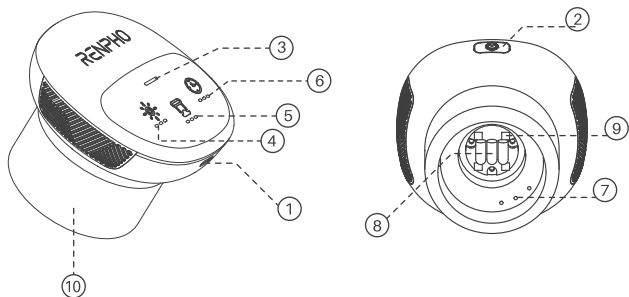

1. USB-C Charging Port

2. Power

3. Power Indicator

4. Light Intensity

5. Suction Intensity

6. Timer

7. Air Inlet

8. LED Light

9. Infrared Heating Light Emitters

10. Interchangeable Cup

USING THE DEVICE

Before Use

1. Charge the Device:

Full charge the device before its first use using the provided cable and a DC 5V charger (not included).

- Charging Indicator: Red flashing light (battery $\leq 99\%$)

Solid blue light (100% charged).

- A full charge (≤ 2.5 hours) provides over 70 minutes of use at maximum suction.

2. Clean & Prep Your Skin:

Wash the skin where the cup will be placed with mild soap and water. Avoid oils, lotions, or creams, as they can break the seal. (Do not use on broken skin, wounds, bruises, varicose veins, moles, or areas of severe inflammation.)

3. Install the Cup

a. Choose a cup size that fits the contour of the body area.



S-Size Cup

For general use on flat surfaces (back, thighs).



L-Size Cup


For larger muscle groups.



Curved Cup

For areas like the spine, neck, or calves.

b. Align the cup with the device's suction port and rotate until it is securely attached.

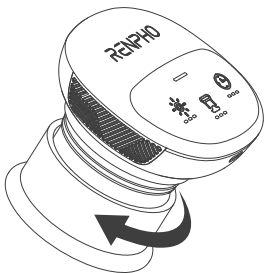

4. Position the Cup

Place the cup firmly against the skin. It must be flat against the skin with no gaps.

Operation

Button	Operation
	• Press and hold to turn On / Off the device, accompanied by a confirmation beep. • When turned on, double press to enter paring mode. • When in paring mode, press to confirm connect.
	• Press to cycle through 3 light intensity levels. • Press and hold to disable light function.
	• Press to cycle through 3 suction levels. Level 1: -30 kPa (1 white LED) Level 2: -45 kPa (2 white LEDs) Level 3: -55 kPa (3 white LEDs) • Press and hold to disable suction function.
	Press to choose a timer: 3 min, 5 min or 8 min when Suction or Infrared function is enabled.

Power Indicator

Status	Indicator
Charging	Flashing Red (Charging) → Solid Blue (Fully charged)
Low Battery	Flashing Red 10 times (Shutdown imminent)
Paring Mode	Flashing Green (Searching) → Solid Green (Connected)
Error/Over-voltage	Flashing Red rapidly

Bruising and Blisters

Bruising after use is normal. This is a common and expected result of cupping. Blisters are a rare but possible side effect. They may occur due to:

Overuse as a beginner (using high intensity or long sessions).

Having skin that is naturally sensitive to suction or heat therapy.

We strongly advise all beginners to start slowly. Use the lowest intensity setting for short sessions (3 minutes per area) to see how your skin reacts. Pay close attention during each session. Stop using the device immediately if you feel strong discomfort or notice any blisters forming.

Acknowledgment of Risk: Cupping therapy can cause temporary skin marks, bruising, or, in rare cases, blisters. The manufacturer is not responsible for skin reactions resulting from misuse or individual skin sensitivity.

Using Multiple Cups (2.4G Wireless Interconnect)

This feature allows you to connect and control multiple Emerge Cups (up to 6) simultaneously as a single network. The actual number depends on how many devices you own.

1. Power On: Turn on all the devices you want to connect.
2. Enter Pairing Mode:
 - a. On EACH individual device, quickly double-press its Power button. Each device's Power indicator will begin flashing green, indicating it is in pairing mode.
 - b. On ONE of the devices (you choose which one), short-press the Power button. Its Power indicator will turn solid green. This device is now broadcasting the network signal.
3. Join Other Devices: For each remaining device (with their green lights still flashing), one by one, short-press their Power button. Once a device successfully joins the network, it will emit a confirmation sound and its Power indicator will turn solid green.

4. After all devices show a solid green light, you can use any device as a controller. Any operation (changing suction intensity, light intensity, or timer) performed on any connected device will be instantly synchronized to all other devices.

To Disconnect (Exit Pairing Mode): To remove a single device or disband the entire network, quickly double-press the Power button on the device(s) you wish to disconnect. The green light will turn off, and revert to blue, indicating it has left the network. Or press and hold the Power button to turn it off.

Note:

- * The device will automatically power off after 2 minutes of inactivity.
- * Do not look directly into the lights.
- * It is not recommended to use the device for more than 10 minutes at a time.

CLEANING & MAINTENANCE

1. After each use, wipe the main unit and cups with a soft, dry cloth. Do not use alcohol, harsh chemicals, abrasive cleaners, or disinfectant wipes.
2. Always turn off the power and disconnect from charging before cleaning.
3. Clean the cups with mild soap and lukewarm water. Rinse and dry thoroughly before storage.
4. Store the device and all components in the included storage bag. Keep in a cool, dry place, away from direct sunlight.
5. Do not submerge the main device unit in water or any other liquid.
6. If the device will be idle for a long period, charge it every 3 months to maintain battery health.

TROUBLESHOOTING GUIDE

Problem	Cause	Solution
The device cannot be turned on.	Battery is weak or depleted.	Charge the device and restart it.
	Power button is not pressed long enough.	Press and hold the Power button for longer duration.
Weak suction.	The air inlet is blocked by debris.	Check the air inlet for possible blockage and clean it.
	The selected suction level is low.	Try a higher suction level (2 or 3).
The cup will not hold suction or loses pressure quickly.	The edge is not in full contact with the skin.	Press the cup firmly and evenly onto the skin, checking for gaps.
The device does not charge.	The charger is faulty.	Try a different 5V charger.
The device stops working suddenly.	The timer is activated.	Restart the device if needed.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not void the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.



TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

+49 40 23969285 (DE ONLY)

Monday-Friday 9:00AM-6:00PM (CET)



Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com (UK&EU)

For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC RF exposure statement:

The device has been evaluated to meet general RF exposure requirement.
The device can be used in portable exposure condition without restriction.

更新日期：	20260112	设计师：	halsey
型 號：	R-WMV01		
尺 寸：	95x135mm		
材 质：	128g铜版纸，双面哑油		
颜 色：	黑白印刷		
工 艺：	/		
装订方式：	折页		