

Mariana ST-A Labels & Markings

Note: All labels provided in this document are SAMPLES only.

1 Mariana ST-A Labels

1.1 Device Labels

1.1.1 Name Plaque

Note: the blue lines will not be printed, those are only identification lines to point out what each item on the label is or for dimensions.

1.1.1.1 AirCurve 11 ST-A (US and Canada)




1.1.1.2 Lumis 11 ST-A (Europe)




Foil product variant badges 'permanently' attached to the side of the product.

1.1.2 Device Serial Number Label – US and Canada

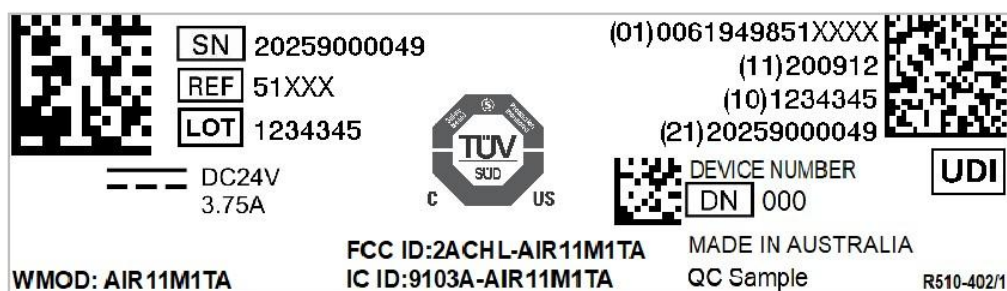

1.1.3 Device Serial Number Label – Europe

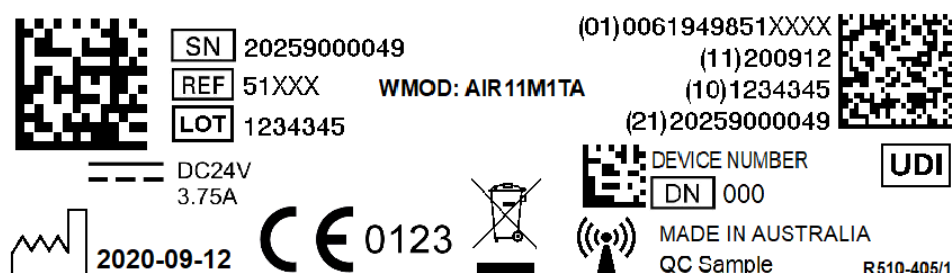

NOTE: The REF will be updated in a future version, it should state "51XXX"

1.1.4 Manufacturer Logo



Branding chrome or colour foiled on front of the product

1.1.5 BF Applied Part Marking



Body Fibrillation Icon moulded into plastic case using an insert on rear of the device, adjacent to the device outlet.

1.1.6 MR Unsafe Marking



White pre-printed label stuck next to the power connector on rear of assembly.


1.1.7 Bottom Case Artwork

 ResMed Pty Ltd 1 Elizabeth Macarthur Drive Bella Vista NSW 2153 Australia See www.resmed.com for other ResMed locations worldwide		EC REP ResMed SAS Parc Technologique de Lyon 292 Allée Jacques Monod 69791 Saint Priest Cedex France	
STORAGE & TRANSPORT RANGES		OPERATING RANGES	
 -25C / 70C -13F / 158F	 95% 5	 2591m 8500ft ALTITUDE LIMITATION	 5C / 35C 41F / 95F
 RTCA/DO-160 Section 21 Category M Compliant & FAA Compliant	MD		IP22
For patent information, see www.resmed.com/ip See symbols glossary at www.resmed.com/symbols		A510-7204/1	

Laser etched on bottom of assembly

ResMed Confidential Proprietary Information, not to be reproduced or made available to third parties without prior consent from ResMed and not to be used in any unauthorised way. © ResMed Ltd 2025

1.1.8 R390-429 Label – FCC & IC NIR M1

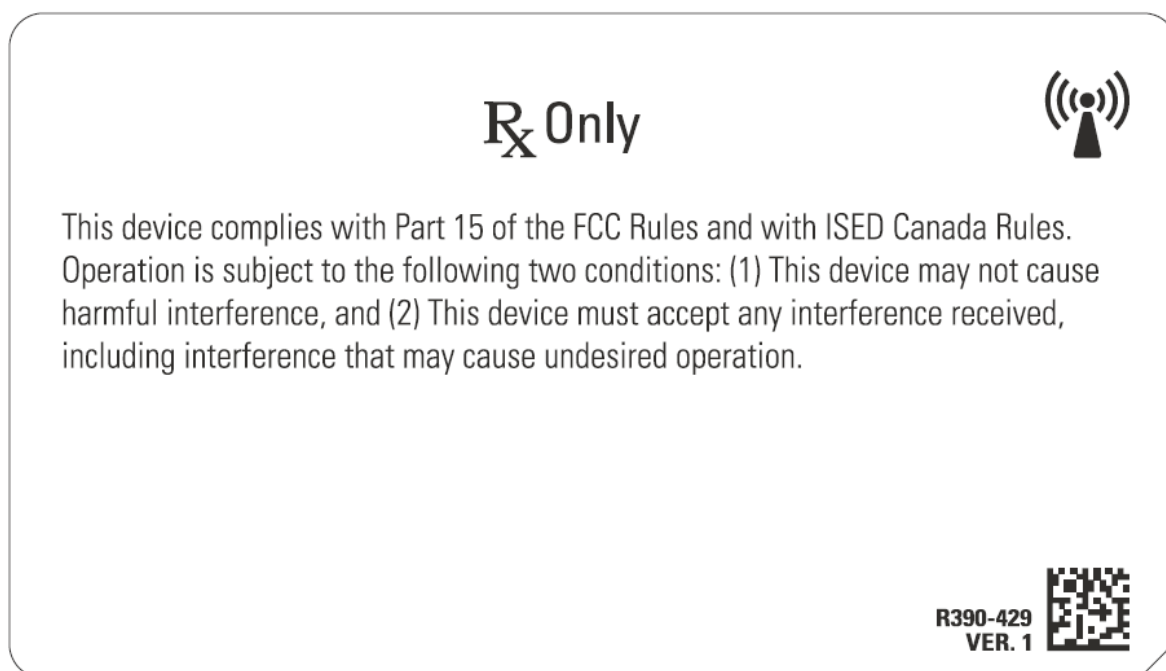


White pre-printed label stuck on bottom of assembly.

1.1.9 Alarm Mute Button Marking



Laser etched on alarm mute button on side of assembly.

1.1.10 Supplemental Oxygen Marking

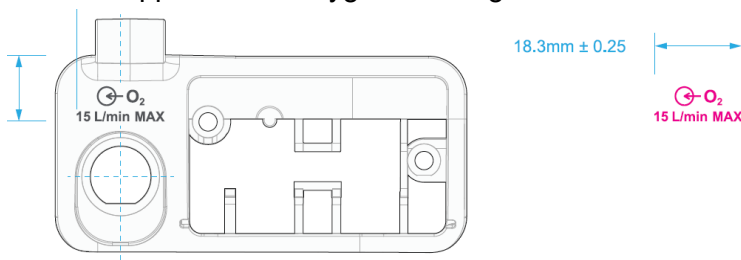

Laser etched above supplemental oxygen port on rear of assembly.

1.2 Device Images

1.2.1 Rear View



1.2.2 Bottom View



1.2.3 Side View – US and Canada Example

