


Set Up

Packaging Station:

- 1 Start the scripClip app.
- 2 Select *Packaging* from the pull-down menu
- 3 Select the color for the Packaging station by clicking on the color square of the scripClip application.
- 4 Position the Barcode scanner and the scripClip Reader as needed.

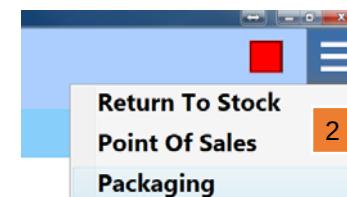

Point of Sale Station:

- 1 Start the scripClip app.
- 2 Select *Point of Sale* from the pull-down menu
- 3 Select the color for the station by clicking on the color square of the scripClip application.
- 4 Position the scripClip Reader.

Return Prescription to Stock:


- 1 Start the scripClip app.
- 2 Select *Return to Stock* from the pull-down menu
- 3 Select the color for the station by clicking on the color square of the scripClip application.
- 4 Position the scripClip Reader.

Packaging the Prescription

Hanging Bag scripClips:

- 1 Using the barcode scanner, scan one by one, up to 10 Rx numbers either from the vial or from the prescription paperwork for a customer.
→ *Scanned Rx numbers appear on the screen.*
- 2 Place the vials and prescription paperwork in an unused Hanging Bag scripClip. Snap the hanging bag closed.
- 3 Hold the Hanging Bag scripClip near the scripClip Reader and press the orange button on scripClip.
→ *The scripClip will flash the station color once.*
- 4 → *scripClip application will display, 'Prescriptions clipped successfully.'*

Place the filled hanging bag scripClip on the will-call rack.


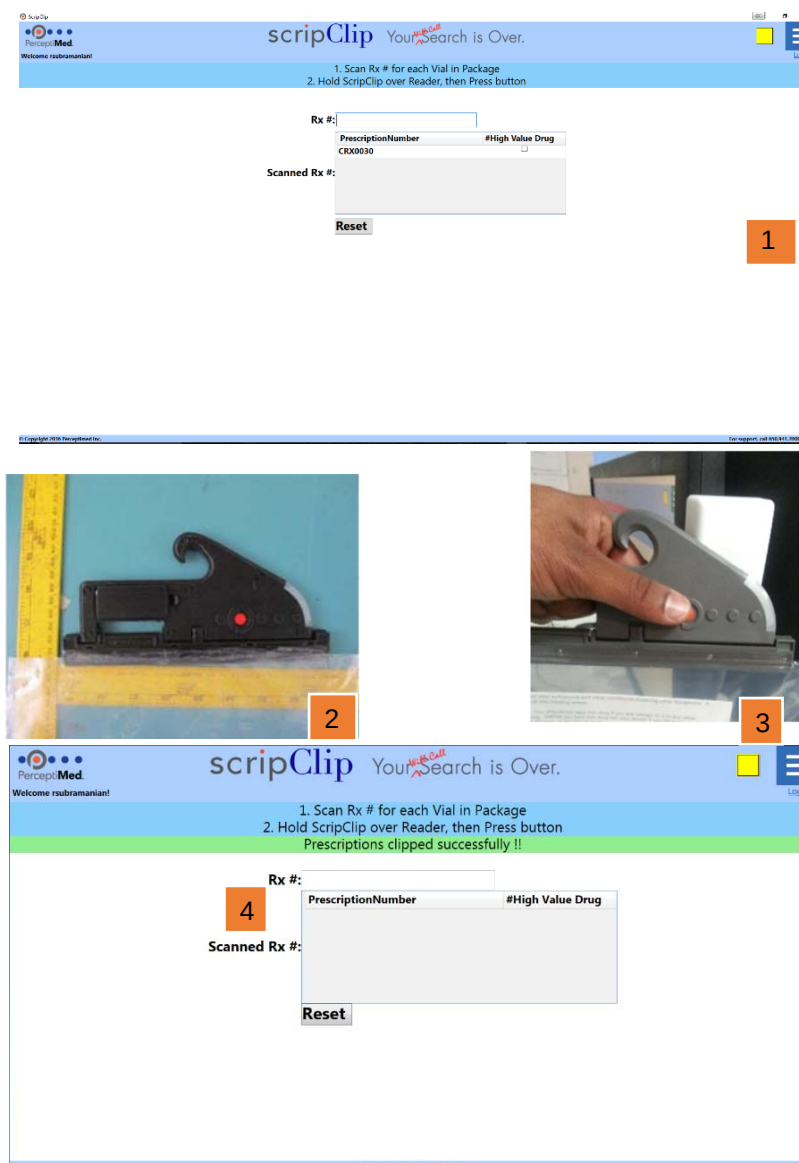


Packaging the Prescription

Clip-on scripClips :

- 1 Using the barcode scanner, scan one by one, up to 10 Rx numbers either from the vial or from the prescription paperwork.
→ Scanned Rx numbers appear on the screen.
- 2 Place the vials and prescription paperwork in the prescription bag.
- 3 Take an unused Clip-on scripClip **3a**, hold near the scripClip Reader **3b** and press the orange button on the scripClip.
→ The scripClip will flash the station color once.
→ scripClip application will display, 'Prescriptions clipped successfully.'
- 4 → scripClip application will display, 'Prescriptions clipped successfully.'
- 5 Attach the scripClip to a closed prescription bag. Place the prescription bag in the will-call area.

The screenshot shows the scripClip application interface. At the top, it says "scripClip Your Search is Over." and "Logged in User (IdentRx)". Below this, there are two instructions: "1. Scan Rx # for each Vial in Package" and "2. Hold ScripClip over Reader, then Press button". In the center, there is a form with the following fields: "Rx #:" (empty), "Scanned Rx #:" (empty), "PrescriptionNumber" (CRX0030), "#PrescriptionNumber" (CRX0030), and "#High Value Drug" (checkbox). Below the form is a "Reset" button. A small orange box with the number "1" is in the bottom right corner.

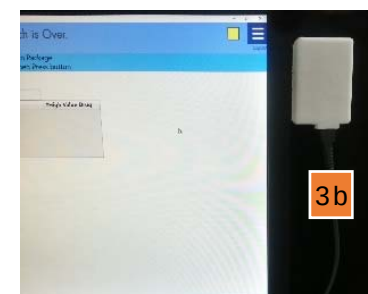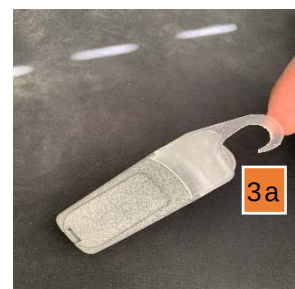

The screenshot shows the scripClip application interface after successful completion. At the top, it says "scripClip Your Search is Over." and "Welcome rsutramanian!". Below this, there are two instructions: "1. Scan Rx # for each Vial in Package" and "2. Hold ScripClip over Reader, then Press button". Below the instructions, it says "Prescriptions clipped successfully!!". In the center, there is a form with the following fields: "Rx #:" (empty), "Scanned Rx #:" (empty), "PrescriptionNumber" (CRX0030), and "#High Value Drug" (checkbox). Below the form is a "Reset" button. A small orange box with the number "3" is in the bottom right corner.



Point of Sale: Locate a Prescription

Locate a prescription:

- 1 Search for the Patient using name or date of birth (or both).
- 2 Press "Locate" to activate the scripClip's flashing light or "Locate with Tone" to activate the scripClip's flashing light and beep.
 - scripClip application shows "Located" and scripClip will flash for 30 seconds. **2a**
 - Retrieve the package with the flashing light that matches the station color.
 - Press the button on the scripClip to stop the flashing light.
 - Hold the Hanging Bag scripClip **2c** near the scripClip reader, and push the orange scripClip button. Or remove the Clip-on scripClip **2d** from the prescription bag, hold the scripClip near the scripClip Reader, and push the scripClip button
 - scripClip will flash the station color once to indicate transaction completed.
 - scripClip application will display 'Package Retrieved.' **2e**
 - Place empty scripClip in storage bin for use on next prescription package.



2c



2d

1

First Name	Last Name	Date of Birth	Phone	Zip	Pkg(s)
John	Miller	2/16/1966	(999)-555-2222	88445	1
John	Miller	8/26/1961	(999)-555-2222	88445	2

2a

Rx IDs	Status
CRX0030	Located in Store
CRX003P	Located in Store

2b

Rx IDs	Status
CRX001R	Package Retrieved




2e

Returning Prescriptions to Stock

Return to stock:

- 1 Select number of days for "Prescriptions older than Days" from pull down menu. Press "Search".
- 2 Scroll to bottom of list and press "Locate Packages" or "Locate Packages with Tone"
- 3 scripClip application shows "Located in Store" for each prescription on list and the scripClips flash for 30 seconds. Retrieve all of the packages with the flashing light that matches the station color.
→ Press the button on the scripClip to stop the flashing light.
- 4 Hold the Hanging Bag scripClip near the scripClip Reader, and push the scripClip button (or) Remove the Clip-on scripClip from the bag, hold the scripClip near the scripClip Reader, and push the scripClip button.
→ scripClip will flash the station color once to indicate transaction completed.
- 5 → scripClip application will display 'Package Retrieved.'
→ Place the empty scripClips in the storage bin for use on next prescription package.
→ Repeat for each packaged to be returned to stock



Select "Packaging" from the pull-down menu in the upper right corner to return to the Point-of-Sale page.

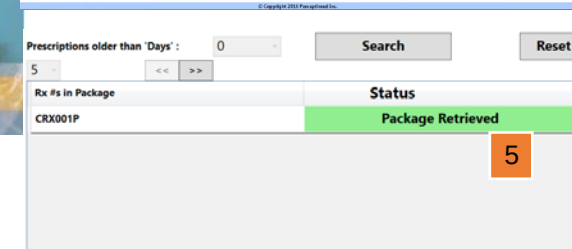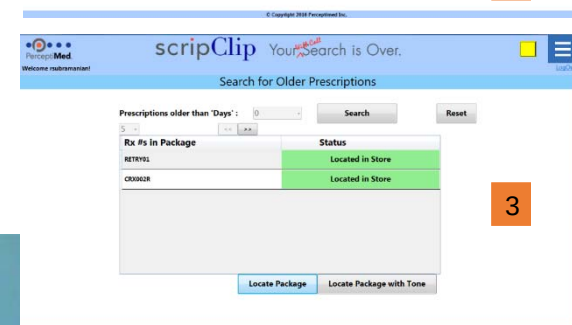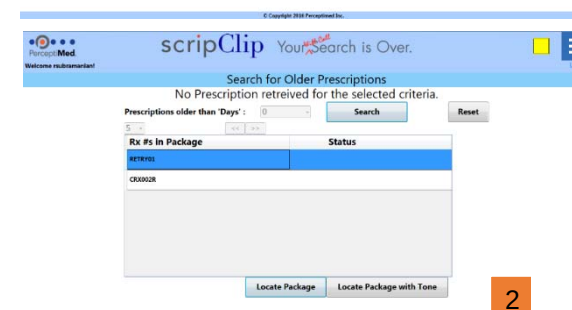

scripClip Tamper Alert (Admin Mode)

Red Alert Message Issued:

- 1 Removing a scripClip from a package without using patient search or the Return To Stock screens causes scripClip to issue an alert. When the 'Alert' message appears, click on the *Alert* button, select the 'Locate' command to activate that scripClip.
- 2 After locating the scripClip, select the Action marked, 'Prescription OK' if the scripClip was accidentally removed or dislodged from the package. If the associated package indicates tampering or is missing, select the Action marked, 'Tampered/Missing Prescription'

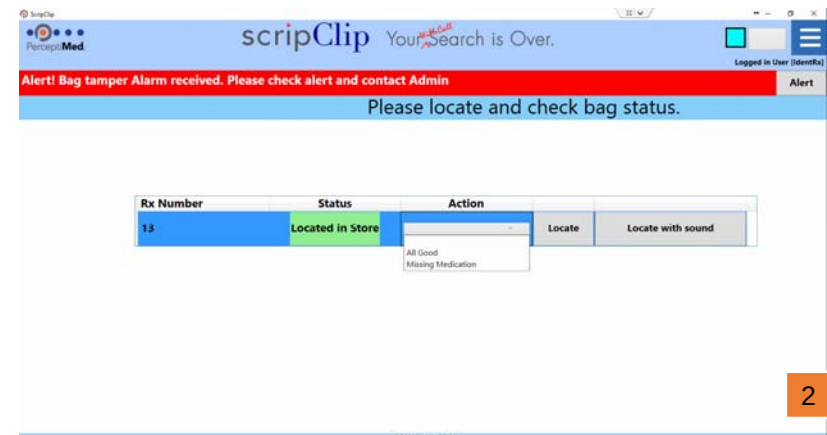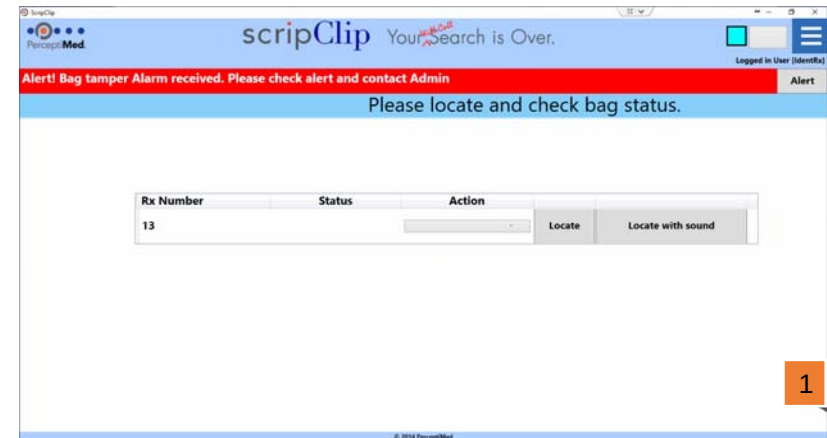

Troubleshooting

Message: **Try again**

Cause: Package not located in store

Solution: Check Rx numbers entered. Press “Locate” or “Locate with Tone” again.

Message: **Wrong package retrieved.**

Cause: scripClip scanned does not match Rx numbers.

Solution: Put scripClip back on package and return to shelf. Press “Locate” again. Retrieve correct package.

Message: **Unable to read Rx numbers.**

Cause: scripClip reader did not receive all Rx numbers from scripClip

Solution: Hold scripClip closer to the scripClip Reader and press the orange button on scripClip again.

Event: **User is unable to find flashing scripClip in store**

Cause: scripClip stops flashing after 30 seconds.

Solution: Press “Locate Prescriptions” again or press “Locate Prescriptions with Tone”. User does not need to re-enter Rx numbers.

Message: **Unable to communicate with the central service**

Cause: Package not located in store

Solution: Call 650 941-7000 ext 5 or send email to

support@perceptiMed.com

Appendix 1: Parts in the scripClip Kit

1 Hanging-Bag scripClip

2 Clip-on scripClip

3 scripClip Reader

4 scripClip Central Computer

5 scripClip Antenna

Quantity of clips and accessories as determined by customer requirements



Information



Battery Replacement:

Do not attempt to replace the battery yourself. Perceptimed will replace battery and verify device as directed by service contract.

Use Restriction:

Do not relocate system and components without a Perceptimed authorized installer. Components are restricted to indoor pharmacy and non-residential application.

Disposal Information:

When this product reaches its end of life, return the device(s) to Perceptimed or dispose as directed by local laws and regulations. The collection and recycling of electronic waste conserves natural resources and protects human health and the environment.

Compliance Information



FCC Compliance Statement:

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is not recommended. If this equipment does cause harmful interference or is interfered with by other services contact Perceptimed or its representative.

Any modifications made to this device and accessories that are not approved by Perceptimed, Inc. may void the authority granted to the user by the FCC to operate this equipment.

Applicable to the Following Devices:

Hanging-Bag scripClip, Model SCLP-HNGBAG-G3. FCC ID:2AJ8D-HNGBAG-G3

Compliance Information



Disclosure of Exposure to Radiofrequency and Electromagnetic Fields:

This equipment complies with the FCC radiation exposure limits set forth per 47 CFR §2.1091 for an uncontrolled environment (ie: both occupational and unrestricted exposure to the public).

The operating power of the scripClip and its associated Reader is deliberately limited. No precautions are required for their use or operation.

The wide area broadcast station and its antennas must be installed and operated such that at least 30cm separation distance is maintained between these components and the operator or the general public. This is normally achieved by overhead ceiling mount of at least 20cm from the level of the floor. These antennas may not be collocated with antennas for WiFi or other applications.

Compliance Information



Industry Canada Compliance Statement:

Applicable to the Following Devices:

Hanging-Bag scripClip - Model SCLP-HNGBAG-G3, FCC ID:2AJ8D-HNGBAG-G3, IC: 23412-HNGBAGG3

Section 8.4 of RSS-GEN

This Device complies with Industry Canada License-exempt RSS standard(s). Operation is subject to the following two conditions: 1) this device may not cause interference, and 2) this device must accept any interference, including interference that may cause undesired operation of the device.

Section 8.3 of RSS-GEN

Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain provided by Perceptimed and approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain is chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that necessary for successful communication.

Antenna types not provided by Perceptimed and approved by Industry Canada are strictly prohibited for use with this device.

Radiation Exposure Statement: This equipment complies with the IC RSS-102 radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20 cm between the radiator and your body.

Informations De Conformité



Industrie Canada Déclaration De Conformité:

Applicable pour les Appareils Suivants:

Hanging-Bag scripClip - Model SCLP-HNGBAG-G3, FCC ID:2AJ8D-HNGBAG-G3, IC: 23412-HNGBAGG3

La Section 8.4 du RSS-GEN

Cet Appareil est conforme à la Licence d'Industrie Canada RSS exemptés de la norme(s). Son fonctionnement est soumis aux deux conditions suivantes: 1) cet appareil ne peut pas causer d'interférences, et 2) cet appareil doit accepter toute interférence, y compris les interférences pouvant causer un fonctionnement indésirable de l'appareil.

La Section 8.4 du RSS-GEN

En vertu de la réglementation d'Industrie Canada, le présent émetteur radio peut fonctionner uniquement à l'aide d'une antenne d'un type et le maximum (ou le) gain fourni par Perceptimed et approuvé pour l'émetteur par Industrie Canada. Pour réduire les risques d'interférences radio à d'autres utilisateurs, le type d'antenne et son gain est choisi que la puissance isotrope rayonnée équivalente (e.i.r.p.) n'est pas supérieure à celle nécessaire pour la réussite de la communication

Types d'antennes ne sont pas fournis par Perceptimed et approuvé par Industrie Canada sont strictement interdits pour une utilisation avec cet appareil.

Énoncé d'exposition aux rayonnements: Cet équipement est conforme aux limites d'exposition aux rayonnements ioniques RSS-102 Pour un environnement incontrôlé. Cet équipement doit être installé et utilisé avec un Distance minimale de 20 cm entre le radiateur et votre corps.