bebird

Enjoy Ear clean,Easy to use

Bebird Smart Visual Ear-Clean Rod A2


bebird

Enjoy Ear clean,Easy to use

Start Bebird A2

1 Scan QR code to download APP

2 Open APP
Click Bebird A2

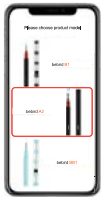


3 Open the cover
Power on automatically



Single mode


1 Connect Product Wi-Fi
"bebird-XXXX"

2 Back to Bebird app
Enter and start to use it

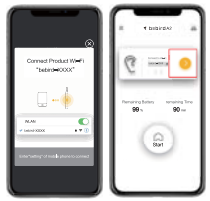

Accessory replacement

Clipped end of ear spoon
It can be replaced once it's pulled out.

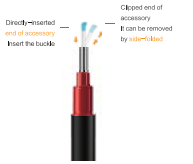

Accessory

1 Q-Elastic Ear Spoon × 4
Optical Fiber PC and Silica Gel Synthetic Europoon

2 Annular Ear Spoon × 2
Dig out the earwax
clean it easily from top to bottom

3 Deep Ear Digging Spoon × 2
Deep into the ear canal
to excavate hidden corners

4 Enjoyable cotton swab × 5
Clean and remove earwax debris

5 Silver Needle Crane Hair stick × 1
Massage the ear with it as gently and relax

6 Silver Needle Crane Hair stick × 1
Massage external auditory canal to stimulate sympathetic nerve

7 Silver Needle Horsetail stick × 1
Massage tympanic membrane to enjoy high purity white noise

8 Alcohol cotton swab × 5
For lens cleaning and ear canal disinfection



Troubleshooting

Trouble	Solution
starting up and no response	Please use the attached charging cable for quick charging The charging time is 1.5hr
Camera does not turn on	Close and open the cover to try a few times
APP quit quickly	Delete the existing APP and re-scan the QR code for installation
Installation of tools in accessory barrel is not smooth	Firstly align the ear-cleaning tool with the lens sleeve, then install it into the lens sleeve
Rod's power indicator is red	Please use the attached charging cable for charging. The charging time is 1.5hr
The display on the app is not reading properly.	Disconnect the phone with product and open cover to restart it,then set phone to reconnect product.
WiFi connection is unstable, sometimes it is disconnected in seconds	Close and open the cover to restart rod and restart app in phone to connect rod

If you still have problems with the above procedures, please contact our Customer Service Representatives directly for further assistance

Hazardous substance

Part name	Hazardous material					
	(Pb)	(Hg)	(Cd)	(Cr(VI))	(PBB)	(PBDE)
Printed Circuit Components	×	○	○	○	○	○
Case	○	○	○	○	○	○
Plastic rubber accessory	○	○	○	○	○	○

This form is formulated in accordance with SJ/T 11364
○: represent the hazardous material ratio is below GB/T26572-2011
×: represent the hazardous material ratio is over GB/T26572-2011

bebird

Product parameters

Brand	bebird
Product Name	Bebird Smart Visual Ear-Clean Rod
Model	A2
Rod weight	29g
Network standard	IEEE 802.11b/g/n
Antenna	Internal FPC antenna
Working frequency	2.4Ghz
Image transmission rate	30fps
image sensor	CMOS
Working temperature	-10 ~ 50℃
Battery	350mAh lithium battery
Battery life	about 90 minutes
Charging time	1.5h
Input current	DC 5V 420mAh
Lens diameter	3.2mm
Optimum focal length	1.5 ~ 2cm
Pixel	5M
Gravity Sensor	3 axis

Certification

Attention

1. When cleaning lens, please wipe it carefully with professional alcohol cotton swab.
2. Before use it, please check the surrounding environment. Do not use the product in the environment where someone runs, in order to avoid strike.
3. This product is not suitable for children under one year old.
4. It is forbidden for children to use this product alone in order to avoid accidental harm.
5. Do not use the endoscope in liquid to avoid damage caused by immersion.
6. Please use our product once a month and charge it,in order to keep the rechargeable polymer Li-battery in good condition.
7. Avoid product exposure to sunshine directly, especially the accessories in the barrel, in order to avoid softening.
8. During using it, the temperature of the product will rise slightly (up to 35 C). Please rest assured that it can be used.
9. In order to ensure the best usage, please keep the child in a sitting position and operate the product horizontally.

Warranty

User info	
User name:	Telephone number:
Address:	
Seller:	Telephone number:

Warning:

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.
If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
-- Reorient or relocate the receiving antenna.
-- Increase the separation between the equipment and receiver.
-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
-- Consult the dealer or an experienced radio/TV technician for help.
NOTE: This device and its antenna(s) must not be co-located or operation in conjunction with any other antenna or transmitter