

Xlink 2.0

Wireless Monitor Bridge

User Manual


Scan or visit to get started:



xchart.com/xlink

Dongle


Front Side

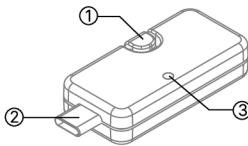

- ① FIND PAIR Button: Press to identify the paired adapter.
- ② USB-C Connector: Used to plug the dongle into your computer.
- ③ ACTIVITY LED: Indicates when the paired adapter is within range.


Adapter

Front Side

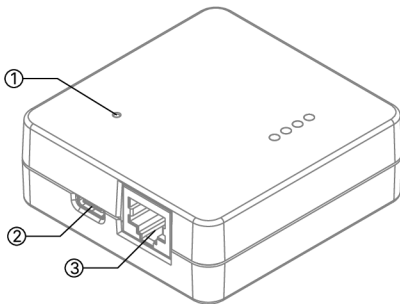

- ① CONFIG Button: Press to enter configuration mode.
- ② POWER Connector: Supplies power to the adapter.
- ③ ETHERNET Connector: Used to connect the adapter to a medical device via Ethernet.

Back Side

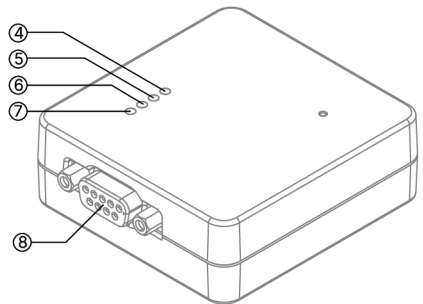

- ④ PAIRED DONGLE PRESSED: Glows if FIND PAIR button pressed.
- ⑤ ETHERNET MODE: Glows when configured for Ethernet mode.
- ⑥ SERIAL MODE: Glows when configured for Serial mode.
- ⑦ DATA ACTIVITY: Blinks when sending or receiving data.
- ⑧ SERIAL Connector: Used to connects the adapter to a medical device via a Serial port.

Xchart.com

Customer Service: support@xchart.com

Xchart, Inc. 304 W Pacific Ave. #210, Spokane, WA 99201

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

This device may not cause harmful interference, and

This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device and an intentional radiator, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

To maintain compliance with FCC's RF Exposure guidelines, This equipment should be installed and operated with minimum distance between 20cm the radiator your body: Use only the supplied antenna.