


Proposed Labeling

This proposed labeling section includes representative copies of the draft product label, draft Instructions for Use, and draft box labeling.

Draft Product Label

Draft product label shown below, to be placed on the back of the CX00 Monitor. The red box indicates the unique section. Smallest label font size: 4pt.

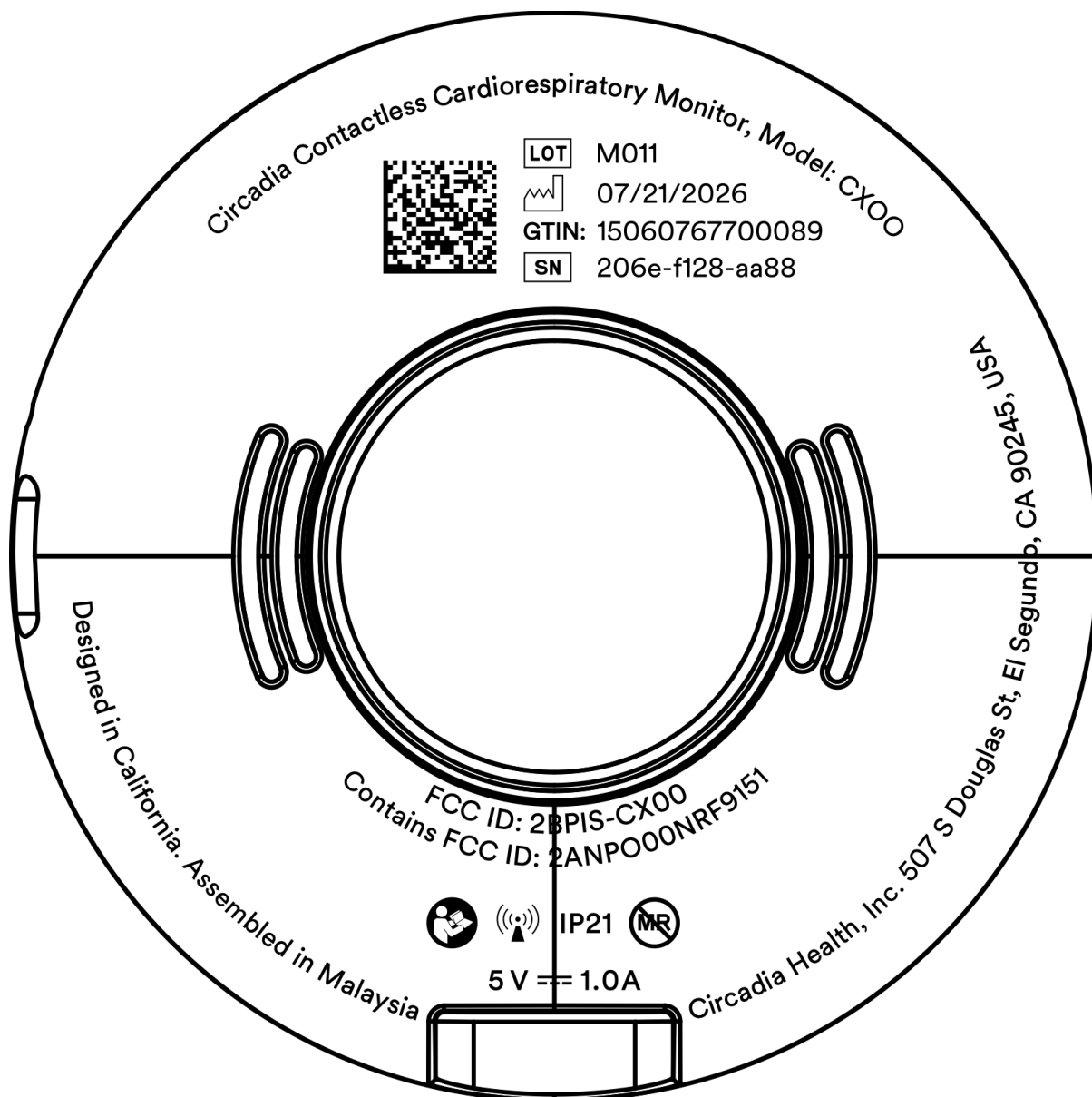

Product labeling will all be laser-engraved directly onto the surface of the enclosure plastic using a 355 nm UV laser. The section below the name of the device and model number is the unique section explained below.

Unique section

On-device labeling to be positioned on the back of the CX00 Monitor, just underneath the name and model number. It will be laser-engraved directly onto the surface of the enclosure plastic using a 355 nm UV laser. Smallest label font size: 4pt.




LOT M011
 **07/21/2026**
GTIN: 15060767700089
SN 206e-f128-aa88

The label will contain:


- GS1 datamatrix, contain all relevant fields (SN#, GTIN, DoM, LOT#)
- LOT: 4-character LOT number
- Date of manufacture (MM/DD/YYYY)
- GTIN: 14-character GS1 GTIN for CX00
- SN: 12-character serial number

Power Adapter Labeling

Power adapter uses the original manufacturer labeling with additional Circadia accessory GTIN and LOT #.

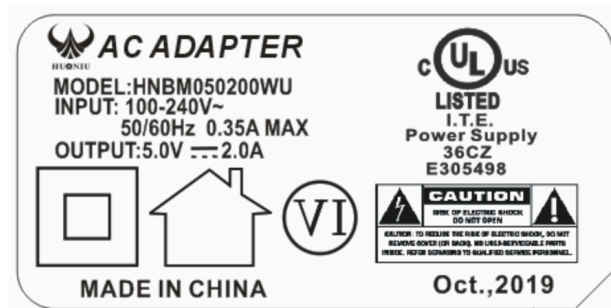

Original manufacturer label on adapter

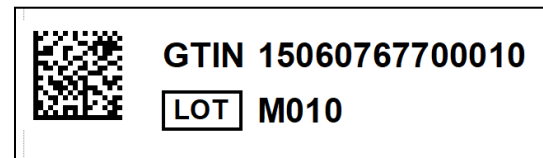

Additional Circadia label on adapter


Draft Box Labeling

Master cartons

The default packaging will be master cartons of 100 x CX00 Monitors.

Each carton will have the following external labeling (dimensions in the example below are for the 100 pcs carton):


- 2 long sides will contain Circadia logo + text
- 1st short side will contain Circadia logo only
- 2nd short side will contain regulatory markings, manufacturer information and contents

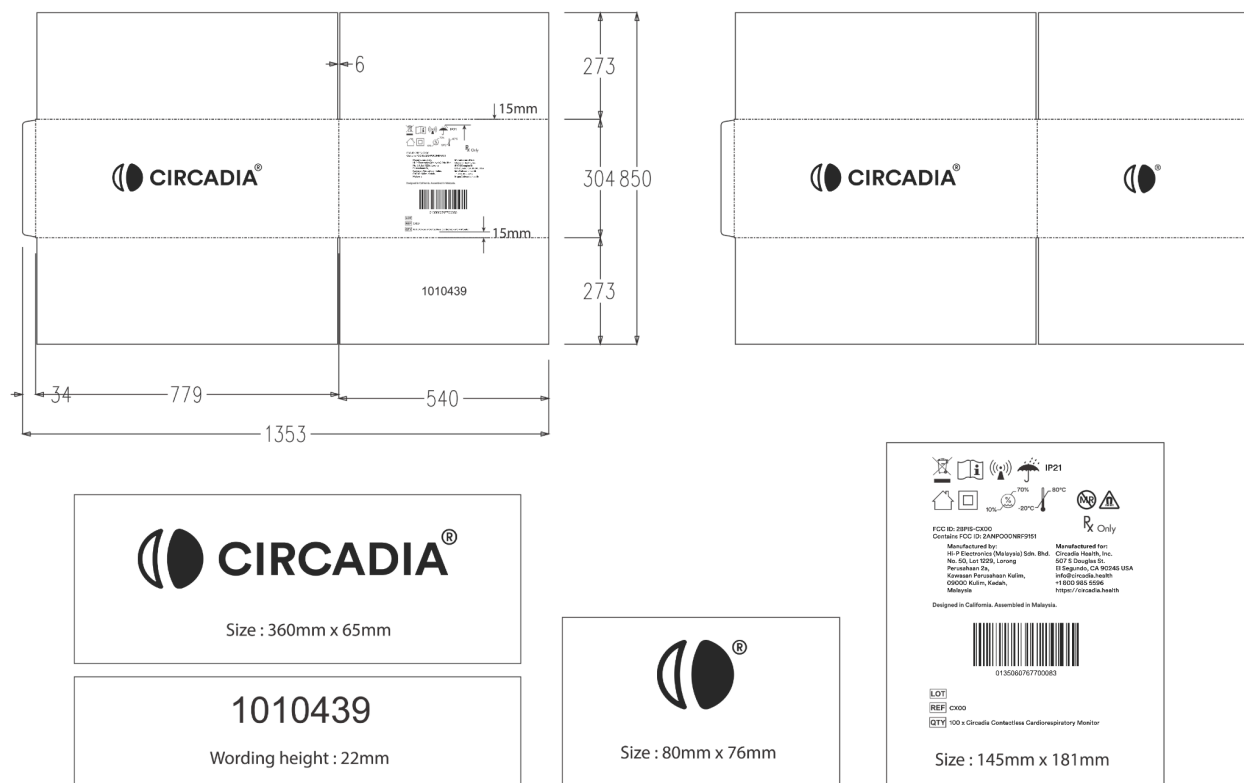

Draft box labeling, for the 100-pcs carton containing 40 CX00 Monitors. Box dimensions and label layout are shown (top), and detailed design of box panels (bottom).

Master Cartons Label - Detailed View

Detailed view of box labeling - 2nd short side of master carton:

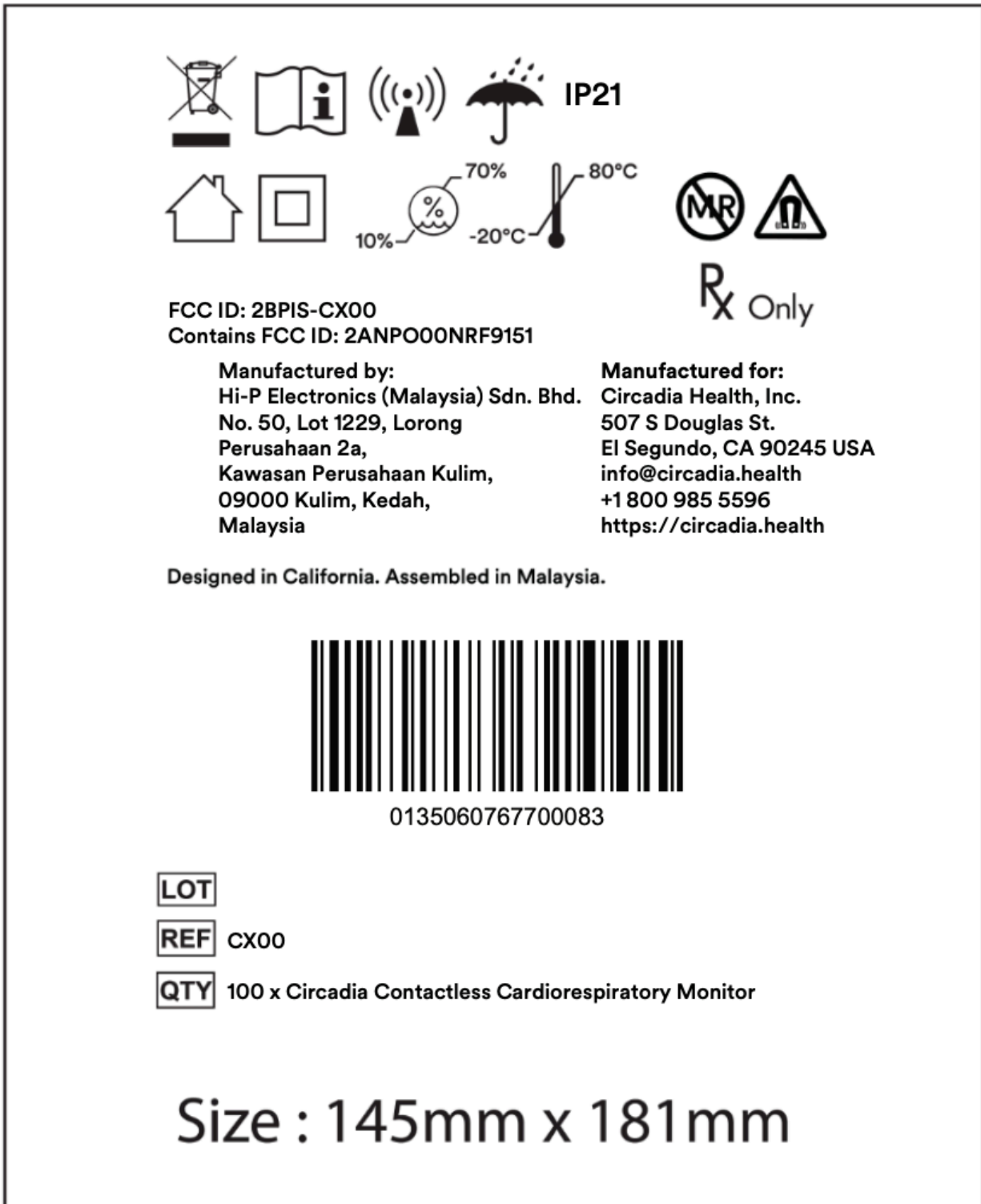

The master carton label contains the following (from top to bottom):

- Regulatory symbols
- Manufacturer details
- GTIN barcode for given size of master carton (different between sizes)
- LOT#
- Product reference
- Quantity of what's inside the box


Ziploc Bags

Individual CX00 Monitors will be packed in bubble-padded Ziploc bags. Ziploc bags will have individual draft labeling as below:


Front view



Back view

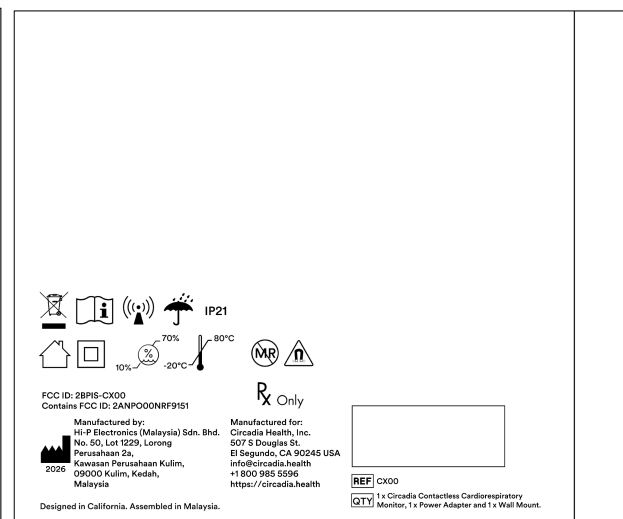

Detailed Ziploc bag labeling:

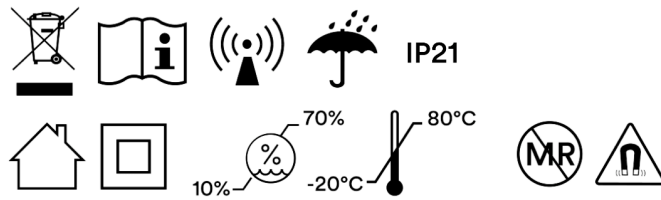

FCC ID: 2BPIS-CX00
Contains FCC ID: 2ANPO00NRF9151

Manufactured by:
Hi-P Electronics (Malaysia) Sdn. Bhd.
No. 50, Lot 1229, Lorong
Perusahaan 2a,
Kawasan Perusahaan Kulim,
09000 Kulim, Kedah,
Malaysia

2026

R_x Only
Manufactured for:
Circadia Health, Inc.
507 S Douglas St.
El Segundo, CA 90245 USA
info@circadia.health
+1 800 985 5596
<https://circadia.health>



REF CX00

QTY 1 x Circadia Contactless Cardiorespiratory
Monitor, 1 x Power Adapter and 1 x Wall Mount.

Designed in California. Assembled in Malaysia.

The position indicated by a red box above in the above label draft will contain a label as below, with fields for:

- GTIN (01)
- Date of manufacture (11)
- Serial number (21)
- LOT number (10)

