

Client

ITS / ETL

CALIBRATION CERTIFICATE

Object(s)

ET3DV6 - SN:1576

Calibration procedure(s)

QA CAL-01.v2
Calibration procedure for dosimetric E-field probes

Calibration date:

August 29, 2003

Condition of the calibrated item

In Tolerance (according to the specific calibration document)

This calibration statement documents traceability of M&TE used in the calibration procedures and conformity of the procedures with the ISO/IEC 17025 international standard.

All calibrations have been conducted in the closed laboratory facility: environment temperature 22 +/- 2 degrees Celsius and humidity < 75%.

Calibration Equipment used (M&TE critical for calibration)

Model Type	ID #	Cal Date (Calibrated by, Certificate No.)	Scheduled Calibration
RF generator HP 8684C	US3642U01700	4-Aug-99 (SPEAG, in house check Aug-02)	In house check: Aug-05
Power sensor E4412A	MY41495277	2-Apr-03 (METAS, No 252-0250)	Apr-04
Power sensor HP 8481A	MY41092180	18-Sep-02 (Agilent, No. 20020918)	Sep-03
Power meter EPM E4419B	GB41293874	2-Apr-03 (METAS, No 252-0250)	Apr-04
Network Analyzer HP 8753E	US37390585	18-Oct-01 (Agilent, No 24BR1033101)	In house check: Oct 03
Fluke Process Calibrator Type 702	SN: 6295803	3-Sep-01 (ELCAL, No 2360)	Sep-03

Calibrated by:

Name

Katja Pokovic

Function

Laboratory Director

Signature

Approved by:

Niels Kuster

Quality Manager

Date issued: September 1, 2003

This calibration certificate is issued as an intermediate solution until the accreditation process (based on ISO/IEC 17025 International Standard) for Calibration Laboratory of Schmid & Partner Engineering AG is completed.

Probe ET3DV6

SN:1576

Manufactured:	April 6, 2001
Last calibration:	February 27, 2002
Recalibrated:	August 29, 2003

Calibrated for DASY Systems

(Note: non-compatible with DASY2 system!)

DASY - Parameters of Probe: ET3DV6 SN:1576

Sensitivity in Free Space

Diode Compression

NormX	1.81 $\mu\text{V}/(\text{V}/\text{m})^2$	DCP X	96	mV
NormY	1.83 $\mu\text{V}/(\text{V}/\text{m})^2$	DCP Y	96	mV
NormZ	1.77 $\mu\text{V}/(\text{V}/\text{m})^2$	DCP Z	96	mV

Sensitivity in Tissue Simulating Liquid

Head 900 MHz $\epsilon_r = 41.5 \pm 5\%$ $\sigma = 0.97 \pm 5\%$ mho/m

Valid for f=800-1000 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	7.0 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	7.0 $\pm 9.5\%$ (k=2)	Alpha	0.24
ConvF Z	7.0 $\pm 9.5\%$ (k=2)	Depth	3.61

Head 1800 MHz $\epsilon_r = 40.0 \pm 5\%$ $\sigma = 1.40 \pm 5\%$ mho/m

Valid for f=1710-1910 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	5.4 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	5.4 $\pm 9.5\%$ (k=2)	Alpha	0.41
ConvF Z	5.4 $\pm 9.5\%$ (k=2)	Depth	2.73

Boundary Effect

Head 900 MHz Typical SAR gradient: 5 % per mm





Probe Tip to Boundary		1 mm	2 mm
SAR _{be} [%]	Without Correction Algorithm	10.2	6.5
SAR _{be} [%]	With Correction Algorithm	0.6	0.6



Head 1800 MHz Typical SAR gradient: 10 % per mm

Probe Tip to Boundary		1 mm	2 mm
SAR _{be} [%]	Without Correction Algorithm	11.4	8.0
SAR _{be} [%]	With Correction Algorithm	0.2	0.2


Sensor Offset

Probe Tip to Sensor Center	2.7	mm
Optical Surface Detection	1.6 $\pm$ 0.2	mm

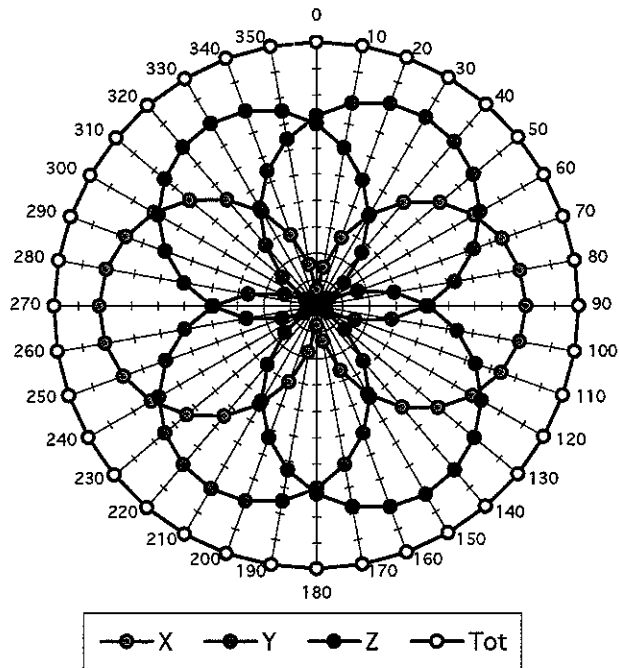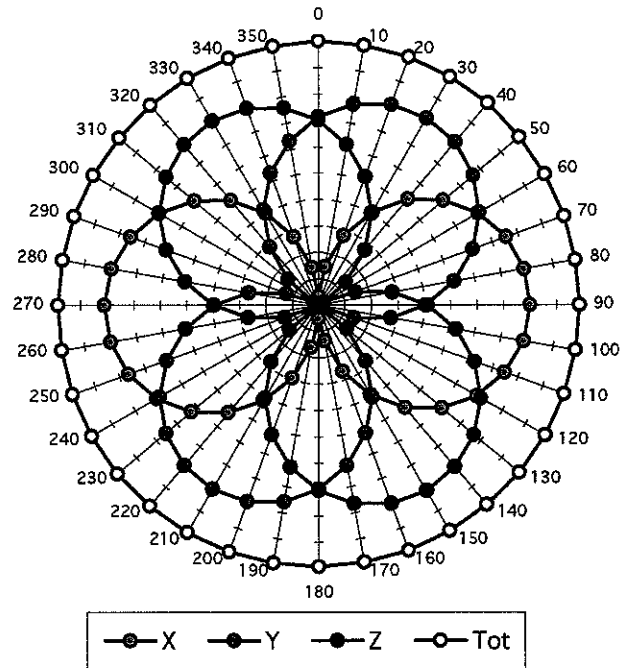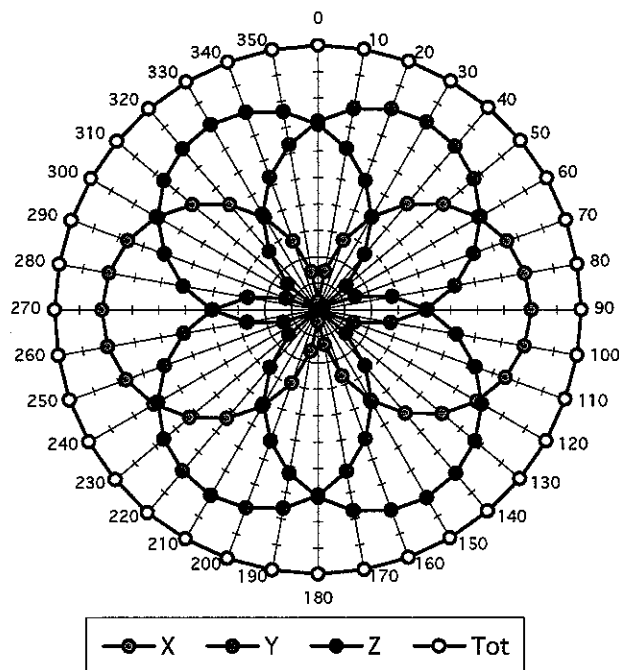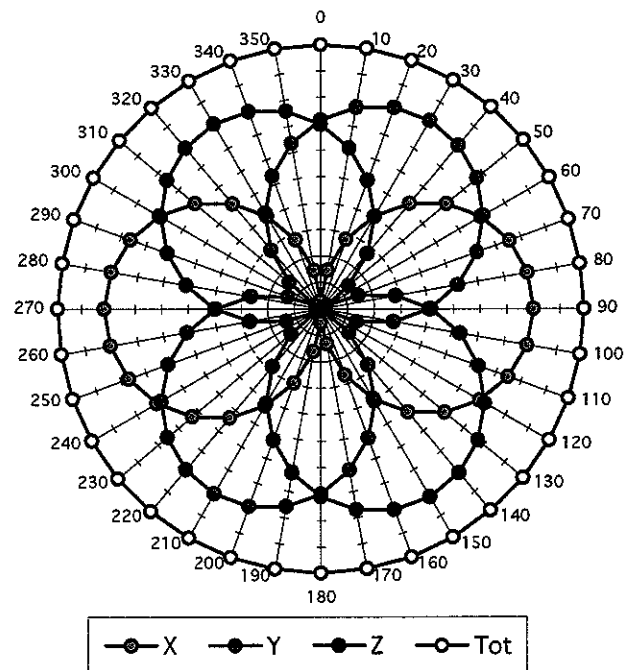Receiving Pattern (ϕ), $\theta = 0^\circ$ **f = 30 MHz, TEM cell if110****f = 100 MHz, TEM cell if110****f = 300 MHz, TEM cell if110****f = 900 MHz, TEM cell if110**

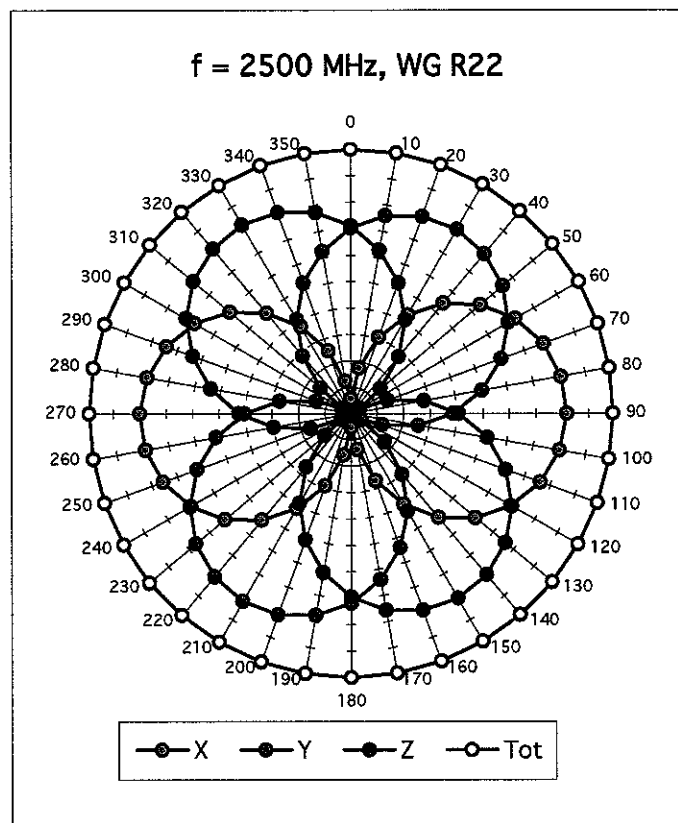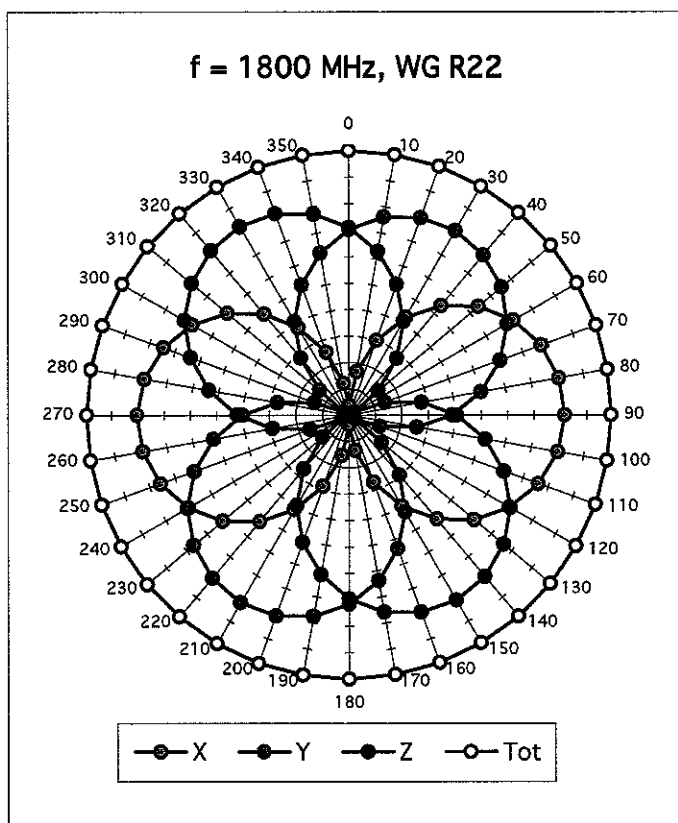


Isotropy Error (ϕ), $\theta = 0^\circ$


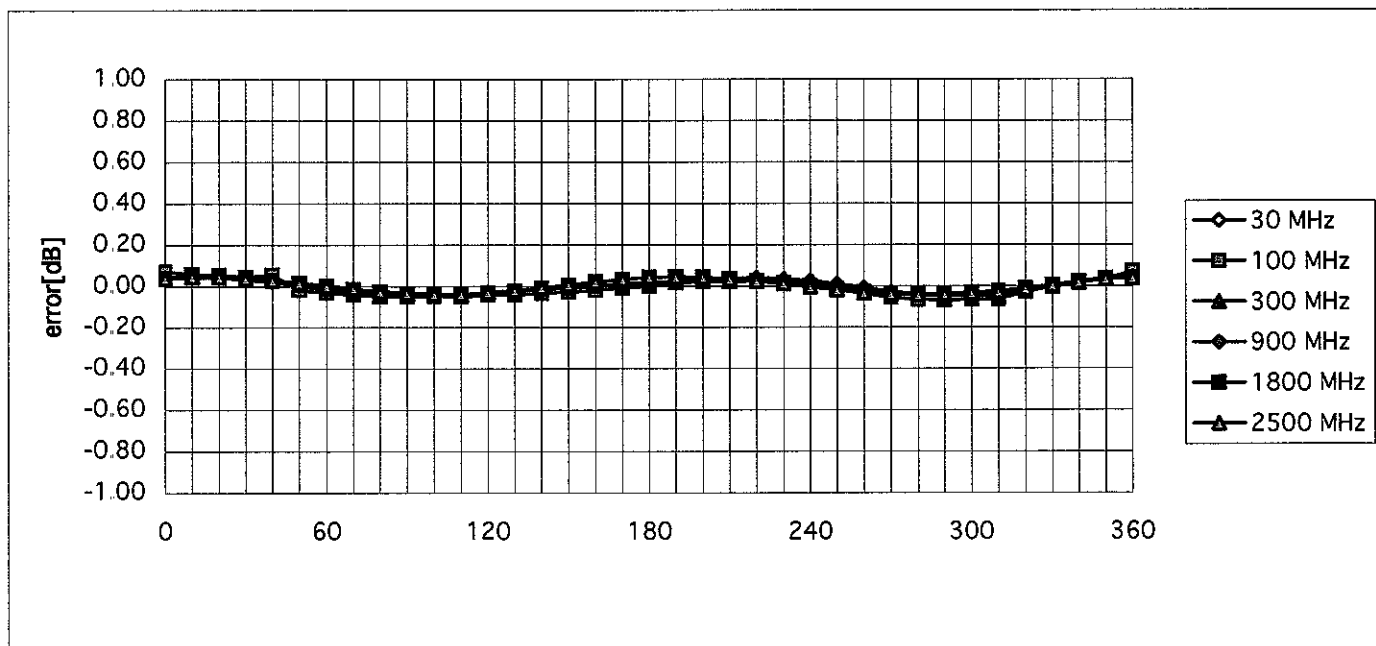


Frequency Response of E-Field



(TEM-Cell:ifi110, Waveguide R22)

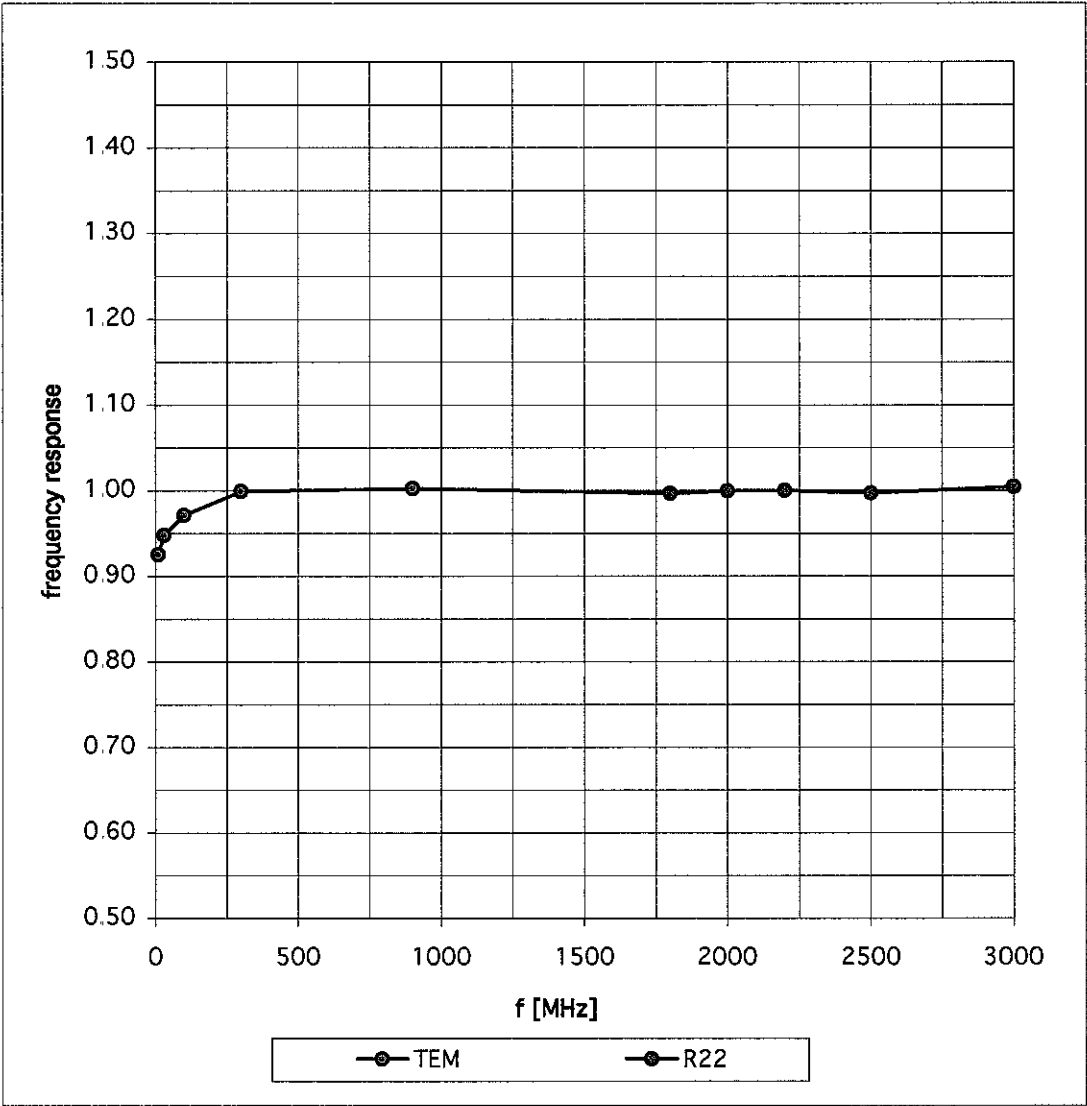

Dynamic Range f(SARhead)
(Waveguide R22)

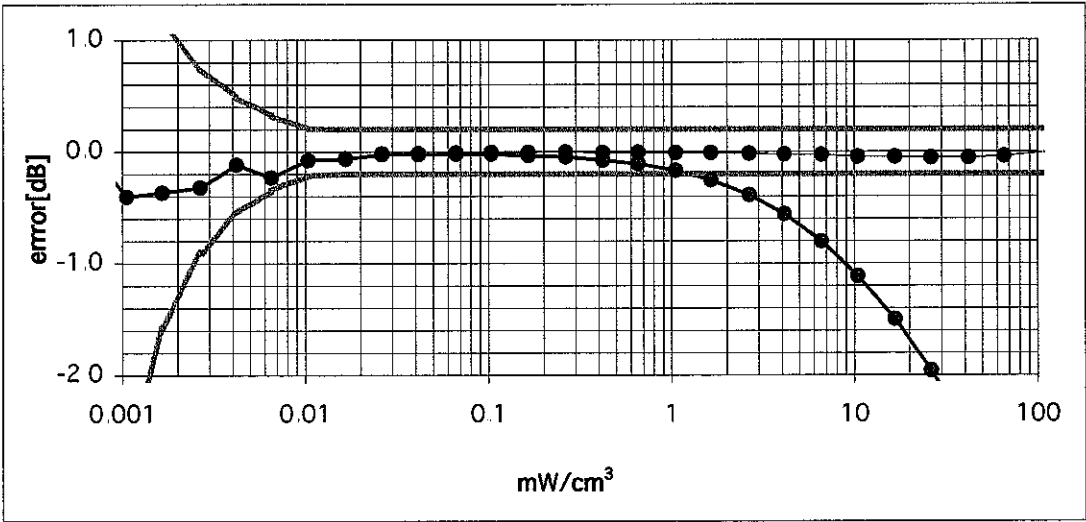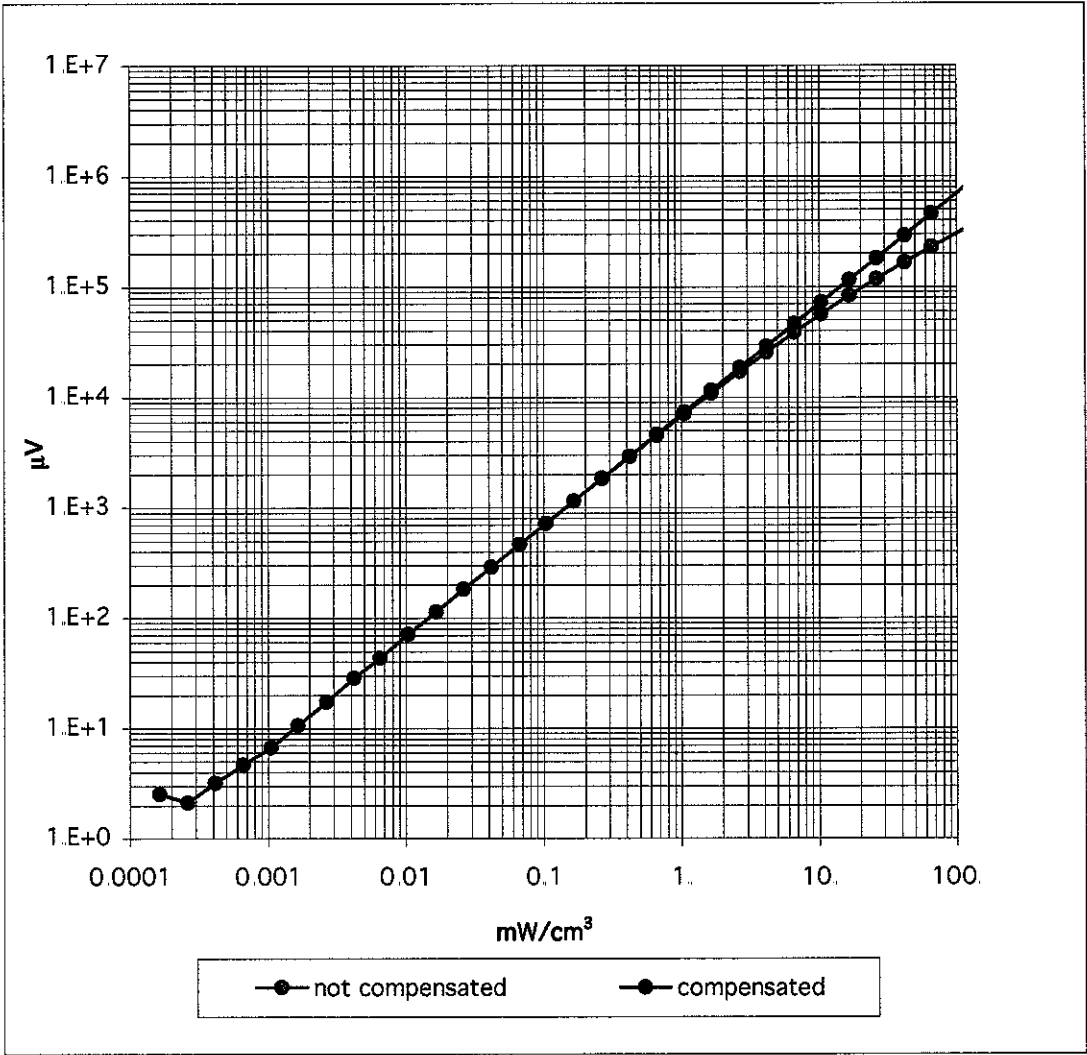

Conversion Factor Assessment

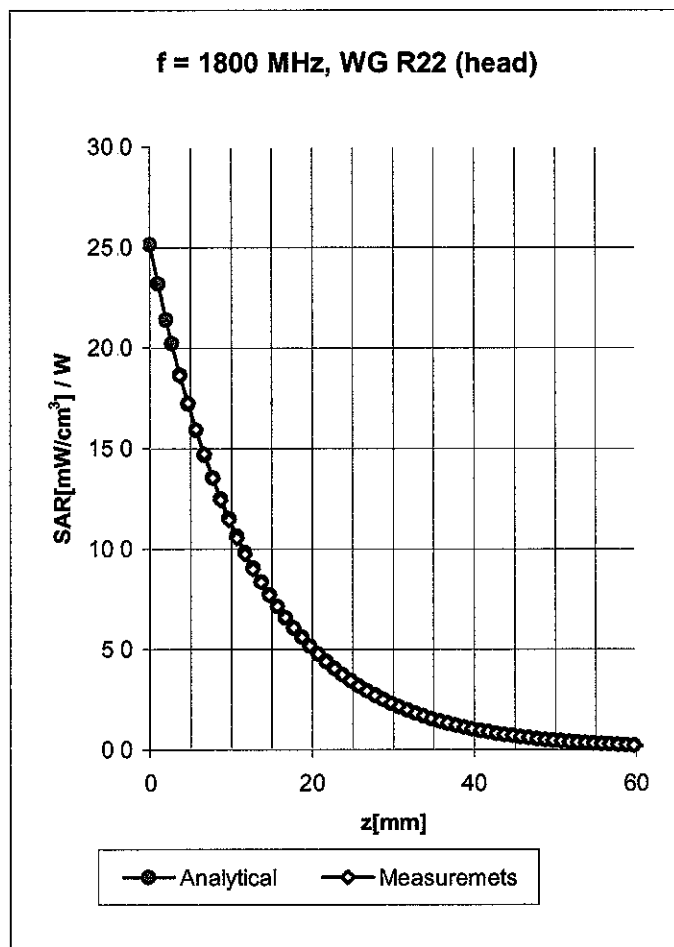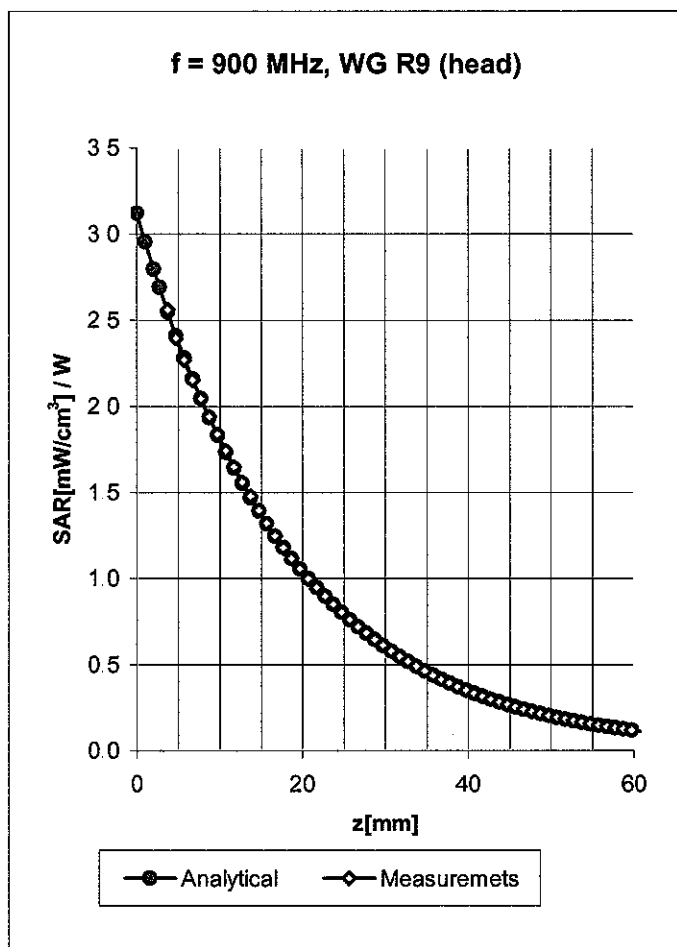

Head 900 MHz $\epsilon_r = 41.5 \pm 5\%$ $\sigma = 0.97 \pm 5\%$ mho/m

Valid for f=800-1000 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X



ConvF X	7.0 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	7.0 $\pm 9.5\%$ (k=2)	Alpha	0.24
ConvF Z	7.0 $\pm 9.5\%$ (k=2)	Depth	3.61

Head 1800 MHz $\epsilon_r = 40.0 \pm 5\%$ $\sigma = 1.40 \pm 5\%$ mho/m

Valid for f=1710-1910 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	5.4 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	5.4 $\pm 9.5\%$ (k=2)	Alpha	0.41
ConvF Z	5.4 $\pm 9.5\%$ (k=2)	Depth	2.73

Conversion Factor Assessment

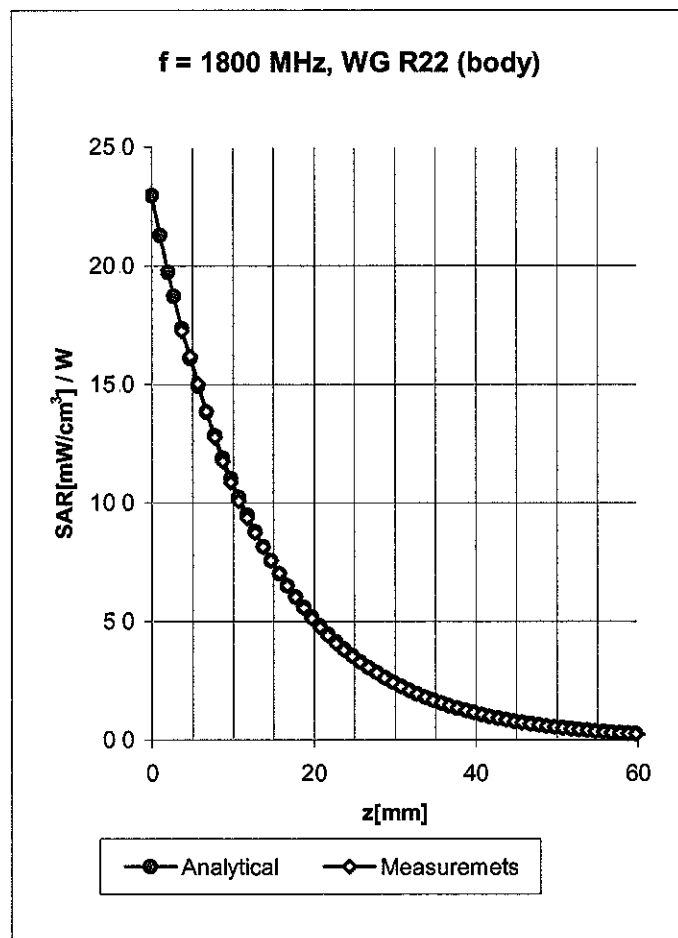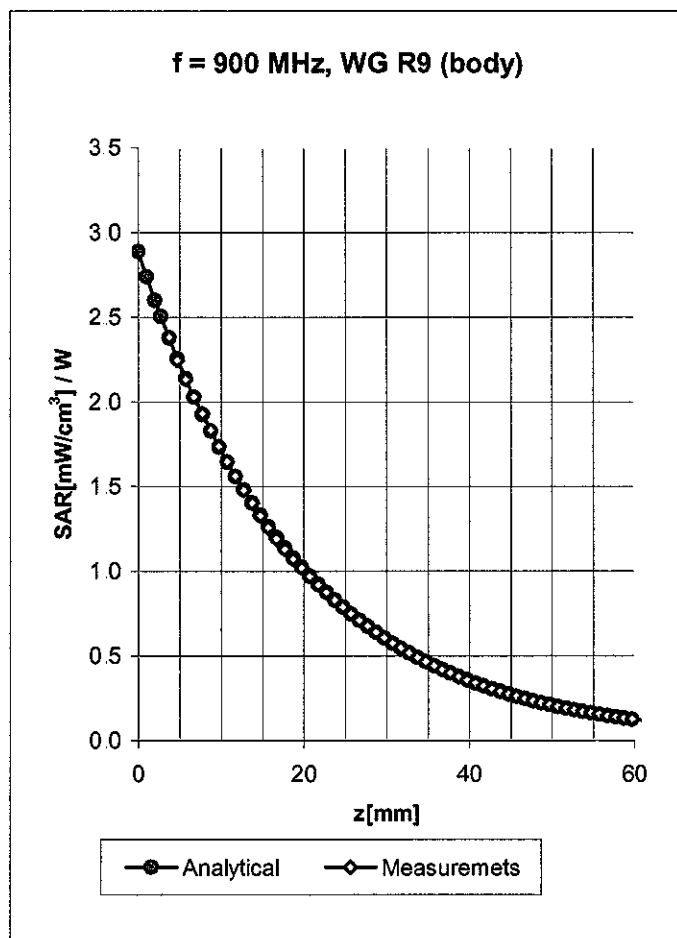

Body 900 MHz $\epsilon_r = 55.0 \pm 5\%$ $\sigma = 1.05 \pm 5\%$ mho/m

Valid for f=800-1000 MHz with Body Tissue Simulating Liquid according to OET 65 Suppl. C



ConvF X	6.7 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	6.7 $\pm 9.5\%$ (k=2)	Alpha	0.30
ConvF Z	6.7 $\pm 9.5\%$ (k=2)	Depth	2.79

Body 1800 MHz $\epsilon_r = 53.3 \pm 5\%$ $\sigma = 1.52 \pm 5\%$ mho/m

Valid for f=1710-1910 MHz with Body Tissue Simulating Liquid according to OET 65 Suppl. C

ConvF X	4.9 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	4.9 $\pm 9.5\%$ (k=2)	Alpha	0.49
ConvF Z	4.9 $\pm 9.5\%$ (k=2)	Depth	2.74

Conversion Factor Assessment

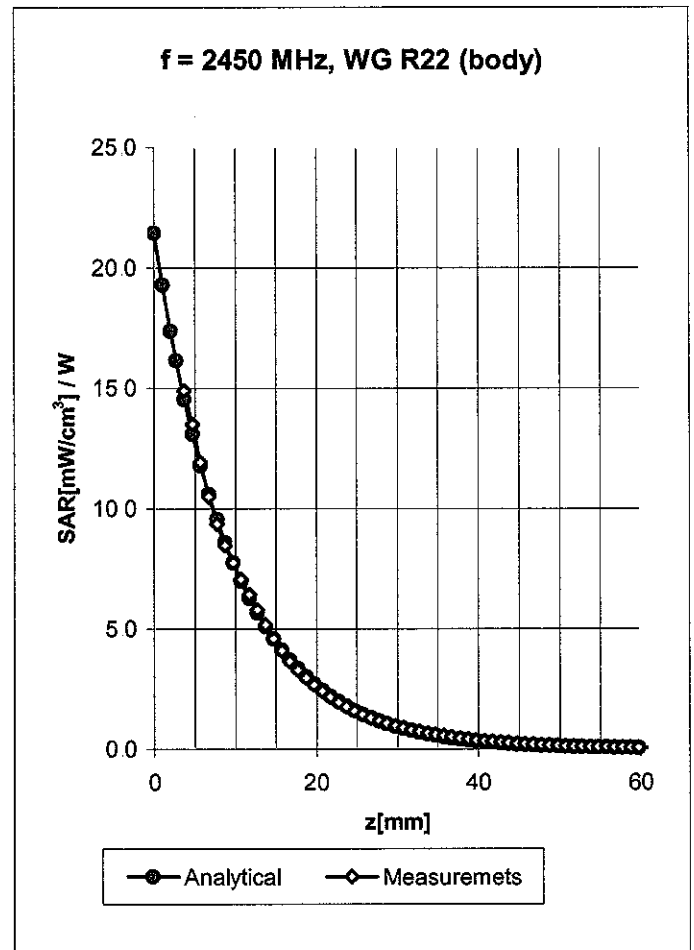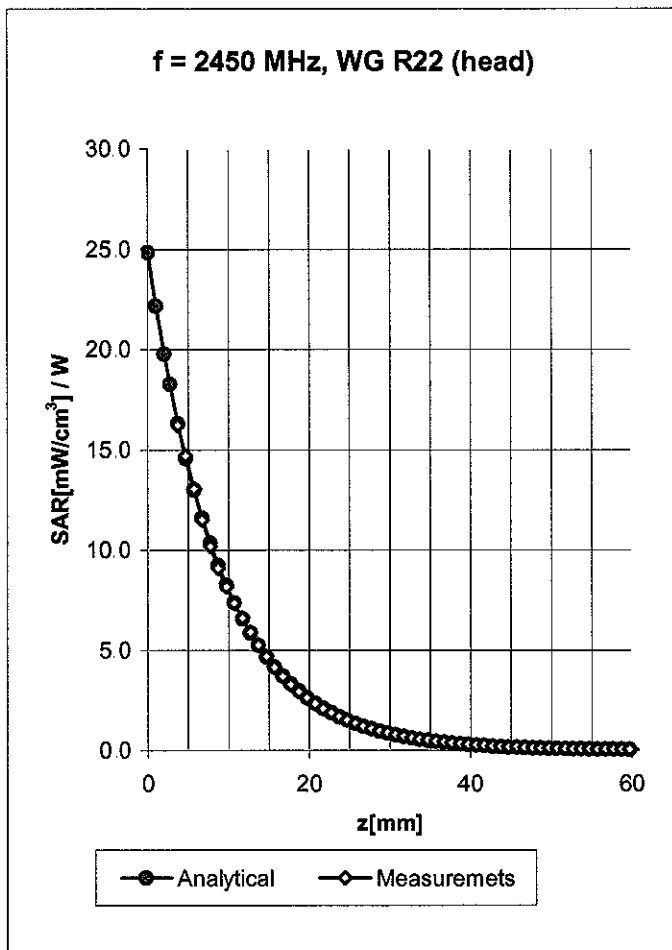

Head 2450 MHz $\epsilon_r = 39.2 \pm 5\%$ $\sigma = 1.80 \pm 5\%$ mho/m

Valid for f=2400-2500 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	4.7 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	4.7 $\pm 9.5\%$ (k=2)	Alpha	0.86
ConvF Z	4.7 $\pm 9.5\%$ (k=2)	Depth	2.00


Body 2450 MHz $\epsilon_r = 52.7 \pm 5\%$ $\sigma = 1.95 \pm 5\%$ mho/m

Valid for f=2400-2500 MHz with Body Tissue Simulating Liquid according to OET 65 Suppl. C

ConvF X	4.3 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	4.3 $\pm 9.5\%$ (k=2)	Alpha	0.90
ConvF Z	4.3 $\pm 9.5\%$ (k=2)	Depth	1.79

Deviation from Isotropy in HSL

Error (θ, ϕ), $f = 900$ MHz

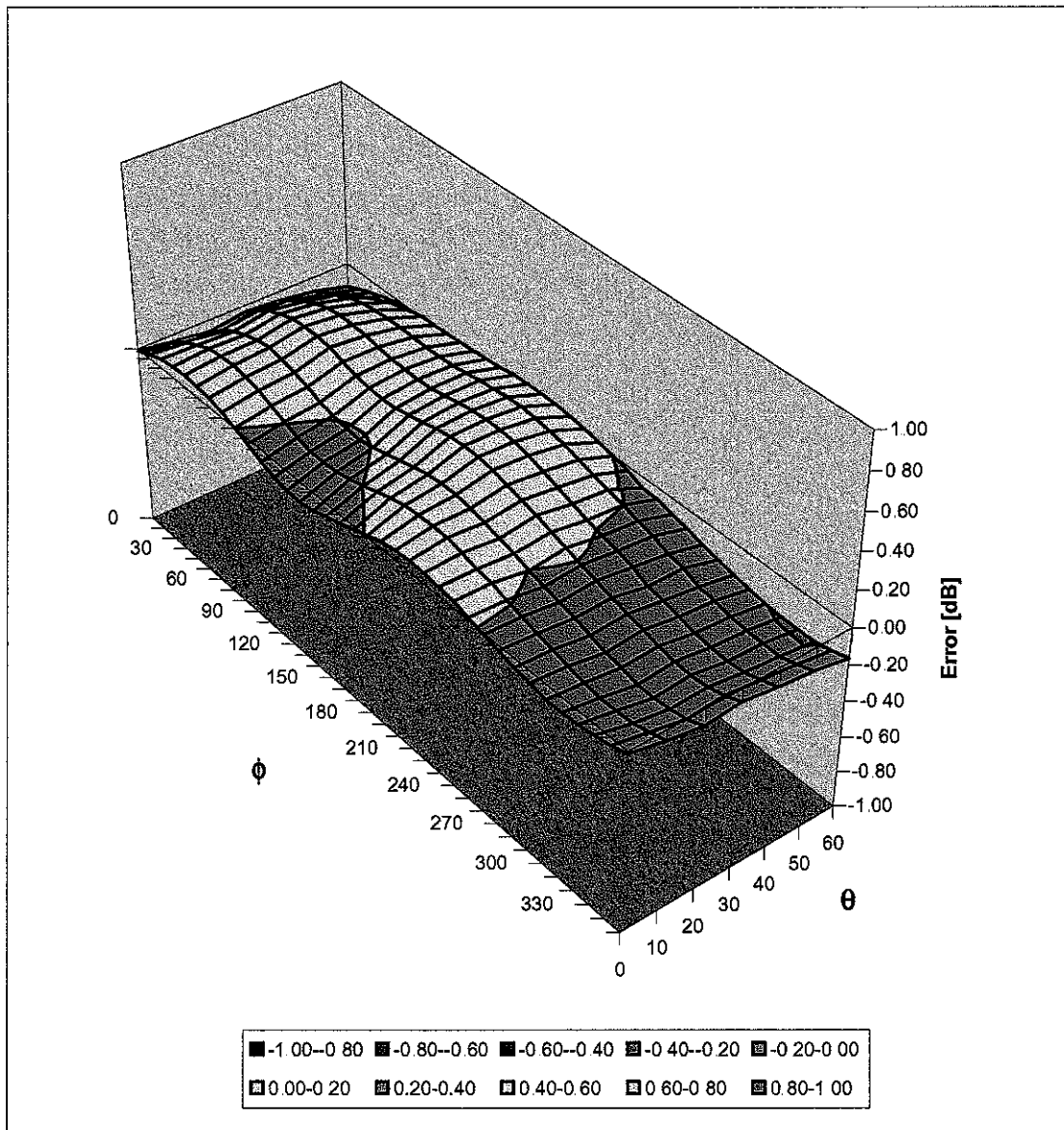

Additional Conversion Factors

for Dosimetric E-Field Probe

Type:

ET3DV6

Serial Number:

1576

Place of Assessment:

Zurich

Date of Assessment:


September 1, 2003

Probe Calibration Date:

August 29, 2003

Schmid & Partner Engineering AG hereby certifies that conversion factor(s) of this probe have been evaluated on the date indicated above. The assessment was performed using the FDTD numerical code SEMCAD of Schmid & Partner Engineering AG. Since the evaluation is coupled with measured conversion factors, it has to be recalculated yearly, i.e., following the re-calibration schedule of the probe. The uncertainty of the numerical assessment is based on the extrapolation from measured value at 900 MHz or at 1800 MHz.

Assessed by:

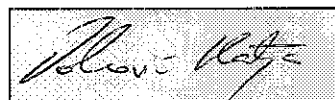

Dosimetric E-Field Probe ET3DV6 SN:1576

Conversion factor ($\pm$ standard deviation)

450 MHz ConvF $7.9 \pm 8\%$

$\epsilon_r = 43.5 \pm 5\%$
 $\sigma = 0.87 \pm 5\% \text{ mho/m}$
(head tissue)

450 MHz ConvF $8.1 \pm 8\%$

$\epsilon_r = 56.7 \pm 5\%$
 $\sigma = 0.94 \pm 5\% \text{ mho/m}$
(body tissue)